

# Operating Manual and Programming Reference

---

## **Model SR785 Dynamic Signal Analyzer**



1290-D Reamwood Avenue  
Sunnyvale, CA 94089 U.S.A.  
Phone: (408) 744-9040 • Fax: (408) 744-9049  
Email: [info@thinkSRS.com](mailto:info@thinkSRS.com) • [www.thinkSRS.com](http://www.thinkSRS.com)

Copyright © 1995, 1996, 1998  
Stanford Research Systems, Inc.  
All Rights Reserved

---

**Certification**

Stanford Research Systems certifies that this product met its published specifications at the time of shipment. Stanford Research Systems further certifies that its calibration measurements are traceable to the United States National Institute of Standards and Technology (NIST).

---

**Warranty**

This Stanford Research Systems product is warranted against defects in materials and workmanship for a period of one (1) year from the date of shipment.

---

**Service**

For warranty service or repair, this product must be returned to a Stanford Research Systems authorized service facility. Contact Stanford Research Systems or an authorized representative before returning this product for repair.

Information in this document is subject to change without notice.

Copyright © Stanford Research Systems, Inc., 1995, 1996, 1998. All rights reserved.

Stanford Research Systems, Inc.  
1290-D Reamwood Avenue  
Sunnyvale, California 94089

Printed in U.S.A.

---

## Safety and Preparation For Use

### WARNING!



Dangerous voltages, capable of causing injury or death, are present in this instrument. Use extreme caution whenever the instrument cover is removed. Do not remove the cover while the unit is plugged into a live outlet.

### Caution



This instrument may be damaged if operated with the LINE VOLTAGE SELECTOR set for the wrong AC line voltage or if the wrong fuse is installed.

### Line Voltage Selection

The SR785 operates from a 100V, 120V, 220V, or 240V nominal AC power source having a line frequency of 50 or 60 Hz. Before connecting the power cord to a power source, verify that the LINE VOLTAGE SELECTOR card, located in the rear panel fuse holder, is set so that the correct AC input voltage value is visible.

Conversion to other AC input voltages requires a change in the fuse holder voltage card position and fuse value. Disconnect the power cord, open the fuse holder cover door and rotate the fuse-pull lever to remove the fuse. Remove the small printed circuit board and select the operating voltage by orienting the printed circuit board so that the desired voltage is visible when pushed firmly into its slot. Rotate the fuse-pull lever back into its normal position and insert the correct fuse into the fuse holder.

### Line Fuse

Verify that the correct line fuse is installed before connecting the line cord. For 100V/120V, use a 1.5 Amp fuse. For 220V/240V, use a 3/4 Amp fuse.

### Line Cord

The SR785 has a detachable, three-wire power cord for connection to the power source and to a protective ground. The exposed metal parts of the instrument are connected to the outlet ground to protect against electrical shock. Always use an outlet which has a properly connected protective ground.

### Service



Do not attempt to service or adjust this instrument unless another person, capable of providing first aid or resuscitation, is present.



Do not install substitute parts or perform any unauthorized modifications to this instrument. Contact the factory for instructions on how to return the instrument for authorized service and adjustment.

### Fan

The fans in the SR785 are required to maintain proper operation. Do not block the vents in the chassis or the unit may not operate properly.



# Contents

Safety and Preparation For Use	i
Contents	iii
Table of Figures	vii
Features	ix
Specifications	xi
<b>Chapter 1 Getting Started</b>	
General Installation	1-2
Front Panel Quick Start	1-3
Things To Watch Out For	1-5
Analyzing a Sine Wave	1-7
Measuring a Frequency Response Function	1-11
Linking (Advanced Operation)	1-15
Triggering and the Time Record	1-19
Octave Analysis	1-23
Capture	1-27
Waterfall Display	1-34
Swept Sine Measurement	1-41
Saving and Recalling	1-47
User Math Functions	1-52
Limit Testing	1-56
Exceedance Statistics	1-60
<b>Chapter 2 Analyzer Basics</b>	
Measurement Groups	2-4
What is an FFT?	2-6
FFT Frequency Spans	2-8
FFT Time Record	2-10
FFT Windowing	2-12
FFT Measurements	2-16
Views	2-21
FFT Averaging	2-24
Real Time Bandwidth and Overlap	2-27
Waterfall Display	2-30
Capture Buffer	2-33
The Source	2-36
Correlation Analysis	2-38
Octave Analysis	2-41
Swept Sine Measurements	2-47
Order Analysis	2-55
Time/Histogram Measurements	2-58
Trace Storage	2-60
User Math Functions	2-61
Signal Inputs	2-65

Input Connections	2-68
Intrinsic Noise Sources	2-70
External Noise Sources	2-71
Curve Fitting and Synthesis	2-75

## Chapter 3 Operation

Overview	3-3
Front Panel Connectors	3-6
Rear Panel Connectors	3-8
Screen Display	3-11
Status Indicators	3-18
Keypad	3-23
Normal and Alternate Keys	3-23
Menu Keys	3-24
Entry Keys	3-25
Control Keys	3-28
Function Keys	3-34
Macros	3-38

## Chapter 4 Menus

<b>Frequency Menus</b>	<b>4-7</b>
FFT Frequency Menu	4-7
Octave Frequency Menu	4-12
Swept Sine Frequency Menu	4-15
Order Frequency Menu	4-19
Time/Histogram Frequency Menu	4-23
<b>Display Setup Menu</b>	<b>4-25</b>
<b>Display Options Menu</b>	<b>4-46</b>
<b>Marker Menu</b>	<b>4-50</b>
Normal Marker Menu	4-54
Harmonic Marker Menu	4-58
Sideband Marker Menu	4-61
Band Marker Menu	4-64
Frequency/Damping Marker Menu	4-66
<b>Source Menus</b>	<b>4-67</b>
Sine Source Menu	4-71
Chirp Source Menu	4-73
Noise Source Menu	4-75
Arbitrary Source Menu	4-78
Swept Sine Source Menu	4-83
<b>Input Menu</b>	<b>4-87</b>
Input Configuration Menu	4-89
Transducer Parameter Menu	4-92
Tachometer Input Menu	4-95
Playback Input Menu	4-97

<b>Trigger Menu</b>	<b>4-100</b>
<b>Average Menus</b>	<b>4-107</b>
FFT, Correlation, and Order Average Menu	4-107
Octave Average Menu	4-115
Swept Sine Average Menu	4-119
<b>User Math Menu</b>	<b>4-122</b>
<b>Window Menu</b>	<b>4-131</b>
<b>Waterfall Menu</b>	<b>4-135</b>
<b>Capture Menu</b>	<b>4-143</b>
<b>Analysis Menu</b>	<b>4-147</b>
Data Table Analysis Menu	4-150
Limit Testing Analysis Menu	4-152
Marker Statistics Analysis Menu	4-156
Exceedance Statistics Analysis Menu	4-159
Curve Fit Menu	4-161
<b>Disk Menu</b>	<b>4-167</b>
Nodal Degree-of-Freedom Menu	4-171
Recall Settings Menu	4-174
Disk Buffers Menu	4-177
Disk Upkeep Menu	4-182
<b>Output Menu</b>	<b>4-184</b>
<b>System Menu</b>	<b>4-191</b>
System Remote Menu	4-194
System Preferences Menu	4-196
System Date/Time Menu	4-199
System Diagnostics Menu	4-200
Edit Macro Menu	4-203

## Chapter 5 Programming

Index of Commands	5-2
Alphabetical List of Commands	5-12
Introduction	5-21
Command Syntax	5-25
Frequency Commands	5-28
Display Setup Commands	5-37
Display Options Commands	5-43
Marker Commands	5-45
Source Commands	5-53
Input Commands	5-60
Trigger Commands	5-66
Average Commands	5-68
User Math Commands	5-74
Window Commands	5-78
Waterfall Commands	5-80
Capture Commands	5-83

Memory Allocation Commands	5-84
Data Table Commands	5-85
Limit Test Commands	5-87
Marker Statistics Commands	5-89
Exceedance Statistics Commands	5-90
Disk Commands	5-94
Output Commands	5-97
System Commands	5-100
Front Panel Commands	5-102
Data Transfer Commands	5-107
Interface Commands	5-117
Nodal Degree-of-Freedom Commands	5-118
Status Reporting Commands	5-120
Status Word Definitions	5-124
Example Program	5-129

## **Chapter 6 File Conversion**

Why File Conversion?	6-2
SR785 File Types	6-3
Supported External File Types	6-4
Using the File Conversion Utility	6-5



## Table of Figures

Figure 2-1 Waterfall Display	2-30
Figure 2-2 Transfer Functions	2-48
Figure 2-3 Capacitive Coupling	2-71
Figure 2-4 Inductive Coupling	2-72
Figure 2-5 Resistive Coupling	2-73
Figure 3-1 Front Panel	3-3
Figure 3-2 Rear Panel	3-8
Figure 3-3 Dual Display Screen	3-11
Figure 3-4 Overlay Display Screen	3-12
Figure 3-5 Vertical Scale Bar	3-13
Figure 3-6 Horizontal Scale Bar	3-14
Figure 3-7 Marker Region	3-15
Figure 3-8 Marker Position Bar	3-15
Figure 3-9 Status Indicator Panel	3-18
Figure 3-10 Front Panel Keypad	3-23



# Features

<b>Measurements</b>	<b>FFT Group</b>		
	FFT (Linear Spectrum)	Time Record	Windowed Time
	Time Capture	Frequency Response	Power Spectrum
	Cross Spectrum	Coherence	Orbit
	User Math Functions		
	<b>Correlation Group</b>		
	Cross Correlation	Autocorrelation	Time Record
	Windowed Time	Time Capture	User Math Functions
	<b>Octave Analysis Group</b>		
	1/1, 1/3, 1/12 Octave	Time Capture	User Math Functions
	$L_{EQ}$	Impulse	
	<b>Swept Sine Group</b>		
Spectrum	Transfer Function	Cross Spectrum	
Normalized Variance	User Math Functions		
<b>Time Histogram Group</b>			
Histogram	Probability Density	Cumulative Density	
Unfiltered Time Record	User Math Functions	Time Capture	
<b>Views</b>	Log Magnitude	Linear Magnitude	Magnitude Squared
	Real Part	Imaginary Part	Phase
	Unwrapped Phase	Nichols Plot	Nyquist Plot
<b>Displays</b>	Single	Dual Upper/Lower	Front/Back
	Waterfall with Skew	Zoom and Pan	
<b>FFT Resolution</b>	100, 200, 400 or 800 FFT lines		
<b>FFT Windows</b>	Hanning	Blackman-Harris	Flattop
	Kaiser	Uniform	Force/Exponential
	User defined		
<b>Correlation Windows</b>	$[-T/2..T/2]$	$[-T/4..T/4]$	$[0..T/2]$
<b>Averaging</b>	RMS	Vector	Peak Hold
	Linear or Exponential	Preview Time Records	Equal Confidence
<b>Analysis</b>	Harmonic	Sideband	Band
	Peak Finding	THD/THD+N	Limit Testing
	Data Table	Exceedance ( $L_N$ )	Statistics
	Waterfall Slice	Curve Fitting	Curve Synthesis
<b>User Math</b>	$+$ , $-$ , $*$ , $/$	Conjugate	Magnitude/Phase
	Real/Imaginary	Sqrt	FFT/Inverse FFT
	$j\omega$	Log/Exp	$d/dx$
	Group Delay	A, B, C Wt	

<b>Source Outputs</b>	Sine Offset Chirp	Two Tone White/Pink Noise Burst Chirp	Swept Sine Burst Noise Arbitrary
<b>Trigger Arming</b>	Auto Arm Time Arm	Manual Arm	RPM Arm
<b>Trigger Sources</b>	Continuous External/TTL	Internal Source	External Manual
<b>Time Capture</b>	Capture time data for later analysis (FFT or Octave). Up to 2 Msamples of data can be saved with standard memory, 8 Msamples with optional memory.		
<b>Storage</b>	3.5", 1.44 Mbytes, DOS formatted disk. Save data, setups and hardcopy.		
<b>Hard Copy</b>	Print to dot matrix or LaserJet/InkJet printers. Plot to HPGL or Postscript plotters. Print/Plot on-line (serial, parallel or IEEE-488) or to disk file. GIF, EPS and PCX graphic formats available for disk output.		
<b>File Translation</b>	Native SR785 binary files can be converted to: ASCII, MATLAB® .MAT Files, Universal File Format, and HP SDF v3 Files. SDF and SR780 files can be converted to native SR785 binary format.		
<b>Interfaces</b>	RS232 serial, Centronics parallel and IEEE-488.		
<b>Help</b>	On screen help system provides Operating Manual and Programming Reference on-line.		

# Specifications

Specifications apply after 30 minutes of warm-up and within 2 hours of last auto-offset. All specifications are with 400 line FFT resolution and anti-alias filters enabled unless stated otherwise.

## Frequency

Range	102.4 kHz or 100 kHz (both displays have the same range).
FFT Spans	195.3 mHz to 102.4 kHz or 191 mHz to 100 kHz. The 2 displays can have different spans and start frequencies.
FFT Resolution	100, 200, 400 or 800 lines
Real Time Bandwidth	102.4 kHz (highest FFT span with continuous data acquisition and averaging on both inputs).
Accuracy	25 ppm from 20° to 40°C

## FFT Dynamic Range

Dynamic Range	-90 dBfs typical, -80 dBfs guaranteed (FFT and Octave), 145 dB (Swept Sine). Includes spurs, harmonic and intermodulation distortion and alias products. Excludes alias responses at extremes of span.
Harmonic Distortion	<-80 dB (Single tone in band).
Intermodulation Distortion	<-80 dB (Two tones in band, each <- 6.02 dBfs).
Spurious	<-80 dBfs
Alias Responses	<-80 dBfs (Single tone outside of span, < 0 dBfs, < 1 Mhz).
Full Span FFT Noise Floor	-100 dBfs typical (Input grounded, Input Range > -30 dBV, Hanning window, 64 RMS averages).
Residual DC Response	< -30 dBfs (FFT with Auto Cal On).

## Amplitude Accuracy

Single Channel	$\pm 0.2$ dB (excluding windowing effects).
Cross Channel	$\pm 0.05$ dB (dc to 102.4 kHz) (Transfer Function measurement, both inputs on the same Input Range, RMS averaged).

## Phase Accuracy

Single Channel	$\pm 3.0$ deg relative to External TTL trigger. (-50 dBfs to 0 dBfs, freq < 10.24 kHz) (Center of frequency bin, DC coupled) For Uniform, Blackman-Harris, Hanning, Flattop and Kaiser windows, phase is relative to a cosine wave at the center of the time record. For Force and Exponential windows, phase is relative to a cosine wave at the beginning of the time record.
Cross Channel	$\pm 0.5$ deg (dc to 51.2 kHz) $\pm 1.0$ deg (dc to 102.4 kHz) (Transfer Function measurement, both inputs on the same Input Range, Vector averaged.)

**Signal Inputs**

Number of Inputs	2
Full Scale Input Range	-50 dBV (3.16 mVpk) to +34 dBV (50 Vpk) in 2 dB steps.
Maximum Input Level	57 Vpk
Input Configuration	Single-ended (A) or True Differential (A-B).
Input Impedance	1 M $\Omega$ + 50 pF
Shield to Chassis	Floating Mode: 1 M $\Omega$ + 0.01 $\mu$ F Grounded Mode: 50 $\Omega$ Shields are always grounded in differential input (A-B).
Maximum Shield Voltage	4 Vpk
AC Coupling	-3 dB rolloff at 0.16 Hz.
CMRR	90 dB at 1 kHz (Input Range < 0 dBV). 80 dB at 1 kHz (Input Range <10 dBV). 50 dB at 1 kHz (Input Range $\geq$ 10 dBV).
ICP Signal Conditioning	Current Source: 4.8 mA Open Circuit Voltage +26 V
A-weight Filter	ANSI Standard S1.4-1983; 10 Hz to 25.6 kHz, Type 0 Tolerance.
Crosstalk	<-145 dB below signal, (Input to Input and Source to Inputs, 50 $\Omega$ receiving input source impedance).
Input Noise	<10 nVrms/ $\sqrt$ Hz (< -160 dBVrms/ $\sqrt$ Hz) above 200 Hz.

**Trigger Input**

Modes	Continuous, Internal, External, or External TTL.
Internal	Level adjustable to $\pm$ 100% of input scale. Positive or Negative slope. Minimum Trigger Amplitude: 5% of input range
External	Level adjustable to $\pm$ 5V in 40 mV steps. Positive or Negative slope. Input Impedance: 1 M $\Omega$ Max Input: $\pm$ 5V Minimum Trigger Amplitude: 100 mV
External TTL	Requires TTL level to trigger (low<0.7V, high>3.0V).
Post-Trigger	Measurement record is delayed up to 100,000 samples after the trigger.
Pre-Trigger	Measurement record starts up to 8000 samples prior to the trigger.

**Tachometer Input**

Modes	External, or External TTL.
External	Level adjustable to $\pm$ 5V in 4 mV steps. or $\pm$ 25V in 20 mV steps. Positive or Negative slope.

**Time Capture**

Mode	Continuous real time data recording to memory.
Maximum Rate	262,144 samples/sec for both inputs. Lower rates may be used for longer capture.

Maximum Capture Length 2M samples standard,  
8M samples optional.

**Octave Analysis**

Standards Conforms to ANSI S1.11-1986, Order 3, Type 1-D.  
 Frequency Range Band centers:  
 Single Channel  
 1/1 Octave 0.125 Hz - 32 kHz  
 1/3 Octaves 0.100 Hz - 40 kHz  
 1/12 Octaves 0.091 Hz - 12.34 kHz  
 Two Channels  
 1/1 Octave 0.125 Hz - 16 kHz  
 1/3 Octaves 0.100 Hz - 20 kHz  
 1/12 Octaves 0.091 Hz - 6.17 kHz  
 Accuracy < 0.2 dB (1 second stable average, single tone at band center).  
 Dynamic Range 80 dB (1/3 Octave, 2 second stable average),  
per ANSI S1.11-1986.  
 Sound Level Exponential time averaged broadband power (L),  
per ANSI S1.4-1983, Type 0.  
Broadband Impulse and Peak power, per  
IEC 651-1979, Type 0.  
Sum of octave bands total power.

**Source Output**

Amplitude Range 0.1 mVpk to 5 Vpk  
 Amplitude Resolution 0.1 mVpk  
 DC Offset <10.0 mV (typical)  
 Output Impedance < 5 Ω; ±100 mA peak output current.

**Sine**

Amplitude Accuracy ±1% of setting, 0 Hz to 102.4 kHz  
 0.1 Vpk to 5.0 Vpk, high impedance load.  
 Offset 0 V to ±5 V, max output ±5 V (ac+dc).  
 Harmonics, Sub-Harmonics and Spurious Signals 0.1 Vpk to 5 Vpk, 0 V offset,  
<-80 dBc (fundamental < 30 kHz),  
<-74 dBc (fundamental > 30 kHz).

**Two Tone**

Amplitude Accuracy ±1% of setting, 0 Hz to 102.4 kHz  
 0.1 Vpk to 5 Vpk, high impedance load.  
 Offset 0 V to ±5 V, max output ±5 V (ac+dc).  
 Harmonics, Sub-Harmonics and Spurious Signals 0.1 Vpk to 2.5 Vpk, 0 V offset,  
< -80 dBc (fundamental < 30 kHz),  
<-74 dBc (fundamental > 30 kHz).

**White Noise**

Time Record Continuous or Burst  
 Bandwidth DC to 102.4 kHz or limited to analysis span.  
 Flatness <0.25 dB pk-pk (typical), <1.0 dB pk-pk (max),

(5000 rms averages).

**Pink Noise**

Time Record Continuous or Burst  
Bandwidth DC to 102.4 kHz  
Flatness <2.0 dB pk-pk, 20 Hz - 20 kHz,  
(measured using averaged 1/3 Octave Analysis).

**Chirp**

Time Record Continuous or Burst  
Output Sine sweep across the FFT span.  
Flatness  $\pm 0.25$  dB pk-pk, Amplitude = 1.0 Vpk.  
Crest Factor 1.94 (5.77 dB)

**Swept Sine**

Auto Functions Source Level, Input Range and Frequency Resolution.  
Dynamic Range 145 dB  
Offset 0 V to  $\pm 5$  V, max output  $\pm 5$  V (ac+dc).

**Arbitrary**

Amplitude Range  $\pm 5$ V  
Record Length 2M samples (playback from Arbitrary Waveform memory or  
a Capture buffer),  
4M and 8M samples optional.  
Variable output sample rate.

**General**

Monitor Monochrome CRT, 800H by 600V resolution.  
Interfaces IEEE-488, RS232 and Printer interfaces standard.  
All instrument functions can be controlled through the  
IEEE-488 and RS232 interfaces. A PC (XT) keyboard input is  
provided for additional flexibility.  
Hardcopy Print to dot matrix and HP LaserJet/InkJet compatible printers.  
Plot to HPGL or Postscript plotters. Print/Plot to RS232 or  
IEEE-488 interfaces or to disk file. Additional file formats  
include GIF, PCX and EPS.  
Disk 3.5 inch DOS compatible format, 1.44 Mbytes capacity.  
Storage of data, setups and hardcopy.  
Preamp Power Power connector for SRS preamplifiers.  
Power 70 Watts, 100/120/220/240 VAC, 50/60 Hz.  
Dimensions 17"W x 8"H x 22"D  
Weight 56 lb.  
Warranty One year parts and labor on materials and workmanship.



## Chapter 1

# Getting Started

These example measurements are designed to acquaint the first time user with the SR785 Dynamic Signal Analyzer. They provide a foundation for understanding how to use the SR785. For a more complete overview of the instrument and its capabilities, refer to the ‘Analyzer Basics’ and ‘Operation’ sections of this manual.

Many of the examples use the test filter enclosed with this manual. The filter is a simple twin-tee 1 kHz passive notch filter. This filter provides an interesting frequency response for these measurements.

## In This Chapter

<b>General Installation</b>	<b>1-2</b>	<b>Analyzing a Sine Wave</b>	<b>1-7</b>
Caution	1-2	<b>Measuring a Frequency Response</b>	
Line Voltage Selection	1-2	<b>Function</b>	<b>1-11</b>
Line Fuse	1-2	<b>Linking (Advanced Operation)</b>	<b>1-15</b>
Line Cord	1-2	<b>Triggering and the Time Record</b>	<b>1-19</b>
Power Switch	1-2	<b>Octave Analysis</b>	<b>1-23</b>
Screen Brightness	1-2	<b>Capture</b>	<b>1-27</b>
Fan	1-2	<b>Waterfall Display</b>	<b>1-34</b>
<b>Front Panel Quick Start</b>	<b>1-3</b>	<b>Swept Sine Measurement</b>	<b>1-41</b>
[Hardkeys]	1-3	<b>Saving and Recalling</b>	<b>1-47</b>
<Softkeys>	1-3	<b>User Math Functions</b>	<b>1-52</b>
[Alt] Keys	1-3	<b>Limit Testing</b>	<b>1-56</b>
Knob	1-3	<b>Exceedance Statistics</b>	<b>1-60</b>
Help	1-4		
<b>Things To Watch Out For</b>	<b>1-5</b>		
Start	1-5		
Live Display	1-5		
Narrow Span	1-5		
Low Detection Frequency	1-5		
Averaging	1-5		
Triggering	1-5		
Scaling and Ranging	1-6		
Local	1-6		
Reset	1-6		

## General Installation

### Caution



This instrument may be damaged if operated with the LINE VOLTAGE SELECTOR set for the wrong AC line voltage or if the wrong fuse is installed.

### Line Voltage Selection

The SR785 operates from a 100V, 120V, 220V, or 240V nominal AC power source having a line frequency of 50 or 60 Hz. Before connecting the power cord to a power source, verify that the LINE VOLTAGE SELECTOR card, located in the rear panel fuse holder, is set so that the correct AC input voltage value is visible.

Conversion to other AC input voltages requires a change in the fuse holder voltage card position and fuse value. Disconnect the power cord, open the fuse holder cover door and rotate the fuse-pull lever to remove the fuse. Remove the small printed circuit board and select the operating voltage by orienting the printed circuit board so that the desired voltage is visible when pushed firmly into its slot. Rotate the fuse-pull lever back into its normal position and insert the correct fuse into the fuse holder.

### Line Fuse

Verify that the correct line fuse is installed before connecting the line cord. For 100V/120V, use a 1.5 Amp fuse. For 220V/240V, use a 3/4 Amp fuse.

### Line Cord

The SR785 has a detachable, three-wire power cord for connection to the power source and to a protective ground. The exposed metal parts of the instrument are connected to the outlet ground to protect against electrical shock. Always use an outlet which has a properly connected protective ground.

### Power Switch

The power switch is on the rear panel. Turn the unit on by depressing the upper half of the power switch. The green power LED on the front panel indicates that the unit is powered.

### Screen Brightness

If the screen is too dark or too bright, adjust the brightness using the Brighter and Dimmer buttons below the softkeys (below right of the display). Do not set the brightness higher than necessary.

### Fan

The fans in the SR785 are required to maintain proper operation. Do not block the vents in the chassis or the unit may not operate properly.

---

## Front Panel Quick Start

There are two types of front panel keys which are referenced in this section. Hardkeys are those keys with labels printed on them. Their function is determined by the label and does not change. Hardkeys are referenced by brackets like this - [Hardkey]. Softkeys are the ten gray keys along the right edge of the screen. Their function is labeled by a menu box displayed on the screen next to the key. Softkey functions change depending upon the menu and instrument configuration. Softkeys are referenced as the <Softkey>.

### [Hardkeys]

The keypad consists of four groups of hardkeys (keys with printed labels).

The ENTRY keys are used to enter numeric parameters which have been highlighted by a softkey. The MENU keys select a menu of softkeys. Pressing a menu key will change the menu boxes which are displayed next to the softkeys. Each menu presents a group of similar or related parameters and functions. The CONTROL keys start and stop data acquisition, toggle the active display and link parameters and functions. These keys are not in a menu since they are used frequently and within any menu. The FUNCTION keys perform common functions such as Auto Scale and Auto Range. These keys can be accessed at any time.

### <Softkeys>

The SR785 has a menu driven user interface. The Menu keys each display a menu of softkeys. The softkeys are at the right of the video display and have different functions depending upon the displayed menu.

There are three types of softkeys - buttons, lists and numeric values. A button performs a function, such as <Full Span>. A list presents a list of choices or options in the entry field (at the top of the screen). Use the knob to make a selection and press [Enter]. <Measurement> is an example of a list. A numeric value presents the current value in the entry field and awaits numeric entry. Enter a new value with the entry keys and press [Enter]. <Start Freq> is an example of a numeric value.

### [Alt] Keys

The [Alt] key is a special key which has no meaning by itself but instead modifies the meaning of another keypress or knob turn. Pressing the [Alt] key toggles the state of the ALT indicator at the top of the screen. Pressing a control key while the ALT indicator is lit selects the "alternate" function labeled underneath the key instead of the normal function. Turning the knob when ALT is lit affects the way the markers move in the waterfall, data tables, and other displays.

### Knob

The knob normally moves the markers within the displays. If a parameter has been highlighted by its softkey, the knob adjusts the parameter. List parameters are most easily modified with the knob. Numeric parameters may also be adjusted with the knob.

Knob list selections are referenced in parenthesis like (Hanning).

## Help



Enter the on screen help system by pressing [Help/Local]. Help on any hardkey or softkey is available simply by pressing the key. Press [1] for the Help Index. Press [0] to exit the help system and return to normal operation.

## Things To Watch Out For

If the analyzer is on but doesn't seem to be taking data, there are a number of things to check.

### Start

Press the [Start/Reset] key to start the measurements. Make sure the Run/Pause indicator at the top of the screen displays 'RUN' instead of 'PAUSE'. Note that in many cases when settings are changed using the menus, the new settings will not take effect until [Start/Reset] is pressed.

### Live Display

If the displays are showing recalled trace data, they are Off-Line and do not display the live measurement data. Set the Display to Live instead of Off-Line (in the [Display Options] menu).

### Narrow Span

If the FFT span is very narrow, the time record is very long (up to 1000's of seconds). Completely new data is available only every time record. Change the Time Record Increment in the [Average] menu) to display overlapped data more often.

### Low Detection Frequency

Swept Sine measurements at very low frequencies ( $\ll 1$  Hz) take a very long time (at least 2 cycles and maybe longer). Do not set the sweep Start to a very low frequency to measure the DC response.

Octave measurements with a very low starting band take a long time to settle before the first measurement is valid. The settling time is related to the bandwidth of the lowest octave band. If the Lowest Band is less than 1 Hz, the settling time can be very long.

### Averaging

Very long averaging times for any measurement may give the appearance that the display does not update. Check the FFT Number Of Averages, the Octave Integration Time or the Swept Sine Integration Time.

When Linear averaging is on, the measurement is paused after the average is completed (unless triggered or waterfall storage is on). Press [Start/Reset] to take another average.

### Triggering

If the analyzer is waiting for a trigger, the Trig Wait indicator at the top of the screen is on.

If the measurement is not meant to be triggered, make sure the Trigger Mode is Auto Arm and the Trigger Source continuous.

If the measurement is meant to be triggered, make sure that the correct Trigger Source is selected and the Trigger Level is appropriate for the trigger signal.

Check that the Trigger Mode is set to Auto Arm. If the Trigger Mode is Manual Arm, then the analyzer will only trigger once and then wait for the next Manual Arm command.

### **Scaling and Ranging**

Check that the inputs are not completely overloaded by using [Auto Range Ch1] and [Auto Range Ch2].

Scale the display to show the entire range of the data with [Auto Scale A] and [Auto Scale B].

### **Local**

Make sure that the analyzer is not in the REMOTE state where the computer interfaces have setup the instrument and locked out the front panel. Press the [Local/Help] key to restore local control.

### **Reset**

If the analyzer still seems to function improperly, turn the power off and turn it back on while holding down the [<-] (backspace) key. This will reset the analyzer into the default configuration. The analyzer should power on running and taking measurements.

# Analyzing a Sine Wave

This measurement investigates the spectrum of a 1.024 kHz sine wave. You will use the SR785 source to provide the sine signal (or you can use a function generator capable of providing a 1.024 kHz sine wave at a level of 100 mV to 1 V, such as the SRS DS345). The actual settings of the generator are not important since you will be using the SR785 to measure and analyze its output.

<p>1. Press [System]</p> <p>Press &lt;Preset&gt;</p> <p>Press [Enter] to confirm Preset.</p>	<p>Display the System menu.</p> <p>Preset returns the unit to its default settings.</p> <p>Preset requires confirmation to prevent accidental reset. Wait until the self tests are completed.</p>
<p>2. Connect the Source Output to the Channel 1 A Input. (Or connect a function generator's output to the Ch1 A Input of the analyzer.)</p>	<p>Setup to analyze the source output.</p> <p>(The input impedance of the analyzer is 1 MΩ. The generator may require a terminator. Many generators have either a 50 Ω or 600 Ω output impedance. Use the appropriate feedthrough termination if necessary. In general, not using a terminator means that the output amplitude will not agree with the generator setting and the distortion may be greater than normal.)</p>
<p>3. Press [Source]</p> <p>Press &lt;Sine&gt;</p> <p>Press &lt;Frequency 1&gt;</p> <p>Press [1] [. ] [0] [2] [4], select (kHz) with the knob, and press [Enter].</p> <p>Press &lt;On&gt;</p> <p>(Or turn on the generator, set the frequency to 1.024 kHz and the amplitude to approximately 1 Vrms.)</p>	<p>Select the Source menu.</p> <p>Choose Sine output.</p> <p>Adjust the output Frequency.</p> <p>Enter 1.024 kHz for the Frequency. Enter the value with the numeric entry keys. Select the units with the knob. Enter the new value with the [Enter] key.</p> <p>Turn the source on. When the instrument is turned on, the source is always off. (Setup the function generator for 1.024 kHz sine output.)</p>

## 1-8 Analyzing a Sine Wave

---

4. Press [Auto Range Ch1]	Let the analyzer automatically set the Input Range to agree with the signal (either from the Source or function generator). Note that the Ch1 Input Range readout at the top of the screen is displayed in inverse when Ch1 Auto Range is on.
5. Press [Freq]  Press <Span>  Use the knob to adjust the Span to 6.4 kHz and press [Enter].	Select the Frequency menu.  Adjust the FFT Span.  Set the Span to display the signal and its first few harmonics.
6. Press [Display Options]  Press <Format>  Select (Single) with the knob and press [Enter].	Select the Display Options menu.  Choose a new Display Format.  Select the desired option from the displayed list and press [Enter]. Single Display Format shows a single large graph.
7. Press [Auto Scale A]  Press [Marker Max]	Automatically scale DisplayA (the active display) to show the entire range of the measurement.  This moves the Marker to the maximum data point in the active display (A). The Marker should now be on the 1.024 kHz signal. The Marker Position shown above the graph displays the frequency and amplitude of the signal.
8. Use the knob to move the Marker around. Take a look at some of the harmonics.	The knob normally adjusts the Marker Position within the active display (DisplayA in this case). If a menu box is highlighted with a softkey, the knob adjusts the selected parameter shown in the entry field at the top of the screen.
9. Let's look at the fundamental only.  Press [Span Down] twice to decrease the Span to 1.6 kHz. The Stop Frequency shown at the bottom right of the graph should read 1.6 kHz.  Press [Marker Max]	You can also use the [Span Up] and [Span Down] keys to adjust the Span.  This isolates the 1.024 kHz fundamental frequency. You may notice that the spectrum takes a noticeable time to settle at this last span. This is because the time record is 250 ms long.  Move the Marker to the peak.



<p>Press [Marker Center]</p>	<p>This sets the span Center Frequency to the Marker Position (for the active display). The signal will be at the center of the span. Further adjustments to the span will keep the center frequency fixed.</p>
<p>10. Let's look at the signal distortion.</p> <p>Press [Freq]</p> <p>Press &lt;Span&gt;</p> <p>Enter [1] [2] [.] [8], select (kHz) with the knob, and press [Enter].</p> <p>Press [Auto Scale A]</p>	<p>Select the Frequency menu.</p> <p>Adjust the Span.</p> <p>You can also use the numeric keypad to enter the span.</p> <p>Enter the 12.8 kHz span numerically. Note that the Center Frequency is no longer 1.024 kHz. This is because a 12.8 kHz span cannot be centered below 6.4 kHz without starting at a negative frequency.</p> <p>Adjust the graph scale and reference to display the entire range of the data. This key can be used at any time.</p>
<p>11. Let's measure some harmonics using the Marker Reference.</p> <p>Press [Marker Max]</p> <p>Press [Marker Ref]</p> <p>Use the knob to move the Marker to the harmonics.</p> <p>Press [Marker Ref]</p>	<p>Move the Marker to the fundamental peak.</p> <p>Set the Marker Offset or Reference to the amplitude of the fundamental. The Marker Position above the graph now reads relative to this offset (~0 dB). This is indicated by the <math>\Delta</math> in front of the Marker Position reading. A small flag shaped symbol is located at the screen location of the reference.</p> <p>The Marker Position shows the distortion peaks relative to the fundamental.</p> <p>Pressing [Marker Ref] again removes the Marker Offset and returns the Marker to absolute readings.</p>
<p>12. Let's have the analyzer measure the distortion.</p> <p>Press [Marker]</p> <p>Press &lt;Mode&gt;</p>	<p>Select the Marker menu.</p> <p>Adjust the Marker Mode.</p>

## 1-10 Analyzing a Sine Wave

---

Select (Harmonic) with the knob and press [Enter].

Press <# Harmonics>

Use the knob to adjust the Number Of Harmonics to 10 and press [Enter].

Press [Marker Max]

Choose the Harmonic Marker for the active display.

Adjust the Number Of Harmonics for analysis.

Enter 10 harmonics.

Move the Marker to the peak (fundamental).

Notice that Harmonic Markers (little triangles) appear on top of all of the harmonic peaks. These indicate which data points are used in the harmonic calculations.

The harmonic calculations are displayed within the menu. THD (total harmonic distortion) is relative to the fundamental. Harmonic power is an absolute measurement of the harmonic power level.

---

This concludes this measurement example. You should have a feeling for the basic operation of the menus, knob and numeric entry, marker movement and some function keys.

---

## Measuring a Frequency Response Function

This example investigates the frequency response of the test filter (enclosed with this manual) using FFT measurements. You will use the SR785 source to provide a broad band chirp and both input channels to measure the input to and output from the device under test.

<p>1. Press [System]</p> <p>Press &lt;Preset&gt;</p> <p>Press [Enter] to confirm Preset.</p>	<p>Display the System menu.</p> <p>Preset returns the unit to its default settings.</p> <p>Preset requires confirmation to prevent accidental reset. Wait until the self tests are completed.</p>
<p>2. Use a BNC Tee to connect the Source Output to the filter input and the Ch1 A Input.</p> <p>Connect the filter output to the Ch2 A Input.</p>	<p>Frequency Response is defined as Ch2 response divided by the Ch1 reference. Thus, Ch1 monitors the filter input (source output) and Ch2 measures the response of the device under test.</p>
<p>3. Press [Source]</p> <p>Press &lt;Chirp&gt;</p> <p>Press &lt;On&gt;</p> <p>Press [Window]</p> <p>Press &lt;Window&gt;</p> <p>Select (Uniform) with the knob and press [Enter].</p>	<p>Select the Source menu.</p> <p>Choose Chirp output. The output is an equal amplitude sine wave at each frequency bin of the FFT spectrum.</p> <p>Turn the source On.</p> <p>Select the Window menu.</p> <p>Adjust the FFT Window function.</p> <p>The Chirp source requires the use of the Uniform window since not all chirp frequency components are present at all points in the time record. The chirp is exactly periodic with the FFT time record and does not 'leak' with the uniform window.</p>
<p>4. Press [Auto Range Ch1]</p> <p>Press [Auto Range Ch2]</p>	<p>Let the analyzer automatically set the Input Ranges to agree with the signals. Note that the Input Range readouts at the top of the screen are displayed in inverse when Auto Range is on.</p>
<p>5. Press [Freq]</p> <p>Press &lt;Span&gt;</p>	<p>Select the Frequency menu.</p> <p>Adjust the FFT Span.</p>

## 1-12 Measuring a Frequency Response Function

---

Use the knob to adjust the Span to 6.4 kHz and press [Enter].	Set the Span to display the filter notch at 1 kHz.  The top display (A) is measuring the filter input and should show a fairly flat spectrum. The bottom display (B) is measuring the filter output and should show a deep notch.  Both displays are measuring absolute signal levels.
6. Press [Display Setup]  Press <Measurement>  Select (Freq. Resp.) with the knob and press [Enter].  Press [Auto Scale A]	Select the Display Setup menu.  Adjust the Measurement of the active display (A).  Choose Frequency Response for the Measurement in DisplayA (top).  Frequency Response is the ratio of the response (Ch2) to the input (Ch1) and is a unitless quantity.  Adjust the scale and reference for DisplayA to show the entire range of the data.
7. Press [Marker]  Press <Width>  Select (Normal) with the knob and press [Enter].  Press <Seeks>  Select (Min) with the knob and press [Enter].  Move the Marker Region with the knob to find the notch frequency and depth. Or press [Marker Min].	Select the Marker menu.  Adjust the Marker Width for DisplayA.  Change to Normal Width (1/2 division).  Adjust what the Marker Seeks within the Marker Region.  Seek the Minimum of the data within the Marker Region.  The Marker Region makes it easy to find narrow peaks and valleys in the graph. The notch should be around 1 kHz and about -60 dB deep.
8. Press [Display Options]  Press <X-Axis>  Press <Link>	Select the Display Options menu.  The graph might look better on a log x axis.  The <Link> key lets you make a choice for both displays. We want both displays on a log x-axis. Note that the linkage indicator at the top of the screen changes from DispA to Link.

<p>Select (Log) with the knob and press [Enter].</p>	<p>Log scale is a common way to display filter response functions. Note that both displays now have log x-axes.</p>
<p>9. Let's show phase response on DisplayB (bottom).</p> <p>Press [Active Display]</p> <p>Press [Display Setup]</p> <p>Press &lt;Measurement&gt;</p> <p>Select (&lt;Freq. Resp.&gt;) with the knob and press [Enter].</p> <p>Press &lt;View&gt;</p> <p>Select (Phase) with the knob and press [Enter].</p> <p>Press [Auto Scale B]</p>	<p>The two displays have separate Measurements.</p> <p>Make DisplayB the active display. The active display has its Marker Position Bar (above the graph) highlighted.</p> <p>Select the Display Setup menu. The setup of DisplayB (the active display) is now shown in the menu.</p> <p>Adjust the Measurement of DisplayB.</p> <p>Choose Frequency Response also.</p> <p>The measured data is a set of complex values which can be viewed in a number of different ways.</p> <p>Choose Phase View to show the phase of the transfer function.</p> <p>Scale DisplayB to show the entire phase transfer function.</p>
<p>10. Let's link the Markers together.</p> <p>Press [Active Display]</p> <p>Press [Marker]</p> <p>Press &lt;Width&gt;</p> <p>Select (Spot) with the knob and press [Enter].</p>	<p>Make DisplayA (top) the active display.</p> <p>Select the Marker menu.</p> <p>Adjust the Marker Width of DisplayA.</p> <p>Change the Marker Width to Spot.</p>
<p>11. Press [Link] and use the knob to move the marker.</p> <p>Press [Enter]</p>	<p>The [Link] key links the two display markers together. This allows simultaneous readout of Transfer Function Magnitude (top) and Phase (bottom).</p> <p>Pressing any key removes the linkage between the markers.</p>

## 1-14 Measuring a Frequency Response Function

---

To permanently link the Markers, go to the Marker menu.

Press [Marker]

Press <Marker>

Select (Link) with the knob and press [Enter].

Move the Marker with the knob.

Select the Marker menu.

Adjust the Marker Type.

Linked Markers move together. Since we changed the DisplayA Marker to Linked, moving the DisplayA Marker moves the DisplayB Marker.

If DisplayB is active, moving its Marker does not move the DisplayA Marker. To do this, change the DisplayB Marker Type to Linked also.

---

This concludes this measurement example. You should have a feeling for the basic operation of two channel measurements and the use of [Active Display].

---

# Linking (Advanced Operation)

This example investigates the test filter (enclosed with this manual) using FFT measurements. You will use the SR785 source to provide a broad band source and both displays to measure the output of the device under test. Display parameter linking and function linking will be explored in greater detail.

<p>1. Press [System]</p> <p>Press &lt;Preset&gt;</p> <p>Press [Enter] to confirm Preset.</p>	<p>Display the System menu.</p> <p>Preset returns the unit to its default settings.</p> <p>Preset requires confirmation to prevent accidental reset. Wait until the self tests are completed.</p>
<p>2. Connect the Source Output to the filter input.</p> <p>Connect the filter output to the Ch2 A Input.</p>	<p>In this example, only the filter output on Channel 2 is required.</p>
<p>3. Press [Source]</p> <p>Press &lt;Noise&gt;</p> <p>Press &lt;On&gt;</p> <p>Press &lt;Type&gt;</p> <p>Select (White) with the knob and press [Enter].</p>	<p>Select the Source menu.</p> <p>Choose Noise output. The output is random noise.</p> <p>Turn the source on.</p> <p>Adjust the Noise Type.</p> <p>This source is White Noise which extends over the entire 0-102.4 kHz frequency range.</p>
<p>4. Press [Auto Range Ch2]</p>	<p>Let the analyzer automatically set the Input Range to agree with the signal. Note that the Input Range readouts at the top of the screen are displayed in inverse when Auto Range is on.</p>
<p>5. Press [Display Setup]</p> <p>Press &lt;Measurement&gt;</p> <p>Select (FFT ch2) with the knob and press [Enter].</p>	<p>Select the Display Setup menu.</p> <p>Change the Measurement of the active display (A).</p> <p>Choose FFT spectrum of Ch2 for the measurement in DisplayA (top). Both displays are independently measuring the filter output spectrum.</p>
<p>6. Press [Input]</p> <p>Press &lt;Analyzer Config&gt;</p>	<p>Select the Input Menu</p> <p>Change the way the input channels are used.</p>

## 1-16 Linking (Advanced Operation)

---

Select (Independent Chan.)	In this example we'll want to set separate frequency spans for the two displays. To do that the analyzer must be in the Independent Channels mode. In the default Dual Channel mode the span is the same for both channels, but dual channel measurements are allowed.
7. Press [Freq]	Select the Frequency menu. The menu shows the frequency parameters for the measurement in DisplayA (active display).
Press <Span>	Highlight the Span. Note that the Link indicator at the top of the screen turns on. This indicates that the highlighted parameter (Span) is linked to both displays. Changing a linked parameter affects both displays at once.
Use the knob to adjust the Span to 12.8 kHz and press [Enter].	Narrow the Span of both displays to show the filter notch at 1 kHz (noisy of course).
Press <Span> again.	Highlight the Span again.
Press [Link]	Pressing [Link] toggles parameter linking off. The Link indicator now shows 'DispA' indicating that this menu box adjusts the span for DisplayA only.
Use the knob to adjust the Span to 3.2 kHz and press [Enter].	Change the Span of DisplayA to 3.2 kHz. The Span of DisplayB remains at 12.8 kHz. The SR785 allows the two displays to have differing Spans and Start frequencies in the Independent Channels mode.  Many parameters affect the displays separately. Linking is a convenient way to adjust the two displays together and keep their settings the same.  The default settings link many measurement parameters, such as frequency and averaging, together as found in many other instruments.
8. Press [Average]	Select the Average menu.
Press <Display Avg>	Select the averaging mode to display. The SR785 always computes all the averaging modes when <Compute Avgs> is set to Yes. The <Display Avg> softkey selects which averaging mode is currently displayed.
Select (RMS) with the knob and press [Enter].	Select RMS Averaging for both displays.



<p>Press &lt;# Avgs&gt;</p> <p>Press [Link]</p> <p>Press [2] [0] and press [Enter].</p>	<p>Change the Number Of Averages for DisplayA.</p> <p>Unlink the Number Of Averages. ‘DispA’ is shown on the Link indicator.</p> <p>Change the Number Of Averages for DisplayA to 20 (instead of 2). DisplayA will average for 10 times as many measurements as DisplayB and be quite a bit smoother.</p>
<p>9. Press [Active Display]</p> <p>Press [Window]</p> <p>Press &lt;Window&gt;</p> <p>Press [Link]</p> <p>Select (Hanning) with the knob and press [Enter].</p>	<p>Let’s change the Window for DisplayB. Make DisplayB the active display.</p> <p>Select the Window menu.</p> <p>The windows are linked by default.</p> <p>Unlink the Window type. ‘DispB’ is shown as the Link indicator.</p> <p>DisplayB is now using the Hanning window while DisplayA is still using the BMH window.</p> <p>So far we have unlinked measurement parameters. Frequency Span, Averaging and Window type affect the actual measurements within the displays. Most analyzers do not allow these measurement parameters to be unlinked.</p> <p>Display parameters, such as Scaling, Views, Units and Marker functions, are usually unlinked. However, linking them can be a convenient way to adjust graph parameters together with a minimum of key presses.</p>
<p>10. Press [Display Setup]</p> <p>Press &lt;Units&gt;</p> <p>Press &lt;Pk Units&gt;</p> <p>Press [Link]</p> <p>Select (rms) with the knob and press [Enter].</p>	<p>Select the Display Setup menu.</p> <p>Enter the Units</p> <p>Select between peak, rms and peak to peak units. Link the Pk Units parameter (default is unlinked).</p> <p>Both displays’ Units become dBVrms with a single parameter entry. The Pk Units remain linked until unlinked with the [Link] key.</p>
<p>11. Press &lt;Return&gt;</p>	<p>Return to the Display Setup Menu.</p>

## 1-18 Linking (Advanced Operation)

---

Press <Y Max>	Let's change the Top Reference of the graphs.
Press [Link]	Link the Y Max parameter (default is unlinked).
Press [-] [1] [0] and press [Enter].	Change the Top Reference of both graphs to -10 dBVrms with a single entry.
Press [Auto Scale A]	Change the scale of DisplayA to center the data.  Another simple way to adjust the scales of both graphs is using the Auto Scale keys.
Press [Link] and then [Auto Scale A]	Pressing [Link] [Auto Scale A] first auto scales DisplayA and then changes the scale of DisplayB to match. This is convenient when you are comparing the two displays.
12. Press [Link] and use the knob to move the Marker. Both Markers move together when linked.	The [Link] key temporarily links the two display Markers together.
Press [Enter]	Pressing any key removes the link between the Markers.  To permanently link the markers, go to the Marker menu and change <Marker> to Link.
Press [Marker Min]	[Marker Min] moves the Marker in the active display (B) to the graph minimum.
Press [Link] then [Marker Min]	Pressing [Link] first moves both Markers to their graph minimums at the same time.  [Link] preceding a function key generally performs the function on both displays at once. [Link] [Auto Scale] matches the active display. [Span Up] and [Span Down] are always linked.
	This concludes this measurement example. You should have a feeling for linking and unlinking and the flexibility of unlinked measurements.

---

# Triggering and the Time Record

This example investigates the trigger and time record. You will need a function generator capable of providing a 100  $\mu$ s wide pulse at 256 Hz with an amplitude of 1V. The output should have a DC level of 0V.

Make sure you have read ‘The FFT Time Record’ in Chapter 2 before trying this exercise.

<p>1. Press [System]</p> <p>Press &lt;Preset&gt;</p> <p>Press [Enter] to confirm Preset.</p>	<p>Display the System menu.</p> <p>Preset returns the unit to its default settings.</p> <p>Preset requires confirmation to prevent accidental reset. Wait until the self tests are completed.</p>
<p>2. Turn on the generator and choose a pulsed output waveform. Set the frequency to 256 Hz, the pulse width to 100 <math>\mu</math>s and the amplitude to 1V. (These settings only need to be approximate.) Make sure that the DC level of the output is near 0V.</p> <p>Connect the generator output to the Ch1 A input of the analyzer.</p>	<p>The input impedance of the analyzer is 1 M<math>\Omega</math>. The generator may require a terminator. Many generators have either a 50 <math>\Omega</math> or 600 <math>\Omega</math> output impedance. Use the appropriate feedthrough termination if necessary. In general, not using a terminator means that the output amplitude will not agree with the generator setting.</p>
<p>3. Press [Input]</p> <p>Press &lt;Input Config&gt;</p> <p>Press &lt;Ch1 Input Range&gt;</p> <p>Press [1.5]. Select (Vpk) with the knob. Now press [Enter]</p>	<p>Select the Input master menu.</p> <p>Select the input configuration submenu</p> <p>Choose an input range that doesn't overload.</p> <p>Set the input range to about 1.5 Vpk. The channel 1 input range indicator at the top of the unit now should read 1.58 Vpk, the closest allowed input range to the value you entered. Adjust the pulse amplitude to that no overloads occur.</p>
<p>4. Press [Active Display]</p> <p>Press [Display Setup]</p> <p>Press &lt;Measurement&gt;</p>	<p>Let's change the Measurement for DisplayB. Make DisplayB the active display.</p> <p>Select the Display Setup menu.</p> <p>Change the Measurement of the active display (B).</p>

## 1-20 Triggering and the Time Record

---

Select (Time1) with the knob and press [Enter].	Choose Time Record of Ch1 for the measurement in DisplayB (bottom). You should see the pulse on the bottom display.
5. Press [Trigger]	Select the Trigger menu.
Press <Trigger Source>	Change the Trigger Source.
Select (Ch 1) with the knob and press [Enter].	Select internal triggering from the Ch 1 input.
Press <Trigger Level>	Adjust the trigger level.
Press [3] [0], select (%) with the knob, and [Enter].	Set the trigger level as a percentage of full scale. Adjust the level for a stable time record in DisplayB.
Press [Auto Scale B]	DisplayB (bottom) should display the pulse waveform at the left edge. In this case, the display shows the signal pulse as a digital oscilloscope would.
6. Press [Window]	Select the Window menu.
	Because the pulse is much shorter than the time record, we need to use the Uniform (or Force) window. The other window functions taper to zero at the start and end of the time record. Always be aware of the effect windowing has on the time record and the FFT.
Press <Window>	Select a new window type for both displays (window type is linked by default).
Select (Uniform) with the knob and press [Enter].	Notice how the spectrum in DisplayA is changed by the Uniform window.
Press [Auto Scale A]	The spectrum in DisplayA is the $\sin x/x$ envelope of a rectangular pulse. The zeroes in the spectrum occur at the harmonics of 1/pulse width (1/100 $\mu$ s or 10 kHz.)
7. Press <Window>	Choose a non-optimum window.
Select (Hanning) with the knob and press [Enter].	Choose the Hanning window. Notice how the spectrum in DisplayA goes away.
Press [Display Setup]	Select the Display Setup menu.

<p>Press &lt;Measurement&gt;</p> <p>Select (WinTime1) with the knob and press [Enter].</p>	<p>Change the Measurement of DisplayB to show the effect of the Hanning window on the time record.</p> <p>The Hanning window is zero at the beginning of the time record and large in the center. This effectively zeroes the signal pulse at the start of the time record leaving nothing in the windowed time record. The FFT operates on this windowed time record and thus the spectrum shows no evidence of the signal pulse.</p>
<p>8. Press [Trigger]</p> <p>Press &lt;Delay1&gt;</p> <p>Press [-] [2], select (ms) with the knob, and [Enter].</p> <p>Press [Auto Scale B]</p>	<p>Select the Trigger menu again.</p> <p>Change the Trigger Delay for the signal on Ch 1.</p> <p>We can get the spectrum back by delaying the time record relative to the trigger so that the pulse is positioned in the center of the time record.</p> <p>A negative delay means that the time record starts before the trigger event. In this case, the time record is about 4 ms long so a delay of -2ms will put the signal pulse in the center of the triggered time record.</p> <p>The trigger delay is specified in time record bins at the current span.</p> <p>Note that the windowed time record in DisplayB shows the signal at the center of the time record.</p> <p>The amplitude of the windowed time record is not the same as the amplitude of the time record itself. This is because the window functions have gain and attenuation at different parts of the time record. The Hanning window is 2.0 at the center so the amplitude of the signal in the windowed time record is twice as large.</p> <p>The Hanning, Flattop, BMH and Kaiser windows are not intended for use with narrow pulse signals. They are used for signals which last the entire time record and normalized as such.</p> <p>The Uniform (and Force) windows have no gain and should be used with pulsed signals such as this.</p>



# Octave Analysis

This example investigates the test filter (enclosed with this manual) using Octave measurements. You will use the SR785 source to provide a broad band source and both displays to measure the output of the device under test.

Refer to ‘Octave Analysis’ in Chapter 2 for more about Octave Analysis measurements.

<p>1. Press [System]</p> <p>Press &lt;Preset&gt;</p> <p>Press [Enter] to confirm Preset.</p>	<p>Display the System menu.</p> <p>Preset returns the unit to its default settings.</p> <p>Preset requires confirmation to prevent accidental reset. Wait until the self tests are completed.</p>
<p>2. Use a BNC Tee to connect the Source Output to the filter input and the Ch1 A Input.</p> <p>Connect the filter output to the Ch2 A Input.</p>	<p>Ch1 measures the source (filter input) and Ch2 measures the filter output.</p>
<p>3. Press [Display Setup]</p> <p>Press &lt;Measure Group&gt;</p> <p>Select (Octave) with the knob and press [Enter].</p> <p>Press [Source], &lt;On&gt;.</p> <p>Press [Link] [Auto Scale A]</p>	<p>Select the Display Setup menu.</p> <p>Change the Measurement Group.</p> <p>Choose the Octave group. Both displays are now making Octave Analysis measurements.</p> <p>Turn the source on.</p> <p>Auto Scale DisplayA and change the scale of DisplayB to match DisplayA.</p> <p>The Octave measurement displays the output of logarithmically spaced bandpass filters. This is not an FFT based measurement. The last bin at the right is a Sound Level measurement and may be calculated independently from the octave bands.</p> <p>Note that even though the source is a single frequency sine wave, the octave display shows a very broad peak. This is because the individual bandpass filters are very broad, 1/3 of an octave in this case.</p>

## 1-24 Octave Analysis

---

<p>4. Press [Source]</p> <p>Press &lt;Noise&gt;</p> <p>Press &lt;Type&gt;</p> <p>Use the knob to select (Pink) and press [Enter].</p> <p>Press [Auto Range Ch1] and [Auto Range Ch2].</p> <p>Press [Auto Scale A] and [Auto Scale B].</p>	<p>Select the Source menu.</p> <p>Choose Noise as the source type. Octave measurements are generally used to measure noise.</p> <p>Change the type of noise.</p> <p>Choose Pink noise. Pink noise rolls off at -3dB per octave. This maintains equal power per octave band and yields a flat octave spectrum.</p> <p>Adjust the input ranges to remove overloads.</p> <p>Notice that the measurement needs to settle after the input range is changed. This is because the measurement is invalid until the input change has propagated through all of the octave band filters. This settling time is related to 1/bandwidth of each filter. Bands which are un-settled are graphed in half intensity. 'Settle' is displayed until all bands in the display are settled.</p> <p>DisplayA (Ch1) shows the flat source spectrum and DisplayB (Ch2) shows the notch filter output.</p>
<p>5. Press [Average]</p> <p>Press &lt;Integration Time&gt;</p> <p>Press [1], select (s) with the knob, and press [Enter].</p>	<p>Select the Average menu. Note that this menu is changed in Octave group.</p> <p>The Integration Time is the averaging time constant. All Octave measurements are rms averaged.</p> <p>Increase the Integration Time to smooth the fluctuations in the spectrum.</p>
<p>6. Press [Freq]</p> <p>Press &lt;Octave Resolution&gt;</p> <p>Use the knob to select (Full) and press [Enter].</p> <p>Press &lt;Octave Resolution&gt;</p>	<p>Select the Frequency menu.</p> <p>Change the number of bands per octave.</p> <p>Choose Full octave bands.</p> <p>Each band represents a full octave with very poor frequency resolution.</p> <p>Change the number of bands per octave again.</p>



<p>Use the knob to select (Twelfth) and press [Enter].</p>	<p>Choose 1/12 octave bands.</p> <p>Each band represents 1/12 of an octave with very good frequency resolution.</p> <p>Note that the measurement requires a long settling time. Narrow bands increases the settling time (1/bandwidth).</p> <p>Choosing narrow bands also increases the number of calculations required and decreases the maximum frequency which can be measured.</p>
<p>7. Press &lt;Octave Resolution&gt;</p> <p>Use the knob to select (Third) and press [Enter].</p> <p>Press &lt;Highest Band&gt;</p> <p>Use the knob to select 20 kHz and press [Enter].</p> <p>Press &lt;Octave Channels&gt;</p> <p>Use the knob to select (1 Channel) and press [Enter].</p> <p>Press &lt;Highest Band&gt;</p> <p>Use the knob to select 40 kHz and press [Enter].</p>	<p>Change the number of bands per octave again.</p> <p>Let's return to 1/3 octaves.</p> <p>Change the highest measured band.</p> <p>Remember, changing the octave resolution has lowered the highest band in the measurement.</p> <p>Set the highest band to 20 kHz.</p> <p>This is the highest allowed band for 2 channel, 1/3 octave analysis.</p> <p>We can change the number of channels which are being measured.</p> <p>Choose single channel octave analysis. This limits the instrument to analyzing only one input, but increases the measurement bandwidth.</p> <p>Change the highest measured band.</p> <p>Single channel analysis has twice the measurement bandwidth of two channel analysis. Thus, to increase the highest measured band, use 1 channel analysis.</p> <p>Set the highest band to 40 kHz.</p> <p>This is the highest allowed band for 1 channel, 1/3 octave analysis.</p>
<p>8. Press [Display Setup]</p>	<p>Select the Display Setup menu.</p>

## 1-26 Octave Analysis

---

Press <Measurement>	Change the Measurement of both displays. Remember, in 1 channel analysis, both displays measure the same input.
Use the knob to select (Oct ch2) and press [Enter].	Both displays now show the filter output on Ch2.
Press <Measurement>	Change the Measurement again.
Use the knob to select (Oct ch1) and press [Enter].	Show the source output on Ch1 on both displays.
9. Press [Input]	Select the Input menu.
Press <Input Config>	Show the Input Configuration submenu.
Press <Ch1 A-Wt Filter>	Choose the Ch1 Input A-Weighting filter.
Use the knob to select (On) and press [Enter].	Turn the A-Wt filter On. The Ch1 Awt indicator at the top of the screen is highlighted.  The hardware A-Wt filter conforms to the ANSI standard and is commonly used in sound measurements. This filter attenuates high and low frequencies according to how people hear and perceive sound.
Press [Link] [Auto Scale A]	The A-Weighted spectrum is a bandpass centered around 2 kHz.
	This concludes this measurement example. You should have a feeling for Octave measurements and how they are setup.

---

# Capture

This example investigates the Capture buffer using FFT measurements. You will use the SR785 to capture a signal and then analyze it from memory.

<p>1. Press [System]</p> <p>Press &lt;Preset&gt;</p> <p>Press [Enter] to confirm Preset.</p>	<p>Display the System menu.</p> <p>Preset returns the unit to its default settings.</p> <p>Preset requires confirmation to prevent accidental reset. Wait until the self tests are completed.</p>
<p>2. Connect the Source Output to the Channel 1 A Input.</p>	<p>Setup to analyze the source output.</p>
<p>3. Press [Source]</p> <p>Press &lt;Sine&gt;</p> <p>Press &lt;Frequency 1&gt;</p> <p>Press [1] [.] [0] [2] [4] select (kHz) with the knob and press [Enter].</p> <p>Press &lt;On&gt;</p>	<p>Select the Source menu.</p> <p>Choose Sine output.</p> <p>Adjust the output Frequency.</p> <p>Enter 1.024 kHz for the Frequency. Enter the value with the numeric entry keys. Select the units with the knob. Enter the new value with the [Enter] key.</p> <p>Turn the source on.</p>
<p>4. Press [Capture]</p> <p>Press &lt;Allocate Memory&gt;</p> <p>Press &lt;Waterfall Memory&gt;</p>	<p>Select the Capture menu.</p> <p>The data memory is allocated between Capture, Waterfall/Order Track, and Arbitrary Waveform storage.</p> <p>Memory is allocated in blocks. Each block can store 2 kPoints. The total number of blocks available is displayed in the &lt;Total Available&gt; menu box. This number is for display only, it cannot be changed from the menu.</p> <p>Larger memory options (up to 4000 blocks) are available.</p> <p>To increase Capture memory, you must first decrease the other allocations so that the sum never exceeds the total available memory.</p>

## 1-28 Capture

---

Press [0] and [Enter].	Decrease the Waterfall allocation to the minimum allowed. Note that the analyzer displays an allocation of 4 blocks. The analyzer always maintains a minimum allocation for each function that uses memory..
Press <Capture Memory>	Select the Capture allocation.
Press [9] [0] [0] and [Enter].	Increase it to 900 blocks. (1,843,200 points).
Press <Confirm Allocation> and <Return>.	You must confirm the new allocation. Changing the memory allocation destroys previously stored data in the memory.
<hr/>	
5. Press <Capture Channels>	Select which inputs to capture.
Use the knob to select (Ch1) and press [Enter].	Choose Ch1 only. In this case, the entire capture buffer is available for Ch1. When both channels are captured, half of the buffer is available for each channel.
Press <Capture Length>	Increase the capture length.
Press [1] [8] [0] [0] and [Enter].	All of the capture allocation (900 blocks) is available. Each block stores 2 kPoints for a total of 1800 kPoints.
<hr/>	
6. Press [Start Capture]	Start the capture. The buffer will take 7.03 seconds to fill. Since the Capture Mode is 1-Shot, the capture stops once the buffer is full. During this time, Capture indicator is highlighted and the Capture Progress indicator shows how much of the buffer has been filled (up to 100%).
	After capture is complete, the Capture indicator shows 'Cap Data' indicating that the Capture buffer contains data.
Press [Active Display]	Make DisplayB (bottom) the active display.
Press [Display Setup]	Select the Display Setup menu.
Press <Measurement>	Change the measurement of DisplayB.
Use the knob to select (Capture1) and press [Enter].	Choose Capture1 to show the contents of the Ch1 Capture buffer.

<p>Press &lt;Zoom&gt;</p> <p>Press [1] [1] and [Enter].</p>	<p>There are far too many points in the buffer to graph each one. The graph shows the envelope of the data in this case.</p> <p>Zoom in to show individual points. 'Expand' below the graph indicates that the graph has been graphically zoomed and does not show all of the data along the X axis.</p> <p>Enter a zoom factor of 11 (<math>2^{11}</math>). The display now shows the signal sine wave clearly.</p>
<p>7. Unplug the signal from the Ch1 input.</p> <p>Press [Input]</p> <p>Press &lt;Input Source&gt;</p> <p>Use the knob to select (Playback) and press [Enter].</p> <p>Press &lt;Playback Config&gt; Press &lt;Playback Length&gt;</p> <p>Press [1] [8] [0] [0] and [Enter].</p> <p>Press &lt;Playback Mode&gt;</p> <p>Use the knob to select (Circular) and press [Enter].</p>	<p>The signal should disappear from the spectrum in DisplayA.</p> <p>Select the Input menu.</p> <p>Change the Input Source to measure from the Capture buffer.</p> <p>Choose Playback instead of the analog inputs.</p> <p>The measurement now takes its input from the data stored in the Ch1 Capture buffer. The signal reappears in the spectrum in DisplayA.</p> <p>Only single channel measurements using Ch1 are allowed in this case (there is no Ch2 data available).</p> <p>Enter the playback configuration submenu. You can choose to playback only a portion of the buffer if desired.</p> <p>Choose the entire buffer by setting the Playback Length equal to the Capture Length.</p> <p>Playback can be 1-Shot (once through the buffer and stop) or Circular (repeat when finished).</p> <p>Choose circular playback. The % indicator at the top of the display shows the current progress through the playback buffer.</p> <p>uffer during playback.</p> <p>DisplayB automatically pans to show the portion of the Capture buffer at the current playback position.</p>

## 1-30 Capture

---

<p>8. Use a BNC TEE to connect the source to both Ch1 and Ch2 A inputs.</p> <p>Press [Input]</p> <p>Press &lt;Input Source&gt;</p> <p>Use the knob to select (Analog) and press [Enter].</p> <p>Press [Display Setup]</p> <p>Press &lt;Measurement&gt;</p> <p>Use the knob to select (FFT ch2) and press [Enter].</p>	<p>Reconnect the analog signal to both Ch1 and Ch2 inputs.</p> <p>Select the Input menu.</p> <p>Change the Input Source.</p> <p>Choose Analog input again. The Capture parameters can not be modified while the measurement input is Playback.</p> <p>Select the Display Setup menu.</p> <p>Change the Measurement for DisplayB.</p> <p>Measure the Ch2 input also.</p>
<p>9. Press [Capture]</p> <p>Press &lt;Capture Channels&gt;</p> <p>Use the knob to select (Ch1+Ch2) and press [Enter].</p> <p>Press &lt;Sampling Rate&gt;</p> <p>Use the knob to decrease the rate to 131.1 kHz and press [Enter].</p> <p>Press [Start Capture]</p>	<p>Select the Capture menu.</p> <p>Change which channels are captured.</p> <p>Choose both channels.</p> <p>The Capture Length is automatically halved to accommodate both channels in the allocated memory.</p> <p>We can increase the capture time by decreasing the Sampling Rate. This decreases the bandwidth of the stored signal.</p> <p>Choose 131.1 kHz as the Sampling Rate. The capture bandwidth is now 51.2 kHz (reduced from 102.4 kHz).</p> <p>During playback from this buffer, the measurement bandwidth will not be allowed to exceed 51.2 kHz.</p> <p>Capture both inputs for 7.03 seconds. Watch for the Capture Progress indicator to reach 100%.</p>
<p>10. Disconnect the signal from both inputs.</p> <p>Press [Input]</p> <p>Press &lt;Input Source&gt;</p>	<p>The signal should disappear from the spectrum in both displays.</p> <p>Select the Input menu.</p> <p>Change the Input Source.</p>

Use the knob to select (Playback) and press [Enter].

Press [Auto Scale A] and [Auto Scale B].

Press [Start/Reset]

Press <Playback Config>  
Press <Playback Speed>

Use the knob to select (Every Time Rec) and press [Enter].

Press [Start/Reset]

Press [Trigger]  
Press <Trigger Source>

Choose Playback from Capture. Since there is captured data for both inputs, both displays start measuring from the capture.

Scale the displays to show the measurements.

Start playback from the beginning of the buffer. Since the default Playback Speed is Normal, the entire playback takes as long as the equivalent real time measurement. In this case, the capture represents 7.73 seconds of data so the playback takes 7.73 seconds as well.

At the current span (51.2 kHz), there are almost 1000 time records of captured data. Not all time records are displayed during Normal playback though all time records contribute to averaged measurements. The display is updated 8 times a second for 7.73 seconds for a total of 62 updates during this playback.

Normal playback is limited to the real time limitations of the equivalent real time analog input measurement.

Enter the Playback Configuration submenu.  
Change the Playback Speed.

Choose Every Time Record. Since the data is stored in memory, we can choose to display every stored time record.

Start the playback at the beginning of the buffer.

In this case, all 1000 time records are measured AND displayed. The display still updates at 8 Hz so playback takes about 125 seconds to complete.

If the time records are overlapped, there may be more than 1000 measurements to display and playback will take even longer.

Every Time Record playback is not limited by real time considerations.

Select the trigger menu.

We'll use manual trigger to show how you can step through each time record in the capture buffer.

## 1-32 Capture

---

<p>Use the knob to select (Manual) and press [Enter] Press &lt;Manual Trigger&gt; several times</p>	<p>Each manual Trigger will step one time record into the capture buffer. Each time record is 1k points long, or represents .11% of the 900k capture buffer. Note that the playback indicator increments by .11% each time &lt;Manual Trigger&gt; is pressed. Return to continuous playback.</p>
<p>Press &lt;Trigger Source&gt;, select (Cont) with the knob and press [Enter].</p>	<p>Return to continuous playback.</p>
<p>11. Press [Freq]</p>	<p>Select the Frequency menu.</p>
<p>Press &lt;Span&gt;</p>	<p>Change the measurement span. The span can not be increased above 51.2 kHz since the captured data is bandwidth limited to 51.2 kHz (because of our capture sampling rate).</p>
<p>Use the knob to select (6.4 kHz) and press [Enter].</p>	<p>Change the span to 6.4 kHz.  Capture playback allows the same captured data to be measured at different spans, windows, averaging, etc. This is useful if the signal is hard to reproduce or occurs infrequently.</p>
<p>Press [Start/Reset]</p>	<p>At this span, the capture buffer only holds 112.48 time records and takes only 14 seconds to playback every time record.</p>
<p>12. Press [Input]</p>	<p>Select the Input menu.</p>
<p>Press &lt;Playback Config&gt; Press &lt;Playback Length&gt;</p>	<p>Enter the playback configuration submenu. It is generally a good idea to playback an exact number of time records. This way, the last record doesn't wrap around and use points from the start of the buffer.  At this span (with no overlap), each time record is 62.5 ms long. The capture sampling rate was 131.1 kHz so 2 kPoints of capture represents 15.625 ms of data. Thus each time record is 8 kPoints of capture long.</p>
<p>Press [8] [9] [6] and [Enter].</p>	<p>We want the Playback Length to be an integer multiple of 8 kPoints. 112 time records uses 896 kPoints and is close to the full capture length.  Change the Playback Length to an exact number of time records.</p>
<p>Press [Start/Reset]</p>	<p>Start the playback again.</p>



---

<p>This concludes this example. Capture and Playback is a way to record a signal and re-analyze it over and over.</p>
---

---

# Waterfall Display

This example demonstrates the use of waterfall displays. Waterfalls are available for FFT and Octave measurements for analog inputs as well as capture playback.

In this example, we will simulate a reverberation measurement measuring the SR785 source. To perform a real measurement, you would use the source to drive a power amplifier and a microphone to receive the signal.

<p>1. Press [System]</p> <p>Press &lt;Preset&gt;</p> <p>Press [Enter] to confirm Preset.</p>	<p>Display the System menu.</p> <p>Preset returns the unit to its default settings.</p> <p>Preset requires confirmation to prevent accidental reset. Wait until the self tests are completed.</p>
<p>2. Connect the Source Output to the Channel 1 A Input.</p>	<p>Setup to analyze the source output.</p>
<p>3. Press [Display Setup]</p> <p>Press &lt;Measure Group&gt;</p> <p>Select (Octave) with the knob and press [Enter].</p>	<p>Select the Display Setup menu.</p> <p>Change the Measurement Group.</p> <p>Choose the Octave group. Both displays are now making Octave Analysis measurements.</p>
<p>4. Press [Freq]</p> <p>Press &lt;Lowest Band&gt;</p> <p>Use the knob to select (50 Hz) and press [Enter].</p>	<p>Select the Frequency menu.</p> <p>Change the lowest band in the display.</p> <p>Select 50 Hz as the lowest band. The lowest band places a limitation on the minimum integration time. Raising this lowest band allows shorter integration times.</p>
<p>5. Press [Source]</p> <p>Press &lt;On&gt;</p> <p>Press &lt;Noise&gt;</p> <p>Press &lt;Type&gt;</p>	<p>Select the Source menu.</p> <p>Turn the source on.</p> <p>Choose Noise as the source type. Octave measurements are generally used to measure noise.</p> <p>Change the type of noise.</p>

---

Use the knob to select (Pink) and press [Enter].	Choose Pink noise. Pink noise rolls off at -3dB per octave. This maintains equal power per octave band and yields a flat octave spectrum.
Press <Burst>	Change the Burst percentage.
Press [9] [6] [Enter]	Using a percentage less than 100% makes the noise source a triggered source. 100% burst outputs noise continuously with 100% duty cycle. Bursts less than 100% will output noise with less than 100% duty cycle and may be triggered.  In this case, the noise will be output for 96% of the source period (100 ms default) or 96 ms.  If the measurement is continuously or source triggered, then the output is noise for 96 ms out of every 100 ms.  If the measurement is externally or manually triggered, the output is 96 ms every trigger with the minimum trigger period equal to 100 ms (the source period).
6. Press [Average]	Select the Average menu. Note that this menu is changed in Octave group.
Press <Averaging Type>	Change the Averaging Type.
Use the knob to select (Linear Time) and press [Enter].	Use Linear Time for best time resolution. Exponential Time averaging takes about 5 integration times to fully respond to a transient. Linear Time averaging responds in a single integration time.
Press <Integration Time>	Change the Integration Time.
Use the knob to select (8 ms) and press [Enter].	Choose the minimum time for the best resolution.
Press <Linear Avg Trig>	Change the Linear Average Trig. This determines how measurements behave when triggering is enabled (not Free Run).

---

## 1-36 Waterfall Display

---

<p>Use the knob to select (Start) and press [Enter].</p>	<p>Choose Start. This means that a new measurement is started as soon as the previous average is complete. In this case, a new measurement is made every 8 ms, regardless of triggering.</p> <p>If we choose One Lin Avg., then a measurement is made only when triggered.</p>
<p>Press &lt;Power Bin&gt;</p>	<p>Change the Power Bin.</p>
<p>Use the knob to select (L) and press [Enter].</p>	<p>Choose L (<math>L_{eq}</math>) as the sound level bin. <math>L_{eq}</math> is a standard broad band sound measurement. When we start the next measurement the result will be displayed as the last bin in the display and labeled 'L'.</p>
<hr/> <p>7. Press [Trigger]</p>	<p>Select the Trigger menu.</p>
<p>Press &lt;Trigger Source&gt;</p>	<p>Change the Trigger Source.</p>
<p>Use the knob to select (Manual) and press [Enter].</p>	<p>Choose Manual trigger. We will start our measurement with a button press. We could also use an external trigger. Triggering on the signal itself requires the use of an external signal source.</p>
<p>Press &lt;Trigd Source Mode&gt;</p>	<p>Change the Triggered Source Mode. This determines whether the source triggers only once (on the first trigger after [Start/Reset]) or on every trigger.</p>
<p>Use the knob to select (One Shot) and press [Enter].</p>	<p>Choose One Shot to trigger the noise burst only once at the start of the measurement.</p>
<hr/> <p>8. Press [Start/Reset]</p>	<p>[Start/Reset] starts the measurement. Since the measurement is triggered (not Free Run), nothing happens until the first trigger is received.</p>
<p>Press &lt;Manual Trigger&gt;</p>	<p>&lt;Manual Trigger&gt; supplies the first trigger. The source outputs a single noise burst (Trigd Source Mode=One Shot). The display starts a continuous stream of octave measurements, each linear averaged for 8 ms and each starting when the previous average is complete (Linear Avg Mode=Continuous).</p>

<p>9. Press [Waterfall]</p> <p>Press &lt;Display&gt;</p> <p>Use the knob to select (Waterfall) and press [Enter].</p> <p>Press &lt;Storage&gt;</p> <p>Press (Save Options) and press [ENTER]</p> <p>Use the knob to select (Active Meas Only) and press [ENTER]</p> <p>Use the knob to select (One Shot) and press [Enter].</p> <p>Press &lt;Total Count&gt;</p> <p>Press [5] [0] [Enter]</p>	<p>Select the Waterfall menu.</p> <p>Change the Display.</p> <p>Choose Waterfall display. This shows measurements scrolling down. Without waterfall storage, this is purely graphical, no data can be read from measurements other than the most recent.</p> <p>Select waterfall Storage.</p> <p>Choose One Shot to fill the waterfall buffer once and stop.</p> <p>Change the number of measurements to store in the waterfall buffer.</p> <p>The Total Count is linked to both displays by default. Entering 50 changes the total count for both displays to 50.</p>
<p>10. Press [Trigger]</p> <p>Press [Start/Reset]</p>	<p>Select the Trigger menu again.</p> <p>[Start/Reset] starts the measurement. Since the measurement is triggered (not Free Run), nothing happens until the first trigger is received.</p>

## 1-38 Waterfall Display

---

Press <Manual Trigger>

<Manual Trigger> supplies the first trigger. The source outputs a single noise burst (Trigd Source Mode=One Shot). The display starts a continuous stream of octave measurements, each linear averaged for 8 ms and each starting when the previous average is complete (Linear Avg Mode=Continuous).

50 averaged measurements are stored in the waterfall buffer starting with the trigger. The first 96 ms (12 measurements) are during the triggered noise burst. The remaining measurements are taken after the noise burst turns off and measure the decay response or reverberation.

The number of records stored in the waterfall is shown in the Vertical Scale Bar as 'wf 50'.

---

11. Press [Alt] and turn the knob clockwise.

Continue until the display does not scroll any further.

[Alt] knob moves the marker in the Z direction (time axis) in the waterfall display. This scrolls the display to show earlier measurement record.

When the keypad and knob are in the alternate mode, the alternate key functions (labeled below each key) are in effect.

The waterfall records are numbered starting with 0 (the most recent measurement) in the back. In this case, we stored 50 records so the earliest record is number 49.

The marker Z position is displayed next to the marker frequency in the Marker Position Bar above the graph. It should read '49' when you have scrolled all the way to the beginning of the buffer.

Notice how the first few records show the spectrum growing at the start of the noise burst. The low frequency bands grow more slowly than the high frequency bands because their filters have longer time constants..

Press [Alt]

Pressing [Alt] removes the keypad and knob from the alternate mode.

Turn the knob clockwise to move the marker along the frequency axis of a single record.

The normal knob function moves the marker along the X axis of a single record.

<p>Move the marker all the way to the right hand edge of the display. The marker position should read 'L:49'.</p>	<p>The last bin in the octave display is the total sound level <math>L_{eq}</math> (as selected by &lt;Power Bin&gt;).</p>
<p>12. Press [Waterfall]</p> <p>Press &lt;More&gt;</p> <p>Press &lt;Slice to Trace&gt;, select (Trace 1) with the knob and press [Enter].</p> <p>Press [Active Display]</p> <p>Press [Alt] [Help/Local], select (Trace 1) with the knob and press [Enter].</p> <p>Press [Auto Scale B]</p> <p>Use the knob to move the marker around in the waterfall slice (DisplayB).</p>	<p>Select the Waterfall menu.</p> <p>Show More of the Waterfall menu.</p> <p>A waterfall slice is the time history of a single X axis location (data at the marker X position from all stored records). The slice data is stored in a trace.</p> <p>Make the bottom display active (DisplayB).</p> <p>Trace to Display is an alternate function. The alternate key functions are labeled below each key (in this case, the [Help/Local] key).</p> <p>Trace to Display recalls trace data to the active display.</p> <p>Auto Scale DisplayB to show the time history of <math>L_{eq}</math>.</p> <p>The points in a slice are numbered and displayed from 0 (oldest) to 49 (newest). Note that this differs from the waterfall display in which the newest record is numbered 0. This is because the slice is a time record with time advancing to the right and it is more natural to number it this way.</p> <p>Note that the first 12 records (96 ms) show a large value for <math>L_{eq}</math> during the noise burst. If this was a real reverberation measurement, the signal would not decay in a single 8 ms measurement but would last for a reverberation time.</p>
<p>13. Press [Active Display]</p> <p>Press [Waterfall]</p> <p>Press &lt;View Count&gt;</p> <p>Press [5] [0] [Enter].</p> <p>Press &lt;Marker Z to&gt;</p>	<p>Make the top display active (DisplayA).</p> <p>Select the Waterfall menu.</p> <p>Change the View Count. This is the number of records which are displayed.</p> <p>Enter 50 to show the entire waterfall buffer.</p> <p>Move the marker to a specific record number.</p>

## 1-40 Waterfall Display

---

Press [0] [Enter].

Enter record 0 (most recent at the back).

Press <Angle>

Change the skew angle of the display.

Use the knob to select (45°) and press [Enter].

Choose 45° to skew the opposite way.

Press <More>

Show More of the Waterfall menu.

Press <Paused Drawing>

Change the waterfall direction.

Use the knob to select (Oldest at Top) and press [Enter].

While the measurement is running, the newest records are added at the top of the waterfall display. When the measurement is done or paused, the waterfall may be drawn with either the newest (Normal) or Oldest at Top.

Since this measurement is measuring a transient decay, it is better to draw the oldest record at the top and the subsequent (smaller amplitude) records in the front.

---

This concludes this example. There are many display parameters in the Waterfall menu which you should familiarize yourself with.

The transient response of any FFT or Octave measurement may be recorded in a waterfall buffer. Using a slice will give a time evolution of a single X axis bin.

---



## Swept Sine Measurement

This example investigates the test filter (enclosed with this manual) using Swept Sine measurements. You will use the SR785 source to provide a sweeping sine source and both inputs to measure the input to and output from the device under test.

<p>1. Press [System]</p> <p>Press &lt;Preset&gt;</p> <p>Press [Enter] to confirm Preset.</p>	<p>Display the System menu.</p> <p>Preset returns the unit to its default settings.</p> <p>Preset requires confirmation to prevent accidental reset. Wait until the self tests are completed.</p>
<p>2. Use a BNC Tee to connect the Source Output to the filter input and the Ch1 A Input.</p> <p>Connect the filter output to the Ch2 A Input.</p> <p>Press [Input], &lt;Input Conifg&gt;</p> <p>Press &lt;Ch1 Input Range&gt;</p> <p>Press [2], select (dBVpk) with the knob, and [Enter].</p>	<p>In this instrument, transfer function is defined as Ch2 response over Ch1 reference. Thus, Ch1 monitors the filter input (source output) and Ch2 measures the response of the device under test.</p> <p>Select the Input Configuration submenu.</p> <p>Adjust the Ch1 input range.</p> <p>Set the range to 2 dBV.</p>
<p>3. Press [Display Setup]</p> <p>Press &lt;Measurements Group&gt;</p> <p>Select (Swept Sine) with the knob and press [Enter].</p>	<p>Select the Display Setup menu.</p> <p>There are six Measurement Groups - FFT, Correlation, Octave, Order Analysis, Swept Sine, and Time/Histogram. The Measurement Group determines which Measurements are available to the displays.</p> <p>Choose the Swept Sine group. The menus now configure swept sine measurements only.</p>
<p>4. Press [Freq]</p> <p>Press &lt;Start&gt;</p> <p>Press [9] [0] [0] select (Hz) with the knob and press [Enter].</p> <p>Press &lt;Stop&gt;</p>	<p>Select the Frequency menu.</p> <p>Adjust the sweep Start Frequency.</p> <p>Enter 900 Hz.</p> <p>Adjust the sweep Stop Frequency.</p>

## 1-42 Swept Sine Measurement

---

Press [1] [.] [1] select (kHz) with the knob and press [Enter].	Enter 1.1 kHz.  The 900 - 1100 Hz sweep covers the filter region of interest.
Press <# of Points>	Adjust the Number Of Points in the sweep.
Press [1] [0] [2] [4] [Enter].	Enter 1024 points.
5. Press [Start/Reset]	Reset and start the sweep. The source is a sine wave whose frequency sweeps from 900 Hz to 1.1 kHz stopping at 1024 discrete frequencies. At each frequency point, the inputs are measured and displayed.
Wait for the sweep to complete at least once.	The small triangular marker moving across the bottom of the graphs indicates the position of the sweep in progress.
Press [Auto Scale A] and [Auto Scale B]	Scale the two displays. The top display is the filter input (source output) and should be fairly flat. The bottom display is the filter output and show the filter notch at 1 kHz.
6. Press [Display Setup]	Select the Display Setup menu.
Press <Measurement>	Change the Measurement of DisplayA (active display).
Select (Freq. Resp.) with the knob and press [Enter].	Choose Frequency Response (filter output divided by filter input).
Press [Display Options]	Select the Display Options menu.
Press <Format>	Change the Display Format.
Select (Single) with the knob and press [Enter].	Choose a Single Display with DisplayA (Frequency Resonse Function) active.
Press [Auto Scale A]	Scale the display to show the Frequency Response.
7. Press [Marker]	Select the Marker menu.
Press <Width>	Change the Marker Region width.
Select (Normal) with the knob and press [Enter].	Choose Normal Width (1/2 division).
Press <Seeks>	Change the Marker Seeks function.

<p>Select (Min) with the knob and press [Enter].</p> <p>Press [Marker Min]</p>	<p>Seek the Minimum data within the Marker Region.</p> <p>Move the Marker to the notch minimum. Read the notch depth and frequency in the marker Position display.</p>
<p>8. Press [Average]</p> <p>Press &lt;Integration Time&gt;</p> <p>Press [4] [0] select (ms) with the knob and press [Enter].</p>	<p>Select the Average menu.</p> <p>Change the Integration Time.</p> <p>Enter 40 ms. The Integration Time is set in increments of 3.9 ms so the entry is rounded to 39 ms.</p> <p>The new estimated sweep time is displayed below the graph.</p> <p>At each frequency point, the inputs measure the amount of signal at the source frequency. This is done by multiplying the input data by the source sine (and cosine) and averaging the results over an integration time. The actual integration time is always rounded up to an exact number of cycles of the source frequency. This rejects signals which are at different frequencies, such as noise and harmonics. Long integration times improve signal to noise while increasing the measurement time.</p> <p>The greater of the Integration Cycles and Integration Time (rounded to the next complete cycle) is used at each frequency. The sweep is at 1 kHz so each cycle is 1 ms. Setting the Integration Time to 40 ms increases the integration time.</p>
<p>9.</p> <p>Press [Auto Range Ch1] and [Auto Range Ch2]</p>	<p>So far, the Input Ranges, Source Level and Sweep Resolution have been constant over the sweep. Let's change these to optimize both the measurement and the measurement time.</p> <p>Change both inputs to Auto Range (the Input Range indicators at the top of the screen are highlighted).</p>

<p>Press [Input], &lt;Input Conifg&gt;</p> <p>Press &lt;Channel&gt;, select (Ch2) with the knob, and press [Enter]. Press &lt;Ch2 Coupling&gt;</p> <p>Select (AC) with the knob and press [Enter].</p>	<p>When Auto Range is on, the Input Range is optimized at each frequency point in the sweep. If the signal is overloaded, the range moves up. If the signal is below half scale, the range moves down. This optimizes the input signal to noise at each point separately and can dramatically improve the S/N of measurements with a large dynamic range. Measurements in excess of 140 dB of dynamic range are possible with swept sine.</p> <p>This notch filter is only -60 dB deep and does not actually require Auto Range for a clean measurement. Note that the Ch2 Input Range changes as the sweep moves through the notch, following the filter output signal.</p> <p>Auto Range increases the sweep time.</p> <p>Select the Input Configuration submenu.</p> <p>Select Ch2 Input Configuration.</p> <p>Change the Input Coupling for Ch2.</p> <p>Choose AC coupling. This eliminates the DC offset from the source and allows Auto Range on Ch2 to cover the entire allowable input range.</p>
<p>10. Press [Freq]</p> <p>Press &lt;Auto Resolution&gt;</p> <p>Select (On) with the knob and press [Enter].</p>	<p>Select the Frequency menu.</p> <p>Change the Auto Resolution mode.</p> <p>Choose Auto Resolution on.</p> <p>If successive points differ by less than the Faster Threshold (on both Ch1 and Ch2), then the sweep starts to skip points. Each successive time this threshold is met, the number of points skipped is increased until the Max # of Skips is reached. This speeds up the sweep in regions where the response is slowly changing.</p>

---

	<p>If successive points differ by more than the Slower Threshold (on either input), then the sweep returns to the earlier point and continues with no skipping. This 'fills' in the region where the response is rapidly changing. The sweep continues from this point, speeding up when allowed and slowing down when required.</p> <p>Note that the sweep progress marker at the bottom of the graph changes speed through the notch.</p> <p>Auto Resolution greatly shortens the measurement time while preserving the resolution where required.</p>
<p>11. Press [Source]</p> <p>Press &lt;Auto Level Ref&gt;</p> <p>Select (Channel 2) with the knob and press [Enter].</p> <p>Press &lt;Ideal Ref&gt;</p> <p>Press [1] [0], select (mV) with the knob, and press [Enter].</p> <p>Press [Active Display]</p>	<p>Select the Source menu.</p> <p>Change the Auto Level Reference.</p> <p>Choose Channel 2 as the Auto Level Reference.</p> <p>Change the Ideal Reference level.</p> <p>Enter 10 mV.</p> <p>Auto Level will try to maintain the Ch2 signal level at the Ideal Reference level by changing the source level at each point of the sweep. The Max Source sets the largest source output allowed.</p> <p>This is useful whenever the transfer function has substantial gain as well as attenuation or if a test requires a constant level within the device under test (usually input or output). In this case, Auto Level is not really required but illustrates its use.</p> <p>Change the active display to DisplayB (which is still measuring the spectrum of Ch2).</p>

---

## 1-46 Swept Sine Measurement

---

Press [Auto Scale B]	Scale the display. On both sides of the notch, you can see Auto Level keeping the Ch2 signal level at 10 mV (-40 dBV). The reference tolerance is 3 dB and is set by the Ref Limits. As the sweep moves into the notch, the source level reaches the Max Source level of 1 V and the Ch2 signal drops to -60 dBV. The spectrum of Ch1 measures the actual source level at each point and the Frequency Response is still calculated correctly.
12.Press [Active Display]	Switch back to DisplayA (Frequency Response).
	This concludes this measurement example. You should have a basic understanding of Swept Sine measurements. The Input Range, Resolution and Source Level optimizations greatly extend the dynamic range of the measurement while minimizing the measurement times.



## 1-48 Saving and Recalling

---

<p>4. Press [Display Setup]</p> <p>Press &lt;Ymid&gt;</p> <p>Press [-] [4] [0] and [Enter].</p> <p>Press [Alt] [Start Capture]</p> <p>Press [Display Ref]</p> <p>Press [Auto Scale A]</p>	<p>Select the Display Setup menu.</p> <p>Change the top reference for the graph.</p> <p>Set the middle reference to -400 dBV.</p> <p>Note that the live data is graphed with the new vertical reference but the reference graph did not move. This allows the reference graph to be offset from the live data (so it is visible).</p> <p>Snap Ref is an alternate function. The alternate key functions are labeled below each key (in this case, the [Start Capture] key).</p> <p>Snap Ref redraws the reference graph with the current display scaling.</p> <p>Pressing [Display Ref] again turns off the reference graph. The reference graph data is lost.</p> <p>To save the reference graph, transfer the reference graph to a trace before turning it off. The trace may be saved to disk.</p> <p>Scale the graph appropriately again.</p>
<p>5. Press [Alt] [Print Screen]</p> <p>Use the knob to select Trace 3 and press [Enter].</p> <p>Press [Span Up] twice to return to full span.</p> <p>Press [Active Display]</p> <p>Press [Alt] [Help/Local]</p>	<p>Display to Trace is an alternate function. This function saves the current data in the active display to a trace.</p> <p>Store the display in Trace 3.</p> <p>There are 5 traces available for data storage. They can store the results of any measurement (other than capture buffer). They can be viewed as complex arrays of data which can be viewed like any other measurement data.</p> <p>Change the live measurement.</p> <p>Make DisplayB (bottom) active.</p> <p>Trace to Display is an alternate function. This function recalls trace data to the active display.</p>



Press [Enter] to select Trace 3.

Only Trace 3 contains data at this time. Only Trace 3 may be selected.

DisplayB now shows the data stored in Trace 3. The display is labeled 'Off-Line' (in its upper left corner) indicating that it is showing static data.

Press [Display Setup]

Select the Display Setup menu.

Note that the Measurement Group and Measurement menu boxes are shown in gray. This indicates that these parameters may not be changed for the active display. In this case, the active display contains stored data so the measurement parameters cannot be changed.

Other measurement parameters, such as averaging and window type, are also not allowed to be changed when the active display is Off-Line.

Press <Ymax>

Change the graph scale of DisplayB.

The view and scale of DisplayB can be changed of course. These parameters simply change the way the stored data is graphed.

Press [8] [0] and [Enter].

Move the graph down.

Press [Active Display]

Make DisplayA (top) the active display. Note that the Measurement Group and Measurement menu boxes are not gray. This is because these parameters may be adjusted for the live measurement in DisplayA.

Press [Active Display]

Make DisplayB (bottom) the active display again.

Press [Display Options]

Select the Display Options menu.

Press <Display>

Make DisplayB live again.

Use the knob to select (Live) and press [Enter].

Choose Live to return the live measurement to DisplayB.

## 1-50 Saving and Recalling

---

<p>6. Put a blank 1.44MB, 3.5" disk into the disk drive.</p> <p>Press [Disk]</p> <p>Press &lt;Disk Upkeep&gt;</p> <p>Press &lt;Format Floppy&gt; and press [Enter] to confirm.</p> <p>Press &lt;Return&gt;</p>	<p>Let's save DisplayA to a disk file.</p> <p>Use a blank disk if possible, otherwise any disk that you don't mind formatting will do. Make sure the write protect tab is off.</p> <p>Select the Disk menu.</p> <p>Choose the Disk Upkeep menu.</p> <p>Make sure that the disk does not contain any information that you want!</p> <p>This function requires a confirmation. Go ahead and confirm. Formatting takes about a minute.</p> <p>Go back to the main Disk menu.</p>
<p>7. Press [Active Display]</p> <p>Press &lt;File Name&gt;</p> <p>Press [Alt]</p> <p>Press [D] [A] [T] [A] [1] [Enter]</p> <p>Press &lt;Display to Disk&gt;</p> <p>Press &lt;Display to Disk&gt;</p>	<p>Make DisplayA (top) the active display again. We are going to save DisplayA to disk.</p> <p>We need a file name.</p> <p>[Alt] lets you enter the letter characters printed below each key. The number and backspace keys function normally.</p> <p>'ALT' is highlighted (at the top of the screen) when the Alternate keys are in use.</p> <p>Enter a file name (any legal DOS file name up to 8 characters). Pressing [Enter] terminates the entry and removes the [Alt].</p> <p>Save the measurement data in DisplayA to disk using the specified file name.</p> <p>The extension .78D is appended automatically.</p> <p>Save the measurement again. Notice that the &lt;File Name&gt; has been changed to DATA2. The SR785 will try to autoincrement any filename containing a number.</p>
<p>8. Press [Active Display]</p> <p>Press &lt;File Name&gt;</p>	<p>Make DisplayB (bottom) the active display.</p> <p>To recall a file, first specify the file name. You can either enter the name or select from the file catalog.</p>

<p>Turn the knob to display the file catalog.</p> <p>Use the knob to select one of the disk files and press [Enter].</p> <p>Press &lt;Disk to Display&gt;</p>	<p>Turning the knob while &lt;File Name&gt; is highlighted displays the file catalog of the current directory.</p> <p>The knob selects a file and scrolls the display. Only the files with the appropriate extension (.78D) are shown. To show all files (*.*), press [Exp].</p> <p>Recall the data in the file to the active display (DisplayB).</p> <p>Once again, DisplayB is 'Off-Line' indicating that it is showing static data.</p>
<p>9. Press [Display Options]</p> <p>Press &lt;Display&gt;</p> <p>Use the knob to select (Live) and press [Enter].</p>	<p>Select the Display Options menu.</p> <p>Make DisplayB live again.</p> <p>Choose Live to return the live measurement to DisplayB.</p>
	<p>This concludes this example. Remember, 'Off-Line' displays are showing stored data, not live measurement results. Many measurement parameters can not be adjusted for an 'Off-Line' display.</p>

## User Math Functions

This example measures the group delay of the test filter (enclosed with this manual) using User Math Functions. You will use the SR785 source to provide a broad band source and both displays to measure the output of the device under test.

<p>1. Press [System]</p> <p>Press &lt;Preset&gt;</p> <p>Press [Enter] to confirm Preset.</p>	<p>Display the System menu.</p> <p>Preset returns the unit to its default settings.</p> <p>Preset requires confirmation to prevent accidental reset. Wait until the self tests are completed.</p>
<p>2. Use a BNC Tee to connect the Source Output to the filter input and the Ch1 A Input.</p> <p>Connect the filter output to the Ch2 A Input.</p>	<p>In this instrument, transfer function is defined as Ch2 response over Ch1 reference. Thus, Ch1 monitors the filter input (source output) and Ch2 measures the response of the device under test.</p>
<p>3. Press [Source]</p> <p>Press &lt;On&gt;</p> <p>Press &lt;Chirp&gt;</p> <p>Press [Window]</p> <p>Press &lt;Window&gt;</p> <p>Select (Uniform) with the knob and press [Enter].</p>	<p>Select the Source menu.</p> <p>Turn on the source.</p> <p>Choose Chirp output. The output is an equal amplitude sine wave at each frequency bin of the FFT spectrum.</p> <p>Select the Window menu.</p> <p>Adjust the FFT Window function.</p> <p>The Chirp source requires the use of the Uniform window since not all chirp frequency components are present at all points in the time record. The chirp is exactly periodic with the FFT time record and does not 'leak' with the uniform window.</p>
<p>4. Press [Auto Range Ch1]</p> <p>Press [Auto Range Ch2]</p>	<p>Let the analyzer automatically set the Input Ranges to agree with the signals. Note that the Input Range readouts at the top of the screen are displayed in inverse when Auto Range is on.</p>
<p>5. Press [Freq]</p> <p>Press &lt;Span&gt;</p>	<p>Select the Frequency menu.</p> <p>Adjust the FFT Span.</p>

<p>Use the knob to adjust the Span to 6.4 kHz and press [Enter].</p>	<p>Set the Span to display the filter notch at 1 kHz.</p> <p>The top display (A) is measuring the filter input and should show a fairly flat spectrum. The bottom display (B) is measuring the filter output and should show a deep notch.</p> <p>Both displays are measuring absolute signal levels.</p>
<p>6. Press [User Math]</p> <p>Press &lt;Function&gt;</p> <p>Use the knob to select (FFTUsrFn2) and press [Enter].</p> <p>Press &lt;Edit Fn&gt;</p> <p>Press &lt;Operands&gt;</p> <p>Use the knob to highlight FFT(2) and press [Enter].</p> <p>Use the knob to highlight / (divide) and press [Enter].</p> <p>Use the knob to highlight FFT(1) and press [Enter].</p>	<p>Select the User Math menu.</p> <p>Choose one of the five user functions available in the FFT measurement group.</p> <p>Select Function2 to edit.</p> <p>Show the Edit Function menu.</p> <p>The display shows the available operands for this function.</p> <p>Use &lt;Operands&gt; and &lt;Operations&gt; to switch the display between operands and operators.</p> <p>Use the knob to highlight the desired box and press [Enter] to insert the selection into the function string at the top of the screen.</p> <p>Choose FFT(2) as the first operand. This is the FFT of Ch2 and is identical to the normal measurement. We are going to define a transfer function (FFT2/FFT1).</p> <p>As soon as the operand is entered, the display switches to show operations.</p> <p>Choose the divide operation next.</p> <p>The display switches back to operands.</p> <p>You can choose another operation instead by pressing &lt;Operations&gt;.</p> <p>Choose FFT(1) as the denominator of the transfer function.</p>

## 1-54 User Math Functions

---

Press <Function String>	This key moves the marker to the function string at the top of the screen. This allows you to delete terms and insert new ones.
Use the knob to move the insertion point to the beginning of the function (highlight FFT2).	Move to the start of the string by highlighting the first term (FFT2).
Press <Operations>	'Ins' above the function string indicates that new terms will be inserted before the highlighted term. 'Rep' indicates that the new term replaces the highlighted one. Use <Insert/Replace> to toggle between insert and replace.
Use the knob to highlight GrpDly and press [Enter].	Insert the Group Delay operator in front of the transfer function.  Closing parentheses are not required (if they are at the very end of the string).
Press <Enter Eq.>	Enter the equation.
<hr/>	
7. Press [Display Setup]	Select the Display Setup menu.
Press <Measurement>	Change the measurement of the active display (DisplayA).
Use the knob to select (FFTUsrFn2) and press [Enter].	Choose FFT User Function2.
Press <View>	Change the view. Group delay is the derivative of the phase with respect to $\omega$ and is a real time.
Use the knob to select (Real Part) and press [Enter].	Choose Real Part. The units automatically switch to linear units.  User Functions are dimensionless quantities. You can choose dB (logarithmic) or units (linear).
Press [Auto Scale A]	Scale the display.  The group delay is the delay time caused by the filter at different frequencies. The data values are seconds of delay. Most points are in the neighborhood of 10 to 600 $\mu$ s.  At the 1 kHz notch, the group delay has a singularity. Remember, the notch filter has a phase discontinuity at the notch frequency.

<p>8. Press [Display Options]</p> <p>Press &lt;d/dx Window&gt;</p> <p>Press [.] [1] [Enter]</p> <p>Press [2] [Enter]</p>	<p>Select the Display Options menu.</p> <p>The group delay is the derivative of the phase. The derivative operation requires an aperture or window.</p> <p>Set the window to 0.1% of the display length. This increases the resolution (makes things narrower) of the graph while decreasing the smoothing of noisy data.</p> <p>Set the window to 2%. This decreases the resolution (makes thing wider) while increasing the smoothing of noisy data.</p>
	<p>This concludes this example. User Functions allow you to define your own measurements starting with the basic SR785 measurements. User Functions can also use stored trace data (for calibrations and normalizations) and user constants.</p>

## Limit Testing

This example is intended to familiarize the user with limit testing. Limit Testing tests the measurement data against a set of defined Limit Segments. When measurement data exceeds a Limit Segment at any point, the test fails. Each display has its own set of Limit Segments.

A Limit Segment is defined as the line between the pair of points (X0,Y0) and (X1,Y1). The segment values between the endpoints are calculated for the displayed span. A segment may be defined as either an Upper or Lower limit. Measurement data which is greater than an Upper limit or less than a Lower limit cause the test to fail.

1. Press [System]  Press <Preset>  Press [Enter] to confirm Preset.	Display the System menu.  Preset returns the unit to its default settings.  Preset requires confirmation to prevent accidental reset. Wait until the self tests are completed.
2. Connect the Source Output to the Channel 1 A Input.  Press [Source], <On> Press [Auto Scale A]	Setup to analyze the source output. The default source is a 10.24 kHz sine.  Turn the source on.  Scale the display to show all of the data.
3. Press [Analysis]  Press <Limit Test>  Press <Edit Limits>	Select the Analysis menu.  Choose Limit Test.  Show the Edit Limits menu. We can't turn on any testing until limits are defined.

---



<p>Press &lt;New Segment&gt;</p>	<p>This function adds a new segment. The new segment has a default position and length.</p>
<p>Press &lt;X0&gt;</p>	<p>The segment is defined by its endpoints, (X0,Y0) and (X1,Y1). These values are specified for the current view and units, in this case, Hz for the x values and dBVpk for the y values.</p>
<p>Press [8] [0] [0] [0] [Enter]</p>	<p>The segments are drawn in half intensity. The arrow markers at the end points point down for upper limits and point up for lower limits. The current segment (whose endpoints are shown and edited in the menu) has two additional arrows at the endpoints.</p>
<p>Press &lt;Y0&gt;</p>	<p>Select X0 first.</p>
<p>Press [-] [1] [0] [Enter]</p>	<p>Enter a value of 8000 Hz.</p>
<p>Press &lt;X1&gt;</p>	<p>Select Y0.</p>
<p>Press [1] [2] [0] [0] [0] [Enter]</p>	<p>Enter a value of -10 dBVpk.</p>
<p>Press &lt;Y1&gt;</p>	<p>Select X1.</p>
<p>Press [-] [1] [0] [Enter]</p>	<p>Enter a value of 12000 Hz.</p>
<p>Press &lt;Return&gt;</p>	<p>Select Y1.</p>
<p>4. Press &lt;Limit Testing&gt;</p>	<p>Enter a value of -10 dBVpk.</p>
<p>Use the knob to select (On) and press [Enter].</p>	<p>The segment should intersect the signal peak at 10.24 kHz.</p>
<p>Press &lt;Edit Limits&gt;</p>	<p>Return to the Limit Test menu.</p>
	<p>Select Limit Testing.</p>
	<p>Turn Limit Testing on. The limit test result is displayed to the left of the graph. In this case 'Fail' should be shown.</p>
	<p>The limit that we drew is an upper limit. Since the data peak exceed this limit, the test fails.</p>
	<p>Go back to the Edit Limits menu.</p>

## 1-58 Limit Testing

---

Press <Shift All>	Select Shift All limit segments. This moves all of the segments together. In this case, there is only one segment.
Press [7] [Enter]	Enter 7 to move the segment up by 7 dBVpk.  The new segment is above the signal peak and the limit test passes.
5. Press <New Segment>	Add another segment.
Press <X0>	Select X0 first.
Press [2] [0] [0] [0] [0] [Enter]	Enter a value of 20000 Hz.
Press <Y0>	Select Y0.
Press [-] [8] [0] [Enter]	Enter a value of -80 dBVpk.
Press <X1>	Select X1.
Press [9] [0] [0] [0] [0] [Enter]	Enter a value of 90000 Hz.
Press <Y1>	Select Y1.
Press [-] [8] [0] [Enter]	Enter a value of -80 dBVpk.  The segment should be above the noise floor. The limit test should still pass.
6. Press <Segment#>	Change the segment which we are editing. The current segment is identified in the order in which they are created.
Press [0] [Enter]	Select segment 0 (the first one). Before editing a segment, make sure that you have chosen the correct one. The current segment is identified in the display by two additional arrows at the endpoints and its endpoint coordinates are displayed in the menu.
Press <Limit Type>	Select the Limit Type for segment 0.
Use the knob to select (Lower) and press [Enter].	Change the limit to a lower limit.  The limit test now fails since data falls below the segment.

---

<p>Press &lt;Return&gt;</p> <p>Press &lt;Limit Beep&gt;</p> <p>Use the knob to select (On) and press [Enter].</p>	<p>Return to the Limit Testing menu.</p> <p>Select Limit Beep.</p> <p>On enables the audible alarm. This alarm alerts you to limit test failures.</p>
<p>7. Press [Display Setup]</p> <p>Press &lt;Measurement&gt;</p> <p>Use the knob to select (Time1) and press [Enter].</p> <p>Press &lt;Measurement&gt;</p> <p>Use the knob to select (FFT ch1) and press [Enter].</p>	<p>Select the Display Setup menu.</p> <p>Change the Measurement.</p> <p>The limit segments are defined for the current measurement, view and units. Changing any of these parameters turns limit testing off.</p> <p>The Limit Testing and Beep are turned off since the limit segments we defined have no meaning for this measurement.</p> <p>Change the Measurement to Time Record Ch1.</p> <p>If we went back to the Limit Testing menu and tried to edit limits now, the previous limit segments would be lost. Each display only has a single set of limits and they are defined for a specific measurement, view and units.</p> <p>Change the Measurement back.</p> <p>Select FFT ch1 again.</p>
<p>8. Press [Analysis]</p> <p>Press &lt;Limit Test&gt;</p> <p>Press &lt;Limit Segments&gt;</p> <p>Use the knob to select (Show) and press [Enter].</p>	<p>Select the Analysis menu.</p> <p>Go to the Limit Testing menu.</p> <p>Select Limit Segments.</p> <p>Show the segments. The two segments defined earlier are still available.</p>
	<p>This concludes this example. Limit testing is a powerful tool for repetitive tests. In an automated test environment, limit segments are usually downloaded from a host computer. The SR785 performs the limit testing in real-time and the results are queried by the host computer.</p>

## Exceedance Statistics

This example is intended to familiarize the user with calculating exceedance centile statistics ( $L_n$ ).  $L_n$  is calculated from measurements stored in the waterfall buffer.

$L_n$  is the amplitude at each bin which is exceeded by  $n\%$  of the records in the waterfall. The SR785 simultaneously calculates  $L_n$  at all frequencies in the spectrum.  $L_n$  is commonly used to characterize environmental noise levels.

<p>1. Press [System]</p> <p>Press &lt;Preset&gt;</p> <p>Press [Enter] to confirm Preset.</p>	<p>Display the System menu.</p> <p>Preset returns the unit to its default settings.</p> <p>Preset requires confirmation to prevent accidental reset. Wait until the self tests are completed.</p>
<p>2. Connect the Source Output to the Channel 1 A Input.</p>	<p>Setup to analyze the source output.</p>
<p>3. Press [Display Setup]</p> <p>Press &lt;Measure Group&gt;</p> <p>Select (Octave) with the knob and press [Enter].</p>	<p>Select the Display Setup menu.</p> <p>Change the Measurement Group.</p> <p>Choose the Octave group. Both displays are now making Octave Analysis measurements.</p>
<p>4. Press [Source]</p> <p>Press &lt;On&gt;</p> <p>Press &lt;Noise&gt;</p> <p>Press &lt;Type&gt;</p> <p>Use the knob to select (Pink) and press [Enter].</p>	<p>Select the Source menu.</p> <p>Turn the Source On</p> <p>Choose Noise as the source type. Octave measurements are generally used to measure noise.</p> <p>Change the type of noise.</p> <p>Choose Pink noise. Pink noise rolls off at -3dB per octave. This maintains equal power per octave band and yields a flat octave spectrum.</p>

---

<p>5. Press [Average]</p> <p>Press &lt;Power Bin&gt;</p> <p>Use the knob to select (L) and press [Enter].</p> <p>Press [Auto Scale A]</p>	<p>Select the Average menu. Note that this menu is changed in Octave group.</p> <p>We will leave the averaging at its default, 100 ms exponential time.</p> <p>Change the Power Bin.</p> <p>Choose L (<math>L_{eq}</math>) as the sound level bin. <math>L_{eq}</math> is a standard broad band sound measurement. The result is displayed as the last bin in the display and is labeled 'L'.</p> <p>Scale DisplayA to show the entire range of the data.</p>
<p>6. Press [Waterfall]</p> <p>Press &lt;Storage&gt;</p> <p>Use the knob to select (One Shot) and press [Enter].</p> <p>Press &lt;Total Count&gt;</p> <p>Press [1] [0] [0] [Enter]</p> <p>Press &lt;Save Option&gt;</p> <p>Use the knob to select (Active Meas Only) and press [Enter]</p>	<p>Select the Waterfall menu. <math>L_n</math> is calculated from measurements stored in the waterfall buffer.</p> <p>Select waterfall Storage.</p> <p>Choose One Shot to fill the waterfall buffer once and stop.</p> <p>Change the number of measurements to store in the waterfall buffer.</p> <p>Enter the Total Count for both displays (100).</p> <p>The Storage Interval is 100 ms. A measurement snapshot is added to the waterfall buffer every Storage Interval so 100 measurements will take 10 s to complete.</p> <p>The SR785 has two waterfall storage options for saving measurements. (All) Means that all possible measurements for the currently selected measurement group will be stored in the waterfall buffer. This option saves more measurements but each trace takes more space in the waterfall buffer, so fewer traces will be available.</p> <p>(Active Meas Only) means that only the currently selected measurement will be saved in the waterfall buffer.</p> <p>To calculate exceedance statistics, we must use the Active Measurement Only save option.</p>

## 1-62 Exceedance Statistics

---

Press [Start/Reset]	<p>Start the measurement. This resets the waterfall buffer. New measurements are added to the waterfall buffer every 100 ms. The number of records stored in the waterfall is shown in the Vertical Scale Bar and increments to 100.</p> <p>Because the waterfall storage is One Shot, the waterfall buffer fills once (100 records). No more records are added after 10 seconds.</p>
7. Press [Analysis]	Select the Analysis menu.
Press <Exceedance Stats>	Select the Exceedance Statistics menu.
Press <Stop Index>	The exceedance is calculated using the records in the waterfall buffer starting with the Start Index (0 is the most recent record) and continuing through the Stop Index (Total Count - 1).
Press [9] [9] [Enter].	Enter 99 to include records 0 through 99 (100 total).
Press [Pause/Cont]	Pause the measurement. Exceedance calculation requires that the active display be paused. This ensures that the waterfall buffer is static and no new records will be added during the calculation.
Press <Calculate Excd>, choose (Trace 1) with the knob and press [Enter].	The exceedance results are stored in a data trace. The trace measurement is the same type as the waterfall measurements. In this case, the trace data for each bin is exceeded by only 1% (Exceedance Pct) of the records stored in the waterfall buffer.
8. Press [Alt] [Link], choose (Trace 1) with the knob and press [Enter].	<p>Trace to Ref is an alternate function. The alternate key functions are labeled below each key (in this case, the [Link] key).</p> <p>Trace to Ref copies the data in a trace to the reference graph of the active display. The reference graph is shown in half intensity.</p> <p>In this case, the reference graph is the level of each octave bin which is exceeded only 1% of the time (<math>L_1</math>).</p> <p>Note that the exceedance is also calculated for the power bin (<math>L_{eq}</math>).</p>

---

<p>Press [Pause/Cont]</p>	<p>Continue the live measurement. Note that the live measurement data rarely exceeds the reference graph.</p>
<p>9. Press [Pause/Cont]</p> <p>Press &lt;Exceedance Pct&gt;</p> <p>Press [9] [9] [Enter].</p> <p>Press &lt;Calculate Excd&gt;, choose (Trace 2) with the knob and press [Enter].</p> <p>Press [Alt] [Link], choose (Trace 2) with the knob and press [Enter].</p>	<p>Pause the live measurement again.</p> <p>Change the Exceedance Percentage.</p> <p>Enter 99%. This level is exceeded 99% of the time.</p> <p>The exceedance results are stored in Trace 2.</p> <p>Trace to Ref is an alternate function. The alternate key functions are labeled below each key (in this case, the [Link] key).</p> <p>Trace to Ref copies the data in a trace to the reference graph of the active display. The reference graph is shown in half intensity.</p> <p>In this case, the reference graph is the level of each octave bin which is exceeded 99% of the time (<math>L_{99}</math>).</p>
	<p>This concludes this example. Exceedance is a common measurement for environmental noise levels, such as airports or highways. Change the octave integration time and waterfall storage interval to optimize the measurement rate. Use a larger waterfall buffer (up to 2000 measurements per display) for long monitoring times.</p>





## Chapter 2

# Analyzer Basics

## In This Chapter

<b>Measurement Groups</b>	<b>2-4</b>	User Function	2-20
Analyzer Configuration	2-4	<b>Views</b>	<b>2-21</b>
<b>What is an FFT?</b>	<b>2-6</b>	Log Magnitude	2-21
Why Look At A Signal's Spectrum?	2-6	Linear Magnitude	2-21
The FFT Analyzer	2-6	Magnitude <sup>2</sup>	2-21
Advantages And Limitations	2-7	Real Part	2-21
<b>FFT Frequency Spans</b>	<b>2-8</b>	Imaginary Part	2-21
Full Span	2-8	Phase	2-21
Spans Less Than Full Span	2-8	Unwrapped Phase	2-22
Baseband Spans	2-9	Nyquist Plot	2-23
Starting the Span Above DC	2-9	Nichols Plot	2-23
<b>FFT Time Record</b>	<b>2-10</b>	<b>FFT Averaging</b>	<b>2-24</b>
Baseband Time Records	2-10	RMS Averaging	2-24
Heterodyned Time Records	2-10	Vector Averaging	2-24
The Time Record Display	2-11	Peak Hold Averaging	2-25
Why Use The Time Record?	2-11	Linear Weighting	2-25
Watch Out For Windowing!	2-11	Waterfall Storage	2-25
<b>FFT Windowing</b>	<b>2-12</b>	Exponential Weighting	2-25
In The Frequency Domain	2-12	<b>Real Time Bandwidth and Overlap</b>	<b>2-27</b>
Uniform	2-12	What is Real Time Bandwidth?	2-27
Hanning	2-13	Averaging Speed	2-27
Flattop	2-13	Overlap Processing	2-27
BMH	2-13	Time Record Increment	2-28
Kaiser	2-14	Settling	2-28
Force-Exponential	2-14	Vector Averaging	2-29
User Defined	2-14	Triggering	2-29
<b>FFT Measurements</b>	<b>2-16</b>	<b>Waterfall Display</b>	<b>2-30</b>
Definitions	2-16	What is a Waterfall?	2-30
FFT Spectrum	2-16	Waterfall Storage	2-30
Power Spectrum	2-16	Waterfall Display	2-31
Time Record	2-17	<b>Capture Buffer</b>	<b>2-33</b>
Windowed Time Record	2-17	Input Sampling	2-33
Orbit	2-18	Capture Fill	2-33
Cross Spectrum	2-18	Capture Playback	2-34
Frequency response	2-18	Capture as the Arbitrary Source	2-35
Coherence	2-19	<b>The Source</b>	<b>2-36</b>
Capture Buffer	2-19	Sine	2-36

## 2-2 Analyzer Basics

---

Two Tone	2-36	Time Record	2-56
Chirp	2-36	Windowed Time Record	2-56
Noise	2-36	RPM Profile	2-57
Arbitrary	2-36	Orbit	2-57
Windowing	2-37	Track	2-57
Source Trigger	2-37	Capture	2-57
External Trigger	2-37	User Function	2-57
<b>Correlation Analysis</b>	<b>2-38</b>	<b>Time/Histogram Measurements</b>	<b>2-58</b>
What is Correlation Analysis?	2-38	Time and Histogram Measurements	2-58
Correlation Measurements	2-38	Measurements in the Time/Histogram	
Capture	2-40	Group	2-58
User Function	2-40	Histogram	2-58
<b>Octave Analysis</b>	<b>2-41</b>	Time Record	2-58
What is Octave Analysis?	2-41	Probability Density Function (PDF)	2-58
Band Center Frequencies	2-41	Cumulative Density Function (CDF)	2-59
Full Octave Bands	2-41	Capture	2-59
1/3 Octave Bands	2-41	User Function	2-59
1/12 Octave Bands	2-41		
Octave Measurements	2-42	<b>Trace Storage</b>	<b>2-60</b>
Octave	2-42	<b>User Math Functions</b>	<b>2-61</b>
Capture	2-42	What is a User Function?	2-61
User Function	2-43	Measurement Groups and Traces	2-61
Octave Averaging	2-43	Operands	2-61
Linear Time	2-43	X Axis	2-62
Exponential Time	2-43	Operations	2-63
Equal Confidence	2-44	User Function Limits	2-64
Peak Hold	2-44		
Sound Level Measurement	2-44	<b>Signal Inputs</b>	<b>2-65</b>
Leq	2-44	Manual Range	2-65
Impulse	2-44	Auto Range	2-65
Peak	2-45	Input Noise	2-65
Total Power	2-45	Input Impedance	2-66
Settling Time	2-46	Anti-aliasing Filter	2-66
		A-Weighting Filter	2-67
		Input Transducer Units	2-67
<b>Swept Sine Measurements</b>	<b>2-47</b>		
Why Use Swept Sine?	2-47	<b>Input Connections</b>	<b>2-68</b>
Swept Sine Measurement Setup	2-49	Single-Ended Connection (A)	2-68
Swept Sine Measurements	2-49	Differential Connection (A-B)	2-68
Spectrum	2-49	Common Mode Signals	2-68
Cross Spectrum	2-50	AC vs DC Coupling	2-69
Frequency response	2-50		
Normalized Variance	2-50	<b>Intrinsic (Random) Noise Sources</b>	<b>2-70</b>
User Function	2-50	Johnson Noise	2-70
Averaging - Settling and Integration	2-50	Shot Noise	2-70
Sweep Frequency and Auto Resolution	2-51	1/f Noise	2-70
Input Auto Ranging	2-52	Total Noise	2-70
Source Auto Level and Ramping	2-53		
<b>Order Analysis</b>	<b>2-55</b>	<b>External Noise Sources</b>	<b>2-71</b>
What is Order Analysis?	2-55	Capacitive Coupling	2-71
Order Spectra and Order Tracking	2-55	Inductive Coupling	2-72
Order Measurement Setup	2-55	Resistive Coupling (Ground Loops)	2-73
Order Measurements	2-56	Microphonics	2-73
Linear Spectrum	2-56	Thermocouple Effects	2-74
Power Spectrum	2-56		

<b>Curve Fitting and Synthesis</b>	<b>2-75</b>
Curve Tables	2-75
Polynomial	2-75
Pole-Zero	2-75
Pole Residue	2-75
Frequency Scale	2-76
Delay	2-76
Trace	2-76
Curve Fitting	2-76
Weighting	2-77
Curve Synthesis	2-77

## Measurement Groups

The SR785 is organized into six Measurement Groups - FFT, Correlation, Octave Analysis, Order Analysis, Time/Histogram, and Swept Sine. Choose the Measurement Group in the [Display Setup] menu.

The SR785 calculates all the measurements in a given group simultaneously, regardless of which measurements are actually being displayed. When a measurement is paused, or finished, you can view the results of **all** measurements in the group without having to take new data.

The Measurement Group determines how the input data is processed. In FFT group, the input data is gathered into time records which are then transformed into spectra. In Correlation group, these spectra are transformed back into the time domain to yield auto and cross correlations. In Time/Histogram group, the time records are processed to give a statistical description of the input signal. In Octave group, the data is passed through a parallel bank of filters and averaged. In Swept Sine group, the input data is processed to determine the spectral content at a single frequency by integrating over an exact number of source cycles. Finally, in the Order Analysis group the input data is processed together with information from a tachometer to give spectra corresponding to the revolutions of a rotating machine, rather than a fixed sampling frequency.

The data processing in each group is governed by the parameters chosen in the menus. For some menus, the choice of the Measurement Group determines which parameters are shown in the menu. The [Frequency], [Average] and [Source] menus have different sets of parameters in each group. For other menus, part or all of the menu is unavailable in certain groups. The [Window], [Capture] and [Waterfall] parameters do not apply to all groups.

Each group has its own set of available measurements. Only the FFT and Order group have time and frequency domain measurements. The Time Histogram and Correlation groups have only time domain measurements, while the Octave and Swept Sine groups have only frequency domain measurements. To select a measurement, choose the Measurement Group, then the Measurement, View and Units.

### Analyzer Configuration

Some measurements made by the SR785 are intrinsically two channel measurements, for instance frequency response, or cross-correlation. Other measurements, such as Autocorrelation or FFT1, only involve a single input channel. The SR785 has a uniquely flexible architecture with regard to processing single and dual channel measurements which is set by the [Input]<Analyzer Coniig> softkey. If two channel measurements are required, this softkey should be set to the default Dual Channel setting. This is the traditional mode in which 2-channel analyzers are operated. All measurements, including two channel measurements are allowed, however, many instrument parameters, such as frequency span, and number of FFT lines must be the same for both channels. Obviously it wouldn't make much sense to take a frequency response measurement where the input and output had different frequency spans.

If two-channel measurements aren't required however, the SR785 allows selection of the Independent Channel configuration. In this mode, each display of the SR785 functions as a completely separate single channel analyzer with an independent set of operating parameters. For instance in the Independent Channel mode display A could be configured to show a broadband measurement of channel 1 with linear averaging while display B could show a narrowband detail of the same input with exponential averaging. The restriction, of course, is that in Independent Channel mode, no two-channel measurements can be computed.

## What is an FFT?

An FFT analyzer takes a time varying input signal, as you would see on an oscilloscope trace, and computes its frequency spectrum. Fourier's basic theorem states that any waveform in the time domain can be represented by the weighted sum of pure sine waves of all frequencies. If the signal in the time domain is periodic, then its spectrum is probably dominated by a single frequency component. The spectrum analyzer represents the time domain signal by its component frequencies.

### Why Look At A Signal's Spectrum?

For one thing, some measurements that are very hard in the time domain are very easy in the frequency domain. Take harmonic distortion. It is hard to quantify the distortion by looking at a good sine wave output from a function generator on an oscilloscope. When the same signal is displayed on a spectrum analyzer, the harmonic frequencies and amplitudes are displayed with amazing clarity. Another example is noise analysis. Looking at an amplifier's output noise on an oscilloscope basically measures just the total noise amplitude. On a spectrum analyzer, the noise as a function of frequency is displayed. It may be that the amplifier has a problem only over certain frequency ranges. In the time domain it would be very hard to tell.

Many of these types of measurements can be done using analog spectrum analyzers. In simple terms, an analog filter is used to isolate frequencies of interest. The filtered signal power is measured to determine the signal strength in certain frequency bands. By tuning the filters and repeating the measurements, a reasonable spectrum can be obtained.

### The FFT Analyzer

An FFT spectrum analyzer works in an entirely different way. The input signal is digitized at a high sampling rate. Nyquist's theorem says that as long as the sampling rate is greater than twice the highest frequency component of the signal, then the sampled data will accurately represent the input signal (in the frequency domain). In the SR785, sampling occurs at 262 kHz. To make sure that Nyquist's theorem is satisfied, the input signal passes through an analog anti-aliasing filter that removes all frequency components above 102.4 kHz. The resulting digital time record is then mathematically transformed into a frequency spectrum using an algorithm known as the Fast Fourier Transform or FFT. The resulting spectrum shows the frequency components of the input signal.

The original digital time record comes from discrete samples taken at the sampling rate. The corresponding FFT yields a spectrum with discrete frequency samples or bins. In fact, the spectrum has half as many frequency bins as there are time points. (Remember Nyquist's theorem.) Suppose that you take 1024 samples at 262 kHz. It takes 3.9 ms to take this time record. The FFT of this record yields 512 frequency points, but over what frequency range? The highest frequency will be determined by the period of 2 time samples or 131 kHz. The lowest frequency is just the period of the entire record or  $1/(3.9 \text{ ms})$  or 256 Hz. The output spectrum thus represents the frequency range from DC to 131 kHz with 512 points spaced every 256 Hz.

**Advantages And Limitations**

The advantage of this technique is its speed. The entire spectrum takes only 3.9 ms to measure. The limitation of this measurement is its resolution. Because the time record is only 3.9 ms long, the frequency resolution is only 256 Hz. Suppose the signal has a frequency component at 380 Hz. The FFT spectrum will detect this signal but place part of it in the 256 Hz bin and part in the 512 Hz bin. One way to measure this signal accurately is to take a time record that is  $1/380$  Hz or 3.846 ms long with 1024 evenly spaced samples. Then the signal would land all in one frequency bin. But this would require changing the sampling rate based upon the signal (which you haven't measured yet). This is not a practical solution. Instead, the way to measure the signal accurately is to lengthen the time record and change the span of the spectrum.

## FFT Frequency Spans

### Full Span

Full span is the widest frequency span corresponding to the fastest available sampling rate. In the SR785, this is DC to 131 kHz using a sampling rate of 262 kHz. Because the signal passes through an anti-aliasing filter at the input, the entire frequency span is not useable. The filter has a flat response from DC to 102.4 kHz and then rolls off steeply from 102.4 kHz to 156 kHz. The range between 102.4 kHz and 131 kHz is therefore not useable and the actual displayed frequency span stops at 102.4 kHz. There is also a frequency bin labeled 0 Hz (or DC). For a time record of 1024 samples (3.9 ms), this bin actually covers the range from 0 Hz to 128 Hz (the lowest measurable frequency) and contains the signal components whose period is longer than the time record (not only DC). So our final displayed spectrum contains 401 frequency bins. The first covers 0 - 128 Hz, the second 128 - 384 Hz, and the 401st covers 102.272 - 102.528 kHz.

### Spans Less Than Full Span

The duration of the time record determines the resolution of the spectrum. What happens if we want a resolution better than 256 Hz? We need to increase the duration of the time record. There are two ways to do this - take more points in each time record or lower the sampling rate. Taking more points is difficult since both the memory and processing requirements increase with the number of points. The longest time record the SR785 can process is 2048 points (800 point FFT).

Instead, we take the approach of lowering the sample rate and making the same number of samples cover a longer time. If we halve the sample rate, this doubles the time record duration and gives us better resolution. However, the sample rate also determines the frequency span. By halving the sample rate, we also halve the frequency span. At a constant number of points in the FFT, we must tradeoff better resolution with narrower frequency spans.

Changing the sample rate of the A-D converter is not practical since that requires changing the analog anti-aliasing filter cutoff frequency. Instead, the incoming data samples (at 262 kHz) are digitally filtered and down-sampled. The advantage is that the digital filter's cutoff frequency can be easily changed. For example, to decrease the sampling rate from 262 kHz to 131 kHz, the incoming data is low-pass filtered to remove any signals above 51.2 kHz. This filter rolls off steeply from 51.2 kHz to 65.6 kHz. Since output of this filter only contains frequencies up to 65.6 kHz, Nyquist only requires a sample rate of 131 kHz and only every other point is kept as part of the time record. The result is a time record of 1024 points sampled at 131 kHz to make up an 7.8 ms record. The FFT processor operates on a constant number of points and the resulting FFT will yield 400 bins from DC to 51.2 kHz. The resolution or linewidth is 128 Hz (1/7.8 ms).

This process of halving the sample rate and span can be repeated by using multiple stages of digital filtering. The SR785 can process a 400 bin spectra with a span of only 195.3 mHz and a time record of 2048 seconds if you have the patience. However, this filtering process only yields baseband measurements (frequency spans which start at DC).



### Baseband Spans

Some points to remember are:

1. The FFT resolution (number of frequency bins in the spectrum) determines the number of points required in the time record.
2. The duration of the time record determines the frequency resolution of the spectrum (spacing of frequency bins in the spectrum).
3. The sampling rate determines the frequency span of the spectrum (Nyquist's theorem). The sampling rate is the number of points in the time record divided by the duration of the time record.

The SR785 allows FFT resolutions of 100, 200, 400 or 800 bins (not counting DC). Changing the resolution does not change the span, instead the time record length is changed.

The various FFT resolutions are summarized below.

FFT Resolution	Frequency Resolution	Time Record
100 bins	Span/100	100/Span
200 bins	Span/200	200/Span
400 bins	Span/400	400/Span
800 bins	Span/800	800/Span

### Starting the Span Above DC

Using digital filtering alone requires that every span start at DC. Frequency shifting is accomplished by heterodyning. Heterodyning is the process of multiplying the incoming signal by a sine (and cosine) wave. The resulting spectrum is shifted by the heterodyne frequency. If we incorporate heterodyning with our digital filtering, we can shift any frequency span so that it starts at DC. The resulting FFT yields a spectrum offset by the heterodyne frequency. When this spectrum is displayed, the frequencies of the X axis are the original frequencies of the signal.

Heterodyning allows the analyzer to compute zoomed spectra (spans that start at frequencies other than DC). The digital processor must filter and heterodyne the input in real time to provide the appropriate filtered and down-sampled time record at all spans and center frequencies.

## FFT Time Record

The FFT operates on time records. A time record is simply a sequence of data samples. The duration of the time record is the FFT resolution/span. There are two types of time records, those corresponding to baseband spans (starting at DC) and those corresponding to zoomed spans (not starting at DC). Zoomed time records are heterodyned (frequency shifted) and do NOT contain the input signal at its original frequencies.

### Baseband Time Records

Baseband time records are very simple to understand. They represent the input signal passed through low-pass filters. At full span, the signal has passed through the analog anti-aliasing filter. The sample rate is 262 kHz. To get the time records for narrower spans, the data is digitally filtered and down-sampled. At a given FFT resolution, each time the span is halved, the bandwidth of the time record is halved and the sampling rate is halved. The length of the time record (in seconds) doubles.

### Heterodyned Time Records

Zoomed time records are more complicated. Heterodyning is a complex operation. The input points are multiplied by  $\cos(\omega t)$  and  $\sin(\omega t)$  to yield a real and an imaginary part.  $\omega$  is  $2\pi$  times the span center frequency. The real and imaginary parts of each point are orthogonal. You can think of the complex time record as two separate records, one real and one imaginary.

The input signal is frequency shifted or heterodyned. This moves signals at the span center to DC and frequencies below span center to negative frequencies. If the span center is at 51.2 kHz, the input range from 0 to 102.4 kHz is shifted to -51.2 kHz to +51.2 kHz. This data is then passed through a low-pass filter which cuts off at  $\pm 51.2$  kHz. This results in a  $\pm 51.2$  kHz (102.4 kHz) useable span centered at 51.2 kHz. The output data only requires a sampling rate of 131 kHz (instead of the original 262 kHz real input rate) so only every other point is saved. Thus, the original 102.4 kHz span is represented by a time record with half as many points and half the sampling rate and the same duration. How can this be?

The complex time record has half as many points as the baseband (real) time record with the same span and resolution. This is because the negative frequency part of the spectrum is kept in the heterodyned case. You can think of the real and imaginary parts of the complex time record as completely independent data streams, each at half of the original sample rate and each with half of the original span. Together, they represent the original span with the original number of samples and the original time record length.

Digital filtering and down-sampling is used to narrow the span of the heterodyned data. This 'zooms' in around the heterodyne frequency (span center). The first digital filter reduces the sample rate by 2 (to 131 kHz) but does not reduce the span.

The second digital filter cuts off at  $\pm 25.6$  kHz and reduces the sample rate by 2 again. The number of points in the time record is NOT halved again (this only happens at the first filter due to the splitting of the real time record into two parts, real and imaginary). The new time record must have twice the original duration and thus, half of the original span. This results in a 51.2 kHz ( $\pm 25.6$  kHz) span centered at 51.2 kHz. The time record

duration is twice the full span time record. The sample rate is 1/4 of the full span baseband sample rate. In comparison, the baseband 0-51.2 kHz span has a sample rate 1/2 of the full span baseband sample rate. This is because the baseband time record is all real and the zoomed time record is complex.

Further filtering and down-sampling reduces the span even further. At each span, the zoomed time record is complex and has half as many points (half the sample rate) as the corresponding real baseband time record.

### **The Time Record Display**

The baseband time record display resembles a digital oscilloscope display. Signals at frequencies above the span have been filtered out. The anti-aliasing filters (both analog and digital) have a steep but finite roll-off at their cutoff frequencies. Signals which are just above the cutoff frequency are outside of the FFT span (not in the displayed spectrum) but appear attenuated in the time record.

Baseband time records are entirely real, they have no imaginary part.

The zoomed time record display does not resemble the original data. The data has been frequency shifted. Signals at the center of the span appear at DC while frequencies at both edges of the span appear as high frequencies. The anti-aliasing filters have a steep but finite roll-off at their cutoff frequencies. Signals which are just outside of the span are not displayed in the FFT but appear frequency shifted and attenuated in the time record.

Zoomed time records are complex, they have both a real and an imaginary part. You can display the magnitude and phase as well as the real or imaginary part. The sampling rate is always half of the equivalent baseband span.

### **Why Use The Time Record?**

The time record display can be useful in determining whether the time record is triggered properly. If the analyzer is triggered and the signal has a large component synchronous with the trigger, then the signal should appear stationary in the time record. If the signal triggers randomly, then the time display will jitter back and forth.

Remember, the time record has a resolution of 1/(sample rate). A triggered time record will always jitter by 1 sample. This jitter is removed in the computation of the phase of the spectrum relative to the trigger.

### **Watch Out For Windowing!**

The SR785 can display both the time record and the windowed time record. Most window functions taper off to zero at the start and end of the time record. If a transient signal occurs at the start of the time record, the corresponding windowed time record and FFT may not show anything because the window function reduces the transient to zero.

# FFT Windowing

A signal which is not exactly periodic within the time record does not fall on an exact frequency bin of the FFT spectrum (integer multiple of the FFT frequency resolution). Its energy is split across multiple adjacent frequency bins. This is true but it's actually worse than that. An FFT spectrum models the time domain as if a time record repeated itself forever. This means the end of the time record is followed by the start of the time record in a circular fashion. If the data is not continuous across the stop to start boundary, the FFT will actually compute the spectrum of the discontinuity and 'leak' energy into all frequencies in the spectrum.

Windows are functions defined over a time record which are periodic in a circular time record. They generally start and end smoothly at zero and are smooth functions in between. When the time record is windowed, its data samples are multiplied by the window function, time point by time point, and the resulting windowed time record is definitely periodic in the circular sense. Windowing eliminates the leakage in the spectrum from signals not exactly periodic with the time record.

### In The Frequency Domain

In the frequency domain, a window acts like a filter. The amplitude of each frequency bin is determined by centering this filter on each bin and measuring how much of the signal falls within the filter.

If the filter is narrow, then only frequencies near the exact bin frequency will contribute to the bin. A narrow filter is called a selective window - it selects a small range of frequencies around each bin. However, since the filter is narrow, frequencies slightly off bin are attenuated and phase shifted. Selective windows are useful for resolving adjacent peaks or improving signal to noise. They should not be used for accurate amplitude measurements (except for signals at exact bin frequencies).

If the filter is wide, then frequencies farther from each exact bin will contribute to the bin amplitude making the signal peaks very wide. However, off bin frequencies are not attenuated. These windows should be used for accurate amplitude measurements rather than good frequency resolution.

Windowing allows the FFT to accurately measure signals at frequencies which are not exact frequency bins. The different types of windows trade off selectivity, amplitude accuracy, and noise floor.

The SR785 offers many types of window functions - Uniform (no windowing), Flattop, Hanning, Blackman-Harris (BMH), Kaiser, Force-Exponential, and User Defined windows.

### Uniform

The Uniform window is actually no window at all. The entire time record is used with equal weighting. A signal will appear in a single frequency bin in the spectrum if its frequency is exactly equal to a frequency bin. (It is exactly periodic within the time record). If its frequency is between bins, it will leak into every bin of the spectrum. These two cases also have a great deal of amplitude variation between them (up to 4 dB).

In general, this window is only useful when looking at transients which do not last the entire time record. The Uniform window may also be used with signals which are exactly periodic in the time record such as a chirp or exact bin sine frequencies.

In the Order Analysis measurement group the uniform window can be used in situations where most of the harmonics of interest will be at integer multiples of the shaft rotation speed. In this case, these harmonics will fall on exact bin frequencies and will be most accurately reproduced with the uniform windows.

### Hanning

The Hanning window is a commonly used window. However, it has an amplitude variation of about 1.5 dB for signals which are not at exact bin frequencies and provides only reasonable selectivity. Its side-lobes are very high and broad for off-bin frequencies. As a result, the Hanning window can limit the performance of the analyzer when looking at signals close together in frequency and very different in amplitude.

The Hanning window is most often used in noise measurements since it has the lowest noise floor.

The Hanning window function is

$$w_i = 1.0 - \cos\left(2\pi \frac{i}{N}\right) \text{ for } i = 0 \dots N-1 \text{ and } N = \text{number of time record points.}$$

### Flattop

The Flattop window has the best amplitude accuracy of any window. The amplitude variation is only about 0.02 dB for signals between exact frequency bins. However, the selectivity is worse. Unlike the other windows, the Flattop window has a very wide pass band and very steep rolloff on either side. Thus, signals appear wide but do not leak across the whole spectrum.

The Flattop window is the best window to use for accurate amplitude measurements.

The Flattop window function is

$$w_i = 1.0 - 1.93 \cdot \cos\left(2\pi \frac{i}{N}\right) + 1.29 \cdot \cos\left(4\pi \frac{i}{N}\right) - 0.388 \cdot \cos\left(6\pi \frac{i}{N}\right) + 0.028 \cdot \cos\left(8\pi \frac{i}{N}\right)$$

for  $i = 0 \dots N-1$  and  $N = \text{number of time record points.}$

### BMH

The BMH window combines good selectivity and reasonable accuracy (about 0.8 dB for signals between exact frequency bins). The BMH window has much lower side-lobes than the Hanning window and very little broadening of non-bin frequencies.

The BMH window is a good window to use for measurements requiring a large dynamic range.

The BMH window function is

## 2-14 FFT Windowing

---

$$w_i = 1.0 - 1.36109 \cdot \cos\left(2\pi \frac{i}{N}\right) + 0.39381 \cdot \cos\left(4\pi \frac{i}{N}\right) - 0.032557 \cdot \cos\left(6\pi \frac{i}{N}\right)$$

for  $i = 0 \dots N-1$  and  $N =$  number of time record points.

### Kaiser

The Kaiser window combines excellent selectivity and reasonable accuracy (about 0.8 dB for signals between exact frequency bins). The Kaiser window has the lowest side-lobes and least broadening for non-bin frequencies. This makes this window the best for selectivity.

The Kaiser window is the best window to use for measurements requiring a large dynamic range.

$$w_i = \frac{I_0 \left[ \pi \cdot \alpha \cdot \sqrt{1 - \left( \frac{i - \frac{N}{2}}{\frac{N}{2}} \right)^2} \right]}{I_0(\pi \cdot \alpha)} \quad \text{for } i = 0 \dots N-1 \text{ and } N = \text{number of time record points.}$$

$\alpha = 0.1R/\pi$ ,  $R = 120.0$  and  $I_0$  is the modified Bessel function of the first kind.

### Force-Exponential

Many impact measurements require a Force window for excitation channel and an Exponential window for the response channel. With the SR785, this is accomplished by selecting the Force/Exponential Window. The Force/Exponential window is actually two windows, either of which can be applied to either input channel.

The Force window is uniform over the beginning of the time record and sets the remainder of the time record to its average value. This method minimizes spectral artifacts caused by zeroing the remainder of the time record. The force length is user specified. This window is used to isolate impulsive signals, such as impact excitations, from noise and other oscillations later in the time record.

The Exponential window attenuates the time record with a decaying exponential time constant. This window is often used in impact testing on the response channel to remove oscillations which last longer than the time record.

When the Force-Exponential window is selected and the analyzer is in the Dual Channel configuration the choice of which window is applied to which input channel is made by the <Channel 1 Window> and <Channel 2 Window> softkeys. When the analyzer is in the Independent Channels configuration, these softkeys govern which window is applied to the measurement in the two displays.

### User Defined

The User window is any function that the user provides. The User window is copied from a stored trace using <Trace to Window>. The trace may contain display data or be loaded from disk or via the computer interface.

Remember, window functions have a great deal of impact on the resulting FFT spectrum. A poorly designed window can result in significant measurement errors.

# FFT Measurements

## Definitions

In the measurement definitions which follow, the following conventions are used. Angle brackets  $\langle \rangle$ , imply that the real and imaginary parts of the quantity within the brackets are averaged with either linear or exponential weighting over some number of time records.  $\text{MAX}()$  implies that the complex quantity inside the parentheses is replaced by the maximum magnitude which that complex value takes over some number of time records. The asterisk (\*) is used to denote the conjugate of a complex quantity.

## FFT Spectrum

The FFT spectrum is the basic measurement of an FFT analyzer. It is simply the FFT of a time record. The FFT spectrum is a complex quantity (it contains phase as well as amplitude information). This is sometimes referred to as the linear spectrum.

The phase of the spectrum is meaningful only if the time record is triggered with a fixed relationship to the input signal. If the signal is repetitive as well (the signal and trigger repeat), then vector averaging can be used to reduce the noise level of the spectrum. The vector averaged spectrum is still a complex quantity.

The precise definition of the FFT1 measurement for all averaging modes is as follows:

### No Average

$$\text{FFT 1} = \text{FFT1}$$

### Vector Average

$$\text{FFT 1} = \langle \text{FFT1} \rangle$$

### RMS Average

$$\text{FFT 1} = \sqrt{\langle \text{FFT1} * \bullet \text{FFT1} \rangle}$$

### Peak Hold Average

$$\text{FFT 1} = \sqrt{\text{MAX}(\text{FFT1} * \bullet \text{FFT1})}$$

## Power Spectrum

The Power Spectrum is derived from the FFT spectrum by multiplying the spectrum by its complex conjugate. The averaged power spectrum is a good approximation to the rms signal and noise amplitudes. The power spectrum is a real quantity and contains no phase information.

The precise definition of the Power Spectrum measurement for all averaging modes is as follows:

### No Average

$$\text{Power Spectrum} = \text{FFT1} * \bullet \text{FFT1}$$

### Vector Average

$$\text{Power Spectrum} = \langle \text{FFT1} \rangle * \bullet \langle \text{FFT1} \rangle$$



**RMS Average**

Power Spectrum =  $\langle \text{FFT1} * \text{FFT1} \rangle$

**Peak Hold Average**

Power Spectrum =  $\text{MAX}(\text{FFT1}) \cdot \text{MAX}(\text{FFT1})$

**Time Record**

A time record is simply a sequence of data samples. The duration of the time record is equal to  $1/\text{FFT}$  resolution.

For baseband spans (spans which start at DC), the time record resembles a digital oscilloscope display. Signals at frequencies above the span have been filtered out. Baseband time records are entirely real, they have no imaginary part.

For zoomed spans (spans which start above DC), the time record display does NOT resemble the original data. The data has been frequency shifted. Signals at the center of the span appear at DC while frequencies at both edges of the span appear as high frequencies. Zoomed time records are complex, they have both a real and an imaginary part. The sampling rate is always half of the equivalent baseband span.

Some points to remember,

1. The time record is not a continuous representation of the input signal. The data is sampled and has a time resolution of  $1/(\text{sample rate})$ . High frequency signals will appear distorted in the time record. However, ALL of the spectral information is preserved by the Nyquist sampling theorem as long as the value of each sample is accurate.
2. Averaging does not affect the time record. Averaging is performed on the FFT spectrum and not on the time data.
3. Amplitude calibration is performed in the frequency domain. Hence, the time record amplitudes are not calibrated.
4. A triggered time record will always jitter by 1 sample. This jitter is removed in the computation of the phase of the spectrum relative to the trigger.

Averaged and calibrated time records can be obtained using a User Math function (using inverse FFT).

**Windowed Time Record**

The FFT operates on windowed time records. The window function is applied to the time record immediately before the FFT. Most window functions taper off to zero at the start and end of the time record. If a transient signal occurs at the start of the time record, the corresponding windowed time record and FFT may not show anything because the window function reduces the transient to zero.

### Orbit

The Orbit measurement is a two channel measurement whose real part is the real part of Time Record Ch1 and whose imaginary part is the real part of Time Record Ch2. Orbit is normally displayed with the Nyquist View (Time2 vertical vs Time1 horizontal).

For baseband spans, the time records are entirely real and the Nyquist view of the Orbit measurement is a Lissajous figure.

### Cross Spectrum

The cross spectrum (sometimes called cross power spectrum) is a two channel measurement defined as

#### No Average

$$\text{Cross Spectrum} = \text{FFT1}^* \cdot \text{FFT2}$$

#### Vector Average

$$\text{Cross Spectrum} = \langle \text{FFT1}^* \rangle \cdot \langle \text{FFT2} \rangle$$

#### RMS Average

$$\text{Cross Spectrum} = \langle \text{FFT1}^* \cdot \text{FFT2} \rangle$$

#### Peak Hold Average

$$\text{Cross Spectrum} = \text{MAX}(\text{FFT2})^* \cdot \sqrt{\langle \text{FFT1}^* \cdot \text{FFT1} \rangle}$$

In these definitions the brackets  $\langle \rangle$  represent averaging the real and imaginary part of the enclosed quantity either linearly or exponentially as described below. The cross spectrum contains both magnitude and phase information. The phase is the relative phase (at each frequency) between the two channels. Vector averaging can be used to eliminate signals which do not have a constant phase relationship between the two channels. In this case, triggering may not be required for vector averaging.

The magnitude is simply the product of the magnitudes of each spectrum. Frequencies where signals are present in both spectra will have large components in the cross spectrum.

### Frequency response

The frequency response (sometimes called transfer function) is a two channel measurement which ratios the spectrum of Ch 2 to the spectrum of Ch 1. Frequency response measures the response of a network or device under test. The reference channel (1) measures the signal at the input to the device and the response channel (2) measures the device output. The result is the complex frequency response of the device.

A broadband source (such as chirp or noise) should be used to measure frequency response.

The definition of frequency response depends on the type of averaging which is displayed.

#### No Average

$$\text{Freq. Response} = \text{FFT2} / \text{FFT1}$$

**Vector Average**

$$\text{Freq. Response} = \langle \text{FFT2} \rangle / \langle \text{FFT1} \rangle$$

**RMS Average**

$$\text{Freq. Response} = (\text{RMSAvg cross spectrum}) / \text{power spectrum 1}$$

$$\text{Freq. Response} = \langle \text{FFT1}^* \cdot \text{FFT2} \rangle / \langle \text{FFT1}^* \cdot \text{FFT1} \rangle$$

**Peak Hold Average**

$$\text{Freq. Response} = \text{MAX}(\text{FFT2}) / \sqrt{\langle \text{FFT1}^* \cdot \text{FFT1} \rangle}$$

Both the RMS averaged and Vector averaged frequency response contain both magnitude and phase information. The phase is the relative phase (at each frequency) between the two channels. The RMS averaged frequency response is computed by taking the ration of the cross spectrum to the input power spectrum, a technique called the tri-spectral average.

**Coherence**

The coherence function is a two channel measurement defined as

$$\text{Coherence} = \text{Mag}^2(\text{RMSAvg CrossSpec}) / (\text{Pwr1} \cdot \text{Pwr2})$$

Averaging is always On and the Averaging Modes are defined by the measurement above. The Type and Number Of Averages are still selected in the [Average] menu. Both channels use the window selected in the [Window] menu.

Coherence is a unitless real quantity and varies from 1.0 (perfect coherence) to 0.0 (no coherence). Coherence measures the percentage of power in the response channel (2) which is caused by (phase coherent with) power in the reference or input channel (1). Ideally, a coherence of 1.0 means that the corresponding frequency response is completely legitimate. All of the response power came from power at the input. If there is noise or other signals generated from within the device under test (which is not related to the input signal), it will result in a coherence of less than 1.0.

**Capture Buffer**

The capture buffer stores sequential time domain input data in memory. See 'Capture' later in this section for more details. The Capture measurement displays the contents of the capture buffer. The FFT and Octave Measurement Groups can use this stored data as input data by choosing Playback as the Input Source in the [Input] menu.

The capture buffer is often very long. To graphically expand a region of the display, use the Pan and Zoom functions in the [Display Setup] menu. The capture buffer display can automatically pan with the capture fill or playback progress through the buffer. During capture fill, if the capture buffer accumulates points faster than they be displayed, some points are not shown. This speeds up the display update so that it keeps up with the real time capture but allows visual aliasing to occur. Once capture is complete, the display is redrawn showing the envelope of all points, eliminating any visual aliasing effects.

To measure from a region of the buffer, set the Playback Start and Length in the [Input] menu. During playback, a marker at the bottom of the graph indicates the current time record position within the buffer.

The capture data is filtered and down-sampled according to the capture Sample Rate. Only baseband data (bandwidth starts at DC) are captured. The capture buffer resembles a digital oscilloscope display. Signals at frequencies above the sample rate/2.56 have been filtered out.

The capture buffer is not a continuous representation of the input signal. The data is sampled and has a time resolution of  $1/(\text{sample rate})$ . High frequency signals will appear distorted in the time record. However, ALL of the spectral information (up to the sampling rate/2.56) is preserved by the Nyquist theorem as long as the value of each sample is accurate.

Amplitude calibration is performed in the frequency domain. Hence, the captured time data amplitudes are not calibrated.

### User Function

User Function displays the results of a user defined math function. User Functions defined within the FFT Measurement Group may include FFT measurements. Use the [User Math] menu to define a math function.

A User Function may not be selected as the measurement if it uses a Trace which is empty or which contains data taken in a different Measurement Group.

See 'User Math' later in this section for more.

---

# Views

Each measurement has an associated view. A view is a way of looking at the complex measurement data.

## Log Magnitude

Log magnitude view graphs the log of the magnitude of the measurement data. The magnitude is  $\sqrt{x^2 + y^2}$  where  $x$  is the real part and  $y$  is the imaginary part. The appearance of a display showing the Log Magnitude view depends on the setting of the <dB Units> softkey. If <dB Units> are on, the graph is shown plotted with a Y axis scaled linearly in dBs. If <dB Units> are off, the Y axis is shown scaled logarithmically in the units of the measurement.

## Linear Magnitude

Linear magnitude view graphs the magnitude of the measurement data. The magnitude is  $\sqrt{x^2 + y^2}$  where  $x$  is the real part and  $y$  is the imaginary part. The Y axis of the display is linear in scaling.

## Magnitude<sup>2</sup>

Magnitude<sup>2</sup> view graphs the magnitude squared of the measurement data. The magnitude squared is  $(x^2 + y^2)$  where  $x$  is the real part and  $y$  is the imaginary part. The Y axis of the display is linear in scaling.

## Real Part

Real Part view graphs the real part of the measurement data. The Y axis of the display is linear in scaling.

## Imaginary Part

Imaginary Part view graphs the imaginary part of the measurement data. The Y axis of the display is linear in scaling.

Real measurement data, such as baseband time record, have zero imaginary part. This view is zero for all points.

Imaginary data arises from the multiplication of the input time data by sine and cosine (heterodyne). Signals in phase with cosine result in real data, signals in phase with sine result in imaginary data. The real and imaginary parts represent data 90 degrees out of phase.

## Phase

Phase view graphs the phase of the measurement data. The phase is a four-quadrant quantity defined as  $\theta = \arctan(y/x)$  where  $y$  is the imaginary part and  $x$  is the real part and  $y$  is positive. If  $y$  is negative, the phase is negated. Phase is graphed linearly from -180 ( $-\pi$ ) to +180 ( $\pi$ ) degrees (radians). To show unwrapped phase, choose the Unwrapped Phase view.

Real measurement data, such as baseband time record, have zero imaginary part. The phase is zero for all points.

Single channel phase is relative to the center of the time record for Uniform, BMH, Hanning, Flattop and Kaiser windows. For Force and Exponential windows, phase is relative to the start of the time record. For User windows, the Window Form is user specified. In general, single channel phase is useful only when the time record is triggered in phase with the signal.

For two channel measurements, the phase is relative between channel 2 and channel 1. Triggering is not always required for meaningful two channel phase measurements.

Phase Suppress sets the phase of small data values to zero. This avoids the messy phase display associated with the noise floor. (Remember, even a small signal has phase.) Set the Phase Suppress threshold in the [Display Options] menu.

### **Watch Out For Triggered Sources**

If a triggered source is selected (Chirp, Burst Chirp, Burst Noise or Arbitrary), the triggered measurement phase is stable only if the input signals are derived from the source output. Turn the source off (or set it to Sine) when making triggered measurements without the source.

### **Watch Out For Phase Errors**

The FFT can be thought of as a set of bandpass filters, each centered on a frequency bin. The signal within each filter shows up as the amplitude of each bin. If a signal's frequency is between exact FFT frequency bins, the filters cause phase errors. Because these filters are very steep and selective, they introduce very large phase shifts for signals not exactly on a frequency bin. Use the SR785 source to generate exact bin frequencies whenever possible.

## **Unwrapped Phase**

Unwrapped Phase view graphs the phase of the measurement data as a continuous function without 'wrapping' around at  $\pm 180$  degrees. This view is generally meaningful only for measurements which have data at every frequency point (such as chirp source or swept sine).

The phase  $\theta$  is calculated for each point as in the 'wrapped' Phase view. The wrapping starts at the left edge of the display. The value of  $\theta \pm (n \times 360 \text{ deg.})$  which is closest to the phase of the previous bin is assigned to each bin ( $n$  is an integer).

Single channel phase is relative to the center of the time record for Uniform, BMH, Hanning, Flattop and Kaiser windows. For Force and Exponential windows, phase is relative to the start of the time record. In general, single channel phase is useful only when the time record is triggered in phase with the signal.

For two channel measurements, the phase is relative between channel 2 and channel 1. Triggering is not always required for meaningful two channel phase measurements.

Phase Suppress sets the phase of small data values to zero. This avoids the messy phase display associated with the noise floor. (Remember, even a small signal has phase.) Set the Phase Suppress threshold in the [Display Options] menu.

### Nyquist Plot

The Nyquist Plot graphs complex measurement data as imaginary part (along the Y axis) versus real part (along the X axis). The graph is scaled linearly. Adjacent frequency (time) points are connected by a line. This view is generally meaningful only for measurements which have data at every frequency point (such as chirp source or swept sine).

The Marker Position Bar shows the real and imaginary parts of each point as well as its frequency (or time). The marker moves sequentially through the frequency (time) points and can be linked to the frequency of the other display.

Real measurement data, such as baseband time record, have zero imaginary part. This view is entirely along the X axis.

### Nichols Plot

The Nichols Plot graphs complex measurement data as log magnitude (along the Y axis) versus unwrapped phase (along the X axis). The Y axis is logarithmic and the X axis is linear. Adjacent frequency (time) points are connected by a line. This view is generally meaningful only for measurements which have data at every frequency point (such as chirp source or swept sine).

The marker readout shows the log magnitude and phase of each point as well as its frequency or time. The marker moves sequentially through the frequency (time) points and can be linked to the frequency of the other display.

Real measurement data, such as baseband time record, have zero phase. This view is entirely along the Y axis.

## FFT Averaging

Averaging successive measurements together improves accuracy and repeatability. For measurements computed from multiple spectra, averaging is performed on measurement results or individual spectra depending upon the measurement definition. Time records are not averaged.

Just as the SR785 computes all measurements in a given measurement group simultaneously, all averaging variants of any measurement are also computed simultaneously. For instance if a frequency response measurement is paused or finished, you can look at the rms averaged, vector averaged, peak hold averaged, or instantaneous (not averaged) frequency response by adjusting the <Display Avg> softkey, there's no need to acquire new data to look at a different averaging type.

You can control whether the SR785 calculates averaged quantities for measurements with the <Compute Avgs> Softkey. If <Compute Avgs> is set to "Yes," the SR785 will compute averages for all measurements. If <Compute Avgs> is "No" the SR785 will not compute any averages and all quantities displayed will be instantaneous values.

### RMS Averaging

RMS averaging involves averaging the result of multiplying a complex quantity by the complex conjugate of another complex quantity. For instance the RMS averaged FFT is defined as:

$$\text{RMSAvg}(\text{FFT1}) = \sqrt{\langle \text{FFT1} * \bullet \text{FFT1} \rangle}$$

The precise definition of what "RMS Averaging" means for each measurement is given in the description of each measurement. Basically, since RMS averaging always involves averaging the "square" of a quantity, RMS averaging reduces fluctuations in the data but does not reduce the actual noise floor (squared values never cancel). With a sufficient number of averages, a very good approximation of the actual noise can be obtained.

Note that the definition given above always yields a real quantity whose phase is zero. This is not true for all RMS averaged quantities computed by the SR785 however. Both the rms averaged frequency response and the rms averaged cross spectrum are complex quantities whose phase is not necessarily zero.

### Vector Averaging

Vector averaging computes the average of the real part (X) and imaginary part (Y) of a measurement according to

$$\text{VecAvg}(X,Y) = (\langle X \rangle, \langle Y \rangle)$$

ie. the Vector average of a complex quantity is the complex quantity formed by the average of its real and imaginary parts independently

Linear averaging computes the equally weighted mean of X and Y over N measurements. Exponential averaging weights new data more than old data and yields a continuous moving average.



Since signed values are combined in the mean, random signals tend to average to zero. This reduces the noise floor since random signals are not phase coherent from measurement to measurement. Only signals with a constant phase have real and imaginary parts which repeat from time record to time record and are preserved. Vector averaging can substantially improve the dynamic range of a measurement as long as the signals of interest have stable phases.

For single channel measurements, vector averaging requires a trigger. The signal of interest MUST be phase synchronous with the trigger to have a stable phase.

For a two channel measurement, the phase is relative between Channel 2 and Channel 1. As long as the signals of interest have stable relative phases, triggering is not required for vector averaging. Triggering is still required to isolate time records which contain the signals of interest.

### Peak Hold Averaging

Peak Hold averaging is similar to rms averaging, in that the rms values of quantities are calculated, but instead of combining the rms values, rather the magnitude of the new data is compared to the magnitude of the averaged data, and if the new magnitude is larger, then the new data becomes the averaged data. This is done on a frequency bin by bin basis. The result is averaged data with the largest magnitudes which occurred over a number of measurements. Peak Hold can compare a fixed number of measurements or run continuously.

### Linear Weighting

Linear weighting combines N (Number Of Averages) measurements with equal weighting in either RMS or Vector averaging.

While Linear averaging is in progress, the number of averages completed is shown in the Horizontal Scale Bar below the graph. When the Number Of Averages has been completed, the measurement stops and 'Done' is displayed to the left of the graph.

### Waterfall Storage

If Waterfall Storage is on, the waterfall buffer only stores the completed linear averages, not each individual measurement. Each time the linear average is done, the result is stored in the waterfall buffer and the average is reset and started over (instead of stopping). Each completed average counts as a single waterfall record.

### Exponential Weighting

Exponential weighting weights new data more than old data. For RMS and Vector averaging, weighting takes place according to

$$\text{Average}_N = (\text{New Data} \cdot 1/N) + (\text{Average}_{N-1}) \cdot (N-1)/N$$

where N is the number of averages.

While Exponential averaging is in progress, the number of averages completed is shown in the Horizontal Scale Bar below the graph. The displayed number stops incrementing at the Number of Averages while the averaging continues.

## 2-26 FFT Averaging

---

Exponential weighting reaches a steady state after approximately  $5N$  measurements. Once in steady state, further changes in the average are detected only if they last for a sufficient number of measurements. Make sure that the number of averages is not so large as to eliminate changes in the data which might be important.

# Real Time Bandwidth and Overlap

## What is Real Time Bandwidth?

The Real Time Bandwidth is the largest frequency span whose corresponding time record exceeds the time it takes to compute the FFT measurement. At this span and narrower spans, it is possible to compute the measurement for every time record without skipping any input data. The spectra are computed in "real time".

At spans larger than the real time bandwidth, the measurement computations take longer than a time record. The analyzer can not keep up and some input data must be ignored. Data is lost between time records while the computations are in progress.

For FFT Resolutions greater than 100 lines, the SR785 can usually compute the measurements in both displays in less than a time record for all spans. This includes the real time digital filtering and heterodyning, the FFT processing, and averaging calculations. The real time bandwidth in this case is 102.4 kHz. Every input sample contributes to a time record and a measurement.

## Averaging Speed

Real time bandwidth has a direct affect on measurement times when Averaging is On. Consider the averaging of full span (102.4 kHz) 400 line FFT's. If the real time bandwidth is only 10 kHz (as in many analyzers), it takes 40 ms to compute the measurement for a 4 ms time record. This means 9 out of every 10 time records are ignored and only 25 full span measurements can be computed in a second. When averaging is on, this usually drops to 10 measurements per second. At this rate it's going to take a couple of minutes to do 500 averages.

The SR785, on the other hand, can make real time measurements at full span (102.4 kHz). This results in 256 measurements per second (on each display!). In fact, this is so fast, that the display can not be updated for each new measurement. The display only updates about 8 times a second. However, when Averaging is On, all of the measurements contribute to the average. The time to complete 500 averages is only a few seconds. (Instead of a few minutes!)

## Overlap Processing

What about narrow spans where the time record is long compared to the processing time? The analyzer computes one measurement per time record and waits and does nothing until the next time record is complete. The update rate is one measurement per time record. With narrow spans, this is quite slow.

With overlap processing, the analyzer does not wait for the next complete time record before computing the next measurement. Instead, it uses data from the previous time record as well as data from the current time record to compute the next measurement. This speeds up the update rate. Remember, most window functions are zero at the start and end of the time record. Thus, the points at the ends of the time record do not contribute much to the FFT. With overlap, these points are "re-used" and appear as middle points in other time records. This is why overlap effectively speeds up averaging and smoothes out window variations.

Typically, time records with 50% overlap provide almost as much noise reduction as non-overlapping time records when RMS averaging is used. When RMS averaging narrow spans, this can reduce the measurement time by 2.

### Time Record Increment

The Time Record Increment is how far the start of each time record is advanced between measurements (as a percentage of the time record length).

If the increment is 100%, the start of the next time record is exactly one time record advanced from the start of the previous time record (no overlap).

If the increment is 25%, then the next time record starts 1/4 of a time record advanced from the start of the previous time record. This is referred to as 75% overlap since the two time records share 75% of a record. The overlap is simply 100% minus the Time Record Increment.

When the Time Record Increment is less than or equal to 100%, the measurement is real time. All time points contribute to one (or more) measurements.

If the increment is 200%, the start of the next time record is advanced from the start of the previous time record by 2 time records leaving a gap of 1 record. This means that the data between the two time records is not measured.

When the Time Record Increment is greater than 100%, then the measurement is not real time and some time points do not contribute to a measurement.

The actual time record increment for the measurements in progress is displayed in the Real Time indicator in the status area (below the Input Ranges). If the indicator shows a value greater than the requested Time Record Increment, it means that the measurement cannot be made with the requested increment but is running with the smallest increment possible. Factors which affect the processor's ability to run real time include the measurement type, averaging and source type.

### Settling

When the frequency span or input signal path (gain, filtering, etc.) is changed, a settling time is required before the FFT measurement is considered settled or valid. The measurement is not settled until the discontinuity in the input data has propagated through the digital filters and a complete new time record has been acquired.

If the time record increment is 100%, unsettled measurements are not displayed. After a change is made which unsettles the measurement, new data is not displayed until the filters are settled and a complete time record has been acquired.

If the measurement is running with a time record increment less than 100% and the measurement is unsettled, unsettled measurements may be displayed. New data is displayed after the filters are settled and a portion of the new time record has been acquired. For example, if the time record increment is 25%, 3 measurements are made before a complete new time record has been acquired. These first 3 measurements have time records which contain data from before AND after the measurement was unsettled. These unsettled measurements are displayed in half intensity indicating that the

measurement is not settled. Once the 4th measurement is complete, the display returns to full intensity since the measurements have completely settled time records.

If averaging is on, changes which unsettle the measurement will restart the average. Unsettled measurements are not included in the new average. Averaging does not start until the measurement is settled. When the time record increment is less than 100%, 'Settle' is displayed (instead of the number of averages) below the graph to indicate that the unsettled measurements are not being averaged.

### **Vector Averaging**

Be sure to select 100% Time Record Increment if you are interested in vector averaged measurement results. Vector averaged results will be incorrect if anything besides 100% Time Record Increment is used because the phase of each time record will not be the same.


### **Triggering**

If the measurement is triggered, then Time Record Increment is ignored. Time records always start with the trigger (with the specified Trigger Delays). The analyzer must use the Continuous Trigger Source to use overlap processing.

# Waterfall Display

## What is a Waterfall?

Waterfall displays show multiple measurement results (records) in a single display. New records are added at the top of the display and older records scroll off the bottom of the display. This gives a time history of the measurement. An example is shown below.

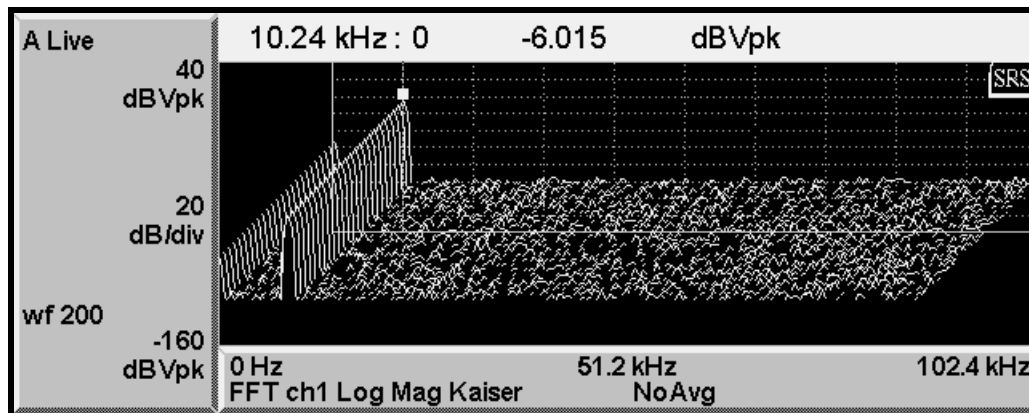

*Figure Chapter 2 -1 Waterfall Display*

Waterfall Storage saves measurement records in waterfall memory. The number of records which may be stored depends upon the allocated memory and the type of measurement. The rate at which records are put into memory is programmable.

Waterfall display WITH storage shows only records which are stored in waterfall memory. While the measurement is running, the display scrolls down and new records are added at the top (back). The display may not scroll fast enough to show every stored record in real time. When the measurement is paused, the display is redrawn with either the newest record in back (normal) or the oldest record in back (reversed). When paused, every stored record which is visible in the display is shown. The marker can access the data in any stored record and scroll the display through the entire waterfall buffer.

Waterfall display WITHOUT storage is simply a graphical display scrolling. The marker is constrained to the most recent record. Older records are only shown graphically, their data has not been saved and may not be accessed in any way. In this case, the display updates and scrolls as quickly as possible.

The [Waterfall] menu configures both waterfall storage and display.

Waterfall displays are not available for Nyquist or Nichols Views or for Swept Sine measurements.

## Waterfall Storage

Waterfall records are stored in waterfall memory. This memory is not retained when the power is off. Use the [Waterfall] <Memory Allocation> menu to allocate memory between the capture buffer, waterfall storage and the arbitrary source waveform. Waterfall memory must be allocated before waterfall storage may be used.

The SR785 can store waterfall measurements in the waterfall buffer in one of two ways depending on the setting of the <Save Option> softkey. If <Save Option> is set to "All," each time the SR785 stores a record it stores enough information to recreate every measurement in the current measurement group. In this mode you can go back after the waterfall storage is complete and look not just at the measurement that was displayed while the waterfall was stored, but at any measurement in the group. Of course, in this mode more waterfall memory is required to store each record.

If <Save Option> is set to "Active Meas Only," only the currently displayed measurement is stored, so you can not go back later and change the measurement. More waterfall records can be stored in this mode however.

The number of measurement records which are stored in waterfall memory is set by the Waterfall Total Count. The maximum Total Count depends upon the allocated memory the type of measurement, and the setting of <Save Option> as described above. The amount of waterfall memory available for each display is always half of the allocated waterfall memory. For example, 500 blocks of allocated waterfall memory is allocated as 250 blocks for each display. For 400 line FFT measurements, This memory can store 250 measurement records from each display in waterfall memory in the "Active Meas. Only" save mode. If <Save Option> is switched to "All" the number of records for each display is reduced to 142.

The rate at which measurement records are added to waterfall memory is set by the Waterfall Skip/Storage Interval. This allows the waterfall memory to hold a longer time history as well as setting a variable storage rate.

In FFT, Order, and Time/Histogram group, for every record added to memory, a skip number of measurements are not stored. For example, a skip of 0 stores every measurement in memory. A skip of 10 stores every 11th measurement in memory (store 1, skip 10). The skipped measurements are still computed and affect (exponential) averaging, they are simply not stored. The elapsed time between stored records ( $1/\text{storage rate}$ ) is simply  $\text{skip}+1$  times the FFT acquisition time times the time record increment. (In the case of triggered measurements,  $\text{skip}+1$  times the trigger period). This allows the time between stored records to be determined exactly.

In Octave group, the Storage Interval is set as a time. In this case, a snapshot is stored to memory every Storage Interval amount of time.

## **Waterfall Display**

To view a waterfall display, set [Waterfall]<Display> to Waterfall. The View Count sets the number of records displayed and the Trace Height sets the percentage of the display height for the Y axis. Scroll Angle and Fast Angles set the angle at which successive records are scrolled in the display. Scrolling at an angle allows changes at a constant frequency to be viewed more easily but takes longer to update the display. Fast Angles limits the choice of Angles to those which scroll faster. The Threshold sets the baseline suppress threshold and Hidden Lines specifies whether portions of a record which are 'hidden' behind other records will be shown.

The Paused Drawing mode determines how the display is drawn when the measurement is paused. Normal displays the newest record at the top (back). Oldest at Top displays

## 2-32 Waterfall Display

---

the oldest record at the top (back). This is reversed from how the display is scrolled while the measurement is running.

A single record can be saved to a trace. A time slice ( history of a single X from all stored records) can also be saved to a trace.



# Capture Buffer

Input Capture samples the analog inputs at a selected sample rate and stores the samples in memory. Storage is continuous and real-time without interruption in the data stream. Triggering only starts the capture, it does not synchronize individual time records.

Captured data can be used as the input for all measurement groups except Swept Sine. The capture buffer can neither be filled nor measured while in the Swept Sine Measurement Group.

Use the [Capture] <Memory Allocation> menu to allocate memory between the capture buffer, waterfall storage and the arbitrary source waveform. Capture memory must be allocated before the capture buffer may be used.

## Input Sampling

The input data is always digitized at the maximum sample rate. The maximum sampling rate is 262.1 kHz when the FFT Base Frequency is 102.4 kHz and 256 kHz when the FFT Base Frequency is 100.0 kHz OR the Measurement Group is Octave. When capturing data in the FFT Measurement Group, it is important to choose the 100.0 kHz Base Frequency if the captured data will be played back in the Octave Group. Only data captured in the Order Measurement Group can be played back in the Order Measurement group because the capture buffer contains special tachometer data in that group that is not present in the other measurement groups.

After the inputs are digitized, they may be filtered and downsampled before being stored in the capture buffer. This allows longer capture times at lower sampling rates. The capture data represents a frequency span from 0 Hz to the sampling rate/2.56. All captured data is baseband. Capturing at less than the maximum sampling rate restricts the playback measurement span to sampling rate/2.56.

In the order measurement group, the sampling rate is determined by the user selected maximum order and maximum rpm.

The Capture Length is selected in 2 kPoint (2048) increments limited by the user memory allocation. A single input or both inputs may be captured. The maximum capture length for a single input is twice the length for both inputs at a given memory allocation.

## Capture Fill

To start capture, press the [Start Capture] key. If the Trigger Mode is Auto Arm and the Trigger Source is Continuous, capture starts immediately. Otherwise, the selected Trigger will be in effect and the next trigger event will start capture. The Capture Progress indicator shows how much of the desired capture length has been filled. If the Capture Mode is 1-Shot, capture stops when the buffer is full. Press [Stop Capture] to halt capture before the buffer is full.

If the Capture Mode is Continuous, once capture is started, it continues indefinitely and fills the capture buffer in a circular fashion. In this case, press [Stop Capture] to halt

capture with the most recently acquired data stored in the buffer. Continuous Capture Mode is not available in the Order measurement group.

While capture is in progress, the displays do not update unless the measurement is a Capture buffer. Use Zoom and Pan to inspect a portion of the buffer. The capture buffer display can automatically pan as the capture progresses through the buffer. During capture, if the capture buffer accumulates points faster than can be displayed, some points are not displayed. This speeds up the display update so that it keeps up with the real time capture data but allows visual aliasing to occur. Once capture is complete, the display is redrawn showing the envelope of all points, eliminating any visual aliasing effects.

### Capture Playback

Captured data can only be played back in a measurement with the same maximum sampling rate. Data digitized at 262.1 kHz can only be played back in FFT Group with the Base Frequency set to 102.4 kHz. Data digitized at 256 kHz can be played back in both Octave and FFT Group as long as the FFT Base Frequency set to 100.0 kHz.

Data that is to be played back in the Order group must be captured in Order group. This is because the order group capture must also capture the tachometer data as well as the signal input.

Choose Playback as the Input Source in any Measurement Group except Swept Sine. Set the Playback Start point within the capture buffer and the Playback Length. When playing back in FFT Group, it is best to set the Playback Length to a multiple of the FFT time record (Acquisition Time).

The Capture Progress indicator shows the playback progress through the buffer (relative to Playback Length).

The capture buffer itself can be displayed as a measurement. Use Zoom and Pan to inspect a portion of the buffer. During playback, the capture buffer display can automatically pan as the playback progresses through the buffer. During playback, a marker at the bottom of the graph indicates the current time record position within the buffer.

The capture data represents a frequency span from 0 Hz to the sampling rate/2.56. The playback measurement cannot exceed this measurement bandwidth. In the Frequency menu, this places a limitation on the maximum FFT span and stop frequency or highest Octave band during playback. The Average, Window and Waterfall menus all operate normally during playback.

The Playback Mode can be either 1-Shot (once through the buffer) or Circular (repeat over and over). In 1-Shot playback, the buffer is played a single time. The playback halts when the end of the buffer is reached. In Circular playback, the playback starts over when the end of the buffer is reached. There is often a discontinuity in the playback measurement when the playback jumps from the end to the start of the buffer. In either mode, press [Start/Reset] to restart playback at the start again.

Normal Playback Speed plays back the capture buffer in real time, i.e. 1 second of capture takes 1 second to play back. Normal Playback Speed looks like the 'live' real time measurement. When the playback is in Octave Group, playback is always Normal Speed.

When playback is in FFT or Time/Histogram Group, the Playback Speed can be either Normal or Every Time Record.

Normal Playback Speed plays back the capture buffer in real time, i.e. 1 second of capture takes 1 second to playback. Normal Speed playback looks like the 'live' real time measurement. Normal playback is limited to the real time limitations of the equivalent real time analog input measurement.

Not all time records are displayed during Normal playback though all time records contribute to averaged measurements. For example, 1 second of capture contains 256 full span FFT time records. Normal Speed playback at full span takes 1 second and updates the display only about 8 times during playback.

Every Time Record playback displays the measurement result for every captured time record. Since the data is stored in memory, Every Time Record playback is not limited by real time considerations. For example, 1 second of capture can contain 256 full span FFT time records. In this case, all 256 time records are measured AND displayed. The display still updates at 8 Hz so playback takes about 32 seconds to complete. If the time records are overlapped, there may be more than 256 measurements to display and playback will take even longer.

When the measurement time records are very long (narrow spans), the time to process and display each measurement is much less than the real time record length. Since the data is already available in the capture buffer, Normal Speed playback means waiting unnecessarily for a real time record to elapse between updates. In this case, Every Time Record playback displays the measurements of all captured time records much faster than Normal Speed real time playback.

When the playback is in Octave Group, playback is always Normal Speed .

### **Capture as the Arbitrary Source**

The contents of the capture buffer can be used as the arbitrary source. This allows a captured waveform to be output as the source. The output sample rate should be the same as the capture sample rate to maintain the correct signal frequencies. If the output sample rate is less than the maximum sampling rate (262 or 256 kHz), then the source output will contain aliases. In this case, be sure to limit the measurement bandwidth to the sampling rate/2.56 to eliminate the effect of these aliases on the measurement.

The amplitude of the arbitrary output depends upon the amplitude of the captured data relative to the Input Range during capture. If the captured data was 100% of the Input Range, then an output amplitude of 100% will be 1 Vpk.

## The Source

The SR785 source provides a variety of test signals which allow the SR785 to measure the response of electronic, mechanical and acoustic devices, without the need for an external generator. In many cases, the SR785 source is better than an external source since it is synchronous with the input sampling.

### Sine

A low distortion sine wave for general purpose gain, distortion and signal/noise measurements. The sine source is synchronous with the FFT, i.e. sine waves can be generated at exact bin frequencies of the FFT. This can eliminate windowing effects in the measured amplitude and phase.

### Two Tone

Two low distortion sine waves can be generated simultaneously for intermodulation distortion tests (IMD). Each tone has independent frequency and amplitude settings.

### Chirp

The Chirp source provides an equal amplitude sine wave at each bin of the displayed spectrum. For a 400 point FFT, the output is the sum of 400 discrete sine waves. The phases of each sine wave are arranged so that they do not add in phase resulting in a low crest factor (ratio of peak to rms). This source is useful for measuring frequency responses quickly without having to make many discrete measurements using a single sine wave.

The Burst Chirp is a sweep over the FFT frequency span in a time less than the time record. Burst Chirp is only available for FFT measurements.

### Noise

Broadband noise is useful for characterizing circuits, mechanical systems or even the audio response of an entire room. White noise provides equal amplitude per root Hz from 0 to 102 kHz, regardless of the measurement span. White noise is useful in electronic applications. Pink noise rolls off at 3 dB/oct providing equal amplitude per octave. Pink noise is preferred in audio applications.

White Noise can be bandlimited to the frequency span of the measurement.

Since the signal is noisy and random, FFT windows are always required when using the Noise source.

Burst Noise is noise output for a fraction of the time record (FFT) or Source Period (Octave).

### Arbitrary

The Arbitrary source plays a waveform stored in memory. The waveform can be either Capture buffer or the Arbitrary Waveform memory. Use the <Memory Allocation> menu to allocate memory between the capture buffer, waterfall storage and the arbitrary source waveform.

The Arbitrary waveform memory can be loaded via the computer interfaces, from a disk file or from a stored trace. The Capture buffer is filled with input samples.

The Arbitrary source can play a portion of memory starting at a specified point. The output sampling rate can also be specified. The full scale amplitude of the source is set as a percentage of 1 Vpk.

### **Windowing**

The Sine, Two Tone, and Chirp sources can be used with or without a window function.

The Sine and Two Tone frequencies can be set at exact bin frequencies of the spectrum, hence they can be exactly periodic in the time record. As long as the signal to noise at the input is high, windowing is not required.

The Chirp waveform consists of many sine waves, each one perfectly periodic in the time record. The Chirp waveform **REQUIRES** a Uniform window to result in a flat spectrum. This is because the individual frequency components do not have a constant amplitude over the time record. Windowing will attenuate certain portions of the spectrum.

Windowing is required when using the Noise source.

### **Source Trigger**

The Sine, Two Tone, Chirp and Arbitrary sources can trigger the FFT time record to measure phase response and/or vector average. Select Source as the Trigger Source in the [Trigger] menu.

For Sine and Two Tone, the source frequencies must be set to a multiple of the linewidth in order for stable time records to be acquired. Random frequencies are not exactly periodic over a time record and do not result in a stable phase even with triggering.

For Burst Chirp and Burst Noise, Source trigger aligns the start of the time record with the start of the burst output.

Source trigger only applies to FFT and Correlation measurements.

### **External Trigger**

For Chirp, Burst Chirp, Burst Noise and Arbitrary, the External trigger can trigger the source output. The source can trigger on only the first trigger or on every trigger after a measurement is started as selected in the [Trigger] menu. The measurements are triggered normally in either case.

# Correlation Analysis

## What is Correlation Analysis?

Correlation is a time domain measurement which reveals similarities and periodicities within a signal (autocorrelation) or between two signals (cross-correlation). Although a time domain measurement, the SR785 computes correlation measurements by transforming input data into the frequency domain, processing it, and retransforming back to the time domain.

## Correlation Measurements

Correlation Analysis is a Measurement Group. Choose the Measurement Group within the [Display Setup] menu. The available measurements are Auto Correlation, Cross Correlation, Time Record, Windowed Time record, Capture Buffer and User Function.

### Auto Correlation

Auto correlation is a single channel measurement. In the time domain, it is a comparison of the signal  $x(t)$  with a time shifted version of itself  $x(t-\tau)$  displayed as a function of  $\tau$ . This is useful for detecting similarities which occur at different times. For example, echoes show up as peaks separated by the echo time. Sine waves appear as sine waves in auto correlation and square waves appear as triangles. Signals which do not repeat or are completely random (such as noise) appear only at  $\tau = 0$ .

The definition of Auto Correlation depends upon the what type of averaging is selected in the [Average] menu.

#### Averaging

Auto Correlation =  $\text{invFFT}(\text{FFTuN} \cdot \text{FFTN}^*)$

#### RMS or Peak Hold Averaging

Auto Correlation =  $\text{invFFT}(\langle \text{FFTuN} \cdot \text{FFTN}^* \rangle)$

#### Vector Averaging On

Auto Correlation =  $\text{invFFT}(\langle \text{FFTuN} \rangle \cdot \langle \text{FFTN}^* \rangle)$

where N is Channel 1 or 2. FFT is a windowed FFT, FFTu is an un-windowed FFT (uniform window) and invFFT is an inverse FFT.

Correlation is a real function and requires a baseband span (real time record). Non-baseband time records do not preserve the original signal frequencies and thus do not yield the correct correlation. A display which is measuring a single channel correlation will have its start frequency set to 0 Hz.

A correlation window is applied to the time record of one FFT in the computation. This is because the FFT models the time domain as a single time record repeating itself over and over. Computing the correlation over a  $\tau$  greater than half of the time record length will result in 'wrap around' where the correlation starts to repeat itself. To avoid this, special windows which zero half of the time record are used. The  $[0..T/2]$  window zeroes the second half of the time record and the  $[-T/4..T/4]$  window zeroes the first and last

quarter of the time record. The  $[-T/2..T/2]$  is a uniform window which should only be used on data which is self windowing (lasts less than half of the time record).

### Cross Correlation

Cross correlation is a two channel measurement. In the time domain, it is a comparison of a signal  $x(t)$  with a time shifted version of another signal  $y(t-\tau)$  displayed as a function of  $\tau$ . This is useful for detecting signals common to both channels but shifted in time.

The definition of Cross Correlation depends upon the displayed average selected in the [Average] menu.

#### No Averaging

Cross Correlation =  $\text{invFFT}(\text{FFT}u2 \cdot \text{FFT}1^*)$

#### Peak Hold or RMS Averaging On

Cross Correlation =  $\text{invFFT}(\langle \text{FFT}u2 \cdot \text{FFT}1^* \rangle)$

#### Vector Averaging On

Cross Correlation =  $\text{invFFT}(\langle \text{FFT}u2 \rangle \cdot \langle \text{FFT}1^* \rangle)$

where FFT1 is the windowed FFT of Channel 1, FFTu2 is the un-windowed FFT (uniform window) of Channel 2 and invFFT is an inverse FFT.

Correlation is a real function and requires a baseband span (real time record). Non-baseband time records do not preserve the original signal frequencies and thus do not yield the correct correlation.

A correlation window is applied to the time record of Ch1 in the computation. This is because the FFT models the time domain as a single time record repeating itself over and over. Computing the correlation over a  $\tau$  greater than half of the time record length will result in 'wrap around' where the correlation starts to repeat itself. To avoid this, special windows which zero half of the time record are used. The  $[0..T/2]$  window zeroes the second half of the time record and the  $[-T/4..T/4]$  window zeroes the first and last quarter of the time record. The  $[-T/2..T/2]$  is a uniform window which should only be used on data which is self windowing (lasts less than half of the time record).

### Time Record

#### Windowed Time Record

These correlation group measurements are similar to their FFT group counterparts. The major differences are:

1. Because correlation is a baseband measurement, the time records are never heterodyned, and are therefore always real.
2. Correlation uses a different set of windows than the FFT measurement group. The Correlation windows, which are described above, are designed to zero half of the time record in order to eliminate "wrap-around" error.

### **Capture**

The capture buffer may be used as a source of data for correlation measurements. See Capture Buffer for more details. The Capture measurement displays the contents of the capture buffer. Correlation measurements can use this stored data as input data by choosing Playback as the Input Source in the [Input] menu.

### **User Function**

User Function displays the results of a user defined math function. User Functions defined within the Correlation Analysis Measurement Group may include correlation measurement results. Use the [User Math] menu to define a math function.



# Octave Analysis

## What is Octave Analysis?

An FFT measurement measures the power within equally divided frequency bins. Octave measurements use logarithmic frequency bands whose widths are proportional to their center frequencies. The bands are arranged in octaves with either 1, 3 or 12 bands per octave (1/1, 1/3 or 1/12 octave analysis). Octave analysis measures spectral power closer to the way people perceive sound, that is, in octaves. 1/3 octave analysis with A-weighting is the most common measurement for acoustics and sound.

To make an octave measurement, the input data passes into a bank of parallel digital filters. The filter center frequencies and shapes are determined by the type of octave analysis, either full (1/1), 1/3 or 1/12 octave, and comply with ANSI S1.11-1986, Order3, Type 1-D. The output from each filter is rms averaged to compute the power and displayed as a bar type graph. This is a real-time measurement of the power within each band and is the only available octave measurement. Since the bands are spaced logarithmically, octave graphs always have a logarithmic X-axis.

In addition to the octave bands, the sound level is also measured and displayed as the last band in the display. Exponential averaged sound power is calculated according to ANSI S1.4-1983, Type 0. Broadband Impulse and Peak measurements are made according to IEC 651-1979, Type 0.

## Band Center Frequencies

Choose the number of bands per octave and the start and stop bands for the measurement within the [Frequency] menu.

The exact band center frequencies are calculated according to the ANSI standard. The displayed frequencies are sometimes rounded to even values for a simpler display. The filter shapes are third-order Butterworth with full, 1/3 or 1/12 octave bandwidth.

### Full Octave Bands

Full octave bands are defined ranging from 125 mHz to 32 kHz. Up to 11 full octaves may be measured at one time. The exact centers are given by

$$\text{Center Frequency} = 1000 \times 2^n \quad \text{for } n = -13 \dots 5$$

### 1/3 Octave Bands

1/3 octave bands are defined ranging from 100 mHz to 40 kHz. Up to 11 octaves (33 bands) may be measured at one time. The exact centers are given by

$$\text{Center Frequency} = 1000 \times 2^{\left(\frac{n-30}{3}\right)} \quad \text{for } n = -10 \dots 46$$

### 1/12 Octave Bands

1/12 octave bands are defined ranging from 100 mHz to 12.34 kHz. Up to 11 octaves (132 bands) may be measured at one time. The exact centers are given by

$$\text{Center Frequency} = 1000 \times 2^{\left(\frac{1}{24}\right)} \times 2^{\left(\frac{n}{12}\right)} \quad \text{for } n = -160 \dots 43$$

### Octave Measurements

Octave Analysis is a Measurement Group. Choose the Measurement Group within the [Display Setup] menu. The available measurements are rms averaged Octave spectrum (power spectrum), Capture Buffer and User Function.

#### Octave

This is the basic measurement of octave analysis. It is simply the rms averaged outputs from the parallel bank of octave filters (also called the power spectrum). The spectrum gives a stable reading of the rms signal amplitudes and noise levels within each band. RMS averaging results in a real spectrum and there is no phase information.

The number of Octave Channels limits the measurement choices of both displays.

If 2 Octave Channels are selected, then both inputs may be measured. The Measurements of both displays are independent. The Highest Band for all Octave Resolutions is decreased by 2 in this case.

If 1 Octave Channel is selected, then only one input may be analyzed. If both displays are making an octave measurement, they must both use the same input and frequency range. Changing the input or frequency range of the active display will change the other display if necessary. Choosing 1 Octave Channel may change the current measurements so that both displays use the same input.

#### Capture

The capture buffer stores sequential time domain data in memory. See 'Capture Buffer' for more details. The Capture measurement displays the contents of the capture buffer. Octave measurements can use this stored data as input data by choosing Playback as the Input Source in the [Input] menu.

The capture buffer is often very long. To graphically expand a region of the display, use the Pan and Zoom functions in the [Display Setup] menu. The capture buffer display can automatically pan with the capture fill or playback progress through the buffer. During capture fill, if the capture accumulates points faster than they can be displayed, some points are not shown. This speeds up the display update so that it keeps up with the real time capture but allows visual aliasing to occur. Once capture is complete, the display is redrawn showing the envelope of all points, eliminating any visual aliasing effects.

To measure from a region of the buffer, set the Playback Start and Stop in the [Input] menu.

The capture data is filtered and down-sampled according to the capture sample rate. Only baseband data (data bandwidth starts at DC) are captured. The capture buffer resembles a digital oscilloscope display. Signals at frequencies above the sample rate/2.56 have been filtered out.

The capture buffer is not a continuous representation of the input signal. The data is sampled and has a time resolution of 1/(sample rate). High frequency signals will appear

distorted in the time record. However, ALL of the spectral information is preserved by the Nyquist sampling theorem as long as the value of each sample is accurate.

Amplitude calibration is performed in the frequency domain. Hence, the captured time data amplitudes are not calibrated.

#### **User Function**

User Function displays the results of a user defined math function. User Functions defined within the Octave Analysis Measurement Group may include octave measurement results. Use the [User Math] menu to define a math function.

A User Function may not be selected as the measurement if it uses a Trace which is empty or which contains data taken in a different Measurement Group.

The number of Octave Channels limits the measurement choices of both displays. If 1 Octave Channel is selected, then only one input may be analyzed. If both displays are making an octave measurement, they must both use the same input. Changing the input of the active display will change the measurement or input of the other display if necessary. User Functions which use both inputs may not be measured. Choosing 1 Octave Channel may change the current measurements so that both displays use the same input.

See 'User Math' later in this section for more.

### **Octave Averaging**

All octave measurements are averaged. There are four types of averaging - Linear Time, Exponential Time, Equal Confidence and Peak Hold.

#### **Linear Time**

The band filter outputs are equally weighted and averaged for an Integration Time.

While Linear averaging is in progress, the integration time completed is shown in the Horizontal Scale Bar below the graph. When the Integration Time has been completed, the measurement stops and 'Done' is displayed below the graph.

If Waterfall Storage is On, the waterfall buffer only stores the completed linear averages, not each individual measurement. Each time the linear average is done, the result is stored in the waterfall buffer and the average is reset and started over (instead of stopping). Each completed average counts as a single waterfall record.

#### **Exponential Time**

New filtered data is weighted more than older data. The exponential 'time constant' is the Integration Time. Averaging continues indefinitely.

While Exponential averaging is in progress, the completed integration time is shown in the Horizontal Scale Bar below the graph. The displayed time stops incrementing at the Integration Time while the averaging continues.

Exponential weighting reaches a steady state after approximately an integration time. Once in steady state, further changes in the average are detected only if they last for a

sufficient number of measurements. Make sure that the integration time is not so large as to eliminate changes in the data which might be important.

### **Equal Confidence**

Equal confidence averaging is exponential averaging where the integration time is set for each band separately. The integration times are set so that there is a 68% probability that the results are within the specified confidence level of the true mean for every band in the measurement. There is a 96% probability that the results are within twice the confidence level (in dB) of the true mean.

The effect of equal confidence averaging is that higher frequency bands have shorter time constants and lower frequency bands have longer time constants. This is very noticeable in the update of the display.

### **Peak Hold**

Peak Hold displays the peak output from each band filter. The peak detection continues indefinitely.

## **Sound Level Measurement**

The last bin in an octave display shows one of many types of sound level measurement.

### **Leq**

Leq exponential time averaged power is computed according to ANSI S1.4-1983. It is computed from real-time low-pass filtered input data, not from the octave band outputs. The bandwidth of Leq is DC-100kHz for 1 octave channel and DC-50kHz for 2 channels.

To measure Leq, set the Averaging Type to Exponential and the Power Bin to L.

When the Averaging Type is Exponential or Linear, the exponential time constant for Leq is the Integration Time. Leq power is always an exponential average.

Standard measurements use either 125 ms (Fast) or 1000 ms (Slow). Note that the lowest octave band places a limit on the averaging time. The standard time constants may not be available if the octave measurement extends to low frequency bands. If the Averaging Type is Equal Confidence, Leq has the same confidence level as the octave measurement. If the Averaging Type is Peak Hold, Leq is not available.

The last bin in the measurement displays the Leq sound level. The last bin is labeled either 'L' (for uniform weighting) or 'LA' (if the input A-Weighting filter is on).

In a User Math function, the L or LA bin is treated the same as the other octave bins, with the exception of the A, B and C weighting operators. These operators do not change the L or LA bin.

### **Impulse**

Broadband Impulse sound level (I) is computed according to IEC 651-1979, Type 0. It is computed from real-time low-pass filtered input data, not from the octave band outputs. The Impulse bandwidth is DC-100kHz for 1 octave channel and DC-50kHz for 2 channels.

To measure I, set the Averaging Type to Exponential, Linear or Equal Confidence and the Power Bin to Impulse.

The time constants for the Impulse measurement are defined by the IEC standard. The Integration Time and Confidence Level are ignored by the Impulse calculation. Impulse power is always an exponential average.

The last bin in the measurement displays the Impulse sound level. The last bin is labeled 'I'.

In a User Math function, the I bin is treated the same as the other octave bins, with the exception of the A, B and C weighting operators. These operators do not change the I bin.

**Peak**

Broadband Peak sound level (P) is computed according to IEC 651-1979, Type 0. It is computed from real-time low-pass filtered input data, not from the octave band outputs. The Peak bandwidth is DC-100kHz for 1 octave channel and DC-50kHz for 2 channels.

To measure P, set the Averaging Type Peak Hold. In this case, the Power Bin must be set to Peak.

The time constants for the Peak measurement are defined by the IEC standard. The Integration Time and Confidence Level are ignored by the Peak calculation. Peak power is always an exponential average.

The last bin in the measurement displays the Peak sound level. The last bin is labeled 'P'.

In a User Math function, the P bin is treated the same as the other octave bins, with the exception of the A, B and C weighting operators. These operators do not change the P bin.

**Total Power**

Total power (T) is computed from the sum of the measured octave bands. The bandwidth is set by the Highest and Lowest Bands of the octave measurement.

To measure Total Power, set the Power Bin to T. Total Power is available for all Averaging Types.

The last bin in the measurement displays the Total Power. The last bin is labeled as follows.

<u>Label</u>	<u>Frequency Weighting</u>
T	uniform weighting
TA	A-Wt (input filter)
Ta	A-Wt (User Math)
Tb	B-Wt (User Math)
Tc	C-Wt (User Math)
Tu	undefined

'Tu' is shown whenever more than one type of weighting exists in the measurement.

In a User Math function, the octave bins are operated on and then summed to provide the T bin. The T bin is computed after the function is evaluated.

### Settling Time

When an octave measurement or input is changed, a settling time is required for the change to propagate through the octave filters and the measurement becomes valid again. This settling time is  $1/(5 \times \text{bandwidth})$  for each band in the measurement.

The bandwidth of an octave filter is given by

$$\text{Bandwidth} = \text{Center Frequency} \times \left[ 2^{\frac{1}{2n}} - \frac{1}{2^{\frac{1}{2n}}} \right] \quad \text{where } n = 1, 3 \text{ or } 12 \text{ (octave)}$$

'Settle' is shown below the graph while settling takes place. Bands which are un-settled are displayed at half intensity in the graph. When Waterfall Storage is selected, records are not added to the waterfall buffer while octavebands are still settling.

# Swept Sine Measurements

A swept sine measurement is basically a sine sweep which steps through a specified sequence of frequency points. At each point, the source maintains a constant frequency and the inputs measure signals only at this frequency. After each measurement point is complete, the source then moves to the next point in the sequence. Unlike the FFT which measures many frequencies at the same time, a swept sine source measures the frequency response a single frequency at a time.


## Why Use Swept Sine?

Frequency response can be measured using the FFT. However, if the frequency response has a large variation within the measurement span, then the FFT may not be the best measurement technique. It's limitation comes from the nature of the chirp source that must be used. The FFT measures the response at all frequencies within the span simultaneously, thus the source must contain energy at all of the measured frequencies. In the time record, the frequency components in the source add up and the peak source amplitude within the time record generally exceeds the amplitude of each frequency component by about 30 dB. Since the input range must be set to accommodate the amplitude peak, each component is measured at -30 dB relative to full scale. This effectively reduces the dynamic range of the measurement by about 30 dB! If the frequency response has a variation from 0 to -100 dB within the measurement span, then each bin of the FFT must measure signals from -30 dBfs to -130 dBfs. Even with a large number of vector averages, this proves difficult, especially with large measurement spans.

Swept sine measurements, on the other hand, can optimize the measurement at *each frequency point*. Since the source is a sine wave, all of the source energy is concentrated at a single frequency, eliminating the 30 dB chirp dynamic range penalty. In addition, if the transfer response drops to -100 dBV, the input range of Channel 2 can auto range to -50 dBV and maintain almost 100 dB of signal to noise. In fact, simply optimizing the input range at each frequency can extend the dynamic range of the measurement to beyond 140 dB!

For frequency responses with both gain and attenuation, the source amplitude can be optimized at each frequency. Reducing the source level at frequencies where there is gain prevents overloads and increasing the amplitude where there is attenuation preserves signal to noise. To optimize the measurement time of sweeps covering orders of magnitude in frequency, the detection bandwidth can be set as a function of frequency. More time can be spent at lower frequencies and less time at higher frequencies. In addition, frequency points can be skipped in regions where the response does not change significantly from point to point. This speeds measurements of narrow response functions.

The figure below illustrates the difference between FFT and swept sine when measuring the frequency response of an elliptic low-pass filter. This filter has a stop band of -80 dB and a zero of about -100 dB. The 400 point FFT measurement was made with 2500 vector averages taking about 10 seconds to complete. The swept sine also took 400 points. Each point was averaged for 16 ms or 10 cycles, whichever was longer. The entire sweep also takes about 12 seconds.

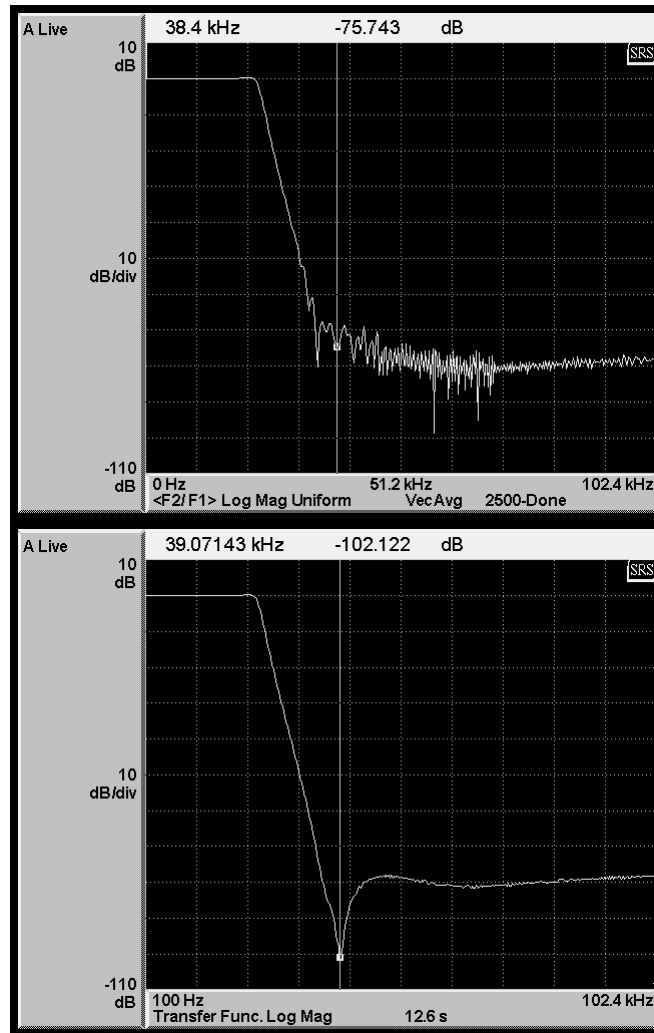

*FFT Measurement*

*Swept Sine Measurement*

*Figure Chapter 2 -2 Frequency responses*

The range of the FFT measurement is limited to about 80 dB. This is because the response signal to Channel 2 contains frequency components within the filter pass band. These components add to create amplitude peaks near 0 dBV (for a 1V chirp amplitude). This requires the input range of Channel 2 to be set near 0 dBV even though each pass band component is only at -30 dBV. The noise floor at full span is about -100 dBV. Thus, the dynamic range of this measurement is roughly 70 dB (pass band to noise floor). Vector averaging in this case improves the dynamic range to about 80 dB.

The swept sine sweep measures each frequency alone, thus optimizing each frequency point individually. This measurement was taken using Auto Range to adjust the input range at each frequency for the best signal to noise. When the frequency is in the stop band, the return signal to Channel 2 is -80 dBV (or less). The input range of Channel 2 is adjusted to -50 dBV for these frequencies eliminating the noise floor limitation. The depth of the zero in the response as well as the shape of the stop band are clearly resolved.



## Swept Sine Measurement Setup

The advantages of the swept sine are apparent in the above comparison. The real usefulness of this technique is that all of the optimizations can be automated in the SR785.

To make a swept sine measurement, select Swept Sine as the Measurement Group. This makes both displays swept sine measurements. Capture and Waterfall are unavailable in the Swept Sine group.

The sweep frequency parameters are set in the [Frequency] menu. The settling and integration time for each point is set in the [Average] menu. The source amplitude is set in the [Source] menu.

Simply press the [Auto Range Ch1] and [Auto Range Ch2] keys to turn on Auto Ranging. Auto Ranging always tracks the input signals during a swept sine measurement.

The signal connections are the same as for an FFT frequency response. The SR785 measures Channel 2 response divided by Channel 1 as the frequency response. The input to the device under test is measured by Channel 1 (typically this is the source output) and the device output is measured by Channel 2.

Swept sine displays differ from FFT displays. The sweep span is not limited to factors of two and the start and stop frequencies are not related to an FFT span. Also, the number of points can range from 10 to 2047. Y axis scaling and X axis zooming are the same as with FFT displays.

Each data point is graphed at its correct frequency and amplitude and connected to adjacent points with a line. The marker only moves to the actual sweep points and the Marker Position Bar shows the actual frequency of each point. If the marker position is displayed with a '?', the point is interpolated from nearby measured points. The actual measurement of these points was skipped due to sweep Auto Resolution.

While a sweep is in progress, a small triangular marker moves across the bottom of the display to show the current sweep point. The current frequency is displayed at the top of the screen.

## Swept Sine Measurements

Only Spectrum, Cross Spectrum, Normalized Variance, Frequency response and User Function measurements are available. Time record, capture and waterfall are not available in this measurement group.

### Spectrum

The swept sine spectrum is simply the measurement of a single channel over a sweep. The spectrum is complex (it contains phase and amplitude information). The phase is relative to the source and is stable but arbitrary. Single channel phase is not generally meaningful.

The spectrum measures the actual signal level at the inputs. If Source Auto Level is On, then the spectrum will tend to be constant. Use Frequency response to remove the effects of a changing source level.

### Cross Spectrum

The swept sine cross spectrum is a two channel measurement defined as

$$\text{Cross Spectrum} = \text{conj}(\text{Spec1}) \cdot \text{Spec2}$$

The cross spectrum contains both magnitude and phase information. The phase is the relative phase (at each frequency) between the two channels.

The magnitude is simply the product of the magnitudes of each spectrum. Frequencies where signal is present in both spectra will have large components in the cross spectrum.

### Frequency response

The swept sine frequency response (sometimes called frequency response) is a two channel measurement defined as

$$\text{Frequency response} = \text{Spec2} / \text{Spec1}$$

The frequency response contains both magnitude and phase information. The phase is the relative phase (at each frequency) between the two channels.

Frequency response measures the response of a network or device under test . The reference channel (1) measures the signal at the input to the device and the response channel (2) measures the device output. The result is the complex frequency response of the device.

### Normalized Variance

Normalized variance is a single channel swept sine measurement similar to the coherence measurement of the FFT measurement group. The normalized variance is an indication of the signal to noise ratio after integration over the specified number of cycles is performed. Normalized variance, like coherence, has a value between 0 and 1. A value of 1 indicates that all noise has been removed from the signal. Values less than one indicate that there is noise remaining. Increasing the Integration Time or Integration Cycles will move the variance closer to 1.

### User Function

User Function displays the results of a user defined math function. User Functions defined within the Swept Sine Measurement Group may include swept sine measurement results. Use the [User Math] menu to define a math function.

A User Function may not be selected as the measurement if it uses a Trace which is empty or which contains data taken in a different Measurement Group.

See 'User Math' later in this section for more.

### Averaging - Settling and Integration

At each frequency point in a swept sine measurement, the inputs measure the amount of signal at the source frequency. This is done by multiplying the input data by  $\sin(\omega t)$  and  $\cos(\omega t)$  and averaging the results over an Integration Time. As with the FFT, real (in-phase) and imaginary (quadrature) signals are measured yielding both magnitude and phase frequency responses.

The Integration Time is always an exact number of cycles at the source frequency. This rejects signals which are at a different frequency, such as noise and harmonics. A long Integration Time results in a narrow detection bandwidth at the source frequency. This improves signal to noise at the cost of longer measurement times. This is similar to the linewidth of an FFT. However, in the FFT the linewidth is exactly related to the frequency span (Span/FFT Resolution). In swept sine, the Integration Time is independent of the frequency span. Thus, wide spans can be measured with narrow linewidths.

The Integration Time is specified in both time and cycles of the source. Times are converted to the next larger exact number of cycles. The larger of the two specifications is used as the Integration Time. A minimum of 1 cycle or 15.6 ms is always measured. To measure each point with the same linewidth, set the Integration Cycles to 1 and the Integration Time to 1/linewidth. To measure each point for a time inversely proportional to the frequency, set the Integration Time to 15.6 ms and the Integration Cycles to the desired number. Remember, the detection bandwidth increases with frequency in this case (the cycles get shorter) which may result in increased detected noise at higher frequencies.

When the integration is complete, the source moves to the next frequency in the sweep. A Settle Time is allowed to pass before any measurement is made at the inputs. This allows the device under test to respond to the frequency change. This can be especially important if the device under test has a high Q. The Settle Time is also specified in both time and cycles of the source. Times are converted to cycles and the larger of the two specified cycles is used.

The Integration and Settle Times are set within the [Average] menu.

### **Sweep Frequency and Auto Resolution**

The span of a swept sine sweep is determined by the Start and Stop frequencies. The entire 102.4 kHz frequency range of the SR785 is available for swept sine measurements. Note that starting a sweep at DC is not possible. In fact, beware of starting at any frequency much less than 1 Hz since the Settle and Integration times are always a minimum of 1 cycle. If the start frequency is very low, the first point can take longer to measure than the entire remainder of the sweep!

The Sweep Number Of Points, or resolution, can be set from 10 to 2047. The points can be in a linear or a logarithmic progression. In many cases it is desirable to sweep over a wide frequency range while still detecting narrow features in the response function. An example might be a filter with many zeroes or a narrow notch. In order to resolve the narrow features, a large number of points must be used in the sweep to improve the frequency resolution. However, a large amount of the sweep time will be spent measuring points between the features of interest. This is where Auto Resolution can save measurement time while preserving resolution.

Auto Resolution is specified by three parameters, the Faster Threshold, the Slower Threshold and the Maximum Step Size. As with all frequency parameters, these are set within the [Frequency] menu.

Auto Resolution examines the measurements of successive frequency points. If the newest measurement is within the Faster Threshold of the previous measurement (for BOTH channels), then the sweep will take larger steps, skipping frequency points. Each successive time this threshold is met, the step size is increased until the Maximum Step Size is reached. This speeds up the sweep in regions where the response is flat (varies less than the Faster Threshold).

If a measurement differs from the previous measurement by more than the Slower Threshold (for EITHER channel), then the sweep returns to the previously measured point and moves to the very next frequency point in the sweep (with no skipping). The sweep continues from this point, speeding up if allowed and slowing down when required. This ‘fills’ in skips in the sweep which vary by more than the Slower Threshold.

Measurements which differ by more than the Faster Threshold (on EITHER channel) but less than the Slower Threshold (on BOTH channels), maintain the present sweep speed. The number of points skipped remains the same in this case.

In order to adjust these parameters appropriately, some knowledge of the frequency response is required. In order to save time, the Faster Threshold must be set to allow some speeding up. Set the Faster Threshold to the desired amplitude resolution remembering that some sacrifice is required to save time. In order to avoid missing the major features in the response, set the Slower Threshold to less than half of the feature size. If the feature is narrow, set the Maximum Step Size to less than half of the feature width (in sweep points).

It is good practice to start by taking a sweep with Auto Resolution Off. This ensures that the span and resolution of the sweep are correct as well as the averaging times. Once these have been established, turn Auto Resolution On to save time on successive measurements. Start by setting the Faster Threshold to something greater than the variations of the flat regions of little interest. Set the Slower Threshold to twice the Faster Threshold and the Maximum Step Size to no more than 5% of the Sweep Number Of Points. Adjust the parameters as necessary to optimize the sweeps where needed.

A well specified Auto Resolution sweep can take 80% off of the sweep time while preserving the essential features of the frequency response. This is very useful in repetitive measurements.

### **Input Auto Ranging**

Input auto ranging has the largest impact on extending the dynamic range of a swept sine measurement. It is also the simplest optimization to setup and use. Simply press the [Auto Range] keys to turn on Auto Ranging. Auto Ranging always tracks the input signals during a swept sine measurement. At each frequency point, each input range is adjusted to keep the signal within -6 dB of full scale (if possible). If the input is overloaded, then the range is increased until the overload is removed.

Auto Ranging increases the measurement times whenever an input range is changed. This is especially true when the signal drops below -6 dBfs. In order to detect this under-range condition, the measurement must be made for the entire integration time. At low

frequencies, this can add a significant amount of sweep time. Overloads can be detected at any time during the integration time and thus ranging up requires less time.

For smooth frequency responses where successive measurement points do not differ by more than 6 dB, auto ranging does not add substantially to the sweep time while increasing the dynamic range dramatically. In some cases, especially at very low frequencies, it might be better to set the range manually to avoid increasing an already long sweep time.

Auto Ranging is set for each channel separately and can be turned on and off during the sweep. The manual input range can be changed at any time during the sweep as well.

### Source Auto Level and Ramping

The source amplitude parameters for a swept sine measurement are set within the [Source] menu. Normally the source amplitude is maintained at a constant level at all frequencies during the sweep. This usually works best for frequency responses which are mostly attenuating and have little or no gain. Simply set the source amplitude such that the device under test is not overloaded anywhere in the sweep.

Source Auto Level will adjust the source amplitude to maintain a constant level, called the Ideal Reference, at the Channel 1 or Channel 2 input. This is useful whenever the frequency response has substantial gain as well as attenuation. Suppose the frequency response of the device under test has +30 dB of gain at one point and -80 dB of loss at another point and overloads at an input or an output of 1.0 V. If a constant source level is used, it must be set to avoid overload during the gain portion of the sweep. This requires a level of -30 dBV. When the sweep reaches the attenuation region, the output signal will drop to -110 dBV. While this is still measurable, it may not be optimum. With Auto Level, the source can try to maintain an Ideal Reference of 1.0 V at the output of the device under test (input of Channel 2) while not exceeding 1.0 V at the device input. Where the gain is +30 dB, the source level is -30 dBV and where the attenuation is -80 dB, the source level increases to the Maximum Source Level (1.0 V in this case). The actual signal output from the device under test varies from 0 dBV to -80 dBV instead of 0 dBV to -110 dBV in the constant level case. Varying the source level narrows the range of the output signals. This avoids overloads when there is gain and increases the output signal to noise when there is attenuation.

Source Auto Level requires input Auto Ranging to be on. This is because the non-reference channel should follow all of the source changes. In the above example, the source changes by 30 dB at Channel 1.

There are five important parameters to set when using Auto Level. The Auto Level Reference Channel determines which input will be maintained at the reference level. This is usually the device under test output (Channel 2). In cases where the device under test is not driven directly from the source output, the device input (Channel 1) may be the reference. The Ideal Reference is the level the source will try to maintain at the Reference Channel. This is determined by the limitations of the device under test or may simply be the desired signal level at which the test is run. The Maximum Source Level is the largest allowed source amplitude. This is limited by the SR785 source output or the device under test input range. The Reference Upper and Lower limits are the allowable tolerances for the Reference Channel. The source level is changed only if the reference

channel measures an input which exceeds the Ideal Reference by more than the Upper Limit or which is less than the Ideal Reference by more than the Lower Limit. In cases where the reference channel input must be limited, the Ideal Reference plus the Upper Limit must be less than or equal to this limit.

Whenever Auto Level is used, the measurement **MUST** be frequency response. The signal connections must be made for a frequency response measurement. This is because source level changes are not normalized in the individual measurements. Only the ratio of Channel 2 to Channel 1 is source level independent.

The Source Ramp Rate is the rate at which the source level changes. If Source Ramping is Off, source level changes are made instantly. If Source Ramping is On, the source level changes are made at the ramp rate. The Settle Time starts after the source reaches the new level. This prevents fast level changes from disturbing the device under test. Of course, this adds to the sweep time.

The source ramps off whenever there are no measurements being made. This occurs before the start of a sweep, at the end of a single sweep or while a sweep is paused (by the user). If Ramping is Off, the source shuts off instantly in these situations.

# Order Analysis

## What is Order Analysis?

In vibration analysis or preventive maintenance of rotating machines, it is often useful to identify spectral components by their relation to the speed of the rotating machine, rather than their absolute frequency. For instance, if we identify a potentially troublesome vibration component at 3 times the rotation speed during a baseline test of a machine at 120 RPM, all we really care about is that the vibration is at "3rd order," not that it's at 6Hz. That way we can compare the results to the 3rd order vibration level in some future test without having to worry about getting the machine to run at exactly 120 RPM again.

In the Order Analysis measurement group the SR785 calculates spectra whose x-axis is in orders, i.e. multiples of the machine speed. To do this, it requires information about the speed and position of the machine's shaft. This is provided by a tachometer input, which provides a signal at fixed positions relative to the shaft rotation. With the SR785, all the traditional measurements of rotating machinery analysis including order maps, rpm profiles, run-up and run-down plots, can be easily made without any additional equipment.

## Order Spectra and Order Tracking

Two types of measurement are typically made with rotating machinery. The first looks at the input amplitude as a function of order. This type of measurement is similar to the standard FFT, except that the x-axis is calibrated in orders instead of in absolute frequency. Typically, a measurement of the order spectrum is first performed to identify orders for further study.

Once we've identified the problem orders, we can move to the second type of measurement, where the amplitude of a specific order, or set of orders is recorded as the machine speed is changed. We may want to ramp the machine from startup to operating speed, and track the behavior of the third order vibration.

The SR785 can perform both types of measurements by selecting a value for the [Freq]<Tracking> softkey. When tracking is turned on the analyzer records the amplitude for each spectral order as a function of RPM. Using the [Freq]<Track Setup> submenu, a 2-dimensional plot of the amplitude of any order of either input channel can be generated. Using RPM trigger arming, the analyzer can be configured so that the tracked data points are acquired at fixed RPM intervals.

The tracking feature uses the SR785's Tracking/Waterfall memory. Thus when tracking is turned on, waterfall storage of spectra is not allowed. If waterfall storage is necessary, turn tracking off.

## Order Measurement Setup

Choose the Maximum RPM, Minimum RPM, Maximum Order, and Delta Order to setup an Order measurement. Maximum RPM affects the available range of Max Order; the higher the selected Maximum RPM, the lower the maximum order that can be selected. The relationship can be summarized as:  $\max(\text{Max RPM}) = 60 \cdot 40,000 / (\text{Max Order})$ . Minimum Order does not affect the measurement at all, it merely sets a default scaling for measurements whose x-axis is calibrated in RPM, such as Order Track.

The combination of Maximum Order and Delta Order also determine the length and resolution of the Order time records. The resolution of the time record, in revolutions, is given by:  $\text{Resolution} = 1/(5.12 \cdot \text{Effective Max Order})$ , where Effective Max Order is equal to  $\text{Delta Order} \cdot 400$ , ( if  $\text{Max Order}/\text{Delta Order} > 200$ ),  $\text{Delta Order} \cdot 200$ , ( if  $100 < \text{Max Order}/\text{Delta Order} \leq 200$ ),  $\text{Delta Order} \cdot 100$ , ( if  $50 < \text{Max Order}/\text{Delta Order} \leq 100$ ), or  $\text{Delta Order} \cdot 100$ , (If  $\text{MaxOrder}/\text{Delta Order} \leq 50$ ).

The length of the time record, in revolutions, is  $1/(\text{Delta Order})$ .

The combination of these parameters also determines the Capture Sampling Rate. While in other measurement groups the Capture Sampling Rate is determined explicitly by the user with the sampling rate softkey. In the order measurement group, the sampling rate is given by:  $\text{Sampling Rate} = 5.12 \cdot \text{Effective Max Order}$ , where the sampling rate is then rounded up to the nearest integer submultiple of 256 kHz.

### Order Measurements

#### Linear Spectrum

The Linear Spectrum in the order measurement group is similar to the FFT measurement in the FFT measurement group except that the x-axis is calibrated in orders instead of absolute frequency. The Linear Spectrum includes both magnitude and phase information. Phase in the Linear Spectrum is measured relative to the tachometer pulse unless an external trigger is used in which case the phase is measured relative to the trigger.

#### Power Spectrum

The Linear Spectrum in the order measurement group is similar to the Power Spectrum measurement in the FFT measurement group except that the x-axis is calibrated in orders instead of absolute frequency. The power spectrum gives the power present in each spectral order, no phase information is included.

#### Time Record

The time record is not similar to the time record measurement in the FFT measurement groups. Time Records in the fft group are measurements of the input signal at fixed time intervals. In the Order measurement group, time records are measurements of the input signal at fixed angular positions around the shaft. Since the analyzer only samples the input at fixed time intervals, the order time record must be calculated by combining the fixed interval sampling with information from the tachometer about where the shaft position was relative to the sampled data.

The x-axis of the order time record is calibrated in revolutions, the resolution and number of resolutions in the time record are determined from the Maximum Order and Delta Order softkeys.

#### Windowed Time Record

The windowed time record shows the Time Record multiplied by the Window function. Remember that in the order measurement group you should use the uniform window if the expected harmonics are integer multiples (orders) of the rotating machine speed.



**RPM Profile**

The RPM profile measurement is a two-dimensional plot showing RPM on the x-axis and time on the y-axis.

**Orbit**

The orbit measurement is a two-dimensional plot of the Channel 2 time record vs. the Channel 1 time record. The orbit measurement is often used in shaft balancing applications.

**Track**

The Order Track measurement is a two dimensional plot of the amplitude of a given order vs. RPM. The Track measurement can only be selected when Tracking is turned on in the [Frequency] menu. The order and input channel displayed by the track measurement is selected with the Track Setup submenu.

The amplitude recorded in the track is the Linear Spectrum, and contains both amplitude and phase information.

Since the tracking feature uses the memory allocated to the Waterfall/Order buffer, the waterfall buffer is not available for waterfall storage when tracking is enable.

**Capture**

The capture buffer may be used as a source of data for order measurements. See Capture Buffer for more details. The Capture measurement displays the contents of the capture buffer. Order measurements can use this stored data as input data by choosing Playback as the Input Source in the [Input] menu.

In the order analysis measurement group, the capture buffer must store information about the tachometer input as well as the signal inputs. Therefore, only capture buffers that were taken while the analyzer is in the order measurement group may be used as a source of data for order measurements during capture playback.

**User Function**

User Function displays the results of a user defined math function. User Functions defined within the Order Measurement Group may include order measurement results. Use the [User Math] menu to define a math function.

A User Function may not be selected as the measurement if it uses a Trace which is empty or which contains data taken in a different Measurement Group.

# Time/Histogram Measurements

## Time and Histogram Measurements

Time and Histogram measurements give a statistical picture of an input signal in the time domain only. Histogram measurements answer questions such as "How often is my signal between 1V and 1.1 V, or "What's the probability that my signal exceeds .5V." Time and Histogram measurements are useful for characterizing noise as well as in many mechanical applications.

## Measurements in the Time/Histogram Group

Time Records, Histograms, Probability Density Functions (PDF), Cumulative Density Functions (CDF), and Capture are the measurements available in the Time/Histogram group. Time/Histogram measurements may be stored in the waterfall buffer and the Capture Buffer may be used as in input.

### Histogram

The Histogram displays the distribution of amplitudes found in the input signal over a time determined by the Histogram Length. Histogram Length may be set in terms of the number of points in the histogram, or the length of time to acquire the histogram

The x-axis of the histogram represents the input amplitude range of the SR785, from negative full scale (plus an overhead factor of 1.25) to positive full scale (plus an overhead factor of 1.25) divided into a number of histogram "bins" given by the Bins softkey. While the histogram is being acquired, the input signal is sampled at intervals given by Sampling Time and the number of samples in each histogram bin is recorded. If Repeat has been selected, the histogram will restart automatically after finishing.

The x-axis unit of the histogram measurement is volts, indicating the voltage corresponding to each histogram bin, while the y-axis unit is counts, indicating how many samples are in the bin. If Engineering Units are turned on, the x-axis will be in the selected Engineering Units.

### Time Record

The Time/Histogram Time Record measurement is similar to the FFT group time record. The key difference is that while the FFT group time record is digitally filtered to eliminate aliases, the Time/Histogram group time record is unfiltered. If the sample time is set to 1s, the analyzer simply samples the input every second. This allows examination of the input signal without the effects of additional filtering. If desired, the input anti-aliasing filters can be turned off to present an even "purer" view of the input signal. The number of samples in a time record is fixed at 1024 in the Time/Histogram group.

### Probability Density Function (PDF)

The probability density function, or PDF, is similar to the histogram, except it is normalized so that the integral of the PDF between any two points is equal to the probability that the amplitude of the input signal is between those two points. Since the integral of the PDF over voltage is a probability, i.e. unitless, the units of the PDF must be inverse volts.

**Cumulative Density Function (CDF)**

The cumulative density function, or CDF, is an integrated version of the PDF. The value of the CDF at a voltage  $V$  is the integral of the PDF from the smallest amplitude up to  $V$ . Thus, the value of the CDF at  $V$  is the probability that the input signal has an amplitude smaller than  $V$ .

**Capture**

The capture buffer stores sequential time domain data in memory. See 'Capture Buffer' for more details. The Capture measurement displays the contents of the capture buffer. Time/Histogram measurements can use this stored data as input data by choosing Playback as the Input Source in the [Input] menu.

The capture buffer is often very long. To graphically expand a region of the display, use the Pan and Zoom functions in the [Display Setup] menu. The capture buffer display can automatically pan with the capture fill or playback progress through the buffer. During capture fill, if the capture accumulates points faster than they can be displayed, some points are not shown. This speeds up the display update so that it keeps up with the real time capture but allows visual aliasing to occur. Once capture is complete, the display is redrawn showing the envelope of all points, eliminating any visual aliasing effects.

To measure from a region of the buffer, set the Playback Start and Stop in the [Input] menu.

The capture buffer is not a continuous representation of the input signal. The data is sampled and has a time resolution of  $1/(\text{sample rate})$ . High frequency signals will appear distorted in the time record. However, ALL of the spectral information is preserved by the Nyquist sampling theorem as long as the value of each sample is accurate.

**User Function**

User Function displays the results of a user defined math function. User Functions defined within the Time/Histogram Measurement Group may include time/histogram measurement results. Use the [User Math] menu to define a math function.

## Trace Storage

Measurement data may be stored in one of 5 Trace buffers. Trace buffers are simply temporary storage for complex measurement results and are shared by all Measurement Groups. Trace data is not retained when the power is turned off.

Display or Reference Display data may be saved to a Trace. Waterfall records or slices may also be saved to a Trace. Trace data is associated with the Measurement Group in which the data was originally taken. A Trace into which an FFT measurement has been saved is associated with the FFT Group.

A stored Trace can be recalled to the active display as Display data. A Trace may be recalled to a display regardless of its associated Measurement Group. When a Trace is recalled as Display data, the Display becomes Off-Line. Since the Off-Line display is showing static trace data, its Measurement, Window, Frequency Span and Averaging may not be modified. Only display related parameters such as View, Scaling and Marker Functions may be changed Off-Line.

A stored Trace can also be recalled to the active display as its Reference Display. In this case the Measurement Group associated with the Trace must be the same as the Measurement Group of the active display (either Live or Off-Line).

A stored Trace may be used as an operand in a User Math Function. In this case, in order to display the User Function, the Trace must contain data AND the Trace must have data stored from within the current Measurement Group. For example, an FFT User Function which uses Trace1 cannot be chosen as a Measurement if Trace1 is empty OR contains octave analysis data.

To preserve Trace data permanently, save the Trace to disk.

A Trace may be recalled from a disk file or loaded from an interface. If a User Function is displayed and it uses a Trace in its equation, then that Trace cannot be changed to data from a different Measurement Group. For example, if an FFT User Function using Trace1 is being displayed, you cannot recall octave data from disk to Trace1.

# User Math Functions

## What is a User Function?

User Functions are used to define measurements which are not already available in each Measurement Group. A User Function is an equation defined by the user in the [User Math] menu. User Functions use Operands such as FFT(1), FFT(2), Traces or Constants, and Operators such as +, -, x, / as well as FFT( ), IFFT( ), Avg( ), Mag( ), Phase( ), etc.

To make a measurement of a defined User Function, select the User Function as the Measurement for the active display. User Functions behave exactly the same as the predefined measurements in general. The display will show the measurement results using real time input data (if the function requires input data). Frequency, Averaging and Displays operate normally for most User Functions. Triggering operates normally within each Measurement Group (a function using an FFT time record is triggered normally).

User Functions are saved and recalled to disk with the instrument Settings.

## Measurement Groups and Traces

Each Measurement Group contains 5 User Functions in addition to its predefined measurements. User Functions can only be measured (displayed) within their own Measurement Group. Each Measurement Group allows different Operands (FFT(1), Oct(1) or Swept Sine Spec(1)) as well as different Operators.

User Functions can also use the 5 stored Traces as operands. In this case, in order to display this User Function, the Trace must contain data AND the Trace must have data stored from within the current Measurement Group. For example, an FFT User Function which uses Trace1 cannot be chosen as a Measurement if Trace1 is empty OR contains octave analysis data.

If a User Function is displayed and it uses a Trace in its equation, then that Trace cannot be changed to data from a different Measurement Group. For example, if an FFT User Function using Trace1 is being displayed, you cannot recall octave data from disk to Trace1.

Changing Measurement Groups and storing Traces may make User Functions unavailable because of the above restrictions. In this case, store Trace data within the current Measurement Group before selecting the User Function as the measurements.

## Operands

Operands are arrays of complex values (real and imaginary). They may be measurement results (such as FFT(1), Time(1) or Oct(1)), a Trace, or a Constant. An array which is real simply has zero for its imaginary parts. The array length of an Operand is determined by the measurement length (number of FFT lines, length of time record, number of octaves, etc.) or Trace length (length of the data which is stored in the Trace). Constants assume the length of the user function.

Operands which are measurement results enclosed in angle brackets, such as <Freq. Resp.> or <Spec(1)>, are exactly the same as the normal measurements. They are

computed using whatever averaging is selected with the <Display Average> softkey. They are computed from the input data (real time Analog or capture Playback) and use the frequency and windowing parameters from the menus. Frequency domain measurements are amplitude calibrated, time records are not. FFT and Time record measurements are also triggered just like the normal measurements.

Operands which are measurement results which are not enclosed in angle brackets, such as FFT1 or Time1 represent unaveraged, instantaneous versions of the measurement.

Operands which contain an explicit averaging type, such as Vec<F1> or PeakHold<F2>) always are averaged according to their indicated type, regardless of the setting of the <Display Average> softkey.

Octave and Swept Sine measurements are always averaged measurements. Use the [Average] menu to set the averaging parameters.

FFT(1) and FFT(2) are the FFT of the Ch1 and Ch2 inputs. These operands use the window chosen in the [Window] menu for the display which is measuring the function.

FFTu(1) and FFTu(2) are un-windowed FFT's of the Ch1 and Ch2 inputs.

Correlation operands such as RMS<Fu1•F1> are used to compute the RMS and Vector averaged auto and cross correlation measurements. The Fu refers to the un-windowed FFT. See the description of the Correlation measurements for exact definitions of the correlation measurements in terms of these operands.

Trace operands are simply the data stored in the Traces. For example, Traces can hold reference data used for normalization or calibration. There are 5 Traces which can be stored. These Traces are shared by all 3 Measurement Groups.

Constants are complex constants which are the same for every array point. Constants such as pi are commonly used in equations. There are 5 user defined constants which are shared by all 3 Measurement Groups.

### **X Axis**

Operands have an X axis type (frequency or time) as well as an X axis (array) length. The X axis type of the User Function determines the X axis labeling of its display.

In general, the length of a User Function is determined by the length of the longest operand array. Operations between operand arrays are performed on a point by point basis, starting at the beginning of each array. Operands of different lengths (400 point FFT1 and 800 point FFT1 stored in a Trace) can be combined in an User Function. The math is performed over the longest operand with the shorter operand array being used in a circular fashion. This usually leads to meaningless results.

Measurement operands (FFT(1), Time(1), etc.) and Traces have an X axis type, either frequency or time domain. Arrays of frequency data start at the lowest frequency. Arrays of time data start at time=0. A frequency domain operand (FFT(1)) can be added to a time domain operand (Time(1)) point by point even though the result is meaningless. The X axis type of a User Function is determined by examining its equation from left to right

and keeping track of domain switching (IFFT and FFT) operations. Operand terms with an incompatible X axis type have no effect on the function's X axis type, though the entire equation is still evaluated.

Constants are simply constant for each point. They also have no effect on the function's X axis type. If a User Function is defined simply in terms of constants (no other operands), then the X axis type defaults to frequency.

If a User Function's X axis type or length is determined by a Trace, then changing the data in the Trace may change the function's X axis type or length as well.

## Operations

Operands are combined with Operations to define a function.

Arithmetic operations (+, -, x, /) combine operand terms on a point by point basis. It is the user's responsibility to ensure that the operand terms have the correct X axis type and lengths in order to produce meaningful results.

View operations ( Mag( ), Mag<sup>2</sup>( ), Phase( ), Real( ), Imag( ) ) simply convert the complex operand array into the desired form. Mag( ) and Mag<sup>2</sup>( ) compute the magnitude ( $\sqrt{x^2 + y^2}$ ) or magnitude squared ( $x^2 + y^2$ ) with a real result. Phase( ) computes the phase ( $\tan^{-1}(y/x)$ ) unwrapped with a real result. Real( ) simply zeroes the imaginary part. Imag( ) zeroes the real part.

Conj( ), Ln( ), Exp( ) and Sqrt( ) are defined as follows,

$$\begin{aligned} \text{Conj}(x + jy) &= x - jy \\ \text{Conj}(x - jy) &= x + jy \\ \text{Ln}(x + jy) &= \ln(r) + j\theta \\ \text{Exp}(x + jy) &= \exp(x) \cdot (\cos(y) + j\sin(y)) \\ \text{Sqrt}(x + jy) &= \sqrt{r} \cdot (\cos(\theta/2) + j\sin(\theta/2)) \\ [X/(1-X)] &= (x + jy)/(1 - x - jy) \end{aligned}$$

where  $r = \sqrt{x^2 + y^2}$  and  $\theta = \tan^{-1}(y/x)$ . Sqrt computes all angles positive from 0. The  $[X/(1-X)]()$  operator is generally used on real arguments only.

The  $j\Omega$  operator is simply  $j\omega$  and zeroes the real part and fills the imaginary part with the value of  $\omega$  ( $2\pi \cdot \text{frequency}$ ) at each point in the array. The frequency is determined by the span of the operand. A measurement operand (FFT(1) or Time(1)) uses the current frequency span. A Trace operand uses the span of the measurement which is stored in the Trace. In the Order Analysis measurement group the current RPM is used to generate the frequency axis.

FFT( ) is the windowed FFT operator. These operations switch the X axis type. The current window chosen in the [Window] menu is used on the operand before the FFT. FFTu( ) is the un-windowed FFT operator. The uniform window is always used. The FFT( ) and FFTu( ) operators do not calibrate their results since their operands are not necessarily input time records. To use a calibrated FFT, use the measurement operand FFT(1) or FFT(2).

IFFT( ) is the inverse complex FFT operation. This operation switches the X axis type.

d/dx( ) is the derivative operator. The derivative is performed with respect to the X-axis bin number, not the x-axis values. For example, to convert d/dx( ) to d/df( ) for a linear frequency display, divide by the frequency spacing of the bins. In the case where the X-axis is logarithmic, the function needs to correct for  $d(\log x)/dx$ . The d/dx( ) operator requires a smoothing aperture or window.

GrpDly( ) is the group delay operator. Group delay is simply  $d\theta/d\omega$ . Use [Display Options]<d/dx Window> to set the aperture. This operator yields valid results for any complex operand which has a frequency x-axis.

AWt( ) is the A-Weighting operator. This filter is a software version of the input A-Weight filter. BWt( ) and CWt( ) are the software B-Weighting and C-Weighting operators. The AWt( ), BWt( ), and CWt( ) operators should only be used on operands with a frequency x-axis. All three filters conform to ANSI Standard S1.4-1983.

### User Function Limits

User Functions and Constants are defined in the [User Math] menu. Only the 5 functions within the current Measurement Group are available for editing in this menu.

Editing a function which is currently being displayed is allowed. In this case, Traces which do not contain compatible measurements are not allowed as operands.

User Functions have a maximum length of 31 operators plus operands. User Functions cannot use another User Function as an operand.



# Signal Inputs

The Input Range on the SR785 varies from a maximum of 34 dBV full scale to a minimum of -50 dBV full scale. A signal which exceeds the current input range will cause the OverLoad indicator to appear at the top of the screen. A signal which exceeds the maximum safe range (35 dBVpk) will turn on the Hi V indicator and set the range to 34 dBV.

The Input Range is displayed in dBV. The maximum and minimum range equivalents are tabulated below.

Max	34 dBVpk 31 dBVrms 50.1 Vpk 35.4 Vrms
Min	-50 dBVpk -53 dBVrms 3.16 mVpk 2.24 mVrms

## Manual Range

The fixed Input Range is specified in the [Input] menu. Signals that exceed the range will overload and become distorted. Signals which fall to a small percentage of the range will lose dynamic range.

## Auto Range

The Input Range can be set to automatically correct for signal variations. Use the [Auto Range Ch1] and [Auto Range Ch2] keys to toggle Auto Ranging On and Off. Specify either Normal or Tracking Auto Range in the [Input] menu.

When Normal Auto Ranging is On and an overload occurs, the Input Range is increased so that the signal no longer overloads. If the signal decreases, the input range is not adjusted. You must take care to ensure that the signal does not fall dramatically after pushing the Input Range to a very (high) insensitive setting.

Tracking Auto Range moves the Input Range up when an overload occurs and down when the signal falls below half full scale.

Auto Range is always Tracking for swept sine measurements.

## Input Noise

Noise measurements in FFT analyzers are usually made using Power Spectral Density (PSD) units. PSD units simply normalize the amplitude of each frequency bin by the square root of the linewidth. A window correction factor is also included to compensate for the bandwidth of the different windows. PSD units have the form  $V_{rms}/\sqrt{Hz}$  or  $dBV_{rms}/\sqrt{Hz}$ .

Most noise sources, including the amplifier input noise, are Gaussian in nature. The amount of noise is proportional to the square root of the bandwidth in which the noise is measured. A noise measurement made using normal units would yield different amplitudes depending upon the linewidth of the FFT span. By using PSD units, the measured noise is independent of linewidth and span.

The input noise of the SR785 inputs is about 10 nVrms/ $\sqrt{\text{Hz}}$ . The input noise of the analog to digital (A/D) converter is about 300 nVrms/ $\sqrt{\text{Hz}}$  (referenced to a full scale of 1 Vpk). What does this mean in terms of the FFT noise floor?

For input ranges below -30 dBV (input gain greater than 30 dB), the measured noise floor is determined by the signal inputs. This is because the signal gain (>30) times the input noise (10 nVrms/ $\sqrt{\text{Hz}}$ ) exceeds the input noise of the A/D converter.

For input ranges above -30 dBV (input gain less than 30 dB), the measured noise floor is determined by the A/D converter. There is not enough gain to amplify the input noise to a level greater than the A/D converter itself. The noise floor in this case is 300 nVrms/ $\sqrt{\text{Hz}}$  relative to 1 Vpk or -130 dBpk/ $\sqrt{\text{Hz}}$  relative to full scale.

All of this assumes that the signal input is driven from a low impedance (and noiseless) source. Remember, resistors have a Johnson noise equal to  $0.13 \times \sqrt{R}$  nVrms/ $\sqrt{\text{Hz}}$ . Even a 50  $\Omega$  resistor has almost 1 nVrms/ $\sqrt{\text{Hz}}$  of noise! A signal source impedance of 6 k $\Omega$  has a Johnson noise equal to the SR785's input noise. To determine the overall noise of multiple noise sources, take the square root of the sum of the squares of the individual noise figures. For example, a 6 k $\Omega$  source impedance has 10 nVrms/ $\sqrt{\text{Hz}}$  of noise. Added to the input noise, the overall noise is  $(10^2 + 10^2)^{1/2}$  or 14 nVrms/ $\sqrt{\text{Hz}}$ .

### Input Impedance

The input impedance of the SR785 is 1 M $\Omega$ . If a higher input impedance is desired, then the SR550 remote preamplifier must be used. The SR550 has an input impedance of 100 M $\Omega$  and is AC coupled above 1 Hz.

### Anti-aliasing Filter

After all of the signal amplification, there is an anti-aliasing filter. This filter is required by the signal digitization process. According to the Nyquist criterion, signals must be sampled at a frequency at least twice the highest signal frequency. In this case, the highest signal frequency is 102.4 kHz and the sampling frequency is 262 kHz so things are ok. However, no signals above 131 kHz can be allowed to reach the A/D converter. These signals would violate the Nyquist criterion and be undersampled. The result of this undersampling is to make these higher frequency signals appear as lower frequencies in the digital data stream. Thus a signal at 175 kHz would appear below 102.4 kHz in the digital data stream and be detected in the spectrum.

To avoid this undersampling, the analog signal is filtered to remove any signals above 160 kHz (when sampling at 262 kHz, signals above 160 kHz will appear below 102 kHz). This filter has a flat pass band from DC to 104 kHz so as not to affect measurements in the operating range of the analyzer. The filter rolls off from 104 kHz to 160 kHz and achieves an attenuation above 160 kHz of at least 90 dB. Amplitude variations and phase shifts due to this filter are calibrated out at the factory and do not affect frequency domain measurements.

For time domain measurements, the anti-aliasing filter limits the rise time of the signal and causes ringing. To eliminate these effects in the time record, the anti-aliasing filter can be removed from the signal path. However, the spectrum of this time record may contain alias frequencies.

### **A-Weighting Filter**

The A-Weighting filter may be inserted after the anti-aliasing filter. The A-Weighting filter simulates the hearing response of the human ear and is often used with Octave Analysis measurements. The input A-Weighting filter conforms to ANSI standard S1.4-1983.

A-, B- and C- Weighting functions are also available as operators in User Math functions.

### **Input Transducer Units**

Transducers (such as accelerometers or microphones) convert a physical quantity (such as acceleration or pressure) into a voltage at the analyzer's input. By assigning Transducer Units to an input, measurements based upon the transducer signal may be displayed in units of the actual physical quantity being measured (such as  $m/s^2$  or Pascals).

Transducer Units require a transducer which is linear over the frequency range of interest.

Measurements of inputs which are assigned units of acceleration, velocity or displacement may be displayed with any of these three units using Transducer Convert in the [Display Options] menu.

# Input Connections

In order to achieve the best accuracy for a given measurement, care must be taken to minimize the various noise sources which can be found in the laboratory. With intrinsic noise (Johnson noise, 1/f noise or input noise), the signal source must be designed with these noise sources in mind. These noise sources are present regardless of the input connections. The effect of noise sources in the laboratory (such as motors, signal generators, etc.) and the problem of differential grounds between the signal source and the analyzer can be minimized by careful input connections.

There are two basic methods for connecting a voltage signal to the analyzer - the single-ended connection is more convenient while the differential connection eliminates spurious pick-up more effectively.

## Single-Ended Connection (A)

In the first method, the analyzer uses the A input in a single-ended mode. The analyzer detects the signal as the voltage between the center and outer conductors of the A input only. The analyzer does not force the shield of the A cable to ground, rather it is internally connected to the analyzer's ground via a resistor. The value of this resistor is selected by the user. Float uses 1 M $\Omega$  and Ground uses 50  $\Omega$ . This avoids ground loop problems between the signal source and the analyzer due to differing ground potentials. The analyzer lets the shield 'quasi-float' in order to sense the source ground. In general, if the source is floating, use a Grounded input. If the source is grounded, use a Floating input. Do not use a Floating input with a floating source since the shield is simply a noise antenna in this case. Noise pickup on the shield alone will appear as noise to the analyzer. Common mode noise, which appears on both the center and shield, is rejected by the common mode rejection (CMR) of the analyzer input, but noise on the shield only is not rejected at all.

## Differential Connection (A-B)

The second method of connection is the differential mode. The analyzer measures the voltage difference between the center conductors of the A and B inputs. Both of the signal connections are shielded from spurious pick-up. Noise pickup on the shields does not translate into signal noise since the shields are ignored. The shields in this case are grounded by 50  $\Omega$  since they are being used as signal shields.

When using two cables, it is important that both cables travel the same path between the signal source and the analyzer. Specifically, there should not be a large loop area enclosed by the two cables. Large loop areas are susceptible to magnetic pickup.

## Common Mode Signals

Common mode signals are those signals which appear equally on both center and shield (A) or both A and B (A-B). With either connection scheme, it is important to minimize both the common mode noise and the common mode signal. If the signal source floats at a non-zero potential, the signal which appears on both the A and B inputs will not be perfectly canceled. The common mode rejection ratio (CMRR) specifies the degree of cancellation. For low frequencies, the CMRR of 90 dB indicates that the common mode signal is canceled to 1 part in 30,000. Even with a CMRR of 90 dB, a 1 V common mode

signal behaves like a 3  $\mu\text{V}$  differential signal! The CMRR decreases by about 6 dB/octave (20 dB/decade) starting at around 1 kHz.

### **AC vs DC Coupling**

The signal input can be either AC or DC coupled. The AC coupling high pass filter passes signals above 160 mHz (0.16 Hz) and attenuates signals at lower frequencies. AC coupling should be used at frequencies above 160 mHz whenever possible. At lower frequencies, DC coupling is required.

For small signals it is common that the DC offset of the signal exceeds the amplitude of the AC signal of interest. In this case, the input range may be limited to accommodate the large DC offset while sacrificing signal to noise in the measurement of the AC signal. If the signal frequency exceeds 0.16 Hz, use AC coupling if possible.

## Intrinsic (Random) Noise Sources

Random noise finds its way into signals in a variety of ways. Good test design can reduce these noise sources and improve the measurement stability and accuracy.

There are a variety of intrinsic noise sources which are present in all electronic signals. These sources are physical in origin.

### Johnson Noise

Every resistor generates a noise voltage across its terminals due to thermal fluctuations in the electron density within the resistor itself. These fluctuations give rise to an open-circuit noise voltage,

$$V_{\text{noise}} (\text{rms}) = (4kTR\Delta f)^{1/2}$$

where  $k$ =Boltzmann's constant ( $1.38 \times 10^{-23}$  J/°K),  $T$  is the temperature in °Kelvin (typically 300°K),  $R$  is the resistance in Ohms, and  $\Delta f$  is the bandwidth in Hz.  $\Delta f$  is the bandwidth of the measurement (usually the FFT linewidth).

### Shot Noise

Electric current has noise due to the finite nature of the charge carriers. There is always some non-uniformity in the electron flow which generates noise in the current. This noise is called shot noise. This can appear as voltage noise when current is passed through a resistor. The shot noise or current noise is given by

$$I_{\text{noise}} (\text{rms}) = (2qI\Delta f)^{1/2}$$

where  $q$  is the electron charge ( $1.6 \times 10^{-19}$  Coulomb),  $I$  is the RMS AC current or DC current depending upon the circuit, and  $\Delta f$  is the bandwidth (usually the FFT linewidth).

### 1/f Noise

Every 10  $\Omega$  resistor, no matter what it is made of, has the same Johnson noise. However, there is excess noise in addition to Johnson noise which arises from fluctuations in resistance due to the current flowing through the resistor. For carbon composition resistors, this is typically 0.1 to 3  $\mu\text{V}$  of rms noise per Volt applied across the resistor. Metal film and wire-wound resistors have about 10 times less noise. This noise has a 1/f spectrum and makes measurements at low frequencies more difficult.

Other sources of 1/f noise include noise found in vacuum tubes and semiconductors.

### Total Noise

All of these noise sources are incoherent. The total random noise is the square root of the sum of the squares of all the incoherent noise sources. Thus, the largest noise source easily dominates all others in determining the noise floor of the measurement.

## External Noise Sources

In addition to the intrinsic noise sources discussed previously, there are a variety of external noise sources within the laboratory.

Most of these noise sources are asynchronous, i.e. they are not related to the signal and do not occur at the signal frequency or its harmonics. Examples include lighting fixtures, motors, cooling units, radios, computer screens, etc. If these noise sources are large, they can determine the input range and hence the noise floor. They can however be removed from the FFT spectrum by using triggering and vector averaging. Since the noise signals are not phase coherent with the trigger and signals of interest, they vector average to zero.

Some noise sources, however, are related to the signal source and, if picked up in the signal, will add or subtract from the actual signal and cause errors in the measurement. A typical source of synchronous noise is a ground loop between the signal source and the analyzer.

Many of these noise sources can be minimized with good laboratory practice and experiment design. There are several ways in which noise sources are coupled into the signal path.

### Capacitive Coupling

An AC voltage from a nearby piece of apparatus can couple to a signal path via a stray capacitance. Although  $C_{\text{stray}}$  may be very small, the coupled noise may still be larger than a weak signal. This is especially damaging if the coupled noise is synchronous with the signal.

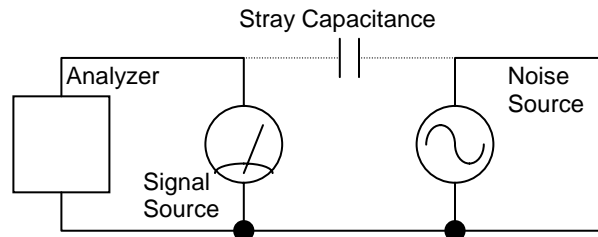


Figure Chapter 2 -3 Capacitive Coupling

We can estimate the noise current caused by a stray capacitance by,

$$i = C_{\text{stray}} dV/dt = \omega C_{\text{stray}} V_{\text{noise}}$$

where  $\omega$  is  $2\pi$  times the noise frequency,  $V_{\text{noise}}$  is the noise amplitude, and  $C_{\text{stray}}$  is the stray capacitance.

For example, if the noise source is a power circuit, then  $f = 60$  Hz and  $V_{\text{noise}} = 120$  V.  $C_{\text{stray}}$  can be estimated using a parallel plate equivalent capacitor. If the capacitance is roughly an area of  $1 \text{ cm}^2$  separated by  $10$  cm, then  $C_{\text{stray}}$  is  $0.009$  pF. The resulting noise current will be  $400$  pA (at  $60$  Hz). This small noise current can be larger than the signal

current. If the noise source is at a higher frequency, the coupled noise will be even greater.


If the noise source is at the signal frequency, then the problem is much worse. Vector averaging rejects noise at other frequencies, but pick-up at the signal frequency appears as signal!

Cures for capacitive noise coupling include:

- 1) Removing or turning off the noise source.
- 2) Keeping the noise source far from the signal source (reducing  $C_{\text{stray}}$ ). Do not bring the signal cables close to the noise source.
- 3) Designing the experiment to measure voltages with low impedance (noise current generates very little voltage).
- 4) Installing capacitive shielding by placing the signal source in a metal box.

### Inductive Coupling

An AC current in a nearby piece of apparatus can couple to the signal path via a magnetic field. A changing current in a nearby circuit gives rise to a changing magnetic field which induces an emf ( $d\Phi_B/dt$ ) in the loop connecting the signal source to the analyzer. This is like a transformer with the analyzer-source loop as the secondary winding.

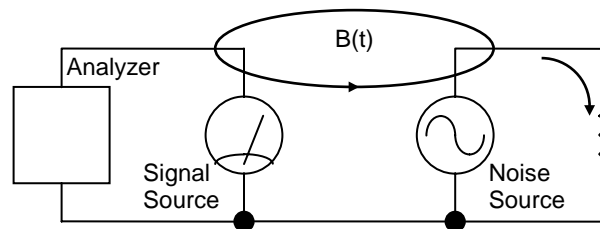

*Figure Chapter 2 -4 Inductive Coupling*

Cures for inductively coupled noise include:

- 1) Removing or turning off the interfering noise source.
- 2) Reduce the area of the pick-up loop by using twisted pairs or coaxial cables, or even twisting the 2 coaxial cables used in differential connections.
- 3) Using magnetic shielding to prevent the magnetic field from crossing the area of the signal path.
- 4) Measuring currents, not voltages, from high impedance sources.



### Resistive Coupling (Ground Loops)

Currents flowing through the ground connections can give rise to noise voltages. This is especially a problem with signal frequency ground currents.

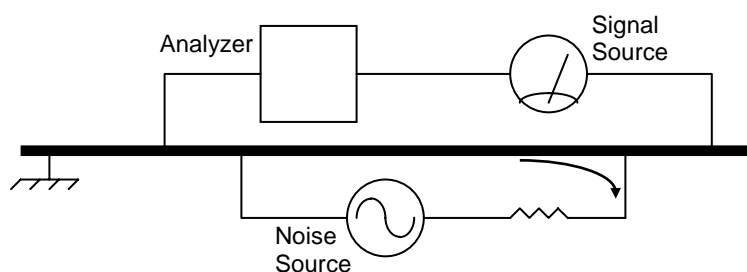


Figure Chapter 2 -5 Resistive Coupling

In this illustration, the analyzer is measuring the signal relative to a ground far from the signal source. The analyzer senses the signal plus the voltage due to the noise source's ground return current passing through the finite resistance of the ground between the analyzer and the source. The analyzer and the source are grounded at different places which, in this case, are at different potentials.

Cures for ground loop problems include:

- 1) Grounding everything to the same physical point.
- 2) Using a heavy ground bus to reduce the resistance of ground connections.
- 3) Removing sources of large ground currents from the ground bus used for small signals.
- 4) Measure the signal relative to the source ground using two cables (A-B). Set the Input Grounding to Float.

### Microphonics

Not all sources of noise are electrical in origin. Mechanical noise can be translated into electrical noise by microphonic effects. Physical changes in the signal source or cables (due to vibrations for example) can result in electrical noise over the entire frequency range of the analyzer.

For example, consider a coaxial cable connecting a signal source to the analyzer. The capacitance of the cable is a function of its geometry. Mechanical vibrations in the cable translate into a capacitance that varies in time, typically at the vibration frequency. Since the cable is governed by  $Q=CV$ , taking the derivative, we have

$$C \, dV/dt + V \, dC/dt = dQ/dt = i$$

Mechanical vibrations in the cable which cause a  $dC/dt$  will give rise to a current in the cable. This current affects the measured signal.

Some ways to minimize microphonic signals are:

- 1) Eliminate mechanical vibrations near the signal path.
- 2) Tie down cables carrying sensitive signals so they do not move.
- 3) Use a low noise cable that is designed to reduce microphonic effects.

### **Thermocouple Effects**

The emf created by junctions between dissimilar metals can give rise to many microvolts of slowly varying potentials. This source of noise is typically at very low frequency since the temperature generally changes slowly. This effect is large on the scale of many signals ( $10^3$ 's of  $\mu\text{V}$ ) and can be a problem for low frequency measurements, especially in the mHz range.

Some ways to minimize thermocouple effects are:

- 1) Hold the temperature of the signal source or detector constant.
- 2) Use a compensation junction, i.e. a second junction in reverse polarity which generates an emf to cancel the thermal potential of the first junction. This second junction should be held at the same temperature as the first junction.

# Curve Fitting and Synthesis

Often, the frequency response measurements made in the FFT and Swept Sine measurement groups need to be compared to a theoretical model of the device's behavior. This comparison is the job of Curve Fitting and Synthesis. In curve fitting, the analyzer extracts the best fit parameters of a linear frequency response function from a measured frequency response function. In Curve Synthesis, the analyzer transforms a set of frequency response parameters into an SR785 frequency response measurement which can be compared with the measured data.

## Curve Tables

Both curve fitting and curve synthesis use the SR785's two Curve Tables. The curve tables allow entry and editing of frequency response parameters in one of three formats: polynomial, pole-zero, and pole residue. Once parameters have been entered into the curve tables the corresponding frequency response function can be synthesized into a trace for comparison with measured data.

### Polynomial

In this format, the curve table represents a frequency response function as the ratio of two polynomials in the complex frequency variable  $s$ .

$$\text{Freq. Re sp.}(s) = \text{Gain} * \frac{(b_n s^n + b_{n-1} s^{n-1} + \dots + b_0)}{a_m s^m + a_{m-1} s^{m-1} + \dots + a_0}$$

The curve tables allow entry of both the numerator and denominator coefficients as well as the order of the numerator and denominator polynomials. The curve table also contains a constant gain factor which multiplies the polynomials.

### Pole-Zero

In the pole-zero format, the numerator and denominator polynomials are factored so that the frequency response curve is described by the ratio of the products of the poles and zeros: To ensure a real impulse response, all complex poles and zeros only occur in complex conjugate pairs.

$$\text{Freq. Re sp.}(s) = \text{Gain} * \frac{(s - z_n)(s - z_{n-1}) \dots (s - z_0)}{(s - p_m)(s - p_{m-1}) \dots (s - p_0)}$$

### Pole Residue

In the pole residue format a partial fraction expansion of the pole-zero form is performed to yield the frequency response as a sum of single pole terms weighted by residues.

$$\text{Freq. Re sp.}(s) = \text{Gain} * \left( \frac{R_m}{(s - p_m)} + \frac{R_{m-1}}{(s - p_{m-1})} + \dots + \frac{R_0}{(s - p_0)} \right)$$

Once again, the residues corresponding to complex conjugate pole pairs are complex conjugates themselves.

If the order of the numerator in the polynomial representation is equal to or greater than the order of the denominator polynomial, the Pole-Residue representation will contain an additional polynomial whose order is the difference between the numerator order and denominator order. This polynomial is added to the equation shown above. The coefficients of this extra polynomial cannot be adjusted by the user; use one of the other representations if this is necessary.

The SR785 will convert curve parameters between any of these three formats. Results of curve fits are always initially presented in pole-zero format although they can always be subsequently converted to any format.

### **Frequency Scale**

Often a theoretical model frequency response function will be specified with respect to a nominal frequency, such as 1 Hz. The curve table contains a frequency scale parameter to allow such a nominal transfer function to be synthesized to any actual frequency. In terms of the pole zero representation, the poles and zeros are all multiplied by the scaling factor before synthesizing the curve. When fitting, after the poles and zeros corresponding to the measured frequency response function are determined, the values are all divided by the Frequency Scale factor. For example, if the curve table contains a pole (s-1), and the frequency scale is 1000, this will create a pole at 1 kHz when synthesizing the table.

### **Delay**

The three formats for a frequency response function discussed above cannot model a pure delay term of the form  $\exp(-sT)$ , where T is the delay time. Because of this, and because measured frequency response data often includes delay, a delay parameter is included in the curve table. On synthesis the analyzer multiplies the Frequency Response function by the delay factor  $\exp(-sT)$ . When fitting, the input data is divided by  $\exp(-sT)$  before the curve parameters are calculated. Note that the curve fitter will not fit an appropriate value of the delay- this value must be supplied by the user.

### **Trace**

The curve table also contains a field which determines which trace will be written to when the curve is synthesized.

## **Curve Fitting**

Before fitting a measured frequency response function several setup parameters must be specified. The fit region determines the region of the active display that will be fit. The number of poles and zeros determine the order of the frequency response function which will be fit to the measured data. The active display must be in the FFT or Swept Sine measurement group to be fit.

When the user presses Start Fit, the analyzer computes in a single pass, the best fit parameters for the active display using the Delay and Frequency Scale parameters for Curve Table 1. These parameters are then placed in curve table 1 in the pole-zero representation. Finally, Curve Table 1 is synthesized and the synthesized trace is placed in the non-active display so it can be compared to the original data in the active display.

**Weighting**

Because the quality of the measured frequency response function is not uniform over the entire region to be fit, it is sometimes necessary to weight the input data when performing the fit.

The analyzer can weight the fit in one of three ways. Uniform weighting weights all parts of the fit region uniformly. Coherence weighting uses the calculated coherence function (or normalized variance in the case of swept sine measurements) to weight the data. Coherence is a measure of the signal to noise ratio of the measured frequency response function. In regions of high signal to noise, the value of the coherence function will be near 1. In regions of low signal to noise, the coherence function will be less than 1. Coherence weighting is not available for Off-Line displays.

User weighting allows the user to specify the weighting function. The analyzer uses whatever data is found in the real part of the specified weighting trace to determine the weighting.

**Curve Synthesis**

The curve table is synthesized when the user presses the Synthesize Table 1 or Synthesize Table 2 softkeys. The analyzer first looks at the active display to determine the frequency range over which the synthesis occurs. Then the SR785 calculates the complex frequency response function corresponding the curve parameters, frequency scale, delay and gain settings contained in the curve table. These values are placed in the trace corresponding to the synthesized curve table. Finally, the active display is taken Off-Line, and the synthesized trace is displayed on the active display.

Synthesis is only allowed when the active display is in the FFT or Swept Sine measurement group.



## Chapter 3

# Operation

## In This Chapter

<b>Overview</b>	<b>3-3</b>		
Power Switch	3-3		
Hardware Reset	3-3		
Software Reset	3-4		
Video Display	3-4		
Keypad	3-4		
Softkeys	3-4		
Knob	3-5		
Disk Drive	3-5		
<b>Front Panel Connectors</b>	<b>3-6</b>		
Ch1 Signal Inputs	3-6		
Ch2 Signal Inputs	3-6		
Trigger Input	3-6		
Source Output	3-6		
PC Keyboard Connector	3-6		
<b>Rear Panel Connectors</b>	<b>3-8</b>		
Power Entry Module	3-8		
IEEE-488 Connector	3-9		
Serial RS232 Connector	3-9		
Parallel Printer Connector	3-9		
Preamp Connector	3-9		
Using SRS Preamps	3-10		
Tachometer Input	3-10		
Start Input	3-10		
<b>Screen Display</b>	<b>3-11</b>		
Displays	3-11		
Data Graph	3-12		
Vertical Scale Bar	3-13		
Horizontal Scale Bar	3-14		
Marker	3-15		
Marker Position Bar	3-15		
Active Display	3-16		
Reference Graphs	3-16		
Menu Display	3-16		
<b>Status Indicators</b>	<b>3-18</b>		
Input Ranges	3-18		
Overloads	3-18		
Input Configurations	3-18		
		AWt	3-18
		ArmWait/Trig-wait/Trig/Acquire	3-18
		Run/Pause/Done	3-19
		Analog/Playback	3-19
		Real Time	3-19
		Sweep Frequency	3-20
		RPM	3-20
		Capture Progress	3-20
		No Cap./Cap. Data	3-20
		Link	3-20
		Enter	3-21
		ALT	3-21
		GPIB/RS232	3-21
		Comm /Error	3-21
		SRQ	3-21
		Local/Remote/LOCK	3-22
		macro	3-22
		<b>Keypad</b>	<b>3-23</b>
		<b>Normal and Alternate Keys</b>	<b>3-23</b>
		<b>Menu Keys</b>	<b>3-24</b>
		<b>Entry Keys</b>	<b>3-25</b>
		Buttons	3-25
		Lists	3-25
		Numeric Values	3-26
		<b>Control Keys</b>	<b>3-28</b>
		[Start/Reset]	3-28
		[Pause/Cont]	3-28
		[Start Capture]	3-28
		[Stop Capture]	3-29
		[Active Display]	3-29
		[Link]	3-30
		[Print Screen]	3-31
		[Help/Local]	3-31
		[Alt] [Control Key]	3-31
		Play Macro ([Alt] [Start/Reset])	3-31
		Macro Rec ([Alt] [Pause/Cont])	3-31
		End Rec ([Alt] [Stop Capture])	3-32
		Snap Ref ([Alt] [Start Capture])	3-32

## 3-2 Operation

---

Ref to Trace ([Alt] [Active Display])	3-32
Trace to Ref ([Alt] [Link])	3-32
Display to Trace ([Alt] [Print Screen])	3-32
Trace to Display ([Alt] [Help/Local])	3-33

### **Function Keys** **3-34**

[Auto Scale A]	3-34
[Auto Scale B]	3-34
[Auto Range Ch1]	3-34
[Auto Range Ch2]	3-35
[Span Up]	3-35
[Span Down]	3-35
[Marker Ref]	3-36
[Display Ref]	3-36
[Marker Center]	3-37
[Marker Max]	3-37
[Marker Min]	3-37
[Show Setup]	3-37

### **Macros** **3-38**

Keypad Macros	3-38
Choosing From Lists	3-38
Menus	3-38



# Overview

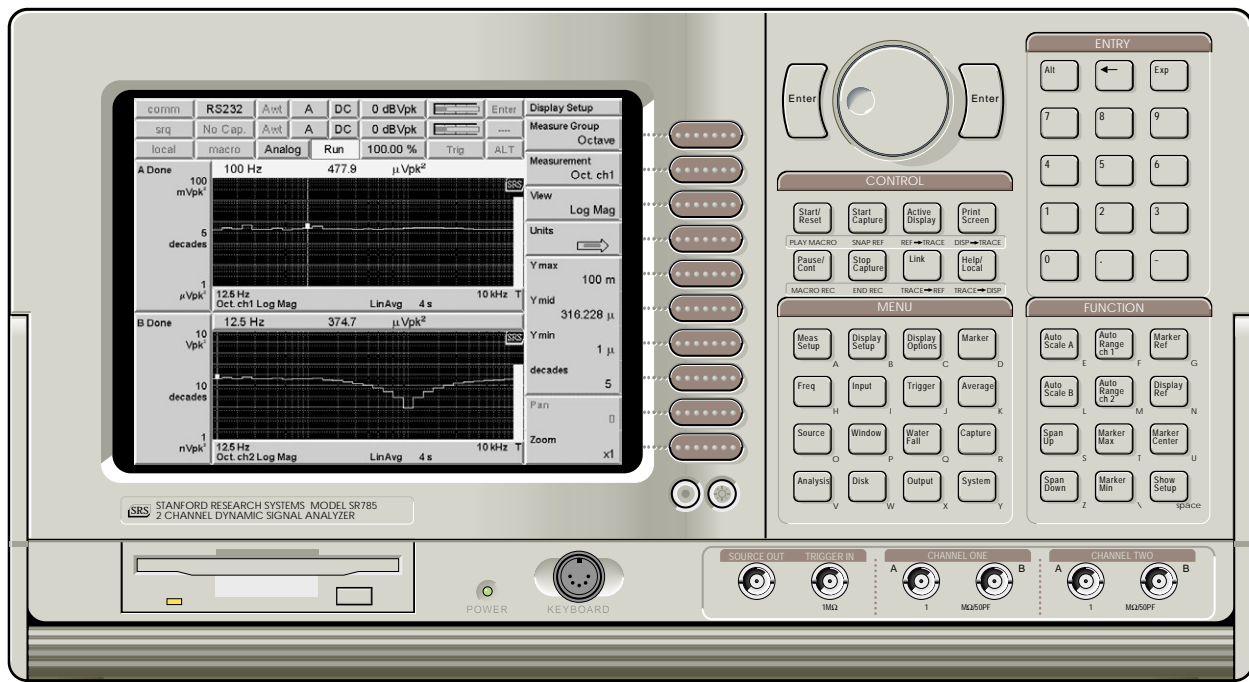


Figure Chapter 3 -1 Front Panel

## Power Switch

The power switch is located on the rear panel. The SR785 is turned on by depressing the upper half of the power switch. The green power LED on the front panel indicates that the unit is powered.

The video display may take a few seconds to warm up and become visible. The brightness is adjusted using the [Brighter] and [Dimmer] buttons below the softkeys. The firmware revision, serial number, memory size and calibration dates of the unit are displayed when the power is turned on.

This power on screen may be displayed at any time using [System]<Show Version>.



## Caution

This instrument may be damaged if operated with the LINE VOLTAGE SELECTOR set for the wrong AC line voltage or if the wrong fuse is installed.

## Hardware Reset

Turn the power on while holding down the backspace [←] key and continue to hold backspace for at least 3 seconds to reset the unit. The unit will perform power on tests and assume the default settings (including the default remote interface settings).

### **Software Reset**

<Preset> in the [System] menu resets the instrument to the default settings. Only the remote interface settings are not changed. All stored data are lost.

This function requires pressing the [Enter] key to confirm.

### **Video Display**

The monochrome video display is the user interface for data display and front panel programming operations. The resolution of the display is 800H by 600V.

The brightness is adjusted using the [Brighter] and [Dimmer] buttons below the softkeys. The contrast is adjusted using [Alt] [Brighter] and [Alt] [Dimmer]. As with most video displays, do not set the brightness higher than necessary. To extend the life of the display, use the Screen Saver in the [System] <Preferences> menu.

A complete description of the screen display follows later in this chapter.

### **Keypad**

The keypad consists of four groups of hardkeys (keys with printed labels).

Hardkeys are referenced in braces like [Active Display] or [Input].

The ENTRY keys are used to enter numeric parameters which have been highlighted by a soft key. The MENU keys select a menu of soft keys. Pressing a menu key will change the menu boxes which are displayed next to the soft keys. Each menu presents a group of similar or related parameters and functions. The CONTROL keys start and stop data acquisition, toggle the active display and link parameters and functions. These keys are not in a menu since they are used frequently and within any menu. The FUNCTION keys perform common functions such as Auto Scale and Auto Range. These keys can be accessed at any time.

A complete description of the keys follows.

### **Softkeys**

The SR785 has a menu driven user interface. The MENU keys each display a menu of softkeys. The softkeys are at the right of the video display and have different functions depending upon the displayed menu.

Softkeys are referenced in brackets like <Span> or <Units>.

There are three types of softkeys - buttons, lists and numeric values. A button performs a function, such as <Full Span>. A list presents a list of choices or options in the entry field (at the top of the screen). Use the knob to make a selection and press [Enter]. <Measurement> is an example of a list. A numeric value presents the current value in the entry field and awaits numeric entry. Enter a new value with the ENTRY keys and press [Enter]. <Start Freq> is an example of a numeric value.

A complete description of the menu softkeys follows in the next chapter.

**Knob**

The knob normally moves the markers within the displays. If a parameter has been highlighted by its softkey, the knob adjusts the parameter. Parameters which have a list of choices are most easily modified with the knob. Numeric entry fields may also be adjusted with the knob.

**Disk Drive**

The 3.5" disk drive is used to store data and instrument settings. Double sided, high density disks (DS/HD) have a capacity of 1.44M bytes and double sided, double density disks (DS/DD) have a capacity of 720k bytes. The disk format is DOS compatible.

Use the [Disk] menu to format a disk or access disk files.

## Front Panel Connectors

### Ch1 Signal Inputs

The Channel 1 A and B inputs are voltage inputs with 1 M $\Omega$ , 50 pF input impedance. The Input Coupling is either DC or AC. The Input Mode may be single-ended (A) or differential (A-B). The A and B connector shields are common and grounded to the chassis by 1 M $\Omega$  + 0.01  $\mu$ F (Float) or 50  $\Omega$  (Ground). The shields should never exceed  $\pm 4$ V.

Do not apply more than  $\pm 50$  V to either input. If the input exceeds  $\pm 57$  V, the input range will be set to 34 dBV to protect the input from damage and the HighV overload status is set. The input range can not be changed while HighV is detected.

For two channel measurements such as Frequency Response, Channel 1 is the reference input to the device under test.

### Ch2 Signal Inputs

The Channel 2 A and B inputs are voltage inputs with 1 M $\Omega$ , 50 pF input impedance. The Input Coupling is either DC or AC. The Input Mode may be single-ended (A) or differential (A-B). The A and B connector shields are common and grounded to the chassis by 1 M $\Omega$  + 0.01  $\mu$ F (Float) or 50  $\Omega$  (Ground). The shields should never exceed  $\pm 4$ V.

Do not apply more than  $\pm 50$  V to either input. If the input exceeds  $\pm 57$  V, the input range will be set to 34 dBV to protect the input from damage and the HighV overload status is set. The input range can not be changed while HighV is detected.

For two channel measurements such as Frequency Response, Channel 2 is the response output from the device under test.

### Trigger Input

When the Trigger Source is External or Ext TTL, the Trigger input triggers the measurement. The input impedance is 1 M $\Omega$  and the minimum pulse width is 25 ns. The Trigger Level is adjustable from -5V to +5V with either Positive or Negative Slope. The minimum trigger amplitude is 100 mV. The trigger input should not exceed  $\pm 5$ V.

### Source Output

The Source can output either Sine, Two Tone, Chirp, Noise or Arbitrary waveforms. The output impedance is less than 5  $\Omega$  and is capable of driving a 50  $\Omega$  load. The output is ground referenced.

### PC Keyboard Connector

An IBM PC or XT compatible keyboard may be attached to the keyboard connector. Most keyboards have a switch on the back to select PC/XT or 8088 mode. The SR785 can be controlled from this keyboard according to the table below.

<u>SR785 keypad</u>	<u>PC Keyboard equivalent</u>
Softkeys	F1-F10
Menu Keys	Letter shown below the Menu key. (e.g. A to show the Frequency Menu).
Function Keys	Letter shown below the Function key. (e.g. E for AutoScale A).
Control Keys	Alt-1 through Alt-4 (top row) Alt-5 through Alt-8 (bottom row)
Number Keys	0 through 9
Knob (slow)	Left and Right arrow
Knob (fast)	Up and Down arrow
Alt	Escape
Enter	Enter or Return

The keyboard should only be connected when the power is off.

# Rear Panel Connectors

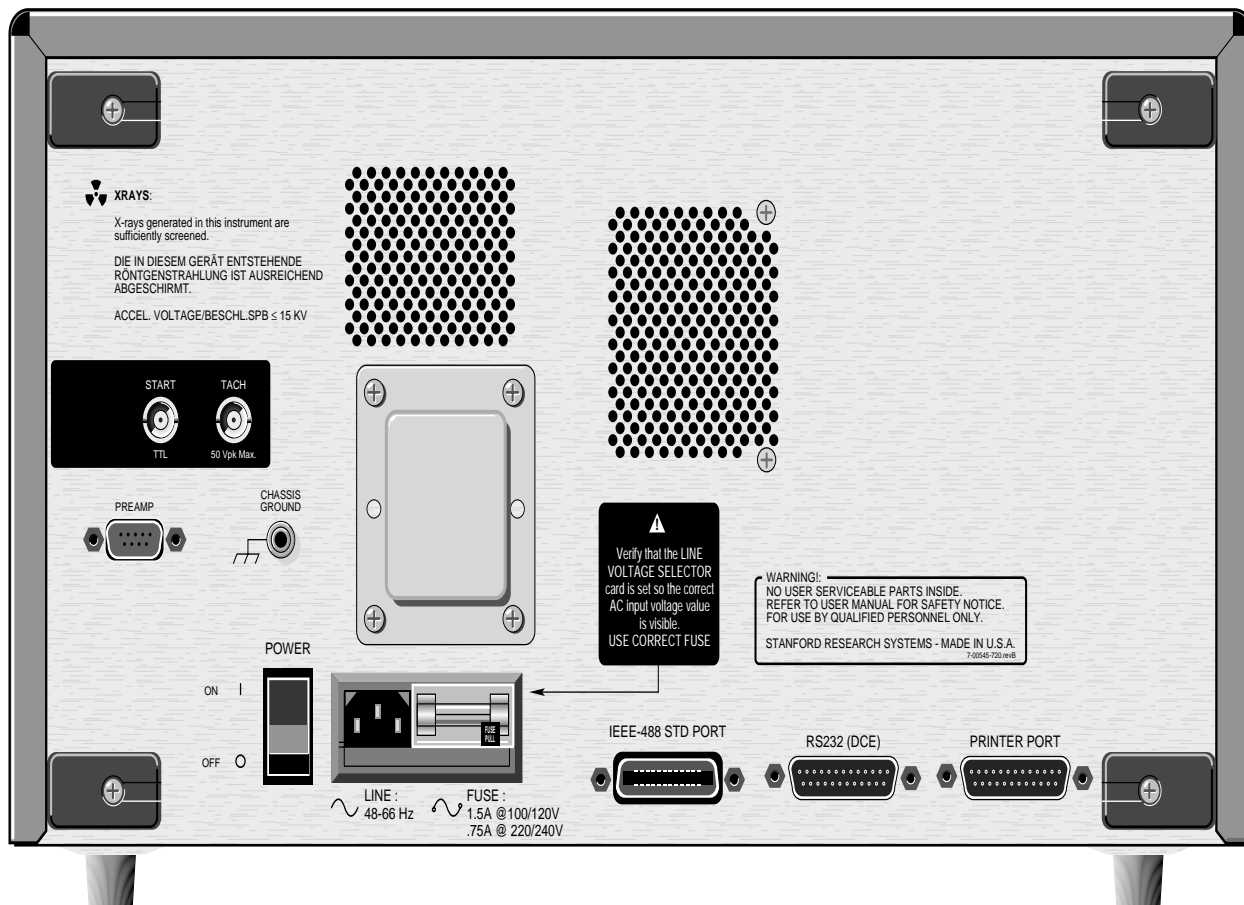


Figure Chapter 3 -2 Rear Panel

## Power Entry Module

The power entry module is used to fuse the AC line, select the line voltage, and block high frequency noise from entering or exiting the instrument. The SR785 has a detachable, three-wire power cord for connection to the power source and to protective ground. The exposed metal parts of the instrument are connected to the outlet ground to protect against electrical shock. Always use an outlet which has a properly connected protective ground.



### Caution

This instrument may be damaged if operated with the LINE VOLTAGE SELECTOR set for the wrong AC line voltage or if the wrong fuse is installed.

Do not attempt to service or adjust this instrument while it is plugged into a live outlet.

### Line Voltage Selection

The SR785 operates from a 100V, 120V, 220V, or 240V nominal AC power source having a line frequency of 50 or 60 Hz. Before connecting the power cord to a power

source, verify that the LINE VOLTAGE SELECTOR card, located in the rear panel fuse holder, is set so that the correct AC input voltage value is visible.

Conversion to other AC input voltages requires a change in the fuse holder voltage card position and fuse value. Disconnect the power cord, open the fuse holder cover door and rotate the fuse-pull lever to remove the fuse. Remove the small printed circuit board and select the operating voltage by orienting it so that the desired voltage is visible when pushed firmly back into its slot. Rotate the fuse-pull lever back into its normal position and insert the correct fuse into the fuse holder.

### Line Fuse

Verify that the correct line fuse is installed before connecting the line cord. For 100V/120V, use a 1.5 Amp fuse. For 220V/240V, use a 3/4 Amp fuse.

### IEEE-488 Connector

The 24 pin IEEE-488 connector allows a host computer to control the SR785 via the IEEE-488 (GPIB) instrument bus. The GPIB Address of the unit is set in the [System] <Remote> menu.

A GPIB plotter with HPGL compatible graphics may be connected to the IEEE-488 port. In this case, the SR785 will control the plotter to generate plots of the display graphs. Select GPIB as the Output Destination and set the Plotter Address and Type in the [Output] menu.

### Serial RS232 Connector

The RS232 interface connector is configured as a DCE (transmit on pin 3, receive on pin 2). The Baud Rate, Parity, and Word Length are set in the [Setup] <Remote> menu. To connect the SR785 to a PC serial adapter, which is usually a DTE, use a straight through serial cable.

A serial plotter or printer may be connected to the RS232 port. Select the Serial interface as the Output Destination and choose the Printer or Plotter Type in the [Output] menu.

### Parallel Printer Connector

Attach a parallel printer to the Printer port. Use a standard printer cable. Select the Centronics printer interface as the Output Destination and choose the Printer Type in the [Output] menu.

### Preamp Connector

This 9 pin 'D' connector provides power and control signals to external preamplifiers such as the SR550 and SR552. The power connections are described below.

<u>Pin</u>	<u>Voltage</u>
1	+20V
2	+5V
6	-20V
7	Signal Ground
8	Ground

### Using SRS Preamps

When using either the SR550 or SR552, connect the power cable (standard 9 pin D connectors) from the preamp to the rear panel Preamp connector. Use BNC cables to connect the A output from the preamp to the A input of either channel on the SR785. The B output from the preamp (preamp ground) may be connected to the B input on the SR785. In this case, use A-B as the Input Mode. Be sure to twist the A and B cables together to reduce noise pick-up.



**The SR550 and SR552 are AC coupled above 1 Hz.** Do not use either preamp to measure signals below 1 Hz!

**The SR785 does not compensate for the gain of the preamp.** Both preamps operate at their highest gain. Measurements made by the SR785 need to be divided by the gain of the preamp. The SR550 has a gain of 10 and the SR552 has a gain of 100.

### Tachometer Input

Attach a tachometer to the rear panel tachometer input to synchronize the SR785 to a rotating machine for order tracking measurements or to trigger a measurement using RPM arming. Parameters relating to the tachometer input are set in the [Input]<Tach Input> Submenu.

### Start Input

This TTL input provides a means of remotely starting a measurement. A TTL rising edge at this input is equivalent to pressing the [Start/Reset]



# Screen Display

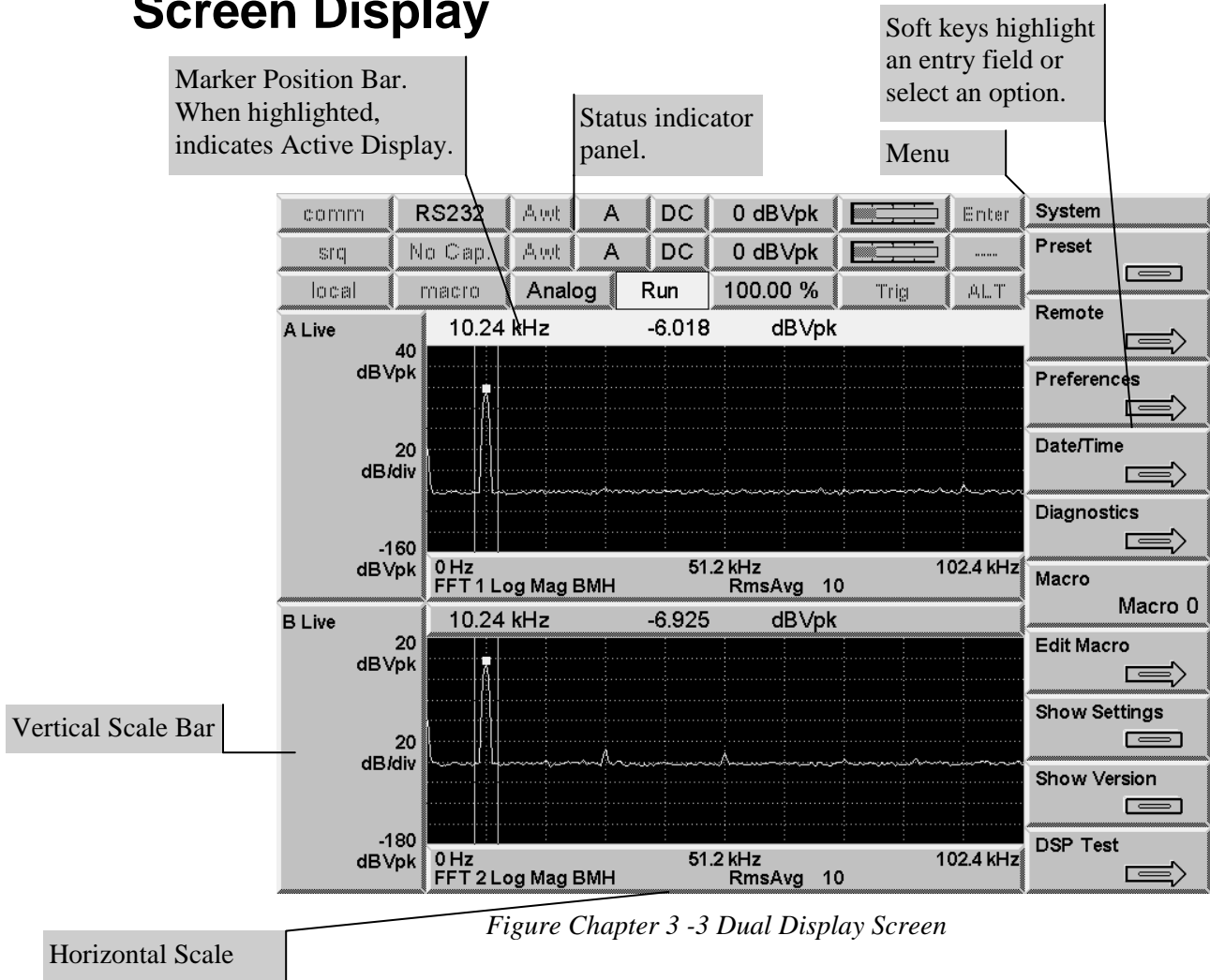


Figure Chapter 3 -3 Dual Display Screen

## Displays

There are two displays, labeled Display A and Display B. Each display consists of a data graph, vertical scale bar, horizontal scale bar and marker position bar.

Each display has an associated measurement. The displays may have different measurements, such as spectrum and time record, different views, such as magnitude and phase, and so on. Many instrument parameters are set independently for each display.

The display format is selected in the [Display Options] menu.

The display shown above is the Dual display format. The [Active Display] key toggles between the two displays. Display A is always on top.

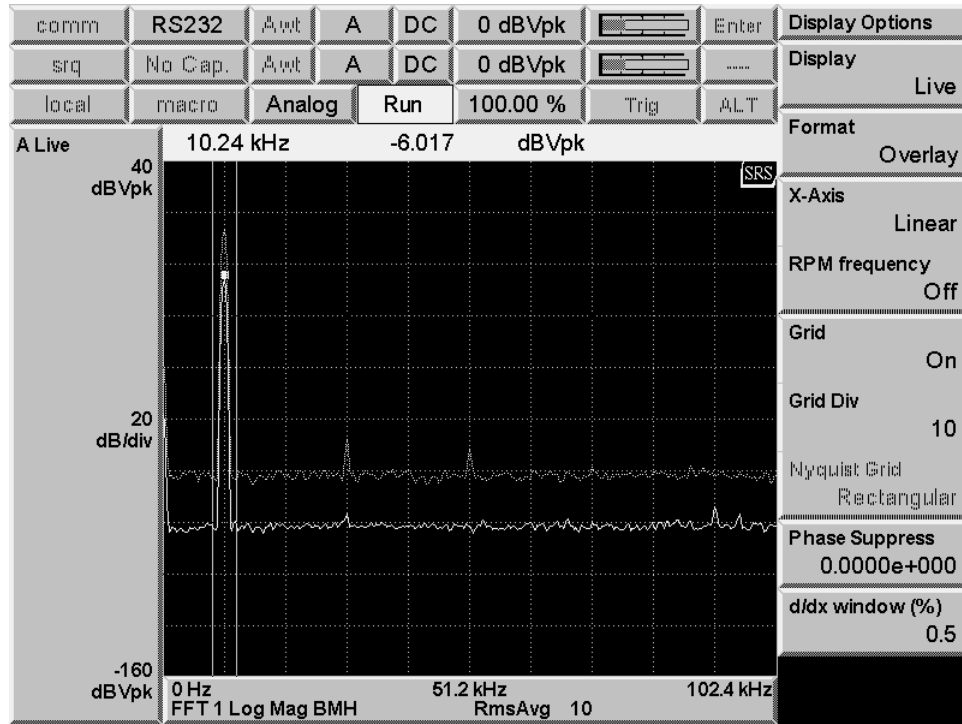


Figure Chapter 3 -4 Overlay Display Screen

The Overlay display format is shown above. The measurement of the inactive display is shown in half-intensity. The Single display format is similar but the inactive display is not shown at all.

### Data Graph

Data is normally graphed with signal amplitude on the Y axis and frequency or time on the X axis. The physical size of the graph remains constant while the vertical and horizontal scales may be changed. The graph area has a dotted grid for reference. There are 10 horizontal divisions and either 8, 10, 12 or 15 vertical divisions. The horizontal span is determined by the measurement span or time. The display normally shows the entire measurement span. The X axis may be expanded (Zoom) and translated (Pan) to display less than the full measurement span. This expansion does not change the measurement span or time, it merely changes the display of the data.

### Half Intensity Data Graph


Unsettled measurement data is graphed in half intensity. Once the measurement becomes settled, the graph returns to full intensity.

FFT and Octave measurements are unsettled when the frequency span or resolution is changed and when the input signal path is changed (gain, filters, etc.). The measurement is not settled again until the discontinuity in the input data has propagated through the digital filters and, for FFT's, a complete new time record has been acquired.

FFT measurements running with a Time Record Increment less than 100% (overlapped time records) display unsettled measurements in half intensity.

Octave measurements show unsettled octave bins in half intensity. Each octave bin is displayed in half intensity until it is settled. The settling time of each bin is related to the  $1/(\text{center frequency})$  and the lowest bins take the longest time to settle.

### Vertical Scale Bar

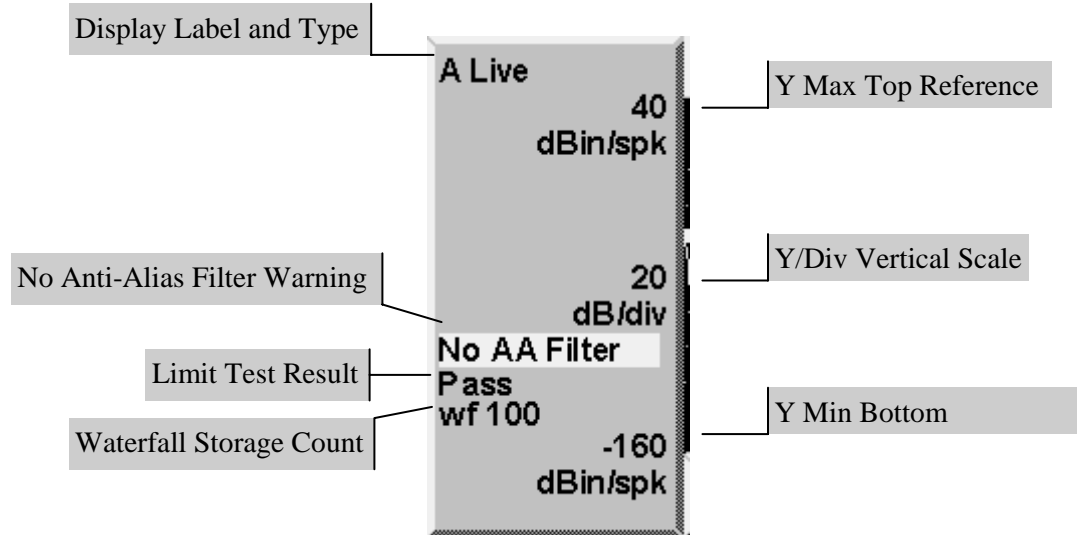

Chapter 3 -5 Vertical Scale Bar

The Vertical Scale Bar is at the left of the data graph for each display.

The display name, either 'A' or 'B' is at the top left corner. The display type (Live, Done, Off-Line or Preview) is next to the display name. 'Live' indicates that the display is showing live measurement results from either the inputs or from capture. 'Done' indicates that the display is showing a completed measurement. An 'Off-Line' display is showing a 'frozen' measurement. A display is Off-Line whenever data is recalled to the display from a Trace or from disk. 'Preview' only occurs when Average Preview is On and indicates that the display is showing a preview time record and is waiting for accept or reject.


The Y Max and Y Min reference values are shown along the right edge. Y Max is the top reference and Y Min is the bottom. The vertical scale per division is shown between Y Max and Y Min. The units and scale/division are selected in the [Display Setup] menu.

The No Anti-Aliasing Filter warning is shown if any measurement input has its anti-aliasing filter off. Measurements made without the anti-aliasing filter may contain alias signals with frequencies above 102.4 kHz.

The Limit Test Result, either 'Pass' or 'Fail' is displayed if limit testing is on for this display.

The Waterfall Storage Count shows how many measurement records are currently stored in the waterfall buffer. This count is not displayed when Waterfall Storage is Off.

## Horizontal Scale Bar

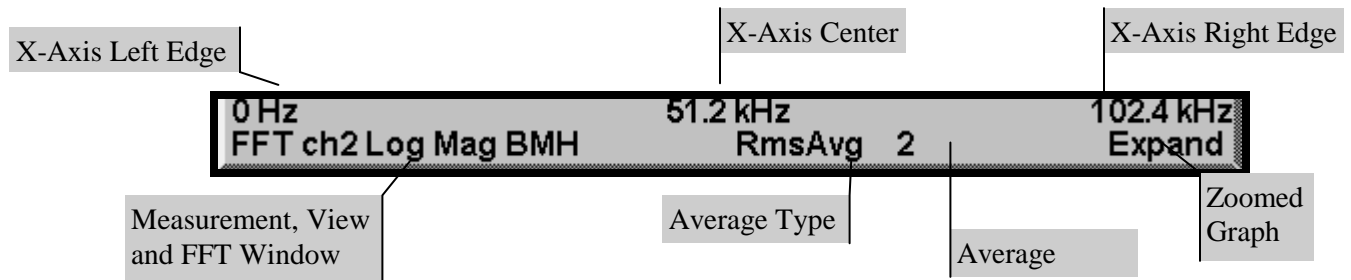

*Figure Chapter 3 -6 Horizontal Scale Bar*

The Horizontal Scale Bar is directly below the data graph.

When displaying spectral measurements with no horizontal expansion, these values are the Start, Center and Stop frequencies of the measurement. The start and stop frequencies can be displayed in RPM instead of Hz using the RPM frequency softkey. When displaying time records, these values are the Start and End of the time record.

Below the X-Axis labels, the measurement, view, window and averaging are shown. These parameters may be different for each display.

The Average Count shows the number of measurements averaged so far (up to the programmed Number of Averages). 'Settle' is shown when the displayed measurement is not settled and is not added to the average. 'Reject' is shown when the Overload Reject is On and a measurement is not included in the average due to input overload.

For Swept Sine measurements, the estimated Sweep Time is displayed instead of the Average Count. This is the time to complete a sweep when no auto functions are selected (no Source Auto Level, no Auto Range, no Auto Resolution). This Sweep Time is simply the sum of the Settle and Integrate times for all points in the sweep.

To horizontally expand a graph, use <Pan> and <Zoom> in the [Display Setup] menu. When the display is expanded on the horizontal axis, the graph labels reflect the displayed data, not the actual measurement span or time. Expanded displays show the 'Expand' indicator.

## Marker

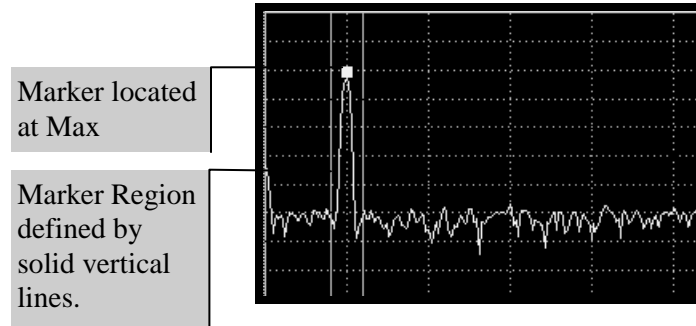


Figure Chapter 3 -7 Marker Region

The Marker Region is the graph region between the two solid vertical marker lines. The Marker Region Width may be set to 1 division (Wide), 1/2 division (Normal), or a single vertical line (Spot). The Marker Region is associated with the display and does not change with graph scaling. The Marker is the small square within the Marker Region. The Marker finds the Minimum, Maximum, or Mean of the data in the Marker Region. When seeking Min or Max, the Marker is located at the position of the minimum or maximum. This allows peaks and valleys in the data to be found easily. When seeking the Mean, the X position of the Marker is at the center of the Marker Region and the Y position is the mean of the data within the region. When a Spot Marker Region is used, the Marker is confined to a single X location. Use the [Marker] menu to configure the Marker for each display.

The Marker position is shown in the Marker Position Bar above the graph.

Use the knob to move the marker in the active display. The knob always moves the marker unless a menu softkey has been pressed and an entry field is displayed at the top of the screen.

To move the marker to the next "peak" in the display, press [Enter] and move the knob. The marker will move to the next peak in the display in the direction of the knob turn. If no peaks are found, the marker will remain in the same position.

## Marker Position Bar

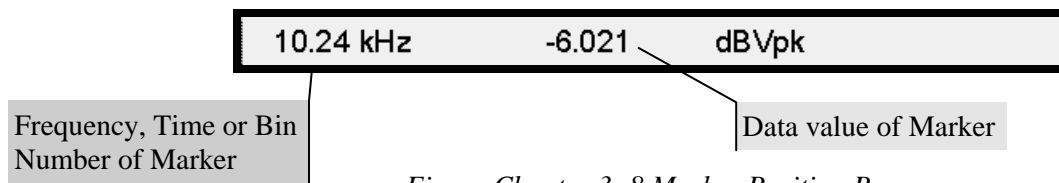


Figure Chapter 3 -8 Marker Position Bar

The Marker Position Bar is displayed above each graph.

The first value is the marker frequency, time or bin number. If the display is a waterfall, the waterfall record number is also shown.

The second value is the marker Y Position or data value. For 2-D displays (Nyquist or Nichols View), X and Y values are both displayed.

If a value is preceded by a  $\Delta$  (delta) symbol, then the value is relative to a reference point or graph. Use the [Marker] menu to configure a relative reading.

If 'N/A' appears instead of a value, measurement data is not available. This often occurs when a measurement with a long acquisition time is changed and new data is not yet available.

If a '?' appears after a value, the displayed value has not actually been measured but is interpolated from actual measured points. This arises in Swept Sine sweeps with Auto Resolution On when some frequency points are not measured.

### Active Display

The active display is indicated by the highlighted Marker Position bar. Use [Active Display] to toggle the active display.

### Reference Graphs

In addition to the data graph, each display can also show a reference graph. The reference graph is a graph of a stored Reference Display and is drawn in half intensity. This allows comparison of live measurement data with a stored Reference Display.

The Reference Display data is copied from the current display data using [Display Ref] or loaded from a data Trace with [Trace to Ref].

Set the Marker Relative to Reference Display (in the [Marker] menu) to read the Marker Position relative to the Reference Display. Press [Display Ref] again to turn off the Reference Display.

Changing the display scale does NOT rescale the Reference Display. This allows the live measurement display to be offset from the Reference Display by changing either Ymax, Ymid or Ymin in the [Display Setup] menu. To graph the Reference Display in the new display scale, press [Alt] [Snap Ref].

### Menu Display

The softkey menu boxes define the functions of the softkeys to the right of the screen. Related functions and parameters are grouped into menus. Pressing each of the Menu keys displays a different set of menu boxes. There are three types of softkeys - buttons, lists and numeric values. A button performs a function, such as <Full Span>. A list presents a list of choices or options in the entry field (at the top of the screen). Use the knob to make a selection and press [Enter]. <Measurement> is an example of a list. A numeric value presents the current value in the entry field and awaits numeric entry. Enter a new value with the entry keys and press [Enter]. <Start Freq> is an example of a numeric value.

Menus which show parameters specific to a single display show the parameters for the active display. Use [Active Display] to change the active display.

A softkey menu box which is shown in gray is not available in the current measurement. Also, the measurement softkeys for an 'Off-Line' display are shown in gray indicating that they may not be changed.

Each menu is described at length in Chapter 4.

## Status Indicators

In addition to the data displays and menu boxes, there are a number of status indicators which are displayed at the top of the screen.

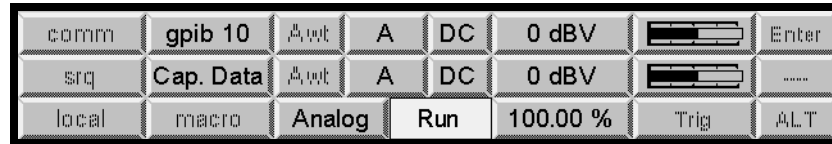


Figure Chapter 3 -9 Status Indicator Panel

### Input Ranges

The Input Ranges for both inputs are displayed. If the range is set manually, the range indicator is on a normal background. If Auto Range is On, then the range is displayed on a highlighted background. The input range is displayed in whatever units it was last entered in with the menu.

### Overloads

The input overload condition of the inputs are displayed to the right of their Input Ranges. The overload indicator is a bar. When the bar is 1 division long (displayed faintly), then the input signal is below half full scale. When the bar is 2 divisions long (displayed in black), the input signal is above half full scale (but not overloaded). When the bar is 4 divisions long (displayed brightly), the input signal exceeds full scale and is overloaded. This also sounds the audible overload alarm.

Once an overload occurs, the outline of the overload bar remains in the indicator. The overload outline indicates that the input has overloaded since the measurement was started. This is important when the measurement is linearly averaged and an overload occurs during the average. Pressing [Start/Reset] restarts the measurement and erases this overload outline.

Whenever an input greater than 57 V is detected at the input, the overload indicator displays 'HighV' on a bright background. The analyzer immediately switches to 34 dBV range to protect the input circuitry. Any attempt to change the input will not be allowed until the input signal is reduced to a safe level.

### Input Configurations

The Input Mode and Coupling for both inputs are displayed. If the coupling is DC or AC, the coupling indicator is on a normal background. If ICP coupling is selected, then ICP is displayed on a highlighted background.

### AWt

The A-Weighting filter status for both inputs are displayed. If the A-Weighting filter is ON, then 'AWt' is displayed on a highlighted background.

### ArmWait/Trig-wait/Trig/Acquire

This indicator reads 'Trig-wait' while waiting for a trigger. Once triggered, the indicator flashes 'Trig' and then displays 'Acquire' while the time record is being acquired. If the



Trigger Mode is Auto Arm, the indicator returns to 'Trig-wait' after the time record is acquired. If the Trigger Mode is Manual Arm, then the indicator reads 'Arm-wait' while waiting for a manual arm command (or keypress).

### **Run/Pause/Done**

'Run' indicates that the measurements of BOTH displays are running. Whenever new data is available for a Live display, the display will update. Off-Line displays do not update.

'Pause' indicates that the measurements of BOTH displays are paused.

'Done' indicates that single shot capture playback is finished.

Use the [Start/Reset] and [Pause/Cont] keys to control the measurements.

### **Analog/Playback**

'Analog' indicates that the measurements of BOTH displays are taking their inputs from the Channel 1 Input or Channel 2 Input.

'Playback' indicates that the measurements of BOTH displays are taking their inputs from the Capture buffers.

### **Real Time**

When the Measurement Group is FFT or Octave, the Real Time indicator shows the actual time record increment for the measurements in progress.

If the displayed increment is 100%, the start of the next time record is exactly one time record advanced from the start of the previous time record. If the increment is 25%, then the next time record starts 1/4 of a time record advanced from the start of the previous time record. This is sometimes referred to as 75% overlap since the two time records share 75% of a record. The overlap is simply 100% minus the Time Record Increment.

If the displayed increment is 200%, the start of the next time record is advanced from the start of the previous time record by 2 time records leaving a gap of 1 record. This means that the data between the two time records is not measured.

When the displayed increment is less than or equal to 100%, the measurement is said to be 'real time'. All time points contribute to one (or more) measurements.

When the displayed increment is greater than 100%, then the measurement is not 'real time' and some time points do not contribute to a measurement.

If the indicator shows a value greater than the requested Time Record Increment, it means that the measurement cannot be made with the requested increment but is running with the smallest increment possible. Factors which affect the processor's ability to run real time include the measurement type, averaging and source type.

See also Real Time Bandwidth and Overlap.

This indicator is replaced by the Capture/Playback Progress indicator when capture (playback) is in progress.

This indicator is replaced by the Sweep Frequency indicator when the Measurement Group is Swept Sine.

This indicator is replaced by the RPM indicator when Show Tach is set to On in the [input]<Tach Input> submenu.

### **Sweep Frequency**

The Sweep Frequency indicator shows the current frequency during Swept Sine sweeps. 'Ramping' is displayed if the source level is currently ramping. This indicator replaces the Realtime indicator when the Measurement Group is Swept Sine.

### **RPM**

The RPM indicator displays the rotational speed of the device under test based on information from the tachometer input. The rotational speed is calculated from the frequency of pulses arriving at the tachometer input and the setting of the Pulses per Revolution softkey. This indicator replaces the Real time indicator when the Measurement Group is Order or when Show Tach is set to On.

### **Capture Progress**

During capture fill, the Capture Progress indicator shows how much of the capture buffer has been filled. During capture playback, this indicator shows the playback progress through the playback portion of the buffer. This indicator replaces the Real time indicator when capture or playback is in progress.

### **No Cap./Cap. Data**

The Capture indicator shows the status of the Capture buffer.

'No Cap.' indicates that the Capture buffer is empty. Input playback from Capture is not allowed. Arbitrary source playback from Capture is also not allowed.

'Capturing' indicates that capture is in progress. The capture progress is shown in the Real Time indicator position.

'Cap. Data' indicates that the Capture buffer contains data for playback.

### **Link**

When an unlinked display parameter is being entered or modified, the Link indicator (next to the top of the menu) reads either 'DispA' or 'DispB', whichever is the active display. Parameter entry or selection modifies only the active display. To change which display will be modified, press [Active Display].

When a linked parameter is highlighted, the Link indicator reads 'Link'. Parameter entry or selection modifies both displays at once.

To link or unlink a parameter, highlight the parameter softkey and press the [Link] key. If the Link indicator is shown in gray, then the parameter linking may not be changed.

The availability of linking and unlinking is affected by the setting of the Analyzer Configuration softkey. For instance, if Analyzer Configuration is set to Independent Channels, Span may be unlinked and the two displays set to have different spans. If Analyzer Config is set to Dual Channel, however, the span softkey will always be linked since Dual Channel measurements, such as frequency response, require both channels to have the same span.

Parameter linking is a convenient way of changing both displays together. For example, Span and Start frequency may be linked while the Measurement is unlinked. This allows the two displays to have different measurements over the same frequency span. Since the Spans are linked, changing the Span does not require separate entries for each display.

### **Enter**

When modifying a parameter, the [Enter] key enters the new parameter value or selection. When selecting from a list of choices with the knob, [Enter] will enter the new selection. When entering a numeric value, [Enter] will enter the new value. Whenever the 'Enter' indicator at the top of the screen is on, [Enter] is required to enter the new choice or entry.

### **ALT**

Indicates that the ALTERNATE keypad is in use. The ALTERNATE keypad uses the alphabetic legends printed below each key. To enter the ALT mode, press the [ALT] key once. Pressing the keys will now enter alphabetic characters into the active entry field. The [0]...[9], [.), [-], [<-] and [ALT] have the same function in the ALTERNATE keypad.

In ALTERNATE mode, the knob moves the waterfall marker along the Z axis.

To return to the normal keypad, press the [ALT] key again.

### **GPIB/RS232**

This indicator shows 'RS232' if the interface responses are directed to the RS232 serial port.

This indicator shows 'GPIB' and the GPIB address if the interface responses are directed to the GPIB port.

Use <Output To> in the [System]<Remote> menu to select the appropriate interface.

### **Comm /Error**

'Comm' is displayed when the computer interface is idle. 'RS232' flashes when there is activity on the RS232 interface. 'GPIB' flashes when there is activity on the GPIB interface.

'ERR' flashes whenever a computer interface error occurs, such as illegal command or out of range parameter is received.

### **SRQ**

This indicator is on whenever a GPIB Service Request is generated by the SR785. SRQ stays on until a serial poll is completed.

### Local/Remote/LOCK

'local' indicates that the unit is under local control (front panel enabled).

'REM' indicates that the unit is under remote interface control. The front panel may be locked out by a computer interface. To return the unit to local control, press the [Help/Local] key. <Override Remote> in the [System]<Remote> menu allows front panel control even in the REMOTE state.

'LOCK' indicates that the front panel is locked out and no front panel keypress will return the unit to the local state. No front panel adjustments may be made. An interface command (LOCL) is required to exit the LOCK state.

### macro

'Record' indicates that a keypad macro is being recorded. Use [Macro Rec] to start recording.

'Play' indicates that a macro is being played. Use [Play Macro] to playback a stored macro.

See 'Macros' later in this chapter for more.

# Keypad

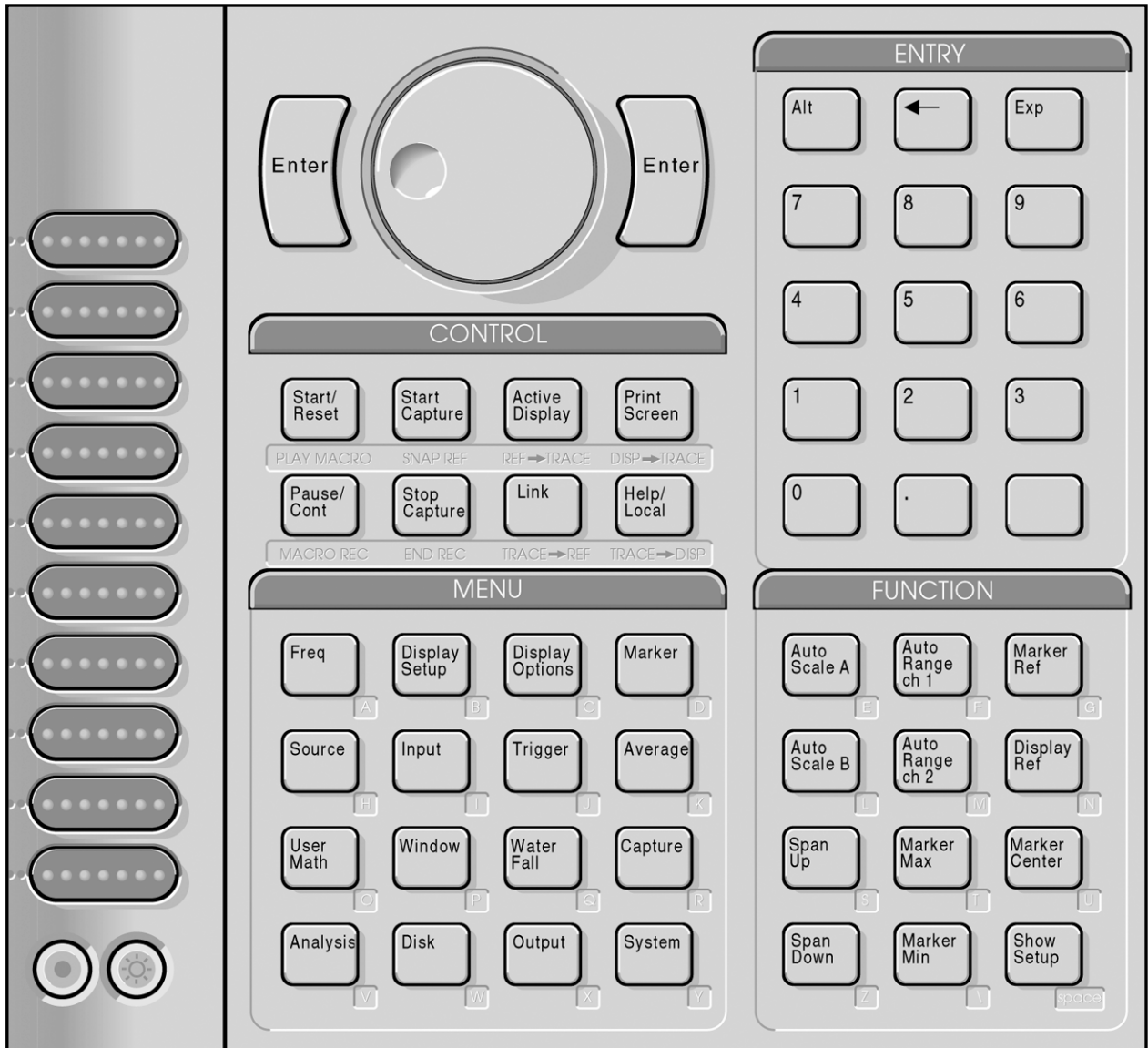


Figure Chapter 3 -10 Front Panel Keypad

## Normal and Alternate Keys

The normal key definitions are printed on each key. In addition, many keys also have alternate definitions printed below them. Press the [Alt] key once to enter the ALT keypad mode. The ALT screen indicator is on indicating that the alternate key definitions are in use. Use the alternate keys to enter alphabetic characters and to access secondary functions. The [0]...[9], [.] , [-], [<-], [Exp] and [Alt] keys have the same definition in both modes. Press [Alt] again to return to the normal keypad if necessary.

## Menu Keys

All operating parameters of the SR785 are grouped into sixteen function menus. The menu keys select a menu of parameters to display next to the ten softkeys. The softkeys either choose a submenu or select a parameter and place it in the entry field at the top of the screen (for numeric entry or knob adjustment). The menus are listed below.

[Freq]	Set the frequency span parameters for measurements.
[Display Setup]	Select the measurement type, view, units and scale.
[Display Options]	Select the display format, grid, user units and log X axis.
[Marker]	Configure the markers.
[Source]	Turns on and configures the source.
[Input]	Configures the signal and tachometer inputs and select the input ranges.
[Trigger]	Setup the trigger.
[Average]	Turn averaging on and select the averaging type.
[User Math]	Define user math functions.
[Window]	Select the measurement window.
[Waterfall]	Setup waterfall displays and storage.
[Capture]	Setup the time domain capture buffer.
[Analysis]	Setup data tables and limit testing.
[Disk]	Store and recall data and settings to and from disk.
[Output]	Configure the hardcopy output.
[System]	Configures the computer interfaces, preferences and clock. Diagnostics menu checks the instrument hardware.

Detailed descriptions of each menu are provided in the next chapter.

## Entry Keys

To modify or examine a parameter, press the menu key which displays the appropriate softkeys. Press a softkey to select (highlight) a parameter or function. There are three types of softkeys or menu boxes - buttons, lists and numeric values.

The active display determines which display's parameters are shown in the menus. Press [Active Display] to toggle between Display A and Display B parameters. Only those parameters which are associated with an individual display have differing values between the displays.

When a display specific parameter is selected, the Link indicator next to the top of the menu is turned on. If the indicator reads 'DispA' or 'DispB', parameter entry or selection modifies only Display A or Display B. If the indicator reads 'Link', both displays are modified at once. To change the Link indicator, use the [Link] key. When 'Link' is activated, the value of the parameter for the non-active display will be made equal to the value of the parameter for the active display.

### Buttons

Buttons are the simplest type of softkey. Pressing a button performs the specified action. Examples of buttons are <Full Span>, <Preferences> and <Display to Disk>. Buttons simply show a graphic 'button' or 'arrow' in their menu box.

### Lists

Parameters with a list of possible selections display their current selection in their softkey menu box. Examples of parameter lists are <View>, <Window> and <Trigger Source>.

To make a new selection, press the softkey to highlight its menu box. The list of possible choices is displayed in the entry field at the top of the screen. For example, to change the View, press [Display Setup] for the menu and then <View> (3rd softkey). The <View> menu box will be highlighted and the View selections will be listed at the top of the screen.

Use the knob to make the desired selection. If a vertical thermometer appears at the right of the entry field, then the selection list is too long to display at one time. The bright bar within the thermometer indicates the position within the complete list. Use the knob to scroll through the entire list. Choices shown in gray are not available in the current configuration. When the correct choice is selected, press [Enter] to change the parameter. The menu box now shows the new selection and becomes un-highlighted.

Many parameters require the [Enter] key to enter the new selection. In these cases, the Enter indicator next to the top of the menu will turn on while the selection is being made. No change is made until [Enter] is pressed. To cancel the operation, press the same softkey again (or another softkey to modify another parameter) instead of [Enter].

Some parameters do not require the [Enter] key to make the new selection valid. <Window> is an example. In this case, moving through the list of windows with the knob changes the Window immediately and displays the selection in the highlighted menu

box. This allows the effect of various choices to be compared while making the selection. Press [Enter], the same softkey, another softkey or a menu key to un-highlight the menu box when finished.

List selection may also be made using the numeric keys. The choices are numbered 0 through 9. Press a number key to make a selection. If [Enter] is not required, the new selection immediately appears in the menu box. If [Enter] is required, the new selection appears when [Enter] is pressed.

If no menu box is highlighted, pressing a numeric entry key automatically selects the most recently modified parameter within the menu and begins numeric entry. This is convenient when a measurement requires a single parameter within a menu (or menus) to be modified repeatedly.

### Numeric Values

Parameters with numeric values display their current value in their menu box. Examples of numeric parameters are <Span>, <# Avgs> and <Ch1 Input Range>. Some numeric values are 'continuous', such as <# Avgs>, while others can only have discrete values, such as <Span>.

To enter a new value, press the softkey to highlight its menu box. The parameter value is displayed in the entry field at the top of the screen. For example, to change the Span, press [Freq] for the menu and then <Span> (1st softkey). The <Span> menu box will be highlighted and the current Span value will be shown at the top of the screen.

The numeric entry keys are used to enter values directly. For example, to set the Span to 12.8 kHz, press [1] [2] [.] [8]. As soon as the [1] is pressed, the entry field shows the new value as it is entered. The list of appropriate units is displayed as well.

The entry field displays the characters as the keys are pressed. The '\_' is the entry point. If an error is made, the backspace [<-] key will erase the last character. Pressing the same softkey again (or another softkey to modify another parameter) will abort the entry operation and leave the value unchanged. Use the knob to select the (kHz) units. When the entry value is correct and the correct units are selected, press [Enter] to change the span to the new value.

Entries may be made in exponential form using the [Exp] key. The entry above may be made by pressing [1] [2] [.] [8] [Exp] [3] and selecting the (Hz) units.

Numeric entry always requires [Enter] to enter the new value.

Some menu boxes have 'Marker' as one of the units choices. If 'Marker' is selected as the units, the analyzer ignores the numeric keys pressed and substitutes the current marker value in the entry field.

If no menu box is highlighted, pressing a numeric entry key automatically selects the most recently modified parameter within the menu and begins numeric entry. This is convenient when a measurement requires a single parameter within a menu (or menus) to be modified repeatedly.



In general, numeric parameters may also be adjusted with the knob. In this case, knob adjustment must be made before any numeric entry is started. [Enter] is not required to enter the new value. Changes are effective immediately while the value is being adjusted. Press [Enter], the same softkey, another softkey or a menu key to un-highlight the menu box when finished. Some entry fields allow only knob adjustment or only numeric entry.

## Control Keys

### **[Start/Reset]**

If the unit is already making measurements, [Start/Reset] resets any averages and starts the measurement over. If the measurement is paused, [Start/Reset] starts the measurement over again. The Run/Pause indicator shows 'Run' when measurements are in progress.

[Start/Reset] also puts into effect any measurement related changes which may have been made since the last measurement was completed. For instance, if a set of averages is completed and the frequency span is changed, the new span does not take effect until the next time [Start/Reset] is pressed.

If the measurement is triggered, [Start/Reset] does not supply a trigger. If the trigger is manual armed, [Start/Reset] does not arm the trigger.

If the displays don't update, see the Operating Hints in Chapter 2 for help.

A TTL rising edge at the rear-panel 'Start' input is equivalent to pressing the [Start/Reset] key.

Command: STRT

### **[Pause/Cont]**

If the unit is already making measurements, [Pause/Cont] halts data acquisition. The measurement is paused where it is and the displays are not updated. The Run/Pause indicator shows 'Pause' when measurements are paused.

If the measurement is paused, [Pause/Cont] resumes the measurement where it was left off. To start over, press [Start/Reset]. The Run/Pause indicator shows 'Run' when measurements are in progress.

For Exponential averaging, [Pause/Cont] simply resumes data acquisition without resetting the average. For Linear averaging, [Pause/Cont] will complete the current average. If the linear average was already complete, [Pause/Cont] will start a new measurement..

For Swept Sine measurements, [Pause/Cont] simply resumes the sweep where it was paused.

If measurement related parameters are changed while the unit is paused, the measurement cannot be resumed. The SR785 will display the message "Can't continue, Measurement has Changed." In this case it is necessary to press [Start/Reset] to restart the measurement.

Command: PAUS and CONT

### **[Start Capture]**

[Start Capture] starts storing data in the capture buffer.

If the Trigger Arming Mode is Auto Arm and the Trigger Source is Continuous, capture starts immediately. Otherwise, the selected Trigger Arming Mode and Source will be in effect and the next trigger event will start capture.

While capture is in progress, the displays do not update. The Capture Progress indicator shows how much of the desired capture length has been completed. If the Capture Mode is 1-Shot, capture stops when the buffer is full. Press [Stop Capture] to halt capture before the buffer is full.

If the Capture Mode is Continuous, once capture is started, it continues indefinitely and fills the capture buffer in a circular fashion. In this case, press [Stop Capture] to halt capture with the most recently acquired data stored in the buffer.

Use the [Capture] <Memory Allocation> menu to allocate memory between the capture buffer, waterfall and order track storage and the arbitrary source waveform. Capture memory must be allocated before the capture buffer may be used.

To playback from the capture buffer, change the Input Source to Playback.

Command: CSTR

### **[Stop Capture]**

[Stop Capture] stops storing data in the capture buffer.

If the Capture Mode is 1-Shot, capture stops when the buffer is full. Press [Stop Capture] to halt capture before the buffer is full.

If the Capture Mode is Continuous, once capture is started, it continues indefinitely and fills the capture buffer in a circular fashion. In this case, press [Stop Capture] to halt capture with the most recently acquired data stored in the buffer.

Turn off Auto Offset during capture to avoid interrupting the input signal.

To playback from the capture buffer, change the Input Source to Playback.

Command: CSTP

### **[Active Display]**

Pressing [Active Display] toggles the active display. In single and overlay display format, the display switches between Display A and Display B. The display label, A or B, is at the upper left of each display. In dual display format, [Active Display] switches which display is active as indicated by the highlighted Marker Position bar above the graph.

The active display determines which display's parameters are displayed in the menus. For example, activating Display A and pressing [Display Setup] allows you to select the measurement for Display A. Pressing [Active Display] once makes Display B active and allows you to select the measurement for Display B using the same menu. Only those parameters which are associated with an individual display have differing values

between the displays. Parameters such as input configuration or source setup are instrument parameters and are not associated with a display.

When a display specific parameter is highlighted (for modification), the Link indicator next to the top of the menu is turned on. If the indicator reads 'DispA' or 'DispB', parameter entry or selection modifies only Display A or Display B. If the indicator reads 'Link', both displays are modified at once. To change the Link indicator, use the [Link] key.

Parameter linking is a convenient way of changing both displays together. For example, Span and Start frequency may be linked while the Measurement is unlinked. This allows the two displays to have different measurements over the same frequency span. Since the Spans are linked, changing the Span does not require separate entries for each display.

Command: ACTD (?) i

### [Link]

The [Link] key is used to link and unlink display specific parameters, to perform a function on both displays at once, or to temporarily link the display markers.

When an unlinked display specific parameter is highlighted (for modification), the Link indicator (next to the top of the menu) reads either 'DispA' or 'DispB', whichever is the active display. Parameter entry or selection modifies only the active display. To change which display will be modified, press [Active Display].

When a linked parameter is highlighted, the Link indicator reads 'Link'. Parameter entry or selection modifies both displays at once.

To link or unlink a parameter, highlight the parameter and press the [Link] key. If the Link indicator is shown in gray, then the parameter linking may not be changed. If linking is available the value of the parameter for the non active display will be made equal to the value of the parameter for the active display.

Parameter linking is a convenient way of changing both displays together. For example, Span and Start frequency may be linked while the Measurement is unlinked. This allows the two displays to have different measurements over the same frequency span. Since the Spans are linked, changing the Span does not require separate entries for each display.

The availability of parameter linking is affected by the setting of the <Analyzer Configuration> softkey. If Analyzer Configuration is set to 'Independent Channels,' parameters such as span may be linked unlinked for the two displays. If <Analyzer Configuration> is set to 'Dual Channel,' these parameters will always be linked. See Chapter 2 for a description of how Analyzer Configuration works.

When a function key, such as [Marker Max], is pressed, the function is performed on the active display. If [Link] is pressed immediately before, then the function is performed on both displays at once. [Link] only affects the next key pressed. If the next key is not a function key then it has no affect. Linking affects each function differently. See the function key definitions for more information.

When no parameter softkey is highlighted, turning the knob adjusts the marker position in the active display. If [Link] is pressed, then the marker of the inactive display is linked to the active display and the knob moves both markers together. Pressing a key will cancel this link. If the key is a function key, its function will be linked.

### **[Print Screen]**

The [Print Screen] hardkey allows you to print or plot the screen from any menu.

Configure printing and plotting in the [Output] menu and choose Print, Plot or Dump Screen as the Hard Copy Button function to configure the [Print Screen] hardkey.

Pressing [<-] (backspace) will abort the printout. No other front panel operations may be performed until printing is completed.

Command: PRNT, PLOT or DUMP

### **[Help/Local]**



[Help/Local] provides on screen help for any key or softkey. Press [Help/Local] to enter the on-line help system. Press [Help/Local] again for more information on the help system.

When a host computer places the unit in the REMOTE state, no keypad or knob input is allowed. To return to front panel operation, press the [Help/Local] key.

### **[Alt] [Control Key]**

Pressing [Alt] and a Control Key performs the function written below each key. These functions are described below.

#### **Play Macro ([Alt] [Start/Reset])**

[Play Macro] plays a recorded keypad macro. Use the knob to select a defined macro and press [Enter]. The Macro indicator at the top of the screen reads 'Play' while the macro is playing.

#### **Macro Rec ([Alt] [Pause/Cont])**

[Macro Rec] starts recording a macro. Use the knob to select which macro (0..9) will be recorded and press [Enter]. The Macro indicator at the top of the screen reads 'Record'. Press the desired sequence of keys. The instrument responds to these keys while the macro is being recorded. Press [End Rec] ([Alt] [Stop Capture]) to stop recording the macro.

Use the numeric keys to choose a parameter from a list (instead of the knob) while recording a macro. The list choices are numbered 0 through 9.

Use [Alt] and the numeric keys to choose the units of a numeric parameter while recording a macro. The units choices are numbered starting with 0 as the first (leftmost) units.

To play a recorded macro, press [Play Macro] ([Alt] [Start/Reset]) and choose the desired macro with the knob.

To edit a recorded macro, use [System] <Edit Macro>.

### **End Rec ([Alt] [Stop Capture])**

[End Rec] stops recording a macro.

To play a recorded macro, press [Play Macro] ([Alt] [Start/Reset]) and choose the desired macro with the knob.

### **Snap Ref ([Alt] [Start Capture])**

[Snap Ref] scales the Reference Display within the active display to the current display scale. Changing the display scale does NOT rescale the Reference Display. This allows the live measurement display to be offset from the Reference Display by changing either Ymax, Ymid or Ymin in the [Display Setup] menu.

Press [Display Ref] to copy the display data into the Reference Display. The Reference Display is graphed in the background.

Command: SNAP d

### **Ref to Trace ([Alt] [Active Display])**

[Ref to Trace] saves the active display's Reference Display (if it is on) to a Trace buffer.

Select a Trace # (1..5) with the knob and press [Enter] to save the Reference Display to the Trace.

A stored trace can be recalled to a Display or Reference Display, used in a User Math Function, saved to disk, or copied to the Arbitrary Waveform buffer.

Command: SVRF d, i

### **Trace to Ref ([Alt] [Link])**

[Trace to Ref] copies Trace data into the active display's Reference Display.

Select a Trace # (1..5) with the knob and press [Enter] to copy the Trace to the Reference Display. Only those Traces which currently have data can be selected. In addition, the Trace data must be compatible with active display measurement. For example, if the Trace is a time record, it cannot be copied into the Reference Display of an FFT measurement.

Command: RCRF d, i

### **Display to Trace ([Alt] [Print Screen])**

[Display to Trace] saves the active display's measurement data to a Trace buffer.

Select a Trace # (1..5) with the knob and press [Enter] to save the display data to the Trace.

A stored trace can be recalled to a Display or Reference Display, used in a User Math Function, saved to disk, or copied to the Arbitrary Waveform buffer.

Command: SVTR d, i

### **Trace to Display ([Alt] [Help/Local])**

[Trace to Display] recalls Trace data to the active display.

Select a Trace # (1..5) with the knob and press [Enter] to recall the Trace data to the active display.

The active display will be set to Off-Line when data is recalled to it. The display will not update with live measurement results until the Display is set back to Live (in the [Display Options] menu).

The measurement, frequency span, averaging and window are recalled with the data and cannot be changed. These menus are shown in gray and reflect the settings or values for the last Live measurement. Parameters which pertain to the display of the data, such as Scale and View, may be changed in the [Display Setup] and [Display Options] menus.

When the Display is returned to Live, the measurement, frequency span, averaging and window all return to the settings in effect before the data was recalled and the live measurement returns to the display.

Command: RCTR d, i

## Function Keys

### [Auto Scale A]

[Auto Scale A] automatically sets the vertical scale of Display A to show the entire range of the data. Horizontal scaling is not affected.

Auto scale only operates on the data which is displayed on the graph. If the graph is expanded, data corresponding to frequency or time bins which are not shown do not figure in the auto scaling calculations.

[Link] [Auto Scale A] auto scales Display A and sets the vertical scale of Display B to match Display A.

Command: ASCL d

### [Auto Scale B]

[Auto Scale B] automatically sets the vertical scale of Display B to show the entire range of the data. Horizontal scaling is not affected.

Auto scale only operates on the data which is displayed on the graph. If the graph is expanded, data corresponding to frequency or time bins which are not shown do not figure in the auto scaling calculations.

[Link] [Auto Scale B] auto scales Display B and sets the vertical scale of Display A to match Display B.

Command: ASCL d

### [Auto Range Ch1]

[Auto Range Ch1] toggles Channel 1 Input Ranging between Manual and Auto. In Manual Ranging, the Input Range is set within the [Input]<Input Config> submenu. In Auto Range, the Input Range is adjusted automatically according to the Ch1 AutoRange Mode (also in the [Input] menu). In Up Only Auto Range, only overloads cause the range to change. In Tracking Auto Range, the Input Range moves up for overloads and down when the signal falls below half scale.

In Swept Sine group, Auto Range is always tracking.

The Ch1 Input Range indicator at the top of the screen will be shown in inverse if Ch1 is Auto Ranging.

Auto Range responds to all frequencies present at the input (except those attenuated by AC coupling), not just those within the measurement span.

[Link] [Auto Range Ch1] toggles the Input Ranging of both channels.

Command: AIRG (?) {i}



**[Auto Range Ch2]**

[Auto Range Ch2] toggles Channel 2 Input Ranging between Manual and Auto. In Manual Ranging, the Input Range is set within the [Input]<Input Config> submenu. In Auto Ranging, the Input Range is adjusted automatically according to the Ch2 AutoRange Mode (also in the [Input] menu). In Up Only Auto Range, only overloads cause the range to change. In Tracking Auto Range, the Input Range moves up for overloads and down when the signal falls below half scale.

In Swept Sine group, Auto Range is always tracking.

The Ch2 Input Range indicator at the top of the screen will be shown in inverse if Ch2 is Auto Ranging.

Auto Range responds to all frequencies present at the input (except those attenuated by AC coupling), not just those within the measurement span.

[Link] [Auto Range Ch2] toggles the Input Ranging of both channels.

Command: A2RG (?) {i}

**[Span Up]**

In FFT group, [Span Up] increases the FFT frequency Span of both displays by a factor of 2. This provides a way of adjusting the spans regardless of which menu is displayed. The Span is adjusted with either a fixed start, Center or End frequency depending upon which frequency was most recently anchored in the [Frequency] menu.

In Correlation group, [Span Up] decreases the time window over which the correlation is computed, but increases the time resolution.

In Octave group, [Span Up] moves the octave measurement span of both displays up by an octave.

[Span Up] affects each display independently if their frequency parameters are unlinked.

[Span Up] has no affect for Swept Sine, Order, or Time/Histogram measurements.

[Link] [Span Up] is the same as [Span Up].

Command: FSPN (?) d {, i}

**[Span Down]**

In FFT group, [Span Down] decreases the FFT frequency Span of both displays by a factor of 2. This provides a way of adjusting the spans regardless of which menu is displayed. The Span is adjusted with either a fixed Start, Center or End frequency depending upon which frequency was most recently anchored in the [Frequency] menu.

In Correlation group, [Span Down] increases the time window over which the correlation is computed, but decreases the time resolution.

In Octave group, [Span Down] moves the octave measurement span of both displays down by an octave.

[Span Down] affects each display independently if their frequency parameters are unlinked.

[Span Down] has no affect for Swept Sine measurements, Order measurements, or Time/Histogram measurements.

[Link] [Span Down] is the same as [Span Down].

Command: FSPN (?) d {, i}

### [Marker Ref]

[Marker Ref] toggles the Marker Rel Mode (in the [Marker] menu) between Off and Relative for the active display. If Marker Rel is Off, [Marker Ref] sets Marker Rel to Relative and sets the X-Rel and Y-Rel Marker Offsets to the current Marker Position. Subsequent marker readings are relative to this offset. The relative Marker Position values are preceded by a  $\Delta$  (delta) symbol. The Marker Offset location within the display is marked by small flag shaped symbols.

If Marker Rel is Relative, [Marker Ref] sets Marker Rel to Off and the Marker Position is absolute.

[Marker Ref] has no affect if the Marker Rel Mode is not Off or Relative.

[Link] [Marker Ref] toggles the Marker Rel Mode of both displays.

Command: MRON (?) d

### [Display Ref]

[Display Ref] toggles the Reference Display within the active display. If the Reference Display is off, [Display Ref] loads the Reference Display with the current display data. The Reference Display is graphed in half intensity. This allows comparison of new data with the stored Reference Display. Set the Marker Relative to Reference Display (in the [Marker] menu) to read the Marker Position relative to the Reference Display. Press [Display Ref] again to turn off the Reference Display.

If display parameters are changed to values incompatible with the Reference Display, the reference display will be turned off.

Changing the display scale does NOT rescale the Reference Display. This allows the live measurement display to be offset from the Reference Display by changing either Ymax, Ymid or Ymin in the [Display Setup] menu. To graph the Reference Display in the new display scale, press [Alt] [Snap Ref].

[Link] [Display Ref] toggles the Reference Display of both displays.

Command: DREF d, i

**[Marker Center]**

[Marker Center] sets the FFT Center frequency to the Marker frequency in the active display. If the new Center frequency and Span conflict with the maximum measurement range, 0 to 102.4(100.0) kHz, then the Span is decreased to the largest span which allows the Marker frequency to be the center.

[Link] [Marker Center] adjusts the span of the active display and sets the Span of the inactive display to match.

Command: MKCN d

**[Marker Max]**

[Marker Max] moves the Marker to the location of the maximum data value within the active display. [Marker Max] only searches the data which is in the display. If the maximum value occurs at more than one location, then the one closest to the left edge is found.

If the Marker Mode (in the [Marker] menu) is Normal, [Marker Max] centers the marker region around the maximum. The Marker Seeks Mode chooses whether the Marker seeks the Maximum, Minimum or Mean of the data within the Marker Region.

For Harmonic or Sideband markers, [Marker Max] moves the Fundamental Marker to the location of the maximum. [Marker Max] has no affect on Band markers.

[Link] [Marker Max] moves the Marker to the location of the maximum within both displays.

Command: MKMX d

**[Marker Min]**

[Marker Min] moves the Marker to the location of the minimum data value within the active display. [Marker Min] only searches the data which is in the display. If the minimum value occurs at more than one location, then the one closest to the left edge is found.

If the Marker Mode (in the [Marker] menu) is Normal, [Marker Min] centers the Marker Region around the minimum. The Marker Seeks Mode chooses whether the marker seeks the Maximum, Minimum or Mean of the data within the Marker Region.

For Harmonic or Sideband markers, [Marker Max] moves the Fundamental Marker to the location of the minimum. [Marker Max] has no affect on Band markers.

[Link] [Marker Min] moves the marker to the location of the minimum within both displays.

Command: MKMN d

**[Show Setup]**

[Show Setup] enters the Help system and displays the measurement setup and system settings. Press [0] to exit Help and return to the measurement displays.

# Macros

## Keypad Macros

A sequence of front panel key presses can be recorded and played as a macro. There are 10 macros which may be defined. Macro definitions are retained at power down. Macros are saved to disk with the instrument settings.

Use [Macro Rec] ([Alt] [Pause/Cont]) to start recording a macro. Use the knob to select which macro (0..9) will be recorded and press [Enter]. Press the desired sequence of keys. The instrument responds to these keys while the macro is being recorded. Press [End Rec] ([Alt] [Stop Capture]) to stop recording the macro.

To play a recorded macro, press [Play Macro] ([Alt] [Start/Reset]), choose the desired macro with the knob and press [Enter].

To edit a recorded macro, use [System] <Edit Macro>.

## Choosing From Lists

The knob is normally used to pick an entry from a list. However, knob movements are recorded as movement relative to the starting position. This leads to ambiguity when using the knob in a macro. Playing the macro will have different results depending upon the starting state of the instrument. Thus, using the knob within a macro to choose a parameter is not recommended.

To choose a parameter from a list, the numeric keys should be used instead of the knob while recording a macro. The list choices are numbered 0 through 9. Press a number key to make a selection and [Enter] if required. For example, to select the Hanning window, press [Window] <Window> [2] [Enter] since Hanning is entry number 2.

If the parameter list has more than 10 choices, FFT Measurement for example, then the knob must be used. In this case, use a numeric key entry to first pick a known selection and then use the knob to select an entry in the list relative to it. For example, to choose Capture1 as the FFT Measurement, press [Display Setup] <Measurement> [0] [Enter] <Measurement> {select Capture1 with the knob} [Enter]. The first entry selects FFT1 with the [0] key regardless of the current measurement. The knob then selects a position in the list relative to FFT1.

## Menus

When recording a macro, remember that the macro is played back EXACTLY as recorded. Thus, do not start a macro with a softkey press since the menu which is displayed at playback may not be the one displayed when the macro was recorded. Always start a macro with a menu key press if you want to change parameters within a menu.

## Chapter 4

# Menus

The SR785 has a menu driven user interface. All operating parameters of the SR785 are grouped into sixteen menus. The Menu keys each display a menu of softkeys. The softkeys at the right of the display change depending upon the displayed menu.

There are three types of softkeys - buttons, lists and numeric values. A button performs a function, such as <Full Span>. A list presents a list of choices or options in the entry field (at the top of the screen). Use the knob to make a selection and press [Enter]. <Measurement> is an example of a list. A numeric value presents the current value in the entry field and awaits numeric entry. Enter a new value with the entry keys and press [Enter]. <Start Freq> is an example of a numeric value.

Some parameters can have different values or settings for each display. In these cases, the menu shows the values or settings for the active display. Press [Active Display] to toggle the active display and show the values or settings for the other display. <Measurement> is an example of such a parameter. Use [Link] to link a parameter to a single display or to both displays at once.



Help for any key or softkey is available on screen by pressing [Help/Local]. Simply press any key for help about that key or function. Press [0] to exit help.

## In This Chapter

<b>Frequency Menus</b>	<b>4-7</b>	<b>Swept Sine Frequency Menu</b>	<b>4-15</b>
<b>FFT Frequency Menu</b>	<b>4-7</b>	Start	4-15
Span	4-7	Stop	4-15
Linewidth	4-8	Repeat	4-16
Acquisition Time	4-8	Type	4-16
Full Span	4-8	Auto Resolution	4-16
FFT Lines	4-9	Number of Points	4-17
Base Frequency	4-9	Maximum Step Size	4-17
Start Frequency	4-9	Faster Threshold	4-17
Center Frequency	4-10	Slower Threshold	4-18
End Frequency	4-10	<b>Order Frequency Menu</b>	<b>4-19</b>
<b>Octave Frequency Menu</b>	<b>4-12</b>	Max RPM	4-19
Highest Band	4-12	Min RPM	4-19
Lowest Band	4-13	Max Order	4-19
Octave Resolution	4-13	Delta Order	4-20
Octave Channels	4-14	Tracking	4-20
		Track Setup	4-20

## 4-2 Menus

---

Track Points	4-20	<b>Normal Marker Menu</b>	<b>4-54</b>
Storage	4-21	Seeks	4-54
Track 1 Order	4-21	Width	4-54
Input Channel (Track 1)	4-21	Rel	4-55
Track 2 Order	4-21	X-Rel	1 4-55
Input Channel (Track 2)	4-22	Y-Rel	4-56
		X-Rel	4-56
		Marker X to	4-57
		Target	4-57
<b>Time/Histogram Frequency Menu</b>	<b>4-23</b>	<b>Harmonic Marker Menu</b>	<b>4-58</b>
Sample Time	4-23	Harmonics	4-58
Histogram Length	4-23	Display	4-58
Bins	4-23	Readout Mode	4-59
Base Frequency	4-24	THD	4-59
Repeat	4-24	Harmonic Power	4-59
<b>Display Setup Menu</b>	<b>4-25</b>	<b>Sideband Marker Menu</b>	<b>4-61</b>
Measurement Group	4-25	Sideband Separation	4-61
Measurement (FFT)	4-26	# Sidebands	4-61
Measurement (Correlation)	4-30	Display	4-62
Measurement (Octave Analysis)	4-32	Readout Mode	4-62
Measurement (Swept Sine)	4-34	Sideband Ratio	4-62
Measurement (Order)	4-35	RMS Sideband Power	4-63
Measurement (Time/Histogram)	4-37	<b>Band Marker Menu</b>	<b>4-64</b>
View	4-38	Modify Band	4-64
Units	4-40	Band Exclude	4-64
Units	4-41	Band Ratio Mode	4-65
dB Units	4-41	Band Power	4-65
Pk Units	4-41	Band Ratio	4-65
PSD Units	4-42	<b>Frequency/Damping Marker Menu</b>	<b>4-66</b>
Phase Units	4-42	Calculate	4-66
dBm Reference Impedance	4-42	Frequency	4-66
Y Max	4-42	Damping Ratio	4-66
Y Mid	4-42	<b>Source Menu</b>	<b>4-67</b>
Y Min	4-43	Source Off	4-67
Y / Div	4-43	Source On	4-67
X Center (Polar)	4-43	Sine	4-67
X / Div (Polar)	4-44	Chirp	4-68
Y Center (Polar)	4-44	Noise	4-68
Y / Div (Polar)	4-44	Arbitrary	4-69
Pan	4-44	<b>Sine Source Menu</b>	<b>4-71</b>
Zoom	4-45	Frequency 1	4-71
		Amplitude 1	4-71
		Offset	4-72
		Tone 2	4-72
		Frequency 2	4-72
		Amplitude 2	4-72
<b>Display Options Menu</b>	<b>4-46</b>	<b>Chirp Source Menu</b>	<b>4-73</b>
Display	4-46	Amplitude	4-73
Format	4-46	Burst	4-73
X-Axis	4-47	Source Display	4-74
RPM Frequency	4-47		
Grid	4-47		
Grid Div	4-48		
Nyquist Grid	4-48		
Phase Suppress	4-48		
d/dx Window	4-49		
<b>Marker Menu</b>	<b>4-50</b>		
Marker	4-50		
Mode	4-50		

<b>Noise Source Menu</b>	<b>4-75</b>	Tach Slope	4-96
Amplitude	4-75	Hold Off Enable	4-96
Type	4-75	Tach Hold Off	4-96
Burst	4-76	Show Tach	4-96
Source Display	4-76		
Source Period	4-77		
<b>Arbitrary Source Menu</b>	<b>4-78</b>	<b>Playback Submenu</b>	<b>4-97</b>
Amplitude	4-78	Playback Start	4-97
Source Play Rate	4-78	Playback Length	4-97
Source	4-79	Set Left Edge	4-97
More	4-79	Set Right Edge	4-98
Arb. Source Start	4-80	Playback Mode	4-98
Arb. Source Length	4-80	Playback Speed	4-98
Set Left Edge	4-80		
Set Right Edge	4-80	<b>Trigger Menu</b>	<b>4-100</b>
Allocate Memory	4-81	Trigger Mode	4-100
Trace to Arb	4-81	Trigger Source	4-101
		Trigger Level	4-102
		Trigger Slope	4-102
		Delay 1	4-102
		Delay B	4-102
<b>Swept Sine Source Menu</b>	<b>4-83</b>	Trigd Source Mode	4-103
Auto Level Reference	4-83	RPM/Time Arm Setup	4-104
Amplitude	4-83	Start RPM	4-104
Ideal Reference	4-84	Start RPM (On/Off)	4-104
Source Ramping	4-84	Delta RPM	4-104
Source Ramp Rate	4-84	Delta RPM Sense	4-105
Reference Upper Limit	4-85	Time Arm Step	4-105
Reference Lower Limit	4-85	Manual Arm	4-105
Maximum Source Level	4-85	Manual Trigger	4-105
Offset	4-86		
<b>Input Menu</b>	<b>4-87</b>	<b>Average Menus</b>	<b>4-107</b>
Input Source	4-87		
Analyzer Configuration	4-87	<b>FFT, Correlation, Order Average Menu</b>	<b>4-107</b>
Auto Offset	4-88	Compute Averages	4-107
		Averaging Type	4-108
<b>Input Configuration Submenu</b>	<b>4-89</b>	Number of Averages	4-109
Channel	4-89	Display Average	4-109
Input Mode	4-89	Time Record Increment	4-111
Grounding	4-89	More	4-112
Coupling	4-90	Overload Reject	4-112
Input Range	4-90	Trigger Average Mode	4-113
AA Filter	4-90	Average Preview	4-113
A-Wt Filter	4-91	Preview Time	4-114
Autoranging	4-91	Accept	4-114
		Reject	4-114
<b>Transducer Parameter Submenu</b>	<b>4-92</b>	<b>Octave Average Menu</b>	<b>4-115</b>
Engineering Units	4-92	Averaging Type	4-115
EU Label	4-93	Integration Time	4-116
EU / Volt	4-93	Confidence	4-117
EUs at Marker	4-93	Power Bin	4-117
User Label	4-94	Linear Average Trigger	4-117
Transducer Convert	4-94		
<b>Tachometer Input Submenu</b>	<b>4-95</b>	<b>Swept Sine Average Menu</b>	<b>4-119</b>
Pulses per Rev	4-95	Settle Time	4-119
Tach Trigger Range	4-95	Settle Cycles	4-119
Tach Level	4-95		

## 4-4 Menus

---

Integration Time	4-120	Total Available	4-145
Integration Cycles	4-120	Capture Memory	4-145
<b>User Math Menu</b>	<b>4-122</b>	Waterfall Memory	4-145
Function	4-122	Arb Memory	4-145
Edit Function	4-122	Confirm Allocation	4-146
Operands	4-123	Clear Allocation	4-146
Operations	4-126	View Header	4-146
Function String	4-128	Auto Pan	4-146
Insert/Replace	4-128	<b>Analysis Menu</b>	<b>4-147</b>
Delete	4-129	Data Table	4-147
Clear Eq.	4-129	Limit Test	4-147
Cancel	4-129	Marker Stats	4-148
Enter Eq.	4-129	Exceedance Stats	4-148
Edit Constant	4-129	Curve Fit	4-148
Constant	4-130	<b>Data Table Analysis Menu</b>	<b>4-150</b>
Real Part	4-130	Data Table	4-150
Imaginary Part	4-130	Insert Line	4-150
Mag	4-130	Modify Line	4-150
Phase(deg)	4-130	Delete Line	4-151
Marker to Mag	4-130	Clear Table	4-151
<b>Window Menu</b>	<b>4-131</b>	<b>Limit Testing Analysis Menu</b>	<b>4-152</b>
Window	4-131	Limit Segments	4-152
Channel 1 Window	4-133	Limit Testing	4-152
Channel 2 Window	4-133	Limit Beep	4-152
Force Length	4-133	Clear Limits	4-153
Expo TC	4-133	Edit Limits	4-153
Trace to Window	4-134	New Segment	4-153
Window to Trace	4-134	Limit Type	4-153
Window Form	4-134	Segment #	4-154
<b>Waterfall Menu</b>	<b>4-135</b>	X0	4-154
Display	4-135	Y0	4-154
Storage	4-136	X1	4-154
Save Option	4-136	Y1	4-154
Total Count	4-137	Delete Segment	4-155
Skip/Storage Interval	4-137	Shift All	4-155
View Count	4-138	<b>Marker Statistics Menu</b>	<b>4-156</b>
Angle	4-138	Marker Stats.	4-156
Marker to Z	4-139	Reset	4-156
Allocate Memory	4-139	Max (Display A)	4-157
More	4-140	Min (Display A)	4-157
Trace Height	4-140	Mean (Display A)	4-157
Fast Angles	4-141	Std. Dev. (Display A)	4-157
Threshold	4-141	Max (Display B)	4-157
Hidden Lines	4-141	Min (Display B)	4-158
Paused Drawing	4-141	Mean (Display B)	4-158
Record to Trace	4-142	Std. Dev. (Display B)	4-158
Slice to Trace	4-142	<b>Exceedance Statistics Menu</b>	<b>4-159</b>
<b>Capture Menu</b>	<b>4-143</b>	Start Index	4-159
Capture Channels	4-143	Stop Index	4-159
Capture Mode	4-143	Exceedance Pct	4-160
Capture Length	4-144	Calculate Excd	4-160
Sampling Rate	4-144		
Allocate Memory	4-145		



<b>Curve Fit Menu</b>	<b>4-161</b>	Load Trace Data (ASCII)	4-178
Start Fit	4-161	Load Trace Data (binary)	4-179
Synthesize Table (1 and 2)	4-161	Buffer	4-179
Table	4-162	Interval	4-180
<b>Fit Setup Submenu</b>	<b>4-163</b>	Disk to Buffer	4-180
Number Poles	4-163	Buffer to Disk	4-181
Number Zeros	4-163	<b>Disk Upkeep Menu</b>	<b>4-182</b>
Weighting	4-163	File Name	4-182
Weighting Trace	4-164	Current Directory	4-182
Set Fit Region	4-164	Make Directory	4-183
<b>Edit Table Submenu</b>	<b>4-165</b>	Del File	4-183
Edit Item	4-165	Del Cur Dir	4-183
Delete Item	4-165	Format Floppy	4-183
Table Format	4-165	<b>Output Menu</b>	<b>4-184</b>
Copy Other Table	4-166	Hard Copy Button	4-184
Clear Table	4-166	Bitmap/Printer	4-185
Enter j	4-166	Bitmap Area	4-185
<b>Disk Menu</b>	<b>4-167</b>	Vector/Plotter	4-185
File Name	4-167	Destination	4-186
Current Directory	4-168	GPIB Control	4-186
Display to Disk	4-168	GPIB Address	4-187
Disk to Display	4-168	Edit A Note	4-187
Settings to Disk	4-169	Text X Position	4-187
Recall Settings	4-169	Text Y Position	4-187
Trace to Disk	4-169	Display	4-188
Disk to Trace	4-169	Visible	4-188
Buffers	4-170	Text String	4-188
Disk Upkeep	4-170	File Start Number	4-188
<b>Nodal Degree-of-Freedom Menu</b>	<b>4-171</b>	Colors	4-189
Reference Name	4-171	Print Bright	4-189
Reference Number	4-171	Print Dim	4-189
Reference Direction	4-172	Print Black	4-189
Response Name	4-172	Print Graph	4-190
Response Number	4-172	Plotter Text Pen	4-190
Response Direction	4-172	Plotter Grid Pen	4-190
Abort Save	4-173	Plotter Trace Pen	4-190
Continue Save	4-173	Plotter Marker Pen	4-190
<b>Recall Settings Menu</b>	<b>4-174</b>	<b>System Menu</b>	<b>4-191</b>
File Name	4-174	Preset	4-191
Current Directory	4-174	Remote	4-191
Measurements	4-175	Preferences	4-191
Sources	4-175	Date/Time	4-192
Analysis	4-175	Diagnostics	4-192
Inputs/Triggers	4-175	Macro	4-192
DRAM settings	4-175	Edit Macro	4-192
General System	4-176	Show Settings	4-193
Macros	4-176	Show Version	4-193
Recall from Disk	4-176	<b>System Remote Menu</b>	<b>4-194</b>
<b>Disk Buffers Menu</b>	<b>4-177</b>	Output To	4-194
File Name	4-177	GPIB Address	4-194
Current Directory	4-177	Override REM?	4-194
		Baud Rate	4-194
		Word Length	4-195

## 4-6 Menus

---

Parity	4-195	String to Printer	4-201
View Queues	4-195	Memory Test	4-201
<b>System Preferences Menu</b>	<b>4-196</b>	System RAM	4-201
Key Click	4-196	System ROM	4-201
Alarms	4-196	Video RAM	4-201
Alarms Volume	4-196	Help ROM	4-201
Done Volume	4-196	Disk Test	4-201
Audible Overload	4-197	Serial Number	4-202
Screen Saver	4-197	Program S/N	4-202
Screen Saver Delay	4-197	Kill AutoOff	4-202
Frequency Format	4-197	<b>Edit Macro Menu</b>	<b>4-203</b>
Node Info Prompt	4-198	Special Keys	4-203
<b>System Date/Time Menu</b>	<b>4-199</b>	Macro String	4-203
Time	4-199	Insert/Replace	4-203
Date	4-199	Inc. Knob Count	4-204
<b>System Diagnostics Menu</b>	<b>4-200</b>	Dec. Knob Count	4-204
Keypad Test	4-200	Delete	4-204
Keyboard Test	4-200	Clear Macro	4-204
Knob Test	4-200	Cancel	4-204
RS232/Printer Test	4-200	Enter Macro	4-204
String to RS232	4-201		

## Frequency Menus

The Measurement Group determines which Frequency menu is displayed. This menu sets the frequency parameters for the measurement. A set of parameters is stored for each group.

### FFT Frequency Menu

When the measurement group is FFT or Correlation, this menu sets the frequency span and location for the active display.

Frequency	
Span	51.2 kHz
Line Width	128 Hz
Acquisition Time	7.81 ms
Full Span	<input type="checkbox"/>
FFT Lines	400
Base Freq.	102.4 kHz
Start Freq.	0 Hz
Center Freq.	25.6 kHz
End Freq.	51.2 kHz

#### Span

Select the FFT Span of the active display. The frequency Span ranges from the FFT Base Frequency (102.4 kHz or 100.0 kHz) to  $2^{-19}$  times the Base Frequency (195.3 mHz or 191 mHz) in factors of 2. A numerically entered value is rounded to the nearest allowable Span. If the new frequency Span would extend below 0 Hz or above the Base Frequency, then the Start and End frequencies will be adjusted.

Changing the Span will change the Linewidth (Span/FFT Resolution) and Acquisition Time (FFT Resolution/Span). The Linewidth and Acquisition time are other ways to change the Span.

The two displays can have different FFT Spans if the Analyzer Configuration is set to Independent Channels. In this configuration, no two channel measurements are allowed (frequency response, cross-spectrum, etc.), but the entry field can be linked to both displays using the [Link] key. If Analyzer Configuration is set to Dual Channel the field is automatically linked to both displays.

Command: FSPN (?) d {, f}

---

### Linewidth

Select the Linewidth of the active display. The Linewidth is defined as the Span divided by the FFT Resolution (100, 200, 400 or 800 lines). A numerically entered value is rounded to the nearest allowable Linewidth.

Changing the Linewidth will change the Span (Linewidth x FFT Resolution) and Acquisition Time (1/Linewidth). The Linewidth and Acquisition Time are other ways to change the Span. If the new frequency Span would extend below 0 Hz or above the FFT Base Frequency, then the Start and End frequencies will be adjusted.

The two displays can have different Linewidths if the Analyzer Configuration is set to Independent Channels. In this configuration, no two channel measurements are allowed (frequency response, cross-spectrum, etc.), but the entry field can be linked or unlinked using the [Link] key. If Analyzer Configuration is set to Dual Channel the field is automatically linked to both displays.

Command: FSPN (?) d {, f}

---

### Acquisition Time

Select the Acquisition Time of the active display. The Acquisition Time is the time record length for the FFT (FFT Resolution/Span). A numerically entered value is rounded to the nearest allowable Acquisition Time.

Changing the Acquisition Time will change the Span (FFT Resolution/Acquisition Time) and Linewidth (1/Acquisition Time). The Linewidth and Acquisition time are other ways to change the Span. If the new frequency Span would extend below 0 Hz or above the FFT Base Frequency, then the Start and End frequencies will be adjusted.

The two displays can have different Acquisition Times if the Analyzer Configuration is set to Independent Channels. In this configuration, no two channel measurements are allowed (frequency response, cross-spectrum, etc.), but the entry field can be linked or unlinked using the [Link] key. If Analyzer Configuration is set to Dual Channel the field is automatically linked to both displays.

Command: FSPN (?) d {, f}

---

### Full Span

Set the FFT Span of both displays to the FFT Base Frequency 102.4 (100.0) kHz. The Linewidth, Acquisition Time, Start frequency and Center frequency are set to 256 (250) Hz, 3.906 (4.00) ms, 0.0 Hz, and 51.2 (50.0) kHz respectively.

## FFT Lines

Select the FFT Resolution of the active display [100, 200, 400, 800 lines].

Changing the FFT Resolution does not change the Span. Instead, the Acquisition Time is changed (FFT Resolution/Span).

Fewer lines means wider linewidths (poorer resolution) but faster measurements. More lines means narrower linewidths (better resolution) but slower measurements.

The various FFT resolutions are summarized below.

FFT <u>Resolution</u>	Frequency <u>Resolution</u>	<u>Time Record</u>
100 lines	Span/100	100/Span
200 lines	Span/200	200/Span
400 lines	Span/400	400/Span
800 lines	Span/800	800/Span

The two displays can have different number of Lines if the Analyzer Configuration is set to Independent Channels. In this configuration, no two channel measurements are allowed (frequency response, cross-spectrum, etc.), but the entry field can be linked or unlinked using the [Link] key. If Analyzer Configuration is set to Dual Channel they field is automatically linked to both displays.

Command: FLIN (?) d {, i}

## Base Frequency

Select the FFT Base Frequency [102.4 kHz, 100.0 kHz].

The Base Frequency sets the Full Span bandwidth for FFT and Correlation measurements. All spans are derived from the base frequency by dividing by powers of 2. 102.4 kHz is commonly used and has an advantage at narrow spans where even bin frequencies of 1 Hz, 0.5 Hz, etc. are available. 100.0 kHz has even bin frequencies of 1000 Hz, 500 Hz, etc. at wide spans.

Changing the FFT Base Frequency affects ALL of the FFT frequency parameters as well as the source frequency.

All FFT measurements on BOTH displays use the same Base Frequency.

Command: FBAS (?) d {, i}

## Start Frequency

Set the Start frequency of the active display FFT span. The Start frequency is the lowest frequency in the measurement span. The knob adjusts the Start frequency in steps equal

## 4-10 FFT Frequency Menu

---

to the Linewidth. If an entered Start frequency is incompatible with the Span, then it will be set to the closest allowed value.

This softkey also anchors the Start frequency. Adjusting the FFT Span will leave the Start frequency unchanged. Enlarging the frequency Span may change the Start frequency if the new measurement would extend below 0 Hz or above 102.4 (100.0) kHz.

If the Measurement Group is Correlation, then the Start frequency is always 0 Hz (baseband).

The two displays can have different Start Frequencies if the Analyzer Configuration is set to Independent Channels. In this configuration, no two channel measurements are allowed (frequency response, cross-spectrum, etc.), but the entry field can be linked or unlinked using the [Link] key. If Analyzer Configuration is set to Dual Channel they field is automatically linked to both displays.

Command: FSTR (?) d {, f}

---

## Center Frequency

Set the Center frequency of the active display FFT span. The Center frequency is the center of the measurement span. The knob adjusts the Center frequency in steps equal to the Linewidth. If an entered Center frequency is incompatible with the Span, then it will be set to the closest allowed value.

This softkey also anchors the Center frequency. Adjusting the FFT Span will leave the Center frequency unchanged. Enlarging the frequency Span may change the Center frequency if the new measurement would extend below 0 Hz or above 102.4 (100.0) kHz.

If the Measurement Group is Correlation, then the Center frequency is always equal to Span/2 (baseband).

The two displays can have different Center Frequencies if the Analyzer Configuration is set to Independent Channels. In this configuration, no two channel measurements are allowed (frequency response, cross-spectrum, etc.), but the entry field can be linked or unlinked using the [Link] key. If Analyzer Configuration is set to Dual Channel they field is automatically linked to both displays.

Command: FCTR (?) d {, f}

---

## End Frequency

Set the End frequency of the active display FFT span. The End frequency is the highest frequency in the measurement span. The knob adjusts the End frequency in steps equal to the Linewidth. If an entered End frequency is incompatible with the Span, then it will be set to the closest allowed value.

This softkey also anchors the End frequency. Adjusting the FFT Span will leave the End frequency unchanged. Enlarging the frequency Span may change the End frequency if the new measurement would extend below 0 Hz or above 102.4 (100.0) kHz.

If the measurement group is Correlation, then the End frequency is always equal to the Span (baseband).

The two displays can have different End Frequencies if the Analyzer Configuration is set to Independent Channels. In this configuration, no two channel measurements are allowed (frequency response, cross-spectrum, etc.), but the entry field can be linked or unlinked using the [Link] key. If Analyzer Configuration is set to Dual Channel they field is automatically linked to both displays.

Command: FEND (?) d {, f}

## Octave Frequency Menu

When the Measurement Group is Octave Analysis, this menu selects the number and spacing of frequency bands for the active display. See ‘Octave Analysis’ in Chapter 2 for a discussion about octave analysis fundamentals.

Frequency
Highest Band 32 kHz
Lowest Band 250 Hz
Octave Resolution Full
Octave Channels 1 Channel

---

### Highest Band

Set the Highest Band for the active display. The limits for the Highest Band are determined by the Octave Resolution and the number of Octave Channels.

#### 1 Octave Channel

- 1/1 Octave: [125 mHz .. 32 kHz]
- 1/3 Octave: [160 mHz .. 40 kHz]
- 1/12 Octave: [190 mHz .. 12.34 kHz]

#### 2 Octave Channels

- 1/1 Octave: [125 mHz .. 16 kHz]
- 1/3 Octave: [160 mHz .. 20 kHz]
- 1/12 Octave: [190 mHz .. 6.17 kHz]

The total power of all measured bands is displayed as a filled in band at the right side of the graph.

The exact band center frequencies are calculated according to the ANSI standard. The displayed frequencies are sometimes rounded to even values for a simpler display.

The measurement may span 1 to 11 octaves. Increasing the Highest Band will also increase the Lowest Band if necessary to keep the measurement span at 11 octaves.

The two displays can have different Highest Bands if the Analyzer Configuration is set to Independent Channels. If Analyzer Configuration is set to Dual Channel they field is automatically linked to both displays.

Command: OHIB (?) d {, f}



## Lowest Band

Set the Lowest Band for the active display. The limits for the Lowest Band are determined by the Octave Resolution.

1/1 Octave: [125 mHz ..]  
1/3 Octave: [100 mHz ..]  
1/12 Octave: [100 mHz ..]

The exact band center frequencies are calculated according to the ANSI standard. The displayed frequencies are the ANSI nominal center frequencies.

The measurement may span 1 to 11 octaves. Decreasing the Lowest Band will also decrease the Highest Band if necessary to keep the measurement span at 11 octaves.

The two displays can have different Lowest Bands if the Analyzer Configuration is set to Independent Channels. If Analyzer Configuration is set to Dual Channel they field is automatically linked to both displays.

Command: OLOB (?) d {, f}

---

## Octave Resolution

Select the Octave Resolution (number of bands per octave) for the active display [Full, Third, Twelfth].

Full octave band centers are calculated by multiplying or dividing 1000 Hz by 2.

Third octave band centers are calculated by multiplying or dividing 1000 Hz by  $2^{1/3}$  or 1.2599.

Twelfth octave band centers are calculated by multiplying or dividing 1000 Hz by  $2^{1/12}$  or 1.0595.

The exact band center frequencies are calculated according to the ANSI standard. The displayed frequencies are sometimes rounded to even values for a simpler display. The filter shapes are third-order Butterworth with full, 1/3 or 1/12 octave bandwidth.

Changing the Octave Resolution will change the Lowest Band and Highest Band to allowed bands of the new resolution.

The two displays can have different Octave Resolution if the Analyzer Configuration is set to Independent Channels. If Analyzer Configuration is set to Dual Channel they field is automatically linked to both displays.

Command: ORES (?) d {, i}

### Octave Channels

Select the number of input channels for the octave measurement. The number of Octave Channels affects the Measurements of both displays.

If 2 Octave Channels are selected, then both inputs may be analyzed. The Measurements of both displays are independent. The Highest Band for all Octave Resolutions is decreased by 2 in this case.

If 1 Octave Channel is selected, then only one input may be analyzed. If both displays are making an octave measurement, they must both use the same input. Changing the input of the active display will change the Measurement or input of the other display if necessary. User Functions which use both inputs may not be measured. Choosing 1 Octave Channel may change the current Measurements so that both displays use the same input.

Command: OCHN (?) d {i}

# Swept Sine Frequency Menu

When the Measurement Group is Swept Sine, this menu sets the frequency sweep parameters. These parameters govern the measurements on both displays. See 'Swept Sine' in Chapter 2 for a discussion about swept sine measurement fundamentals.

Frequency	
Start	1 kHz
Stop	20 kHz
Repeat	Continuous
Type	Linear
Auto Resolution	On
# of Points	200
Max Step Size	8
Faster Threshold	< 0.1 dB
Slower Threshold	> 0.5 dB

---

## Start

Set the Start frequency [1 mHz .. 102.4 kHz]. Note that measurements at frequencies less than 1 Hz take a significant amount of time.

The Start frequency is usually less than the Stop frequency. If the Start frequency is greater than the Stop frequency, then the sweep proceeds downward from the Start frequency.

If the Start frequency is changed during a sweep, the sweep will be reset.

Command: SSTR d (,?) {f}

---

## Stop

Set the Stop frequency [1 mHz .. 102.4 kHz]. Note that measurements at frequencies less than 1 Hz take a significant amount of time.

The Stop frequency is usually greater than the Start frequency. If the Stop frequency is less than the Start frequency, then the sweep is downward from the Start frequency.

If the Stop frequency is changed during a sweep, the sweep will be reset.

Command: SSTP (?) d {, f}

### Repeat

Select the Sweep Repeat Mode [Single Shot, Continuous].

In Single Shot mode, the measurement is paused at the completion of the sweep and the source ramps off.

In Continuous mode, the measurement is repeated at the completion of each sweep. The source moves immediately to the Start frequency and begins the sweep again. Data from the previous sweep is preserved on the display and can be examined with the cursor.

Command: SRPT (?) d {, i}

---

### Type

Select the Sweep Type [Linear, Log].

Linear sweep computes the measurement points in a linear progression from the start to the stop frequency. Choosing linear sweep also sets the X axis to linear.

Log sweep uses a logarithmic progression of points from the start to the stop frequency. Choosing log sweep also sets the X axis to logarithmic.

Command: SSTY (?) d {, i}

---

### Auto Resolution

Select Auto Resolution [Off, On].

If Auto Resolution is Off, all points in the sweep are measured.

Auto Resolution On allows the sweep to skip points if sequential points do not vary by more than a specified amount. Auto Resolution is specified by three parameters, the Faster Threshold, the Slower Threshold and the Maximum Step Size. These parameters can only be adjusted when Auto Resolution is On.

Auto Resolution examines the results of successive measurements. If the newest measurement is within the Faster Threshold of the previous measurement (for BOTH channels), then the sweep will take larger steps, skipping frequency points. Each successive time this threshold is met, the step size is increased until the Maximum Step Size is reached. This speeds up the sweep in regions where the response is flat (varies less than the Faster Threshold).

If a measurement differs from the previous measurement by more than the Slower Threshold (for EITHER channel), then the sweep returns to the previously measured point and moves to the very next frequency point in the sweep (with no skipping). The sweep continues from this point, speeding up if allowed and slowing down when required. This 'fills' in skips in the sweep which vary by more than the Slower Threshold.

Measurements which differ by more than the Faster Threshold (on EITHER channel) but less than the Slower Threshold (on BOTH channels), maintain the present sweep speed. The number of points skipped remains the same in this case.

If the marker position is displayed with a '?', the point has not actually been measured but is interpolated from actual measured points. The measurement of these points was skipped due to sweep Auto Resolution.

Command: SARS (?) d {, i}

---

## Number of Points

Set the Number Of Points [10..2048]. The points are in a Linear or a Logarithmic progression as set by the Sweep Type.

If the Number Of Points is changed during a sweep, the sweep will be reset.

A sweep with a large number of points can detect narrow features. A sweep with a small number of points will take less time. In order to save time while maintaining a high resolution, use Auto Resolution.

Command: SNPS (?) d {, i}

---

## Maximum Step Size

Set the Maximum Step Size for Auto Resolution frequency sweeps [2..256].

When Auto Resolution is On, each successive time the Faster Threshold condition is met (on BOTH channels), the frequency step size is increased until the Maximum Step Size is reached. This sets the maximum speed at which the sweep will continue until the Slower Threshold is exceeded.

Generally, this number should not exceed 5% of the Number Of Points in the sweep.

Command: SSKP (?) d {, i}

---

## Faster Threshold

Set the Faster Threshold [0.01..3 dB].

Auto Resolution examines the results of successive measurements. If the newest measurement is within the Faster Threshold of the previous measurement for BOTH channels, then the sweep will take larger steps, skipping frequency points. Each successive time this threshold is met, the step size is increased until the Maximum Step Size is reached. This speeds up the sweep in regions where the response is flat (varies less than the Faster Threshold).

Measurements which differ by more than the Faster Threshold on EITHER channel, but less than the Slower Threshold on BOTH channels, maintain the present sweep speed.

## 4-18 Swept Sine Frequency Menu

---

Generally, the Faster Threshold should be set to the desired amplitude resolution. In order for skipping to occur, this should be something greater than the variations in the flattest regions of the measurement.

Command: SFST (?) d {, i}

---

### Slower Threshold

Set the Slower Threshold [0.05..6 dB].

If Auto Resolution is skipping points and a measurement differs from the previous measurement by more than the Slower Threshold for EITHER channel, then the sweep returns to the previously measured point and moves to the very next frequency point in the sweep (with no skipping). The sweep continues from this point, speeding up if allowed and slowing down when required. This 'fills' in skips in the sweep which vary by more than the Slower Threshold.

Measurements which differ by more than the Faster Threshold on EITHER channel, but less than the Slower Threshold on BOTH channels, maintain the present sweep speed.

Generally, the Slower Threshold should be set to less than half of the smallest desired feature size (relative to the region before the feature). A good rule is to set the Slower Threshold at twice the Faster Threshold.

Command: SSLO (?) d {, i}

## Order Frequency Menu

When the Measurement Group is Order, this menu sets parameters related to which orders the SR785 will compute and track. These parameters govern the measurements on both displays. See ‘Order Analysis’ in Chapter 2 for a discussion about order measurement fundamentals.

Frequency	
Max RPM	6e+004
Min RPM	30
Max Order	10
Delta Order	0.1
Tracking	Off
Track Setup	

---

### Max RPM

Set the maximum rotational speed allowed during order measurements. The maximum rpm can take on a value between that specified by Min RPM and a value calculated from the Max Order by the formula:  $\max(\text{Max RPM}) = 60 * 40,000 / \text{Max Order}$ . Selecting a lower Max RPM will allow calculation of higher orders.

Command: ORMX (?) d {, f}

---

### Min RPM

This parameter does not affect the way order measurements are calculated by the SR785. It merely provides a default lower value by which to scale measurements whose x-axis is RPM, such as RPM profiles and Tracked Orders.

Command: ORMN (?) d {, f}

---

### Max Order

Specifies the maximum order, i.e. multiple of the rotational speed, calculated by the SR785. Allowed values are between 3 and 200. Choosing a higher Max Order will limit the allowed range of Max RPM.

Command: OMAX (?) d {, f}

### Delta Order

Sets the resolution for order spectra. The value may range between Max Order divided by 400, and 1.0. Choosing a smaller value for Delta Order (higher resolution) increases the time it takes to acquire an order spectrum.

Command: ODLT (?) d {, f}

---

### Tracking

Turns order tracking on and off. When tracking is On, the Track1 and Track2 measurements become available in the order measurement group. The Track measurements are plots of the amplitude of one particular order as a function of RPM.

Because track measurements use the memory allocated to the waterfall buffer, you cannot use waterfall storage while Tracking is on.

Command: OTRK (?) d {, i}

---

### Track Setup

Displays the track setup submenu. This submenu displays parameters which define the Track1 and Track2 measurements.

The screenshot shows a vertical menu titled "Track Setup". The menu items are: "# Track Pts" with a value of 400; "Storage" with a value of Continuous; "Track 1 Order" with a value of 1; "Input Channel" with a value of Ch 1; "Track 2 Order" with a value of 1; "Input Channel" with a value of Ch 2; and "Return" with a right arrow icon.

---

### Track Points

Sets the number of points for track measurements. When the specified number of points have been acquired the measurement will either stop, or wrap around depending on the setting of the track storage softkey.

A track measurement is a two-dimensional plot of the amplitude of a particular order vs. rpm with time as the parameter. Individual points in the measurement can be taken continuously, or at specified time or rpm intervals using Time arming and RPM Arming.

Command: ONPT (?) d {, i}



**Storage**

Determines the storage mode for track measurements. If Storage is set to One Shot, the track measurement will end after the specified number of points have been acquired. If Storage is set to continuous, the measurement will not stop, but the new data will overwrite the oldest data.

Command: OSTO (?) d {, i}

---

**Track 1 Order**

Specifies the order associated with the Track 1 measurement. If the entered value is not an integer multiple of delta order, the SR785 will display the closest available order to the entered value.

Regardless of the order specified, the SR785 always tracks all orders specified by the max order and delta order softkeys. The Track 1 Order softkey merely selects which order will be displayed. Once a track measurement is done or paused, any order track can be displayed without retaking any data.

Command: O1TK (?) d {, f}

---

**Input Channel (Track 1)**

Specifies the input channel associated with the Track 1 measurement.

Regardless of the input channel specified, the SR785 always tracks all orders specified by the max order and delta order softkeys for both input channels. The Input Channel (Track 1) softkey merely selects which channel will be displayed. Once a track measurement is done or paused, data for each channel can be displayed without retaking any data.

Command: O1BN (?) d {, i}

---

**Track 2 Order**

Specifies the order associated with the Track 2 measurement. If the entered value is not an integer multiple of delta order, the SR785 will display the closest available order to the entered value.

Regardless of the order specified, the SR785 always tracks all orders specified by the max order and delta order softkeys. The Track 2 Order softkey merely selects which order will be displayed. Once a track measurement is done or paused, any order track can be displayed without retaking any data.

Command: O2TK (?) d {, f}

---

### Input Channel (Track 2)

Specifies the input channel associated with the Track 2 measurement.

Regardless of the input channel specified, the SR785 always tracks all orders specified by the max order and delta order softkeys for both input channels. The Input Channel (Track 2) softkey merely selects which channel will be displayed. Once a track measurement is done or paused, data for each channel can be displayed without retaking any data.

Command: O2BN (?) d {, i}

## Time/Histogram Frequency Menu

When the Measurement Group is Time/Histogram, this menu sets parameters related to the resolution and number of points in the histograms. These parameters govern the measurements on both displays. See 'Time/Histogram Analysis' in Chapter 2 for a discussion about time/histogram measurement fundamentals.

Frequency	
Sample Time	3.81 $\mu$ s
Histogram Length	3000 samples
Bins	256
Base Freq.	102.4 kHz
Repeat	On

---

### Sample Time

Sets the sampling interval for time records and histograms in the Time/Histogram measurement group. The allowed sampling times depend on the setting of the base frequency softkey. If the base frequency is 102.4 kHz the allowed sampling times range from 3.81 $\mu$ s to 2s in powers of 2. If the base frequency is 100.0 kHz the allowed sampling times range from 3.91 $\mu$ s to 2.048s in powers of 2.

Command: SPAN (?) d {, i}

---

### Histogram Length

Sets the number of points in each histogram. The number of points can be specified in samples, records, (each record is 1024 samples), or in seconds. When the specified number of samples, records, or seconds, has been reached, the histogram will update.

If repeat is on the SR785 will begin accumulating a new histogram at the end of the prior histogram.

Command: HLEN (?) d {, f}

---

### Bins

Sets the number of histogram 'bins' along the x-axis. Values from 4 to 1024 in power of 2 increments are allowed. The x-axis display for histogram measurements ranges from -1.25 times the input range to +1.25 times the input range. The resolution of the histogram is the therefore 2.5 times the input range divided by the number of bins.

Command: HBIN (?) d {, i}

## 4-24 Time/Histogram Frequency Menu

---

### Base Frequency

This softkey sets the allowed values of the sampling time. If set to 100 kHz the fundamental sampling rate of the instrument is 256 kHz and the allowed sampling times will be multiples of the reciprocal of this frequency (3.91 $\mu$ s). If set to 102.4 kHz the the fundamental sampling rate of the instrument is 262.144 kHz and the allowed sampling times will be multiples of the reciprocal of this frequency (3.814 $\mu$ s).

Command: FBAS (?) d {, i}

---


### Repeat

When repeat is set 'off', the SR785 computes a single histogram and stops. If repeat is on, the analyzer will automatically acquire new histograms continually as each one is finished. This can be useful for analyzing signals which change slowly with time.

Command: HRPT (?) d {, i}\\$

# Display Setup Menu

The Display Setup menu allows the user to choose the Measurement, View and Display Scaling for the active display.

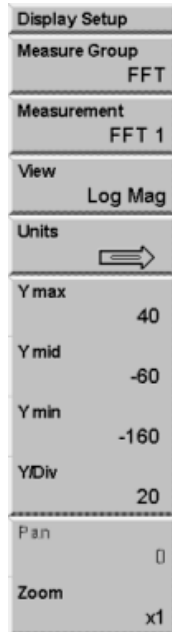

## Measurement Group

Select the Measurement Group for both displays [FFT, Correlation, Octave Analysis, Swept Sine, Order, Time/Histogram]. To choose a specific measurement, choose the Measurement Group, then select a Measurement.

The FFT Measurement Group uses the Fast Fourier Transform (FFT) to compute the frequency spectrum of finite time records.

The Correlation Measurement Group calculates the time domain auto and cross correlation of signals.

The Octave Analysis Measurement Group measures the power within logarithmically spaced frequency bands covering up to 11 octaves.

The Swept Sine Measurement Group uses the sine source to sweep a frequency range making optimized measurements at discrete frequency points along the way.

The Order Measurement Group makes spectrum and time measurements which are synchronized to the rotation of a rotating machine.

The Time/Histogram Measurement Group calculates histograms and statistical parameters of the input signal.

The choice of Measurement Group changes the [Frequency], [Average] and [Source] menus.

Command: MGRP(?) d {, i}

---

### Measurement (FFT)

Select the Measurement of the active display when the Measurement Group is FFT.

Each Measurement has an associated View. Changing the Measurement changes the View to the View last used with the new Measurement.

The two displays generally have different Measurements. This entry field can be linked to both displays using the [Link] key.

Measurements which are listed as '2 channel' may only be selected if <Analyzer Config> is set to 'Dual Channel.'

Command: MEAS(?) d {, i}

The following Measurements are available in the FFT Measurement Group:

#### FFT (Linear Spectrum)

The FFT spectrum is the basic measurement of an FFT analyzer. It is simply the FFT of a time record. The spectrum is a complex quantity (it contains phase as well as amplitude information). This is sometimes referred to as the linear spectrum.

The phase of the spectrum is meaningful only if the time record is triggered with a fixed relationship to the input signal. If the signal is periodic as well (the signal and trigger repeat), then vector averaging can be used to reduce the noise level of the spectrum. The vector averaged spectrum is still a complex quantity

The precise definition of the FFT1 measurement for all averaging modes is as follows:

#### No Average

$$\text{FFT 1} = \text{FFT1}$$

#### Vector Average

$$\text{FFT 1} = \langle \text{FFT1} \rangle$$

#### RMS Average

$$\text{FFT 1} = \sqrt{\langle \text{FFT1} * \bullet \text{FFT1} \rangle}$$

#### Peak Hold Average

$$\text{FFT 1} = \sqrt{\text{MAX}(\text{FFT1} * \bullet \text{FFT1})}$$

#### Power Spectrum

The Power Spectrum is derived from the FFT spectrum by multiplying the spectrum by its complex conjugate. The averaged power spectrum is a good approximation to the rms

signal and noise amplitudes. The power spectrum is a real quantity and contains no phase information.

The precise definition of the Power Spectrum 1 measurement for all averaging modes is as follows:

**No Average**

$$\text{Power Spectrum1} = \text{FFT1} * \text{FFT1}$$

**Vector Average**

$$\text{Power Spectrum1} = \langle \text{FFT1} \rangle * \langle \text{FFT1} \rangle$$

**RMS Average**

$$\text{Power Spectrum1} = \langle \text{FFT1} * \text{FFT1} \rangle$$

**Peak Hold Average**

$$\text{Power Spectrum1} = \text{MAX}(\text{FFT1}) * \text{MAX}(\text{FFT1})$$

**Time Record**

A time record is simply a sequence of data samples. The length of the time record is the FFT resolution/span.

For baseband spans (spans which start at DC), the time record resembles a digital oscilloscope display. Signals at frequencies above the span have been filtered out. Baseband time records are entirely real, they have no imaginary part.

For zoomed spans (spans which start above DC), the time record does NOT resemble the original data. The data has been frequency shifted. Signals at the center of the span appear at DC while frequencies at both edges of the span appear as high frequencies. Zoomed time records are complex, they have both a real and an imaginary part. The sampling rate is always half of the equivalent baseband span.

Remember, the time record is not a continuous representation of the input signal. The data is sampled and has a time resolution of  $1/(\text{sample rate})$ . High frequency signals will appear distorted in the time record. However, ALL of the spectral information is preserved by the Nyquist sampling theorem as long as the value of each sample is accurate. A triggered time record will always jitter by 1 sample. This jitter is removed in the computation of the phase of the spectrum relative to the trigger.

Averaging does not affect the time record. Averaging is performed on the FFT spectrum and not on the time data.

Amplitude calibration is also performed in the frequency domain. Hence, the time data amplitudes are not calibrated.

**Windowed Time Record**

The FFT operates on windowed time records. The window function is applied to the time record immediately before the FFT. Most window functions taper off to zero at the start and end of the time record. If a transient signal occurs at the start of the time record, the corresponding windowed time record and FFT may not show anything because the window function reduces the transient to zero.

### Orbit

The Orbit measurement is a two channel measurement whose real part is the real part of Time Record Ch1 and whose imaginary part is the real part of Time Record Ch2. Orbit is normally displayed with the Nyquist View (Time2 vertical vs Time1 horizontal).

For baseband spans, the time records are entirely real and the Nyquist view of the Orbit measurement is a Lissajous figure.

### Cross Spectrum

The cross spectrum (sometimes called cross power spectrum) is a two channel measurement defined as

#### No Average

$$\text{Cross Spectrum} = \text{FFT1}^* \cdot \text{FFT2}$$

#### Vector Average

$$\text{Cross Spectrum} = \langle \text{FFT1}^* \rangle \cdot \langle \text{FFT2} \rangle$$

#### RMS Average

$$\text{Cross Spectrum} = \langle \text{FFT1}^* \cdot \text{FFT2} \rangle$$

#### Peak Hold Average

$$\text{Cross Spectrum} = \text{MAX}(\text{FFT2})^* \cdot \sqrt{\langle \text{FFT1}^* \cdot \text{FFT1} \rangle}$$

The cross spectrum contains both magnitude and phase information. The phase is the relative phase (at each frequency) between the two channels. Vector averaging can be used to eliminate signals which do not have a constant phase relationship between the two channels. In this case, triggering is not required for vector averaging.

The magnitude is simply the product of the magnitudes of each spectrum. Frequencies where signal is present in both spectra will have large components in the cross spectrum.

### Frequency Response

The Frequency Response (sometimes called transfer function) is a two channel measurement which ratios the spectrum of Ch 2 to the spectrum of Ch 1. Frequency Response measures the response of a network or device under test. The reference channel (1) measures the signal at the input to the device and the response channel (2) measures the device output. The result is the complex Frequency Response of the device.

The definition of frequency response depends on the type of averaging which is displayed.

#### No Average

$$\text{Freq. Response} = \text{FFT2} / \text{FFT1}$$

#### Vector Average

$$\text{Freq. Response} = \langle \text{FFT2} \rangle / \langle \text{FFT1} \rangle$$

#### RMS Average

$$\text{Freq. Response} = (\text{RMSAvg cross spectrum}) / \text{power spectrum 1}$$



$$\text{Freq. Response} = \langle \text{FFT1}^* \cdot \text{FFT2} \rangle / \langle \text{FFT1}^* \cdot \text{FFT1} \rangle$$

### Peak Hold Average

$$\text{Freq. Response} = \text{MAX}(\text{FFT2}) / \sqrt{\langle \text{FFT1}^* \cdot \text{FFT1} \rangle}$$

Both the RMS averaged and Vector averaged frequency response contains both magnitude and phase information. The phase is the relative phase (at each frequency) between the two channels. The RMS averaged frequency response is computed by taking the ration of the cross spectrum to the input power spectrum, a technique called the tri-spectral average.

### Coherence

The coherence function is a two channel measurement defined as

$$\text{Coherence} = \text{Mag}^2(\text{RMSAvg}(\text{CrossSpec})) / (\text{Pwr1} \cdot \text{Pwr2})$$

Averaging is always On and the Averaging Modes are defined by the measurement above. The Type and Number Of Averages are still selected in the [Average] menu.

Coherence is a unitless real quantity and varies from 1.0 (perfect coherence) to 0.0 (no coherence). Coherence measures the percentage of power in the response channel (2) which is caused by (phase coherent with) power in the reference or input channel (1). Ideally, a coherence of 1.0 means that the corresponding Frequency Response is completely legitimate. All of the response power came from power at the input. If there is noise or other signals generated from within the device under test (which is not related to the input signal), it will result in a coherence of less than 1.0.

### Capture

The capture buffer stores sequential time domain data in memory. See 'Capture' in Chapter 2 for more details. The capture measurement displays the contents of the capture buffer. FFT measurements can use this stored data as input data by choosing Playback as the Input Source in the [Input] menu.

The capture buffer often stores many time records. To graphically expand a region of the display, use the Pan and Zoom functions in the [Display Setup] menu. The capture buffer display will automatically pan as the capture fill and playback progress through the buffer. During capture fill, if the capture buffer contains more points than can be displayed, points are skipped. This speeds up the display update so that it keeps up with the real time capture but allows visual aliasing to occur. Once capture is complete, the display is redrawn showing the envelope of all points, eliminating any visual aliasing effects.

To measure from a region of the buffer, set the Playback Start and Length in the [Input] menu.

The capture data is filtered and down-sampled according to the capture Sample Rate. Only baseband data (data bandwidth starts at DC) are captured. The capture buffer resembles a digital oscilloscope display. Signals at frequencies above 1/2.56 times the sample rate have been filtered out.

The capture buffer is not a continuous representation of the input signal. The data is sampled and has a time resolution of  $1/(\text{sample rate})$ . High frequency signals will appear distorted in the time record. However, ALL of the spectral information (up to  $1/2.56$  times the Sample Rate) is preserved by the Nyquist theorem as long as the value of each sample is accurate.

Amplitude calibration is performed in the frequency domain. Hence, the captured time data amplitudes are not calibrated.

### User Functions

A User Function displays the results of a user defined math function. User Functions defined within the FFT Measurement Group may include FFT measurement results. Use the [User Math] menu to define a function.

A User Function may not be selected as the measurement if it uses a Trace which is empty or which contains data taken in a different Measurement Group.

See 'User Math Functions' in Chapter 2 for more.

---

## Measurement (Correlation)

Select the Measurement of the active display when the Measurement Group is Correlation.

Each Measurement has an associated View. Changing the Measurement changes the View to the View last used with the new Measurement.

Command: MEAS(?) d {, i}

The following Measurements are available in the Correlation Measurement Group:

### Auto Correlation

Auto correlation is a single channel measurement. In the time domain, it is a comparison of a signal  $x(t)$  with a time shifted version of itself  $x(t-\tau)$  displayed as a function of  $\tau$ . This is useful for detecting similarities which occur at different times. For example, echoes show up as peaks separated by the echo time. Sine waves appear as sine waves and square waves appear as triangles. Signals which do not repeat or are completely random for all time appear only at  $\tau = 0$ .

### Averaging

Auto Correlation =  $\text{invFFT}(\text{FFTuN} \cdot \text{FFTN}^*)$

### RMS or Peak Hold Averaging

Auto Correlation =  $\text{invFFT}(\langle \text{FFTuN} \cdot \text{FFTN}^* \rangle)$

### Vector Averaging On

Auto Correlation =  $\text{invFFT}(\langle \text{FFTuN} \rangle \cdot \langle \text{FFTN}^* \rangle)$

where N is Channel 1 or 2. FFT is a windowed FFT, FFTu is an un-windowed FFT (uniform window) and invFFT is an inverse FFT.

Correlation is a real function and requires a baseband span (real time record). Non-baseband time records do not preserve the original signal frequencies and thus do not yield the correct correlation. The Start Frequency is always set to 0 in the Correlation Measurement Group.

A window is applied to the time record of the first FFT in the computation. This is because the FFT models the time domain as a single time record repeating itself over and over. Computing the correlation over a  $\tau$  greater than half of the time record length will result in 'wrap around' errors where data starts to repeat itself. To avoid this, special windows which zero half of the time record are used. The  $[0..T/2]$  window zeroes the second half of the time record and the  $[-T/4..T/4]$  window zeroes the first and last quarter of the time record. The  $[-T/2..T/2]$  is a uniform window which should only be used on data which is self windowing (lasts less than half of the time record).

### Cross Correlation

Cross correlation is a two channel measurement which is only available when the [Input]<Analyzer Conig> softkey is set to 'Dual Channel.' In the time domain, it is a comparison of a signal  $x(t)$  with a time shifted version of another signal  $y(t-\tau)$  displayed as a function of  $\tau$ . This is useful for measuring time delays between two common signals.

The definition of Cross Correlation depends upon the Average Mode selected in the [Average] menu.

#### Averaging

Cross Correlation =  $\text{invFFT}(\text{FFT}_{u2} \cdot \text{FFT}_{1*})$

#### Peak Hold or RMS Averaging On

Cross Correlation =  $\text{invFFT}(\langle \text{FFT}_{u2} \cdot \text{FFT}_{1*} \rangle)$

#### Vector Averaging On

Cross Correlation =  $\text{invFFT}(\langle \text{FFT}_{u2} \rangle \cdot \langle \text{FFT}_{1*} \rangle)$

where FFT1 is the windowed FFT of Channel 1 and FFTu2 is the un-windowed FFT (uniform window) of Channel 2.

Correlation is a real function and requires a baseband span (real time record). Non-baseband time records do not preserve the original signal frequencies and thus do not yield the correct correlation. A display which is measuring a single channel correlation will have its start frequency set to 0 Hz. A display which measures cross correlation will set the start frequency of both displays to 0 Hz.

A correlation window is applied to the time record of the first FFT in the computation. This is because the FFT models the time domain as a single time record repeating itself over and over. Computing the correlation over a  $\tau$  greater than half of the time record length will result in 'wrap around' errors where data starts to repeat itself. To avoid this, special windows which zero half of the time record are used. The  $[0..T/2]$  window zeroes the second half of the time record and the  $[-T/4..T/4]$  window zeroes the first and last quarter of the time record. The  $[-T/2..T/2]$  is a uniform window which should only be used on data which is self windowing (lasts less than half of the time record).

### Capture

The capture buffer stores sequential time domain data in memory. See 'Capture' in Chapter 2 for more details. The capture measurement displays the contents of the capture buffer. FFT measurements can use this stored data as input data by choosing Playback as the Input Source in the [Input] menu.

The capture buffer often stores many time records. To graphically expand a region of the display, use the Pan and Zoom functions in the [Display Setup] menu. The capture buffer display will automatically pan as the capture fill and playback progress through the buffer. During capture fill, if the capture buffer contains more points than can be displayed, points are skipped. This speeds up the display update so that it keeps up with the real time capture but allows visual aliasing to occur. Once capture is complete, the display is redrawn showing the envelope of all points, eliminating any visual aliasing effects.

To measure from a region of the buffer, set the Playback Start and Length in the [Input]< menu.

The capture data is filtered and down-sampled according to the capture Sample Rate. Only baseband data (data bandwidth starts at DC) are captured. The capture buffer resembles a digital oscilloscope display. Signals at frequencies above 1/2.56 times the sample rate have been filtered out.

The capture buffer is not a continuous representation of the input signal. The data is sampled and has a time resolution of 1/(sample rate). High frequency signals will appear distorted in the time record. However, ALL of the spectral information (up to 1/2.56 times the Sample Rate) is preserved by the Nyquist theorem as long as the value of each sample is accurate.

Amplitude calibration is performed in the frequency domain. Hence, the captured time data amplitudes are not calibrated.

### User Functions

A User Function displays the results of a user defined math function. User Functions defined within the Correlation Measurement Group may include correlation measurement results. Use the [User Math] menu to define a function.

A User Function may not be selected as the measurement if it uses a Trace which is empty or which contains data taken in a different Measurement Group.

See 'User Math Functions' in Chapter 2 for more.

---

### Measurement (Octave Analysis)

Select the Measurement of the active display when the Measurement Group is Octave Analysis.

Each Measurement has an associated View. Changing the Measurement changes the View to the View last used with the new Measurement.

Command: MEAS(?) d {, i}

The following Measurements are available in the Octave Analysis Measurement Group:

### **Octave Spectrum**

The octave spectrum is the basic measurement of octave analysis. It is simply the rms averaged outputs from the parallel bank of octave filters (also called the power spectrum). The spectrum gives a stable reading of the rms signal amplitudes and noise levels within each band. RMS averaging results in a real spectrum and there is no phase information.

### **Capture**

The capture buffer stores sequential time domain data in memory. See 'Capture' in Chapter 2 for more details. The Capture measurement displays the contents of the capture buffer. Octave measurements can use this stored data as input data by choosing Playback as the Input Source in the [Input] menu.

The capture buffer is often very long. To graphically expand a region of the display, use the Pan and Zoom functions in the [Display Setup] menu. The capture buffer display will automatically pan as the capture fill and playback progress through the buffer. During capture fill, if the capture buffer contains more points than can be displayed, points are skipped. This speeds up the display update so that it keeps up with the real time capture but allows visual aliasing to occur. Once capture is complete, the display is redrawn showing the envelope of all points, eliminating any visual aliasing effects.

To measure from a region of the buffer, set the Playback Start and Length in the [Input] menu.

The capture data is filtered and down-sampled according to the capture sample rate. Only baseband data (data bandwidth starts at DC) are captured. The capture buffer resembles a digital oscilloscope display. Signals at frequencies above  $1/2.56$  times the sample rate have been filtered out.

The capture buffer is not a continuous representation of the input signal. The data is sampled and has a time resolution of  $1/(\text{sample rate})$ . High frequency signals will appear distorted in the time record. However, ALL of the spectral information (up to  $1/2.56$  times the Sample Rate) is preserved by the Nyquist theorem as long as the value of each sample is accurate.

Amplitude calibration is performed in the frequency domain. Hence, the captured time data amplitudes are not calibrated.

### **User Functions**

A User Function displays the results of a user defined math function. User Functions defined within the Octave Analysis Measurement Group may include octave measurement results. Use the [User Math] menu to define a function.

A User Function may not be selected as the measurement if it uses a Trace which is empty or which contains data taken in a different Measurement Group.

See 'User Math Functions' in Chapter 2 for more.

---

### Measurement (Swept Sine)

Select the Measurement of the active display when the Measurement Group is Swept Sine.

Each Measurement has an associated View. Changing the Measurement changes the View to the View last used with the new Measurement.

Command: MEAS(?) d {, i}

The following Measurements are available in the Swept Sine Measurement Group:

#### Spectrum

The swept sine spectrum is simply the measurement of a single channel over a sweep. The spectrum is complex (it contains phase and amplitude information). The phase is relative to the source and is stable but arbitrary. Single channel phase is not generally meaningful.

The spectrum measures the actual signal level at the inputs. If source auto level is On, then the spectrum will tend to be constant. Use Frequency Response to remove the effects of a changing source level.

#### Normalized Variance

Normalized variance is a measure of the signal to noise ratio of a swept sine measurement after the signal has been integrated for a specified time, or a specified number of cycles. Variance values near 1 indicate that the final signal to noise ratio is high, while values near zero indicate that the final signal to noise ratio is poor.

#### Cross Spectrum

The swept sine cross spectrum is a two channel measurement defined as

$$\text{Cross Spectrum} = \text{conj}(\text{Spec1}) \cdot \text{Spec2}$$

The cross spectrum contains both magnitude and phase information. The phase is the relative phase (at each frequency) between the two channels.

The magnitude is simply the product of the magnitudes of each spectrum. Frequencies where signal is present in both spectra will have large components in the cross spectrum.

#### Frequency Response

The swept sine Frequency Response (sometimes called frequency response) is a two channel measurement defined as

$$\text{Frequency Response} = \text{Spec2} / \text{Spec1}$$

The Frequency Response contains both magnitude and phase information. The phase is the relative phase (at each frequency) between the two channels.

Frequency Response measures the response of a network or device under test. The reference channel (1) measures the signal at the input to the device and the response channel (2) measures the device output. The result is the complex Frequency Response of the device.

### **User Functions**

A User Function displays the results of a user defined math function. User Functions defined within the Swept Sine Measurement Group may include swept sine measurement results. Use the [User Math] menu to define a function.

A User Function may not be selected as the measurement if it uses a Trace which is empty or which contains data taken in a different Measurement Group.

See 'User Math Functions' in Chapter 2 for more.

---

## **Measurement (Order)**

### **Linear Spectrum**

The Order Linear Spectrum displays the complex spectrum of the input signal as a function of order, e.g. multiples of the rotation frequency. This complex measurement contains both magnitude and phase information.

The phase of the order linear spectrum is relative to the point on the machine shaft corresponding to the tachometer pulse. In the case of several tach pulses per revolution the choice of which tach pulse is arbitrary, but is consistent from record to record.

### **Power Spectrum**

The Order Power Spectrum displays the power in the input signal as a function of order (multiple of the rotation frequency.) It is derived from the Linear Spectrum by multiplying each point by its complex conjugate.

### **Time Record**

The order time record displays the amplitude of the input signal at equiangular sampling points around the machine shaft. Note that strictly speaking, this is not a '**time** record:' the time at which the shaft reached any particular angular depends on the rotational speed of the shaft and is communicated to the analyzer via the tachometer input.

The number of time record points in each shaft revolution is equal to 5.12 times the setting of the Max Order softkey.

In the order measurement group, the displayed time record will be averaged if averaging is on.

### **Windowed Time Record**

The order group Windowed Time Record shows the Time Record multiplied by the currently selected Window Function. Remember that window functions prevent events which occur at the beginning and end of the time record from appearing in spectral measurements. If this happens, use a uniform window, or adjust the trigger to move the event to the center of the time record.

### Track

The Order Track measurement is a two dimensional plot of the amplitude of a given order vs. RPM. The Track measurement can only be selected when Tracking is turned on in the [Frequency] menu. The order and input channel displayed by the track measurement is selected with the Track Setup submenu.

The amplitude recorded in the track is the Linear Spectrum, and contains both amplitude and phase information.

Since the tracking feature uses the memory allocated to the Waterfall/Order buffer, the waterfall buffer is not available for waterfall storage when tracking is enable.

### RPM Profile

The RPM profile measurement is a two-dimensional plot showing RPM on the x-axis and time on the y-axis.

### Orbit

The order Orbit measurement is a two channel measurement whose real part is the real part of Time Record Ch1 and whose imaginary part is the real part of Time Record Ch2. Orbit is normally displayed with the Nyquist View (Time2 vertical vs Time1 horizontal).

Orbit measurements are often used in the order measurement group for shaft balancing.

### Capture

The capture buffer stores sequential time domain data in memory. See 'Capture' in Chapter 2 for more details. The capture measurement displays the contents of the capture buffer. Order measurements can use this stored data as input data by choosing Playback as the Input Source in the [Input] menu. In the order measurement group, the capture buffer stores tachometer information as well as input data. Therefore, only a capture buffer which was taken while the instrument was in the order group can be used as a source of playback data for order measurements.

Note that in order group, the data in the capture buffer will not look like the data in the order time record. The capture buffer contains data sampled at fixed time intervals. The time record contains data which is processed to yield samples at equal positions around the shaft.

To graphically expand a region of the display, use the Pan and Zoom functions in the [Display Setup] menu. The capture buffer display will automatically pan as the capture fill and playback progress through the buffer. During capture fill, if the capture buffer contains more points than can be displayed, points are skipped. This speeds up the display update so that it keeps up with the real time capture but allows visual aliasing to occur. Once capture is complete, the display is redrawn showing the envelope of all points, eliminating any visual aliasing effects.

### User Functions

A User Function displays the results of a user defined math function. User Functions defined within the Order Measurement Group may include order measurement results. Use the [User Math] menu to define a function.



A User Function may not be selected as the measurement if it uses a Trace which is empty or which contains data taken in a different Measurement Group.

See 'User Math Functions' in Chapter 2 for more.

---

## Measurement (Time/Histogram)

### Histogram

This is the basic measurement in the Time/Histogram group. The display shows a histogram of the number of samples found in some length of the input signal versus the amplitude of the sample. The amplitude covers the entire input range of the analyzer.

### PDF

The PDF, or probability density function, is a normalized version of the histogram.. The value of the PDF for a given voltage  $V$  multiplied by a small voltage interval,  $\Delta V$ , is the probability that a sample of the signal will fall between  $V$  and  $V+\Delta V$ .

### CDF

The CDF, or cumulative density function, is the integral of the probability density function. The CDF at a given voltage is the probability that a sample of the input signal will be below that voltage.

### Time

The time record in the time/histogram measurement group shows the value of the input signal sampled at intervals specified by <Sample Time> softkey. Unlike the FFT time record, the histogram time record is not filtered as the sampling frequency is decreased and may contain aliased signals. While this is important for the frequency domain measurements in the FFT measurement group, it is not a consideration in time/histogram group where the goal is to represent the original signal as closely as possible.

### Capture

The capture buffer stores sequential time domain data in memory. See 'Capture' in Chapter 2 for more details. The capture measurement displays the contents of the capture buffer. Time/Histogram measurements can use this stored data as input data by choosing Playback as the Input Source in the [Input] menu.

The capture buffer often stores many time records. To graphically expand a region of the display, use the Pan and Zoom functions in the [Display Setup] menu. The capture buffer display will automatically pan as the capture fill and playback progress through the buffer. During capture fill, if the capture buffer contains more points than can be displayed, points are skipped. This speeds up the display update so that it keeps up with the real time capture but allows visual aliasing to occur. Once capture is complete, the display is redrawn showing the envelope of all points, eliminating any visual aliasing effects.

### User Functions

A User Function displays the results of a user defined math function. User Functions defined within the Histogram Measurement Group may include histogram measurement results. Use the [User Math] menu to define a function.

A User Function may not be selected as the measurement if it uses a Trace which is empty or which contains data taken in a different Measurement Group.

See 'User Math Functions' in Chapter 2 for more.

---

### View

Select the View of the active display Measurement. A View is a way of looking at complex measurement data.

Each display has its own Measurement and each Measurement has an associated View. Changing Measurements changes the View to the View associated with the new Measurement. The two displays generally have different Measurements each with their own associated Views. This entry field can be linked to the Measurement on both displays using the [Link] key.

Command: VIEW (?) d {, i}

The following Views are available:

#### Log Magnitude

Log magnitude view graphs the log of the magnitude of the measurement data. The magnitude is simply  $\sqrt{x^2 + y^2}$  where x is the real part and y is the imaginary part. The Y axis of the display is logarithmic. If <dB Units> in the <Units> submenu are 'Off,' the Y-axis will display a logarithmic grid. If <dB Units> are 'On,' the Y-axis will have a grid which is linear in decibels. The shape of the graph is the same in the two cases.

#### Linear Magnitude

Linear magnitude view graphs the magnitude of the measurement data. The magnitude is simply  $\sqrt{x^2 + y^2}$  where x is the real part and y is the imaginary part. The Y axis of the display is linear in scaling.

#### Magnitude<sup>2</sup>

Magnitude<sup>2</sup> view graphs the magnitude squared of the measurement data. The magnitude squared is simply  $(x^2 + y^2)$  where x is the real part and y is the imaginary part. The Y axis of the display is linear in scaling.

#### Real Part

Real Part view graphs the real part of the measurement data. The Y axis of the display is linear in scaling.

#### Imaginary Part

Imag Part view graphs the imaginary part of the measurement data. The Y axis of the display is linear in scaling.

Real measurement data, such as baseband time record, have zero imaginary part. This view is zero for all points.

Imaginary data arises from the multiplication of the time data by sine and cosine (heterodyne and FFT). Signals in phase with cosine result in real data, signals in phase with sine result in imaginary data. The real and imaginary parts represent data 90 degrees out of phase.

### **Phase**

Phase view graphs the phase of the measurement data. The phase is a four-quadrant quantity defined as  $\theta = \arctan(y/x)$  where  $y$  is the imaginary part and  $x$  is the real part and  $x$  is positive. If  $x$  is negative, an additional 180 deg. is added to the phase. Phase is graphed linearly from -180 ( $-\pi$ ) to +180 ( $\pi$ ) degrees (radians). To show unwrapped phase, choose the Unwrapped Phase view.

Real measurement data, such as baseband time record, have zero imaginary part. The phase is zero for all points.

Single channel phase is relative to the center of the time record for Uniform, BMH, Hanning, Flattop and Kaiser windows. For Force and Exponential windows, phase is relative to the start of the time record. In general, single channel phase is useful only when the time record is triggered in phase with the signal.

For two channel measurements, the phase is channel 2 relative to channel 1. Triggering is not generally required for meaningful two channel phase measurements.

Phase suppress sets the phase of small data values to zero. This avoids the messy phase display associated with the noise floor. (Remember, even a small signal has phase.) Set the phase suppress threshold in the [Display Options] menu.

The FFT can be thought of as a set of bandpass filters, each centered on a frequency bin. The signal within each filter shows up as the amplitude of each bin. If a signal's frequency is between bins, the filters cause phase errors. Because these filters are very steep and selective, they introduce very large phase shifts for signals not exactly on a frequency bin. Use the SR785 source to generate exact bin frequencies whenever possible.

### **Unwrapped Phase**

Unwrapped Phase view graphs the phase of the measurement data as a continuous function without 'wrapping' around at  $\pm 180$  degrees. This view is generally meaningful only for measurements which have data at every frequency point (such as chirp source or swept sine).

The phase  $\theta$  is calculated for each point as in the 'wrapped' Phase view. The value of  $\theta \pm (n \times 360 \text{ deg.})$  which is closest to the phase of the previous point is assigned to each point ( $n$  is an integer). The wrapping starts at the left edge of the display.

Single channel phase is relative to the center of the time record for Uniform, BMH, Hanning, Flattop and Kaiser windows. For Force and Exponential windows, phase is relative to the start of the time record. In general, single channel phase is useful only when the time record is triggered in phase with the signal.

For two channel measurements, the phase is channel 2 relative to channel 1. Triggering is not generally required for meaningful two channel phase measurements.

Phase suppress sets the phase of small data values to zero. This avoids the messy phase display associated with the noise floor. (Remember, even a small signal has phase.) Set the phase suppress threshold in the [Display Options] menu.

### **Nyquist Plot**

The Nyquist plot graphs complex measurement data as imaginary part (along the Y axis) versus real part (along the X axis). The graph is scaled linearly. Adjacent frequency points are connected by a line. This view is generally meaningful only for measurements which have data at every frequency point (such as chirp source or swept sine).

The marker readout shows the real and imaginary parts of each point as well as its frequency. The marker moves sequentially through the frequency points and can be linked to the frequency of the other display.

Real measurement data, such as baseband time record, have zero imaginary part. This view is entirely along the X axis.

### **Nichols Plot**

The Nichols plot graphs complex measurement data as log magnitude (along the Y axis) versus unwrapped phase (along the X axis). The Y axis is logarithmic and the X axis is linear. Adjacent frequency points are connected by a line. This view is generally meaningful only for measurements which have data at every frequency point (such as chirp source or swept sine).

The marker readout shows the log magnitude and phase of each point as well as its frequency. The marker moves sequentially through the frequency points and can be linked to the frequency of the other display.


Real measurement data, such as baseband time record, have zero phase. This view is entirely along the Y axis.

---

## **Units**

Display the Units submenu. Press <Return> for the main [Display Setup] menu.

Choices in the units submenu are limited to those which are appropriate for the current Measurement and View. For instance, dB units may not be selected when the view is Real Part. Settings in the units submenu are associated with the Measurement and View. Changing Measurements or Views may change the settings within the Units submenu

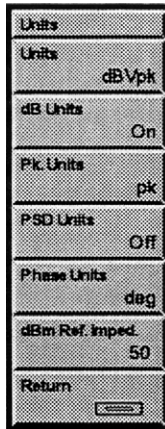


---

## Units

This display-only menu box shows the units of the active display. The units will depend upon the measurement, the view, the Transducer Parameter settings, and the settings of the <dB Units>, <Pk. Units>, <PSD Units>, <Phase Units>, and <dBm Ref. Imped> softkeys.

---

## dB Units

Determines whether a Log Mag view will be shown in dB (decibel) units. If off, no decibel conversion will be performed. If set to 'On,' the measurement results will be converted to decibels relative to the underlying units of the display. If the view is Log Mag and dB Units is off, the graph will be shown with a logarithmic y-axis. If dB units is on, the graph will be shown with a y-axis linear in dbs.

Setting this softkey to 'dBm,' will convert measurement results to decibels relative to 1 mW. The voltage to power conversion is done by assuming that  $\text{power} = \text{volts} \cdot \text{volts} / \text{impedance}$  where impedance is the setting of the <dBm Ref. Imped> softkey.

dBspl will convert measurement results to sound pressure level units. dBspl measures sound pressure in decibels relative to a sound pressure of 20  $\mu$ Pascals. dBspl is only allowed if the Engineering Units for the corresponding input are set to Pascals.

Command: UNDB (?) d {, i}

---

## Pk Units

Determines whether the active display units will be peak, peak to peak, or rms. This softkey is inactive if the underlying measurement and view is unitless, e.g. frequency response.

Command: UNPK (?) d {, i}

### PSD Units

Power Spectral Density (PSD) is the energy in each frequency bin normalized to the noise bandwidth of the bin. The noise bandwidth is calculated from the FFT Linewidth and the Window function. This allows comparison of results taken at different spans and with different Windows. If <PSD Units> are on, the units of the display will be  $V/\sqrt{\text{Hz}}$ , or  $V^2/\text{Hz}$ .

Command: PSDU (?) d {, i}

---

### Phase Units

Selects radians or degrees as units for the phase and unwrapped phase views.

Command: UNPH (?) d {, i}

---

### dBm Reference Impedance

Set the dBm Reference Impedance for both displays.

The dBm Reference Impedance is used to calculate dBm units. dBm is simply  $10 \cdot \log(\text{power}/1\text{mW})$  where  $\text{power} = V_{\text{rms}}^2/\text{Reference Impedance}$ .

Command: DBMR (?) {f}

---

### Y Max

Set the Y Maximum (top reference) of the active display. The top reference is the Y value of the top of the graph. The top reference has the same units as the display.

When Y Max is changed, the value of Y-Min is kept fixed and Y-mid and Y/Div are adjusted.

The two displays can have different top reference values. This entry field can be linked to both displays using the [Link] key.

Command: YMAX (?) d {, x}

---

### Y Mid

Set the Y Midpoint (center reference) of the active display. The center reference is the Y value of the center of the graph. The center reference has the same units as the display.

The two displays can have different center reference values. This entry field can be linked to both displays using the [Link] key.

When Y Mid is changed, the value of Y/Div is kept fixed and Y-Max and Y-Min are adjusted.

---

Command: YMID (?) d {, x}

---

## Y Min

Set the Y Minimum (bottom reference) of the active display. The bottom reference is the Y value of the bottom of the graph. The bottom reference has the same units as the display.

The two displays can have different bottom reference values. This entry field can be linked to both displays using the [Link] key.

When Y Min is changed, the value of Y-Max is kept fixed and Y-mid and Y/Div are adjusted.

Command: YMIN (?) d {, x}

---

## Y / Div

Set the Y/Division scale of the active display. This value is the vertical scale of the graph.

The knob adjusts the vertical scale in a 1-2-5-10 sequence. Knob adjustments set the Y Midpoint (center reference) to the current Marker X Position. The Y Maximum and Y Minimum are adjusted to give the effect of vertical zooming. Entering a scale numerically leaves the Y Min (bottom reference) unchanged. The Y Max and Y Mid are adjusted.

If the View is 'Log Mag' and <dB Units> are set 'Off,' then the graph will be logarithmically scaled. In this case the box will display the total number of decades along the Y-axis.

The two displays can have different vertical scales. This entry field can be linked to both displays using the [Link] key.

Command: YDIV (?) d {, x}

---

## X Center (Polar)

Set the X Center of the active display when the View is Nyquist or Nichols plot. The X Center is the X axis value at the center of the graph. The X Center has the same units as the display.

The two displays can have different X Center values. This entry field can be linked to both displays using the [Link] key.

Command: XCEN (?) d {, x}

### X / Div (Polar)

Set the X/Division scale of the active display when the View is Nyquist or Nichols plot. This value is the horizontal scale of the graph. The knob adjusts the scale in a 1-2-5-10 sequence. Changing the scale does not change the X Center location.

The X axis is scaled linearly and the scale is expressed in the display units.

The two displays can have different scales. This entry field can be linked to both displays using the [Link] key.

Command: XDIV (?) d {, x}

---

### Y Center (Polar)

Set the Y Center of the active display when the View is Nyquist or Nichols plot. The Y Center is the Y axis value at the center of the graph. The Y Center has the same units as the display.

The two displays can have different Y Center values. This entry field can be linked to both displays using the [Link] key.

Command: YCEN (?) d {, x}

---

### Y / Div (Polar)

Set the Y/Division scale of the active display when the view is Nyquist or Nichols plot. This value is the vertical scale of the graph. The knob adjusts the scale in a 1-2-5-10 sequence. Changing the scale does not change the Y Center location.

If the Y axis is scaled linearly, the vertical scale is expressed in the display units. If the Y axis is logarithmic, the vertical scale is always dB/division.

The two displays can have different scales. This entry field can be linked to both displays using the [Link] key.

Command: Y2DV (?) d {, x}

---

### Pan

Set the Pan (translation) of the active display. This value is the left most displayed bin when the X axis is Zoomed (expanded).

The two displays can have different translations when expanded. This entry field can be linked to both displays using the [Link] key.

Command: XPAN (?) d {, i}



## Zoom

Select the Zoom factor (X axis expand) of the active display. 'Expand' is shown in the display whenever it is expanded.

Expanding the X axis is a convenient way of examining closely spaced details in a display without changing the measurement. Expansion is about the Marker Position whenever possible. Use Pan to translate a zoomed display.

Zoom is not allowed when the X axis is logarithmic or when the View is polar (Nyquist or Nichols plot).

The two displays can have different Zoom values when expanded. This entry field can be linked to both displays using the [Link] key.

Command: XZOM (?) d {,i}

## Display Options Menu

The Display Options menu configures the active display format.

Display Options	
Display	Live
Format	Dual
X-Axis	Linear
RPM frequency	Off
Grid	On
Grid Div	10
Nyquist Grid	Rectangular
Phase Suppress	0.0000e+000
d/dx window (%)	0.5

---

### Display

Select the Display Update of the active display [Live, Off-Line].

A Live display updates continuously to show live measurement results (as long as the measurement is not paused).

An Off-Line display does not update. The display shows a ‘frozen’ measurement result or trace. A display becomes Off-Line whenever a trace is recalled to the display from disk or memory. Since an Off-Line display is showing old data, its Measurement, Window, Frequency Span and Averaging may not be modified. Only display related parameters such as View, Scaling and Marker Functions may be changed Off-Line.

The two displays can have different Display Update Modes. This entry field can be linked to both displays using the [Link] key.

A single display may be taken Off-Line while the other display is still Live. This allows comparison of live results with a previous result. This is unlike the [Pause/Cont] key which pauses ALL measurements.

Command: DISP (?) d {, i}

---

### Format

Select the screen Display Format [Single, Dual, or Overlay].

In Single format, the active display is shown. Use the [Active Display] key to switch between DisplayA and DisplayB.

In Dual format, DisplayA is always on top. The [Active Display] key switches the active display from top (A) to bottom (B).

In the Overlay format, the inactive display is shown in half intensity video on the same graph as the active display. The axes markers and menus still reflect the values for the active display.

Command: DFMT (?) {i}

---

## **X-Axis**

Select the X Axis Scale Type for the active display [Linear, Logarithmic].

Logarithmic scaling is only available for FFT measurements with frequency as the X axis. Displaying a measurement with a time axis (Time Record, Capture) will change the scaling to linear (and disable this softkey).

Octave Analysis is always displayed on a logarithmic X axis. The X axis type for Swept Sine measurements is set by the Sweep Type. This softkey is disabled for the Octave Analysis, Swept Sine, Order, and Time/Histogram Measurement Groups.

The logarithmic scale does not affect the frequency bins of the FFT measurement. The resolution is still linear. The points are simply graphed with a logarithmic axis. DC bins are not shown.

The two displays can be scaled differently. This entry field can be linked to both displays using the [Link] key.

Command: XAXS (?) d {, i}

---

## **RPM Frequency**

Select the Frequency Axis Units for the active display [Off, On].

RPM frequency is only available for FFT measurements with frequency as the X axis. If RPM Frequency is On, frequency values will be converted to RPM. If RPM Frequency is Off, the values are displayed in Hz.

Command: RPMF(?) {i}

---

## **Grid**

Select the Grid for the active display [On, Off].

The grid is the set of dotted lines on the display which mark each scale division.

## 4-48 Display Options Menu

---

The two displays have their own Grid. This entry field can be linked to both displays using the [Link] key.

Command: GRID (?) d {, i}

---

### Grid Div

Select the number of Grid Divisions for the active display [8, 10, 12, 15].

Changing the number of Grid Divisions changes the vertical scaling (Y/div) and leaves the display references unchanged. If the view is polar, the horizontal scaling (X/div) is also changed. Change the scaling in the [Display Setup] menu.

The two displays have their own number of Grid Divisions. This entry field can be linked to both displays using the [Link] key.

Command: GDIV (?) d {, i}

---

### Nyquist Grid

Select the type of grid for 2-dimensional (Nyquist View) graphs.

If this softkey is set to rectangular, the Nyquist View is displayed with a rectangular grid. In this case, X-Center, Y-Center, Y/div, and X/div are all individually adjustable.

If the grid is set to polar, a polar grid is displayed. In this case, X center and Y center are fixed at 0. Y/div is adjustable and gives the number of Y units in each radial spoke of the grid.

Command: GDIV (?) d {, i}

---

### Phase Suppress

Set the Phase Suppress Threshold for the active display [0.. ]. When the display View is Phase or Unwrapped Phase, the phase of data points with a squared magnitude less than the threshold are displayed as zero. This avoids the messy phase display associated with the noise floor. (Remember, even a small signal has phase.)

To display the phase of all points, set the threshold to 0.0.

To suppress the phase of signals below -60 dBV (1.0E-3 V), set the threshold to 1.0E-6 V<sup>2</sup> (magnitude squared).

The two displays have their own Phase Suppress Threshold. This entry field can be linked to both displays using the [Link] key.

Command: PHSL (?) d {, x}

### **d/dx Window**

Set the d/dx Window for the active display as a percentage of the display length [0.0005..10 ].

The percentage is converted to a number of display bins (minimum of 1) and the aperture is 2 times this number plus 1. The aperture is used in the calculation of d/dx and group delay in user math functions. Larger windows result in smoother graphs with less x-axis resolution.


Derivative calculations are not available in Octave Analysis measurements and this menu box cannot be changed.

The two displays have their own d/dx Window. The same function may viewed in the two displays with different d/dx windows. This entry field can be linked to both displays using the [Link] key.

Command: DDXW (?) d {, x}

## Marker Menu

The Marker menu configures the display markers. Marker functions such as Harmonic Distortion, Sideband Power or Band Analysis are available in this menu.

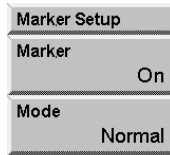

---

### Marker

Select the Marker Tracking of the active display [Off, On, Track, Link].

Off turns the Marker off. It may be desirable to do this before printing the display.

On turns the Marker on. Use the knob (whenever no entry is pending) to move the Marker Region.

Track automatically moves the Marker position to the maximum point in the display. For Normal Marker, minimum tracking is also available by choosing Minimum for the Marker Seeks selection.

Link moves the Marker Region in the other display with the Marker Region in this display. The marker setup of the other display still configures its marker display.

Each display has its own Marker Tracking. This entry field can be linked to both displays using the [Link] key.

Command: MRKR (?) d {, i}

---

### Mode

Select the Marker Mode of the active display [Normal, Harmonic, Sideband, Band]. Choosing the Marker Mode also configures the rest of this menu.

FFT and Order measurement groups can use any Marker Mode.

Octave measurement group can use Normal and Band modes only.

Correlation, Swept Sine, and Time/Histogram measurement groups can only use Normal mode.

Each display has its own Marker Mode. This entry field can be linked to both displays using the [Link] key.

Command: MKMD (?) d {, i}

### Normal

The Normal Marker Mode defines a Marker Region between two solid vertical lines. The Width may be set to a single point (Spot) or 1/2 (Normal) or 1 (wide) grid division. The Marker is a small square which Seeks the Max, Min or Mean of the data within the Marker Region.

When seeking Max or Min, the Marker is located at the max or min data point. This allows peaks and valleys in the display to be easily read by the Marker Position display (above the graph). When the Marker Seeks the Mean, the X position is the center of the Marker Region and the Y position is the Mean of the data within the region.

When a Spot Marker is used, the Marker is confined to a single X axis location.

Use [Marker Max] or [Marker Min] to move the Marker to the maximum or minimum of the graph.

The Marker Position display above the graph shows the X position (frequency or time) and the Y value (amplitude) of the Marker.

The knob moves the Marker Region of the active display whenever there is no pending entry. When the display is a Waterfall, use [Alt] knob to move the marker from record to record (along the Z axis) and scroll the display through the waterfall memory. Use the knob (without [Alt]) to move the marker along the X axis within a record.

### Harmonic

The Harmonic Marker Mode defines a Fundamental Marker with a solid vertical line. In addition to the Fundamental Marker, a number of harmonics are identified by small triangular Harmonic Markers. The Marker Position Display can show the position of the Fundamental or of a single Harmonic.

The Total Harmonic Distortion (THD) is calculated within the [Marker] menu (in dB and percent). Only those harmonics within the measurement span which are identified by Harmonic Markers contribute to the THD. The Harmonic Power (in Vrms) is also shown. To measure THD+N, use Band Marker.

Harmonic Marker is only available for FFT and Order measurements in the frequency domain.

The knob moves the Fundamental Marker of the active display whenever there is no pending entry. Use [Marker Max] to move the Fundamental Marker to the peak of the graph. The harmonic identification works best if the fundamental frequency is a bin frequency. If the fundamental is off-bin, then it is best to use the Flattop window to broaden all of the peaks.

Harmonic Marker is not available for Correlation, Octave, Swept Sine, or Time/Histogram measurement groups.

### Sideband

The Sideband Marker Mode defines a Fundamental (carrier) Marker with a solid vertical line. In addition to the Fundamental Marker, a number of sidebands (discrete

frequencies) are identified by small triangular Sideband Markers. The Marker Position Display can show the position of the Fundamental or of a single Sideband.

The Sideband Power is calculated within the [Marker] menu (in display units and dB). Only those sidebands within the measurement span which are identified by the Sideband Markers contribute to the Sideband Power.

Sideband Marker is only available for FFT measurements in the frequency domain.

The knob moves the Fundamental Marker of the active display whenever there is no pending entry. Use [Marker Max] to move the Fundamental Marker to the peak of the graph. The sideband identification works best if the fundamental frequency is a bin frequency. If the fundamental is off-bin, then it is best to use the Flattop window to broaden all of the peaks.

Sideband Marker is not available for Correlation, Octave or Swept Sine, or Time/Histogram measurement groups.

### **Band**

The Band Marker Mode defines two frequency bands in the display. The Upper (2) band is delimited by dashed vertical markers and the Lower (1) band with solid vertical markers. The pairs of vertical markers are connected by shaded bands at the bottom of the graph. These shaded bands define the bins within the Upper and Lower bands.

Use <Modify Band> to modify either the Upper or Lower band. Use the knob to move selected band. Use [Alt] knob to re-size the selected band.

The two Band Powers are calculated within the [Marker] menu (in Vrms). Only those frequencies within the two band limits contribute to the Band Powers.

Band Exclude can be used to remove a small band from a larger band. For example, exclude '2from1' excludes those bins in the Upper band from the Lower band. The Lower Band Power only includes those bins in the Lower band which are not also in the Upper band.

The Band Ratio is the ratio of the Upper Band Power to either the Lower Band Power or the sum of the two band powers. The ratio is displayed in % and dB.

Band Ratio can be used to measurement THD+N. In this case, use the Lower (1) band to define a small region around the fundamental and the Upper (2) band to define the total bandwidth of interest. Exclude '1from2' so the Upper (2) band does not measure the fundamental. Choose the  $2/(1+2)$  ratio mode to measure THD+N (power with fundamental notched out/total power including fundamental).

Band Marker is only available for frequency domain FFT, Order, and Octave measurements.

Band Marker is not available for the Swept Sine, Correlation, or Time/Histogram measurement groups.



### **Frequency/Damping**

The Frequency/Damping cursor fits the resonant frequency and damping ratio for a single degree of freedom in a frequency response curve. The fit area is moved using the knob. Use [Alt] knob to re-size the fit region. The Frequency/Damping cursor should only be used with frequency response data..

The Frequency/Damping Cursor is only available in the FFT and Swept Sine measurement groups.

## Normal Marker Menu

This menu is displayed when the active display Marker Mode is Normal.

Marker Setup	
Marker	On
Mode	Normal
Seeks	Max
Width	Spot
Rel	Off
X-Rel	0
Y-Rel	0
X-Rel	Off
Marker X to	<input type="text"/>
Target	0

---

### Seeks

Select the Normal Marker Seeks mode of the active display [Max, Min, Mean].

The Marker searches the data points within the Marker Region for the Maximum or Minimum data value or calculates the Mean of the region. The results are shown in the Marker Position display above the graph. When Seeking Max or Min, the Marker is located at the max or min data point. When the Seeking the Mean, the X position is the center of the Marker Region and the Y position is the Mean of the data within the region.

Each display has its own Marker Seeks Mode. This entry field can be linked to both displays using the [Link] key.

Command: MSEK (?) d {, i}

---

### Width

Select the Normal Marker Width of the active display [Spot, Normal, Wide].

The Normal Marker Mode defines a Marker Region between two solid vertical lines. The Width may be set to a single point (Spot) or 1/2 (Normal) or 1 (Wide) grid division. The Marker is the small square which Seeks the Max, Min or Mean of the data within the Marker Region.

Each display has its own Normal Marker Width. This entry field can be linked to both displays using the [Link] key.

Command: MWID (?) d {, i}

---

## Rel

Select the Normal Marker Rel Mode of the active display [Off, Relative, Reference Display, Other Display].

The Marker Position display above the graph shows the X position (frequency or time) and the Y value (amplitude) of the Marker. A relative reading is preceded by a  $\Delta$  (delta) symbol.

The X position may be shown absolute or relative in all cases as selected by <X-Rel>.

When Rel is Off the absolute marker position is displayed.

Relative calculates the marker position relative to the Marker Offset (X-Rel, Y-Rel). The Marker Offset location (if it is within the display) is marked by a small flag shaped symbol.

Reference Display calculates the marker position relative to the Reference Display. The X position is absolute and the Y position is offset by the Reference Display value at the X position. Ref. Disp. is only available if the Reference Display is On.

Other Display calculates the marker position relative to the other display's marker position. The marker offset is just the marker position of the other display. This is generally valid only if the two displays have a similar measurement.

Pressing [Marker Ref] changes from Off to Rel. and sets the Marker Offset to the current Marker Position. [Marker Ref] again changes from Rel. back to Off. [Marker Ref] does not change the other modes.

Each display has its own Normal Marker Rel Mode. This entry field can be linked to both displays using the [Link] key.

Command: MREL (?) d {, i}

---

## X-Rel

Set the X Offset for the Normal Marker of the active display.

When the Marker Rel Mode is Relative, the marker position is calculated relative to the Marker Offset (X-Rel, Y-Rel). The Marker Offset location (if it is within the display) is marked by a small flag shaped symbol.

The Marker Offset has no associated units. The relative calculation simply subtracts the Marker Offset (in the display units) from the absolute marker position. Changing the display units invalidates the Marker Offset. Re-enter the Marker Offset or use [Marker Ref] to reset the Marker Offset in the new units.

## 4-56 Normal Marker Menu

---

Pressing [Marker Ref] changes from Off to Rel. and sets the Marker Offset to the current Marker Position. [Marker Ref] again changes from Rel. back to Off.

Each display has its own Normal Marker X Offset. This entry field can be linked to both displays using the [Link] key.

Command: MROX (?) d {, x}

---

### Y-Rel

Set the Y Offset for the Normal Marker of the active display. Enter a new offset numerically or adjust it with the knob.

When the Marker Rel Mode is Relative, the marker position is calculated relative to the Marker Offset (X-Rel, Y-Rel). The Marker Offset location (if it is within the display) is marked by a small flag shaped symbol.

The Marker Offset has no associated units. The relative calculation simply subtracts the Marker Offset (in the display units) from the absolute marker position. Changing the display units invalidates the Marker Offset. Re-enter the Marker Offset or use [Marker Ref] to reset the Marker Offset in the new units.

Pressing [Marker Ref] changes from Off to Rel. and sets the Marker Offset to the current Marker Position. [Marker Ref] again changes from Rel. back to Off.

Each display has its own Normal Marker Y Offset. This entry field can be linked to both displays using the [Link] key.

Command: MROY (?) d {, x}

---

### X-Rel

Set the X Rel Mode for the Normal Marker of the active display [Off, On].

When the Marker Rel Mode is Relative to the Marker Offset or Relative to the Other Display, the marker Y value is shown relative to the offset or other display.

Choose On to show the marker X value relative to the offset or other display as well.

Choose Off to show the absolute marker X value. This is useful when making relative amplitude measurements while still needing to read the marker X position in absolute units.

Each display has its own X Rel Mode. This entry field can be linked to both displays using the [Link] key.

Command: MXRL (?) d {, i}

**Marker X to**

Sends the marker to the X position specified by the Target softkey. The target position is entered in the current units of the x-axis. If the exact x value specified is not available, the marker will be sent to the closest available X value.

---

**Target**

Specifies the X position the marker will be sent to with the Marker X to button. The target is entered in the current units of the x-axis. If the exact x value specified is not available, the marker will be sent to the closest available X value.

## Harmonic Marker Menu

This menu is displayed when the active display Marker Mode is Harmonic.

Marker Setup	
Marker	On
Mode	Harmonic
# Harmonics	3
Display	Fundamental
Readout Mode	Absolute
THD	0.00737 %
	-82.6 dB
RMS Harm. Power	18.44 $\mu$ Vrms

### Harmonics

Set the Number Of Harmonics for the active display [1.. 100]. The second harmonic (2xfundamental) is identified as #1, etc.

The Harmonic Marker identifies the fundamental frequency with the Fundamental Marker (solid vertical line) and the number of harmonics with small triangular Harmonic Markers. The Marker Position Display can read the Fundamental or a single Harmonic.

Only those harmonics within the measurement span which are identified by the Harmonic Markers contribute to the THD calculations. The harmonic identification works best if the fundamental frequency is a bin frequency. If the fundamental is off-bin, then it is best to use the Flattop window to broaden all of the peaks.

Each display has its own Number Of Harmonics. This entry field can be linked to both displays using the [Link] key.

Command: HRMN (?) d {, i}

### Display

Select the Harmonic and Sideband Display for the active display. The Marker Position can display the Fundamental or any identified Harmonic or Sideband.

For Harmonic Marker, the minimum is 0 (Fundamental) and the maximum is the Number Of Harmonics. The second harmonic (2xfundamental) is identified as #1, etc.

For Sideband Marker, the value ranges from 0 (Fundamental) to plus and minus the Number Of Sidebands. The first lower sideband is identified as #-1, the first upper sideband is #1, etc.

The small square marker will move to the selected harmonic or sideband. If the selection is not within the measurement span, the marker will move to the highest or lowest measured frequency (which will be reported as the Marker Position).

Each display has its own Harmonic and Sideband Display. This entry field can be linked to both displays using the [Link] key.

Command: HDSP (?) d {, i}

---

**Readout Mode**

Select the Readout Mode for the Harmonic and Sideband Marker of the active display [Absolute, Relative]. The Marker Position Display shows the amplitude of the Fundamental or any identified Harmonic or Sideband in either Absolute units or Relative to the fundamental.

Relative Mode reports the amplitude of the selected harmonic or sideband relative to the fundamental. The marker amplitude units are dBFundamental. The frequency is always absolute.

Each display has its own Readout Mode. This entry field can be linked to both displays using the [Link] key.

Command: HRDO (?) d {, i}

---

**THD**

This menu box displays the Total Harmonic Distortion (THD) of the active display. The THD squared is the harmonic power (sum of the squared magnitudes of the harmonics) divided by the sum of the fundamental power plus the harmonic power.

Only those harmonics within the measurement span which are identified by Harmonic Markers contribute to the THD.

The THD is displayed in dB ( $20\log(\text{THD})$ ) and percent.

To measure THD+N, use the Band Marker.

Command: HTHD ? d , i

---

**Harmonic Power**

This menu box displays the Total Harmonic Power of the active display. The harmonic power is the sum of the squared magnitudes of the harmonics identified with a Harmonic

## 4-60 Harmonic Marker Menu

---

Marker. The result is shown in Vrms, or dBVrms depending on the setting of dB Units. To convert a value in Vrms to power, square the result.

Command: HPWR ? d



# Sideband Marker Menu

This menu is displayed when the active display Marker Mode is Sideband.

Marker Setup	
Marker	On
Mode	Sideband
Sideband Sep.	4.096 kHz
# Sidebands	1
Display	Fundamental
Readout Mode	Absolute
Sideband Ratio	-98.3 dB
RMS SB Power	3.038 $\mu$ Vrms

## Sideband Separation

Set the Sideband Separation for the active display.

The Sideband Marker identifies the carrier frequency with the Fundamental Marker (solid vertical line) and the number of sidebands with small triangular Sideband Markers. The Sideband Separation determines the location of these Sideband Markers relative to the Fundamental.

Only those sidebands within the measurement span which are identified by the Sideband Markers contribute to the Sideband Power calculations.

Each display has its own Sideband Separation. This entry field can be linked to both displays using the [Link] key.

Command: SSEP (?) d {, f}

## # Sidebands

Select the Number Of Sidebands for the active display [1..50]. The first upper sideband (fundamental+separation) is identified as #1, the first lower sideband (fundamental-separation) is #-1, etc.

The Sideband Marker identifies the carrier frequency with the Fundamental Marker (solid vertical line) and the number of sidebands with small triangular Sideband Markers. The Sideband Separation determines the location of these Sideband Markers relative to the fundamental.

## 4-62 Sideband Marker Menu

---

Only those sidebands within the measurement span which are identified by the Sideband Markers contribute to the Sideband Power calculations. The sideband identification works best if the fundamental frequency is a bin frequency. If the fundamental is off-bin, then it is best to use the Flattop window to broaden all of the peaks.

Each display has its own Number Of Sidebands. This entry field can be linked to both displays using the [Link] key.

Command: SDBN (?) d {, i}

---

### Display

Select the Harmonic and Sideband Display for the active display. The marker display can read the position of the fundamental or any identified harmonic or sideband.

For sideband marker, the value ranges from 0 (fundamental) to plus and minus the Number of Sidebands. The first lower sideband is identified as #-1, the first upper sideband is #1, etc.

The small square marker will move to the selected harmonic or sideband. If the selection is not within the measurement span, the marker will move to the highest or lowest measured frequency (which will be reported as the marker position).

Each display has its own harmonic/sideband display. This entry field can be linked to both displays using the [Link] key.

Command: HDSP (?) d {, i}

---

### Readout Mode

Select the Readout Mode for the Harmonic and Sideband Marker of the active display [Absolute, Relative]. The Marker Position Display shows the amplitude of the Fundamental or any identified Harmonic or Sideband in either Absolute units or Relative to the fundamental.

Relative Mode reports the amplitude of the selected harmonic or sideband relative to the fundamental. The marker amplitude units are dBFundamental. The frequency is always absolute.

Each display has its own Readout Mode. This entry field can be linked to both displays using the [Link] key.

Command: HRDO (?) d {, i}

---

### Sideband Ratio

This menu box displays the Sideband Ratio of the active display. The sideband power is the sum of the squared magnitudes of the sidebands. Only those sidebands within the

measurement span which are identified by Sideband Markers contribute to the calculation.

The Sideband Ratio in dB is  $10\log(\text{sideband power}/(\text{fundamental power}+\text{sideband power}))$ .

Command: SPWR ? d , i

---

**RMS Sideband Power**

This menu box displays the Sideband Power of the active display. The Sideband Power is the sum of the squared magnitudes of the sidebands. Only those sidebands within the measurement span which are identified by Sideband Markers contribute to the calculation. The result is shown in  $V_{rms}$ , or  $dBV_{rms}$  depending on the setting of dB Units. To convert a value in  $V_{rms}$  to power, square the result

Command: SPWR ? d , i

## Band Marker Menu

This menu is displayed when the active display Marker Mode is Band.

Marker Setup	
Marker	On
Mode	Band
Modify Band	Upper <sup>2</sup>
Band Exclude	none
Band Ratio Mode	<sup>2</sup> / <sub>1</sub>
Band Power	2.958 $\mu$ Vrms
	38.56 $\mu$ Vrms
Band Ratio	7.67 %
	-22.3 dB

### Modify Band

Select the Upper (2) or Lower (1) Band Marker to modify. Use the knob to move the selected band. Use [Alt] knob to re-size the selected band.

The Marker Position display above the graph shows the left and right limits of the selected band and the Band Power. The Band Power is the sum of the squared magnitudes of all frequencies within the band (shown in Vrms). The label '1' or '2' preceding the Marker Position display denotes Lower or Upper Band Marker.

Command: BMKL (?) d, i {, j} or BMKR (?) d, i {, j}

### Band Exclude

Select Band Exclusion for the Band Marker in the active display [none, 1from2, 2from1].

None means that the two bands are independent and include all bins within their left and right limits.

1from2 means the bins in band 1 (Lower) are excluded from the band 2 (Upper) power. The region where the two bands overlap is not included in band 2 (Upper).

2from1 means the bins in band 2 (Upper) are excluded from the band 1 (Lower) power. The region where the two bands overlap is not included in band 1 (Lower).

For example, to measure THD+N, use the Lower (1) band to define a small region around the fundamental and the Upper (2) band to define the total bandwidth of interest.

Exclude '1from2' so the Upper (2) band does not measure the fundamental. Choose the 2/(1+2) ratio mode to measure THD+N (power with fundamental notched out/total power including fundamental).

Command: BEXC (?) d {, i}

---

### Band Ratio Mode

Select the Band Ratio Mode for the Band Marker in the active display [2/1, 2/(1+2)]. The calculated ratio is displayed in the <Band Ratio> menu box.

2/1 specifies that the Band Ratio is computed from the ratio of the Upper (2) band power to the Lower (1) band power.

2/(1+2) specifies that the Band Ratio is computed from the ratio of the Upper (2) band power to the sum of the Lower (1) and Upper (2) band powers.

For example, to measure THD+N, use the Lower (1) band to define a small region around the fundamental and the Upper (2) band to define the total bandwidth of interest. Exclude '1from2' so the Upper (2) band does not measure the fundamental. Choose the 2/(1+2) ratio mode to measure THD+N (power with fundamental notched out/total power including fundamental).

Command: BRAM (?) d {, i}

---

### Band Power

This menu box displays the Upper (top) and Lower (bottom) Band Powers of the active display. The Band Power is the sum of the squared magnitudes of all bins within each Band Marker. The result is shown in Vrms, or dBVrms depending on the setting of dB Units. To convert a value in Vrms to power, square the result

Command: BPWR (?) d, i

---

### Band Ratio

This menu box displays the band ratio, either 2/1 or 2/(1+2) as selected by the Band Ratio Mode, for the active display. The result is shown in % and dB.

Command: BRAT (?) d, i

---

## Frequency/Damping Marker Menu

This menu is displayed when the active display Marker Mode is Frequency/Damping.



---

### Calculate

Calculate the resonant frequency and damping parameters for the data in the cursor region. This calculation will only yield meaningful results for frequency response data.

Command: MMCA? d

---

### Frequency

Displays the resonant frequency of the data in the cursor area. The resonant frequency is calculated by fitting a single-pole frequency response function to the data. The resonant frequency is the imaginary part of this pole.

Command: MMCA? d

---

### Damping Ratio

Displays the damping ratio of the data in the cursor area. The damping ratio frequency is calculated by fitting a single-pole frequency response function to the data. The damping ratio is negative the ratio of the real part of the pole to the magnitude of the pole.

Command: MMCA? d


## Source Menu

The Source menu selects and configures the source waveform. The Measurement Group determines the available source types. Choosing a source type selects the source waveform and displays the appropriate setup menu.



### Watch Out For Triggered Sources

If a triggered source is selected (Chirp, Burst Chirp, Burst Noise or Arbitrary), the triggered measurement phase is stable only if the input signals are derived from the source output. Turn the source off (or set it to Sine) when making triggered measurements without the source.

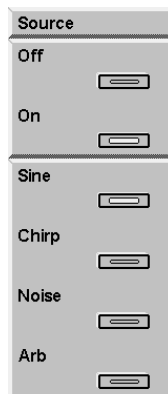


---

### Source Off

Turn off the source output. The output is held at 0 V.

If a triggered source is selected (Chirp, Burst Chirp, Burst Noise or Arbitrary), the phase of the measurement is relative to the source, not the trigger.

Command: SRCO (?) {i}

---

### Source On

Turn on the source output. The output is determined by the selected source type.

Command: SRCO (?) {i}

---

### Sine

Set the source output to Sine and display the Sine Setup menu.

The output is the sum of two tones (sine waves) plus the offset. To generate a single tone, set the amplitude of one of the tones to zero. The frequencies should be exact bin frequencies of the FFT. This eliminates windowing effects in the measured amplitude and phase.

The sine source is not triggered. The output is always continuous.

Command: STYP (?) {i}

---

### Chirp

Set the source output to Chirp and display the Chirp Setup menu.

The output is an equal amplitude sine wave at each frequency bin of the FFT spectrum. This source is useful for measuring transfer functions quickly. Use Source Trigger for stable time records and phase. The Uniform window should be used since individual frequencies are not present at all times in the time record. The chirp is exactly periodic with the FFT time record so the Uniform window does not cause 'leakage' problems.

This source is only available for the FFT and Correlation Measurement Groups. The chirp frequency span and burst time record is linked to either DisplayA or DisplayB. If the other display has a different span, the chirp span will not be appropriate for that display.

The chirp source is triggered by External triggers. Set the Trigd Source Mode to trigger once or with every trigger. With the continuous Trigger, the chirp repeats over and over.



#### Trigger Caution

If the Chirp source is selected, the triggered FFT measurement phase is stable only if the input signals are derived from the triggered source output. Turn the source off (or set it to Sine) when making triggered measurements of external signals (not the source).

Command: STYP (?) {i}

---

### Noise

Set the source output to Noise and display the Noise Setup menu.

The output is broadband noise. This source is useful for estimating transfer functions quickly. White noise provides equal amplitude per  $\sqrt{\text{Hz}}$ . White noise is useful in electronic applications. Pink noise rolls off at 3 dB/oct providing equal energy per octave. Pink noise is preferred for octave analysis.

Burst noise is selected whenever the Burst percentage is less than 100%. The burst noise source is triggered by External triggers. Set the Trigd Source Mode to trigger once or with every trigger. With Free Run Trigger Mode, the burst repeats over and over.

In FFT group, the noise bandwidth and burst time record is linked to either DisplayA or DisplayB. If the other display has a different span, the source may not be appropriate for that display.

Since the signal is random, windows are always required when making FFT measurements using the noise source.



In Octave group, the burst period is not linked to a display but is specified in time. The bandwidth is always full bandwidth.

**Trigger Caution**

If the Noise source is selected and the Burst is less than 100%, the triggered FFT measurement phase is stable only if the input signals are derived from the triggered source output. Turn the source off (or set it to Sine) when making triggered measurements of external signals (not the source).

Command: STYP (?) {i}

---

**Arbitrary**

Set the source output to Arbitrary and display the Arbitrary Setup menu.

The Arbitrary source plays a waveform stored in memory. The waveform can be either Capture buffer or the Arbitrary Waveform memory. Use the <Memory Allocation> menu to allocate memory between the capture buffer, waterfall storage and the arbitrary source waveform. The Capture buffers may not be selected as the Arbitrary source if they do not contain data.

The Arbitrary waveform memory can be loaded via the computer interfaces, from a disk file or from a stored trace. The Capture buffer is filled by capturing an input signal.

The Arbitrary source can play a portion of memory starting at a specified point. The output sampling rate can also be specified. The amplitude of the source is set relative to 1 V (up to a maximum of 5 V).

**Triggering**

The Arbitrary source is triggered by External triggers. Set the Trigd Source Mode to trigger once (1-Shot) or with every trigger (Continuous).

In the 1-Shot case, the source only triggers on the first trigger after [Start/Reset]. This same trigger starts the measurement (and waterfall). This is useful when the source length is very long compared to the measurement time. This allows a sequence of many measurements to be triggered at the start of a single playback of a long source.

In the Continuous case, the source triggers on every trigger, with a minimum trigger period set by the Arbitrary Source Length. This is useful when the source length is the same as the FFT time record length.

Source trigger outputs the source continuously over and over. The FFT time records are synchronized to the source start.

Continuous Trigger outputs the source continuously over and over.

Do not use Ch1 or Ch2 input trigger since the output will not start until a trigger is received.




### **Trigger Caution**

If the Arbitrary source is selected, the triggered FFT measurement phase is stable only if the input signals are derived from the triggered source output. Turn the source off (or set it to Sine) when making triggered measurements of external signals (not the source).

Command: STYP (?) {i}

# Sine Source Menu

The Sine Source menu is used to set the frequency, offset, and amplitude of the sine source.

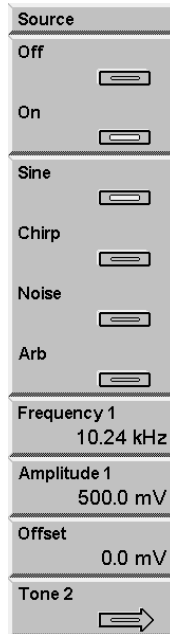

## Frequency 1

Set the Frequency of Tone 1. The sine output is the sum of two tones (sine waves). When the knob is used to adjust the frequency, the resolution is equal to the Linewidth of the active display FFT Span. The knob always sets the frequency to an exact multiple of the Linewidth. For octave analysis, the Linewidth of the most recent FFT Spans are used (with an FFT Base of 100.0 kHz).

The keypad allows arbitrary frequencies to be entered. Remember, the output is periodic over the FFT time record only if the frequency is an exact multiple of the Linewidth. Source Trigger will not result in a stable phase for non-periodic frequencies.

The sine source is not triggered. The output is always continuous.

Command: S1FR (?) {f}

## Amplitude 1

Set the peak Amplitude of Tone 1 [0..5V]. The amplitude resolution is 0.1 mV. The sine output is the sum of two tones (sine waves) and the offset (constant). To generate a single tone, set the amplitude of one of the tones to zero. Note that the sum of the amplitudes of Tone 1, Tone 2 and the absolute value of the offset cannot exceed 5 V.

Command: S1AM (?) {x}

## 4-72 Chirp Source Menu

---

### Offset

Set the DC Offset of the Sine Source [-5V..5V]. The offset resolution is 0.1 mV. The sine output is the sum of two tones (sine waves) and the offset (constant). Note that the sum of the amplitudes of Tone 1, Tone 2 and the absolute value of the offset cannot exceed 5 V.


Using large offsets with small tone amplitudes will degrade the distortion performance of the sine source.

Command: SOFF (?) {x}

---

### Tone 2

Display the Tone 2 Settings menu. Press <Return> for the main [Source] menu.

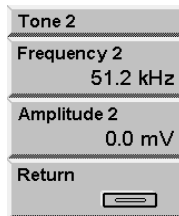

### Frequency 2

Set the Frequency of Tone 2. The sine output is the sum of two tones (sine waves). When the knob is used to adjust the frequency, the resolution is equal to the Linewidth of the active display FFT Span. The knob always sets the frequency to an exact multiple of the Linewidth. For octave analysis, the Linewidth of the most recent FFT Spans are used (with an FFT Base of 100.0 kHz).

The keypad allows arbitrary frequencies to be entered. Remember, the output is periodic over the FFT time record only if the frequency is an exact multiple of the Linewidth. Source Trigger will not result in a stable phase for non-periodic frequencies.

The sine source is not triggered. The output is always continuous.

Command: S2FR (?) {f}

---

### Amplitude 2

Set the Amplitude of Tone 2 [0..5V]. The amplitude resolution is 0.1 mV. The sine output is the sum of two tones (sine waves) and the DC offset (constant). To generate a single tone, set the amplitude of one of the tones to zero. Note that the sum of the amplitudes of Tone 1, Tone 2 and the absolute value of the offset cannot exceed 5 V.

Command: S2AM (?) {x}

---

## Chirp Source Menu

The Chirp Source menu is used to configure the chirp source.

<b>Source</b>	
Off	<input type="checkbox"/>
On	<input type="checkbox"/>
<b>Sine</b>	
Sine	<input type="checkbox"/>
Chirp	<input type="checkbox"/>
Noise	<input type="checkbox"/>
Arb	<input type="checkbox"/>
<b>Amplitude</b>	1000.0 mV
<b>Burst</b>	100 %
<b>Source&lt;-&gt;Display</b>	Display A

### Amplitude

Set the peak Amplitude of the chirp. The peak output level is approximate due to the ripple in the source output reconstruction filter.

Changing the FFT Resolution changes the amplitudes of the individual frequency components relative to the peak. If the FFT Resolution is 400 lines, the amplitude of each frequency component is about -30 dB relative to the peak amplitude. If the individual frequency components were perfectly random, each component would be  $1/\sqrt{400}$  (-26 dB) relative to the peak. However, the chirp waveform is identical from time record to time record and each component has a fixed phase relative to all the other components. This worsens the crest factor by a few dB and reduces the amplitude of each component for a specified peak amplitude. The input dynamic range of the measurement is reduced when using the chirp source.

Command: CAMP (?) {x}

### Burst

Set the Burst Percentage of the chirp [0.1%..100%]. The chirp waveform is output over a percentage of the FFT time record of the display selected as the Source Display. Burst is only available for the FFT Measurement Group.

For a continuous output, use 100% burst with continuous Trigger Source or Source Trigger. Source Trigger will synchronize the FFT time record with the chirp waveform.

## 4-74 Chirp Source Menu

---

With External Trigger, the chirp waveform is triggered along with the FFT time record. Set the Trigd Source Mode to Continuous to trigger the source with every time record.

Do not use Ch1 or Ch2 input trigger since the output will not start until a trigger is received.

Command: CBUR (?) {x}

---

### Source Display

Select the Source Display [DisplayA, DisplayB]. The Source Display only applies for the FFT Measurement Group.

The Source Display is the display which determines the burst period (FFT time record) and the bandwidth for bandlimited noise and chirp (FFT Span). There is a single Source Display for both Chirp and Noise outputs. Changing this entry for one type of source also changes it for the other.

Select which display sets the time record and frequency span of the chirp. If the other display has a different span, the chirp will not be appropriate for that span.

Burst Chirp is a percentage of the Source Display FFT time record length.

Command: CSRC (?) {i}

# Noise Source Menu

The Noise Source menu is used to configure the noise source.

<b>Source</b>	
Off	<input type="checkbox"/>
On	<input type="checkbox"/>
<b>Sine</b>	
Sine	<input type="checkbox"/>
Chirp	<input type="checkbox"/>
Noise	<input type="checkbox"/>
Arb	<input type="checkbox"/>
<b>Amplitude</b> 500.0 mV	
<b>Type</b> BL White	<b>Type</b> White
<b>Burst</b> 100 %	<b>Burst</b> 100.000 %
<b>Source&lt;-&gt;Display</b> Display A	<b>Source Period</b> 100 ms

## Amplitude

Set the peak Amplitude of the noise output. The source output will overshoot this amplitude by as much as 100% a small percentage of the time. Because of the nature of noise, the peak amplitude is not perfectly defined.

Command: NAMP (?) {x}

## Type

Select the Type of noise output. [White, Pink, Bandlimited White or Pink].

### White Noise

White noise provides equal noise density ( $V/\sqrt{\text{Hz}}$ ) from 0 to above 102.4 kHz (regardless of the measurement bandwidth of the displays). The spectrum of white noise appears flat in an FFT spectrum. The power spectral density (PSD) remains constant at all spans.

### Pink Noise

Pink noise rolls off at 3 dB/oct providing equal energy per octave and extends beyond 102 kHz. The spectrum of pink noise appears flat in octave analysis.

### **BandLimited White**

Bandlimiting restricts the noise bandwidth to the measurement span of the Source Display. Bandlimited noise is available only for the FFT and CorrelationMeasurement Group.

Command: NTYP (?) {i}

---

### **Burst**

Set the Burst Percentage of the noise [0.1%..100%].

In FFT group, the noise waveform is output for a percentage of the FFT time record of the display selected as the Source Display.

In Octave group, the burst period is not linked to a display but is specified in time by the Source Period.

For a continuous output, use 100% burst.

#### **Triggering**

When the burst percentage is less than 100%, the burst noise source is triggered by External triggers. Set the Trigd Source Mode to trigger once (1-Shot) or with every trigger (Continuous). In the 1-Shot case, the burst source only triggers on the first trigger after [Start/Reset]. This same trigger starts the measurement (and waterfall). In the Continuous case, the source triggers on every trigger, with a minimum trigger period set by the time record of the Source Display (FFT, Correlation, Time/Histogram) or the Source Period (Octave, Order).

Source trigger outputs a noise burst every source period or FFT time record. The FFT time records are synchronized to the burst.

Free Run Trigger Mode outputs a noise burst every source period.

Do not use Ch1 or Ch2 input trigger since the output will not start until a trigger is received.

Command: NBUR (?) {x}

---

### **Source Display**

Select the Source Display [DisplayA, DisplayB]. The Source Display only applies for the FFT and Correlation Measurement Groups. The Source Display is the display which determines the burst time record (FFT time record) and the bandwidth for bandlimited noise and chirp (FFT Span). There is a single Source Display for both Chirp and Noise outputs. Changing this entry for one type of source also changes it for the other.

Select which display sets the bandwidth for bandlimited noise. Burst Noise is a percentage of the Source Display FFT time record length.

---



The noise output is continuous for Octave Analysis.

Command: CSRC (?) {i}

---

## Source Period

Set the Burst Period for noise in Octave group [4 ms..1ks].

In Octave and Order group, the burst period is not linked to a display but is specified in time by the Source Period. The burst output is present for a Burst percentage of the Source Period.

### Triggering

When the burst percentage is less than 100%, the burst noise source is triggered by External triggers. Set the Trigd Source Mode to trigger once (1-Shot) or with every trigger (Continuous). In the 1-Shot case, the burst source only triggers on the first trigger after [Start/Reset]. This same trigger starts the measurement (and waterfall). In the Continuous case, the source triggers on every trigger, with a minimum trigger period set by the Source Period.

Source trigger outputs a noise burst every source period.


Free Run Trigger Mode outputs a noise burst every source period.

Do not use Ch1 or Ch2 input trigger since the output will not start until a trigger is received.

Command: NPER (?) {x}

## Arbitrary Source Menu

The Arbitrary Source menu is used to configure the arbitrary source.

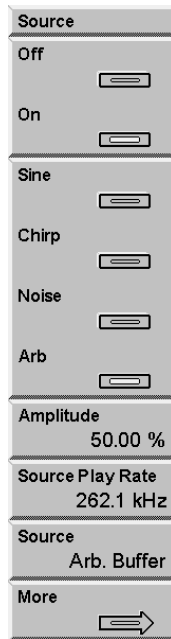


---

### Amplitude

Set the Amplitude of the arbitrary source relative to 1V [0%..500%]. The maximum output is 5V.

If the Arbitrary Waveform has been copied from a Trace, the data is normalized so that the point with the largest absolute value is output at  $\pm 1V$  when the Amplitude is 100%.

When the Arbitrary Waveform is loaded from the computer interface, the data ranges from -1.0 to +1.0. A value of 1.0 is output at 1V when the Amplitude is 100%.

When playback is from a Capture buffer, the input range corresponds to full scale. For example, if a 0.1 Vpk (-20 dBV) signal is captured on an Input Range of -20 dBV (full scale at the input), it will play back at 1V when the Amplitude is 100%. To duplicate the actual signal amplitude, set the Amplitude to 10%.

Command: AAMP (?) {i}

---

### Source Play Rate

Select the Arbitrary Source Play Rate. The Play Rate can be 1, 1/2, 1/4, 1/8, ... times the maximum sampling rate. The maximum sampling rate is 262.1 kHz when the FFT Base Frequency is 102.4 kHz and 256 kHz when the FFT Base Frequency is 100.0 kHz OR the Measurement Group is Octave.

When the Arbitrary source is a Capture buffer, the Play Rate defaults to the capture Sampling Rate. The Play Rate should be equal to the capture Sampling Rate to preserve the frequency content of the buffer.

When the Arbitrary source is the Arbitrary waveform memory, the Play Rate should be set according to the rate used to calculate the waveform samples.

When the Arbitrary source is copied from a time record trace, the Play Rate should be set to the sampling rate at which the time record was acquired.

**Caution**

In both cases, if the Play Rate is LESS than the maximum sampling rate, the output will contain alias frequency components above 1/2 the sampling rate. In general, the frequency span of any measurement using this source should not exceed 1/2 times the Play Rate to avoid measuring these alias frequencies.

Command: ARAT (?) {i}

**Source**


Select Arbitrary waveform playback from the Arbitrary Waveform memory or from a Capture buffer. Use the <Memory Allocation> menu to allocate memory between the capture buffer, waterfall storage and the arbitrary source waveform.

The Arbitrary waveform memory can be loaded from a stored trace, via the computer interfaces or from a disk file. The Capture buffer is filled by capturing an input signal.

Command: ASRC (?) {i}

**More**

Display the Arbitrary Source Settings menu. Press <Return> for the main [Source] menu.

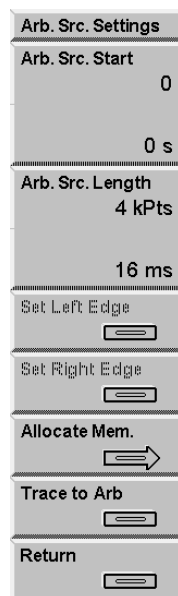

---

### Arb. Source Start

Select the Arbitrary Source Start position within the selected source buffer (either Arbitrary Waveform memory or Capture buffer). The start position is specified as a bin number or individual point position. During Arbitrary playback, the points are separated by a time equal to 1/Play Rate. The corresponding start time relative to the start of source buffer is displayed as well.

Command: ASTR (?) {i}

---

### Arb. Source Length

Select the Arbitrary Source Length in 2 kPoint (2048 points) increments. The Arbitrary Source Length (in points) plus the Arbitrary Source Start cannot exceed the length of the Arbitrary source buffer, either the Arbitrary waveform memory or the Capture Length. The corresponding source time length is displayed as well.

Command: ALEN (?) {i}

---

### Set Left Edge

The Arbitrary Source playback must be from a Capture Buffer and the Measurement of the active display must be the Capture Buffer for this key to be active.

Choose the correct Capture Buffer as the Measurement and use Zoom and Pan to show the region of interest. Press <Set Left Edge> to set the Arbitrary Source Start to the marker position. The Arbitrary Source Length is not changed.

If the requested Source Start plus the Source Length (in points) exceeds the Capture Length of the buffer, the Source Length is reduced.

---

### Set Right Edge

The Arbitrary Source playback must be from a Capture Buffer and the Measurement of the active display must be the Capture Buffer for this key to be active.

Choose the correct Capture Buffer as the Measurement and use Zoom and Pan to show the region of interest. Press <Set Right Edge> to set the Arbitrary Source Length to the marker position minus the Source Start. The Source Start is not changed.

The Source Length is always set to a 2 kPoint (2048 points) increment.

If the marker position is to the left of the Source Start, then the Source Length is set to the minimum (2 kPts).

---

**Allocate Memory**

Display the Memory Allocation menu.

**Total Available**

Displays the total memory available for storage, either 980, 2004 or 4052 blocks, depending upon the amount of installed memory. Each block is 2 kPoints (2048 points).

The total of the Capture, Waterfall and Arbitrary memory allocations cannot exceed the Total Available memory. It may be necessary to decrease one allocation in order to increase another.

Command: MMEM ?

**Capture Memory**

Allocates memory blocks for the capture buffer. The allocated Capture Memory sets the limit for the Capture Length. Capture memory must be allocated before the capture buffer may be used.

**Waterfall Memory**

Allocates memory blocks for waterfall storage. Waterfall memory must be allocated before waterfall displays may be used.

**Arb. Memory**

Allocates memory blocks for the Arbitrary Source waveform. The Arbitrary Source can also use the capture buffer as the output waveform.

**Clear Allocation**

Clears the existing memory allocations in this menu.

**Confirm Allocation**

Places the memory allocations in this menu into effect. Exiting this menu without pressing <Confirm Allocation> will cancel any adjustments made in this menu.

Command: MALC i, j, k

---

**Trace to Arb**

Copy the real part of a stored trace into the Arbitrary Waveform memory. Select a Trace [1..5] and press [Enter]. The Arb Length is set to 2 kpoints and the Arbitrary Waveform Source is set to the Arbitrary Buffer.

Only traces which contain FFT measurements can be copied to the Arbitrary memory. In almost all cases, the trace should contain a time record.

If the length of the trace is less than 2 kPoints (2048 points), then the trace is repeated until 2 kPoints is reached. This can cause discontinuities if the trace is not continuous from its end to its beginning. If this source is measured with an FFT time record equal to the trace length (in time), then windowing should solve this problem.

## 4-82 Arbitrary Source Menu

---

Only the real part of the trace is used. Baseband time records are completely real and can be reproduced as an Arbitrary waveform. Non-baseband time records or complex traces do not reproduce well.

The amplitude of the Arbitrary Waveform is normalized to the maximum value in the Trace. The largest trace value is considered full scale (100% amplitude = 1V).

Be sure to set the Play Rate to the appropriate value to reproduce the frequencies in the trace correctly.

Command: TARB i

## Swept Sine Source Menu

The Source menu sets the source amplitude parameters when the Measurement Group is Swept Sine. See 'Swept Sine' in Chapter 2 for a discussion about swept sine measurement fundamentals.

Source
Auto Level Ref Channel 1
Ideal Ref 100.0 mV
Source Ramping Off
Source Ramp Rate 1 V/S
Ref Upper Limit 3.0 dB
Ref Lower Limit -3.0 dB
Max Source 1000.0 mV
Offset 0.0 mV

---

### Auto Level Reference

Select the Auto Level Reference mode [Off, Channel 1, Channel 2].

Auto Level Off maintains the source amplitude at a constant level at all frequencies during the swept sine sweep. This usually works best for transfer functions which are fairly flat.

Auto Level Channel 1 or 2 will adjust the source level to maintain a constant level, called the Ideal Reference, at the Channel 1 or Channel 2 input. Generally, this is useful for frequency response functions with substantial gain as well as attenuation. Input Auto Ranging is turned on for both inputs when Auto Level is on. Whenever Auto Level is used, the measurement **MUST** be frequency response.

The Auto Level Reference Channel determines which input will be maintained at the Ideal Reference Level. This is usually the device under test output (Channel 2). In some cases where the device under test is not driven directly from the source output, the device input (Channel 1) may be the reference.

Command: SSAL (?) {i}

---

### Amplitude

Set the sine Amplitude [0 mV .. 5000 mV]. This parameter is adjustable only if Auto Level Reference is Off.

The swept sine source turns off whenever there are no measurements being made. This is before the sweep is started, at the end of a single sweep or while a sweep is paused. If

## 4-84 Swept Sine Source Menu

---

Source Ramping is Off, the source shuts off instantly in these situations. If Source Ramping is On, the source will ramp off. When the sweep is started or re-started, the source will ramp on for the first measurement.

The amplitude may be changed at any time during a sweep.

Command: SSAM (?) {x}

---

### Ideal Reference

Set the Ideal Reference [0 mV .. 5000 mV]. This parameter is adjustable only if Auto Level Reference is set to Channel 1 or Channel 2.

The Ideal Reference is the signal level that the source maintains at the Reference Channel to within the Reference Limits. This is determined by the limitations of the device under test or may simply be the desired signal level for the test. If the Ideal Reference requires a source amplitude greater than the Maximum Source Level, then the Reference Channel signal will fall below the Ideal Reference. In cases where the required source amplitude is very small, the resolution of the source amplitude may result in a Reference Channel signal which is not within tolerance of the Ideal Reference.

The Ideal Reference may be changed at any time during a sweep.

Command: SSRF (?) {x}

---

### Source Ramping

Select the Source Ramping Mode [Off, On].

If Source Ramping is Off, source level changes are made instantly.

If Source Ramping is On, source level changes are made at the Source Ramp Rate. Settling starts after the source amplitude reaches the desired level. This prevents fast level changes from disturbing the device under test.

The swept sine source turns off whenever there are no measurements being made. This is before the sweep is started, at the end of a single sweep or while a sweep is paused. If Source Ramping is Off, the source shuts off instantly in these situations. If Source Ramping is On, the source will ramp off. When the sweep is started or re-started, the source will ramp on for the first measurement.

Command: SRMP (?) {i}

---

### Source Ramp Rate

Set the Source Ramp Rate [0.001 V/s .. 500 V/s].

The Source Ramp Rate is the rate at which the source amplitude changes when Source Ramping is On. If Source Ramping is Off, source amplitude changes are made instantly.



Note that very slow ramp rates can add significantly to the sweep time, especially if Auto Level is on and many changes to the source amplitude are required during the sweep. Also, the source ramps up when the sweep is started and ramps off when paused. This causes a delay before the first point is measured.

Command: SRAT (?) {x}

---

**Reference Upper Limit**

Set the Reference Upper Limit [0.1 dB .. 30.0 dB]. This parameter is adjustable only if Auto Level Reference is set to Channel 1 or Channel 2.

When Auto Level is on, the Ideal Reference is the signal level the source will try to maintain at the Reference Channel. The Reference Upper and Lower Limits are the allowable tolerances for the Reference Channel. The source level is changed only if the reference channel measures an input which exceeds the Ideal Reference by more than the Upper Limit or which is less than the Ideal Reference by more than the Lower Limit. In cases where the reference channel input must be limited, the Ideal Reference plus the Upper Limit must be less than or equal to this limit. If the required source amplitude is very small, the resolution of the source amplitude may result in a Reference Channel signal which is not within tolerance of the Ideal Reference.

The Reference Upper Limit may be changed during a sweep.

Command: SSUL (?) {x} or SSLL (?) {x}

---

**Reference Lower Limit**

Set the Reference Lower Limit [-0.1 dB .. -30.0 dB]. This parameter is adjustable only if Auto Level Reference is set to Channel 1 or Channel 2.

When Auto Level is on, the Ideal Reference is the signal level the source will try to maintain at the Reference Channel. The Reference Upper and Lower Limits are the allowable tolerances for the Reference Channel. The source level is changed only if the reference channel measures an input which exceeds the Ideal Reference by more than the Upper Limit or which is less than the Ideal Reference by more than the Lower Limit. If a source amplitude greater than the Maximum Source Level is required, then the Reference Channel signal may fall below the Reference Lower Limit.

The Reference Lower Limit may be changed during a sweep.

Command: SSLL (?) {x}

---

**Maximum Source Level**

Set the Maximum Source Level [0 mV .. 5000 mV]. This parameter is adjustable only if Auto Level Reference is set to Channel 1 or Channel 2.

## 4-86 Swept Sine Source Menu

---

The Maximum Source Level is the largest allowed source amplitude. This is limited by the SR785 source output or the device under test input range. If the Reference Lower Limit requires a source amplitude greater than the Maximum Source Level, then the source amplitude is set to the Maximum Source Level. In this case, the Reference Channel may fall below the Reference Lower Limit.

The Maximum Source Level may be changed at any time during a sweep.

Command: SMAX (?) {x}

---


### Offset

Set the DC Offset of the Swept Sine Source [-5V..5V]. The offset resolution is 0.1 mV. Note that the sum of the offset and Amplitude (if Auto Level Reference is off) or Maximum Source (if Auto Level Reference is on) cannot exceed 5 V. Also, using large offsets with small tone amplitudes will degrade the distortion performance of the sine source.

Command: SOFF (?) {x}

# Input Menu

The Input menu configures the analog signal inputs, the tachometer input, and sets the parameters for capture playback.

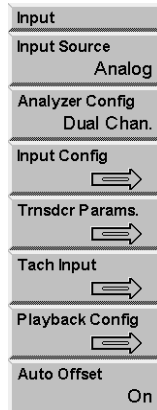

## Input Source

Select the measurement Input Source [Analog, Playback].

Analog uses the Ch1 and Ch2 front panel inputs as the source for all measurements. This menu adjusts the input configurations for both inputs.

Playback uses the data stored in the capture buffer as the input for all measurements. This selection is not valid until capture data has been acquired. This menu sets the playback parameters. Capture is not available for swept sine measurements.

Both displays use the same Input Source for their measurements. The Analog/Playback indicator shows the input source at the top of the screen.

Command: ISRC (?) {i}

## Analyzer Configuration

Select the analyzer configuration [Independent Channels, Dual Channel].

Independent Channels allows the two displays to have different frequency parameters. For instance, Display A could be set to full span while Display B showed a zoomed high resolution display of a single spectral feature. In Independent Channels mode, only single channel measurements are allowed. Dual channel measurements, such as frequency response, or orbit, cannot be selected.

Dual Channel allows selection of all measurements, both single and dual channel. However, in Dual Channel mode the two displays always have the same frequency parameters.

## 4-88 Input Menu

---

When the measurement group is Swept Sine or Order, the Analyzer Configuration is always set to Dual Channel.

Command: LINK (?) {i}

---

### Auto Offset


Enable Auto Offset input calibration [Off, On].

When Auto Offset is On, the analyzer will periodically perform an input offset calibration. This involves internally grounding both inputs and correcting their DC offsets. This only takes a second. However, the measurements will need to settle again after an offset calibration. To avoid disrupting a lengthy measurement, turn Auto Offset Off.

A warning message is displayed 15 seconds before an Auto Offset calibration is performed. Pressing [<-] (backspace) will abort the scheduled calibration.

Command: IAOM (?) {i}

## Input Configuration Submenu

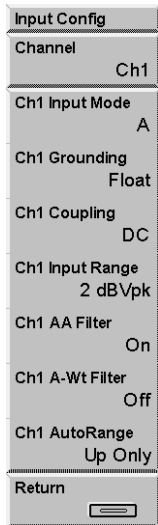

### Channel

Selects which input channels will be affected by the input configuration submenu [Ch1, Ch2, Both]. When Ch1 is selected, items on the menu will display and change parameters for the Ch1 input. When Ch2 is selected, items on the menu will display and change parameters for the Ch2 input. When Both is selected, the menu will display parameters for the Ch1 input but any changes made will affect both inputs.

### Input Mode

Select the Input Mode for the selected input[A (single-ended), A-B (differential)]. The Input indicator shows the current mode at the top of the screen.

In general, when looking at very small signals, connect A to the signal source and B to the signal ground and use A-B. In this case, make sure that the two input cables do not encompass any loop area (twist them together or run them side by side).

Command: I1MD (?) {i}

Command: I2MD (?) {i}

### Grounding

Select the Input Grounding for the selected input shields [Float, Ground].

The shields of A and B are connected.

Float connects the shields to chassis ground through  $1\text{ M}\Omega + 0.01\ \mu\text{F}$ . This allows the shield of a single ended input to 'float'.

## 4-90 Input Configuration Menu

---

Ground connects the shields to chassis ground with 50  $\Omega$ . In this mode, do not exceed 3 V on the shields.

The impedance between signal input and chassis ground is always 1 M $\Omega$ .

Command: I1GD (?) {i}

Command: I2GD (?) {i}

---

### Coupling

Select the Input Coupling for the selected input [DC, AC, ICP]. The Input indicator shows the current coupling at the top of the screen.

The 3 dB bandwidth of the AC coupling is 0.16 Hz.

ICP coupling connects a 5 mA current source (26 VDC open circuit) to the center conductor of the A input connector. This supply powers ICP accelerometers. The signal is AC coupled from the center conductor. Set the Input Mode to A (not A-B).

Command: I1CP (?) {i}

Command: I2CP (?) {i}

---

### Input Range

Select the Input Range for the selected input [-50 dBV .. +34 dBV]. The Input Range is the full scale signal input just before overload.

The actual underlying Input Range of the SR785 varies by 2 dB steps from -50 dBVpk to +34 dBVpk. The input range can be set in units of dBVpk, dBVrms, dBVpp, Vpk, Vrms, dBEUpk, dBEUrms, dBEUpp, EUpk, EURms, or EUpp. If a value is entered in one of these units the closest actual hardware input range is selected.

The Input Range indicator shows the current range at the top of the screen.

Pressing this softkey turns off Auto Range for the selected input channels

If the input signal exceeds 35 dBV, the input range is automatically set to 34 dBV. The range may not be changed while this HighV condition exists.

Command: I1RG (?) {i}

Command: I2RG (?) {i}

---

### AA Filter

Select the Anti-Aliasing Filter for selected input [Off, On].

The anti-aliasing filter should generally be left On. Frequency domain measurements may have spurious alias signals if the filter is Off.

For time domain measurements involving square waves or impulses, removing the filter can remove ringing and improve the accuracy of the time record.

The specifications apply only if the filter is On.

Command: I1AF (?) {i}

Command: I2AF (?) {i}

---

### **A-Wt Filter**

Select the A-Weighting Filter for the selected input [Off, On]. The AWt indicator at the top of the screen shows 'AWt' in bright whenever the A-Weighting filter is On.

The A-Weighting filter simulates the hearing response of the human ear and is often used with Octave Analysis measurements. The input A-Weighting filter conforms to ANSI standard S1.4-1983.

A-, B- and C- Weighting functions are also available as operators in User Math functions.

Command: I1AW (?) {i}

Command: I2AW (?) {i}

---

### **Autoranging**

Select the AutoRange Mode for the selected input [Up Only, Tracking]. This mode only applies when Ch1 Auto Range is On.

In Up Only, only overload causes the Input Range to change. The Input Range only moves up.

In Tracking Mode, the Input Range moves up for overloads and down when the signal falls below half scale. Some signals, such as low frequency noise, can cause the Input Range to oscillate. Do not use Tracking Mode in these cases.

Auto Range responds to all frequencies present at the input (except those attenuated by AC coupling or the anti-aliasing filter), not just those within the measurement span.

Command: I1AR (?) {i}

Command: I2AR (?) {i}

## Transducer Parameter Submenu

Display the Transducer Parameters menu. Press <Return> or [Input] for the main [Input] menu.

Transducers (such as accelerometers or microphones) convert a physical quantity (such as acceleration or pressure) into a voltage at the analyzer's input. By assigning Transducer Units to an input, measurements based upon the transducer signal may be displayed in units of the actual physical quantity being measured (such as  $m/s^2$  or Pascals).

Transducer Units require a transducer which is linear over the frequency range of interest.

Measurements of inputs which are assigned units of acceleration, velocity or displacement may be displayed with any of these three units using Transducer Convert in the [Input]<Transducer Parameter> submenu.

Transducer Units	
Channel	Ch1
Ch1 Eng. Units	On
Ch1 EU Label	$m/s^2$
Ch1 EU / Volt	1 EU/V
Ch1 EU at Marker	<input type="checkbox"/>
Ch1 User Label	EU
Transducer Convert	
	$m/s^2$
Return	<input type="checkbox"/>

## Engineering Units

Select Engineering Units for the selected input [On, Off].

Measurements which require input data from an input which has Engineering Units On will be displayed with units of EU Labels instead of Volts.

The two inputs may have different Engineering Units. For instance, if Ch2 has engineering units On and EU Label set to Pascals, and Ch1 has engineering units off, then the units of the frequency response measurement will be Pascals/Volt.

Note that this softkey applies to an input channel, unlike most softkeys which apply to a specific display.

Command: EU1M (?) {i}  
EU2M (?) {i}



## EU Label

Choose the Engineering Units Label for the selected input. The displayed units change from Volts to the EU Label when Engineering Units are On. Choose USER for the user defined User Label.

Input units of acceleration ( $m/s^2$ ,  $in/s^2$ , or  $g$ ), velocity ( $m/s$  or  $in/s$ ), and displacement ( $m$  or  $in$ ) may be converted in the frequency domain using Transducer Convert. For example, the FFT of an input configured for an accelerometer with Engineering Units of  $m/s^2$  may be displayed in displacement (meters) using Transducer Convert.

The two inputs may have different EU Labels. For instance, if Ch2 has EU Label set to  $m/s^2$ , and Ch1 has EU Label set to dyne, then the units of the frequency response measurement will be  $m/s^2/dyne$ .

Note that this softkey applies to an input channel, unlike most softkeys which apply to a specific display.

Command: EU1L (?) {i}

Command: EU2L (?) {i}

---

## EU / Volt

Set the Engineering Units per Volt scale factor for the selected input. This factor is used to scale the display units from Volts to Engineering Units when Engineering Units are On. The displayed units change from Volts to the EU Label.

The EU scale factor may be entered in EU/V, V/EU, or dB(EU/V).

Note that this softkey applies to an input channel, unlike most softkeys which apply to a specific display.

Command: EU1V (?) {x}

EU2V (?) {x}

---

## EUs at Marker

Calculate the EU/Volt scale factor for either the selected input based upon the marker reading in the active display. Enter the desired calibrated value for the marker reading (in the display units).

The active display measurement must be a single channel measurement and Engineering Units must be On.

This feature is normally used for the input channel which is being measured in the active display. In this case, the active display marker reading changes to the calibrated value.

## 4-94 Transducer Parameter Menu

---

However, you can use this feature for the other input channel. In this case, the active display is not affected but other measurements based upon the other input are affected.

---

### User Label

Enter the User Label for the selected input Engineering Units. Use the [ALT] keys (for alphabetic entry) to enter a user defined name for the engineering units. The displayed units change from Volts to the User Label when Engineering Units are On and USER is chosen for the EU Label.

Note that this softkey applies to an input channel, unlike most softkeys which apply to a specific display.

Command: EU1U (?) {s}  
EU2U (?) {s}

---

### Transducer Convert

Select the Transducer Conversion for the selected input. Transducer conversion only has an effect for measurements of inputs with Transducer Units of acceleration, velocity, or displacement.

Transducers (such as accelerometers) convert a physical quantity, (such as acceleration) into a voltage at the analyzer's input. By assigning Transducer Units to an input, measurements based upon the transducer signal (acceleration) may be displayed in units of the actual physical quantity being measured (velocity or displacement). Transducer conversion is accomplished in the frequency domain by multiplying and dividing by  $j\omega$  which is equivalent to differentiating and integrating in the time domain.

The two input channels may have different transducer conversions.

Command: TD1C (?) {s}  
EU2U (?) {s}

## Tachometer Input Submenu

This submenu sets parameters for the tachometer input. The tachometer input is used in the Order Measurement Group to synchronize the analyzer to a rotating machine. The tachometer input is also used by the trigger in the RPM Arming trigger mode.

Tach Input
Pulses per Rev 1
Tach Trig Rng ± 25 V
Tach Level 10.00 V
Tach Slope Rising
Hold Off Enable Off
Tach Hold Off 2 μs
Show Tach Off
Return 

---

### Pulses per Rev

Sets the number of tach pulses in each revolution of the device under test. The allowed range is 1 to 2048. Non-integer values are allowed.

Command: TAPR(?) {f}

---

### Tach Trigger Range

Specify the range over which the tachometer input can be set to trigger [5V, 25V, TTL]. If TTL is selected the tachometer will trigger on a TTL rising edge. If 5V is selected, the tachometer trigger level can be set from -5V to +5V with 40 mV resolution. If 25V is selected, the threshold resolution is 200 mV.

TTL Trigger range will yield the best phase accuracy.

Command: TARG(?) {i}

---

### Tach Level

Sets the tachometer trigger level when the trigger range is set to 5V or 25V. The resolution of the trigger level is 40 mV at the 5V range and 200 mV on the 25V range. The tachometer will trigger when the input passes through the specified level in the direction specified by the tach slope softkey.

If the tachometer signal is noisy near the specified tach trigger level, use the hold off to increase the noise immunity of the tach input.

## 4-96 Tachometer Input Menu

---

Command: TALV(?) {f}

---

### Tach Slope

Sets the slope with which the tachometer signal must pass through the tachometer level to trigger the tachometer [Rising, Falling].

Command: TASL(?) {i}

---

### Hold Off Enable

Enables or disables tachometer holdoff. When tachometer holdoff is enabled, the tachometer is prevented from triggering after a tach pulse has been received for a time specified by the Tach Hold Off.

If the tachometer signal is not noisy in the region near the tachometer trigger threshold, it should not be necessary to use the holdoff. If the signal is noisy, i.e. passes through the threshold several times for a given tachometer pulse, set the holdoff time long enough to prevent retriggering for each tach pulse.

Command: TAHO(?) {i}

---

### Tach Hold Off

Specify the interval during which the tachometer input is prevented from triggering when tachometer hold off is enabled.

For a noisy tachometer signal, the tachometer holdoff time should be set long enough so that the tachometer triggers only once for each pulse on the tachometer input.

Command: TAHD(?) {f}

---

### Show Tach

Specifies whether the RPM corresponding to the tachometer input should be displayed in the RPM indicator location.

When the measurement group is Order, Show Tach is always On. Likewise, if the trigger mode is RPM Arm, Show Tach is always on.

If Show Tach is off, the indicator is displayed in the same position as the RPM indicator.

Command: TASH(?) {i}

## Playback Submenu

This menu sets the parameters for capture playback.




---

### Playback Start

Select the Capture Playback Start position within the capture buffer. The start position is specified as a bin number or individual point position. The bins are separated by a time equal to  $1/\text{Capture Sampling Rate}$ . The corresponding start time relative to the start of the buffer is displayed as well.

Both channels playback starting at the Playback Start.

Command: ISTR (?) {i}

---

### Playback Length

Select the Capture Playback Length in 2 kPoint (2048 points) increments. The Playback Length (in points) plus the Playback Start cannot exceed the Capture Length of the buffer. The corresponding playback time is displayed as well.

Both channels playback the same Playback Length.

Command: ILEN (?) {i}

---

### Set Left Edge

The active display must be a Capture Buffer measurement for this key to be active.

## 4-98 Playback Input Menu

---

Choose a Capture Buffer as the Measurement and use Zoom and Pan to show the region of interest. Press <Set Left Edge> to set the Capture Playback Start to the marker position. The Playback Length is not changed.

If the requested Start position plus the Playback Length (in points) exceeds the Capture Length of the buffer, the Playback Length is reduced.

---

### Set Right Edge

The active display must be a Capture Buffer measurement for this key to be active.

Choose a Capture Buffer as the Measurement and use Zoom and Pan to show the region of interest. Press <Set Right Edge> to set the Capture Playback Length to the marker position minus the Playback Start. The Playback Start is not changed.

If the marker position is to the left of the Playback Start, then the Playback Length is set to the minimum (2 kPts).

---

### Playback Mode

Select the Capture Playback Mode [1-Shot, Circular].

In 1-Shot playback, the buffer is played a single time. The playback halts when the end of the buffer is reached.

In Circular playback, the playback starts over when the end of the buffer is reached. There is often a discontinuity in the playback measurement when the playback jumps from the end to the start of the buffer.

In either mode, press [Start/Reset] to restart playback at the start again.

The Capture Progress indicator at the top of the screen shows the playback progress through the buffer.

Both channels playback with the same Playback Mode.

Command: IMOD (?) {i}

---

### Playback Speed

Select the Capture Playback Speed [Normal, Every Time Record].

Normal Playback Speed plays back the capture buffer in real time, i.e. 1 second of capture takes 1 second to playback. Normal Speed playback looks like the 'live' real time measurement. Normal playback is limited to the real time limitations of the equivalent real time analog input measurement.

Not all time records are displayed during Normal playback though all time records contribute to averaged measurements. For example, 1 second of capture contains 256 full

span FFT time records. Normal Speed playback at full span takes 1 second and updates the display only about 8 times during playback.

Every Time Record playback displays the measurement result for every captured time record. Since the data is stored in memory, Every Time Record playback is not limited by real time considerations. For example, 1 second of capture contains 256 full span FFT time records. In this case, all 256 time records are measured AND displayed. The display still updates at 8 Hz so playback takes about 32 seconds to complete. If the time records are overlapped, there may be more than 256 measurements to display and playback will take even longer.

When the measurement time records are very long (narrow spans), the time to process and display each measurement is much less than the real time record length. Since the data is already available in the capture buffer, Normal Speed playback means waiting unnecessarily for a real time record to elapse between updates. In this case, Every Time Record playback displays the measurements of all captured time records much faster than Normal Speed real time playback.

When the playback is in Octave Group, playback is always Normal Speed.

The Capture Progress indicator at the top of the screen shows the playback progress through the buffer.

Both channels playback with the same Playback Speed.

Command: `ISPD (?) {i}`

# Trigger Menu

The Trigger menu configures the Trigger Mode, Source, Level and Delay.



### Watch Out For Triggered Sources

If a triggered source is selected (Chirp, Burst Chirp, Burst Noise or Arbitrary), the triggered FFT measurement phase is stable only if the input signals are derived from the source output. Turn the source off (or set it to Sine) when making triggered measurements without the source.

Trigger
Trigger Mode Auto Arm
Trigger Source Ch 1
Trigger Level 7.09 %
Trigger Slope Rising
Delay 1 0 s
Delay 2 0 s
Trigd. Source Mode Continuous
RPM/Time Setup 
Manual Arm 
Manual Trigger 

---

## Trigger Mode

Select the Trigger Arming Mode [Auto Arm, Manual Arm, RPM Arm, Time Arm]. The trigger must be armed before a trigger event will be recognized.

Auto Arm arms the trigger automatically whenever possible, without waiting for an arming event. The combination of Auto Arm Trigger Mode and the Continuous Trigger Source is often known as 'Free Run' triggering since the analyzer will take data as quickly as possible.

In Manual Arm, the trigger is armed only after <Manual Arm> is pressed, or the TARM interface command is received. Once the trigger is armed, the next trigger event will trigger the measurement.

In RPM Arming, the trigger is armed when the RPM passes through a start threshold and at fixed rpm increments thereafter.

In Time Arming, the trigger is armed after a specified amount of time has elapsed. The time interval is specified with the Time Arm Step softkey.



Triggers are ignored unless measurements are started (with [Start/Reset]) and the trigger is armed. Triggers are ignored while paused.

For FFT, Correlation, Time/Histogram, and Order measurements, the trigger triggers individual time records.

For octave and swept sine measurements, the first trigger after arming starts the measurement and subsequent triggers are ignored. Use [Start/Reset] to reset the measurement and start again.

For FFT, Correlation, Octave, and Time/Histogram modes, the trigger can be configured to trigger a group of linear averages instead of individual time records with the Trigger Avg. Mode Softkey.

Command: TMOD (?) {i}

---

## Trigger Source

Select the Trigger Source [Cont, Ch1, Ch2, External, Ext TTL, Source, Manual].

Continuous trigger automatically generates a trigger as soon as the trigger is armed. The combination of Auto Arm Trigger Mode and the Continuous Trigger Source is often known as 'Free Run' triggering since the analyzer will take data as quickly as possible.

Ch1 or Ch2 internally trigger on the Ch1 or Ch2 input signal. The trigger is detected after the anti-aliasing filter (if On).

External and Ext TTL trigger on the front panel Trigger Input. External has a variable Trigger Level and Ext TTL triggers on TTL signals.

Manual triggers with <Manual Trigger> or an interface command. Use Auto Arm Trigger Mode with Manual Trigger since unwanted triggers are not a problem.

External, Ext TTL and Manual also trigger the burst Chirp, Noise, and arbitrary source waveforms.

Source triggers synchronously with the source waveform. The source runs continuously and Source Trigger synchronizes the time record with the source waveform. Source Trigger is not valid for swept sine measurements.

For FFT measurements, the phase of signals which are synchronous with the trigger is stable. Vector averaging preserves these signals while attenuating random signals, thus increasing the dynamic range of the measurement.

For Order measurements, the phase is measured relative to the trigger rather than the tachometer pulse if a non continuous trigger source is used.

The Trigger Source applies to both displays.

Command: TSRC (?) {i}

## 4-102 Trigger Menu

---

### Trigger Level

Set the Trigger Level [-100% .. 99%]. The Trigger Level applies to Ch1, Ch2 and External Trigger Sources. The Trigger Level is specified as a percentage of the Input Range for Ch1 or Ch2 and as a percentage of 5 V for External trigger.

The trigger detector requires a minimum signal amplitude of 4% of the Input Range (200 mV for External). For Ch1 or Ch2 internal trigger, the signal must exceed -28 dBfs in order to trigger. Internal trigger is detected after the anti-aliasing filter (if On).

Command: TLVL (?) {i}

---

### Trigger Slope

Select the Trigger Slope [Rising, Falling]. The Trigger Slope applies to Ch1, Ch2, External and Ext TTL Trigger Sources.

Command: TSLP (?) {i}

---

### Delay 1

Set the Trigger Delay for Ch1, or Display A. depending on the Analyzer Configuration. The Trigger Delay applies only to FFT and Correlation measurements.

When the Analyzer Configuration is Dual Channel, the delay setting affects the Ch1 input. When the Analyzer Configuration is Independent Channels, the delay setting affects the measurement displayed in Display A.

The triggered time record does not have to start with the trigger event. The time record can start before the trigger (negative delay) or after the trigger (positive delay). A delay of 0 starts the time record at the trigger. The Trigger Delay does not apply if the Trigger Source is Continuous.

The trigger delay can range from -8,000 to +100,000 times the underlying FFT sample rate.

When a large trigger delay is used, the display may update more slowly since the acquisition time for each record (time record length plus trigger delay) is noticeably long.

Command: TDLA (?) {i}

---

### Delay B

Set the Trigger Delay for Ch1, or Display A. depending on the Analyzer Configuration. The Trigger Delay applies only to FFT and Correlation measurements.

When the Analyzer Configuration is Dual Channel, the delay setting affects the Ch2 input. When the Analyzer Configuration is Independent Channels, the delay setting affects the measurement displayed in Display B.

The triggered time record does not have to start with the trigger event. The time record can start before the trigger (negative delay) or after the trigger (positive delay). A delay of 0 starts the time record at the trigger. The Trigger Delay does not apply if the Trigger Source is Continuous.

The trigger delay can range from -8,000 to +100,000 times the underlying FFT sample rate.

When a large trigger delay is used, the display may update more slowly since the acquisition time for each record (time record length plus trigger delay) is noticeably long.

Command: TDLB (?) {i}

---

## Trigd Source Mode

Set the trigger mode for the source [1-Shot, Continuous]. The Chirp, Burst Chirp, Burst Noise and Arbitrary sources are all triggered sources. This box is unavailable if the source is not a triggered source.

When the trigger is set to Auto-Arm, Continuous, the source is output continuously and this menu box is ignored.

### Triggering

The triggered sources are triggered by External triggers.

In the 1-Shot case, the source only triggers on the first trigger after [Start/Reset]. This same trigger starts the measurement (and waterfall). This is useful when the source length is very long compared to the measurement time or many measurements are required to determine the response to a single source output. This allows a sequence of many measurements to be triggered at the start of a single triggered source output.

In the Continuous case, the source triggers on every trigger, with a minimum trigger period set by the source length (FFT time record, noise source period or arbitrary source length). This is useful when the source length is the same as the FFT time record length.

Do not use Ch1 or Ch2 input trigger since the output will not start until a trigger is received.

Source trigger outputs the source continuously over and over. FFT time records are synchronized to the source start.



### Watch Out For Triggered Sources

If a triggered source is selected (Chirp, Burst Chirp, Burst Noise or Arbitrary), the triggered FFT measurement phase is stable only if the input signals are derived from the

## 4-104 Trigger Menu

---

source output. Turn the source off (or set it to Sine) when making triggered measurements without the source.

Command: STMD (?) {i}

---

### RPM/Time Arm Setup

Select the RPM/Time Arm menu to view and adjust parameters relating to RPM and Time arming.

RPM/Time Arm	
Start RPM	50
Start RPM	<input type="radio"/> Off
Delta RPM	10
Delta RPM Sense	Abs. Change
Time Arm Step	100 ms
Return	<input type="button" value="Return"/>

---

#### Start RPM

Specify the threshold RPM for RPM arming. Start RPM is only active if Start RPM (On/Off) is On.

For the trigger to initially arm, the RPM must exceed the threshold RPM (if Delta RPM Sense is Increasing), or must be below the threshold (if Delta RPM Sense is Decreasing.). If Delta RPM sense is set to Abs. Change, the Start RPM is ignored.

Command: TRSR (?) {f}

---

#### Start RPM (On/Off)

Activates the Start RPM threshold. If Start RPM is set On, the trigger will not be armed until the RPM passes the threshold. If Start RPM is Off, the threshold setting is ignored.

If Start RPM is On, the RPM threshold must be met for both the RPM arming and Time Arming trigger modes.

Command: TRSM (?) {i}

---

#### Delta RPM

Set the RPM increments for which the trigger will be armed. Once the Start RPM condition, if any, has been met. The trigger will arm each time the rpm changes by the specified amount.

If Delta RPM Sense is set to 'Incr. RPM' the trigger will arm only on an increasing RPM change. If it is set to 'Decr. RPM' the trigger will arm only on a decreasing RPM change. If set to 'Abs. Change' the trigger will arm on any change by this amount.

Command: TRDR (?) {f}

---

**Delta RPM Sense**

Sets the sign of the rpm change which arms the trigger in the RPM Arming trigger mode.

The trigger can be set to arm on positive, negative, or any change in RPM of an amount set by the Delta RPM softkey.

Command: TRDM (?) {i}

---

**Time Arm Step**

Sets the time interval for the Time Arming trigger mode. Once the Start RPM condition, if any, has been met, the trigger will arm repetitively with this time interval.

Command: TIAS (?) {f}

---

**Manual Arm**

Manually arm the trigger. This function only applies if the Trigger Mode is Manual Arm.

Manual Arm enables triggering after [Start/Reset] only after the trigger is armed using <Manual Arm> or from an interface command. Once the trigger is armed, the next trigger event will trigger the measurement.

For FFT, correlation, and Time Histogram measurements each time record requires the trigger to be re-armed. For octave and swept sine measurements, the first trigger after arming starts the measurement and subsequent triggers are ignored.

Command: TARM

---

**Manual Trigger**

Trigger the measurements (and burst Chirp & Noise sources) regardless of Trigger Source.

Manual Trigger Source requires <Manual Trigger> or an interface command to trigger. Use Auto Arm Trigger Mode with Manual Trigger Source since unwanted triggers are not a problem.

## 4-106 Trigger Menu

---

For FFT, correlation, and Time Histogram measurements each time record requires the trigger to be re-armed. For octave and swept sine measurements, the first trigger after arming starts the measurement and subsequent triggers are ignored.




Command: TMAN

## Average Menus

The Measurement Group determines which Average menu is displayed. This menu sets the averaging parameters for the measurement. A set of parameters is stored for each group. No averaging parameters are available in the Time/Histogram measurement group.

### FFT, Correlation, Order Average Menu

When the Measurement Group is FFT, Correlation, or Order, this menu sets the averaging parameters for the active display. If the Analyzer Configuration is Dual Channel, the averaging parameters are the same for both displays. Averaging successive measurements together improves accuracy and repeatability.

Averaging	
Compute Avgs	Yes
Averaging Type	Linear / Fix Len
# Avgs	10
Display Avg	Vector
Time Record Incr.	100.00%
More	
Avg Preview	Manual
Preview time	2 s
Accept	
Reject	

### Compute Averages

Turns the computation of averages on and off. If Compute Averages is off, no averaged quantities will be computed or displayed. Although this results in a slight improvement in the speed of some measurements, for normal operation Compute Averages should be left On.

Some FFT measurements, such as coherence, have averaging as part of their definition. If one of these measurements is selected, Compute Averages will be forced On.

Command: FAVG (?) d {, i}

### Averaging Type

Select the Averaging Type (weighting). Choose Linear/Fixed Length or Exponential/Continuous.

The [Start/Reset] key resets the current average and starts a new average. The [Pause/Cont] key pauses the average in progress. Pressing [Pause/Cont] again will continue the average from where it was paused.

If the Analyzer Configuration is 'Independent Channels,' each display can have its own Averaging Type. This entry field can be linked to both displays by using the [Link] key.

Command: FAVT (?) d {, i}

### Linear/Fixed Length Weighting

Linear weighting combines N (Number Of Averages) measurements with equal weighting in either RMS or Vector averaging.

While Linear averaging is in progress, the number of averages completed is shown in the Horizontal Scale Bar below the graph. When the Number Of Averages has been completed, the measurement stops and 'Done' is displayed along the y-axis of the graph. If the 'Done' alarm is enabled, an audible alert is also sounded.

When used with peak hold averaging, this weighting mode is known as Fixed Length averaging. This means that a fixed number, N, of measurements are compared to determining the peak.

### Waterfall Storage

If Waterfall Storage is on, the waterfall buffer only stores the completed linear averages, not each individual measurement. Each time the linear average is done, the result is stored in the waterfall buffer and the average is reset and started over (instead of stopping). Each completed average counts as a single waterfall record.

### Exponential/Continuous Weighting

Exponential weighting weights new data more than old data. For RMS and Vector averaging, weighting takes place according to the formula,

$$\text{Average}[N] = (\text{New Data})(1/N) + (\text{Average}[N-1])(N-1)/N$$

where N is the Number Of Averages.

While Exponential averaging is in progress, the number of averages completed is shown in the Horizontal Scale Bar below the graph. The displayed number stops incrementing at the Number of Averages while the averaging continues.

Exponential weighting reaches a steady state after approximately 5N measurements. Once in steady state, further changes in the average are detected only if they last for a sufficient number of measurements. Make sure that the number of averages is not so large as to eliminate changes in the data which might be important.



When used with peak hold averaging, this weighting mode is known as Continuous averaging. This means that new measurements are continually being compared to the current maximum to determine the new peak.

---

## Number of Averages

Set the Number Of Averages for the active display [2..32767].

The Number Of Averages specifies the number of measurements for Linear/Fixed Length averages or the weighting of new data in Exponential/Continuous averages.

The [Start/Reset] key resets the current average and starts a new average. The [Pause/Cont] key pauses the average in progress. Pressing [Pause/Cont] again will continue the average from where it was paused.

Changing the Number Of Averages does not start averaging.

If the Analyzer Configuration is 'Independent Channels,' each display can have its own Number of Averages. This entry field can be linked to both displays by using the [Link] key.

Command: FAVN (?) d {, i}

---

## Display Average

Select the type of averaging for the measurement displayed on the active display [None, RMS, Vector, Peak Hold]. Note that when Compute Averages is On, the SR785 always computes the non-averaged, rms averaged, vector averaged, and peak-hold averaged version of all measurements simultaneously. The Display Average softkey only determines which version of the measurement is currently displayed. You can change Display Average to view a different average of the current measurement without retaking any data.

Some measurements, such as coherence, have their own unique averaging mode. In these cases, the Display Average softkey is ignored.

The [Start/Reset] key resets the current average and starts a new average. The [Pause/Cont] key pauses the average in progress. Pressing [Pause/Cont] again will continue the average from where it was paused.

Each display has its own Display Average. This entry field can be linked to both displays by using the [Link] key.

Command: FAVM (?) d {, i}

## RMS Averaging

For a simple FFT measurement, the definition of the RMS averaged measurement is:

$$\text{RMS Avg (FFT)} = \langle \text{conj(FFT)} \cdot \text{FFT} \rangle = \langle X^2 + Y^2 \rangle$$

Where X and Y are the real and imaginary parts of the instantaneous FFT measurements and the  $\langle \rangle$  angle brackets indicate either linearly weighted or exponentially weighted averaging over N FFT records. For other measurements, the definition of the RMS averaged measurement depends upon the measurement. In the above example, the RMS averaged FFT is real, and has zero phase. In general however, rms averaged measurements have real and imaginary parts and have non-zero phase.

Linear averaging computes the equally weighted mean over N measurements. Exponential averaging weights new data more than old data and yields a continuous moving average.

RMS averaging reduces fluctuations in the data but does not reduce the actual noise floor (squared values never cancel). With a sufficient number of averages, a very good approximation of the actual noise floor can be measured.

### Vector Averaging

For a simple FFT measurement, the definition of the Vector average measurement is:

$$\text{Vector Avg(FFT)} = \langle X \rangle + j \langle Y \rangle$$

Where X and Y are the real and imaginary parts of the instantaneous FFT measurements and the  $\langle \rangle$  angle brackets indicate either linearly weighted or exponentially weighted averaging over N FFT records. For other measurements, the definition of the vector averaged measurement depends upon the measurement.

Vector averaging is often called time averaging, since vector averaging in the frequency domain yields the same value as taking the FFT of the averaged input time record.

Linear averaging computes the equally weighted mean of X and Y over N measurements. Exponential averaging weights new data more than old data and yields a continuous moving average.

Since signed values are combined in the mean, random signals tend to average to zero. This reduces the noise floor since random signals are not phase coherent from measurement to measurement. Signals with a constant phase have real and imaginary parts which repeat from time record to time record and are preserved. Vector averaging can substantially improve the dynamic range of a measurement as long as the signals of interest have stable phases.

For single channel measurements, vector averaging requires a trigger. The signal of interest MUST be phase synchronous with the trigger to have a stable phase.

For a two channel measurement, the phase is relative between Channel 2 and Channel 1. As long as the signals of interest have stable relative phases, triggering is not required for vector averaging. Triggering is still required to isolate time records which contain the signals of interest.

The Time Record Increment should be set to 100% when vector averaged measurements are being used.

### Peak Hold Averaging

Peak hold averaging is similar to RMS averaging, in that the RMS measurement quantities are calculated. However, instead of averaging the RMS measurements together, in peak-hold averaging the new data is compared to the old data and the maximum value is kept.

In Continuous Peak-Hold Averaging, the new data is continually compared with the current maximum. In Fixed-Length PeakHold Averaging, a fixed number of records is examined for the maximum.

---

### Time Record Increment

Set the Time Record Increment for the active display [0% .. 300%]. The Time Record Increment is how far the start of each time record is advanced between measurements (as a percentage of the time record length).

If the increment is 100%, the start of the next time record is exactly one time record advanced from the start of the previous time record.

If the increment is 25%, then the next time record starts 1/4 of a time record advanced from the start of the previous time record. This is sometimes referred to as 75% overlap since the two time records share 75% of a record. The overlap is simply 100% minus the Time Record Increment.

When the Time Record Increment is less than or equal to 100%, the measurement is said to be 'real time'. All time points contribute to one (or more) measurements.

If the increment is 200%, the start of the next time record is advanced from the start of the previous time record by 2 time records leaving a gap of 1 record. This means that the data between the two time records is not measured.

When the Time Record Increment is greater than 100%, then the measurement is not 'real time' and some time points do not contribute to a measurement.

The actual time record increment for the measurements in progress is displayed in the Real Time indicator in the status area (below the Input Ranges). If the indicator shows a value greater than the requested Time Record Increment, it means that the measurement cannot be made with the requested increment but is running with the smallest increment possible. Factors which affect the processor's ability to run real time include the measurement type, averaging and source type.

Be sure to set the Time Record Increment to 100% when vector averaged measurements are being used.

The Time Record Increment is ignored whenever the measurement is triggered. In this case, the time records start with the trigger.

### Settling

When the frequency span or input signal path (gain, filtering, etc.) is changed, a settling time is required before the FFT measurement is considered settled or valid. The

## 4-112 Average Menus

---

measurement is not settled until the discontinuity in the input data has propagated through the digital filters and a complete new time record has been acquired.

If the time record increment is 100%, unsettled measurements are not displayed. After a change is made which unsettles the measurement, new data is not displayed until the filters are settled and a complete time record has been acquired.

If the measurement is running with a time record increment less than 100% and the measurement is unsettled, unsettled measurements may be displayed. New data is displayed after the filters are settled and a portion of the new time record has been acquired. For example, if the time record increment is 25%, 3 measurements are made before a complete new time record has been acquired. These first 3 measurements have time records which contain data from before AND after the measurement was unsettled. These unsettled measurements are displayed in half intensity indicating that the measurement is not settled. Once the 4th measurement is complete, the display returns to full intensity since the measurements have completely settled time records.

If averaging is on, changes which unsettle the measurement will restart the average. Unsettled measurements are not included in the new average. Averaging does not start until the measurement is settled. When the time record increment is less than 100%, 'Settle' is displayed (instead of the number of averages) below the graph to indicate that the unsettled measurements are not being averaged.

See also Real Time Bandwidth and Overlap.

If the Analyzer Configuration is 'Independent Channels,' each display can have its own Time Record Increment. This entry field can be linked to both displays by using the [Link] key.

Command: FOVL (?) d {, x}

---

### More

This submenu selects additional averaging parameters.



---

### Overload Reject

Select Overload Reject for both displays [Off, On].

When Overload Reject is On, time records which contain points which were overloaded at the input are not included in the measurement average. Whenever a time record is rejected, the 'Reject' indicator below the graph (next to the average type and number) turns on briefly. This has no effect if Averaging is Off.

When Overload Reject is Off, all time records are part of the measurement average. Overloads can greatly disturb the average so be certain that the correct Input Range is used.

Command: FREJ (?) d {, i}

---

### Trigger Average Mode

Selects the triggering mode for FFT, Correlation, and Time/Histogram measurements [Time Record, Start Average].

Normally, when Trigger Average Mode is set to Time Record each trigger triggers a single time record in each of these measurement groups. So if the Averaging Type is set to Linear and the Number of Averages is set to 50, 50 triggers will be required to complete the average.

When Trigger Average Mode is set to Start Average, the first trigger will trigger the first time record, and subsequent time records will be acquired as quickly as possible until the linear average is complete. In the above example, only one trigger would be required to acquire all 50 time records and complete the average.

Command: TAVM (?) d {, i}

---

### Average Preview

Select Average Preview for both displays [Off, Manual or Timed].

When averaging is on, Average Preview allows each individual measurement to be accepted (added to the average) or rejected (not added to the average) based upon the input time records. This is useful for rejecting bad time records from corrupting an averaged measurement.

When Average Preview is Manual, each input time record is displayed after it is acquired. Choose <Accept> or <Reject> to either use or ignore the measurement based upon the displayed time record. After the time record is accepted or rejected, the display reverts back the measurement display until the next time record is available (usually triggered).

When Average Preview is Timed, the input time record is automatically accepted after the Preview Time (unless rejected first).

The preview time records are displayed in place of the actual measurements. For single channel measurements, the appropriate channel input time record is shown in each display. For 2-channel measurements, the DisplayA shows Ch1 time record and DisplayB shows Ch2 time record. In this case, use the dual display format to view both time records.

While the preview time records are displayed, the displays are labeled 'Preview' instead of 'Live' (in the Vertical Scale Bar). While previewing the time record, the measurement

## 4-114 Average Menus

---

can be changed so as to preview any other measurement in the current measurement group. Subsequent time records will return to the newly selected measurement rather than the time record for preview.

Command: PAVO (?) d {, i}

---

### Preview Time

Select the Preview Time for both displays [0.5 s .. 1 ks].

When Average Preview is Timed, the input time record is automatically accepted after the Preview Time (unless rejected first).

Command: PAVT (?) d {, x}

---

### Accept

Accept the displayed preview time record and add the measurement to the average.

Average Preview allows each individual measurement to be accepted (added to the average) or rejected (not added to the average) based upon the input time records. This is useful for rejecting bad time records from corrupting an averaged measurement.

Command: PAVA

---

### Reject

Reject the displayed preview time record and do not add the measurement to the average.

Average Preview allows each individual measurement to be accepted (added to the average) or rejected (not added to the average) based upon the input time records. This is useful for rejecting bad time records from corrupting an averaged measurement.

Command: PAVR

# Octave Average Menu

When the Measurement Group is Octave Analysis, this menu sets the averaging parameters for the active display.

Averaging	
Averaging Type	Exp. Time
Integration Time	100 ms
Confidence	0.125 dB
Power Bin	Total
Linear Avg. Trig	One Lin. Avg.

---

## Averaging Type

Select the Averaging Type for the active display [Linear Time, Exponential Time, Peak Hold, Equal Confidence].

Choosing Peak Hold also sets the Power Bin to Peak.

Octave measurements are always rms averaged to measure the power in each band.

The [Start/Reset] key resets the current averages and starts the measurement over. The [Pause/Cont] key pauses the measurement. Pressing [Pause/Cont] again resets the averages and starts the measurement over.

The Averaging Type of each display can be selected separately. This entry field can be linked to both displays by using the [Link] key.

Command: OTYP (?) d {, i}

### Linear Time

The band filter outputs are equally weighted and averaged for an Integration Time.

While Linear averaging is in progress, the integration time completed is shown in the Horizontal Scale Bar below the graph. When the Integration Time has been completed, the measurement stops and 'Done' is displayed below the graph. If the 'Done' alarm is enabled, an audible alert is also sounded.

### Waterfall Storage

If Waterfall Storage is On, the waterfall buffer only stores the completed linear averages, not each individual measurement.

When the measurement is free running (Auto Arm, Continuous Trigger), each time the linear average is done, the result is stored in the waterfall buffer and the average is reset and started over (instead of stopping). Each completed average counts as a single waterfall record.

When the measurement is triggered, each linear average may be triggered as well. Set the Linear Avg Mode to choose triggered averages or automatically reset averages.

### **Exponential Time**

New filtered data is weighted more than older data. The exponential 'time constant' is the Integration Time. Averaging continues indefinitely.

While Exponential averaging is in progress, the completed integration time is shown in the Horizontal Scale Bar below the graph. The displayed time stops incrementing at the Integration Time while the averaging continues.

Exponential weighting reaches a steady state after approximately an integration time. Once in steady state, further changes in the average are detected only if they last for a sufficient number of measurements. Make sure that the integration time is not so large as to eliminate changes in the data which might be important.

### **Equal Confidence**

Equal Confidence averaging is exponential averaging where the integration time is set for each band separately. The integration times are set so that there is a 68% probability that the results are within the specified Confidence Level of the true mean for every band in the measurement. There is a 96% probability that the results are within twice the Confidence Level (in dB) of the true mean.

The effect of Equal Confidence averaging is that higher frequency bands have shorter time constants and lower frequency bands have longer time constants. This is very noticeable in the update of the display.

Averaging continues indefinitely.

### **Peak Hold**

Peak Hold displays the peak output from each band filter. The peak detection continues indefinitely.

Choosing Peak Hold also sets the Power Bin to Peak.

---

## **Integration Time**

Set the Integration Time for the active display [4 ms .. 1 ks] in 4 ms increments. The Integration Time affects Linear Time, Exponential Time, and Peak Hold averaging.

The minimum Integration Time is determined by the Lowest Band in the Octave display. 4 ms Integration Time is only available if the Lowest Band is 100 Hz or greater. Decreasing the Lowest Band below 100 Hz increases the minimum value of the Integration Time.

The [Start/Reset] key resets the current averages and starts the measurement over. The [Pause/Cont] key pauses the measurement. Pressing [Pause/Cont] again resets the averages and starts the measurement over.



Each display has its own Integration Time. This entry field can be linked to both displays by using the [Link] key.

Command: OTIM (?) d {, x}

---

**Confidence**

Set the Confidence Level for the active display [0.125 dB .. 2 dB]. The Confidence Level only affects Equal Confidence averaging.

The [Start/Reset] key resets the current averages and starts the measurement over. The [Pause/Cont] key pauses the measurement. Pressing [Pause/Cont] again resets the averages and starts the measurement over.

Each display has its own Confidence Level. This entry field can be linked to both displays by using the [Link] key.

Command: OCNF (?) d {, i}

---

**Power Bin**

Choose the type of broadband power measurement for the active display [Total, Impulse, L, Peak]. If the Averaging Type is Peak Hold, only Peak power may be selected. If the Averaging Type is not Peak Hold, Total, Impulse and L may be selected.

In addition to the octave bands, a broadband sound level is measured and displayed as the last band in the octave graph. Exponential averaged sound power (Leq) is calculated according to ANSI S1.4-1983, Type 0. Broadband Impulse and Peak measurements are made according to IEC 651-1979, Type 0.

Exponential averaged power (Leq), Impulse (I) and Peak (P) power are computed from real-time input data. These are broadband measurements and are not computed from the octave bands. The band-limited Total (T) power is computed from the octave bands.

Each display has its own Power Bin. This entry field can be linked to both displays by using the [Link] key.

Command: OIMP (?) d {, i}

---

**Linear Average Trigger**

Choose how linear averaged octave measurements behave when triggered [One Lin. Avg., Start]. The Linear Avg Mode only applies for linear averaging when the Trigger Mode and Trigger Source are not Auto Arm and Continuous, respectively.

When the measurement is triggered the Linear Average Mode determines how the measurement restarts.

One Linear Average means that each trigger starts a new measurement.

## 4-118 Octave Average Menu

---

Start means that the next measurement starts immediately after the previous measurement finishes. Only the first measurement after [Start/Reset] is started with a trigger. This mode is useful for measuring reverberation where the first measurement is triggered with the noise burst (use the Burst Noise source with 1-Shot Trigd Source Mode) and subsequent measurements follow continuously. All complete averages are stored to a waterfall to examine the sound level decay.

Both displays always have the same Linear Average Trigger.

Command: OLAT (?) d {, i}

## Swept Sine Average Menu

When the Measurement Group is Swept Sine, this menu sets the settle and integration times. These parameters govern the measurement of both channels. See ‘Swept Sine’ in Chapter 2 for a discussion about swept sine measurement fundamentals.

Averaging	
Settle Time	20 ms
Settle Cycles	5
Integration Time	100 ms
Integration Cycles	10

---

### Settle Time

Set the Settle Time in increments of 3.906 ms [7.8125 ms .. 1 ks]

At each frequency point, a settling time is allowed to pass before any measurement is made. This allows the device under test to respond to the frequency change. The Settle Time is rounded up to the next multiple of 3.906 ms. The actual settling time is the larger of the Settle Time and the Settle Cycles (rounded to the next multiple of 3.906 ms). The settling time is always a minimum of 1 cycle or 7.8125 ms.

Changes made to the Settle Time during a sweep take effect immediately.

The estimated sweep time is displayed in the Horizontal Scale Bar. This time is simply the sum of the Settle and Integrate times for all points in the sweep. Auto functions (Source Auto Level, Auto Range, Auto Resolution) will change the actual sweep time.

Command: SSTM (?) d {, x}

---

### Settle Cycles

Set the number of Settle Cycles [1..32767].

At each frequency point, a settling time is allowed to pass before any measurement is made. This allows the device under test to respond to the frequency change. The actual settling time is the larger of the Settle Time and the Settle Cycles (rounded to the next multiple of 3.906 ms). The settling time is always a minimum of 1 cycle or 7.8125 ms.

Changes made to the number of Settle Cycles during a sweep take effect immediately.

The estimated sweep time is displayed in the Horizontal Scale Bar. This time is simply the sum of the Settle and Integrate times for all points in the sweep. Auto functions (Source Auto Level, Auto Range, Auto Resolution) will change the actual sweep time.

## 4-120 Swept Sine Average Menu

---

Command: SSCY (?) d {, i}

---

### Integration Time

Set the Integration Time in increments of 3.906 ms [15.625 ms .. 1 ks].

At each frequency point, the inputs measure the signal at the source frequency. This is done by multiplying the input data by the source sine (and cosine) wave and averaging the results over an integration time. The actual integration time is always an exact number of cycles at the source frequency. This rejects signals which are at a different frequency, such as noise and harmonics. A long integration time results in a narrow detection bandwidth at the source frequency. This improves signal to noise at the cost of longer measurement times .

The Integration Time is converted to next largest number of exact cycles. The actual number of integration cycles is the larger of the Integration Time (in cycles) and the Integration Cycles. The integration time is always a minimum of 1 cycle or 15.625 ms.

To measure each point for the same amount of time (constant detection bandwidth), set the Integration Cycles to 1 and the Integration Time to the desired time. The Integration Time should be greater than 1 cycle of the lowest frequency in the sweep.

Changes made to the Integration Time during a sweep take effect immediately.

The estimated sweep time is displayed in the Horizontal Scale Bar. This time is simply the sum of the Settle and Integrate times for all points in the sweep. Auto functions (Source Auto Level, Auto Range, Auto Resolution) will change the actual sweep time.

Command: SITM (?) d {, x}

---

### Integration Cycles

Set the number of Integration Cycles [1..32767].

At each frequency point, the inputs measure the signal at the source frequency. This is done by multiplying the input data by the source sine (and cosine) wave and averaging the results over an integration time. The actual integration time is always an exact number of cycles at the source frequency. This rejects signals which are at a different frequency, such as noise and harmonics. A long integration time results in a narrow detection bandwidth at the source frequency. This improves signal to noise at the cost of longer measurement times .

The integration time is an exact number of cycles. The actual number of integration cycles is the larger of the Integration Time (in cycles) and the Integration Cycles. The integration time is always a minimum of 1 cycle or 15.625 ms.

To measure each point for a time inversely proportional to the frequency, set the Integration Time to 15.625 ms (minimum) and the Integration Cycles to the desired

number. Remember, the detection bandwidth increases with frequency in this case (the cycles get shorter) which may result in increased detected noise at higher frequencies.

Changes made to the number of Integration Cycles during a sweep take effect immediately.

The estimated sweep time is displayed in the Horizontal Scale Bar. This time is simply the sum of the Settle and Integrate times for all points in the sweep. Auto functions (Source Auto Level, Auto Range, Auto Resolution) will change the actual sweep time.


Command: SSCY (?) d {, i}

## User Math Menu

The User Math menu sets up user defined measurements and math calculations.

Use a User Function to define measurements which are not already available in each Measurement Group. To make a measurement of a defined User Function, select the User Function as the Measurement for the active display.

See 'User Math Functions' in Chapter 2 for more.

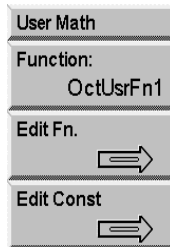

---

### Function

To edit a User Functions, use the knob to choose a function and press [Enter]. Only the 5 User Functions within the current Measurement Group may be edited. The current definition of the functions is displayed as the knob moves through the list.

---

### Edit Function

Display the Edit Function menu to edit the selected User Function. Press <Cancel> to exit back to the main [User Math] menu .

To enter an equation, use <Operands> and <Operations> and the knob to select operands and operations. [-] (Backspace) deletes the term before the cursor.

To edit an equation, use <Func String> and the knob to move the flashing cursor within the equation. <Delete> will delete the highlighted term. <Operands> and <Operations> will either insert or replace at the cursor. Use <Insert/Replace> to switch between insert and replace mode. <Clear Eq> will delete the entire equation.

<Cancel> will discard any changes made in this menu and exit this menu.

<Enter Eq> will enter the new function equation and exit this menu. The new User Function definition is used immediately if it is being displayed.

- Command: USRF (?) i { , j , k , l , m , ... } (FFT Group)
- Command: USRR (?) i { , j , k , l , m , ... } (Correlation Group)
- Command: USRO (?) i { , j , k , l , m , ... } (Octave Group)
- Command: USRS (?) i { , j , k , l , m , ... } (Swept Sine Group)
- Command: USRT (?) i { , j , k , l , m , ... } (Order Group)
- Command: USRH (?) i { , j , k , l , m , ... } (Time/Histogram Group)



### Operands

Use the knob to pick one of the displayed Operands and press [Enter] to place it in the equation at the cursor location. The display then automatically switches to the Operations display. To enter another operand instead, press <Operands> again.

The available measurement operands depend upon the current Measurement Group.

### FFT Measurement Group

Time(1)	Time(2)	FFT(1)	FFT(2)	< Pwr(F1) >
< Pwr(F2) >	Coherence()	< Cross Spec >	< Freq Resp >	
Vec< F1 >	Vec< F2 >	RMS< F1 >	RMS< F2 >	
PeakHold<F1>	PeakHold<F2>	Rms< xSpec >	j	
Const1	Const2	Const3	Const4	Const5
Trace1	Trace2	Trace3	Trace4	Trace5

**Correlation Measurement Group**

Time(1)	Time(2)	Auto_Corr(1)	Auto_Corr(2)	X_Corr()
FFT(1)	FFT(2)	FFTu(1)	FFTu(2)	Vec<Fu1><F1>
Vec<Fu2><F2>	Vec<Fu1><F2>	Rms<Fu1 F1>	Rms<Fu2 F2>	Rms<Fu2 F1>
Const1	Const2	Const3	Const4	Const5
Trace1	Trace2	Trace3	Trace4	Trace5
j				

**Octave Measurement Group**

Oct(1)	Oct(2)	j		
Const1	Const2	Const3	Const4	Const5
Trace1	Trace2	Trace3	Trace4	Trace5

**Swept Sine Measurement Group**

Spec(1)	Spec(2)			
Norm_Var(1)	Norm_Var(2)	Cross_Spec()	Freq_Resp()	j
Const1	Const2	Const3	Const4	Const5
Trace1	Trace2	Trace3	Trace4	Trace5

**Order Measurement Group**

Time(1)	Time(2)	<Spec(1)>	<Spec(2)>	
<Pwr(1)>	<Pwr(2)>	Vec<T1>	Vec<T2>	j
Const1	Const2	Const3	Const4	Const5
Trace1	Trace2	Trace3	Trace4	Trace5



**Time/Histogram Measurement Group**

Hist(1)	Hist(2)	PDF(1)	PDF(2)	
CDF(1)	CDF(2)	Time(1)	Time(2)	
Const1	Const2	Const3	Const4	Const5
Trace1	Trace2	Trace3	Trace4	Trace5
j				

Operands are arrays of complex values (real and imaginary). They may be measurement results (such as FFT(1), Time(1) or Oct(1)), a Trace, or a Constant. An array which is real simply has zero for its imaginary parts. The array length of an Operand is determined by the measurement length (number of FFT lines, length of time record, number of octaves, etc.) or Trace length (length of the data which is stored in the Trace). Constants assume the length of the user function.

Operands which are measurement results enclosed in angle brackets, such as <Freq. Resp.> or <Spec(1)>, are exactly the same as the normal measurements. They are computed using whatever averaging is selected with the <Display Average> softkey. They are computed from the input data (real time Analog or capture Playback) and use the frequency and windowing parameters from the menus. Frequency domain measurements are amplitude calibrated, time records are not. FFT and Time record measurements are also triggered just like the normal measurements.

Operands which are measurement results which are not enclosed in angle brackets, such as FFT1 or Time1 represent unaveraged, instantaneous versions of the measurement.

Operands which contain an explicit averaging type, such as Vec<F1> or PeakHold<F2> always are averaged according to their indicated type, regardless of the setting of the <Display Average> softkey.

Octave and Swept Sine measurements are always averaged measurements. Use the [Average] menu to set the averaging parameters.

FFT(1) and FFT(2) are the FFT of the Ch1 and Ch2 inputs. These operands use the window chosen in the [Window] menu for the display which is measuring the function.

FFTu(1) and FFТу(2) are un-windowed FFT's of the Ch1 and Ch2 inputs.

Correlation operands such as RMS<Fu1 F1> are used to compute the RMS and Vector averaged auto and cross correlation measurements. See chapter 2 for a description of how to compute all the averaged versions of all the predefined measurements.

Trace operands are simply the data stored in the Traces. For example, Traces can hold reference data used for normalization or calibration. There are 5 Traces which can be stored. These Traces are shared by all 3 Measurement Groups.

## 4-126 User Math Menu

---

Constants are complex constants which are the same for every array point. Constants such as pi are commonly used in equations. There are 5 user defined constants which are shared by all 3 Measurement Groups.

---

### Operations

Use the knob to pick one of the displayed Operations and press [Enter] to place it in the equation at the cursor location. The display then automatically switches to the Operands display. To enter another operation instead, press <Operations> again.

The available operations depend upon the current Measurement Group.

#### FFT Measurement Group

+	-	x	/	(
)	Mag()	Phase()	Real()	Imag()
Mag <sup>2</sup> ()	Sqrt()	Conj()	Ln()	Exp()
FFT()	FFTu()	IFFT()	jOmega()	d/dx()
GrpDly()	AWt()	BWt()	CWt()	[X / (1-X)]()

#### Correlation Measurement Group

+	-	x	/	(
)	Mag()	Phase()	Real()	Imag()
Mag <sup>2</sup> ()	Sqrt()	Conj()	Ln()	Exp()
FFT()	FFTu()	IFFT()	jOmega()	d/dx()
GrpDly()	AWt()	BWt()	CWt()	[X / (1-X)]()

#### Octave Measurement Group

+	-	x	/	(
)	Mag()	Phase()	Real()	Imag()
Mag <sup>2</sup> ()	Sqrt()	Conj()	Ln()	Exp()
FFT()	FFTu()	IFFT()	d/dx()	
AWt()	BWt()	CWt()	[X / (1-X)]()	

**Swept Sine Measurement Group**

+	-	x	/	(
)	Mag(	Phase(	Real(	Imag(
Mag <sup>2</sup> (	Sqrt(	Conj(	Ln(	Exp(
jOmega(	d/dx(	GrpDly(	AWt(	BWt(
CWt(	[X / (1-X)](			

**Order Measurement Group**

+	-	x	/	(
)	Mag(	Phase(	Real(	Imag(
Mag <sup>2</sup> (	Sqrt(	Conj(	Ln(	Exp(
FFT(	FFTu(	IFFT(	d/dx(	
AWt(	BWt(	CWt(	[X / (1-X)](	

**Time/Histogram Measurement Group**

+	-	x	/	(
)	Mag(	Phase(	Real(	Imag(
Mag <sup>2</sup> (	Sqrt(	Conj(	Ln(	Exp(

Arithmetic operations (+, -, x, /) combine operand terms on a point by point basis. It is the user's responsibility to ensure that the operand terms have the correct X axis type and lengths in order to produce meaningful results.

View operations ( Mag( ), Mag<sup>2</sup>( ), Phase( ), Real ( ), Imag( ) ) simply convert the complex operand array into the desired form. Mag( ) and Mag<sup>2</sup>( ) compute the magnitude ( $\sqrt{x^2 + y^2}$ ) or magnitude squared ( $x^2 + y^2$ ) with a real result. Phase( ) computes the phase ( $\tan^{-1}( y/x )$ ) unwrapped with a real result. Real( ) simply zeroes the imaginary part. Imag( ) zeroes the real part.

Conj( ), Ln( ), Exp( ) and Sqrt( ) are defined as follows,

$$\begin{aligned}
 \text{Conj}(x + jy) &= x - jy \\
 \text{Conj}(x - jy) &= x + jy \\
 \text{Ln}(x + jy) &= \ln(r) + j\theta \\
 \text{Exp}(x + jy) &= \exp(x) \cdot (\cos(y) + j\sin(y)) \\
 \text{Sqrt}(x + jy) &= \sqrt{r} \cdot (\cos(\theta/2) + j\sin(\theta/2)) \\
 [X/(1-X)] &= (x + jy)/(1 - x - jy)
 \end{aligned}$$

where  $r = \sqrt{x^2 + y^2}$  and  $\theta = \tan^{-1}(y/x)$ . Sqrt computes all angles positive from 0. The  $[X/(1-X)]()$  operator is generally used on real arguments only.

The  $j\Omega$  operator is simply  $j\omega$  and zeroes the real part and fills the imaginary part with the value of  $\omega$  ( $2\pi \cdot$  frequency) at each point in the array. The frequency is determined by the span of the operand. A measurement operand (FFT(1) or Time(1)) uses the current frequency span. A Trace operand uses the span of the measurement which is stored in the Trace.

FFT() is the windowed FFT operator. These operations switch the X axis type. The current window chosen in the [Window] menu is used on the operand before the FFT. FFTu() is the un-windowed FFT operator. The uniform window is always used. The FFT() and FFTu() operators do not calibrate their results since their operands are not necessarily input time records. To use a calibrated FFT, use the measurement operand FFT(1) or FFT(2).

IFFT() is the inverse complex FFT operation. This operation switches the X axis type.

d/dx() is the derivative operator. The derivative is performed with respect to the X-axis bin number, not the x-axis values. For example, to convert d/dx() to d/df() for a linear frequency display, divide by the frequency spacing of the bins. In the case where the X-axis is logarithmic, the function needs to correct for  $d(\log x)/dx$ . The d/dx() operator requires a smoothing aperture or window.

GrpDly() is the group delay operator. Group delay is simply  $d\theta/d\omega$ . Use [Display Options]<d/dx Window> to set the aperture. This operator yields valid results for any complex operand which has a frequency x-axis.

AWt() is the A-Weighting operator. This filter is a software version of the input A-Weight filter. BWt() and CWt() are the software B-Weighting and C-Weighting operators. The AWt(), BWt(), and CWt() operators should only be used on operands with a frequency x-axis. All three filters conform to ANSI Standard S1.4-1983.

---

### Function String

Move the cursor to the function equation window at the top of the screen. Use the knob to move the highlighted cursor within the equation. <Delete> will delete the highlighted term. <Operands> and <Operations> will either insert or replace at the cursor. Use <Insert/Replace> to switch between insert and replace mode.

[<-] (Backspace) deletes the term before the cursor.

---

### Insert/Replace

Toggle between insert and replace mode while editing a User Function equation.

If editing in insert mode, 'Ins' appears in the upper right corner of the edit window. If editing in replace mode, 'Rep' appears.

When editing a User Function equation, selecting <Operands> or <Operations> will insert or replace before the cursor location.

---

**Delete**

Delete the highlighted term in a User Function equation. To move the cursor in a User Function equation, use <Func String> and the knob.

---

**Clear Eq.**

Clear the entire equation definition. If a cleared User Function is entered, it will be undefined and unavailable for display. If a display is currently measuring this function, then it may not be cleared. In this case, change the measurement and then clear the equation.

---

**Cancel**

Discard any changes made in this menu and exit this menu.

---

**Enter Eq.**

Enter the displayed equation as the new User Function and exit this menu. The new User Function definition is used immediately if it is being displayed.

If the new function uses an incompatible or empty Trace it may not be displayed as a measurement. If a display is already measuring this function, then it may not be entered. In this case, change the measurement and then change the equation.

---

**Edit Constant**

Display the Edit Constant menu. Press <Return> or [User Math] for the main [User Math] menu .

Constants may be used as operands in a User Function. There are 5 constants which may be defined and are shared by all Measurement Groups.

Command: USRC (?) i { , x, y }

Edit Const	
Constant	Const1
Real Part	1
Imag Part	0
Mag	1
Phase (deg)	0
Marker->Mag	<input type="button" value="=&gt;"/>
Return	<input type="button" value="=&gt;"/>

---

### Constant

Use the knob to pick one of the five Constants to edit and press [Enter].

---

### Real Part

Enter a new value for the real part of the selected Constant. The imaginary part is left unchanged and the magnitude and phase are updated.

Command: USRC (?) i {, x, y}

---

### Imaginary Part

Enter a new value for the imaginary part of the selected Constant. The real part is left unchanged and the magnitude and phase are updated.

Command: USRC (?) i {, x, y}

---

### Mag

Enter a new value for the magnitude of the selected Constant. The phase is left unchanged and the real and imaginary parts are updated.

---

### Phase(deg)

Enter a new value for the phase (in degrees) of the selected Constant. The magnitude is left unchanged and the real and imaginary parts are updated.


---

### Marker to Mag

Copy the marker reading (exactly as shown in the Marker Position Bar) of the active display to the magnitude of the selected Constant. No unit translation takes place. Make sure that the units of the active display are correct before using this feature.

# Window Menu

The Window menu allows the user to choose the window function for FFT measurements. See ‘Windowing’ in Chapter 2 for more information.

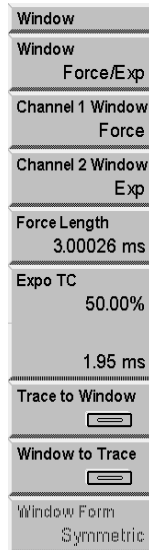

## Window

Select the Window Function for the active display.

If the Analyzer Configuration is set to Independent Channels, the two displays can have different Windows. If the Analyzer Configuration is set to Dual Channel, both displays use the same Window.

The following windows are available for FFT and Order measurements.

Command: FWIN (?) d { , i }

### Uniform

The time record is used with no windowing. This window provides amplitude accuracy only for exact bin frequencies and very poor frequency selectivity. This window is a poor choice for continuous signals. In general, this window is only useful when looking at transients which do not last the entire time record. The uniform window may also be used with signals which are exactly periodic in the time record such as a chirp.

The uniform window is also useful in the order analysis measurement group when the harmonics of interest are integer multiples of the shaft rotation speed. This guarantees that these order fall on exact bin frequencies.

### Hanning

The Hanning window is the most commonly used window. However, it has an amplitude variation of about 1.5 dB for off-bin signals and provides only reasonable selectivity. Its side-lobes are very high and broad for off-bin frequencies. As a result, the Hanning

window can limit the performance of the analyzer when looking at signals close together in frequency and very different in amplitude. The Hanning window is most often used in noise measurements since it has the lowest noise floor.

### **Flattop**

The Flattop window has the best amplitude accuracy of any window. Its off-bin amplitude variation is about 0.02 dB. However, the selectivity is worse. Unlike the other windows, the Flattop window has a very wide pass band and very steep rolloff on either side. Thus, signals appear wide but do not leak across the whole spectrum. The Flattop window is the best window to use for accurate amplitude measurements.

### **BMH**

The BMH window has reasonable off-bin amplitude accuracy (about 0.8 dB) and much lower side-lobes than the Hanning window. Thus, it has very good selectivity and very little broadening of non-bin frequencies. The BMH window is a good window to use for measurements requiring a large dynamic range.

### **Kaiser**

The Kaiser window has the lowest side-lobes and least broadening for non-bin frequencies. This makes this window the best for selectivity. The off-bin amplitude variation is the same as the BMH (about 0.8 dB). The Kaiser window is the best window to use for measurements requiring a large dynamic range.

### **Force-Exponential**

The Force window is uniform over the beginning of the time record and zero over the remainder. The force length is user specified. This window is used to isolate impulsive signals, such as impact excitations, from noise and other oscillations later in the time record.

The Exponential window attenuates the time record with a decaying exponential time constant. This window is often used in impact testing on the response channel to remove oscillations which last longer than the time record.

Many impact measurements require the Force window for the excitation channel and the Exponential window for the response channel. Use the <Channel 1 Window> and <Channel 2 Window> softkeys to select which channel to apply the force and exponential windows to.

### **User Defined**

The User window is any function that the user provides. The User window is copied from a stored trace using <Trace to Window>. The trace may contain stored data or may be loaded from disk or via the computer interface.

Remember, window functions have a great deal of impact on the resulting FFT spectrum. A poorly designed window can result in significant measurement errors.



### **Correlation Windows:**

When the measurement group is Correlation, the available windows are:

#### **[-T/2..T/2]**

This is a uniform window which should only be used on data which is self windowing (lasts less than half of the time record).

#### **[0..T/2]**

This window zeroes the second half of the time record.

#### **[-T/4..T/4]**

This window zeroes the first and last quarter of the time record.

Octave, Swept Sine, and Time/Histogram measurements do not use windows.

---

### **Channel 1 Window**

Selects whether the Force or Exponential window is applied to the Channel 1 input. This softkey is only active if the selected window is Force/Exponential.

Command: W1FE (?) {, i}

---

### **Channel 2 Window**

Selects whether the Force or Exponential window is applied to the Channel 2 input. This softkey is only active if the selected window is Force/Exponential.

Command: W2FE (?) {, i}

---

### **Force Length**

Set the Force Window Length. The maximum value that can be entered is the length of the FFT time record. Points in the time record up to the Force Length are unmodified. Points in the time record past the Force Length are zeroed.

Choose Force/Exponential as the window before setting the Force Length.

Command: FWFL (?) d {, i}

---

### **Expo TC**

Set the Expo Window Time Constant as a percentage of the FFT time record [5%..1000%]. This is the point where the window function reaches 1/e.

Choose Force/Exponential as the window before setting the Expo Time Constant.

Command: FWTC (?) d {, i}

### Trace to Window

Copy a stored data trace (real part only) to the User window. Use the knob to pick a trace which already contains data and press [Enter]. To see the effects of the new window the User window must be selected.

The trace data is interpolated if necessary in order to provide the window function with the correct number of points. The trace length is mapped to the entire window length. The User window may be used with time records of any length.

Command: TRWI d, i

---

### Window to Trace

Copy the current window of the active display to the real part of a stored data trace. Use the knob to pick a trace and press [Enter].

<Window to Trace> always stores a trace with 2048 points. This corresponds to the longest possible time record (800 line baseband). Using this trace in User Math requires attention to the lengths of other operands in any user function.

Command: WITR d, i

---

### Window Form

Select the Window Form of the User window [Non-Symmetric, Symmetric]. This softkey only affects the User window.



Use Symmetric if the User window is centered or symmetric in the time record. Phase will be measured relative to the center of the time record. Uniform, Hanning, Flattop, BMH and Kaiser windows are all symmetric.

Use Non-Symmetric if the User window is non-symmetric. Phase will be measured relative to the start of the time record. Force and Exponential windows are non-symmetric.

Command: WSYM (?) d {, i}

# Waterfall Menu

The Waterfall menu configures the waterfall displays. Waterfall displays are not available for Nyquist or Nichols Views or Swept Sine measurements.

Waterfall	
Display	Normal
Storage	Off
Save Option	All
Total Count	139
Skip	30
View Count	10
Angle	-30°
Marker Z to	10
Allocate Mem.	
More...	

## Display

Select the Normal or Waterfall Display for the active display. Waterfall Display only applies to FFT and Octave measurements. Waterfall displays are not available for Nyquist or Nichols Views or Swept Sine measurements.

Normal Display only shows the results of a single measurement.

Waterfall Display WITHOUT storage is simply a graphical display scrolling. The marker is constrained to the most recent record (at the back of the waterfall). Older records are only shown graphically, their data has not been saved and may not be accessed in any way. In this case, the display updates and scrolls as quickly as possible.

Waterfall display WITH storage shows only those records which are stored in waterfall memory. While the measurement is running, the display scrolls down and the newest record is added to the top (back). While running, the display may not update fast enough to show every stored record in real time. Pause the measurement to redraw the display with either the newest record at the top (normal) or the oldest record at the top (reversed). When paused, the display shows every stored record which is visible. The number of records stored in waterfall memory is displayed in the Vertical Scale Bar when storage is on.

A waterfall display may not be taken Off-Line.

### Moving the Marker

The marker may be moved to a record other than 0 only if storage is on.

## 4-136 Waterfall Menu

---

Use [Alt] knob to move the marker from record to record (along the Z axis) and scroll the display through the waterfall buffer. Use <Cursor Z to> to move to a specified record. The marker position along the Z axis is displayed in the Marker Position Bar next to the marker frequency or time. The Z axis position is the record number starting from 0 (most recently stored at the back of the waterfall).<Record to Trace> saves the record selected by the marker to a Trace.

Use the knob without [Alt] to move the marker along the X axis within a record. <Slice to Trace> saves a 'slice' to a Trace. A slice is the history of a single X position (data at the marker X position from all stored records).

The Waterfall Display for each display can be selected separately. This entry field can be linked to both displays by using the [Link] key.

Command: WDSP (?) d {, i}

---

### Storage

Select the type of Waterfall Storage [Off, Continuous, One Shot]. The two displays always have the same Waterfall Storage mode. Waterfall storage is not available for Swept Sine measurements.

Storage Off does not store any measurements to waterfall memory and changes the Waterfall Display to a graphical scroll only.

Continuous storage adds measurement records to the waterfall memory until [Pause/Cont] is pressed. The last Total Count number of records are saved in memory and may be displayed.

One shot stores the Total Count number of measurement records and pauses the waterfall storage. The first Total Count number of records are saved in memory. Pause the measurements to scroll through the waterfall memory.

Waterfall memory must be allocated before waterfall storage may be used.

#### Linear Averaging

If Waterfall Storage is On, the waterfall buffer only stores the completed linear averages, not each individual measurement. Each time the linear average is done, the result is stored in the waterfall buffer and the average is reset and started over (instead of stopping). Each completed average counts as a single waterfall record.

Command: WSTO (?) d {, i}

---

### Save Option

Select whether the analyzer should store all the measurements in the current measurement group, or only the currently displayed measurements.

When All is selected the SR785 stores all the measurements in the current measurement group to the waterfall buffer. This allows later viewing of any measurement, but uses up more space in the waterfall buffer.

When Active Measurement Only is selected, the SR785 stores only the currently selected measurement for each display. This takes less space in the waterfall buffer and allows a higher value for Total Count, but only allows the currently selected measurements to be displayed.

Command: WFSB (?) d{, i}

---

**Total Count**

Set the Total Count (number of records) to store in the waterfall memory for the active display.

The maximum Total Count depends upon the allocated memory and the setting of Save Option>. For example, if 500 blocks of allocated waterfall memory are allocated, and the <Save Option> is set to Active Measurement each display will store 250 400 line FFT measurements. If <Save Option> is changed to Save All the buffer will be able to hold only 142 measurements. Changing frequency resolution may reduce the Total Count since more buffer memory is required to hold measurements with higher resolution. .

The Total Count is limited to 2048 records (per display) regardless of the amount of installed memory.

The rate at which measurement records are added to waterfall memory is set by the Waterfall Skip. This allows the waterfall memory to hold a longer time history as well as setting a variable storage rate.

Command: WTOT (?) d {, i}

---

**Skip/Storage Interval**

Set the Waterfall Skip/Storage Interval for the active display. The rate at which measurement records are added to waterfall memory is set by the Skip/Storage Interval. This allows the waterfall memory to hold a longer time history as well as setting a variable storage rate.

The interval has no effect if Storage is Off. The Waterfall Display updates at the normal rate for the display measurement and span.

In the FFT, Correlation, Order, and Time/Histogram groups, for every measurement added to memory, a skip number of measurements are not stored. For example, a skip of 0 stores every measurement in memory. A skip of 10 stores every 11th measurement in memory (store 1, skip 10). The skipped measurements are still computed and affect (exponential) averaging, they are simply not stored. The elapsed time between stored records (1/storage rate) is simply skip+1 times the FFT acquisition time times the time

## 4-138 Waterfall Menu

---

record increment. (In the case of triggered measurements, skip+1 times the trigger period). This allows the time between stored records to be determined exactly.

In Octave group, the Storage Interval is set as a time. In this case, a snapshot is stored to memory every Storage Interval amount of time (4 ms minimum with 4 ms resolution).

The Skip/Storage Interval has a slightly different effect if the measurement is linear averaged. In this case, the result of every Nth complete linear average is stored in the waterfall buffer and the average is reset and started over again. Change the Number of Averages (FFT, et. al.) or Integration Time (Octave) to change the storage rate.

Depending upon the storage rate, the waterfall display may not scroll fast enough to show every record being saved. In this case, the displays are redrawn with all visible records shown when the measurement is paused.

Command: WFSK (?) d {, i} or WOSK (?) d {, x}

---

### View Count

Set the Waterfall View Count for the active display [2..2048]. The View Count is the number of records shown within the display. This number is independent of the Total Count (number of records stored). The View Count is the Z axis scaling of the waterfall display and should be less than the Total Count in order to use the display area effectively.

Setting the View Count to a large number (>40 or so) will cause the display to skip some records. In this case, every nth record is drawn in order to show the desired View Count range of the waterfall buffer. A large View Count is a good way to display an overview of the entire buffer in order to identify a region of interest. Remember, View Count is strictly a display parameter, it does not affect the actual storage of waterfall records.

The View Count and the Trace Height determine the available scroll Angles.

Each display has its own Waterfall View Count. This entry field can be linked to both displays by using the [Link] key.

Command: WVCT (?) d {, i}

---

### Angle

Set the requested Waterfall Scroll Angle for the active display.

Waterfall displays show multiple measurement results (records) in a single display. New records are added at the top of the display and older records scroll off the bottom of the display. When the Scroll Angle is 0, the older records scroll straight down. When the Scroll Angle is negative, the older records shift left as they scroll down. When the Scroll Angle is positive, they shift right as they scroll down.

The View Count and the Trace Height determine the available scroll angles. When Fast Angles is On, the available scroll angles are limited to those which scroll faster (due to the graphics architecture). This menu box sets the requested scroll angle. The display scrolls at the allowed angle closest to the requested angle. When adjusted with the knob, only available angles are selected.

Each display has its own Waterfall Scroll Angle. This entry field can be linked to both displays by using the [Link] key.

Command: WANG (?) d {, i}

---

**Marker to Z**

Move the marker to a specified record within the waterfall buffer in the active display. Storage must be on to move to records other than the most recent (record 0). Enter a record number (starting from 0 as the most recent) and press [Enter]. The waterfall display will re-draw to include the specified record number.

[Alt] knob also moves the marker from record to record (along the Z axis) and scrolls the display through the waterfall buffer. <Marker to Z> is a faster way to move to an arbitrary record.

The marker position along the Z axis is displayed in the Marker Position Bar next to the marker frequency or time. The Z axis position is the record number starting from 0 (most recently stored at the back of the waterfall).<Record to Trace> saves the record selected by the marker to a Trace.

Command: MWFL d, i

---

**Allocate Memory**

Display the Memory Allocation menu.

**Total Available**

Displays the total memory available for storage, either 980, 2004 or 4022 blocks, depending upon the amount of installed memory. Each block is 2 kPoints (2048 points).

The total of the Capture, Waterfall and Arbitrary memory allocations cannot exceed the Total Available memory. It may be necessary to decrease one allocation in order to increase another.

Command: MMEM ?

**Capture Memory**

Allocates memory blocks for the capture buffer. The allocated Capture Memory sets the limit for the Capture Length. Capture memory must be allocated before the capture buffer may be used.

## 4-140 Waterfall Menu

---

### Waterfall/ Order Memory

Allocates memory blocks for waterfall storage and tracking orders. Waterfall memory must be allocated before waterfall displays may be used.

### Arb. Memory

Allocates memory blocks for the Arbitrary Source waveform. The Arbitrary Source can also use the capture buffer as the output waveform.

### Clear Allocation

Clears the existing memory allocations in this menu.

### Confirm Allocation

Places the memory allocations in this menu into effect. Exiting this menu without pressing <Confirm Allocation> will cancel any adjustments made in this menu.

Command: MALC i, j, k

---

## More

Display the More Waterfall menu. Press <Return> or [Waterfall] for the main Waterfall] menu.

Waterfall	
Trace Height	70%
Fast Angles	Off
Threshold	0%
Hidden Lines	Invisible
Paused Drawing	Normal
Record to Trace	<input type="checkbox"/>
Slice to Trace	<input type="checkbox"/>
Return	<input type="checkbox"/>

---

## Trace Height

Set the Waterfall Trace Height for the active display [20%..80%]. The Trace Height specifies the percentage of the display height occupied by the Y axis of a single record.

The View Count and the Trace Height determine the available scroll Angles.

Each display has its own Waterfall Trace Height. This entry field can be linked to both displays by using the [Link] key.

Command: WHIT (?) d {, i}



### Fast Angles

Select the Waterfall Fast Angles for the active display [On, Off].

The View Count and the Trace Height determine the available scroll Angles. When Fast Angles is On, the available Scroll Angles are limited to those which update faster (due to the graphics architecture).

Each display has its own Fast Angles mode. This entry field can be linked to both displays by using the [Link] key.

Command: WFST (?) d {, i}

---

### Threshold

Set the Waterfall Threshold for the active display [0%..99%]. The Threshold is specified as a percentage of the Y axis (Trace Height). Only those portions of each record which exceed the Threshold are displayed. This removes baseline noise from the waterfall display.

Note that the marker still moves along the actual data in each record, even if it is below the threshold. In this case, the marker will be located where there is no displayed data. The Marker Position Bar displays the actual data value in all cases.

Each display has its own Waterfall Threshold. This entry field can be linked to both displays by using the [Link] key.

Command: WTHR (?) d {, i}

---

### Hidden Lines

Set the Waterfall Hidden Line Mode for the active display [Invisible, Visible]. In a waterfall display, as the older records scroll downward, they are considered to move towards the 'front' of the display. New records which appear at the top are 'behind' the older records. Invisible Hidden Lines will remove the portions of the new record which are 'behind' older records. This cleans up a waterfall display by removing intersecting lines. If records at the front of the waterfall are obscuring details hidden behind them, reverse the waterfall direction.

Each display has its own Waterfall Hidden Line Mode. This entry field can be linked to both displays by using the [Link] key.

Command: WHID (?) d {, i}

---

### Paused Drawing

Set the waterfall display direction for a paused measurement [Normal, Oldest at Top].

While the measurement is running, the waterfall display scrolls down and the newest record is added to the top (back) of the waterfall.

## 4-142 Waterfall Menu

---

When the measurement is paused, the waterfall is redrawn with either the newest record at the top (Normal) or reversed (Oldest at Top). If records at the front of the waterfall are obscuring details 'hidden' behind them, reverse the waterfall direction.

Each display has its own Paused Drawing Mode. This entry field can be linked to both displays by using the [Link] key.

Command: WREV (?) d {, i}

---

### Record to Trace

Save a single record or measurement from the waterfall to a Trace. <Record to Trace> saves the record selected by the marker to a Trace. Choose a Trace with the knob and press [Enter] to complete the operation.

A stored trace can be recalled to a Display or Reference Display, used in a User Math Function, saved to disk, or copied to the Arbitrary Waveform buffer.

When the measurement is paused, use [Alt] knob to move the marker from record to record (in the Z axis) and scroll the display through the waterfall memory. Use the knob without [Alt] to move the marker along the X axis within a record.

Command: WTRC (?) d, i, j

---

### Slice to Trace

Save a time 'slice' from the waterfall to a Trace. A slice is the history of a single X position (data at the marker X position from all stored records). Choose a Trace with the knob and press [Enter] to complete the operation.


The stored slice contains as many points as were stored in the Waterfall buffer (up to Total Count). The points in a slice trace are numbered and displayed from 0 (data from oldest record) to n (data from newest record). Note that this differs from the waterfall display in which the newest record is numbered 0. A recalled slice resemble a time record (X axis is time advancing left to right) with the X axis labeled as waterfall record number (instead of time).

A stored slice trace can be recalled to a Display or Reference Display, used in a User Math Function or saved to disk. A slice can not be copied to the Arbitrary Waveform buffer.

Command: WSLC (?) d, i, j

# Capture Menu

The Capture menu configures the Capture Buffer. See 'Capture Buffer' in Chapter 2 for more information.

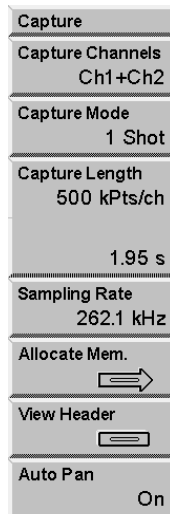

## Capture Channels

Select which inputs to capture [Ch1, Ch2, Ch1+Ch2]. Choose either input or both inputs.

Capturing a single input limits playback to measurements of a single input.

Capturing both inputs allows playback of all measurements but the Capture Length will be half as long as for a single input.

When the measurement group is order. Capture Channels is always set to Ch1+Ch2.

Turn off Auto Offset during capture to avoid interrupting the input signal.

Command: CCHN (?) {i}

## Capture Mode

Select the Capture Mode [1 Shot, Continuous].

Press [Start Capture] to begin storing data in the capture buffer. If the Capture Mode is 1 Shot, capture stops when the buffer is full. Press [Stop Capture] to halt capture before the buffer is full.

If the Capture Mode is Continuous, once capture is started, it continues indefinitely and fills the capture buffer in a circular fashion. In this case, press [Stop Capture] to halt capture with the most recent data in the buffer.

Continuous capture is not allowed in the order measurement group.

## 4-144 Capture Menu

---

Command: CMOD (?) {i}

---

### Capture Length

Set the Capture Length in 2 kPoint (2048 points) increments.

Memory is allocated in blocks. Each block is 2 kPoints. Thus 500 blocks of capture memory allows a maximum Capture Length of 1000 kPoints for a single input and a maximum Capture Length of 500 kPoints for two inputs.

The Capture Length can be less than the allocated memory.

Command: CLEN (?) {i}

---

### Sampling Rate

Select the capture Sampling Rate. The Sampling Rate can be 1, 1/2, 1/4, 1/8, ... times the maximum sampling rate.

The inputs are always digitized at the maximum sampling rate. The maximum sampling rate is 262.1 kHz when the FFT Base Frequency is 102.4 kHz and 256 kHz when the FFT Base Frequency is 100.0 kHz OR the Measurement Group is Octave.

If the selected capture Sampling Rate is the maximum rate, then the digitized input data is stored in the capture buffer. If the Sampling Rate is less than the maximum rate, the input data is filtered and downsampled to the desired rate and then stored in the capture buffer. This allows the capture buffer to hold a longer time period of data at a lower bandwidth. The captured data represents a frequency span from 0 Hz to 1/2.56 times the Sampling Rate. Capturing at less than the maximum sampling rate restricts the playback measurement span to 1/2.56 times the sampling rate.

In the order measurement group, the sampling rate cannot be directly set by the user. Instead the analyzer calculates the sampling rate based on the Max Order and Delta Order. Based on the computed sample rate and the capture length, the analyzer will display the length of the capture buffer in seconds.

Command: CRAT (?) {i}

---

## Allocate Memory

Display the Memory Allocation menu.

Memory Allocation
Total Available 1000 Blks
Capture Memory. 500 Blks
Waterfall Memory 500 Blks
Arb. Memory 0 Blks
Confirm Allocation <input type="button" value="=&gt;"/>
Clear Allocation <input type="button" value="=&gt;"/>
Return <input type="button" value="=&gt;"/>

---

### Total Available

Total Available displays the total memory available for storage, either 980, 2004 or 4054 blocks, depending upon the amount of installed memory. Each block is 2 kPoints (2048 points).

The total of the Capture, Waterfall and Arbitrary memory allocations cannot exceed the Total Available memory. It may be necessary to decrease one allocation in order to increase another.

Command: MMEM ?

---

### Capture Memory

Capture Memory allocates memory blocks for the capture buffer. The allocated Capture Memory sets the limit for the Capture Length.

Command: MALC i, j, k

---

### Waterfall Memory

Waterfall Memory allocates memory blocks for waterfall storage.

Command: MALC i, j, k

---

### Arb Memory

Arb. Memory allocates memory blocks for the Arbitrary Source waveform. The Arbitrary Source can also use the capture buffer as the output waveform.

Command: MALC i, j, k

---

### Confirm Allocation

Confirm Allocation places the new memory allocations into effect. Exiting this menu without pressing <Confirm Allocation> will cancel any adjustments made in this menu.

---

### Clear Allocation

Clear Allocation clears the existing memory allocations in this menu. The new allocations do not take effect unless <Confirm Allocation> is pressed.

---

### View Header

View information about the capture buffer data.

---

### Auto Pan


Set Auto Pan On or Off.

Auto Pan On automatically pans zoomed capture buffer displays to keep up with capture fill and playback. The measurement must be a capture buffer and the display needs to be Zoomed in order for this to have any effect. During capture fill, this keeps the capture display showing the new points being added to the buffer. During playback, this keeps the capture display showing the points which are currently being measured.

Off leaves the display Pan to that set in the [Display Setup] menu.

Command: CPAN (?) {i}

# Analysis Menu

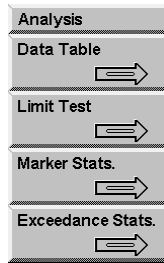

---

## Data Table

Display the Data Table menu for the active display.

The Data Table reports the Y data values for user specified X coordinates. For example, the entries could be a set of harmonic frequencies which need monitoring. Each display has its own Data Table though only the table associated with the active display is shown at any time. If the Display Format is Dual, the active display's Data Table is in the inactive display's location. If the Display Format is Single or Overlay, the Format is first changed to Dual.

To generate a report of the measurement, the active display Data Table may be printed. Choose ASCII Dump for the [Output] <Hard Copy Button> and use [Print Screen].

Data Table x coordinates are saved to disk with the instrument settings.

A Data Table with coordinates defined in the frequency domain will revert to bin locations if the measurement is changed to the time domain (and vice versa). Data Tables are generally only applicable for the domain (time or frequency) and span in which they were originally defined.

---

## Limit Test

Display the Limit Test menu for the active display.

Limit Testing tests the measurement results against a set of defined Limit Segments. When display data exceeds a Limit Segment at any point, the test fails. Each display has its own set of Limit Segments. Limit Segments are saved to disk with the instrument settings.

A Limit Segment is defined as the line between the pair of points (X0,Y0) and (X1,Y1). The segment values between the endpoints are calculated for the displayed span. A segment may be defined as either an Upper or Lower limit. Measurement data which is greater than an Upper limit or less than a Lower limit cause the test to fail.

Limit Segments are defined for the current View. Changing the View invalidates the Limit Segments and limit testing is not available in the new View. Either return to the

## 4-148 Analysis Menu

---

original View or define new segments in the new View. Limit testing is not available for Nyquist or Nichols views.

Measurement data is compared with the limits only over the range of X values for which Limit Segments have been defined. Segments do not have to cover the entire display span or connect with each other. Only segments or the portions of segments within the measurement span are tested. Limits outside the measurement span are ignored.

The test result is shown in the display (to the left of the graph). An audible alarm may be sounded when a test fails.

---

### Marker Stats

Display the Marker Statistics menu.

The Marker Statistics menu displays various statistics about the marker Y values including Mean, Max, Min and Standard Deviation.

When Marker Statistics are On, the various statistical quantities are updated whenever new data is available on the display. The accumulation of the statistical quantities continues even when this menu is not displayed. Return to this menu to see the latest statistical values.

---

### Exceedance Stats

Display the Exceedance Statistics menu.

This menu is used to calculate  $L_n$  (exceedance centile) from Octave or FFT measurements stored in the waterfall buffer. The calculation of  $L_n$  requires that the active display have records stored in the waterfall buffer. The display must be paused or done with a 1-shot waterfall. Waterfall display is not required to be on. The save option must be set to Active Measurement Only to calculate exceedance.

$L_n$  is the amplitude at each bin which is exceeded by n% of the records in the waterfall.  $L_1$  is a large amplitude exceeded only 1% of the time.  $L_{99}$  is a small amplitude exceeded by 99% of the measurements.

The range of records within the waterfall buffer which are analyzed is specified by the Start and Stop index. The Exceedance Pct sets the centile (1-99).

Press <Calculate Excd> to start the calculation. The result is stored in a data trace and has the same measurement type as the waterfall measurements. To view the result, recall the trace to a display or reference graph.

---

### Curve Fit

Display the Curve Fit/Synthesis menu.




Curve fitting is used to identify a mathematical model with closely approximates a measured frequency response. Curve synthesis generates a SR785 frequency response from a mathematical model that can be compared with actual data.

The SR785 does curve fitting and synthesis in two curve tables. Each table contains a model of frequency response function containing up to 20 poles and zeros. The model can be expressed in pole-zero format, pole residue format, or polynomial format. Either table can be synthesized into a trace. When fitting, the best fit model is calculated and placed in curve table 1.

# Data Table Analysis Menu

The Data Table menu is used to edit and display the data table for the active display.

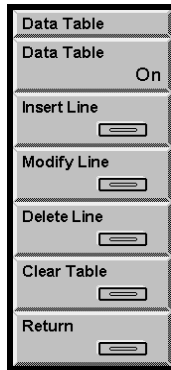

---

### Data Table

Turn the Data Table display for the active display On or Off. The Data Table is shown in the inactive display's location. If the Display Format is Single, the Format is first changed to Dual. Turning off the Data Table does not change the Display Format back.

Command: DTBL (?) d {, i}

---

### Insert Line

Insert a new line in the Data Table for the active display. This key has no effect unless the Data Table is turned On.

Move the marker to the desired frequency (or time) and press this key to insert a new line in the Data Table. The new line will be inserted after the highlighted entry. To select an entry in the table, use the backspace key [ $\leftarrow$ ] or press [Alt] and turn the knob. Press [Alt] again to return the keypad to normal mode.

Command: DINS d, i, x

---

### Modify Line

Modify the x value for an existing line in the Data Table for the active display. This key has no effect unless the Data Table is turned On.

Move the marker to the desired frequency (or time) and press this key to change the x value for the selected line in the Data Table. To select an entry in the table, use the backspace key [ $\leftarrow$ ] or press [Alt] and turn the knob. Press [Alt] again to return the keypad to normal mode.

Command: DMOD d, i, x

### Delete Line

Delete the highlighted line in the Data Table display for the active display. This key has no effect unless the Data Table is turned On. To select an entry in the table, use the backspace key [`<-`] or press [`Alt`] and turn the knob. Press [`Alt`] again to return the keypad to normal mode.

You can not delete all of the entries in the table.

Command: DDLT d, i

---


### Clear Table

Clear the Data Table for the active display. The table is left with a single entry for the first bin in the display.

Command: DCLR d

# Limit Testing Analysis Menu

The Limit Testing menu is used to edit the limits and display the limit test results for the active display.

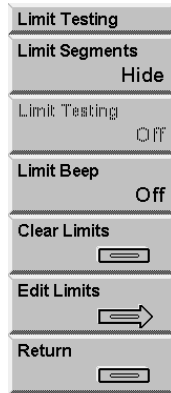

---

## Limit Segments

Display the Limit Segments in the active display [Hide, Show].

Hide turns off the Limit Segments in the display. Limit testing may still be performed.

Show displays the Limit Segments in the display. Choosing <Edit Limits> also sets this to Show.

Command: LSON (?) d {, i}

---

## Limit Testing

Turn Limit Testing for the active display On or Off.

Off turns off limit testing. The limit segments are not affected.

On turns on limit testing. Test results [Pass, Fail] are shown in the display to the left of the graph. Limit testing may not be turned on until there is at least one defined limit segment.

Command: LTST (?) d {, i}

---

## Limit Beep

Turn the audible Limit Alarm for the active display On or Off.

Off turns off the alarm.

On turns on the alarm. Limit tests which fail will sound an audible alarm.

Command: LALM (?) d {, i}

---

### Clear Limits

Clear the Limit Segments for the active display. All limit segments are cleared and Limit Testing is turned Off.

Command: LCLR d

---

### Edit Limits

Display the Edit Limits menu. Press <Return> for the <Limit Test> menu. Choosing this menu also Shows Limit Segments.

The image shows a vertical menu titled "Edit Limits". The menu items are: "New Segment" with a right-pointing arrow icon; "Limit Type" with "Upper" selected; "Segment #", "X0", "Y0", "X1", and "Y1" each with a small square icon; "Delete Segment" with a right-pointing arrow icon; "Shift All" with a small square icon; and "Return" with a right-pointing arrow icon.

---

### New Segment

Add a new segment to the Limit Test for the active display. The new segment starts at the end of the last existing segment. This simplifies the building of a continuous limit table. The new segment is also the Current Segment.

Command: LMAX (?) d {, i}

---

### Limit Type

Select the Limit Type for the Current Segment [Upper, Lower]. Each segment is identified by small triangular arrow markers at the segment end points. These markers are above the segment and point downwards for Upper limits. They are below the segment and point upwards for Lower limits.

## 4-154 Limit Testing Analysis Menu

---

Command: LSEG (?) d, i {, j, x0, y0, x1, y1}

---

### Segment #

Select the Current Segment for editing [0..n]. The limit Type and Start (X0, Y0) and End (X1, Y1) points for this segment are edited in this menu. The Current Segment is identified by two triangular arrows located at the segment endpoints. These arrows are above the Upper or below the Lower Limit Type markers.

Use <Delete Segment> to remove this segment. The numbering of the remaining segments may change as a result. Use the Current Segment markers to identify the correct segment before editing.

---

### X0

Enter or adjust the X coordinate for the Current Segment Start point. The X coordinates are assumed to be in the display units (usually Hz).

Use Zoom and Pan to expand the display when drawing small segments.

Command: LSEG (?) d, i {, j, x0, y0, x1, y1}

---

### Y0

Enter or adjust the Y coordinate for the Current Segment Start point. The Y coordinates are assumed to be in the display units.

Use Zoom and Pan to expand the display when drawing small segments.

Command: LSEG (?) d, i {, j, x0, y0, x1, y1}

---

### X1

Enter or adjust the X coordinate for the Current Segment End point. The X coordinates are assumed to be in the display units (usually Hz).

Use Zoom and Pan to expand the display when drawing small segments.

Command: LSEG (?) d, i {, j, x0, y0, x1, y1}

---

### Y1

Enter or adjust the Y coordinate for the Current Segment End point. The Y coordinates are assumed to be in the display units.

Use Zoom and Pan to expand the display when drawing small segments.

Command: LSEG (?) d, i {, j, x0, y0, x1, y1}

---

### **Delete Segment**

Delete the Current Segment. The numbering of the remaining segments may change as a result. Use the Current Segment markers to identify the correct segment before editing.

Command: LDLT d, i

---


### **Shift All**

Shift all of the Limit Segments up or down together. Enter a value or use the knob to shift by 1/10's of a division.

Command: LSFT d, x

## Marker Statistics Menu

The Marker Statistics menu displays various statistics about the marker Y values.

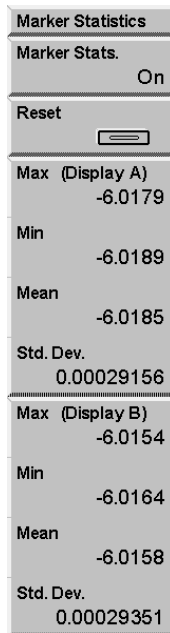

---

### Marker Stats.

Turn the accumulation of marker statistics On or Off.

When Marker Statistics are On, the various statistical quantities are updated whenever new data is available on the display. The accumulation of the statistical quantities continues even when this menu is not displayed. Return to this menu to see the latest statistical values.

When Marker Statistics are Off, the various quantities are frozen. New data is not included in the statistics.

Command: MSAO (?) d {, i}

---

### Reset

Reset and start the accumulation of marker statistics. The various quantities are reset to zero and the accumulation of marker statistics begins.

Use this key whenever the marker position is changed to avoid mixing data from different marker positions.

Command: MSRS



---

### **Max (Display A)**

This menu box displays the maximum value of the marker Y value for display A since Marker Stats was turned On or Reset.

This max value is updated whenever new data is available for Display A and does not update if Display A is not visible.

Command: MSAA ?

---

### **Min (Display A)**

This menu box displays the minimum value of the marker Y value for display A since Marker Stats was turned On or Reset.

This min value is updated whenever new data is available for Display A and does not update if Display A is not visible.

Command: MSIA ?

---

### **Mean (Display A)**

This menu box displays the mean value of the marker Y value for display A since Marker Stats was turned On or Reset.

This mean value is updated whenever new data is available for Display A and does not update if Display A is not visible.

Command: MSEA ?

---

### **Std. Dev. (Display A)**

This menu box displays the standard deviation of the marker Y value for display A since Marker Stats was turned On or Reset.

This standard deviation is updated whenever new data is available for Display A and does not update if Display A is not visible.

Command: MSSA ?

---

### **Max (Display B)**

This menu box displays the maximum value of the marker Y value for display B since Marker Stats was turned On or Reset.

---

## 4-158 Marker Statistics Menu

---

This max value is updated whenever new data is available for Display B and does not update if Display B is not visible.

Command: MSAB ?

---

### Min (Display B)

This menu box displays the minimum value of the marker Y value for display B since Marker Stats was turned On or Reset.

This min value is updated whenever new data is available for Display B and does not update if Display B is not visible.

Command: MSIB ?

---

### Mean (Display B)

This menu box displays the mean value of the marker Y value for display B since Marker Stats was turned On or Reset.

This mean value is updated whenever new data is available for Display B and does not update if Display B is not visible.

Command: MSEB ?

---

### Std. Dev. (Display B)


This menu box displays the standard deviation of the marker Y value for display B since Marker Stats was turned On or Reset.

This standard deviation is updated whenever new data is available for Display B and does not update if Display B is not visible.

Command: MSSB ?

## Exceedance Statistics Menu

The Exceedance Statistics menu calculates exceedance centile statistics from Octave and FFT measurements stored in the waterfall buffer. Exceedance calculation is only available when the waterfall buffer has been saved using the Active Measurement Only Save Option.

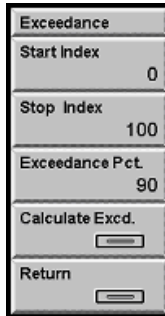

The screenshot shows a vertical menu with the following fields and buttons:

Exceedance	
Start Index	0
Stop Index	100
Exceedance Pct.	90
Calculate Excd.	<input type="button" value="=&gt;"/>
Return	<input type="button" value="=&gt;"/>

---

### Start Index

Set the starting record in the waterfall buffer for the calculation of  $L_n$  (exceedance centile). The analysis starts at the start index and includes all records through the stop index.

Index 0 is the most recent record, larger indexes are older.

The total number of records currently stored and available in the waterfall buffer is displayed in the Vertical Scale Bar. The Start Index should not exceed the total number of records in the buffer.

Both displays use the same Start Index.

Command: ESTR (?) {i}

---

### Stop Index

Set the stop record in the waterfall buffer for the calculation of  $L_n$  (exceedance centile). The analysis starts at the start index and includes all records through the stop index.

Index 0 is the most recent record, larger indexes are older.

The total number of records currently stored and available in the waterfall buffer is displayed in the Vertical Scale Bar. The Stop Index should not exceed the total number of records in the buffer.

Both displays use the same Stop Index.

Command: ESTP (?) {i}

### Exceedance Pct

Set the Exceedance Centile ( $L_n$ ).

$L_n$  is the amplitude at each bin which is exceeded by  $n\%$  of the records in the waterfall.  $L_1$  is a large amplitude exceeded only 1% of the time.  $L_{99}$  is a small amplitude exceeded by 99% of the measurements.

The range of records within the waterfall buffer which are analyzed is specified by the Start and Stop index.

Both displays use the same Exceedance Centile.

Command: EPCT (?) {i}

---

### Calculate Excd

Start the exceedance centile calculation for the active display. Choose a Trace to store the result and press [Enter] to proceed.


The display must be paused or done with a 1-shot waterfall. Waterfall display is not required to be on. The save option must be set to Active Measurement Only.

The result is stored in a data trace and has the same measurement type as the waterfall measurements. To view the result, recall the trace to a display or reference graph.

Command: EXCE d, i

## Curve Fit Menu

The Curve Fit Menu calculates best fit curve parameters from a measured frequency response function. Synthesis of user entered curve parameters is also provided.

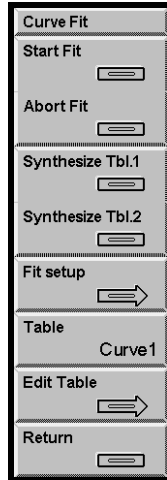

### Start Fit

When start fit is pressed the analyzer begins the curve fitting process. The SR785 will determine a fit to the data in the active display, which must be part of the FFT or swept-sine measurement group. When determining the fit, the analyzer uses the fit region and weighting specified in the fit setup submenu.

After the fit is determined, the fit parameters are placed in curve Table 1. The curve in table 1 is then synthesized into whichever trace is specified by the table. The non-active display is taken off-line, and the trace with the fit data is displayed.

The fit parameters include the frequency scaling specified in curve table 1. For instance, if the frequency scaling is set to 1000, a pole at 1 kHz will show up as a term (s-1) in the pole/zero display.

Command: EFIT d

### Synthesize Table (1 and 2)

Synthesizes the curve in table 1 or 2 into a trace. The frequency range and resolution over which the synthesis is performed is determined by the active display, which must be part of the FFT or swept-sine measurement group.

When the synthesis is complete the data in the active display is replaced by the trace containing the synthesis.

## 4-162 Curve Fit Menu

---

The synthesis uses the frequency scaling specified in the curve table being synthesized. For instance, if the frequency scaling is set to 1000, the term (s-1) in the curve table will produce a pole at 1 kHz in the synthesized frequency response.


Command: ESYN i,d

---

### Table

Turns the curve table display on and off. When the curve table display is on you can use the edit table menu to modify poles, zeros, residues, polynomial coefficients, overall gain, delay, or frequency scaling.

## Fit Setup Submenu

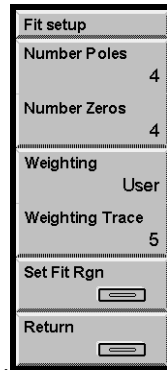


---

### Number Poles

Specifies the number of poles to use in the fit [0...20]. The analyzer will compute the best fit to the data in the active display using the specified number of poles.

Each complex pole is actually a member of a complex conjugate pair and counts as 2 poles.

Command: ENPL (?) {i}

---

### Number Zeros

Specifies the number of zeros to use in the fit [0...20]. The analyzer will compute the best fit to the data in the active display using the specified number of zeros.

Each complex zero is actually a member of a complex conjugate pair and counts as 2 poles.

Command: ENZE (?) {i}

---

### Weighting

Specifies the type of weighting to use when performing the fit [Auto, User].

When Auto weighting is selected, the analyzer computes an appropriate weighting function for the type of data being fit. Live FFT data displays will be weighted with the Coherence function. Live Swept Sine data are weighted with the normalized variance of channel 2. Off line displays, including synthesized traces and recalled data, are weighted uniformly.

When User weighting is selected the analyzer uses data in one of the 5 traces to weight the fit. The trace used is selected with the Weighting Trace softkey. The length of the trace being used must match the length of the display being fit. Only the real part of the trace data is used

## 4-164 Fit Setup Menu

---

Command: EWTU (?) {i}

---

### **Weighting Trace**

Selects which trace [1,5] will be used to weight the fit when User Weighting is selected.

Command: EWTT (?) {i}

---

### **Set Fit Region**


Selects the region of the active display which will be fit. When the Fit Setup menu is entered, a band at the bottom of the graph displays the fit region. The position of the band can be adjusted using the knob. The width of the band can be adjusted by pressing [Alt] and turning the knob.

When the fit parameters are calculated the fit curve will be synthesized only inside the fit region. Outside of the fit region, the data in the fit trace will be zero.

Command: ERNG (?) d, {i,j}



## Edit Table Submenu

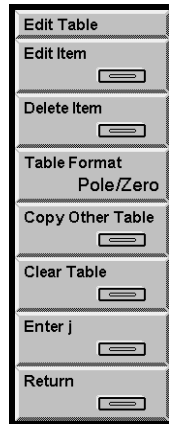


---

### Edit Item

Edits the highlighted item in the curve table. When a curve table is displayed, turning the knob moves the highlighted item. Items available for editing include the curve parameters (poles, zeros, residues, and polynomial coefficients), the frequency scale, delay, gain, or trace number.

Curve parameters may be edited or added using the Edit Item softkey. To edit the item, use the knob to move the highlighted region to the parameter to be edited. Pressing Edit Item brings up the edit window and allows the value to be changed.

To add a new curve parameter, use the knob to move the highlighted region to the blank field at the bottom of the appropriate column in the table. Pressing Edit Items allows entry of a new curve parameter in that field.

Command: EPOL (?) d, i, {f, g} (Poles)  
 Command: EZER (?) d, i, {f, g} (Zeros)  
 Command: ERES (?) d, i, {f, g} (Residues)  
 Command: EPLY (?) d, i, j, {f} (Polynomial Coefficient)  
 Command: EGAN (?) d, {f} (Gain)  
 Command: EFSC (?) d, {f} (Frequency Scale)  
 Command: EDLY (?) d, {f} (Delay)  
 Command: ETRC (?) d, {i} (Trace)

---

### Delete Item

Deletes the highlighted item in the curve table. This option is only available when the highlighted item is a pole, zero, residue, or polynomial coefficient.

---

### Table Format

Changes the representation of the curve in the curve fit table [polynomial, pole zero, or pole residue].

In the polynomial format, the curve is described as the ratio of polynomials in  $s$  multiplied by an overall gain factor..

$$\text{Freq. Re sp.}(s) = \text{Gain} * \frac{(b_n s^n + b_{n-1} s^{n-1} + \dots + b_0)}{a_m s^m + a_{m-1} s^{m-1} + \dots + a_0}$$

In the pole-zero format, the numerator and denominator polynomials are factored so that the frequency response curve is described by the ratio of the products of the poles and zeros: To ensure a real impulse response, all complex poles and zeros only occur in complex conjugate pairs.

$$\text{Freq. Re sp.}(s) = \text{Gain} * \frac{(s - z_n)(s - z_{n-1}) \dots (s - z_0)}{(s - p_m)(s - p_{m-1}) \dots (s - p_0)}$$

In the pole residue format a partial fraction expansion of the pole-zero form is performed to yield the frequency response as a sum of single pole terms weighted by residues.

$$\text{Freq. Re sp.}(s) = \text{Gain} * \left( \frac{R_m}{(s - p_m)} + \frac{R_{m-1}}{(s - p_{m-1})} + \dots + \frac{R_0}{(s - p_0)} \right)$$

Command: EITM? i, j

---

### Copy Other Table

Copies the curve from the other curve fit table into the table currently being edited. Since this option will erase, the curve currently being edited, the analyzer will prompt for a confirmation before copying the table data.

---

### Clear Table

Clears all data from the curve table currently being edited. Since this option will erase the curve currently being edited, the analyzer will prompt for a confirmation before clearing the table data.

Command: ECLR? i

---


### Enter j

Allows the entry of complex poles, zeros, and residues into the curve tables. To enter a complex value, enter the real part, then press the <Enter j> Softkey. A 'j' will appear in the edit window. Now enter the imaginary part.

Complex poles, zeros, or residues only occur as part of complex conjugate pairs. Thus, if the string 3 j 7 is entered for a pole, for instance, it corresponds to the term  $(s - (3 + j7))(s - (3 - j7))$  in the frequency response.

## Disk Menu

The Disk menu saves and recalls data and settings to and from disk. Utility functions are provided to manage directories, erase files and format disks. Files are saved as DOS files on either 720k or 1.4M 3.5 inch disks.

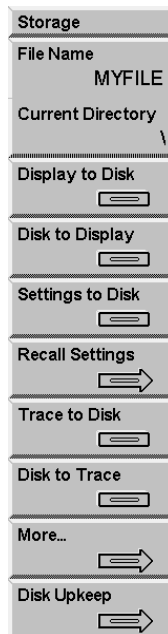


---

### File Name

Enter the save and recall Current File Name.

Turn the knob to bring up the file catalog display listing all files in the Current Directory with the extension .78? (SR785 files). Press [Exp] to display all files in the directory (\*.\*). Choose a file name with the knob and press [Enter] to make it the Current File Name. This file name is used for saving and recalling displays, traces and settings. Remember, saving to this file will write over the existing file!

New file names are entered using the [Alt] key and the letters associated with each key. Press [Alt] again to return to normal keypad operation. DOS file name conventions must be followed. File names are 8 characters or less with an extension (.ext) of up to 3 characters. Default extensions are automatically supplied if no extension is specified. The default extensions should be used since the directory display assumes an extension appropriate for the current menu.

If the filename specified ends in a sequence of digits, the filename will be incremented any time a save-to-disk operation is performed. For instance if SRS001.78D is shown in the menubox, it will be incremented to SRS002.78D after a save operation is performed.

## 4-168 Disk Menu

---

Command: FNAME (?) {s}

Command: FREE ?

Command: FXST ? s

Command: FRST ?

Command: FNXT ?

---

### Current Directory

Enter the Current Disk Directory.

Turning the knob will bring up the directory tree display which lists all of the sub-directories on the disk. Choose a directory with the knob and press [Enter] to make it the Current Directory.

A directory may be entered using the [Alt] key and the letters associated with each key. Press [Alt] again to return to normal keypad operation. An error results if the entered directory does not exist.

New directories are created with <Make Directory>.

Command: FDIR (?) {s}

---

### Display to Disk

Save the active display data to the Current File in the Current Directory. If the Current File has no specified extension, the default extension .78D is used.

The measurement, frequency span, averaging and window are stored with the data in order to preserve the integrity of the data when recalled.

Command: FSAV d

---

### Disk to Display

Recall data from the Current File in the Current Directory to the active display. The active display will be set to Off-Line when data is recalled to it. The display will not update with live measurement results until the Display is set back to Live (in the [Display Options] menu).

The measurement, frequency span, averaging and window are recalled with the data and cannot be changed. These menus are shown in gray and reflect the settings or values for the last Live measurement. Parameters which pertain to the display of the data, such as scale and view, may be changed in the [Display Setup] and [Display Options] menus.

When the Display is returned to Live, the measurement, frequency span, averaging and window all return to the settings in effect before the data was recalled and the live measurement returns to the display.

Command: FRCL d

---

### Settings to Disk

Save the instrument settings to the Current File in the Current Directory. If the Current File has no specified extension, the default extension .78S is used.

The entire instrument setup is stored to the file. Portions of the setup may be excluded from the recall.

Command: SSAV

---

### Recall Settings

Display the Recall Settings menu. Choose the portions of the setup to recall from the Current File in the Current Directory. The recalled settings become effective immediately.

Press [Disk] again to return to the main Disk menu.

Command: SRCL i

---

### Trace to Disk

Save Trace data to the Current File in the Current Directory. If the Current File has no specified extension, the default extension .78D is used.

Select a Trace # (1..5) with the knob and press [Enter] to save it to disk. Only traces which currently have data in them may be selected. To save a display to a Trace, use [Display to Trace] ([Alt] [Print Screen]).

Command: TSAV i

---

### Disk to Trace

Recall data from the Current File in the Current Directory to a Trace. If the Current File has no specified extension, the default extension .78D is used.

Select a Trace # (1..5) with the knob and press [Enter] to recall from disk to the Trace. The recalled data will replace any data presently stored in the Trace.

A stored trace can be recalled to a Display or Reference Display, used in a User Math Function or saved to disk.

Command: TRCL i

## 4-170 Disk Menu

---

### Buffers

Display the Disk Buffers menu. Press <Return> or [Disk] for the main Disk menu. The Disk Buffers menu loads traces from disk files. It also loads and saves Capture, Arbitrary Source and Waterfall data buffers to disk.

### Disk Upkeep

Display the Disk Upkeep menu. Press <Return> or [Disk] for the main Disk menu. The Disk Upkeep menu allows files to be deleted, directories to be created and removed and disks to be formatted.

## Nodal Degree-of-Freedom Menu

The nodal DOF menu appears after <Display to Disk> or <Trace to Disk> is selected if the <Node Info?> softkey in the [System]<Preferences> menu is set to "on." The nodal DOF menu allows entry of descriptive parameters relating to the nodal DOFs which will be saved with the file and available to external programs which convert the files to formats usable by modal analysis programs.

<b>Node Information</b>	
Reference Name	NODE1
Reference Num	1
Ref. Direction	+X
<hr/>	
Response Name	NODE2
Response Num	2
Resp. Direction	+Y
<hr/>	
Abort Save	<input type="button" value="=&gt;"/>
Continue Save	<input type="button" value="=&gt;"/>

---

### Reference Name

Enter a descriptive name for the reference node. For a frequency response function, the reference node is the node where the stimulus is applied. The name may be up to six characters long.

The SR780 does not use the nodal degree-of-freedom information itself. The information is stored with the .78D file and is available for use by external programs which convert the SR780 files into formats suitable for modal analysis programs.

Command: RFNA (?) {s}

---

### Reference Number

Enter the identification for the reference node. For a frequency response function, the reference node is the node where the stimulus is applied. The number may be any integer.

The SR780 does not use the nodal degree-of-freedom information itself. The information is stored with the .78D file and is available for use by external programs which convert the SR780 files into formats suitable for modal analysis programs.

Command: RFNU (?) {i}

### Reference Direction

Enter the direction associated with the reference node, i.e. the direction of the stimulus. Direction can be specified either along the  $\pm X$ ,  $Y$ , or  $Z$  axes for linear stimulus, or as  $\pm \theta_x, \theta_y$ , or  $\theta_z$  direction for angular stimulus.

The SR780 does not use the nodal degree-of-freedom information itself. The information is stored with the .78D file and is available for use by external programs which convert the SR780 files into formats suitable for modal analysis programs.

Command: RFDR (?) {i}

---

### Response Name

Enter a descriptive name for the response node. For a frequency response function, the response node is the node where the response is measured, e.g. the node where the accelerometer or other detector is placed. The name may be up to six characters long.

The SR780 does not use the nodal degree-of-freedom information itself. The information is stored with the .78D file and is available for use by external programs which convert the SR780 files into formats suitable for modal analysis programs.

Command: RSNA (?) {s}

---

### Response Number

Enter the identification for the reference node. For a frequency response function, the response node is the node where the response is measured, e.g. the node where the accelerometer or other detector is placed. The number may be any integer.

The SR780 does not use the nodal degree-of-freedom information itself. The information is stored with the .78D file and is available for use by external programs which convert the SR780 files into formats suitable for modal analysis programs.

Command: RSNU (?) {i}

---

### Response Direction

Enter the direction associated with the response node, i.e. the axis of the accelerometer or other measuring device. Direction can be specified either along the  $\pm X$ ,  $Y$ , or  $Z$  axes for linear stimulus, or as  $\pm \theta_x, \theta_y$ , or  $\theta_z$  for angular measurements.

The SR780 does not use the nodal degree-of-freedom information itself. The information is stored with the .78D file and is available for use by external programs which convert the SR780 files into formats suitable for modal analysis programs.

Command: RSDR (?) {i}

---



**Abort Save**


Press to abort the selected disk operation while in the nodal DOF information menu.

---

**Continue Save**

Press when finished entering nodal DOF information to continue the selected disk operation.

## Recall Settings Menu

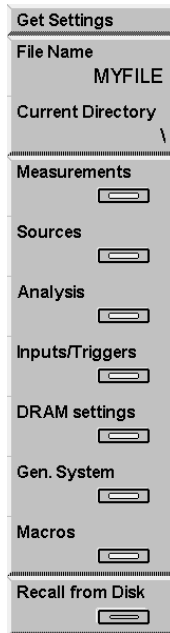

---

### File Name

Enter the save and recall Current File Name.

Turn the knob to bring up the file catalog display listing all files in the Current Directory with the extension .78? (SR785 files). Press [Exp] to display all files in the directory (\*.\*). Choose a file name with the knob and press [Enter] to make it the Current File Name. This file name is used for saving and recalling displays, traces and settings. Remember, saving to this file will write over the existing file!

New file names are entered using the [Alt] key and the letters associated with each key. Press [Alt] again to return to normal keypad operation. DOS file name conventions must be followed. File names are 8 characters or less with an extension (.ext) of up to 3 characters. Default extensions are automatically supplied if no extension is specified. The default extensions should be used since the directory display assumes an extension appropriate for the current menu.

Command: FNAM (?) {s}  
Command: FREE ?  
Command: FXST ? s  
Command: FRST ?  
Command: FNXT ?

---

### Current Directory

Enter the Current Disk Directory.

Turning the knob will bring up the directory tree display which lists all of the sub-directories on the disk. Choose a directory with the knob and press [Enter] to make it the Current Directory.

A directory may be entered using the [Alt] key and the letters associated with each key. Press [Alt] again to return to normal keypad operation. An error results if the entered directory does not exist.

New directories are created with <Make Directory>.

Command: FDIR (?) {s}

---

### Measurements

Include the [Frequency], [Display Setup], [Display Options], [Marker], [Average], [Window], [User Math], [Waterfall] and [Capture] menu settings in the recall. This may invalidate current [Analysis] menu settings.

Command: SRCL i

---

### Sources

Include the [Source] menu settings in the recall.

Command: SRCL i

---

### Analysis

Include the [Analysis] menu settings in the recall. Do not recall Data Tables and Limits unless Measurements are also recalled. This is because the analysis parameters are defined for specific frequency Spans and Views.

Command: SRCL i

---

### Inputs/Triggers

Include the [Input] and [Trigger] menu settings in the recall.

Command: SRCL i

---

### DRAM settings

Include the memory allocation (Capture, Arbitrary Source, or Waterfall) settings in the recall.

Command: SRCL i

## 4-176 Recall Settings Menu

---

### General System

Include the [System] and [Output] menu settings in the recall.

Command: SRCL i

---

### Macros

Include stored macros in the recall.

Command: SRCL i

---


### Recall from Disk

Recall the selected instrument settings from the Current File in the Current Directory. The recalled settings become effective immediately.

Command: SRCL i

---

# Disk Buffers Menu

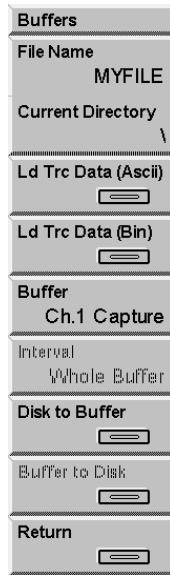


---

## File Name

Enter the save and recall Current File Name.

Turn the knob to bring up the file catalog display listing all files in the Current Directory with the extension .78? (SR785 files). Press [Exp] to display all files in the directory (\*.\*). Choose a file name with the knob and press [Enter] to make it the Current File Name. This file name is used for saving and recalling displays, traces and settings. Remember, saving to this file will write over the existing file!

New file names are entered using the [Alt] key and the letters associated with each key. Press [Alt] again to return to normal keypad operation. DOS file name conventions must be followed. File names are 8 characters or less with an extension (.ext) of up to 3 characters. Default extensions are automatically supplied if no extension is specified. The default extensions should be used since the directory display assumes an extension appropriate for the current menu.

Command: FNAM (?) {s}  
 Command: FREE ?  
 Command: FXST ? s  
 Command: FRST ?  
 Command: FNXT ?

---

## Current Directory

Enter the Current Disk Directory.

## 4-178 Disk Buffers Menu

---

Turning the knob will bring up the directory tree display which lists all of the sub-directories on the disk. Choose a directory with the knob and press [Enter] to make it the Current Directory.

A directory may be entered using the [Alt] key and the letters associated with each key. Press [Alt] again to return to normal keypad operation. An error results if the entered directory does not exist.

New directories are created with <Make Directory>.

Command: FDIR (?) {s}

---

### Load Trace Data (ASCII)

Load ASCII data from the Current File in the Current Directory into an existing Trace. This is a way to import calibration data from a file into a trace. The trace can then be used in a user function to calibrate live measurements.

The target trace must already store data of the type and length to be loaded from disk. The file only contains data points, no measurement information is recalled. The recalled data simply takes the place of the stored trace data and assumes the trace's measurement parameters (Measurement, Averaging, Window, Frequency span) and length.

The ASCII file consists of N pairs of ASCII floating point values. Each pair of values represents the real and imaginary parts of a single complex data point. Thus, the file contains 2N floating point numbers. The values are linearly scaled (usually with units of Vpk).

The first value in the file is N, followed by the real part of the first point, then the imaginary part of the first point, etc. Each value is delimited by a comma, white space, tab or carriage return. Traces which are entirely real should store 0.0 as the imaginary part of every point.

The value of N depends upon the type and length of the target trace. N should be picked at least as large as the number of points displayed in the trace. If N is specified too small, the trace will be zero padded. If N is too large, the extra points will be ignored.

For example, the ASCII file for a 400 line FFT might be

```
400
0.000, 0.500
1.000, 1.500
2.000, 2.500
...
399.000, 399.500
```

where (0.000, 0.500) is the first complex data point and (399.000, 399.500) is the 400th data point.

To read this file into a trace, a 400 line FFT must first be stored in the trace. This can be done using <Display to Trace> or <Disk to Trace>. Loading the ASCII file into the trace then replaces the trace data with the data from the disk file.

Command: TASC ? i, n

---

### Load Trace Data (binary)

Load binary data from the Current File in the Current Directory into a Trace. This is a way to import calibration data from a file into a trace. The trace can then be used in a user function to calibrate live measurements.

The target trace must already store data of the type and length to be loaded from disk. The file only contains data points, no measurement information is recalled. The recalled data simply takes the place of the stored trace data and assumes the trace's measurement parameters (Measurement, Averaging, Window, Frequency span) and length.

The binary file contains N complex points. Each complex data point is two floating point numbers, representing the real and imaginary parts. Thus, the file contains 2N floating point numbers. The values are linearly scaled (usually with units of Vpk).

The first value in the file is N (4-byte long integer), followed by the real part of the first point (4-byte IEEE float), then the imaginary part of the first point (4 byte IEEE float), etc. There are no delimiters in a binary file. The data is packed. Traces which are entirely real should store the value 0.0 as the imaginary part of every point.

The binary file has the following form where [?] represents a single byte.

```
[N lsb] [N] [N] [N msb] [X0 lsb] [X0] [X0] [X0 msb] [Y0 lsb] [Y0] [Y0] [Y0 msb]
[X1 lsb] [X1] [X1] [X1 msb] [Y1 lsb] ... [Y(N-1) msb]
```

where (X0, Y0) is the first complex data point and (XN-1, YN-1) is the Nth data point. Each value is 4-byte IEEE float.

The value of N depends upon the type and length of the target trace. N should be picked at least as large as the number of points displayed in the trace. If N is specified too small, the trace will be zero padded. If N is too large, the extra points will be ignored.

Command: TLOD ? i, n

---

### Buffer

Choose which data buffer (Ch1 capture, Ch2 capture, Arbitrary Waveform or Waterfall) to save to or load from disk. Buffers which have no allocation may not be chosen. Only the capture buffers specified by <Capture Channels> may be chosen.

Ch1 capture and Ch2 capture are saved and recalled separately. Both can be recalled into a buffer configured for Ch1+Ch2capture channels.

## 4-180 Disk Buffers Menu

---

The waterfall buffers for both displays are saved together (all stored waterfall records for both displays are saved). If one display has waterfall storage off, then no data is saved for that display.

Capture and Arbitrary files are compatible. Data saved from a capture buffer can be loaded into the arbitrary waveform buffer (and vice versa).

---

### Interval

Choose to save the entire data buffer or only the playback portion (Capture or Arbitrary Waveform). All stored waterfall records for both displays are saved.

This only affects the saving of data to disk. When a buffer is loaded from disk, the buffer is configured to be exactly the length of the recalled data.

---

### Disk to Buffer

Load binary data from the Current File in the Current Directory into the selected data buffer. This function only recalls files made using <Buffer to Disk>. The selected buffer (Capture, Arbitrary Waveform or Waterfall) must have sufficient memory allocated before loading from disk.

#### Load Arbitrary Waveform

A saved Arbitrary or Capture file may be loaded into the Arbitrary Waveform buffer. The Arb. Src. Length will be set to the length of the disk data.

#### Load A Capture Buffer

A saved Arbitrary or Capture file may be loaded into a Capture buffer.

If the disk data has the same Capture Length and Sampling Rate as shown in the current [Capture] menu, the data is simply loaded into the selected channel's buffer. In the case of Ch1+Ch2 capture, loading one channel does not disturb the contents of the other channel's buffer.

If either the length or sampling rate is different from the current [Capture] menu, they are changed to the length and sampling rate of the disk data. In the case of Ch1+Ch2 capture, loading Ch1 (Ch2) capture will zero the Ch2 (Ch1) capture buffer if the capture length or sampling rate is changed by this operation.

To recall both capture channels, configure <Capture Channels> for Ch1+Ch2. Both files must contain the same length buffers with the same sampling rate. Use <Buffer> to choose either Ch1 or Ch2 and load the appropriate disk file to the buffer. This will set the capture length and sampling rate according to the recalled disk file. Then use <Buffer> to choose the other channel and load the other disk file to the buffer. Since the existing Capture Length and Sampling Rate are already compatible, this operation does not zero the first capture buffer.

#### Load the Waterfall Buffer

A saved Waterfall file may be loaded into the Waterfall buffer.



Loading a Waterfall file loads data into the waterfall buffer *and* recalls the measurement parameters which were in effect when the waterfall buffer was saved. The recalled parameters include items in the [Frequency], [Display Setup], [Display Options], [Average], [Window], [User Math], and [Waterfall] menus.

If the file contains records from both displays, data is recalled for both displays. The Waterfall Storage is set to Recalled. Recalled simply means that the waterfall displays contain data which did not come from the two inputs but rather from disk.

Command: APUT ? i  
Command: CPUT ? i, j  
Command: WPUT ?  
Command: TPUT ? i

---

**Buffer to Disk**

Save binary data from the selected data buffer to the Current File in the Current Directory. Use <Interval> to select the playback portion or entire buffer.

Command: AGET ? i  
Command: CGET ? i  
Command: WGET ?  
Command: TGET ? i

## Disk Upkeep Menu



---

### File Name

Enter the File Name for Disk Upkeep. Keeping this file name separate from the save/recall file name makes it harder to erase a file accidentally.

Turn the knob to bring up the file catalog display listing all files in the Current Directory with the extension .78? (SR785 files). Press [Exp] to display all files in the directory (\*.\*). Choose a file name with the knob and press [Enter] to make it the Upkeep File Name. This file may be erased with <Del File>.

New file names are entered using the [Alt] key and the letters associated with each key. Press [Alt] again to return to normal keypad operation. DOS file name conventions must be followed. File names are 8 characters or less with an extension (.ext) of up to 3 characters.

Command: DNAM (?) {s}

---

### Current Directory

Enter the Current Disk Directory.

Turning the knob will bring up the directory tree display which lists all of the sub-directories on the disk. Choose a directory with the knob and press [Enter] to make it the Current Directory.

A directory may be entered using the [Alt] key and the letters associated with each key. Press [Alt] again to return to normal keypad operation. An error results if the entered directory does not exist.

New directories are created with <Make Directory>.

Command: FDIR (?) {s}

---

### **Make Directory**

Make a new directory on the disk. Enter a directory name with the [Alt] key. The new directory will be created in the Current Directory.

Command: MDIR s

---

### **Del File**

Delete the Disk Upkeep File from the Current Directory. Make sure that the file and directory are correct before pressing this key!

Command: DELF

---

### **Del Cur Dir**

Delete the Current Directory from disk. Make sure that the directory is correct before pressing this key! The directory must be empty before it can be deleted.

Command: DELD

---

### **Format Floppy**

Format the disk. The disk capacity is 720k for DS/DD disks and 1.44M for DS/HD disks. A disk must be formatted before files can be stored on it. Formatting a disk involves erasing all information from the disk and rewriting the directory. **Formatting a disk destroys all data presently on disk!**

---

# Output Menu

The Output menu prints or plots the screen display. Printer and plotter configurations are set in this menu.



---

## Hard Copy Button

Assign either Bitmap/Print, Vector/Plotter ASCII Dump to the [Print Screen] hardkey.

The [Print Screen] hardkey allows you to print, plot or dump the screen from any menu.

Command: POUT

### Bitmap/Print

Prints the screen using the selected Bitmap/Printer and Bitmap Area to the selected Destination (Interface or Disk).

Other front panel operations are disabled until printing is completed. Pressing backspace will abort the operation.

Command: PRNT

### Vector/Plot

Plot the screen using the selected Vector Plotter to the selected Destination (Interface or Disk).

Other front panel operations are disabled until plotting is completed. Pressing backspace will abort the operation.

Command: PLOT

### ASCII Dump

Dump the data in the active display in ASCII format to the selected Destination (Interface or Disk). The data is written f,y,x[cr] with the data for each point on a single line. f is the frequency or time value, y is the y-axis value and x is the x-axis value (if necessary).

Other front panel operations are disabled until the dump is completed. Pressing backspace will abort the operation.

Command: DUMP

---

### Bitmap/Printer

Select the Bitmap/Printer type for the Bitmap/Print operation [Epson FX, HP PCL, HP Small PCL, PCX 2 bit, GIF, PCX 8 bit].

Epson FX is used for Epson compatible dot matrix graphics printers.

HP PCL and HP Small PCL are used for HP LaserJet or InkJet compatible printers. HP PCL is a full page printout and HP Small PCL is half page.

The PCX and GIF types are bit mapped file formats. Choose one of these if you want to print to a disk file and import the image into a PC application.

The Destination for all Bitmap/Printer types can be either an on-line Interface or a Disk file.

Command: PRTP (?) {i}

---

### Bitmap Area

Select the portion of the screen for the Bitmap/Print operation [Graphs, Menu, Indicators, All]. Vector/Plot only plots the displayed graphs.

Graphs prints the graph areas. This is a single or dual graph.

Menu prints the menu of softkeys.

Indicators prints the status indicators at the top of the screen.

All prints the entire screen image including the menu softkeys and status area.

Command: PSCR (?) {i}

---

### Vector/Plotter

Select the Vector/Plotter type for the Vector/Plot operation [HPGL, Postscript]. Vector/Plot only plots the displayed graphs.

## 4-186 Output Menu

---

HPGL is used for HPGL compatible plotters.

Postscript is used for Postscript compatible plotters or printers. Postscript output to a disk file is a convenient way to import the screen plot to a PC application.

The Destination for all Plotter Types can be either an on-line Interface or a Disk file.

Command: PLTP (?) {i}

---

### Destination

This softkey selects the Print, Plot or ASCII Dump Destination [Disk File, Centronics, Serial, GPIB].

Choose Disk File to output to a disk file. All Bitmap/Printer and Vector/Plotter types can save to Disk File. This is convenient for using a printer or plotter which is not connected to this unit or to import the screen image into a PC application. Files are written to the Current Directory specified in the [Disk] menu. Files are named SCRNX.XXX.EXT where XXXX is a 4 digit number which automatically increments starting at the File Start Number. The .EXT extensions are

.TXT	ASCII Dump
.EFX	Epson
.PCL	HP Laser
.GIF	GIF
.PCX	PCX
.HPG	HPGL
.EPS	Postscript

Centronics selects the parallel printer interface. Most dot matrix printers use this.

Serial selects the RS-232 serial interface. Some plotters and printers use this.

GPIB selects the IEEE-488 interface. Most HPGL plotters use this.

Command: PDST (?) {i}

---

### GPIB Control

Choose the GPIB Controller [Host or SR785].

Choose SR785 if there are no other controllers attached to the GPIB interface. This is typically the case where only a plotter is attached to the SR785.

Choose Host if there is host computer which is in control. This is the case where both the SR785 and the plotter are controlled by a host computer. The host is responsible for issuing the PLOT command to the SR785, and then making the plotter a listener and the SR785 a talker. The plotter commands will then be transferred from the SR785 to the

plotter. The host can periodically untalk the SR785 and serial poll the SR785 (via the Instrument Status) to determine when the plot is finished.

Command: PCIC (?) {i}

---

### GPIB Address

Enter the GPIB Address for a GPIB plotter or printer [0..30].

If an HPGL plotter is used with the GPIB interface, this GPIB Address must agree with the actual address of the plotter.

Command: PLTA (?) {i}

---

### Edit A Note

Select a Note to edit [0..7] and press [Enter] to display the Edit Note menu. A Note is a text annotation within a display graph. Notes are visible only within this menu and may be printed or plotted with the graph using [Print Screen].

Press <Return> for the [Output] menu.

Command: NOTE i, j {, k, l, m, s}

Edit A Note	
Text X position	50
Text Y position	50
Display	0
Visible	Yes
Text String	ABC
Return	<input type="button" value="Return"/>

---

### Text X Position

Adjust the horizontal position of the Note with the knob.

Command: NOTE i, j {, k, l, m, s}

---

### Text Y Position

Adjust the vertical position of the Note with the knob.

Command: NOTE i, j {, k, l, m, s}

## 4-188 Output Menu

---

---

### Display

Select which display the Note appears in [0=DisplayA, 1=DisplayB].

Command: NOTE i, j {, k, l, m, s}

---

### Visible

Make the Note visible [Yes (visible), No (invisible)].

Command: NOTE i, j {, k, l, m, s}

---

### Text String

Enter the Note text. Use the [Alt] key to enter the note text.

Command: NOTE i, j {, k, l, m, s}

---

## File Start Number

Set the Starting File Number for printing, plotting or dumping to Disk [XXXX (up to 4 digit integer)].

All Bitmap/Printer and Vector/Plotter types can save to Disk. This is convenient for using a printer which is not connected to this unit or to import the screen image into a PC application. Files are written to the Current Directory specified in the [Disk] menu. Files are named SCRNXXXX.EXT where XXXX is a 4 digit number which automatically increments starting at the File Start Number.


Command: PFIL (?) {i}



---

## Colors

Display the Hardcopy Colors menu. Press <Return> or [Output] for the main [Output] menu.

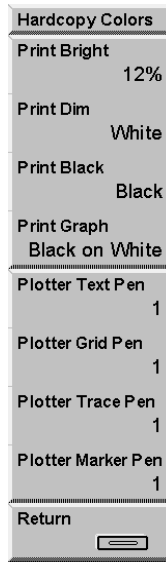

### Print Color Selections

Printing is a bit mapped operation. Different areas of the screen can be printed with different print densities giving a range of gray between black and white (white paper). These colors only affect Epson (FX and 24 pin) and HP (Laser and Laser small) printer types.

---

#### Print Bright

Select the print density of highlighted areas (e.g. active display marker display) [White (none on white paper), 6%, 12%, 25%, 50%, (100%) Black].

Command: PBRI (?) {i}

---

#### Print Dim

Select the print density of the normal background (e.g. X and Y axis graph labels) [White (none on white paper), 6%, 12%, 25%, 50%, (100%) Black].

Command: PDIM (?) {i}

---

#### Print Black

Select the print density for black text [Black, White (none on white paper)]. Use White only if the Print Bright and Print Dim are both black or gray.

Command: PBLK (?) {i}

---

### **Print Graph**

Select the graph print mode [Black on White, White on Black]. The graph is normally printed as Black (data and grid) on White (background). This uses the least printer ink. Choosing White (data and grid) on Black (background) prints a large black area on the paper.

Command: PGRF (?) {i}

### **Plot Color Selections**

Many plotters have a multipen carousel. In this case, different items on the graph may be plotted using different color pens. When using a single pen plotter, all items are plotted using the one pen regardless of the pen definitions. Only the graph display areas are plotted.

---

### **Plotter Text Pen**

Assign a plotter pen number to the text labels on the graphs [1 to 8].

Command: PLTX (?) {i}

---

### **Plotter Grid Pen**

Assign a plotter pen number to the graph grid [1 to 8].

Command: PLGD (?) {i}

---

### **Plotter Trace Pen**

Assign a plotter pen number to the graph data trace [1 to 8].

Command: PLTR (?) {i}

---


### **Plotter Marker Pen**

Assign a plotter pen number to the graph marker [1 to 8]. Turn the marker off to avoid plotting it.

Command: PLMK (?) {i}

# System Menu

The System menu configures the interfaces, clock/calendar and system preferences and edits recorded macros. The Diagnostics menu accesses various hardware tests.

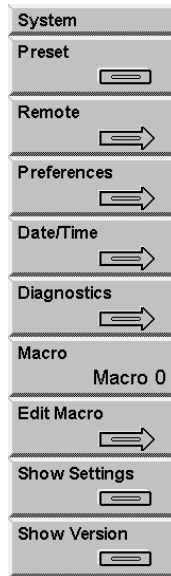

---

## Preset

Reset the instrument to the default settings. Only the remote interface settings are not changed. All stored data are lost.

This function requires pressing the [Enter] key to confirm.

To completely reset the instrument (including remote interface settings), power the unit on with the backspace key held down.

Command: \*RST

---

## Remote

Display the [System] <Remote> interface menu. Remote interface parameters should not be altered while the computer interface is active. Press <Return> or [System] for the [System] menu.

---

## Preferences

Display the [System] <Preferences> menu. Press <Return> or [System] for the [System] menu.

## 4-192 System Menu

---

### Date/Time

Display the [System] <Date/Time> menu. Press <Return> or [System] for the [System] menu.

---

### Diagnostics

Display the [System] <Diagnostics> test menu. Press <Return> or [System] for the [System] menu.

---

### Macro

Use the knob to pick one of the ten Macros to edit and press [Enter]. All macros may be edited including those which have not been recorded.

See 'Macros' in Chapter 3 for more.

---

### Edit Macro

Display the Edit Macro menu to edit the selected macro. Press <Cancel> to exit back to the main [System] menu.

It is generally easier to record a macro than it is to enter it in this menu. Use this menu to edit existing macros. See 'Macros' in Chapter 3 for more.

To enter a macro string, enter Menu, Function, Control and Number keys by simply pressing the desired key (these keys do not have their normal effect in this menu). Use <Special Keys> and the knob to select <Softkeys> and [Enter]. [-] (Backspace) deletes the term before the cursor.

To edit a macro string, use <Macro String> and the knob to move the flashing cursor within the string. <Delete> will delete the highlighted term. Menu, Function, Control and Number keys and <Special Keys> will either insert or replace at the cursor. Use <Insert/Replace> to switch between insert and replace mode. <Clear Macro> will delete the entire string.

To increase (CW) or decrease (CCW) a knob event, use <Inc. Knob Count> and <Dec. Knob Count>.

<Cancel> will discard any changes made in this menu and exit this menu.

<Enter Macro> will enter the new macro and exit this menu. Pressing [Enter] while the cursor is in the macro string display also enters the macro and exits this menu.

---

**Show Settings**

Show the instrument setup. This key enters the Help system and displays the measurement setup and system settings. Press [0] to exit Help.

---

**Show Version**

Show the power on screen. This screen shows the version number as well as the results of the power on tests. The size of the installed data memory is also displayed (2 Ms, 4 Ms or 8Ms).

## System Remote Menu

Remote
Output To GPIB
GPIB Address 10
Override REM? Yes
Baud Rate 9600 bd
Word Length 8 bits
Parity None
View Qs <input type="checkbox"/>
Return <input type="checkbox"/>

---

### Output To

Select the Output Interface [RS232, GPIB]. All interface query responses are directed to the selected Output Interface regardless of which interface received the query. Use the OUTX command to select the correct interface at the beginning of every program.

Command: OUTX (?) {i}

---

### GPIB Address

Enter or adjust the instrument GPIB Address [0..30].

---

### Override REM?

Select the Override Remote mode [Yes, No]. In general, every GPIB interface command will put the unit into the REMOTE state with the front panel locked out. To defeat this feature, set this entry to Yes. In this mode, the front panel is not locked out when the unit is in the REMOTE state.

To return from REMOTE to LOCAL (front panel enabled), press [Help/Local].

Command: OVRM (?) {i}

---

### Baud Rate

Select the RS232 (Serial) interface Baud Rate [300, 1200, 2400, 4800, 9600, 19.2k]. Most PC's use 9600 as a default.

### Word Length

Select the RS232 (Serial) interface Word Length [7 bits, 8 bits]. 8 bits is standard.

---

### Parity

Select the RS232 (Serial) interface Parity [None, Even, Odd]. None is standard.

---

### View Queues


Display the Remote Interface Queue. The SR785 buffers the characters which are received or transmitted over the interface. The Queue display shows the interface history. The Queue display may slow down communications and should be displayed only when testing or debugging a host program.

The upper half of the screen is the Receive Queue. These are the most recent characters which have been received by the SR785. Commands which have already been executed are shown in normal text. Commands which have not yet been executed are shown with a bright background. Command errors are shown in inverse text.

The lower half of the screen is the Transmit Queue. These are the most recent characters which have been placed in output buffer. Characters which have already been sent over the interface are shown in normal text. Characters which are waiting to be sent are shown with a bright background.

Press any key to exit from this display.

## System Preferences Menu

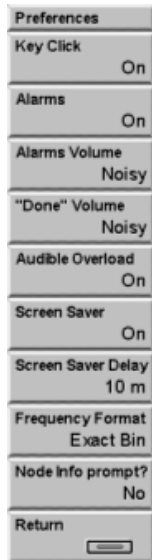

---

### Key Click

Turn the Key Click on or off.

Command: KCLK (?) {i}

---

### Alarms

Turn the system Alarm messages On or Off.

Command: ALRM (?) {i}

---

### Alarms Volume

If Alarms Volume is set to Quiet, alarm messages have no accompanying sound. If Alarms are set to Noisy, alarm messages are accompanied by an audible alarm.

Set system Alarms to Off for no alarm message at all.

Command: ALRT (?) {i}

---

### Done Volume

If Done Volume is set to Noisy, an audible alarm is sounded when a measurement is done or completed. For example, when linear averaging is complete, an alarm is sounded.

Command: ADON (?) {i}



**Audible Overload**

If Audible Overload is set to On, an audible alarm is sounded whenever an input is overloaded.

Command: AOVL (?) {i}

---

**Screen Saver**

Turns the Screen Saver On or Off. Using the Screen Saver extends the life of the video display and prevents screen 'burn'.

When the Screen Saver is On and the Screen Saver Delay passes without a key press, knob turn or interface command, the screen will be blanked and random patterns are displayed. Press any key to return the screen to normal.

Command: \SAVR (?) {i}

---

**Screen Saver Delay**

Set the Screen Saver Delay in minutes [1..59].

When the Screen Saver is on and the Screen Saver Delay passes without a key press, knob turn or interface command, the screen will be blanked and random patterns displayed. Press any key to return the screen to normal.

Command: SDLY (?) {i}

---

**Frequency Format**

Select the Frequency Format [Exact Bin, Rounded].

Exact Bin displays frequencies exactly. This can make the frequency display have many more digits than required to distinguish adjacent frequency bins. However, this format is useful when determining the exact frequency of a bin.

Rounded shortens the frequency displays to show only those digits which are required to distinguish adjacent frequency bins. This simplifies and shortens the frequency displays but does not show the exact frequency of a bin.

For example, for a 400 line FFT on 100.0 kHz base with a span of 97.66 Hz, the first frequency bin is shown as 244.140625 mHz (Exact Bin) or 244.1 mHz (Rounded). The second bin is 488.28125 mHz (Exact Bin) and 488.3 mHz (Rounded).

Command: FFMT (?) {i}

---

### Node Info Prompt


Turns the node information prompt when saving files on or off.

Modal analysis programs require that a given spectrum, cross spectrum, or frequency response function (transfer function) be associated with information regarding the nodal degrees of freedom (DOF). When Node Info Prompt is set on, the SR785 allows entry of DOF parameters when saving data using Trace to Disk or Display to Disk. When the <Trace to Disk> or <Display to Disk> is initiated using a remote command, the node info prompt will not be displayed, however the information entered in the [Node Info] menu will still be saved with the disk file as long as the <Node Info Prompt?> softkey is on.

The nodal DOF information is saved with the disk file and is used when translating SR785 files to Universal File Format (UFF) and other file formats used by popular modal analysis programs.

Command: SVNI (?) {i}

## System Date/Time Menu

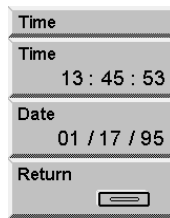

---

### Time

Set the System Time [hh:mm:ss]. The System Time is entered as hours:minutes:seconds and all entries are 6 digits. 24 hour format is used. [Enter] sets the new time or press this softkey again to abort the entry.

Command: TIME (?) {i}


---

### Date

Set the System Date [mm:dd:yy]. The System Date is entered as month:day:year and all entries are 6 digits. [Enter] sets the new date or press this softkey again to abort the entry.

Command: DATE (?) {i}

## System Diagnostics Menu

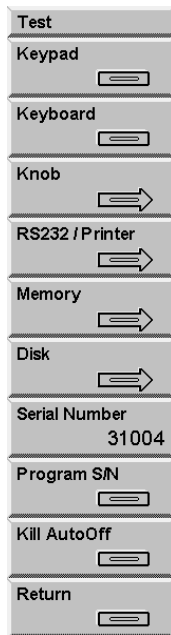

---

### Keypad Test

Display the Keypad Test screen. The test screen is a map of the keypad. Press each key until all keys are marked. When all of the keys are marked, the test is completed successfully. Turn the knob to abort this test and return to the <Diagnostics> menu.

---

### Keyboard Test

Display the Keyboard Test screen. Characters typed on an attached PC keyboard (PC/XT or 8088 mode required) are displayed on the test screen. If the displayed characters are accurate, then the keyboard interface is functioning correctly. If not, check that the keyboard is in the correct mode. Many keyboards have a switch on the bottom to select PC/XT (8088) or AT (80286) mode. Turn the knob to abort this test and return to the <Diagnostics> menu.

---

### Knob Test

Display the Knob Test screen. A circle with a marker is displayed. Select one of the 4 speeds displayed in the menu [Snail, Human Baby, Adult, Lunatic]. Turning the knob will move the marker around the circle, verifying knob action and direction. Press <Return> for the <Diagnostics> menu.

---

### RS232/Printer Test

Display the RS232/Printer Test screen.

A loop back adapter is required to complete the RS232 test. This adapter is simply a mating connector with pins 2 and 3 connected so characters transmitted by the interface will be received as well. Press <Begin> to start the test. Press <Return> for the <Diagnostics> menu.

---

### **String to RS232**

Send an ASCII test string to a serial printer connected to the RS232 port.

---

### **String to Printer**

Send an ASCII test string to a parallel printer connected to the Printer port.

---

## **Memory Test**

Display the <Memory Test> menu. Select a test and press <Begin>.

---

### **System RAM**

Display the System RAM Test screen. Press <Begin> to test the main CPU program RAM. Press <Return> for the <Memory Test> menu.

---

### **System ROM**

Display the System ROM Test screen. Press <Begin> to test the main CPU program ROM. Press <Return> for the <Memory Test> menu.

---

### **Video RAM**

Display the Video Ram Test screen. Press <Begin> to test the video display memory. Press <Return> for the <Memory Test> menu.

---

### **Help ROM**

Display the Help ROM Test screen. Press <Begin> to test the help system memory. Press <Return> for the <Memory Test> menu.

---

## **Disk Test**

Display the Disk Drive Test screen.

**Continuing with this test will destroy any data on the disk currently in the drive!**

Remove the data disk and insert a scratch disk. This test will check the controller, format the disk and perform a read/write check on the disk. The entire test takes about 2 minutes. Press <Begin> to start the test. When the test reaches the Disk Changed Sensor Out phase, remove the disk and then insert it again.

## 4-202 System Diagnostics Menu

---

Press <Return> for the <Diagnostics> menu.

---

### Serial Number

Reset the internally stored serial number after replacing the DSP board. The serial number which is reported on the power-on screen and in the GPIB identification string is stored in the calibration ROM on the DSP board.

**Do not use this function unless the DSP board has been replaced!**

After replacing the DSP board, enter the unit serial number (shown on the rear panel) and use <Program S/N> to store the serial number permanently.

Press <Return> for the <Diagnostics> menu.

---

### Program S/N

Reset the internally stored serial number after replacing the DSP board. The serial number which is reported on the power-on screen and in the GPIB identification string is stored in the calibration ROM on the DSP board.

**Do not use this function unless the DSP board has been replaced!**

After replacing the DSP board, use <Serial Number> to enter the unit serial number (shown on the rear panel) and use <Program S/N> to store the serial number permanently.

Press <Return> for the <Diagnostics> menu.

---


### Kill AutoOff

Clear the internal input offset calibrations and disable Auto Offset.

**Do not use this function unless you are making input offset or CMR adjustments to an analog input board. See the Service Manual for more information.**

The unit must be turned off and back on after using this function to restore the input calibrations!

# Edit Macro Menu

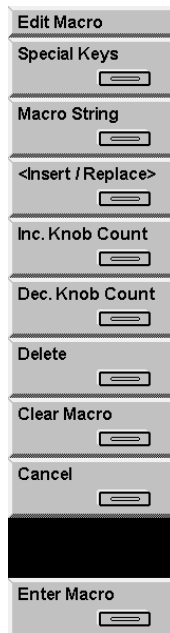


---

## Special Keys

Use the knob to pick a <Softkey> or other special key from the display and press [Enter] to place it in the macro string at the cursor location.

To place [Enter] into the macro string, use <Special Keys>. Pressing the [Enter] key while the cursor is in the macro string displays the macro and exits this menu.

To place [-] (backspace) into the macro string, use <Special Keys>. The [-] key deletes the term before the cursor.

---

## Macro String

Move the cursor to the macro string window at the top of the screen. Use the knob to move the highlighted cursor within the string. <Delete> will delete the highlighted term. Menu, Function, Control and Number keys and <Special Keys> will either insert or replace at the cursor. Use <Insert/Replace> to switch between insert and replace mode.

[-] (Backspace) deletes the term before the cursor.

---

## Insert/Replace

Toggle between insert and replace mode while editing a macro string.

If editing in insert mode, 'Ins' appears in the upper right corner of the edit window. If editing in replace mode, 'Rep' appears.

## 4-204 Edit Macro Menu

---

When editing a macro string, Menu, Function, Control and Number keys and <Special Keys> will insert or replace before the cursor.

---

### **Inc. Knob Count**

To edit a knob event, move the cursor to the macro string window at the top of the screen and move to the desired knob event (displayed as <Knob n>). <Inc. Knob Count> will increase the knob count by 1 which is the same as rotating clockwise one choice in a parameter list.

---

### **Dec. Knob Count**

To edit a knob event, move the cursor to the macro string window at the top of the screen and move to the desired knob event (displayed as <Knob n>). <Dec. Knob Count> will decrease the knob count by 1 which is the same as rotating counterclockwise one choice in a parameter list.

---

### **Delete**

Delete the highlighted term in a macro string. To move the cursor in a macro string, use <Macro String> and the knob.

---

### **Clear Macro**

Clear the entire macro string.

---

### **Cancel**

Discard any changes made in this menu and exit this menu.

---

### **Enter Macro**

Enter the displayed string as the new macro and exit this menu. Use [Play Macro] to playback the macro string.

---



## Chapter 5

# Remote Programming

## In This Chapter

<b>Index of Commands</b>	<b>5-2</b>	Arbitrary Source Commands	5-57
<b>Alphabetical List of Commands</b>	<b>5-12</b>	Swept Sine Source Commands	5-58
<b>Introduction</b>	<b>5-21</b>	Input Commands	5-60
Communicating With GPIB	5-21	Transducer Parameter Commands	5-62
Communicating With RS232	5-21	Tachometer Input Commands	5-64
Screen Indicators And Queues	5-21	Input Playback Commands	5-65
Command Format	5-22	Trigger Commands	5-66
Command Synchronization	5-22	Average Commands (FFT, Correlation, and Order)	5-68
Data Synchronization	5-23	Average Commands (Octave)	5-71
Command Context Errors	5-24	Average Commands (Swept Sine)	5-73
Example Program	5-24	User Math Commands	5-74
<b>Command Syntax</b>	<b>5-25</b>	Window Commands	5-78
Commands and Units	5-25	Waterfall Commands	5-80
Help	5-26	Capture Commands	5-83
Things to Remember!	5-27	Memory Allocation Commands	5-84
Output Interface (RS232 or GPIB)	5-27	Data Table Commands	5-85
Active Display and Linking	5-27	Limit Test Commands	5-87
Analyzer Configuration and Linking	5-27	Marker Statistics Commands	5-89
Command Context Errors	5-27	Exceedance Statistics Commands	5-90
Frequency Commands (FFT and Correlation)	5-28	Curve Fit Commands	5-91
Frequency Commands (Octave)	5-30	Disk Commands	5-94
Frequency Commands (Swept Sine)	5-32	Output Commands	5-97
Frequency Commands (Order Analysis)	5-34	System Commands	5-100
Frequency Commands (Time/Histogram)	5-36	Front Panel Commands	5-102
Display Setup Commands	5-37	Data Transfer Commands	5-107
Display Options Commands	5-43	Interface Commands	5-117
Marker Commands	5-45	Nodal Degree-of-Freedom Commands	5-118
Marker Commands (Normal)	5-47	Status Reporting Commands	5-120
Marker Commands (Harmonic)	5-49	<b>Status Word Definitions</b>	<b>5-124</b>
Marker Commands (Sideband)	5-50	<b>Example Program</b>	<b>5-129</b>
Marker Commands (Band)	5-51		
Marker Commands (Frequency/ Damping)	5-52		
Source Commands	5-53		
Sine Source Commands	5-54		
Chirp Source Commands	5-55		
Noise Source Commands	5-56		

# Index of Commands

## Variables

d	display (0=DisplayA, 1=DisplayB, 2=Both)
i, j, k, l, m	integers
x, y	real numbers
f	frequency in Hz
s	text string

---

## Frequency (FFT and Correlation Measurement Group)

FSPN (?) d {, f}	5-28	FFT Frequency Span
FLIN (?) d {, i}	5-28	FFT Resolution
FBAS (?) d {, i}	5-28	FFT Base Frequency
FSTR (?) d {, f}	5-28	FFT Start Frequency
FCTR (?) d {, f}	5-28	FFT Center Frequency
FEND (?) d {, f}	5-29	FFT End Frequency
UNST d	5-29	Unsettle Measurement

---

## Frequency (Octave Measurement Group)

OHIB (?) d {, f}	5-30	Octave Highest Band
OLOB (?) d {, f}	5-30	Octave Lowest Band
ORES (?) d {, i}	5-30	Octave Resolution
OCHN (?) d {, i}	5-30	Octave Channels
UNST d	5-30	Unsettle Measurement

---

## Frequency (Swept Sine Measurement Group)

SSTR (?) d {, f}	5-32	Swept Sine Start Frequency
SSTP (?) d {, f}	5-32	Swept Sine Stop Frequency
SSFR ?	5-32	Swept Sine Progress?
SRPT (?) d {, i}	5-32	Swept Sine Repeat Mode
SSTY (?) d {, i}	5-32	Swept Sine Sweep Type
SARS (?) d {, i}	5-32	Swept Sine Auto Resolution
SNPS (?) d {, i}	5-33	Swept Sine Number of Points
SSKP (?) d {, i}	5-33	Swept Sine Max Step Size
SFST (?) d {, x}	5-33	Swept Sine Faster Threshold
SSLO (?) d {, x}	5-33	Swept Sine Slower Threshold

---

## Frequency (Order Measurement Group)

ORMX (?) d {, f}	5-34	Max RPM
ORMN (?) d {, f}	5-34	Min RPM
OMAX (?) d {, f}	5-34	Max Order
ODLT (?) d {, f}	5-34	Delta Order
OTRK (?) d {, i}	5-34	Order Tracking
ONPT (?) d {, i}	5-34	Track Points
OSTO (?) d {, i}	5-35	Track Storage Mode
O1TK (?) d {, f}	5-35	Track 1 Order
O1BN (?) d {, i}	5-35	Track 1 BNC
O2TK (?) d {, f}	5-35	Track 2 Order
O2BN (?) d {, i}	5-35	Track 2 BNC

---

## Frequency (Time/Histogram Measurement Group)

FSPN (?) d {, i}	5-28	Sampling Time
------------------	------	---------------

---

---

HLEN (?) d {, f}	5-36	Histogram Length
HBIN (?) d {, i}	5-36	Histogram Bins
FBAS (?) d {, i}	5-28	Base Frequency
HRPT (?) d {, i}	5-36	Histogram Repeat

---

### Display Setup

MGRP (?) d {, i}	5-37	Measurement Group
MEAS (?) d {, i}	5-37	Display Measurement
VIEW (?) d {, i}	5-39	Display View
UNIT? d	5-40	Display Units
UNDB (?) d {, i}	5-40	dB Units
UNPK (?) d {, i}	5-40	Peak Units
PSDU (?) d {, i}	5-40	PSD Units On/Off
UNPH (?) d {, i}	5-40	Phase Units
DBMR (?) {x}	5-40	dBm Reference Impedance
YMAX (?) d {, x}	5-40	Y Maximum
YMID (?) d {, x}	5-40	Y Midpoint
YMIN (?) d {, x}	5-41	Y Minimum
YDIV (?) d {, x}	5-41	Y/Division
YMDX d x,y	5-41	Ymax + Y/Division
XCEN (?) d {, x}	5-41	X Center (Polar)
XDIV (?) d {, x}	5-41	X/Division (Polar)
YCEN (?) d {, x}	5-41	Y Center (Polar)
Y2DV (?) d {, x}	5-41	Y/Division (Polar)
XPAN (?) d {, i}	5-41	X Pan
XZOM (?) d {, i}	5-42	X Zoom

---

### Display Options

DISP (?) d {, i}	5-43	Display Update Mode
DFMT (?) {i}	5-43	Display Format
ACTD (?) i	5-43	Active Display
XAXS (?) d {, i}	5-43	X Axis Scale Type
RPMF (?) {i}	5-43	RPM Frequency
GRID (?) d {, i}	5-43	Grid On/Off
GDIV (?) d {, i}	5-43	Number of Grid Divisions
GPOL (?) d {, i}	5-43	Nyquist Grid Type
PHSL (?) d {, x}	5-44	Phase Suppress Threshold
DDXW (?) d {, x}	5-44	d/dx Window

---

### Marker

MRKR (?) d {, i}	5-45	Marker Tracking
MKMD (?) d {, i}	5-45	Marker Mode
MBIN d, i	5-45	Marker Move to bin i
MWFL d, i	5-45	Marker Move Waterfall to record i
MRKX ? d	5-45	Query the Marker X Position
MRKB ? d	5-45	Query the Marker Bin
MRKY ? d	5-45	Query the Marker Y Position
MRKZ ? d	5-46	Marker Z Read
MKMX d	5-46	Move the Marker to the Maximum
MKMN d	5-46	Move the Marker to the Minimum
MKCN d	5-46	Center of FFT Span to Marker

## 5-4 Index of Commands

---

### Normal Marker

MSEK (?) d {, i}	5-47	Normal Marker Seeks Mode
MWID (?) d {, i}	5-47	Normal Marker Width
MREL (?) d {, i}	5-47	Normal Marker Relative Mode
MROX (?) d {, x}	5-47	Normal Marker X Offset
MROY (?) d {, x}	5-47	Normal Marker Y Offset
MXRL (?) d {, i}	5-47	Normal Marker X Rel Mode
MRON (?) d	5-47	Set Marker to Ref

---

### Harmonic Marker

HRMN (?) d {, i}	5-49	Number of Harmonics
HDSP (?) d {, i}	5-49	Harmonic/Sideband Display
HRDO (?) d {, i}	5-49	Harmonic/Sideband Readout
HTHD ? d, i	5-49	Query the Total Harmonic Distortion
HPWR ? d	5-49	Query the Harmonic Power

---

### Sideband Marker

SSEP (?) d {, f}	5-50	Sideband Separation
SDBN (?) d {, i}	5-50	Number of Sidebands
HDSP (?) d {, i}	5-50	Harmonic/Sideband Display
HRDO (?) d {, i}	5-50	Harmonic/Sideband Readout
SPWR ? d, i	5-50	Query the Sideband Power

---

### Band Marker

BMKL (?) d, i {, j}	5-51	Band i Left Bin j
BMKR (?) d, i {, j}	5-51	Band i Right Bin j
BEXC ? d {, i}	5-51	Band Exclude
BRAM ? d {, i}	5-51	Band Ratio Mode
BPWR ? d, i	5-51	Band Power
BRAT ? d, i	5-51	Band Ratio

---

### Frequency Damping Marker

MMCA ? d	5-52	Frequency/Damping Results
----------	------	---------------------------

---

### Source

SRCO (?) {i}	5-53	Source On/Off
STYP (?) {i}	5-53	Source Type

---

### Sine Source

S1FR (?) {f}	5-54	Sine Frequency 1
S1AM (?) {x}	5-54	Sine Amplitude 1
S2FR (?) {f}	5-54	Sine Frequency 2
S2AM (?) {x}	5-54	Sine Amplitude 2
SOFF (?) {x}	5-54	Sine Offset

---

### Chirp Source

CAMP (?) {x}	5-55	Chirp Amplitude
CBUR (?) {x}	5-55	Chirp Burst Percentage
CSRC (?) {i}	5-55	Source Display

---

### Noise Source

NAMP (?) {x}	5-56	Noise Amplitude
NTYP (?) {i}	5-56	Noise Type

---

---

NBUR (?) {x}	5-56	Noise Burst Percentage
NPER (?) {x}	5-56	Noise Source Period
CSRC (?) {i}	5-56	Source Display

---

**Arbitrary Source**

AAMP (?) {i}	5-57	Arbitrary Source Amplitude
ARAT (?) {i}	5-57	Arbitrary Source Rate
ASRC (?) {i}	5-57	Arbitrary Source Buffer
ASTR (?) {i}	5-57	Arbitrary Source Start
ALEN (?) {i}	5-57	Arbitrary Source Length
TARB i	5-57	Trace to Arb
ALOD ? n	5-110	Download Arbitrary Binary
AGET ? i	5-112	Upload Arbitrary Buffer
APUT ? i	5-113	Download Arbitrary Buffer

---

**Swept Sine Source (Swept Sine Measurement Group)**

SSAL (?) {i}	5-58	Swept Sine Auto Level
SSAM (?) {x}	5-58	Swept Sine Amplitude
SSRF (?) {x}	5-58	Swept Sine Ideal Reference
SRMP (?) {i}	5-58	Swept Sine Source Ramping Off/On
SRAT (?) {x}	5-58	Swept Sine Ramp Rate
SSUL (?) {x}	5-58	Swept Sine Upper Limit
SSLL (?) {x}	5-59	Swept Sine Lower Limit
SMAX (?) {x}	5-59	Swept Sine Max Level
SOFF (?) {x}	5-59	Swept Sine Offset

---

**Inputs**

ISRC (?) {i}	5-60	Input Source
LINK (?) {i}	5-60	Analyzer Configuration
I1MD (?) {i}	5-60	Ch1 Input Mode
I1GD (?) {i}	5-60	Ch1 Input Grounding
I1CP (?) {i}	5-60	Ch1 Input Coupling
I1RG (?) {x}	5-60	Ch1 Input Range
A1RG (?) {i}	5-60	Ch1 AutoRange Off/On
I1AR (?) {i}	5-60	Ch1 AutoRange Mode
I1AF (?) {i}	5-60	Ch1 Anti-Alias Filter
I1AW (?) {i}	5-61	Ch1 A-Weight Filter
I2MD (?) {i}	5-61	Ch2 Input Mode
I2GD (?) {i}	5-61	Ch2 Input Grounding
I2CP (?) {i}	5-61	Ch2 Input Coupling
I2RG (?) {x}	5-61	Ch2 Input Range
A2RG (?) {i}	5-61	Ch2 AutoRange Off/On
I2AR (?) {i}	5-61	Ch2 AutoRange Mode
I2AF (?) {i}	5-61	Ch2 Anti-Alias Filter
I2AW (?) {i}	5-61	Ch2 A-Weight Filter
IAOM (?) {i}	5-61	Auto Offset
INPC ? i	5-122	Input Ovld Read

---

**Transducer Parameters**

EU1M (?) d {,i}	5-62	Ch1 Engineering Units Off/On
EU1L (?) d {,i}	5-62	Ch1 EU Label
EU1V (?) d {,x}	5-62	Ch1 EU/Volt
EU1U (?) d {,s}	5-62	Ch1 User Label
TD1C (?) d {,i}	5-62	Ch1 Transducer Convert

---

## 5-6 Index of Commands

---

EU2M (?) {i}	5-62	Ch2 Engineering Units Off/On
EU2L (?) {i}	5-63	Ch2 EU Label
EU2V (?) {x}	5-63	Ch2 EU/Volt
EU2U (?) {s}	5-63	Ch2 User Label
TD2C (?) d {,i}	5-63	Ch2 Transducer Convert

---

### Input Playback

ISTR (?) {i}	5-65	Capture Playback Start
ILEN (?) {i}	5-65	Capture Playback Length
IMOD (?) {i}	5-65	Capture Playback Mode
ISPD (?) {i}	5-65	Capture Playback Speed

---

### Trigger

TMOD (?) {i}	5-66	Trigger Arming Mode
TSRC (?) {i}	5-66	Trigger Source
TLVL (?) {x}	5-66	Trigger Level
TSLP (?) {i}	5-66	Trigger Slope
TDLA (?) {x}	5-66	Trigger Delay A
TDLB (?) {x}	5-66	Trigger Delay B
TARM	5-66	Manual Trigger Arm
STMD (?) {i}	5-66	Triggered Source Mode
TMAN	5-66	Manual Trigger
TRSR (?) {x}	5-67	Starting RPM
TRSM (?) {i}	5-67	Start RPM Mode
TRDR (?) {x}	5-67	Delta RPM
TRDM (?) {i}	5-67	Delta RPM Mode
TIAS (?) {x}	5-67	Time Arm Increment

---

### Average (FFT, Correlation, and Order)

FAVG (?) d {, i}	5-68	Average On
FAVM (?) d {, i}	5-68	Display Average
FAVT (?) d {, i}	5-68	Average Type
FAVN (?) d {, i}	5-68	Average Number
NAVG ? d	5-68	Averages Completed?
FOVL (?) d {, x}	5-68	FFT Time Record Increment
FREJ (?) d {, i}	5-69	FFT Overload Reject
TAVM (?) { i}	5-69	Trigger Average Mode
PAVO (?) d {, i}	5-69	Average Preview
PAVT (?) d {, x}	5-69	Preview Time
PAVA	5-69	Accept Preview
PAVR	5-69	Reject Preview

---

**Average (Octave Measurement Group)**

OTYP (?) d {, i}	5-71	Octave Average Type
OTIM (?) d {, x}	5-71	Octave Average Time
NAVG ? d	5-71	Octave Averages Completed?
OCNF (?) d {, i}	5-71	Octave Confidence Level
OIMP (?) d {, i}	5-71	Octave Power Bin
OLAT (?) d {, i}	5-72	Octave Linear Average Mode

**Average (Swept Sine Measurement Group)**

SSTM (?) d {, x}	5-73	Swept Sine Settle Time
SSCY (?) d {, i}	5-73	Swept Sine Settle Cycles
SITM (?) d {, x}	5-73	Swept Sine Int. Time
SICY (?) d {, i}	5-73	Swept Sine Int. Cycles

**User Math Functions**

USRF (?) i {, j, k, l ...}	5-74	FFT User Function i
USRR (?) i {, j, k, l ...}	5-74	Correlation User Function i
USRO (?) i {, j, k, l ...}	5-74	Octave User Function i
USRS (?) i {, j, k, l ...}	5-74	Swept Sine User Function i
USRT (?) i {, j, k, l ...}	5-74	Order User Function i
USRH (?) i {, j, k, l ...}	5-74	Time/Histogram User Function i
USRC (?) i {, x, y}	5-77	User Constant i

**Window**

FWIN (?) d {, i}	5-78	Window
FWFL (?) d {, i}	5-78	Force Length
FWTC (?) d {, i}	5-79	Expo TC
TRWI i	5-79	Trace i to User Window
WTR d, i	5-79	Window to Trace i
WSYM (?) d {, i}	5-79	Window Form

**Waterfall**

WDSP (?) d {, i}	5-80	Waterfall Display
WSTO (?) d {, i}	5-80	Waterfall Storage
WFSB (?) d {, i}	5-80	Waterfall Save Option
WAVA ? d	5-80	Waterfall Records Stored?
WTOT (?) d {, i}	5-80	Waterfall Total Count
WFSK (?) d {, i}	5-80	Waterfall FFT Skip
WOSK (?) d {, x}	5-80	Waterfall Octave Skip
WVCT (?) d {, i}	5-81	Waterfall View Count
WHIT (?) d {, i}	5-81	Waterfall Trace Height
WANG (?) d {, i}	5-81	Waterfall Angle
WFST (?) d {, i}	5-81	Waterfall Fast Angles
WTHR (?) d {, i}	5-81	Waterfall Threshold
WHID (?) d {, i}	5-81	Waterfall Hidden Lines
WREV (?) d {, i}	5-81	Waterfall Paused Drawing
WTRC d, i, j	5-81	Waterfall Record to Trace
WSLC d, i, j	5-82	Waterfall Slice to Trace
WGET ?	5-114	Upload Waterfall Buffer
WPUT ?	5-115	Download Waterfall Buffer

**Capture**

CCHN (?) {i}	5-83	Capture Channels
CMOD (?) {i}	5-83	Capture Mode

## 5-8 Index of Commands

---

CLEN (?) {i}	5-83	Capture Length
CRAT (?) {i}	5-83	Capture Rate
CPAN (?) {i}	5-83	Capture Auto Pan
CSTR	5-83	Capture Start
CSTP	5-83	Capture Stop
CGET ? i, j	5-113	Upload Capture Buffer
CPUT ? i, j	5-32	Download Capture Buffer

---

### Memory

MMEM ?	5-84	Total Memory Available
MALC i, j, k	5-84	Allocate Memory

---

### Data Table

DTBL (?) d {, i}	5-85	Data Table On
DMAX ? d {, i}	5-85	Data Table Length
DINS d, i, j	5-85	Data Table Insert bin j at line i
DDLT d, i	5-85	Data Table Delete
DCLR d	5-86	Data Table Clear
DTRD ? d {, i}	5-86	Data Table Query

---

### Limit Testing

LSO (?) d {, i}	5-87	Show Limit Segments
LTST (?) d {, i}	5-87	Limit Testing On
LALM (?) d {, i}	5-87	Limit Beep
LFAL ? d	5-87	Limit Test Fail?
LCLR d	5-87	Clear Limits
LMAX (?) d {, i}	5-87	Last Limit Segment #
LSEG (?) d, i {, j, x0, y0, x1, y1}	5-87	Set Limit Segment i
LDLT d, i	5-88	Delete Limit Segment
LSFT d, x	5-88	Shift Limit Segments

---

### Marker Statistics

MSAO (?) d {, i}	5-89	Marker Statistics On
MSRS	5-89	Reset Marker Statistics
MSAA ?	5-89	Marker Max A Read
MSIA ?	5-89	Marker Min A Read
MSEA ?	5-89	Marker Mean A Read
MSSA ?	5-89	Marker Std Dev A Read
MSAB ?	5-89	Marker Max B Read
MSIB ?	5-89	Marker Min B Read
MSEB ?	5-89	Marker Mean B Read
MSSB ?	5-89	Marker Std Dev B Read

---

### Exceedance Statistics

ESTR (?) {i}	5-90	Exceed Start Index
ESTP (?) {i}	5-90	Exceed Stop Index
EPCT (?) {i}	5-90	Exceed Centile
EXCE d, i	5-90	Calculate Exceedance

---

### Curve Fit

EFIT d	5-91	Start Fit
ESYN i, d	5-91	Synthesis
ENPL (?) {, i}	5-91	Number Poles
ENZE (?) {, i}	5-91	Number Zeros

---



---

EWTU (?) {,i}	5-91	Weighting
EWTT (?) {,i}	5-91	Weighting Trace
ERNG (?) d {,i,j}	5-91	Fit Range
EDLY (?) i {,x}	5-91	Table Delay
EFSC (?) i {,x}	5-91	Table Frequency Scale
EFSC (?) i {,x}	5-91	Table Frequency Scale
EGAN (?) i {,x}	5-92	Table Gain
ETRC (?) i {,j}	5-92	Table Trace
ETRC (?) i j	5-92	Table Item Query
EPLY (?) i,j,k{x}	5-92	Polynomials
EPOL (?) i,j{x,y}	5-92	Poles
EZER (?) i,j{x,y}	5-92	Zeros
ERES (?) i,j{x,y}	5-92	Residues

---

**Disk**

FNAM (?) {s}	5-94	Save/Recall File Name
FDIR (?) {s}	5-94	Save/Recall Directory
MDIR s	5-94	Make Directory
FXST ? s	5-94	File Exist?
FREE ?	5-94	Disk Free Space?
FRST ?	5-95	Reset Disk Catalog
FNXT ?	5-95	Read Disk Catalog
FSAV d	5-94	Display to Disk
FRCL d	5-95	Disk to Display
SSAV	5-95	Settings to Disk
SRCL i	5-95	Recall Settings
TSAV i	5-96	Trace i to Disk
TRCL i	5-96	Disk to Trace i
DNAM (?) {s}	5-96	Upkeep File Name
DELF	5-96	Delete File
DELD	5-96	Delete Directory

---

**Output**

POUT (?) {i}	5-97	Print Screen hardkey
PDST (?) {i}	5-97	Print/Plot Destination
PFIL (?) {i}	5-97	File Start Number
DUMP	5-97	Dump Display Data
PRNT	5-97	Print Screen
PRTP (?) {i}	5-97	Bitmap/Printer Type
PSCR (?) {i}	5-97	Bitmap Area
PBRI (?) {i}	5-98	Print Bright
PDIM (?) {i}	5-98	Print Dim
PBLK (?) {i}	5-98	Print Black
PGRF (?) {i}	5-98	Print Graph
PLOT	5-98	Plot Screen
PLTP (?) {i}	5-98	Vector/Plotter Type
PLTA (?) {i}	5-98	Plotter GPIB Address
PCIC (?) {i}	5-98	GPIB Control
PLTX (?) {i}	5-98	Plotter Text Pen
PLGD (?) {i}	5-98	Plotter Grid Pen
PLTR (?) {i}	5-98	Plotter Trace Pen
PLMK (?) {i}	5-98	Plotter Marker Pen
NOTE i, j {, k, l, m, s}	5-99	Display Note

---

## 5-10 Index of Commands

---

### System

OUTX (?) {i}	5-100	Output Remote Interface
OVRM (?) {i}	5-100	Override Remote
KCLK (?) {i}	5-100	Key Click
ALRM (?) {i}	5-100	Alarms On
ALRT (?) {i}	5-100	Alarms Volume
ADON (?) {i}	5-100	Done Volume
AOVL (?) {i}	5-100	Audible Overload
SAVR (?) {i}	5-100	Screen Saver On
SDLY (?) {i}	5-100	Screen Saver Delay
FFMT (?) {i}	5-100	Frequency Format
TIME (?) {i, j, k}	5-100	Time
DATE (?) {i, j, k}	5-101	Date

---

### Front Panel

ACTD (?) i	5-102	Active Display
STRT	5-102	Start/Reset
PAUS	5-102	Pause
CONT	5-102	Continue
UNST d	5-102	Unsettle Measurement
CSTR	5-83	Capture Start
CSTP	5-83	Capture Stop
SVTR d, i	5-102	Display d to Trace i
RCTR d, i	5-103	Trace i to Display d
SVRF d, i	5-103	Reference to Trace
RCRF d, i	5-103	Trace to Reference
DUMP	5-103	Dump Display Data
PRNT	5-103	Print Screen
PLOT	5-103	Plot Screen
ASCL d	5-103	AutoScale
A1RG (?) {i}	5-103	Ch1 AutoRange Off/On
A2RG (?) {i}	5-103	Ch2 AutoRange Off/On
MRON (?) d	5-103	Set Marker to Ref
MKMX d	5-104	Move the Marker to the Maximum
MKMN d	5-104	Move the Marker to the Minimum
MKCN d	5-104	Center of FFT Span to Marker
DREF d, i	5-104	Display Ref
SNAP d	5-104	Snap Reference
SBRI (?) {i}	5-100	Screen Brightness
SCON (?) {i}	5-104	Screen Contrast
KEYP i	5-104	Key Press
KNOB i	5-106	Knob
TONE i, j	5-106	Tone
PLAY i	5-106	Play Tune

---

### Traces

SVTR d, i	5-102	Display d to Trace i
RCTR d, i	5-103	Trace i to Display d
TSAV i	5-96	Trace i to Disk
TRCL i	5-96	Disk to Trace i
TLOD ? i, n	5-109	Download Trace i Binary
TASC ? i, n	5-109	Download Trace i Ascii
TGET ? i	5-111	Upload Trace i Buffer
TPUT ? i	5-112	Download Trace i Buffer

---

**Data Transfer**


---

DSPN ? d	5-107	Display d Length
DSPY ? d {, j}	5-107	Read Display d (bin j) ASCII
REFY ? d, j	5-107	Read Ref Display d bin j
DSPB ? d {, j}	5-107	Read Display d Binary
DBIN ? d, j	5-108	Read Display d Bin Freq or Time
FBIN ? d, x	5-109	Read Display d Bin Number
TLOD ? i, n	5-109	Download Trace i Binary
TASC ? i, n	5-109	Download Trace i Ascii
ALOD ? n	5-110	Download Arbitrary Binary
TGET ? i	5-111	Upload Trace i Buffer
TPUT ? i	5-112	Download Trace i Buffer
AGET ? i	5-112	Upload Arbitrary Buffer
APUT ? i	5-113	Download Arbitrary Buffer
CGET ? i, j	5-113	Upload Capture Buffer
CPUT ? i, j	5-32	Download Capture Buffer
WGET ?	5-114	Upload Waterfall Buffer
WPUT ?	5-115	Download Waterfall Buffer

---

**Interface**

*RST	5-117	Reset
*IDN ?	5-117	Device Identification?
LOCL (?) {i}	5-117	Local/Remote
OVRM (?) {i}	5-117	Override Remote

---

**Status**

*CLS	5-120	Clear All Status Registers
*PSC (?) {i}	5-120	Power On Status Clear
*SRE (?) {i} {, j}	5-120	Serial Poll Status Enable
*STB ? {i}	5-120	Serial Poll Status Read
*ESE (?) {i} {, j}	5-120	Standard Event Status Enable
*ESR ? {i}	5-121	Standard Event Status Read
ERRE (?) {i} {, j}	5-121	Error Status Enable
ERRS ? {i}	5-121	Error Status Read
INSE (?) {i} {, j}	5-121	Instrument Status Enable
INST ? {i}	5-122	Instrument Status Read
DSPE (?) {i} {, j}	5-122	Display Status Enable
DSPS ? {i}	5-122	Display Status Read
INPE (?) {i} {, j}	5-122	Input Status Enable
INPS ? {i}	5-122	Input Status Read
INPC ? i	5-122	Input Ovlid Read

---

**Nodal Degree-of-Freedom Information**

SVNI(?) {i}	5-118	Save Nodal Information
RFNA(?) {s}	5-118	Reference Node Name
RFNU(?) {i}	5-118	Reference Node Number
RFDR(?) {i}	5-118	Reference Node Direction
RSNA(?) {s}	5-119	Response Node Name
RSNU(?) {i}	5-119	Response Node Number
RSDR(?) {i}	5-118	Response Node Direction

---

## Alphabetical List of Commands

---

*		
*CLS	5-120	Clear All Status Registers
*ESE (?) {i} {,j}	5-120	Standard Event Status Enable
*ESR ? {i}	5-121	Standard Event Status Read
*IDN ?	5-117	Device Identification?
*PSC (?) {i}	5-120	Power On Status Clear
*RST	5-117	Reset
*SRE (?) {i} {,j}	5-120	Serial Poll Status Enable
*STB ? {i}	5-120	Serial Poll Status Read

---

### A

A1RG (?) {i}	5-60	Ch1 AutoRange Off/On
A2RG (?) {i}	5-61	Ch2 AutoRange Off/On
AAMP (?) {i}	5-57	Arbitrary Source Amplitude
ACTD (?) i	5-102	Active Display
ADON (?) {i}	5-100	Done Volume
AGET ? i	5-112	Upload Arbitrary Buffer
ALEN (?) {i}	5-57	Arbitrary Source Length
ALOD ? n	5-110	Download Arbitrary Binary
ALRM (?) {i}	5-100	Alarms On
ALRT (?) {i}	5-100	Alarms Volume
AOVL (?) {i}	5-100	Audible Overload
APUT ? i	5-113	Download Arbitrary Buffer
ARAT (?) {i}	5-57	Arbitrary Source Rate
ASCL d	5-103	AutoScale
ASRC (?) {i}	5-57	Arbitrary Source Buffer
ASTR (?) {i}	5-57	Arbitrary Source Start

---

### B

BEXC ? d {, i}	5-51	Band Exclude
BMKL (?) d, i {, j}	5-51	Band i Left Bin j
BMKR (?) d, i {, j}	5-51	Band i Right Bin j
BPWR ? d, i	5-51	Band Power
BRAM ? d {, i}	5-51	Band Ratio Mode
BRAT ? d, i	5-51	Band Ratio

---

### C

CAMP (?) {x}	5-55	Chirp Amplitude
CBUR (?) {x}	5-55	Chirp Burst Percentage
CCHN (?) {i}	5-83	Capture Channels
CGET ? i, j	5-113	Upload Capture Buffer
CLEN (?) {i}	5-83	Capture Length
CMOD (?) {i}	5-83	Capture Mode
CONT	5-102	Continue
CPAN (?) {i}	5-83	Capture Auto Pan
CPUT ? i, j	5-32	Download Capture Buffer
CRAT (?) {i}	5-83	Capture Rate
CSRC (?) {i}	5-55	Source Display
CSTP	5-83	Capture Stop
CSTR	5-83	Capture Start

---

**D**


---

DBIN ? d, j	5-108	Read Display d Bin Freq or Time
DBMR (?) {x}	5-40	dBm Reference Impedance
DCLR d	5-86	Data Table Clear
DDLT d, i	5-85	Data Table Delete
DDXW (?) d {, x}	5-44	d/dx Window
DELD	5-96	Delete Directory
DELF	5-96	Delete File
DFMT (?) {i}	5-43	Display Format
DINS d, i, j	5-85	Data Table Insert bin j at line i
DISP (?) d {, i}	5-43	Display Update Mode
DMAX ? d {, i}	5-85	Data Table Length
DNAM (?) {s}	5-96	Upkeep File Name
DREF d, i	5-104	Display Ref
DSPB ? d {, j}	5-107	Read Display d Binary
DSPE (?) {i} {, j}	5-122	Display Status Enable
DSPN ? d	5-107	Display d Length
DSPS ? {i}	5-122	Display Status Read
DSPY ? d {, j}	5-107	Read Display d (bin j) ASCII
DTBL (?) d {, i}	5-85	Data Table On
DTRD ? d {, i}	5-86	Data Table Query
DUMP	5-103	Dump Display Data
DUMP	5-97	Dump Display Data

---

**E**

EDLY (?) i {,x}	5-91	Table Delay
EFIT d	5-91	Start Fit
EFSC (?) i {,x}	5-91	Table Frequency Scale
EGAN (?) i {,x}	5-92	Table Gain
ETRC (?) i, j	5-92	Table Item Query
ENPL (?) {,i}	5-91	Number Poles
ENZE (?) {,i}	5-91	Number Zeros
EPCT (?) {i}	5-90	Exceed Centile
EPLY (?) i,j,k{,x}	5-92	Polynomials
EPOL (?) i,j{,x,y}	5-92	Poles
ERES (?) i,j{,x,y}	5-92	Residues
ERNG (?) d {,i,j}	5-91	Fit Range
ERRE (?) {i} {, j}	5-121	Error Status Enable
ERRS ? {i}	5-121	Error Status Read
ESTP (?) {i}	5-90	Exceed Stop Index
ESTR (?) {i}	5-90	Exceed Start Index
ESYN i,d	5-91	Synthesis
ETRC (?) i {,j}	5-92	Table Trace
EU1L (?) d {,i}	5-62	Ch1 EU Label
EU1M (?) d {,i}	5-62	Ch1 Engineering Units Off/On
EU1U (?) d {,s}	5-62	Ch1 User Label
EU1V (?) d {,x}	5-62	Ch1 EU/Volt
EU2L (?) {i}	5-63	Ch2 EU Label
EU2M (?) {i}	5-62	Ch2 Engineering Units Off/On
EU2U (?) {s}	5-63	Ch2 User Label
EU2V (?) {x}	5-63	Ch2 EU/Volt
EWTT (?) {,i}	5-91	Weighting Trace
EWTU (?) {,i}	5-91	Weighting
EXCE d, i	5-90	Calculate Exceedance

---

## 5-14 Index of Commands

---

EZER (?) i,j{x,y}	5-92	Zeros
<hr/>		
<b>F</b>		
FAVG (?) d {, i}	5-68	Average On
FAVM (?) d {, i}	5-68	Display Average
FAVN (?) d {, i}	5-68	Average Number
FAVT (?) d {, i}	5-68	Average Type
FBAS (?) d {, i}	5-28	Base Frequency
FBAS (?) d {, i}	5-28	FFT Base Frequency
FBIN ? d, x	5-109	Read Display d Bin Number
FCTR (?) d {, f}	5-28	FFT Center Frequency
FDIR (?) {s}	5-94	Save/Recall Directory
FEND (?) d {, f}	5-29	FFT End Frequency
FFMT (?) {i}	5-100	Frequency Format
FLIN (?) d {, i}	5-28	FFT Resolution
FNAM (?) {s}	5-94	Save/Recall File Name
FNXT ?	5-95	Read Disk Catalog
FOVL (?) d {, x}	5-68	FFT Time Record Increment
FRCL d	5-95	Disk to Display
FREE ?	5-94	Disk Free Space?
FREJ (?) d {, i}	5-69	FFT Overload Reject
FRST ?	5-95	Reset Disk Catalog
FSAV d	5-94	Display to Disk
FSPN (?) d {, f}	5-28	FFT Frequency Span
FSPN (?) d {, i}	5-28	Sampling Time
FSTR (?) d {, f}	5-28	FFT Start Frequency
FWFL (?) d {, i}	5-78	Force Length
FWIN (?) d {, i}	5-78	Window
FWTC (?) d {, i}	5-79	Expo TC
FXST ? s	5-94	File Exist?
<hr/>		
<b>G</b>		
GDIV (?) d {, i}	5-43	Number of Grid Divisions
GPOL (?) d {, i}	5-43	Nyquist Grid Type
GRID (?) d {, i}	5-43	Grid On/Off
<hr/>		
<b>H</b>		
HBIN (?) d {, i}	5-36	Histogram Bins
HDSP (?) d {, i}	5-50	Harmonic/Sideband Display
HLEN (?) d {, f}	5-36	Histogram Length
HPWR ? d	5-49	Query the Harmonic Power
HRDO (?) d {, i}	5-50	Harmonic/Sideband Readout
HRMN (?) d {, i}	5-49	Number of Harmonics
HRPT (?) d {, i}	5-36	Histogram Repeat
HTHD ? d, i	5-49	Query the Total Harmonic Distortion
<hr/>		
<b>I</b>		
I1AF (?) {i}	5-60	Ch1 Anti-Alias Filter
I1AR (?) {i}	5-60	Ch1 AutoRange Mode
I1AW (?) {i}	5-61	Ch1 A-Weight Filter
I1CP (?) {i}	5-60	Ch1 Input Coupling
I1GD (?) {i}	5-60	Ch1 Input Grounding
I1MD (?) {i}	5-60	Ch1 Input Mode
I1RG (?) {x}	5-60	Ch1 Input Range

---

---

I2AF (?) {i}	5-61	Ch2 Anti-Alias Filter
I2AR (?) {i}	5-61	Ch2 AutoRange Mode
I2AW (?) {i}	5-61	Ch2 A-Weight Filter
I2CP (?) {i}	5-61	Ch2 Input Coupling
I2GD (?) {i}	5-61	Ch2 Input Grounding
I2MD (?) {i}	5-61	Ch2 Input Mode
I2RG (?) {x}	5-61	Ch2 Input Range
IAOM (?) {i}	5-61	Auto Offset
ILEN (?) {i}	5-65	Capture Playback Length
IMOD (?) {i}	5-65	Capture Playback Mode
INPC ? i	5-122	Input Ovld Read
INPC ? i	5-122	Input Ovld Read
INPE (?) {i} {, j}	5-122	Input Status Enable
INPS ? {i}	5-122	Input Status Read
INSE (?) {i} {, j}	5-121	Instrument Status Enable
INST ? {i}	5-122	Instrument Status Read
ISPD (?) {i}	5-65	Capture Playback Speed
ISRC (?) {i}	5-60	Input Source
ISTR (?) {i}	5-65	Capture Playback Start

---

**K**

KCLK (?) {i}	5-100	Key Click
KEYP i	5-104	Key Press
KNOB i	5-106	Knob

---

**L**

LALM (?) d {, i}	5-87	Limit Beep
LCLR d	5-87	Clear Limits
LDLT d, i	5-88	Delete Limit Segment
LFAL ? d	5-87	Limit Test Fail?
LINK (?) {i}	5-60	Analyzer Configuration
LMAX (?) d {, i}	5-87	Last Limit Segment #
LOCL (?) {i}	5-117	Local/Remote
LSEG (?) d, i {, j, x0, y0, x1, y1}	5-87	Set Limit Segment i
LSFT d, x	5-88	Shift Limit Segments
LSON (?) d {, i}	5-87	Show Limit Segments
LTST (?) d {, i}	5-87	Limit Testing On

---

**M**

MALC i, j, k	5-84	Allocate Memory
MBIN d, i	5-45	Marker Move to bin i
MDIR s	5-94	Make Directory
MEAS (?) d {, i}	5-37	Display Measurement
MGRP (?) d {, i}	5-37	Measurement Group
MKCN d	5-104	Center of FFT Span to Marker
MKCN d	5-46	Center of FFT Span to Marker
MKMD (?) d {, i}	5-45	Marker Mode
MKMN d	5-104	Move the Marker to the Minimum
MMCA ? d	5-52	Frequency/Damping Results
MMEM ?	5-84	Total Memory Available
MREL (?) d {, i}	5-47	Normal Marker Relative Mode
MRKB ? d	5-45	Query the Marker Bin
MRKR (?) d {, i}	5-45	Marker Tracking
MRKX ? d	5-45	Query the Marker X Position

---

## 5-16 Index of Commands

---

MRKY ? d	5-45	Query the Marker Y Position
MRKZ ? d	5-46	Marker Z Read
MRON (?) d	5-103	Set Marker to Ref
MROX (?) d {, x}	5-47	Normal Marker X Offset
MROY (?) d {, x}	5-47	Normal Marker Y Offset
MSAA ?	5-89	Marker Max A Read
MSAB ?	5-89	Marker Max B Read
MSAO (?) d {, i}	5-89	Marker Statistics On
MSEA ?	5-89	Marker Mean A Read
MSEB ?	5-89	Marker Mean B Read
MSEK (?) d {, i}	5-47	Normal Marker Seeks Mode
MSIA ?	5-89	Marker Min A Read
MSIB ?	5-89	Marker Min B Read
MSRS	5-89	Reset Marker Statistics
MSSA ?	5-89	Marker Std Dev A Read
MSSB ?	5-89	Marker Std Dev B Read
MWFL d, i	5-45	Marker Move Waterfall to record i
MWID (?) d {, i}	5-47	Normal Marker Width
MXRL (?) d {, i}	5-47	Normal Marker X Rel Mode

---

### N

NAMP (?) {x}	5-56	Noise Amplitude
NAVG ? d	5-68	Averages Completed?
NAVG ? d	5-71	Octave Averages Completed?
NBUR (?) {x}	5-56	Noise Burst Percentage
NOTE i, j {, k, l, m, s}	5-99	Display Note
NPER (?) {x}	5-56	Noise Source Period
NTYP (?) {i}	5-56	Noise Type

---

### O

O1BN (?) d {, i}	5-35	Track 1 BNC
O1TK (?) d {, f}	5-35	Track 1 Order
O2BN (?) d {, i}	5-35	Track 2 BNC
O2TK (?) d {, f}	5-35	Track 2 Order
OCHN (?) d {, i}	5-30	Octave Channels
OCNF (?) d {, i}	5-71	Octave Confidence Level
ODLT (?) d {, f}	5-34	Delta Order
OHIB (?) d {, f}	5-30	Octave Highest Band
OIMP (?) d {, i}	5-71	Octave Power Bin
OLAT (?) d {, i}	5-72	Octave Linear Average Mode
OLOB (?) d {, f}	5-30	Octave Lowest Band
OMAX (?) d {, f}	5-34	Max Order
ONPT (?) d {, i}	5-34	Track Points
ORES (?) d {, i}	5-30	Octave Resolution
ORMN (?) d {, f}	5-34	Min RPM
ORMX (?) d {, f}	5-34	Max RPM
OSTO (?) d {, i}	5-35	Track Storage Mode
OTIM (?) d {, x}	5-71	Octave Average Time
OTRK (?) d {, i}	5-34	Order Tracking
OTYP (?) d {, i}	5-71	Octave Average Type
OUTX (?) {i}	5-100	Output Remote Interface
OVRM (?) {i}	5-100	Override Remote

---



**P**


---

PAUS	5-102	Pause
PAVA	5-69	Accept Preview
PAVO (?) d {, i}	5-69	Average Preview
PAVR	5-69	Reject Preview
PAVT (?) d {, x}	5-69	Preview Time
PBLK (?) {i}	5-98	Print Black
PBRI (?) {i}	5-98	Print Bright
PCIC (?) {i}	5-98	GPIB Control
PDIM (?) {i}	5-98	Print Dim
PDST (?) {i}	5-97	Print/Plot Destination
PFIL (?) {i}	5-97	File Start Number
PGRF (?) {i}	5-98	Print Graph
PHSL (?) d {, x}	5-44	Phase Suppress Threshold
PLAY i	5-106	Play Tune
PLGD (?) {i}	5-98	Plotter Grid Pen
PLMK (?) {i}	5-98	Plotter Marker Pen
PLOT	5-103	Plot Screen
PLOT	5-98	Plot Screen
PLTA (?) {i}	5-98	Plotter GPIB Address
PLTP (?) {i}	5-98	Vector/Plotter Type
PLTR (?) {i}	5-98	Plotter Trace Pen
PLTX (?) {i}	5-98	Plotter Text Pen
POUT (?) {i}	5-97	Print Screen hardkey
PRNT	5-103	Print Screen
PRNT	5-97	Print Screen
P RTP (?) {i}	5-97	Bitmap/Printer Type
PSCR (?) {i}	5-97	Bitmap Area
PSDU (?) d {, i}	5-40	PSD Units On/Off

---

**R**

R CRF d, i	5-103	Trace to Reference
R CTR d, i	5-103	Trace i to Display d
REFY ? d, j	5-107	Read Ref Display d bin j
RFDR(?) {i}	5-118	Reference Node Direction
RFNA(?) {s}	5-118	Reference Node Name
RFNU(?) {i}	5-118	Reference Node Number
RPMF (?) {i}	5-43	RPM Frequency
RSDR(?) {i}	5-118	Response Node Direction
RSNA(?) {s}	5-119	Response Node Name
RSNU(?) {i}	5-119	Response Node Number

---

**S**

S1AM (?) {x}	5-54	Sine Amplitude 1
S1FR (?) {f}	5-54	Sine Frequency 1
S2AM (?) {x}	5-54	Sine Amplitude 2
S2FR (?) {f}	5-54	Sine Frequency 2
SARS (?) d {, i}	5-32	Swept Sine Auto Resolution
SAVR (?) {i}	5-100	Screen Saver On
SBRI (?) {i}	5-100	Screen Brightness
SCON (?) {i}	5-104	Screen Contrast
SDBN (?) d {, i}	5-50	Number of Sidebands
SDLY (?) {i}	5-100	Screen Saver Delay
SFST (?) d {, x}	5-33	Swept Sine Faster Threshold

---

## 5-18 Index of Commands

---

SICY (?) d {, i}	5-73	Swept Sine Int. Cycles
SITM (?) d {, x}	5-73	Swept Sine Int. Time
SMAX (?) {x}	5-59	Swept Sine Max Level
SNAP d	5-104	Snap Reference
SNPS (?) d {, i}	5-33	Swept Sine Number of Points
SOFF (?) {x}	5-54	Sine Offset
SOFF (?) {x}	5-59	Swept Sine Offset
SPWR ? d, i	5-50	Query the Sideband Power
SRAT (?) {x}	5-58	Swept Sine Ramp Rate
SRCL i	5-95	Recall Settings
SRCO (?) {i}	5-53	Source On/Off
SRMP (?) {i}	5-58	Swept Sine Source Ramping Off/On
SRPT (?) d {, i}	5-32	Swept Sine Repeat Mode
SSAL (?) {i}	5-58	Swept Sine Auto Level
SSAM (?) {x}	5-58	Swept Sine Amplitude
SSAV	5-95	Settings to Disk
SSCY (?) d {, i}	5-73	Swept Sine Settle Cycles
SSEP (?) d {, f}	5-50	Sideband Separation
SSFR ?	5-32	Swept Sine Progress?
SSKP (?) d {, i}	5-33	Swept Sine Max Step Size
SSLL (?) {x}	5-59	Swept Sine Lower Limit
SSLO (?) d {, x}	5-33	Swept Sine Slower Threshold
SSRF (?) {x}	5-58	Swept Sine Ideal Reference
SSTM (?) d {, x}	5-73	Swept Sine Settle Time
SSTP (?) d {, f}	5-32	Swept Sine Stop Frequency
SSTR (?) d {, f}	5-32	Swept Sine Start Frequency
SSTY (?) d {, i}	5-32	Swept Sine Sweep Type
SSUL (?) {x}	5-58	Swept Sine Upper Limit
STMD (?) {i}	5-66	Triggered Source Mode
STRT	5-102	Start/Reset
STYP (?) {i}	5-53	Source Type
SVNI(?) {i}	5-118	Save Nodal Information
SVRF d, i	5-103	Reference to Trace
SVTR d, i	5-102	Display d to Trace i

---

### T

TARB i	5-57	Trace to Arb
TARM	5-66	Manual Trigger Arm
TASC ? i, n	5-109	Download Trace i Ascii
TASC ? i, n	5-109	Download Trace i Ascii
TAVM (?) { i}	5-69	Trigger Average Mode
TD1C (?) d {,i}	5-62	Ch1 Transducer Convert
TD2C (?) d {,i}	5-63	Ch2 Transducer Convert
TDLA (?) {x}	5-66	Trigger Delay A
TDLB (?) {x}	5-66	Trigger Delay B
TGET ? i	5-111	Upload Trace i Buffer
TGET ? i	5-111	Upload Trace i Buffer
TIAS (?) {x}	5-67	Time Arm Increment
TIME (?) {i, j, k}	5-100	Time
TLOD ? i, n	5-109	Download Trace i Binary
TLOD ? i, n	5-109	Download Trace i Binary
TLVL (?) {x}	5-66	Trigger Level
TMAN	5-66	Manual Trigger
TMOD (?) {i}	5-66	Trigger Arming Mode

---

TONE i, j	5-106	Tone
TPUT ? i	5-112	Download Trace i Buffer
TPUT ? i	5-112	Download Trace i Buffer
TRCL i	5-96	Disk to Trace i
TRCL i	5-96	Disk to Trace i
TRDM (?) {i}	5-67	Delta RPM Mode
TRDR (?) {x}	5-67	Delta RPM
TRSM (?) {i}	5-67	Start RPM Mode
TRSR (?) {x}	5-67	Starting RPM
TRWI i	5-79	Trace i to User Window
TSAV i	5-96	Trace i to Disk
TSLP (?) {i}	5-66	Trigger Slope
TSRC (?) {i}	5-66	Trigger Source

---

**U**

UNDB (?) d {, i}	5-40	dB Units
UNIT? d	5-40	Display Units
UNPH (?) d {, i}	5-40	Phase Units
UNPK (?) d {, i}	5-40	Peak Units
UNST d	5-30	Unsettle Measurement
USRC (?) i {, x, y}	5-77	User Constant i
USRF (?) i {, j, k, l ...}	5-74	FFT User Function i
USRH (?) i {, j, k, l ...}	5-74	Time/Histogram User Function i
USRO (?) i {, j, k, l ...}	5-74	Octave User Function i
USRR (?) i {, j, k, l ...}	5-74	Correlation User Function i
USRS (?) i {, j, k, l ...}	5-74	Swept Sine User Function i
USRT (?) i {, j, k, l ...}	5-74	Order User Function i

---

**V**

VIEW (?) d {, i}	5-39	Display View
------------------	------	--------------

---

**W**

WANG (?) d {, i}	5-81	Waterfall Angle
WAVA ? d	5-80	Waterfall Records Stored?
WDSP (?) d {, i}	5-80	Waterfall Display
WFSB (?) d {, i}	5-80	Waterfall Save Option
WFSK (?) d {, i}	5-80	Waterfall FFT Skip
WFST (?) d {, i}	5-81	Waterfall Fast Angles
WGET ?	5-114	Upload Waterfall Buffer
WGET ?	5-114	Upload Waterfall Buffer
WHID (?) d {, i}	5-81	Waterfall Hidden Lines
WHIT (?) d {, i}	5-81	Waterfall Trace Height
WITR d, i	5-79	Window to Trace i
WOSK (?) d {, x}	5-80	Waterfall Octave Skip
WPUT ?	5-115	Download Waterfall Buffer
WPUT ?	5-115	Download Waterfall Buffer
WREV (?) d {, i}	5-81	Waterfall Paused Drawing
WSLC d, i, j	5-82	Waterfall Slice to Trace
WSTO (?) d {, i}	5-80	Waterfall Storage
WSYM (?) d {, i}	5-79	Window Form
WTHR (?) d {, i}	5-81	Waterfall Threshold
WTOT (?) d {, i}	5-80	Waterfall Total Count
WTRC d, i, j	5-81	Waterfall Record to Trace
WVCT (?) d {, i}	5-81	Waterfall View Count

---

## 5-20 Index of Commands

---

### X

XAXS (?) d {, i}	5-43	X Axis Scale Type
XCEN (?) d {, x}	5-41	X Center (Polar)
XDIV (?) d {, x}	5-41	X/Division (Polar)
XPAN (?) d {, i}	5-41	X Pan
XZOM (?) d {, i}	5-42	X Zoom

---

### Y

Y2DV (?) d {, x}	5-41	Y/Division (Polar)
YCEN (?) d {, x}	5-41	Y Center (Polar)
YDIV (?) d {, x}	5-41	Y/Division
YMAX (?) d {, x}	5-40	Y Maximum
YMDX d x,y	5-41	Ymax + Y/Division
YMID (?) d {, x}	5-40	Y Midpoint
YMIN (?) d {, x}	5-41	Y Minimum

## Introduction

The SR785 Dynamic Signal Analyzer may be remotely programmed via either the RS232 or GPIB (IEEE-488) interfaces. Any computer supporting one of these interfaces may be used to program the SR785. Both interfaces are receiving at all times, however, the SR785 will send responses only to the Output Interface specified in the [System] <Remote> menu.



---

Use the OUTX command at the beginning of every program to direct the SR785 responses to the correct interface.

---

### Communicating With GPIB

The SR785 supports the IEEE-488.1 (1978) interface standard. It also supports the required common commands of the IEEE-488.2 (1987) standard. Before attempting to communicate with the SR785 over the GPIB interface, the SR785's Device Address must be set in the [System] <Remote> menu.

### Communicating With RS232

The SR785 is configured as a DCE (transmit on pin 3, receive on pin 2) device and supports CTS/DTR hardware handshaking. The CTS signal (pin 5) is an output indicating that the SR785 is ready, while the DTR signal (pin 20) is an input that is used to control the SR785's data transmission. If desired, the handshake pins may be ignored and a simple 3 wire interface (pins 2,3 and 7) may be used. The RS232 interface Baud Rate, Word Length, and Parity must be set in the [System] <Remote> menu.

### Screen Indicators And Queues

To assist in programming, the SR785 has 4 interface indicators which are displayed at the top of the screen.

The RS232/gpib indicator shows 'RS232' if the interface responses are directed to the RS232 serial port and 'gpib' and the address if the interface responses are directed to the GPIB port.

The comm indicator flashes 'RS232' when there is activity on the RS232 interface and 'GPIB' when there is activity on the GPIB interface. 'ERR' flashes whenever a computer interface error occurs, such as illegal command or out of range parameter is received.

The REM indicator is on whenever the SR785 is in a remote state (front panel locked out).

The SRQ indicator is on when the SR785 generates a service request. SRQ stays on until a GPIB serial poll is completed.

To help find program errors, the SR785 can display the interface buffers on the screen. This screen is accessed by <View Qs> in the [System] <Remote> menu. The last 256 characters received and transmitted by the SR785 are displayed.

### Command Format

Communications with the SR785 uses ASCII characters. Commands may be in either UPPER or lower case. A command to the SR785 consists of a four character command mnemonic with optional ?, arguments if necessary, and a command terminator. The command, arguments and terminator may be separated by spaces. The terminator must be a linefeed <lf> or carriage return <cr> on RS232, or a linefeed <lf> or EOI on GPIB. No command processing occurs until a terminator is received. Commands function identically on GPIB and RS232 whenever possible. Command mnemonics beginning with an asterisk '\*' are IEEE-488.2 (1987) defined common commands. These commands also function identically on RS232. Commands may require one or more parameters. Multiple parameters are separated by commas (,).

Multiple commands may be sent on one command line by separating them with semicolons (;).

There is no need to wait between commands. The SR785 has a 256 character input buffer and processes commands in the order received. If the buffer fills up, the SR785 will hold off handshaking on the GPIB and attempt to hold off handshaking on RS232. Similarly, the SR785 has a 256 character output buffer to store output until the host computer is ready to receive it. If either buffer overflows, both buffers are cleared and an error reported.

The present value of a particular parameter may be determined by querying the SR785 for its value. A query is formed by appending a question mark '?' to the command mnemonic and omitting the desired parameter from the command. Values returned by the SR785 are sent as a string of ASCII characters terminated by a carriage return <cr> on RS232 and by a line-feed <lf> on GPIB. If multiple query commands are sent on one command line (separated by semicolons, of course), the answers will be returned individually, each with a terminator.

#### Examples of Commands

TSLP 0 <lf>	Set the Trigger Slope to Positive
FCTR 1, 10E3 <lf>	Set the Center Frequency of DisplayB to 10000 Hz (10 kHz)
FCTR ? 1	Query the Center Frequency of DisplayB
*IDN ? <lf>	Query the Device Identification String
STRT <lf>	Start the measurement (same as [Start/Reset] key)

### Command Synchronization

IFC (Interface Ready, bit 7) in the Serial Poll status signals that the SR785 is ready to receive and execute a command. When a command is received, this bit is cleared, indicating that command execution is in progress. No other commands will be processed until this command is completed. Commands received during this time are stored in the buffer to be processed later. Only GPIB serial polling will generate a response while a command is in progress. When all pending commands have executed, the IFC bit is set again. By checking IFC with serial polls, a host computer can ensure that all previously sent commands have finished before sending a new command.

Since most commands execute very quickly, the host computer does not need to continually check the IFC bit. Commands may be sent one after another and they will be processed immediately.

However, some commands, such as file and print/plot commands and data transfer operations, may require a long time to execute. In addition, the host program may need to check that these operations executed without error. In these cases, after the command is sent, the status should be queried.

When using the GPIB interface, serial polling may be used to check the IFC bit in the Serial Poll status while an operation is in progress. After the IFC bit becomes set, signaling the completion of the command, then the IERR bit may be checked to verify successful completion of the command.

If the RS232 interface is used, or serial polling is not available, then the \*STB? query command may be used to read the Serial Poll status word. However, \*STB NEVER returns the IFC bit set (since \*STB is itself a command).

Since the SR785 processes one command at a time, status queries will not be processed until the previous operation is finished. Thus a response to a status query in itself signals that the previous command is finished. The query response may then be checked for various errors.

For example, the command line FSAV 0; ERRS ? <lf> will save the DisplayA data to disk and return the Error Status word when finished. The Disk Error bit may be checked to make sure that the FSAV (Display to Disk) command completed without error. Since the FSAV command may take a long time to execute, it is important that the host computer interface does not time out while waiting for the response to the ERRS ? query. In the case where the host interface times out before the ERRS ? response, the host program must wait before sending the ERRS ? query.

Interface commands, like their keyboard equivalents, may not take effect immediately if the analyzer is paused, or finished with a set of linear averages. If the display is paused or Done, and a command is sent to change a frequency parameter, the menu will reflect the new parameter, querying the parameter with a command will return the new parameter, but *the new parameter does not take effect until [Start] is pressed or the STRT command is sent*. You should always send STRT after a series of setup commands to ensure that the new parameters are in effect.

## Data Synchronization

Changing measurement parameters often invalidates the display data. For example, changing FFT spans or resolution requires some time before new data taken with the new span or resolution is available. Another example is changing the FFT average mode while averaging is on. In this case, the average is re-started and some time is required before a new average is completed.

For these types of commands, simply waiting for IFC to be set does not ensure that the display data reflects the operation just completed. Before querying or saving the display data, the Display Status word needs to be queried until either NEW, AVG or STL is set

indicating new data is available, average completed or new settled data is available. This will ensure that the display data has been updated since the previous operation.

Remember, these status bits are only reset when read. The status word needs to be cleared before waiting for the desired bits to become set.

For example, to change the span and wait until new data is available, a program might use the following procedure.

send (“\*CLS; FSPN 0,6400; STRT”) *\*CLS clears all status words, FSPN changes the span of DisplayA to 6.4 kHz. STRT starts the measurement with the new span*

query (“DSPS?”) *DSPS? reads the Display Status*  
until NEWA and STLA are set *NEWA&STLA indicates new settled data is available at the new span.*

### Command Context Errors

Commands which perform menu softkey functions are not allowed unless the corresponding softkey function is allowed. If the softkey is not available in the current context, either shown in gray or not displayed at all, then the corresponding set command is not allowed.

In general, softkeys are disabled for two reasons - the function is specifically not allowed in the current measurement context, or the function is only applicable for certain measurement groups.

For example, FFT <Start Frequency> is grayed when the measurement is full span. The FSTR set command is not allowed in this situation. However, the FSTR? query command is still allowed.

An example of group specific parameter are the commands which set Swept Sine parameters. These set commands are not allowed unless the measurement group is Swept Sine (just like the menu interface).

### Example Program

An example program is included at the end of this chapter. This program is a good reference for writing your own programs to control the SR785.



## Command Syntax

The four letter mnemonic (shown in CAPS) in each command sequence specifies the command. The rest of the sequence consists of parameters. Parameters shown in { } are not always required. Generally, parameters in { } are required to set a value in the SR785. Multiple parameters are separated by commas. Multiple commands may be sent on one command line by separating them with semicolons (;).

The present value of a parameter may be determined by sending a query command. Commands that may be queried have a question mark in parentheses (?) after the mnemonic. Commands that may ONLY be queried have a ? after the mnemonic. Commands that MAY NOT be queried have no ?. A query is formed by including the question mark ? after the command mnemonic and omitting the queried parameter from the command. The query parameters shown in { } are NOT sent with a query. The query returns the value of these parameters. Values are returned as a string of ASCII characters (unless otherwise noted).

Do NOT send ( ) or { } as part of the command.

For example, the command sequence FSTR (?) d {, f} is used as follows.

```
FSTR 0, 1.024E3      Set the Start Frequency of DispA to 1024 Hz.
FSTR ? 0           Query the Start Frequency of DispA.
```

Variables are defined as follows.

d	display (0=DispA, 1=DispB, 2=Both Displays) d=2 is not allowed for queries.
i, j, k, l, m	integers
x, y	real numbers
f	frequency in Hz
s	text string

Commands to set values which may be different for each display require the d parameter. These values must be queried separately for each display (d=2 is not allowed for queries) even if the parameter is linked.

All numeric variables may be expressed in integer, floating point or exponential formats (i.e., the number five can be either 5, 5.0, or 0.5E1). Strings are sent as a sequence of ASCII characters.

### Commands and Units

Some parameters display a choice of units in the entry field when a numeric value is entered. These same units can be sent in the corresponding command. For instance, when a numeric value is entered for the span, the following choice of units appears:

kHz    Hz    mHz     $\mu$  Hz

Likewise, to set the span to 51.2 kHz one could send any one of these five equivalent commands:

```
FSPN 2, 51200
FSPN 2, 51.2 kHz
FSPN 2, 51200 Hz
FSPN 2, 51200000 mHz
FSPN 2,5.12E10 uHz
```

Actually, it is only necessary to send sufficient characters to uniquely identify the units from the list. Thus FSPN 2, 51.2 k would set the span to 51.2 kHz as well.

Regardless of the units used to set the parameter, the analyzer almost always returns the result of queries in one unit. So if the span were queried after any of the above commands were sent, the analyzer would return 51200.

It is never necessary to send a unit string with a set command. Each command has a "default unit" which it will use if no units are sent with the value. The default units are also the units used when the analyzer responds to queries.

Some commands must identify the units in which queries are responded to. Examples would be the source amplitude and input range commands which can return values in dB and linear units. In these cases, the response to the query is in the form: f,j where f is the actual floating point return value, and j is the index into the list of units which uniquely identifies the units of f. These commands are identified in the individual command description.

In the detailed command descriptions, lists of allowed units are indicated in angle brackets <>. Within the list, the default units are listed surrounded by square brackets []. For instance, for the above command, the allowed unit list would be written as:

```
<kHz, [Hz], mHz, μHz>
```

### Help

The detailed command list is available on screen by pressing [Help/Local] to enter the help system. Press [4] to show the command list. Commands are also cross referenced in the help about each key.





---

## Things to Remember!

---

### Output Interface (RS232 or GPIB)

All responses are directed to the interface selected by <Output To> in the [System] <Remote> menu, *regardless* of which interface received the query. Use the OUTX command to select the correct interface at the beginning of every program.

### Active Display and Linking

Set commands which require a 'd' parameter specify DisplayA (0), DisplayB (1) or Both Displays (2) *without* regard to the active display or any parameter linking. Active display and linking are screen features which do *not* apply to interface commands.

### Analyzer Configuration and Linking

When the Analyzer Configuration is set to Dual Channel, many commands that have a 'd' (display) parameter require d be set to '2' (both channels). This reflects the fact that the two displays cannot have different frequency parameters in the dual channel configuration.

### Command Context Errors

Commands which perform softkey functions are not allowed unless the corresponding softkey function is allowed. If the softkey is not available, either shown in gray or not displayed at all, then the corresponding set command is not allowed.

Be careful to send commands in the correct order to avoid context errors. A good practice is to send the commands in the same order as programming the instrument using the softkeys.

## Frequency Commands (FFT and Correlation)

---

### **FSPN (?) d {, f} <kHz, [Hz], mHz, $\mu$ Hz>**

The FSPN command sets (queries) the FFT frequency Span of display d. The parameter f is a real number representing the span in the selected units. The available spans depend upon the FFT Base Frequency.

The set command changes the FFT Span to the available span closest to f.

The query command always returns the span in Hz.

This command is valid only when the Measurement Group is FFT or Correlation. The set command requires display d to be Live.

---

### **FLIN (?) d {, i}**

The FLIN command sets (queries) the Resolution of display d when the measurement group is FFT or Correlation. The parameter i selects 100 (0), 200 (1), 400 (2) or 800 (3) lines. FLIN does not change the Span.

This command is valid only when the Measurement Group is FFT or Correlation. The set command requires display d to be Live.

---

### **FBAS (?) d {, i}**

The FBAS command sets (queries) the FFT Base Frequency of display d. The parameter i selects 100.0 kHz (0) or 102.4 kHz (1).

To set the Base Frequency, d must be 2 (both displays). There is only one Base Frequency for both displays. Changing the Base Frequency will affect all of the FFT and source frequency parameters.

This command is valid only when the Measurement Group is FFT, Correlation, or Time/Histogram. The set command requires at least one display to be Live.

---

### **FSTR (?) d {, f} <kHz, [Hz], mHz, $\mu$ Hz>**

The FSTR command sets (queries) the FFT Start frequency of display d. The parameter f is a frequency (real number of the specified units). Values of f which would cause the span to exceed the 0 to 102.4 (100.0) kHz range cause an error.

This command is valid only when the Measurement Group is FFT. The set command requires display d to be Live.

Since correlation is always a baseband measurement, this command is not valid when the measurement group is correlation.

---

### **FCTR (?) d {, f} <kHz, [Hz], mHz, $\mu$ Hz>**

The FCTR command sets (queries) the FFT Center frequency of display d. The parameter f is a frequency (real number of the specified units). Values of f which would cause the span to exceed the 0 to 102.4 (100.0) kHz range cause an error.

---

This command is valid only when the Measurement Group is FFT. The set command requires display d to be Live.

Since correlation is always a baseband measurement, this command is not valid when the measurement group is correlation.

---

**FEND (?) d {, f} <kHz, [Hz], mHz,  $\mu$ Hz>**

The FEND command sets (queries) the FFT End frequency of display d. The parameter f is a frequency (real number of the specified units). Values of f which would cause the span to exceed the 0 to 102.4 (100.0) kHz range cause an error.

This command is valid only when the Measurement Group is FFT. The set command requires display d to be Live.

Since correlation is always a baseband measurement, this command is not valid when the measurement group is correlation.

---

**UNST d**

The UNST command unsettles the measurement of display d. The measurement is not actually perturbed by the UNST command. The settling status is set to unsettled and the full settling time of the measurement is required before the status returns to settled.

The measurement is unsettled by changing any one of several measurement parameters. For example, changing the input range or FFT span will unsettle the measurement. If the signal comes from an external source and is changed in such a way as to require the measurement to settle, it is convenient to use the UNST command and wait for settling to finish.

This command is valid only when the Measurement Group is FFT or Octave.

## Frequency Commands (Octave)

---

**OHIB (?) d {, f} <kHz, [Hz], mHz,  $\mu$ Hz>**

The OHIB command sets (queries) the Highest Band in the Octave measurement of display d. The parameter f is a frequency (real number in the specified units). The Highest Band is specified to the nearest octave.

For 1 Channel octave analysis, the set command requires d=2.

This command is valid only when the Measurement Group is Octave. The set command requires display d to be Live.

---

**OLOB (?) d {, f} <kHz, [Hz], mHz,  $\mu$ Hz>**

The OLOB command sets (queries) the Lowest Band in the Octave measurement of display d. The parameter f is a frequency (real number in the specified units). The Lowest Band is specified to the nearest octave.

For 1 Channel octave analysis, the set command requires d=2.

This command is valid only when the Measurement Group is Octave. The set command requires display d to be Live.

---

**ORES (?) d {, i}**

The ORES command sets (queries) the Octave Resolution of display d. The parameter i selects 1/1 (0), 1/3 (1) or 1/12 (2) octave analysis. Changing the resolution may change the values of the Highest and Lowest Band for display d.

For 1 Channel octave analysis, the set command requires d=2.

This command is valid only when the Measurement Group is Octave. The set command requires display d to be Live.

---

**OCHN (?) d {, i}**

The OCHN command sets (queries) the Octave Channels. The parameter i selects 1 (0) or 2 (1) channel octave analysis.

This command is valid only when the Measurement Group is Octave. The set command requires both displays to be Live.

---

**UNST d**

The UNST command unsettles the measurement of display d. The measurement is not actually perturbed by the UNST command. The settling status is set to unsettled and the full settling time of the measurement is required before the status returns to settled.

The measurement is unsettled by changing any one of several measurement parameters. For example, changing the input range or FFT span will unsettle the measurement. If the signal comes from an external source and is changed in such

---

a way as to require the measurement to settle, it is convenient to use the UNST command and wait for settling to finish.

This command is valid only when the Measurement Group is FFT, Correlation or Octave.

## Frequency Commands (Swept Sine)

---

**SSTR (?) d {, f} <kHz, [Hz], mHz,  $\mu$ Hz>**

The SSTR command sets (queries) the swept sine Start frequency of display d. The parameter f is a frequency (real number of specified units). The range of f is 0.001 to 102.4E3. The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Swept Sine. The set command requires a display to be Live.

---

**SSTP (?) d {, f} <kHz, [Hz], mHz,  $\mu$ Hz>**

The SSTP command sets (queries) the swept sine Stop frequency of display d. The parameter f is a frequency (real number of specified units). The range of f is 0.001 to 102.4E3. The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Swept Sine. The set command requires a display to be Live.

---

**SSFR ?**

The SSFR command queries which sweep point was the most recently measured. While a sweep is in progress, SSFR? returns values from 0 to the Number Of Points - 1.

After a sweep is started, SSFR? will not return 0 until the first point is measured. This may take some time if Source Ramping is On or if the Settling Time is very long.

This command is valid only when the Measurement Group is Swept Sine.

---

**SRPT (?) d {, i}**

The SRPT command sets (queries) the swept sine Repeat Mode of display d. The parameter i selects Single Shot (0) or Continuous (1). The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Swept Sine. The set command requires a display to be Live.

---

**SSTY (?) d {, i}**

The SSTY command sets (queries) the swept sine Sweep Type of display d. The parameter i selects Linear (0) or Logarithmic (1). The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Swept Sine. The set command requires a display to be Live.

---

**SARS (?) d {, i}**

The SARS command sets (queries) the swept sine Auto Resolution Mode of display d. The parameter i selects Off (0) or On (1). The set command requires d=2 (both displays).

---



This command is valid only when the Measurement Group is Swept Sine. The set command requires a display to be Live.

---

**SNPS (?) d {, i}**

The SNPS command sets (queries) the swept sine Number Of Points of display d. The parameter i is a number of points from 10 to 2047. The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Swept Sine. The set command requires a display to be Live.

---

**SSKP (?) d {, i}**

The SSKP command sets (queries) the swept sine Maximum Number Of Skips of display d. The parameter i is a maximum number of skips from 2 to 200. The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Swept Sine. The set command requires a display to be Live.

---

**SFST (?) d {, x}**

The SFST command sets (queries) the swept sine Faster Threshold of display d. The parameter x is a level from 0.01 to 3.0 (dB). The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Swept Sine. The set command requires a display to be Live.

---

**SSLO (?) d {, x}**

The SSLO command sets (queries) the swept sine Slower Threshold of display d. The parameter x is a level from 0.05 to 6.0 (dB). The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Swept Sine. The set command requires a display to be Live.

---

## Frequency Commands (Order Analysis)

---

### **ORMX (?) d {, x}**

The ORMX command sets (queries) the maximum RPM for display d. The parameter x is the maximum rpm (a real number). The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Order. The set command requires a display to be live.

---

### **ORMN (?) d {, x}**

The ORMN command sets (queries) the minimum RPM for display d. The parameter x is the minimum rpm (a real number). The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Order. The set command requires both displays to be Live.

---

### **OMAX (?) d {, x}**

The OMAX command sets (queries) the maximum order for display d. The parameter x is the maximum order (a real number). The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Order. The set command requires a display to be live.

---

### **ODLT (?) d {, x}**

The ODLT command sets (queries) delta order for display d. The parameter x is delta order (a real number). The range of f is between Max Order/400 and 1. The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Order. The set command requires a display to be live.

---

### **OTRK (?) d {, i}**

The OTRK command sets (queries) tracking for display d. The parameter i selects Off (0) or On (1). The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Order. The set command requires a display to be Live.

---

### **ONPT (?) d {, i}**

The ONPT command sets (queries) number of points in a track measurement for display d. The parameter i is the number of points. The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Order. The set command requires a display to be Live.

---

**OSTO (?) d {, i}**

The OSTO command sets (queries) the track storage mode for display d. The parameter i selects Continuous (0) or One Shot (1). The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Order. The set command requires a display to be Live.

---

**O1TK (?) d {, x}**

The O1TK command sets (queries) the order associated with the track 1 measurement for display d. The parameter x is the order (a real number). The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Order. The set command requires a display to be Live.

---

**O1BN (?) d {, i}**

The O1BN command sets (queries) the input channel associated with the track 1 measurement for display d. The parameter i sets the input channel to Channel One (0) or Channel 2 (1). The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Order. The set command requires a display to be Live.

---

**O2TK (?) d {, x}**

The O2TK command sets (queries) the order associated with the track 2 measurement for display d. The parameter x is the order (a real number). The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Order. The set command requires a display to be Live.

---

**O2BN (?) d {, i}**

The O2BN command sets (queries) the input channel associated with the track 2 measurement for display d. The parameter i sets the input channel to Channel One (0) or Channel 2 (1). The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Order. The set command requires a display to be Live.

---

## Frequency Commands (Time/Histogram)

---

**FSPN (?) d {, f} <kHz, [Hz], mHz,  $\mu$ Hz>**

The FSPN command sets (queries) the sampling time for display d. If the base frequency is set to 102.4 kHz the sampling time ranges from 3.81 $\mu$ S to 2S in power of 2 increments. If the base frequency is set to 100 kHz the sampling time ranges from 3.91 $\mu$ S to 2.048S in power of 2 increments. The parameter f is a frequency (real number of specified units) related to the sampling time by the equation  $f = 25 / (64 \cdot t)$ .

The set command requires d=2 (both displays). The set command requires a display to be Live.

---

**HLEN (?) d {, x} <samples, records,  $\mu$ S, mS, [S]>**

The HLEN command sets (queries) the histogram length of display d. The parameter x is the histogram length in the specified units

This command is valid only when the Measurement Group is Time/Histogram. The set command requires the display d to be Live.

---

**HBIN (?) d {, i}**

The HBIN command sets (queries) the number of histogram bins for display d. The parameter i is related to the number of histogram bins by the equation  $\text{Bins} = 2^{*(i+2)}$ . The range of i is from 0 to 8.

This command is valid only when the Measurement Group is Time/Histogram. The set command requires the display d to be Live.

---

**FBAS (?) d {, i}**

The FBAS command sets (queries) the Base Frequency of display d. The parameter i selects 100.0 kHz (0) or 102.4 kHz (1).

To set the Base Frequency, d must be 2 (both displays). There is only one Base Frequency for both displays. Changing the Base Frequency will affect all of the FFT and source frequency parameters.

The set command requires at least one display to be Live.

---

**HRPT (?) d {, i}**

The HRPT command sets (queries) the histogram repeat mode for display d. The parameter i sets the repeat mode to Off (0) or On (1).

This command is valid only when the Measurement Group is Time/Histogram. The set command requires the display d = 2 (both displays)..

---

# Display Setup Commands

---

## MGRP (?) d {, i}

The MGRP command sets (queries) the Measurement Group of display d. The parameter i selects FFT (0), Correlation (1), Octave (2), Swept Sine (3), Order (4), ord Time/Histogram (5)..

Changing the Measurement Group changes the Frequency, Display Setup and Average menus. In addition, the Source and Capture menus may also change. The Measurement Group specific settings in these menus change to those last used with the new group. Parameters within these menus may only be changed with commands which are valid within the new group.

The set command requires d=2 (both displays).

---

## MEAS (?) d {, i}

The MEAS command sets (queries) the Measurement of display d. The parameter i selects the measurement from the list below. Only those measurements available in the current Measurement Group are allowed.

Each measurement has an associated view. Changing the Measurement changes the View to the view last used with the new Measurement.

The set command requires display d to be Live.

i	FFT Group <u>Measurement</u>
0	FFT 1
1	FFT 2
2	Power Spectrum 1
3	Power Spectrum 2
4	Time 1
5	Time 2
6	Windowed Time 1
7	Windowed Time 2
8	Orbit
9	Coherence
10	Cross Spectrum
11	Frequency Response
12	Capture Buffer 1
13	Capture Buffer 2
14	FFT User Function 1
15	FFT User Function 2
16	FFT User Function 3
17	FFT User Function 4
18	FFT User Function 5

## 5-38 Display Setup Commands

---

	Correlation Group
<u>i</u>	<u>Measurement</u>
19	Auto Correlation 1
20	Auto Correlation 2
21	Cross Correlation
22	Time 1
23	Time 2
24	Windowed Time 1
25	Windowed Time 2
26	Capture Buffer 1
27	Capture Buffer 2
28	Correlation Function 1
29	Correlation Function 2
30	Correlation Function 3
31	Correlation Function 4
32	Correlation Function 5
	Octave Group
<u>i</u>	<u>Measurement</u>
33	Octave 1
34	Octave 2
35	Capture 1
36	Capture 2
37	Octave User Function 1
38	Octave User Function 2
39	Octave User Function 3
40	Octave User Function 4
41	Octave User Function 5
	Swept Sine Group
<u>i</u>	<u>Measurement</u>
42	Spectrum 1
43	Spectrum 2
44	Normalized Variance 1
45	Normalized Variance 2
46	Cross Spectrum
47	Frequency Response
48	Swept Sine User Function 1
49	Swept Sine User Function 2
50	Swept Sine User Function 3
51	Swept Sine User Function 4
52	Swept Sine User Function 5
	Order Group
<u>i</u>	<u>Measurement</u>
53	Linear Spectrum 1
54	Linear Spectrum 2
55	Power Spectrum 1
56	Power Spectrum 2
57	Time 1

58	Time 2
59	Windowed Time 1
60	Windowed Time 2
61	RPM Profile
62	Orbit
63	Track 1
64	Track 2
65	Capture Buffer 1
66	Capture Buffer 2
67	Order User Function 1
68	Order User Function 2
69	Order User Function 3
70	Order User Function 4
71	Order User Function 5
	Time/Histogram Group
i	<u>Measurement</u>
72	Histogram 1
73	Histogram 2
74	PDF 1
75	PDF 2
76	CDF 1
77	CDF 2
78	Time 1
79	Time 2
80	Capture Buffer 1
81	Capture Buffer 2
82	Histogram User Function 1
83	Histogram User Function 2
84	Histogram User Function 3
85	Histogram User Function 4
86	Histogram User Function 5

A measurement should be setup by first choosing the Measurement Group, then the Measurement, then the View, then the Units. Finally, the display scale and references should be set.

---

**VIEW (?) d {, i}**

The VIEW command sets (queries) the View of display d. The parameter i selects the view from the list below.

Each view has associated units. Changing the View changes the units parameters to the values last used with the new view.

## 5-40 Display Setup Commands

---

<u>i</u>	<u>View</u>
0	Log Magnitude
1	Linear Magnitude
2	Magnitude Squared
3	Real Part
4	Imaginary Part
5	Phase
6	Unwrapped Phase
7	Nyquist
8	Nichols

A measurement should be setup by first choosing the Measurement Group, then the Measurement, then the View, then the Units. Finally, the display scale and references should be set.

---

### UNIT? d

The UNIT query command queries the Units of display d. The unit string depends on the measurement, the view, and the settings of dB Units, pkUnits, PSD Units and Phase Units.

---

### UNDB (?) d {, i}

The UNDB command sets (queries) the dBUnits setting for display d. The parameter i sets dBUnits to Off (0), On(1), dBm (2), and dBspl (3).

---

### UNPK (?) d {, i}

The UNPK command sets (queries) the pk units setting for display d. The parameter i sets pkUnits to Off (0), pk (1), rms (2), and pp (3).

---

### PSDU (?) d {, i}

The UNPK command sets (queries) psd units setting for display d. The parameter i sets psd Units to Off (0), or On (1).

---

### UNPH (?) d {, i}

The UNPH command sets (queries) phase units setting for display d. The parameter i sets phase Units to Degrees (0), or Radians (1).

---

### DBMR (?) {x}

The DBMR command sets (queries) the dBm Reference Impedance (for BOTH displays). The parameter x is a real value of Ohms. This affects the calculation of dBm units in both displays.

---

### YMAX (?) d {, x}

The YMAX command sets (queries) the Y Maximum (top reference) of display d. The parameter x is a real number in the display units. This command is not valid when the View is Nichols or Nyquist.

---

### YMID (?) d {, x}

The YMID command sets (queries) the Y Midpoint (center reference) of display d. The parameter x is a real number in the display units. This command is not valid when the View is Nichols or Nyquist.



**YMIN (?) d {, x}**

The YMIN command sets (queries) the Y Minimum (bottom reference) of display d. The parameter x is a real number in the display units. This command is not valid when the View is Nichols or Nyquist.

---

**YDIV (?) d {, x}**

The YDIV command sets (queries) the Y/Division scale of display d. The parameter x is a real number in the display units. This command is valid for all Views.

Changing the scale of a display may change the Y Max and Y Mid values. Always use a YMAX, YMID or YMIN command after the YDIV command.

---

**YMDX d {, x, y}**

The YMDX command sets the YMax value and Y/Div value simultaneously. The parameter x is the value of Ymax. The parameter y is the value of Y/div. This command is not valid when the View is Nichols or Nyquist..

---

**XCEN (?) d {, x}**

The XCEN command sets (queries) the X Center of display d. The parameter x is a real number in the display units. This command is only valid when the View is Nichols or Nyquist.

---

**XDIV (?) d {, x}**

The XDIV command sets (queries) the polar X/Division scale of display d. The parameter x is a real number in the display units. This command is only valid when the View is Nichols or Nyquist.

Changing the scale of a display may change the X Center value. Always use the XCEN command after the XDIV command.

---

**YCEN (?) d {, x}**

The YCEN command sets (queries) the Y Center of display d. The parameter x is a real number in the display units. This command is only valid when the View is Nichols or Nyquist.

---

**Y2DV (?) d {, x}**

The Y2DV command sets (queries) the polar Y/Division scale of display d. The parameter x is a real number in the display units. This command is only valid when the View is Nichols or Nyquist.

Changing the scale of a display may change the Y Center value. Always use the YCEN command after the Y2DV command.

---

**XPAN (?) d {, i}**

The XPAN command sets (queries) the Pan (translation) of display d. The parameter i is the left most displayed bin. This command is only valid when the display is Zoomed (expanded).

---

## 5-42 Display Setup Commands

---

### **XZOM (?) d {, i}**

The XZOM command sets (queries) the Zoom factor (X axis expand) of display d. The parameter i is the zoom factor (1-5). This command is not valid when the X axis is logarithmic or when the View is Nyquist or Nichols.

---

# Display Options Commands

---

**DISP (?) d {, i}**

The DISP command sets (queries) the Display Update Mode of display d. The parameter i selects Off-Line (0) or Live (1).

Many parameters (Frequency, Measurement, etc.) are not adjustable for an Off-Line display. Commands which set these parameters for an Off-Line display are not valid.

---

**DFMT (?) {i}**

The DFMT command sets (queries) the Display Format. The parameter i selects Single (0), Dual (1), or Overlay (2).

Use the ACTD command to select the active display in the Single format. In the Single format, the inactive display is still accessible via commands.

---

**ACTD (?) i**

The ACTD command sets the Active Display to display i. The parameter i selects DisplayA (0) or DisplayB (1).

---

**XAXS (?) d {, i}**

The XAXS command sets (queries) the X Axis Scale Type of display d. The parameter i selects Linear (0) or Logarithmic (1).

This command is valid only when display d is an FFT measurement with frequency as the X axis. FFT measurements with a time X axis (Time, WinTime, Capture or User Time Function) always have a linear scale.

---

**RPMF (?) {i}**

The RPMF command sets (queries) the RPM frequency scaling for both displays. The parameter i selects Hz (0) or RPM (1) as the units for frequency.

---

**GRID (?) d {, i}**

The GRID command sets (queries) the Grid On/Off Mode of display d. The parameter i selects Off (0) or On (1).

---

**GDIV (?) d {, i}**

The GDIV command sets (queries) the number of Grid Divisions of display d. The parameter i selects 8 (0), 10 (1), 12 (2) or 15 (3).

Changing the Grid Divisions changes the vertical scaling (Y/div) and horizontal scaling (X/div) (Nyquist and Nichols views).

---

**GPOL (?) d {, i}**

The GPOL command sets (queries) the Nyquist Grid type for display d. The parameter i selects Rectangular Grid (0) or Polar Grid (1).

---

## 5-44 Display Options Commands

---

---

### **PHSL (?) d {, x}**

The PHSL command sets (queries) the Phase Suppress Threshold of display d. The parameter x is a real magnitude squared value. This affects the calculation of phase for display d.

---

### **DDXW (?) d {, x}**

The DDXW command sets (queries) the d/dx Window of display d. The parameter x is a percentage of the display width. This affects the calculation of d/dx and group delay for user math functions in display d.

---

# Marker Commands

---

**MRKR (?) d {, i}**

The MRKR command sets (queries) the Marker Tracking of display d. The parameter i selects Off (0), On (1), Track (2) or Link (3).

---

**MKMD (?) d {, i}**

The MKMD command sets (queries) the Marker Mode of display d. The parameter i selects Normal (0), Harmonic (1), Sideband (2), Band (3), or Frequency/Damping (4).

All modes are allowed in the FFT and Order measurement groups.

Harmonic, Sideband, and Frequency Damping are not allowed in Octave measurement group.

Harmonic, Sideband, Band, and Frequency/Damping are not allowed in the Swept Sine, Time/Histogram, or Correlation measurement group.

---

**MBIN d, i**

The MBIN command moves the Marker or Marker Region of display d to bin i. Bin 0 is the left most bin in the display. The marker of each display must be moved separately (d=2 is not allowed).

A Spot marker will move to bin i. A marker region will be centered on bin i (if possible). The Normal marker will still seek the Max, Min or Mean within the region.

This command is not valid if the Marker of display d is Off.

---

**MWFL d, i**

The MWFL command moves the Waterfall Marker of display d to record i (in the Z axis). Record 0 is the most recent record in the display. The marker of each display must be moved separately (d=2 is not allowed).

This command is valid only if display d is a Waterfall and the measurement is paused. This command is not valid if the Marker of display d is Off.

---

**MRKX ? d**

The MRKX ? command queries the time, frequency or bin number of the marker position in display d. The returned value is always the absolute position of the marker, even when the on-screen marker is relative. The returned value has the same units as the first value shown in the Marker Position Bar.

This command is not valid if the Marker of display d is Off.

---

**MRKB ? d**

The MRKB ? command queries the marker bin number of display d. The value returned is the bin number of the marker. Bin 0 is the left most bin in the display.

---

This command is not valid if the Marker of display d is Off.

---

### **MRKY ? d**

The MRKY ? command queries the vertical axis data value at the marker position of display d. The returned value is always the absolute position of the marker, even when the on-screen marker is relative.

This command is not valid if the Marker of display d is Off.

---

### **MRKZ ? d**

The MRKZ ? command queries the Z position at the marker position of display d.

If the display is a 2-D display (Nyquist or Nichols view), MRKZ? queries the horizontal axis data value at the marker position. The returned value is always the absolute position of the marker, even when the on-screen marker is relative. MRKY? queries the position along the vertical axis and MRKX? queries the time or frequency.

If the display is a waterfall, MRKZ? queries the Z-axis record number of the marker. Lower record numbers are more recent.

This command is not valid if the Marker of display d is Off.

---

### **MKMX d**

The MKMX command performs Marker to Max on display d. (Same as [Marker Max] key). The marker of each display must be moved separately (d=2 is not allowed).

This command is not valid if the Marker of display d is Off.

---

### **MKMN d**

The MKMN command performs Marker to Min on display d. (Same as [Marker Min] key). The marker of each display must be moved separately (d=2 is not allowed).

This command is not valid if the Marker of display d is Off.

---

### **MKCN d**

The MKCN command performs Marker to Center on display d. (Same as [Marker Center] key). The center frequency of the FFT span is set to the marker frequency on display d. The span is decreased if necessary. The center of each display must be moved separately (d=2 is not allowed).

This command is only valid when the Measurement Group is FFT and display d is Live. This command is not valid if the Marker of display d is Off.

---

## Marker Commands (Normal)

---

**MSEK (?) d {, i}**

The MSEK command sets (queries) the Normal Marker Seeks Mode of display d. The parameter i selects Max (0), Min (1) or Mean (2).

This command is only valid if the Marker Mode for display d is Normal.

---

**MWID (?) d {, i}**

The MWID command sets (queries) the Normal Marker Width of display d. The parameter i selects Spot (0), 1/2 division (1) or 1 division (2).

This command is only valid if the Marker Mode for display d is Normal.

---

**MREL (?) d {, i}**

The MREL command sets (queries) the Normal Marker Relative Mode of display d. The parameter i selects Off (0), Relative to Offset (1), Relative to Reference Display (2) or Relative to Other Display (3).

This command is only valid if the Marker Mode for display d is Normal.

---

**MROX (?) d {, x}**

The MROX command sets (queries) the X Offset for the Normal Marker of display d. The parameter x is a real number in display units.

This command is only valid if the Marker Mode for display d is Normal.

---

**MROY (?) d {, x}**

The MROY command sets (queries) the Y Offset for the Normal Marker of display d. The parameter x is a real number in display units.

This command is only valid if the Marker Mode for display d is Normal.

---

**MXRL (?) d {, i}**

The MXRL command sets (queries) the X Rel Mode for the Normal Marker of display d. The parameter i selects Absolute (Off) or Relative (On).

This command is only valid if the Marker Mode for display d is Normal.

---

**MRON (?) d**

The MRON command sets toggles the Normal Marker between Marker Rel Off and Marker Rel to Offset. This is similar to the [Marker Ref] key. MRON operates on a single display (d=2 is not allowed).

If Marker Rel is Off, MRON d sets the Normal Marker offsets (X and Y) to the current marker position and sets the Marker to Relative to Offset (relative marker readings).

---

## 5-48 Marker Commands

---

If Marker Rel is Relative to Offset, MRON d sets the Marker Rel to Off (absolute marker readings).

MRON? d returns 0 if Marker Rel is Off and 1 if Marker is Rel to Offset.

This command is only valid if the Marker Mode for display d is Normal.



---

## Marker Commands (Harmonic)

---

**HRMN (?) d {, i}**

The HRMN command sets (queries) the Number of Harmonics for display d. The parameter i is a number of harmonics from 1 to 100.

This command is only valid if the Marker Mode for display d is Harmonic.

---

**HDSP (?) d {, i}**

The HDSP command sets (queries) the Harmonic and Sideband Marker Display for display d.

For Harmonic Marker, i=0 selects the Fundamental, i=1 selects 2xFundamental, etc. The maximum value of i is the Number of Harmonics.

For Sideband Marker, i=0 selects the Fundamental, i=-1 selects the first lower sideband, i=1 selects the first upper sideband. The maximum value of i plus or minus the Number of Sidebands. i<0 is valid only for Sideband marker.

This command is only valid if the Marker Mode for display d is Harmonic or Sideband.

---

**HRDO (?) d {, i}**

The HRDO command sets (queries) the Harmonic and Sideband Marker Readout for display d. The parameter i selects Absolute (0) or Relative to Fundamental (1).

This command is only valid if the Marker Mode for display d is Harmonic or Sideband.

---

**HTHD ? d, i**

The HTHD ? command queries the Total Harmonic Distortion for display d. The parameter i selects percent (0) or dB Relative to Fundamental (1) and is required.

This command is only valid if the Marker Mode for display d is Harmonic.

---

**HPWR ? d**

The HPWR ? command queries the Harmonic Power for display d. The returned value is a real value of Vrms or dBVrms depending on the setting of dB Units>.

This command is only valid if the Marker Mode for display d is Harmonic.

---

## Marker Commands (Sideband)

---

**SSEP (?) d {, f} <kHz, [Hz], mHz,  $\mu$ Hz>**

The SSEP command sets (queries) the Sideband Separation for display d. The parameter f is the sideband separation in Hz.

This command is only valid if the Marker Mode for display d is Sideband.

---

**SDBN (?) d {, i}**

The SDBN command sets (queries) the Number of Sidebands for display d. The parameter i is a number of sidebands from 1 to 50.

This command is only valid if the Marker Mode for display d is Sideband.

---

**HDSP (?) d {, i}**

The HDSP command sets (queries) the Harmonic and Sideband Marker Display for display d.

For Sideband Marker, i=0 selects the Fundamental, i=-1 selects the first lower sideband, i=1 selects the first upper sideband. The maximum value of i plus or minus the Number of Sidebands. i<0 is valid only for Sideband marker.

This command is only valid if the Marker Mode for display d is Harmonic or Sideband.

---

**HRDO (?) d {, i}**

The HRDO command sets (queries) the Harmonic and Sideband Marker Readout for display d. The parameter i selects Absolute (0) or Relative to Fundamental (1).

This command is only valid if the Marker Mode for display d is Harmonic or Sideband.

---

**SPWR ? d, i**

The SPWR ? command queries the Sideband Power for display d. The parameter i selects dB Relative to Fundamental (0) or Vrms or dBVrms (1) and is required.

This command is only valid if the Marker Mode for display d is Sideband.

---

---

## Marker Commands (Band)

---

**BMKL (?) d, i {, j}**

The BMKL command sets (queries) the Band Left Bin for display d. The parameter i selects Upper (0) or Lower (1) band and j is the bin number (0 is the most left) and cannot exceed number of bins in the display.

This command is only valid if the Marker Mode for display d is Band.

---

**BMKR (?) d, i {, j}**

The BMKR command sets (queries) the Band Right Bin for display d. The parameter i selects Upper (0) or Lower (1) band and j is the bin number (0 is the most left) and cannot exceed number of bins in the display.

This command is only valid if the Marker Mode for display d is Band.

---

**BEXC ? d {, i}**

The BEXC command queries the Band Exclusion for display d. The parameter i selects none (0), 1from2 (1) or 2from1 (2).

This command is only valid if the Marker Mode for display d is Band.

---

**BRAM ? d {, i}**

The BRAM command queries the Band Ratio Mode for display d. The parameter i selects 2/1 (0) or 2/(1+2) (1).

This command is only valid if the Marker Mode for display d is Band.

---

**BPWR ? d, i**

The BPWR command queries the Band Power for display d. The parameter i selects Upper (0) or Lower (1) band power.

This command is only valid if the Marker Mode for display d is Band.

---

**BRAT ? d, i**

The BRAT command queries the Band Ratio for display d. The parameter i selects percent (0) or dB (1) band ratio.

This command is only valid if the Marker Mode for display d is Band.

---

## Marker Commands (Frequency/Damping)

---

### MMCA ? d

The MMCA command queries the Frequency/Damping marker result. The query returns two values in the form f, g where f is the resonant frequency and g is the damping factor.

This command is only valid if the Marker Mode for display d is Frequency/Damping.

## Source Commands

---

**SRCO (?) {i}**

The SRCO command sets (queries) the Source On or Off. The parameter *i* selects Off (0) or On (1).

This command is not valid when the Measurement Group is Swept Sine.

---

**STYP (?) {i}**

The STYP command sets (queries) the Source Type. The parameter *i* selects Sine (0), Chirp (1), Noise (2) or Arbitrary (3).

When the Measurement Group is Swept Sine, the Source Type may not be changed and this command is not valid.

## Sine Source Commands

---

**S1FR (?) {f} <kHz, [Hz], mHz,  $\mu$ Hz>**

The S1FR command sets (queries) the Frequency of Sine Tone 1. The parameter f is a frequency in the specified units. The query command always returns a value in Hz.

This command is valid only when the Source Type is Sine.

---

**S1AM (?) {x} <mVpk, mVpp, mVrms, [Vpk], Vrms, dBVpk, dBVpp, dBVrms>**

The S1AM command sets (queries) the Amplitude of Sine Tone 1. The parameter x is a real number of the specified units. The query returns two numbers of the form y,i where y is a real number and i is an index indicating the units.

This command is valid only when the Source Type is Sine.

---

**SOFF (?) {x} <mV, [V]>**

The SOFF command sets (queries) the Offset of the Sine Source. The parameter x is a real number of the specified units. The query returns two numbers of the form y,i where y is a real number and i is an index indicating the units.

This command is valid only when the Source Type is Sine.

---

**S2FR (?) {f} <kHz, [Hz], mHz,  $\mu$ Hz>**

The S2FR command sets (queries) the Frequency of Sine Tone 2. The parameter f is a frequency in the specified units. The query command always returns a value in Hz.

This command is valid only when the Source Type is Sine.

---

**S2AM (?) {x} <mVpk, mVpp, mVrms, [Vpk], Vrms, dBVpk, dBVpp, dBVrms>**

The S2AM command sets (queries) the Amplitude of Sine Tone 2. The parameter x is the amplitude in the specified units. The query returns two numbers of the form y,i where y is a real number and i is an index indicating the units.

This command is valid only when the Source Type is Sine.

---

## Chirp Source Commands

---

**CAMP (?) {x} <mV, [V], dBVpk>**

The CAMP command sets (queries) the Chirp Amplitude. The parameter x is the amplitude in the specified units. The query returns two numbers of the form y,i where y is a real number and i is an index indicating the units.

This command is valid only when the Source Type is Chirp.

---

**CBUR (?) {x}**

The CBUR command sets (queries) the Chirp Burst Percentage. The parameter x is a percentage from 1 to 100.

This command is valid only when the Source Type is Chirp and the Measurement Group is FFT or Correlation.

---

**CSRC (?) {i}**

The CSRC command sets (queries) the Source Display. The parameter i selects DisplayA (0) or DisplayB (1)

This command is valid only when the Source Type is Chirp or Noise and the Measurement Group is FFT or Correlation.

---

## Noise Source Commands

---

**NAMP (?) {x} <mV, [V], dBVpk>**

The NAMP command sets (queries) the Noise Amplitude. The parameter x is the amplitude in the specified units. The query returns two numbers of the form y,i where y is a real number and i is an index indicating the units.

This command is valid only when the Source Type is Noise.

---

**NTYP (?) {i}**

The NTYP command sets (queries) the Noise Type. The parameter i selects Bandlimited White (0), White (1) or Pink (2).

This command is valid only when the Source Type is Noise.

---

**NBUR (?) {x}**

The NBUR command sets (queries) the Noise Burst Percentage. The parameter x is a percentage from 1 to 100.

This command is valid only when the Source Type is Noise and the Measurement Group is FFT or Time/Histogram.

---

**NPER (?) {x} <ks, [s], ms>**

The NPER command sets (queries) the Noise Burst Source Period. The parameter x is a time from 4 ms to 1 ks.

This command is valid only when the Source Type is Noise and the Measurement Group is Octave or Order.

---

**CSRC (?) {i}**

The CSRC command sets (queries) the Source Display. The parameter i selects DisplayA (0) or DisplayB (1)

This command is valid only when the Source Type is Chirp or Noise and the Measurement Group is FFT.

---



---

# Arbitrary Source Commands

---

**AAMP (?) {x}**

The AAMP command sets (queries) the Arbitrary Source Amplitude. The parameter x is a percentage of full scale (1V).

This command is valid only when the Source Type is Arbitrary.

---

**ARAT (?) {i}**

The ARAT command sets (queries) the Arbitrary Source Playback Rate. The parameter i [0 .. 19] selects a rate equal to the maximum rate divided by 2 to the ith power. Thus i=0 selects the maximum rate and i=1 selects half the maximum rate.

This command is valid only when the Source Type is Arbitrary.

---

**ASRC (?) {i}**

The ASRC command sets (queries) the Arbitrary Source Buffer. The parameter i selects Arbitrary Waveform memory (0), Ch1 Capture (1) or Ch2 Capture (2). The selected buffer must contain data otherwise an error occurs.

This command is valid only when the Source Type is Arbitrary.

---

**ASTR (?) {i}**

The ASTR command sets (queries) the Arbitrary Source Start point. The parameter i is a point in the Arbitrary Waveform (either in Arbitrary memory or a Capture buffer).

This command is valid only when the Source Type is Arbitrary.

---

**ALEN (?) {i}**

The ALEN command sets (queries) the Arbitrary Source Length. The parameter i is the length in kPts in the Arbitrary Waveform (either in Arbitrary memory or a Capture buffer). i must be an even number from 2 to the length of the waveform.

This command is valid only when the Source Type is Arbitrary.

---

**TARB i**

The TARB command copies Trace i to the Arbitrary Waveform buffer. The Arbitrary Length is changed to 2 kPoints and the Arbitrary Source is change to Arb. Buffer. Trace i must contain FFT measurement data (usually a time record).

This command is valid only when the Source Type is Arbitrary.

---

## Swept Sine Source Commands

---

### **SSAL (?) {i}**

The SSAL command sets (queries) the Auto Level Reference. The parameter *i* selects Off (0), Channel 1 (1) or Channel 2 (2).

This command is valid only when the Measurement Group is Swept Sine.

---

### **SSAM (?) {x} <mV, [V], dBVpk>**

The SSAM command sets (queries) the Swept Sine Amplitude. The parameter *x* is the amplitude in the specified units. The query returns two numbers of the form *y,i* where *y* is a real number and *i* is an index indicating the units.

This command is valid only when the Measurement Group is Swept Sine and Auto Level Reference is Off.

---

### **SSRF (?) {x} <mV, [V], dBVpk>**

The SSRF command sets (queries) the Swept Sine Ideal Reference. The parameter *x* is the ideal reference in the specified units. The query returns two numbers of the form *y,i* where *y* is a real number and *i* is an index indicating the units.

This command is valid only when the Measurement Group is Swept Sine and Auto Level Reference is Ch1 or Ch2.

---

### **SRMP (?) {i}**

The SRMP command sets (queries) the Swept Sine Source Ramping. The parameter *i* selects Off (0) or On (1).

This command is valid only when the Measurement Group is Swept Sine.

---

### **SRAT (?) {x}**

The SRAT command sets (queries) the Swept Sine Source Ramp Rate. The parameter *x* is the ramp rate in V/s.

This command is valid only when the Measurement Group is Swept Sine.

---

### **SSUL (?) {x}**

The SSUL command sets (queries) the Swept Sine Reference Upper Limit. The parameter *x* is a ratio in dB from 0.1 to 30.0 dB.

This command is valid only when the Measurement Group is Swept Sine and Auto Level Reference is Ch1 or Ch2.

---

**SSL (*x*)**

The SSL command sets (queries) the Swept Sine Reference Lower Limit. The parameter *x* is a ratio in dB from -0.1 to -30.0 dB.

This command is valid only when the Measurement Group is Swept Sine and Auto Level Reference is Ch1 or Ch2.

---

**SMA (*x*)** <mV, [V], dBVpk>

The SMA command sets (queries) the Swept Sine Maximum Level. The parameter *x* is the maximum level in the specified units. The query returns two numbers of the form *y,i* where *y* is a real number and *i* is an index indicating the units.

This command is valid only when the Measurement Group is Swept Sine and Auto Level Reference is Ch1 or Ch2.

---

**SOF (*x*)** <mV, [V]>

The SOF command sets (queries) the Offset of the Swept Sine Source. The parameter *x* is a real number of the specified units. The query returns two numbers of the form *y,i* where *y* is a real number and *i* is an index indicating the units.

This command is valid in the swept sine measurement group.

# Input Commands

---

**ISRC (?) {i}**

The ISRC command sets (queries) the Input Source. The parameter *i* selects Analog (0) or Capture (1).

Capture is not valid when the Measurement Group is Swept Sine.

---

**LINK (?) {i}**

The LINK command sets (queries) the Analyzer Configuration. The parameter *i* selects Independent Channels (0) or Dual Channels (1).

---

**I1MD (?) {i}**

The I1MD command sets (queries) the Ch1 Input Mode. The parameter *i* selects A (single-ended) (0) or A-B (differential) (1).

---

**I1GD (?) {i}**

The I1GD command sets (queries) the Ch1 Input Grounding. The parameter *i* selects Float (0) or Ground (1).

---

**I1CP (?) {i}**

The I1CP command sets (queries) the Ch1 Input Coupling. The parameter *i* selects DC (0), AC (1) or ICP (2).

---

**I1RG (?) {x} <[dBVpk], dBVpp, dBVrms, Vpk, Vpp, Vrms, dBEUpk, dBEUpp, dBEUrms, EUpk, EUpp,EUrms>**

The I1RG command sets (queries) the Ch1 Input Range. The parameter *x* is the full scale input range in the specified units. If the exact input range specified is not available, the closest available input range will be used.

The query command returns a value of the form *x,j* where *x* is the value and *j* is the index into the list of units.

If Ch1 AutoRange is On, the I1RG *i* command will turn Ch1 AutoRange Off and the set the Ch1 Range to *f*.

---

**A1RG (?) {i}**

The A1RG command sets (queries) the Ch1 AutoRanging Off/On. The parameter *i* selects Off (Manual) (0) or On (AutoRanging) (1). If *i*=1 and Ch1 AutoRange is already On, a new AutoRange is performed.

---

**I1AR (?) {i}**

The I1AR command sets (queries) the Ch1 AutoRange Mode. The parameter *i* selects Up Only (0) or Tracking (1).

---

**I1AF (?) {i}**

The I1AF command sets (queries) the Ch1 Anti-Aliasing Filter Off/On. The parameter *i* selects Off (0) or On (1).

---

---

<b>I1AW (?) {i}</b>	The I1AW command sets (queries) the Ch1 A-Weighting Filter Off/On. The parameter i selects Off (0) or On (1).
<b>I2MD (?) {i}</b>	The I2MD command sets (queries) the Ch2 Input Mode. The parameter i selects A (single-ended) (0) or A-B (differential) (1).
<b>I2GD (?) {i}</b>	The I2GD command sets (queries) the Ch2 Input Grounding. The parameter i selects Float (0) or Ground (1).
<b>I2CP (?) {i}</b>	The I2CP command sets (queries) the Ch2 Input Coupling. The parameter i selects DC (0), AC (1) or ICP (2).
<b>I2RG (?) {x} &lt;[dBVpk], dBVpp, dBVrms, Vpk, Vpp, Vrms, dBEUpk, dBEUpp, dBEUrms, EUpk, EUpp,EUrms&gt;</b>	<p>The I2RG command sets (queries) the Ch2 Input Range. The parameter x is the full scale input range in the specified units. If the exact input range specified is not available, the closest available input range will be used.</p> <p>The query command returns a value of the form x,j where x is the value and j is the index into the list of units.</p> <p>If Ch2 AutoRange is On, the I2RG i command will turn Ch2 AutoRange Off and the set the Ch2 Range to f.</p>
<b>I2AR (?) {i}</b>	The I2AR command sets (queries) the Ch2 AutoRange Mode. The parameter i selects Normal (0) or Tracking (1).
<b>A2RG (?) {i}</b>	The A2RG command sets (queries) the Ch2 AutoRanging Off/On. The parameter i selects Off (Manual) (0) or On (AutoRanging) (1). If i=1 and Ch2 AutoRange is already On, a new AutoRange is performed.
<b>I2AF (?) {i}</b>	The I2AF command sets (queries) the Ch2 Anti-Aliasing Filter Off/On. The parameter i selects Off (0) or On (1).
<b>I2AW (?) {i}</b>	The I2AW command sets (queries) the Ch2 A-Weighting Filter Off/On. The parameter i selects Off (0) or On (1).
<b>IAOM (?) {i}</b>	The IAOM command sets (queries) the Input Auto Offset Off/On. The parameter i selects Off (0) or On (1).

---

## Transducer Parameter Commands

---

**EU1M (?) d {,i}**

The EU1M command sets (queries) the Ch1 Engineering Units Mode for display d. The parameter i selects Off (0) or On (1).

The set command requires d = 2 (both displays).

---

**EU1L (?) d {,i}**

The EU1L command sets (queries) the Ch1 Engineering Units Label for display d. The set command requires d = 2 (both displays).

The parameter i selects the unit label.

<u>i</u>	<u>Label</u>	<u>i</u>	<u>Label</u>
0	m/s <sup>2</sup>	8	kg
1	m/s	9	lbs
2	m	10	N
3	in/s <sup>2</sup>	11	dyne
4	in/s	12	Pas
5	in	13	bar
6	mil	14	USER
7	g		

---

**EU1V (?) d {,x}**

The EU1V command sets (queries) the Ch1 Engineering Units per Volt scale for display d. The set command requires d = 2 (both displays). The parameter x is real number of EU/Volt (EU's per Volt).

---

**EU1U (?) d {,s}**

The EU1U command sets (queries) the Ch1 User Label for display d. The set command requires d = 2 (both displays). The string s is the user label.

---

**TD1C (?) d {,i}**

The TD2C command sets (queries) the Ch1 Transducer Conversion. for display d. The set command requires d = 2 (both displays). The parameter i selects the units to be converted to:

<u>i</u>	<u>Label</u>
0	m/s <sup>2</sup>
1	m/s
2	m
3	in/s <sup>2</sup>
4	in/s
5	in
6	mil
7	g

---

**EU2M (?) d {,i}**

The EU2M command sets (queries) the Ch2 Engineering Units Mode for display d. The set command requires d = 2 (both displays).. The parameter i selects Off (0) or On (1).

**EU2L (?) d {,i}**

The EU2L command sets (queries) the Ch2 Engineering Units Label for display d. The set command requires d = 2 (both displays). The parameter i selects the unit label.

<u>i</u>	<u>Label</u>	<u>i</u>	<u>Label</u>
0	m/s <sup>2</sup>	8	kg
1	m/s	9	lbs
2	m	10	N
3	in/s <sup>2</sup>	11	dyne
4	in/s	12	Pas
5	in	13	bar
6	mil	14	USER
7	g		

**EU2V (?) d {,x}**

The EU2V command sets (queries) the Ch2 Engineering Units per Volt scale for display d. The set command requires d = 2 (both displays). The parameter x is real number of EU/Volt (EU's per Volt).

**EU2U (?) d {,s}**

The EU2U command sets (queries) the Ch2 User Label for display d. The set command requires d = 2 (both displays). The string s is the user label.

**TD2C (?) d {,i}**

The TD2C command sets (queries) the Ch2 Transducer Conversion. for display d. The set command requires d = 2 (both displays). The parameter i selects the units to be converted to:

<u>i</u>	<u>Label</u>
0	m/s <sup>2</sup>
1	m/s
2	m
3	in/s <sup>2</sup>
4	in/s
5	in
6	mil
7	g

## Tachometer Input Commands

---

**TAPR (?) {x}**

The TAPR command sets (queries) the number of tach pulses per revolution. The parameter x is the number of pulses per revolution.

---

**TARG (?) {i}**

The TARG command sets (queries) the Tachometer Trigger Range. The parameter i selects 5V (0), 25V(1), or TTL(2).

---

**TALV (?) {x}**

The TALV command sets (queries) the Tachometer Trigger Level. The parameter x is the tach trigger level in volts.

---

**TASL (?) {i}**

The TARG command sets (queries) the Tachometer Trigger Slope. The parameter i selects Rising (0) or Falling (1).

---

**TAHO (?) {i}**

The TAHO command sets (queries) the Tachometer Holdoff Enable. The parameter i selects Off (0) or On (1).

---

**TAHD (?) {x}**

The TAHD command sets (queries) the Tachometer Holdoff Time. The parameter x is the holdoff time in seconds.

---

**TASH (?) {i}**

The TASH command sets (queries) the Show Tach status. The parameter i selects Off (0) or On (1).

---



## Input Playback Commands

---

**ISTR (?) {i}**

The ISTR command sets (queries) the Capture Playback Start point. The parameter *i* is a point in the Capture buffer.

This command is valid only when the Input Source is Playback.

---

**ILEN (?) {i}**

The ILEN command sets (queries) the Capture Playback Length. The parameter *i* is the length in kPts in the Capture buffer. *i* must be an even number from 2 to the length of the waveform.

This command is valid only when the Input Source is Playback.

---

**IMOD (?) {i}**

The IMOD command sets (queries) the Capture Playback Mode. The parameter *i* selects 1-Shot (0) or Circular (1).

This command is valid only when the Input Source is Playback.

---

**ISPD (?) {i}**

The ISPD command sets (queries) the Capture Playback Speed. The parameter *i* selects Normal (0) or Every Frame (1).

This command is valid only when the Input Source is Playback.

---

# Trigger Commands

---

**TMOD (?) {i}**

The TMOD command sets (queries) the Trigger Arming Mode. The parameter i selects , Auto Arm (0), Manual Arm (1), RPM Arm (2), or Time Arm(3).

---

**TSRC (?) {i}**

The TSRC command sets (queries) the Trigger Source. The parameter i selects Continuous (0), Ch1 (1), Ch2 (2), External (3), External TTL (4), Source (5) or Manual (6).

---

**TLVL (?) {x} <[%], V, mV, ChU>**

The TSRC command sets (queries) the Trigger Level. The parameter x is the trigger level in the specified units. The query command returns two values in the form: f,j where f is the trigger level and j is the index into the list of units.

The ChU Unit are equivalent to EUs if EUs are turned on for the active display, otherwise ChU are equivalent to Volts.

---

**TSLP (?) {i}**

The TSRC command sets (queries) the Trigger Slope. The parameter i selects Rising (0) or Falling (1).

---

**TDLA (?) {x} <[s], ms,  $\mu$ s>**

The TDLA command sets (queries) the Trigger Delay 1. The parameter x is the delay in the specified units.

This command is valid only when the Measurement Group is FFT or Correlation.

---

**TDLB (?) {x} <[s], ms,  $\mu$ s>**

The TDLB command sets (queries) the Trigger Delay 2. The parameter x is the delay in the specified units.

This command is valid only when the Measurement Group is FFT or Correlation.

---

**TARM**

The TARM command Manually Arms the trigger. This command may not be queried.

---

**STMD (?) {i}**

The STMD command sets (queries) the Triggered Source Mode. The parameter i selects 1-Shot (0) or Continuous (1).

---

**TMAN**

The TMAN command Manually Triggers if armed. This command may not be queried.

---

**TRSR (?) {x}**

The TRSR command sets (queries) the RPM Arming Start RPM. The parameter x is the start RPM.

---

**TRSM (?) {i}**

The TRSM command sets (queries) the RPM Arming Start RPM mode. The parameter i selects Start RPM Off (0) or On (1).

---

**TRDR (?) {x}**

The TRDR command sets (queries) the RPM Arming Delta RPM. The parameter x is delta rpm.

---

**TRDM (?) {i}**

The TRDM command sets (queries) the RPM Arming Delta RPM sense. The parameter i selects Absolute Change (0), Increasing RPM (1), or Decreasing RPM (2)..

---

**TIAS (?) {x}**

The TIAS command sets (queries) the Time Arming Step. The parameter x is the time increment in seconds..

---

## Average Commands (FFT, Correlation, and Order)

---

### **FAVG (?) d {, i}**

The FAVG command sets (queries) the Compute Average Off/On for display d. The parameter i selects Off (0) or On (1).

This command not valid only when the Measurement Group is Octave, Swept Sine, or Time/Histogram. The set command requires display d to be Live.

---

### **FAVM (?) d {, i}**

The FAVM command sets (queries) the type of averaging shown on display d. The parameter i selects None (0), Vector (1), RMS (2) or Peak Hold (3).

This command is valid only when the Measurement Group is FFT, Correlation or Order. The set command requires display d to be Live.

---

### **FAVT (?) d {, i}**

The FAVT command sets (queries) the FFT Averaging Type for display d. The parameter i selects Linear/Fixed Length (0), or Exponential/Continuous (1).

This command is valid only when the Measurement Group is FFT, Correlation, or Order. The set command requires display d to be Live.

---

### **FAVN (?) d {, i}**

The FAVN command sets (queries) the Number of Averages for display d. The parameter i is a number of averages from 2 to 32767.

This command is valid only when the Measurement Group is FFT, Correlation, or Order. The set command requires display d to be Live.

---

### **NAVG ? d**

The NAVG command queries the Number of Averages completed for display d.

For linear averaging, the returned value is less than or equal to the FFT Number of Averages.

For exponential averaging, the returned value is the actual number of averages completed and eventually exceeds the FFT Number of Averages.

This command is valid in the FFT, Correlation, Order, and Octave Measurement Groups.

---

### **FOVL (?) d {, x}**

The FOVL command sets (queries) the FFT Time Record Increment for display d. The parameter x is a percentage up to 300.

This command is valid only when the Measurement Group is FFT. The set command requires display d to be Live.

---

**FREJ (?) d {, i}**

The FREJ command sets (queries) the Overload Reject. The parameter i selects Off (0) or On (1). The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is FFT or Correlation. The set command requires a display to be Live.

---

**TAVM (?) {i}**

The TAVM command sets (queries) the Trigger Average Mode. The parameter i selects Time Records (0) or Averages (1).

This command is valid only when the Measurement Group is FFT, Time/Histogram, or Correlation. The set command requires a d=2 (both displays).

---

**PAVO (?) d {, i}**

The PAVO command sets (queries) the Average Preview. The parameter i selects Off (0), Manual (1) or Timed (2). The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is FFT. The set command requires a display to be Live.

---

**PAVT (?) d {, x}**

The PAVT command sets (queries) the Preview Time. The parameter x is a real number of seconds. The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is FFT. The set command requires a display to be Live.

---

**PAVA**

The PAVA command accepts the previewed time record. The averages of both displays is updated to include the previewed time record.

Poll the NEWA and NEWB (New Data) status bits in the Display Status word to determine when each time record has been acquired. After accepting or rejecting the time record, the display reverts back to showing the actual measurement. This change does not set NEWA or NEWB but occurs upon receipt of the accept or reject command (or after a Preview Time).

This command is valid only when the Measurement Group is FFT. This command has no effect unless the previewed time records are displayed.

---

**PAVR**

The PAVR command rejects the previewed time record. The averages of both displays ignores the previewed time record.

Poll the NEWA and NEWB (New Data) status bits in the Display Status word to determine when each time record has been acquired. After accepting or rejecting the time record, the display reverts back to showing the actual measurement.

---

## 5-70 Average Commands

---

This change does not set NEWA or NEWB but occurs upon receipt of the accept or reject command (or after a Preview Time).

This command is valid only when the Measurement Group is FFT. This command has no effect unless the previewed time records are displayed.

---

## Average Commands (Octave)

---

**OTYP (?) d {, i}**

The OTYP command sets (queries) the Octave Averaging Type for display d. The parameter i selects Linear Time (0), Exponential Time (1), Peak Hold (2) or Equal Confidence (3).

For 1 Channel octave analysis, the set command requires d=2.

This command is valid only when the Measurement Group is Octave. The set command requires display d to be Live.

---

**OTIM (?) d {, x} <ks, [s], ms>**

The OTIM command sets (queries) the Octave Averaging Time for display d. The parameter x is the averaging time from 0.004 to 1000 seconds.

For 1 Channel octave analysis, the set command requires d=2.

This command is valid only when the Measurement Group is Octave. The set command requires display d to be Live.

---

**NAVG ? d**

The NAVG command queries the amount of averaging completed for display d.

The returned value is the amount of averaging completed in 4 ms increments. For example, a completed linear average of 1 second returns 250.

For exponential averaging, the returned value is the actual averaging completed and eventually exceeds the Octave Average Time.

This command is valid in FFT and Octave Measurement Group.

---

**OCNF (?) d {, i}**

The OCNF command sets (queries) the Octave Averaging Confidence Level for display d. The parameter i selects 0.125 dB (0), 0.25 dB (1), 0.5 dB (2), 1.0 dB (3) or 2.0 dB (4).

For 1 Channel octave analysis, the set command requires d=2.

This command is valid only when the Measurement Group is Octave. The set command requires display d to be Live.

---

**OIMP (?) d {, i}**

The OIMP command sets (queries) the Octave Power Bin for display d. The parameter i selects Total (0), Impulse (1), L (2) or Peak (3).

Peak power is not allowed if the Averaging Type is not Peak Hold.

Peak power is automatically selected if the Averaging Type is Peak Hold. Use the OTYP command to set Peak Hold averaging.

---

## 5-72 Average Commands

---

For 1 Channel octave analysis, the set command requires d=2.

This command is valid only when the Measurement Group is Octave. The set command requires display d to be Live.

---

### **OLAT (?) d {, i}**

The OLAT command sets (queries) the Octave Linear Average Mode for display d. The parameter i selects Triggered (0) or Continuous (1).

The set command requires d=2 (both displays have the same mode).

This command is valid only when the Measurement Group is Octave. The set command requires a display to be Live.



## Average Commands (Swept Sine)

---

**SSTM (?) d {, x} <ks, [s], ms>**

The SSTM command sets (queries) the Swept Sine Settle Time for display d. The parameter x is a settle time from 0.008 to 1000 seconds. The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Swept Sine. The set command requires a display to be Live.

---

**SSCY (?) d {, i}**

The SSCY command sets (queries) the Swept Sine Settle Cycles for display d. The parameter i is a number of cycles from 1 to 32767 seconds. The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Swept Sine. The set command requires a display to be Live.

---

**SITM (?) d {, x} <ks, [s], ms>**

The SITM command sets (queries) the Swept Sine Integration Time for display d. The parameter x is an integration time from 0.016 to 1000 seconds. The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Swept Sine. The set command requires a display to be Live.

---

**SICY (?) d {, i}**

The SICY command sets (queries) the Swept Sine Integration Cycles for display d. The parameter i is a number of cycles from 1 to 32767 seconds. The set command requires d=2 (both displays).

This command is valid only when the Measurement Group is Swept Sine. The set command requires a display to be Live.

---

# User Math Commands

---

USRF (?) i {j, k, l, ...}  
USRR (?) i {j, k, l, ...}  
USRO (?) i {j, k, l, ...}  
USRS (?) i {j, k, l, ...}  
USRT (?) i {j, k, l, ...}  
USRH (?) i {j, k, l, ...}

The USRF, USRR, USRO and USRS, USRT and USRH commands define (query) the user functions for the FFT, Correlation, Octave, Swept Sine, and Time/Histogram Measurement Groups.

The parameter i selects a function number from 1 to 5.

The parameters j, k, l, ... are integers which are tokens for the available operands and operations. The tokens are listed below. The letters 'F', 'C', 'O' and 'S', 'Or', and 'H' indicate in which measurement groups each token is valid (FFT, Correlation, Octave, Swept Sine, Order, and Histogram).

For example, to define FFT User Function 3 as FFT(1)+FFT(2), the command USRF 3, 10, -1, 20 is used. If the function cannot be entered, a command execution error occurs.

<u>Operand/Operation</u>	<u>Token</u>	<u>Valid in Groups</u>
Mag(	-101	F,C,O,S,Or,H
Conj(	-102	F,C,O,S,Or,H
Real(	-103	F,C,O,S,Or,H
Imag(	-104	F,C,O,S,Or,H
Ln(	-105	F,C,O,S,Or,H
Exp(	-106	F,C,O,S,Or,H
FFT(	-107	F,C,Or
FFTu(	-108	F,C,Or
IFFT(	-109	F,C,Or
Sqrt(	-110	F,C,O,S,Or
jOmega(	-113	F,C,S
Phase(	-114	F,C,O,S,Or,H
Mag <sup>2</sup> (	-115	F,C,O,S,Or,H
AWt(	-116	F,C,O,S,Or
BWt(	-117	F,C,O,S,Or
CWt(	-118	F,C,O,S,Or
d/dx(	-119	F,C,O,S,Or
GrpDly(	-120	F,C,O,S
[X/(1-X)](	-121	F,C,O,S,Or

+	-1	F,C,O,S,Or,H
-	-2	F,C,O,S,Or,H
x	-3	F,C,O,S,Or,H
/	-5	F,C,O,S,Or,H
(	-6	F,C,O,S,Or,H
)	-7	F,C,O,S,Or,H

FFT Group Operands

Time(1)	20
Time(2)	21
FFT(1)	22
FFT(2)	23
<Pwr(1)>	14
<Pwr(2)>	15
<Coherence>	70
<CrossSpec>	71
<Freq.Resp.>	72
Vec<F1>	24
Vec<F2>	25
RMS<F1>	26
RMS<F2>	27
PeakHold<F1>	28
PeakHold<F2>	29
RMS<CrossSpec>	67

Correlation Group Operands

Time(1)	20
Time(2)	21
Auto_Corr(1)	64
Auto_Corr(2)	65
X_Corr	73

FFT(1)	22
FFT(2)	23
FFTu(1)	30
FFTu(2)	31
Vec<Fu1><F1>	32
Vec<Fu2><F2>	33
Rms<Fu1F1>	34
Rms<Fu2F2>	35
Rms<Fu2><F1>	69
Octave Operands	
Oct(1)	36
Oct(2)	37
Swept Sine Operands	
Spec(1)	38
Spec(2)	39
Norm Var(1)	40
Norm Var(2)	41
Cross Spectrum()	74
Freq. Resp.()	75
Order Operands	
Time(1)	42
Time(2)	43
<Spec(1)>	60
<Spec(2)>	61
<Pwr(1)>	62
<Pwr(2)>	63
Vec<Time(1)>	44
Vec<Time(2)>	45
Time/Histogram Operands	
Histo(1)	48
Histo(2)	49
PDF(1)	50
PDF(2)	51
CDF(1)	52
CDF(2)	53
Time(1)	20
Time(2)	21
Operands for All Groups	
Trace1	1
Trace2	2
Trace3	3
Trace4	4
Trace5	5
Const1	6

Const2	7
Const3	8
Const4	9
Const5	10
j	11

---

**USRC (?) i {, x, y}**

The USRC command sets (queries) the real and imaginary parts of User Constant i. The parameter i selects a User Constant from 1 to 5.

The parameters x and y are floating point values for the real and imaginary parts.

# Window Commands

---

**FWIN (?) d {, i}**

The FWIN command sets (queries) the FFT Window for display d. The parameter i selects a window below.

<u>i</u>	<u>Window</u>
0	Uniform
1	Flattop
2	Hanning
3	BMH
4	Kaiser
5	Force/Exponential
6	User
7	[-T/2..T/2] (Correlation only)
8	[0..T/2] (Correlation only)
9	[-T/4..T/4] (Correlation only)

The Correlation windows may be selected only if display d is a live Correlation measurement.

This command is valid only when the Measurement Group is FFT, Correlation, or Order. The set command requires display d to be Live.

---

**W1FE (?) {i}**

The W1FE command sets (queries) the window which will be applied to Channel 1 when the main window selection is Force/Exponential. If the Analyzer Configuration is set to Independent Channels the command sets (queries) the window which will be applied to the measurement in Display A. The parameter i selects Force Window (0) or Exponential Window (1).

This command is valid only when the Measurement Group is FFT or Order. The set command requires display d to be Live and the Force Exponential Window to be selected.

---

**W2FE (?) {i}**

The W2FE command sets (queries) the window which will be applied to Channel 2 when the main window selection is Force/Exponential. If the Analyzer Configuration is set to Independent Channels the command sets (queries) the window which will be applied to the measurement in Display B. The parameter i selects Force Window (0) or Exponential Window (1).

This command is valid only when the Measurement Group is FFT or Order. The set command requires display d to be Live and the Force Exponential Window to be selected.

---

**FWFL (?) d {, x} <[s], ms,  $\mu$ s>**

The FWFL command sets (queries) the FFT Force Window Length for display d. The parameter x is the duration of the force window in the specified units.

---

This command is valid only when the Measurement Group is FFT or Order. The set command requires display d to be Live and the Force Exponential window to be selected..

---

**FWTC (?) d {, x}**

The FWTC command sets (queries) the FFT Expo Window Time Constant for display d. The parameter x is a percentage of the time record.

This command is valid only when the Measurement Group is FFT or Order. The set command requires display d to be Live and the Force Exponential window to be selected..

---

**TRWI i**

The TRWI command copies the real part of stored Trace i to the User window. Trace i must contain stored data.

This command is valid only when the Measurement Group is FFT or Order.

---

**WITR d, i**

The WITR copies the window of display d to the real part of Trace i. Trace i will have a length of 2048 points.

This command is valid only when the Measurement Group is FFT or Order.

---

**WSYM (?) d {, i}**

The WSYM sets (queries) User Window Form for display d. The parameter i selects Non-Symmetric (0) or Symmetric (1).

This command is valid only when the Measurement Group is FFT or Order and the window of display d is User.

## Waterfall Commands

---

**WDSP (?) d {, i}**

The WDSP command sets (queries) the Waterfall Display Mode for display d. The parameter i selects Normal (0) or Waterfall (1).

This command is not valid when the Measurement Group is Swept Sine.

---

**WSTO (?) d {, i}**

The WSTO command sets (queries) the Waterfall Storage for display d. The parameter i selects Off (0), Continuous (1), or One Shot (2). The set command requires d=2 (both displays).

This command is not valid when the Measurement Group is Swept Sine.

---

**WFSB (?) d {, i}**

The WFSB command sets (queries) the Waterfall Save Option for display d. The parameter i selects Save All Measurements (0) or Save Active Measurement Only (1).

This command is not valid when the Measurement Group is Swept Sine.

---

**WTOT (?) d {, i}**

The WTOT command sets (queries) the Waterfall Total Count for display d. The parameter i is a number of records. If the allocated memory is too small for the requested number of records, an error occurs.

This command is not valid when the Measurement Group is Swept Sine.

---

**WAVA ? d**

The WAVA command queries the number of records available in the waterfall buffer for display d.

Records are not stored in the buffer unless Waterfall Storage is on (Continuous or One Shot).

This command is not valid when the Measurement Group is Swept Sine.

---

**WFSK (?) d {, i}**

The WFSK command sets (queries) the Waterfall Skip count for FFT measurements for display d. The parameter i is a number of records.

This command is valid only when the Measurement Group is FFT, Correlation, Order or Time/Histogram.

---

**WOSK (?) d {, x} <ks, [s], ms>**

The WOSK command sets (queries) the Waterfall Skip count for Octave measurements for display d. The parameter x is a time from .008 to 1000 s.

---



This command is valid only when the Measurement Group is Octave.

---

**WVCT (?) d {, i}**

The WVCT command sets (queries) the Waterfall View Count for display d. The parameter i is a number of records.

This command is not valid when the Measurement Group is Swept Sine.

---

**WHIT (?) d {, i}**

The WHIT command sets (queries) the Waterfall Trace Height for display d. The parameter i is a percentage [20..80] of the total display height.

This command is not valid when the Measurement Group is Swept Sine.

---

**WANG (?) d {, i}**

The WANG command sets (queries) the requested Waterfall Angle for display d. The parameter i is a signed integer number of degrees [-75..+75]. The display will scroll as close to this angle as possible.

This command is not valid when the Measurement Group is Swept Sine.

---

**WFST (?) d {, i}**

The WFST command sets (queries) the Waterfall Fast Angles for display d. The parameter i selects Off (0) or (On).

This command is not valid when the Measurement Group is Swept Sine.

---

**WTHR (?) d {, i}**

The WTHR command sets (queries) the Waterfall Threshold for display d. The parameter i is a percentage of the full scale trace height.

This command is not valid when the Measurement Group is Swept Sine.

---

**WHID (?) d {, i}**

The WHID command sets (queries) the Waterfall Hidden Lines for display d. The parameter i selects Invisible (0) or Visible (1).

This command is not valid when the Measurement Group is Swept Sine.

---

**WREV (?) d {, i}**

The WREV command sets (queries) the Waterfall Paused Drawing for display d. The parameter i selects Normal (Newest at Top) (0) or Reversed (Oldest at Top) (1).

This command is not valid when the Measurement Group is Swept Sine.

---

**WTRC d, i, j**

The WTRC command saves waterfall record j [0 is most recent, 1 is next, etc.] from display d to Trace i [1..5].

## 5-82 Waterfall Commands

---

This command is valid only when the measurement is paused with waterfall storage on. If there is no record *j*, then an error occurs.

---

### **WSLC *d, i, j***

The WSLC command saves the waterfall time slice of bin *j* [0 is left most on x axis] from display *d* to Trace *i* [1..5].

This command is valid only when the measurement is paused with waterfall storage on. If there is no bin *j*, then an error occurs.

---

# Capture Commands

---

**CCHN (?) {i}**

The CCHN command sets (queries) the Capture Channels. The parameter *i* selects Ch1 (0), Ch2 (1) or Ch1+Ch2 (2).

This command is invalid if the source is arbitrary playback from the capture buffer or if the measurement group is Swept Sine.

---

**CMOD (?) {i}**

The CMOD command sets (queries) the Capture Mode. The parameter *i* selects 1-Shot (0) or Continuous (1).

This command is invalid if the source is arbitrary playback from the capture buffer or if the measurement group is Swept Sine.

---

**CLEN (?) {i}**

The CLEN command sets (queries) the Capture Length. The parameter *i* is a length in kPts per channel. *i* must be an even number from 2 to the maximum allocated memory.

This command is invalid if the source is arbitrary playback from the capture buffer or if the measurement group is Swept Sine.

---

**CRAT (?) {i}**

The CRAT command sets (queries) the Capture Rate. The parameter *i* [0 .. 19] selects a rate equal to the maximum rate divided by 2 to the *i*th power. Thus *i*=0 selects the maximum rate and *i*=1 selects half the maximum rate.

This command is invalid if the source is arbitrary playback from the capture buffer or if the measurement group is Swept Sine or Order.

---

**CPAN (?) {i}**

The CPAN command sets (queries) the Capture Auto Pan. The parameter *i* selects Off (0) or On (1).

---

**CSTR**

Same as [Start Capture] key. Capture memory must already be allocated. The Input Source cannot be Playback and the Source cannot be Arbitrary playback from Capture.

---

**CSTP**

Same as [Stop Capture] key. This command has no effect if Capture is not in progress.

---

## Memory Allocation Commands

---

### **MMEM ?**

The MMEM command queries the Total Available Memory for Capture, Waterfall/Order Track and Arbitrary Waveform. MMEM ? returns the total number of memory blocks (2 kPts/block) in the instrument.

---

### **MALC i, j, k**

The MALC command sets the Memory Allocation to *i* blocks for Capture, *j* blocks for Waterfall/Order Track and *k* blocks for Arbitrary Waveform. The sum of  $i+j+k$  cannot exceed the Total Available Memory.

The MALC command automatically confirms the allocation.

---

# Data Table Commands

---

**DTBL (?) d {, i}**

The DTBL command sets (queries) the Data Table for display d. The parameter i selects Off (0) or On (1).

If the Data Table is On for both displays, use the ACTD command to select the active display. Only the Data Table of the active display may be edited or queried.

---

**DMAX (?) d {, i}**

The DMAX command sets (queries) the Data Table Length for display d. The parameter i is a last line number in the table up to 199.

The DMAX ? command queries the last line number in the table.

The DMAX d, i command sets the last line number in the table to i. If i is greater than the current last line number, new entries are created up to line i. These new entries use the current marker position as their X position. If the value of i is less than the current last line number, an error is returned. To remove entries from the table, use the DDLT command.

This command is valid only if the Data Table for display d is On and display d is the active display. Use the ACTD command to select the active display. Only the Data Table of the active display may be edited or queried.

---

**DINS d, i {, j}**

The DINS command sets (queries) the X-Axis bin number for line number i in the Data Table for display d. The parameter i selects the line number. If the value of i is greater than the last line number, an error is returned. The parameter j is the X axis bin number of line i. Bin 0 is the leftmost bin in the display. If there is no bin j in the display, an error occurs.

This command is valid only if the Data Table for display d is On and display d is the active display. Use the ACTD command to select the active display. Only the Data Table of the active display may be edited or queried.

---

**DDLT d, i**

The DDLT command deletes line number i. The value of i may not exceed the last line number in the table. The remaining lines are renumbered sequentially. The Data Table can not be erased completely. The last remaining line may not be deleted.

This command is valid only if the Data Table for display d is On and display d is the active display. Use the ACTD command to select the active display. Only the Data Table of the active display may be edited or queried.

### DCLR d

The DCLR command clears the Data Table for display d. The Data Table can not be erased completely. The last remaining line may not be deleted.

This command is valid only if the Data Table for display d is On and display d is the active display. Use the ACTD command to select the active display. Only the Data Table of the active display may be edited or queried.

---

### DTRD ? d {, i}

The DTRD ? command queries the Data Table measurement results for display d.

The DTRD ? d command queries the entire table. The data is returned in the form X0,Y0,X1,Y1,X2,Y2 ... Xn,Yn[lf] where Xn,Yn are the last entries in the table (index n). Any table entry not within the current measurement span returns the value -1.0E-34 for both X and Y.

The DTRD ? d, i command queries the X and Y values for line i only. The data is returned X, Y [lf]. If the table entry is not within the current measurement span, the values -1.0E-34 are returned for both X and Y.

This command is valid only if the Data Table for display d is On and display d is the active display. Use the ACTD command to select the active display. Only the Data Table of the active display may be edited or queried.

---

## Limit Test Commands

---

**LSON (?) d {, i}**

The LSON command sets (queries) the Limit Segments Display for display d. The parameter i selects Hide (0) or Show (1).

---

**LTST (?) d {, i}**

The LTST command sets (queries) the Limit Testing for display d. The parameter i selects Off (0) or On (1).

---

**LALM (?) d {, i}**

The LALM command sets (queries) the Limit Fail Alarm for display d. The parameter i selects Off (0) or On (1).

---

**LFAL ? d**

The LFAL ? command queries the result of the most recently completed limit test of display d. If the test passed, 0 is returned. If the test failed, 1 is returned.

---

**LCLR d**

The LCLR command clears the Limit Segments for display d. All defined limit segments for display d are erased. The limits for each display are cleared separately (d=2 is not allowed).

This command should always be sent before defining any limit segments. Limit segments may not be defined or edited after the Measurement or Measurement Group is changed until LCLR clears the existing segments.

---

**LMAX (?) d {, i}**

The LMAX command sets (queries) the Last Limit Segment Number for display d. The parameter i is a last limit segment number from 0 to 199. The set command is not valid for d=2 (both displays).

The LMAX ? command queries the last limit segment number.

The LMAX d, i command sets the last limit segment number to i. If i is greater than the current last limit segment number, new entries are created up to line i. These entries simply follow the last limit segment. If the value of i is less than the last line number, an error is returned. To remove entries from the table, use the LDLT command.

This command also sets Limit Segments to Show.

---

**LSEG (?) d, i {, j, x0, y0, x1, y1}**

The LSEG command sets (queries) the Endpoints and Type of limit segment i for display d. The parameter i selects the limit segment number from 0 to the last segment. If i exceeds the last segment number (as set by LMAX), an error is reported.

---

## 5-88 Limit Test Commands

---

The segments are defined separately for each display. The set command is not valid for d=2 (both displays).

The parameter j selects Upper (0) or Lower (1) limit. The parameters x0 and y0 are the coordinates of the limit start (X0,Y0). The parameters x1 and y1 are the coordinates of the limit end (X1,Y1). x0, y0, x1 and y1 are real numbers in the display units.

The LSEG ? d, i command queries the endpoints and type of limit segment i. The data is returned j, f1, f2, y1, y2[lf].

The LSEG d, i, j, x0, y0, x1, y1 command sets the endpoints of segment i to x0, y0 and x1, y1 and the type to Upper (j=0) or Lower (j=1).

This command also sets Limit Segments to Show.

---

### **LDLT d, i**

The LDLT command deletes Limit Segment Number i for display d. The value of i may not exceed the last limit segment number. The remaining segments are renumbered sequentially. This command is not valid for d=2 (both displays).

This command also sets Limit Segments to Show.

---

### **LSFT d, x**

The LSFT command shifts all Limit Segments for display d. The parameter x is real number. The Y coordinates of all segments are shifted by x (in display units). This command is not valid for d=2 (both displays).

This command also sets Limit Segments to Show.



---

# Marker Statistics Commands

---

**MSAO (?) d {, i}**

The MSAO command sets (queries) the Marker Statistics On or Off for display d. The parameter i selects Off (0) or On (1). The set command requires d=2 (both displays).

When Marker Statistics are On, the various statistical quantities are updated whenever new data is available on the display. The accumulation of the statistical quantities continues even when the Marker Statistics menu is not displayed.

---

**MSRS**

The MSRS command resets and starts the accumulation of marker statistics. The various quantities are reset to zero and the accumulation of marker statistics begins.

---

**MSAA ?**

The MSAA ? command queries Max for the Display A Marker.

---

**MSIA ?**

The MSIA ? command queries Min for the Display A Marker.

---

**MSEA ?**

The MSEA ? command queries Mean for the Display A Marker.

---

**MSSA ?**

The MSSA ? command queries Standard Deviation for the Display A Marker.

---

**MSAB ?**

The MSAB ? command queries Max for the Display B Marker.

---

**MSIB ?**

The MSIB ? command queries Min for the Display B Marker.

---

**MSEB ?**

The MSEB ? command queries Mean for the Display B Marker.

---

**MSSB ?**

The MSSB ? command queries Standard Deviation for the Display B Marker.

---

## Exceedance Statistics Commands

---

### **ESTR (?) {i}**

The ESTR command sets (queries) the Exceedance Start Index (for both displays). The parameter *i* is a waterfall index [0..].

This command is not valid when the measurement group is swept-sine.

---

### **ESTP (?) {i}**

The ESTP command sets (queries) the Exceedance Stop Index (for both displays). The parameter *i* is a waterfall index [0..]. Remember, the records are numbered starting from 0. The last record is the number of stored records minus one. The exceedance calculation can not proceed if the stop index is greater than the number of records stored in the waterfall buffer minus one.

This command is not valid when the measurement group is swept-sine.

---

### **EPCT (?) {i}**

The EPCT command sets (queries) the Exceedance Centile (for both displays). The parameter *i* is a centile from 1 to 99.

This command is not valid when the measurement group is swept-sine.

---

### **EXCE d, i**

The EXCE command starts the exceedance centile calculation for display *d*. The parameter *i* selects a Trace [1..5] to store the result.

Display *d* must be paused or done with waterfall storage on and records stored in the waterfall buffer using the active measurement only save option. Waterfall display is not required to be on.

The result is stored in a data trace and has the same measurement type as the waterfall measurements.

This command is not valid when the measurement group is swept-sine.

---

---

# Curve Fit Commands

---

**EFIT d**

The EFIT command starts the curve fit calculation for display d. The fit parameters will be stored in curve table 1.

Display d must be a non-waterfall FFT or Swept-Sine measurement with a frequency x-axis.

---

**ESYN i,d**

The ESYN command synthesizes the parameters in one of the two curve tables and displays the result. The parameter i specifies curve table 1 (0) or curve table 2 (1). The parameter d specifies the display on which the result will appear: Display A (0) or Display B(1).

---

**ENPL (?) {,j}**

The ENPL command sets (queries) the number of poles in the fitting function. The parameter j is the number of poles, between 0 and 20.

---

**ENZE (?) {,j}**

The ENZE command sets (queries) the number of zeros in the fitting function. The parameter j is the number of zeros, between 0 and 20.

---

**EWTT (?) {,j}**

The EWTT command sets (queries) the trace which will be used if user weighting is selected. The parameter j is the trace number (1-5).

---

**EWTT (?) {,j}**

The EWTT command sets (queries) the trace which will be used if user weighting is selected. The parameter j is the trace number (1-5).

---

**ERNG (?) d {,i,j}**

The ERNG command sets (queries) the fit range. The parameter d is the display number. The parameters i and j are the starting and stopping bin numbers for the fit range.

---

**EDLY (?) i {,x}**

The EDLY command sets (queries) the delay value associated with one of the curve tables. The parameter i selects curve table 1 (0) or curve table 2 (1) The parameter x is the delay in seconds.

---

**EFSC (?) i {,x}**

The EFSC command sets (queries) the frequency scaling for one of the curve tables. The parameter i selects curve table 1 (0) or curve table 2 (1) The parameter x is the frequency scaling value.

---

## 5-92 Curve Fit Commands

---

---

### **EGAN (?) i {,x}**

The EGAN command sets (queries) the gain for one of the curve tables. The parameter i selects curve table 1 (0) or curve table 2 (1) The parameter x is the gain value.

---

### **ETRC (?) i {,j}**

The ETRC command sets (queries) the trace associated with one of the curve tables. The parameter i selects curve table 1 (0) or curve table 2 (1) The parameter j is the trace number, from 1 to 5.

---

### **EITM ? i ,j**

The EITM command sets queries the number of items associated with a specific table format for one of the curve tables. The parameter i selects curve table 1 (0) or curve table 2 (1) The parameter j selects the type of item as Zeros (0), Poles (1), Residues (2), Numerator Polynomial (3), or Denominator Polynomial (4). This command will switch the table format to the specified representation before returning the number of items.

---

### **EPLY(?) i, j, k {,x}**

The EPLY command sets (queries) the polynomial coefficients for one of the curve tables. The parameter i selects curve table 1 (0) or curve table 2 (1). The parameter j selects the numerator polynomial (1), or the denominator polynomial (0). The parameter k selects which polynomial coefficient will be set (queried.) from 0 up to the number of items in the list. The parameter x is the actual coefficient.

This command will change the specified table format to polynomial.

---

### **EPOL(?) i, j {,x [,y]}**

The EPOL command sets (queries) the poles for one of the curve tables. The parameter i selects curve table 1 (0) or curve table 2 (1). The parameter j selects which pole will be set (queried.) from 0 up to the number of items in the list. The parameter x is the real part of the pole. The optional parameter y is the imaginary part (if any) of the pole

This command will change the specified table format to pole/zero.

---

### **EZER(?) i, j {,x [,y]}**

The EZER command sets (queries) the zeros for one of the curve tables. The parameter i selects curve table 1 (0) or curve table 2 (1). The parameter j selects which zero will be set (queried.) from 0 up to the number of items in the list. The parameter x is the real part of the zero. The optional parameter y is the imaginary part (if any) of the zero

This command will change the specified table format to pole/zero.

---

### **ERES(?) i, j {,x [,y]}**

The ERES command sets (queries) the residues for one of the curve tables. The parameter i selects curve table 1 (0) or curve table 2 (1). The parameter j selects which residue will be set (queried.) from 0 up to the number of items in the list.

---

The parameter  $x$  is the real part of the residue. The optional parameter  $y$  is the imaginary part (if any) of the residue

This command will change the specified table format to pole/residue.

---

**ECLR  $i$**

The ECLR command clears one of the curve tables. The parameter  $i$  selects curve table 1 (0) or curve table 2 (1)

## Disk Commands

---

### **FNAM (?) {s}**

The FNAM command sets (queries) the Save/Recall File Name. All save and recall disk operations use the name specified by this command. Be sure to use the FNAM s command before any file operation commands.

For example, the “FNAM MYDATA.DAT” command sets the file name to MYDATA.DAT. DOS file name conventions must be followed, i.e. file names are 8 characters or less with an optional extension of up to 3 characters. If the extension is omitted, the SR785 uses default extensions.

File access is to the current directory.

---

### **FDIR (?) {s}**

The FDIR command sets (queries) the Save/Recall File Directory. All save and recall disk operations use the directory specified by this command. Be sure to use the FDIR s command before any file operation commands.

‘FDIR DIRNAME’ looks within the current directory for a sub-directory called DIRNAME. If it exists, the current directory is changed to DIRNAME. FDIR moves down a directory path one level at a time. Use multiple FDIR commands to move down a directory path.

‘FDIR ..’ moves up one directory towards the root. ‘FDIR \’ sets the current directory to the root.

Only the last directory in the path is displayed in the menu box.

---

### **MDIR s**

The MDIR command creates directory s within the current directory.

For example, the “MDIR TODAY” command creates the sub-directory TODAY within the current directory. DOS file conventions must be followed. Use FDIR to set the current directory.

---

### **FXST ? s**

The ‘FXST ? FNAME’ query returns a 1 if the file FNAME exists in the current directory. 0 is returned if there is no file (or subdirectory) named FNAME in the current directory.

Wild cards (\* and ?) are not allowed in the file name. The volume label and sub-directories are treated the same as actual files.

---

### **FREE ?**

The FREE ? query returns the number of bytes available on the disk.

---

**FRST ?**

The FRST ? query returns the first entry in the current directory. The entire entry string is returned, "filename.ext, NNN, m/d/y, h:m:s" where NNN=file size, m/d/y=creation date and h:m:s=creation time.

Use FRST? to reset the catalog file pointer and then repeat FNXT? to query the rest of the directory entries.

**FNXT ?**

The FNXT ? query returns the next entry in the current directory. The entire entry string is returned, "filename.ext, NNN, m/d/y, h:m:s" where NNN=file size, m/d/y=creation date and h:m:s=creation time. FNXT? returns the string "\*" to signify that there are no more entries.

Use FRST? to reset the catalog file pointer and then repeat FNXT? to query the rest of the directory entries until "\*" is returned.

**FSAV d**

The FSAV command saves display d data to disk. The file name is specified by FNAM and the directory is specified by FDIR.

**FRCL d**

The FRCL command recalls data from disk to display d. Display d will be set Off-Line. The file name is specified by FNAM and the directory is specified by FDIR.

**SSAV**

The SSAV command saves the instrument settings to disk. The file name is specified by FNAM and the directory is specified by FDIR.

**SRCL i**

The SRCL command recalls the instrument settings from disk. The individual bits of integer i determine which parameter groups are recalled.

<u>bit</u>	<u>Parameter Group</u>
0	Measurements
1	Sources
2	Analysis
3	Inputs/Triggers
4	DRAM Settings
5	General System
6	Macros

If a bit is 0, the corresponding group IS NOT recalled. If a bit is 1, the corresponding group IS recalled. To determine the value of i, start with i=0. For each bit which is 1, add 2 raised to the bit number. For example, SRCL 127 recalls all groups. SRCL 17 recalls only the Measurements ( $2^0=1$ ) and DRAM Settings ( $2^4=16$ ).

**TSAV i**

The TSAV command saves Trace i data to disk. The file name is specified by FNAME and the directory is specified by FDIR.

---

**TRCL i**

The TRCL command recalls data from disk to Trace i. The file name is specified by FNAME and the directory is specified by FDIR.

---

**DNAM (?) {s}**

The DNAM command sets (queries) the Upkeep File Name. The Delete File function will destroy this file! Be sure to use the DNAM s command before any file delete command.

For example, the “DNAM MYDATA.DAT” command sets the file name to MYDATA.DAT. DOS file name conventions must be followed, i.e. file names are 8 characters or less with an optional extension of up to 3 characters. The extension is NOT automatically supplied.

File access is to the current directory.

---

**DEL F**

The DEL F command deletes the file specified by the Upkeep File Name from the current directory. The file name is specified by DNAM and the directory is specified by FDIR.

---

**DEL D**

The DEL D command deletes the current directory. The directory must be empty otherwise no action is taken. If the directory is deleted, the current directory is changed to the directory one level closer to the root. The current directory is specified by FDIR.

---



# Output Commands

---

## POUT (?) {i}

The POUT command sets (queries) the [Print Screen] hardkey assignment. The parameter *i* selects Bitmap/Print (0), Vector/Plot (1) or ASCII Dump (2).

---

## PDST (?) {i}

The PDST command sets (queries) the Print/Plot/Dump Destination. The parameter *i* selects Disk File (0), Centronics Printer Port (1), RS232 Serial Port (2) or GPIB Port (3).

---

## PFIL (?) {i}

The PFIL command sets (queries) the Disk File Start Number. Print/Plot/Dump to disk automatically numbers the disk file names starting with the File Start Number. The parameter *i* is an integer up to 4 digits.

The PFIL ? command queries the next file number to be used.

---

## DUMP

The DUMP command dumps the data in the active display in ASCII format to the selected Destination (Interface or Disk). The data is written *f,y,x[cr]* with the data for each point on a single line. *f* is the frequency or time value, *y* is the y-axis value and *x* is the x-axis value (if necessary)

---

## PRNT

The PRNT command prints the screen using the selected Bitmap/Printer type and Destination. All other operations are disabled until printing is completed.

---

## PRTP (?) {i}

The PRTP command sets (queries) the Bitmap/Printer type. The parameter *i* selects a type below.

<i>i</i>	<u>Printer Type</u>
0	Epson FX
1	HP PCL
2	HP Small PCL
3	PCX 2 bit (file)
4	GIF (file)
5	PCX 8 bit (file)

---

## PSCR (?) {i}

The PSCR command sets (queries) what portion of the screen to print. The parameter *i* selects Graphs Only (0), Menu Only (1), Status Indicators Only (2) or All (3).

This only affects printing. Plotting only plots the display graphs.

## 5-98 Output Commands

---

---

**PBRI (?) {i}**

The PBRI command sets (queries) the print density of highlighted areas. The parameter *i* selects White (0), 6% (1), 12% (2), 25% (3), 50% (4) or Black (5).

---

**PDIM (?) {i}**

The PDIM command sets (queries) the print density of normal areas. The parameter *i* selects White (0), 6% (1), 12% (2), 25% (3), 50% (4) or Black (5).

---

**PBLK (?) {i}**

The PBLK command sets (queries) the print density of black text. The parameter *i* selects Black (0) or White (1).

---

**PGRF (?) {i}**

The PGRF command sets (queries) the graph print mode. The parameter *i* selects Black on White (0) or White on Black (1).

---

**PLOT**

The PLOT command plots the graph displays using the selected Vector/Plotter type and Destination. All other operations are disabled until plotting is completed.

---

**PLTP (?) {i}**

The PLTP command sets (queries) the Vector/Plotter type. The parameter *i* selects HPGL (0) or PostScript (1).

---

**PLTA (?) {i}**

The PLTA command sets (queries) the Plotter GPIB Address. The parameter *i* ranges from 0 to 30 and should agree with the address of the plotter in use.

---

**PCIC (?) {i}**

The PCIC command sets (queries) the GPIB Control mode. The parameter *i* selects Host (0) or SR785 (1).

---

**PLTX (?) {i}**

The PLTX command sets (queries) the Plotter Text Pen Number. The parameter *i* is a pen number from 1 to 8.

---

**PLGD (?) {i}**

The PLGD command sets (queries) the Plotter Grid Pen Number. The parameter *i* is a pen number from 1 to 8.

---

**PLTR (?) {i}**

The PLTR command sets (queries) the Plotter Trace Pen Number. The parameter *i* is a pen number from 1 to 8.

---

**PLMK (?) {i}**

The PLMK command sets (queries) the Plotter Marker Pen Number. The parameter *i* is a pen number from 1 to 8.

---

**NOTE i, j {, k, l, m, s}**

The NOTE command defines a Display Note. The parameter i selects a Note from 0 to 9. The parameter j makes the note Not Visible (0) or Visible (1). To turn a Note On and Off, use the NOTE i, j command.

The parameter k puts the note in Display A (0) or Display B (1). The parameters l and m are the Text X and Text Y position from 0 to 100. The string s is the note text. The parameters k, l, m and s must always be sent together.

# System Commands

---

**OUTX (?) {i}**

The OUTX command sets (queries) the Output Interface. The parameter i selects GPIB (0) or RS232 (1).

The OUTX i command should always be sent at the start of any program to direct query responses to the correct interface.

---

**OVRM (?) {i}**

The OVRM command sets (queries) the GPIB Override Remote mode. The parameter i selects No (0) or Yes (1). When Override Remote is set to Yes, the front panel is not locked out when the unit is in the REMOTE state.

---

**KCLK (?) {i}**

The KCLK command sets (queries) the Key Click Off/On. The parameter i selects Off (0) or On (1).

---

**ALRM (?) {i}**

The ALRM command sets (queries) the Alarms Off/On. The parameter i selects Off (0) or On (1).

---

**ALRT (?) {i}**

The ALRT command sets (queries) the Alarms Volume. The parameter i selects Quiet (0) or Noisy (1).

---

**ADON (?) {i}**

The ADON command sets (queries) the Done Volume. The parameter i selects Quiet (0) or Noisy (1).

---

**AOVL (?) {i}**

The AOVL command sets (queries) the Audible Overload. The parameter i selects Off (0) or On (1).

---

**SAVR (?) {i}**

The SAVR command sets (queries) the Screen Saver Off/On. The parameter i selects Off (0) or On (1).

---

**SDLY (?) {i}**

The SDLY command sets (queries) the Screen Saver Delay. The parameter i is a delay in minutes from 1 to 59.

---

**FFMT (?) {i}**

The FFMT command sets (queries) the Frequency Format. The parameter i selects Exact Bin (0) or Rounded (1).

---

**TIME (?) {i, j, k}**

The TIME command sets (queries) the System Time.

---

The TIME ? command queries the time. A string of the form “HH:MM:SS” is returned.

The TIME i, j, k command sets the time to i [0..23] hours, j [0..59] minutes and k [0..59] seconds. The parameters i, j and k must all be sent.

---

### **DATE (?) {i, j, k}**

The DATE command sets (queries) the System Date.

The DATE ? command queries the date. A string of the form “MM/DD/YY” is returned.

The DATE i, j, k command sets the date to i [1..12] month, j [1..31] day and k [0..99] year. The parameters i, j and k must all be sent.

## Front Panel Commands

---

### **ACTD (?) i**

The ACTD command sets the Active Display to display i. The parameter i selects DisplayA (0) or DisplayB (1).

Use the ACTD command to select the active display in the Single format. In the Single format, the inactive display is still accessible via commands.

---

### **STRT**

The STRT command starts the measurement. Any average in progress is reset and started over. If the measurement is paused, STRT starts the measurement over.

This function is the same as pressing the [Start/Reset] key.

---

### **PAUS**

If the measurement is already in progress, PAUS pauses the measurement.

If the measurement is paused, PAUS has no effect.

This function is similar to pressing the [Pause/Cont] key.

---

### **CONT**

If the measurement is paused, CONT continues the measurement.

If the measurement is running, CONT has no effect.

This function is similar to pressing the [Pause/Cont] key.

---

### **UNST d**

The UNST command unsettles the measurement of display d. The measurement is not actually perturbed by the UNST command. The settling status is set to unsettled and the full settling time of the measurement is required before the status returns to settled.

The measurement is unsettled by changing any one of several measurement parameters. For example, changing the input range or FFT span will unsettle the measurement. If the signal comes from an external source and is changed in such a way as to require the measurement to settle, it is convenient to use the UNST command and wait for settling to finish.

This command is valid only when the Measurement Group is FFT or Octave.

---

### **SVTR d, i**

The SVTR command saves display d to Trace i.

---

---

<b>RCTR d, i</b>	The RCTR command recalls Trace i to display d. If Trace i does not have data, then an error occurs.
------------------	---

---

<b>SVRF d, i</b>	The SVRF command saves the Reference Display of display d to Trace i.
------------------	---

---

<b>RCRF d, i</b>	The RCRF command recalls Trace i to the Reference Display of display d. If Trace i does not have data or is not compatible with the active display measurement, an error occurs. For example, if Trace i is a time record, it cannot be recalled into the Reference Display of an FFT measurement.
------------------	--

---

<b>PRNT</b>	The PRNT command prints the screen using the selected Printer Type and Destination. All other operations are disabled until printing is completed.
-------------	--

---

<b>PLOT</b>	The PLOT command plots the screen using the selected Plotter Type and Destination. All other operations are disabled until plotting is completed.
-------------	---

---

<b>DUMP</b>	The DUMP command dumps the data in the active display in ASCII format to the selected Destination (Interface or Disk). The data is written f,y,x[cr] with the data for each point on a single line. f is the frequency or time value, y is the y-axis value and x is the x-axis value (if necessary)
-------------	--

---

<b>ASCL d</b>	The ASCL command AutoScales display d. This is the same as pressing the [Auto Scale A] or [Auto Scale B] keys. The displays are AutoScaled separately (d=2 is not allowed).
---------------	---

---

<b>A1RG (?) {i}</b>	The A1RG command sets (queries) the Ch1 AutoRanging Off/On. The parameter i selects Off (Manual) (0) or On (AutoRanging) (1). If i=1 and Ch1 AutoRange is already On, a new AutoRange is performed.
---------------------	---

---

<b>A2RG (?) {i}</b>	The A2RG command sets (queries) the Ch2 AutoRanging Off/On. The parameter i selects Off (Manual) (0) or On (AutoRanging) (1). If i=1 and Ch2 AutoRange is already On, a new AutoRange is performed.
---------------------	---

---

<b>MRON d</b>	The MRON command sets the Normal Marker offsets (X and Y) to the current marker position and sets the Marker Relative to Offset (if Marker Rel was Off). This is similar to the [Marker Ref] key.  This command is only valid if the Marker Mode for display d is Normal.
---------------	---

---

## 5-104 Front Panel Commands

---

---

### **MKMX d**

The MKMX command performs Marker to Max on display d. (Same as [Marker Max] key).

This command is not valid if the Marker of display d is Off.

---

### **MKMN d**

The MKMN command performs Marker to Min on display d. (Same as [Marker Min] key).

This command is not valid if the Marker of display d is Off.

---

### **MKCN d**

The MKCN command performs Marker to Center on display d. (Same as [Marker Center] key). The center frequency of the FFT span is set to the marker frequency on display d. The span is decreased if necessary.

This command is only valid when the Measurement Group is FFT and display d is Live. This command is not valid if the Marker of display d is Off.

---

### **DREF d, i**

The DREF command sets and clears the Display Reference for display d. If i is 0, the Display Reference is turned Off. If i is 1, the current data becomes the Display Reference and the Display Reference is turned On. DREF d,1 only sets the Display Reference if the Display Reference is currently off.

This is similar to the [Display Ref] key.

---

### **SNAP d**

The SNAP command rescales the Reference Display of display d to the current display scale. (Same as [Snap Ref] key). This command has no effect if there is no Reference Display in display d.

---

### **SBRI (?) {i}**

The SBRI command sets (queries) the Screen Brightness. The parameter i is a brightness level from 150 (dimmiest) to 255 (brightest).

---

### **SCON (?) {i}**

The SCON command sets (queries) the Screen Contrast. The parameter i is a contrast level from 0 (no contrast) to 90 (most contrast).

---



**KEYP i**

The KEYP command performs the same function as pressing key i in the table below.

<u>key</u>	<u>i</u>
<Softkey 1> (top)	0
<Softkey 2>	12
<Softkey 3>	11
<Softkey 4>	10
<Softkey 5>	9
<Softkey 6>	8
<Softkey 7>	16
<Softkey 8>	24
<Softkey 9>	32
<Softkey 10> (bottom)	40
[Freq]	33
[Display Setup]	34
[Display Options]	35
[Marker]	36
[Source]	41
[Input]	42
[Trigger]	43
[Average]	44
[User Math]	49
[Window]	50
[Waterfall]	51
[Capture]	52
[Analysis]	57
[Disk]	58
[Output]	59
[System]	60
[AutoScale A]	37
[AutoScale B]	45
[Span Up]	53
[Span Down]	61
[AutoRange Ch1]	38
[AutoRange Ch2]	46
[Marker Max]	54
[Marker Min]	62
[Marker Ref]	39
[Display Ref]	47
[Marker Center]	55
[Show Setup]	63
[Start/Reset]	17
[Pause/Cont]	25
[Start Capture]	18
[Stop Capture]	26
[Active Display]	19
[Link]	27
[Print Screen]	20

## 5-106 Front Panel Commands

---

[Help/Local]	28
[Alt]	7
[Backspace]	6
[Exp]	5
[0]	29
[1]	21
[2]	22
[3]	23
[4]	13
[5]	14
[6]	15
[7]	2
[8]	3
[9]	4
[.]	30
[-]	31
[Enter]	1
[Brighter]	56
[Dimmer]	48

---

### **KNOB i**

The KNOB command simulates turning the knob. The parameter *i* selects counter-clockwise (0) or clockwise (1).

---

### **TONE i, j**

The TONE command makes an audible tone. The parameter *i* is the duration (in 5 ms increments) and *j* selects a note from 0 to 66.

Sequential TONE commands need to be separated by a pause which is at least as long as the preceding tone.

---

### **PLAY i**

The PLAY command plays one of the SR785's pre-programmed sounds. The parameter *i* selects a sound from 0 to 6 (most recent TONE).

For example, PLAY 2 sounds an alert.

---

# Data Transfer Commands

---

**DSPN ? d**

The DSPN ? command queries the length of display d. The parameter d selects Display A (0) or Display B (1).

The returned value is the length (number of points or bins) of display d. The points are numbered from 0 to (length-1).

Always use the DSPN ? command to determine the display length before reading all of the data from a display.

---

**DSPY ? d {, j}**

The DSPY ? command queries the data in display d. The parameter d selects Display A (0) or Display B (1).

The DSPY ? d, j command queries the data value of bin j only. The bins are numbered from 0 to length-1. The value is returned as an ASCII real number.

The DSPY ? d command queries all of the data in display d. Do not serial poll for IFC (bit 7) in the Serial Poll status after this command is sent. IFC will NOT be set until the transfer is complete. Send the command and then make the host computer a listener and the SR785 a talker. Data is returned continuously starting with bin 0 and ending with bin #(length-1). The data points are ASCII real numbers, separated by commas, and the last data point is followed by a terminator. If the DSPY? d command is used with the RS232 interface, the host computer interface must be fast enough to keep up with the data stream.

The DSPY ? command returns data from the last available display update. To keep the data static, pause the measurement or take the display Off-Line.

The returned data depends upon the display View and Units. The data values are the same as if they were read with the marker. In the case of a 2-D view (Nyquist or Nichols), two values per bin are returned (separated by a comma). The values are in the same order as shown in the Marker Position Bar. In this case, there are twice as many points returned.

To read the data in a Trace, recall the Trace to a display and then read the display.

When the display d is a waterfall display, this command queries the value for the trace corresponding to the waterfall marker.

---

**REFY ? d, j**

The REFY ? command queries the data in Reference Graph of display d. The parameter d selects Display A (0) or Display B (1). An error occurs if the display has no Reference Graph.

## 5-108 Data Transfer Commands

---

The REFY ? d, j command queries the Reference Graph data value of bin j only. The bins are numbered from 0 to length-1. The value is returned as an ASCII real number.

The returned data depends upon the display View and Units. The data values are the same as if they were read with a marker.

---

### DSPB ? d {, j}

The DSPB? d command returns the data in display d in binary format. The parameter d selects Display A (0) or Display B (1). This command is only available with the GPIB interface.

The DSPB ? d, j command returns the data value of bin j only. The bins are numbered from 0 to length-1. The value is returned as a 4 byte IEEE float.

The DSPB ? d command returns all of the data in display d in binary format.

Do not serial poll for IFC (bit 7) in the Serial Poll status after DSPB ? is sent. IFC will NOT be set until the transfer is complete. Send the command and then make the host computer a listener and the SR785 a talker.

The DSPB? d returns 4 bytes per bin starting with bin 0 and continuing to bin #(length-1). There is no separation between data points. No line-feed follows the last data point, instead, EOI is asserted with the last byte. The 4 bytes are an IEEE floating point number.

The returned data depends upon the display View and Units. The data values are the same as if they were read with the marker. In the case of a 2-D view (Nyquist or Nichols), two values per bin are returned (8 bytes per bin). The first value is the Y axis value and the second is the X axis value. In this case, there are twice as many points returned.

The host interface must be capable of binary transfer, i.e. accepting line feeds and carriage returns as data rather than terminators. In addition, the host program must read exactly the correct number of bytes (4xlength or 8xlength). While a binary transfer is in progress, the analyzer will not respond to any other queries and the displays will not update.

If the host program does not start reading the points within 1 second, or pauses for 1 second while reading, the binary transfer will be aborted.

To read the data in a Trace, recall the Trace to a display and then read the display.

---

### DBIN ? d, j

The DBIN? command queries the frequency or time of bin j in display d. The parameter d selects Display A (0) or Display B (1). The parameter j is bin from 0 to length-1.

The returned value is either a frequency (spectra), a time (time record or correlation) or a band center frequency (octave).

---

**FBIN ? d, x**

The FBIN ? command queries the bin number of x (frequency or time) in display d. The parameter d selects Display A (0) or Display B (1). The parameter x is a floating point value of Hz or seconds.

The returned value is an integer bin number.

---

**TLOD ? i, n**

The TLOD command loads binary data into Trace i. The parameter i selects Trace 1 through 5. This command is valid only via the GPIB interface.

The target Trace i must already store data of the measurement type and length to be loaded from the interface. The loaded data simply takes the place of the stored trace data and assumes the measurement parameters (Measurement, Averaging, Window, Frequency span, etc.) and length of Trace i.

TLOD loads the actual complex data values in the trace, NOT the data displayed in a particular view. Remember, views are simply different ways of showing the actual complex data.

A trace contains n complex points. The value of n depends upon the type and length of the target trace. N should be chosen equal to or less than the length of the target trace. If the length is less than that of the target trace the remainder of the trace will be filled with zeros.

The download sequence is as follows:

**Host**

Send TLOD ? i, n. Do NOT wait for IFC to be set in the Serial Poll status.

**SR785**

Checks to make sure that Trace i has length n. Returns 1 (4-byte binary long int) if OK to begin data transfer. A return of 0 indicates that n is too large for Trace i.

**Host**

On receipt of 1 (4-byte binary long int), executes a binary transfer to the SR785 of 2n 4-byte IEEE floats. The order is real part of point 0, imag part of point 0, real part of point 1, imag part of point 1, etc. Each 4-byte float is transmitted least significant byte first. A total of 8n bytes are transferred. Assert EOI with the final byte of the transfer.

Serial Poll until IFC (bit 7) is set in the Serial Poll status before sending another command.

---

**TASC ? i, n**

The TASC command loads ASCII data into Trace i. The parameter i selects Trace 1 through 5. This command is valid only via the GPIB interface.

## 5-110 Data Transfer Commands

---

The target Trace *i* must already store data of the measurement type and length to be loaded from the interface. The loaded data simply takes the place of the stored trace data and assumes the measurement parameters (Measurement, Averaging, Window, Frequency span, etc.) and length of Trace *i*.

TASC loads the actual complex data values in the trace, NOT the data displayed in a particular view. Remember, views are simply different ways of showing the actual complex data.

A trace contains *n* complex points. The value of *n* depends upon the type and length of the target trace. *N* should be chosen equal to or less than the length of the target trace. If the length is less than that of the target trace the remainder of the trace will be filled with zeros.

Each complex point is two floating point numbers, representing the real and imaginary parts. Thus, the TASC command loads  $2n$  floating point numbers (in ASCII format).

The downloaded values are assumed to have the same units as the target trace. Depending upon the trace measurement, the values are assumed to be either unitless, Volts or Volts<sup>2</sup>. See TLOD for more.

The download sequence is as follows:

### **Host**

Send TASC ? *i*, *n*. Do NOT wait for IFC to be set in the Serial Poll status.

### **SR785**

Checks to make sure that Trace *i* has length *n*. Returns 1 (4-byte binary long int) if OK to begin data transfer. A return of 0 indicates that *n* is too large for Trace *i*.

### **Host**

On receipt of 1 (4-byte binary long int), executes a transfer to the SR785 of  $2n$  ASCII floats. The order is real part of point 0, imag part of point 0, real part of point 1, imag part of point 1, etc. Each ASCII float is delimited by a comma, space, tab or carriage return (NOT EOI). Assert EOI with the final byte of the transfer.

Serial Poll until IFC (bit 7) is set in the Serial Poll status before sending another command.

---

### **ALOD ? *n***

The ALOD command loads binary data into the Arbitrary Waveform buffer. This command is valid only via the GPIB interface.

The parameter *n* specifies the number of points to be loaded and must be a multiple of 2048. The Arbitrary Waveform buffer must already have sufficient memory allocated. After the load is complete, the Arbitrary Length will be *n* points.

Each point is a 4-byte IEEE float number and range from -1.0 to +1.0. A value of 1.0 is output at 100% amplitude as 1 Volt. The data should not exceed this range.

The download sequence is as follows:

**Host**

Send ALOD ? n. Do NOT wait for IFC to be set in the Serial Poll status.

**SR785**

Checks to make sure that the Arbitrary buffer has sufficient length. Returns 1 (4-byte binary long int) if OK to begin data transfer. A return of 0 indicates that n is too large for the buffer.

**Host**

On receipt of 1 (4-byte binary long int), executes a binary transfer to the SR785 of n 4-byte IEEE floats. The order is point 0, point 1, etc. Each 4-byte float is transmitted least significant byte first. A total of 4n bytes are transferred. Assert EOI with the final byte of the transfer.

Serial Poll until IFC (bit 7) is set in the Serial Poll status before sending another command.

---

### TGET ? i

The TGET command uploads a Trace buffer to the host computer. This command is valid only via the GPIB interface. The parameter i selects a Trace (1-5). Trace i must contain data.

The uploaded data should be saved in its entirety by the host computer. The saved data can be downloaded back to the SR785 at a later time using TPUT. The TGET and TPUT commands allow a host computer to save and reload a trace buffer without using disks.

TGET uploads trace information as well as trace data. TPUT restores the trace type, length and other attributes as well as the data.

The upload sequence is as follows:

**Host**

Send TGET ? i. Do NOT wait for IFC to be set in the Serial Poll status.

**SR785**

Returns n (4-byte binary long int) which is the number of bytes needed to transfer the Trace buffer.

**Host**

On receipt of n (4-byte binary long int), executes a binary read from the SR785 of n bytes. Expect EOI with the final byte of the transfer.

Serial Poll until IFC (bit 7) is set in the Serial Poll status before sending another command.

## 5-112 Data Transfer Commands

---

### **TPUT ? i**

The TPUT command downloads Trace buffer data to the SR785. This command is valid only via the GPIB interface. The parameter *i* selects a Trace (1-5). Trace *i* does not need to already contain data.

The download data must be data which was originally uploaded by TGET in its entirety. The TGET and TPUT commands allow a host computer to save and reload a Trace buffer without using disks.

The download sequence is as follows:

#### **Host**

Send TPUT ? *i*. Do NOT wait for IFC to be set in the Serial Poll status.

#### **SR785**

Returns 1 (4-byte binary long int) when OK to begin binary transfer to the SR785.

#### **Host**

On receipt of 1 (4-byte binary long int), executes a binary transfer to the SR785 of *n* bytes (as uploaded using TGET). Asserts EOI with the final byte of the transfer.

#### **SR785**

Receives *n* bytes and the EOI. If an error is detected during transfer, the Trace buffer reverts to the empty state and a command execution error occurs.

Serial Poll until IFC (bit 7) is set in the Serial Poll status before sending another command.

---

### **AGET ? i**

The AGET command uploads the Arbitrary Waveform buffer to the host computer. This command is valid only via the GPIB interface. The parameter *i* selects upload the whole buffer (0) or just the playback portion (1).

The uploaded data should be saved in its entirety by the host computer. The saved data can be downloaded back to the SR785 at a later time using APUT. The AGET and APUT commands allow a host computer to save and reload the Arbitrary buffer without using disks.

The upload sequence is as follows:

#### **Host**

Send AGET ? *i*. Do NOT wait for IFC to be set in the Serial Poll status.

#### **SR785**

Returns *n* (4-byte binary long int) which is the number of bytes needed to transfer the Arbitrary buffer.



**Host**

On receipt of n (4-byte binary long int), executes a binary read from the SR785 of n bytes. Expect EOI with the final byte of the transfer.

Serial Poll until IFC (bit 7) is set in the Serial Poll status before sending another command.

---

**APUT ?**

The APUT command downloads Arbitrary Waveform buffer data to the SR785. This command is valid only via the GPIB interface. After the download is complete, the Arbitrary Length will be set to the length of the download.

The download data must be data which was originally uploaded by AGET in its entirety. The AGET and APUT commands allow a host computer to save and reload the Arbitrary buffer without using disks.

The download sequence is as follows:

**Host**

Send APUT ?. Do NOT wait for IFC to be set in the Serial Poll status.

**SR785**

Returns 1 (4-byte binary long int) when OK to begin binary transfer to the SR785. A return of 0 indicates that there is no Arbitrary Waveform memory allocated.

**Host**

On receipt of 1 (4-byte binary long int), executes a binary transfer to the SR785 of n bytes (as uploaded using AGET). Asserts EOI with the final byte of the transfer.

**SR785**

Receives n bytes and the EOI. If an error is detected during transfer, the Arbitrary buffer reverts to the empty state and a command execution error occurs. This may occur if there is insufficient Arbitrary memory allocated.

Serial Poll until IFC (bit 7) is set in the Serial Poll status before sending another command.

---

**CGET ? i, j**

The CGET command uploads a Capture buffer to the host computer. This command is valid only via the GPIB interface. The parameter i selects Ch1 (0) or Ch2 (1). The parameter j selects upload the whole buffer (0) or just the playback portion (1).

The uploaded data should be saved in its entirety by the host computer. The saved data can be downloaded back to the SR785 at a later time using CPUT. The CGET and CPUT commands allow a host computer to save and reload the Capture buffer without using disks.

The upload sequence is as follows:

---

## 5-114 Data Transfer Commands

---

**Host**

Send CGET ? i, j. Do NOT wait for IFC to be set in the Serial Poll status.

**SR785**

Returns n (4-byte binary long int) which is the number of bytes needed to transfer the Capture buffer.

**Host**

On receipt of n (4-byte binary long int), executes a binary read from the SR785 of n bytes. Expect EOI with the final byte of the transfer.

Serial Poll until IFC (bit 7) is set in the Serial Poll status before sending another command.

---

**CPUT ? i**

The CPUT command downloads Capture buffer data to the SR785. This command is valid only via the GPIB interface. The parameter i selects Ch1 (0) or Ch2 (1). After the download is complete, the Capture Length will be set to the length of the download.

The download data must be data which was originally uploaded by CGET in its entirety. The CGET and CPUT commands allow a host computer to save and reload the Capture buffer without using disks.

The download sequence is as follows:

**Host**

Send CPUT ? i. Do NOT wait for IFC to be set in the Serial Poll status.

**SR785**

Returns 1 (4-byte binary long int) when OK to begin binary transfer to the SR785. A return of 0 indicates that there is no Capture memory allocated.

**Host**

On receipt of 1 (4-byte binary long int), executes a binary transfer to the SR785 of n bytes (as uploaded using CGET). Asserts EOI with the final byte of the transfer.

**SR785**

Receives n bytes and the EOI. If an error is detected during transfer, the Capture buffer reverts to the empty state and a command execution error occurs. This may occur if there is insufficient Capture memory allocated.

Serial Poll until IFC (bit 7) is set in the Serial Poll status before sending another command.

---

**WGET ?**

The WGET command uploads the Waterfall buffer to the host computer. This command is valid only via the GPIB interface.

The uploaded data should be saved in its entirety by the host computer. The saved data can be downloaded back to the SR785 at a later time using WPUT. The WGET and WPUT commands allow a host computer to save and reload the Waterfall buffer without using disks.

The upload sequence is as follows:

**Host**

Send WGET ?. Do NOT wait for IFC to be set in the Serial Poll status.

**SR785**

Returns n (4-byte binary long int) which is the number of bytes needed to transfer the Waterfall buffer.

**Host**

On receipt of n (4-byte binary long int), executes a binary read from the SR785 of n bytes. Expect EOI with the final byte of the transfer.

Serial Poll until IFC (bit 7) is set in the Serial Poll status before sending another command.

---

**WPUT ?**

The WPUT command downloads Waterfall buffer data to the SR785. This command is valid only via the GPIB interface.

The download data must be data which was originally uploaded by WGET in its entirety. The WGET and WPUT commands allow a host computer to save and reload the Waterfall buffer without using disks.

The download sequence is as follows:

**Host**

Send WPUT ?. Do NOT wait for IFC to be set in the Serial Poll status.

**SR785**

Returns 1 (4-byte binary long int) when OK to begin binary transfer to the SR785. A return of 0 indicates that there is no Waterfall memory allocated.

**Host**

On receipt of 1 (4-byte binary long int), executes a binary transfer to the SR785 of n bytes (as uploaded using WGET). Asserts EOI with the final byte of the transfer.

**SR785**

Receives n bytes and the EOI. If an error is detected during transfer, the Waterfall buffer reverts to the empty state and a command execution error occurs. This may occur if there is insufficient Waterfall memory allocated.

Serial Poll until IFC (bit 7) is set in the Serial Poll status before sending another command.

## 5-116 Data Transfer Commands

---

Loading the Waterfall buffer loads data into the waterfall buffer *and* recalls the measurement parameters which were in effect when the waterfall buffer was saved. The recalled parameters include items in the [Frequency], [Display Setup], [Display Options], [Average], [Window], [User Math], and [Waterfall] menus.

# Interface Commands

---

## **\*RST**

The \*RST command resets the SR785 to its default configurations. The communications setup is not changed. All other modes and settings are set to their default conditions and values. This command takes some time to complete.

Do not send other commands on the same command line after \*RST. Do not send another command until \*RST is complete. Either wait a fixed amount of time (approximately 12 seconds) or use serial polling to wait for IFC to be set again in the Serial Poll status before sending more commands.

---

## **\*IDN ?**

The \*IDN ? query returns the SR785's device identification string. This string is in the format "Stanford\_Research\_Systems,SR785, s/n00001, ver007". In this example, the serial number is 00001 and the firmware version is 007.

---

## **LOCL (?) {i}**

The LOCL command sets the local/remote state. The parameter *i* selects LOCAL (0), REMOTE (1) or LOCAL LOCKOUT (2). The states duplicate the GPIB local/remote states. In the LOCAL state, both command execution and keyboard input are allowed. In the REMOTE state command execution is allowed but the keyboard and knob are locked out except for the [Help/Local] key which returns the SR785 to the LOCAL state. In the LOCAL LOCKOUT state all front panel operation is locked out, including the [Help/Local] key.

---

## **OVRM (?) {i}**

The OVRM command sets (queries) the GPIB Override Remote mode. The parameter *i* selects No (0) or Yes (1). When Override Remote is set to Yes, the front panel is not locked out when the unit is in the REMOTE state.

---

## Nodal Degree-of-Freedom Commands

---

**SVNI (?) {i}**

The SVNI command determines (queries) whether nodal degree-of-freedom information will be stored with SR785 data (.78D) files. If this softkey is set to "on" nodal DOF information will be saved with the files and will be available for use by external programs which convert SR785 files into formats used by modal analysis programs. If set to "off" no nodal DOF information will be saved and the external programs must obtain this information from the user.

---

**RFDR (?) {i}**

The RFDR command sets (queries) the reference node direction. The parameter i specifies the direction according to the following table:

<u>i</u>	<u>Node Direction</u>
0	-θ Z
1	-θ Y
2	-θ X
3	-Z
4	-Y
5	-X
6	Scalar (no direction)
7	X
8	Y
9	Z
10	+θ X
11	+θ Y
12	+θ Z

---

**RFNA (?) {s}**

The RFNA command sets (queries) the name of the reference node. The name may be up to 6 characters long.

---

**RFNU (?) {i}**

The RFNU command sets (queries) the number of the reference node.

---

**RSDR (?) {i}**

The RSDR command sets (queries) the response node direction. The parameter i specifies the direction according to the following table:

---

<u>i</u>	<u>Node Direction</u>
0	-θ Z
1	-θ Y
2	-θ X
3	-Z
4	-Y
5	-X
6	Scalar (no direction)
7	X
8	Y
9	Z
10	+θ X
11	+θ Y
12	+θ Z

---

**RSNA (?) {s}**

The RSNA command sets (queries) the name of the reference node. The name may be up to 6 characters long.

---

**RSNU (?) {i}**

The RSNU command sets (queries) the number of the reference node.

## Status Reporting Commands

The Status Word definitions follow this section.

---

### **\*CLS**

The \*CLS command clears all status registers. The enable registers are not changed.

---

### **\*PSC (?) {i}**

The \*PSC command sets the value of the power-on status clear bit.

If i=1 the power-on status clear bit is set and all status registers and enable registers are cleared on power up.

If i=0 the bit is cleared and all enable registers are stored at power down. The status registers are cleared and the enable registers are restored to their stored values on power up. This allows a service request to be generated at power up.

---

### **\*SRE (?) {i} {, j}**

The \*SRE i command sets the Serial Poll enable register to the decimal value i (0-255). The \*SRE i, j command sets bit i (0-7) to j (0 or 1).

The \*SRE? command queries the value (0-255) of the serial poll enable register. The \*SRE? i command queries the value (0 or 1) of bit i (0-7).

When a bit becomes set in BOTH the Serial Poll status word AND the Serial Poll enable register, an SRQ (GPIB service request) is generated. The SRQ is cleared by performing a serial poll. The bit in the Serial Poll status word which caused the SRQ must be cleared before this bit can cause another SRQ. To clear this bit, the condition which causes it to be set in the Serial Poll status word needs to be cleared. For the INST, DISP, INPT, IERR or ESB bits, this is accomplished by clearing the enabled status bits in the Instrument, Display, Input, Error or Standard Event status words (by reading them).

---

### **\*STB ? {i}**

The \*STB? command queries the value of the Serial Poll status word. The value is returned as a decimal number from 0 to 255. The \*STB? i command queries the value (0 or 1) of bit i (0-7).

The value of bit 6 (SRQ) when read using \*STB? returns 1 if a bit is set in BOTH the Serial Poll status word AND the Serial Poll enable register. This is independent of serial polling and SRQ's. Bit 6 is the SRQ bit only when serial polled.

\*STB? has no effect on the value of the Serial Poll status word. To clear a bit in the Serial Poll status, the condition which causes it to be set must be cleared. For the INST, DISP, INPT, IERR or ESB bits, this is accomplished by clearing the enabled status bits in the Instrument, Display, Input, Error or Standard Event status words (by reading them).



**\*ESE (?) {i} {, j}**

The \*ESE i command sets the Standard Event enable register to the decimal value i (0-255). The \*ESE i, j command sets bit i (0-7) to j (0 or 1).

The \*ESE? command queries the value (0-255) of the status word enable register. The \*ESE? i command queries the value (0 or 1) of bit i (0-7).

When a bit becomes set in BOTH the Standard Event status word AND the Standard Event enable register, bit 6 (ESB) of the Serial Poll status word is set. This causes an SRQ if bit 6 in the Serial Poll enable register is set. To clear a bit in the Standard Event status word, use \*ESR?.

---

**\*ESR ? {i}**

The \*ESR? command queries the value of the Standard Event status word. The value is returned as a decimal number from 0 to 255. The \*ESR? i command queries the value (0 or 1) of bit i (0-7).

\*ESR? clears the entire word while \*ESR? i clears just bit i.

---

**ERRE (?) {i} {, j}**

The ERRE i command sets the Error enable register to the decimal value i (0-65535). The ERRE i, j command sets bit i (0-15) to j (0 or 1).

The ERRE? command queries the value (0-65535) of the error status enable register. The ERRE? i command queries the value (0 or 1) of bit i (0-15).

When a bit becomes set in BOTH the Error status word AND the Error status enable register, bit 3 (IERR) of the Serial Poll status word is set. This causes an SRQ if bit 3 in the Serial Poll enable register is set. To clear a bit in the Error status word, use ERRS?.

---

**ERRS ? {i}**

The ERRS? command queries the value of the Error status word. The value is returned as a decimal number from 0 to 65535. The ERRS? i command queries the value (0 or 1) of bit i (0-15).

ERRS? clears the entire word while ERRS? i clears just bit i.

---

**INSE (?) {i} {, j}**

The INSE i command sets the Instrument status enable register to the decimal value i (0-255). The INSE i, j command sets bit i (0-7) to j (0 or 1).

The INSE? command queries the value of the Instrument status enable register. The INSE? i command queries the value (0 or 1) of bit i (0-7).

When a bit becomes set in BOTH the Instrument status word AND the Instrument status enable register, bit 0 (INST) of the Serial Poll status word is set. This causes an SRQ if bit 0 in the Serial Poll enable register is set. To clear a bit in the Instrument status word, use INST?.

---

## 5-122 Status Reporting Commands

---

---

### INST ? {i}

The INST? command queries the value of the Instrument status word. The value is returned as a decimal number from 0 to 255. The INST? i command queries the value (0 or 1) of bit i (0-7).

INST? clears the entire word while INST? i clears just bit i.

---

### DSPE (?) {i} {, j}

The DSPE i command sets the Display status enable register to the decimal value i (0-65535). The DSPE i, j command sets bit i (0-15) to j (0 or 1).

The DSPE? command queries the value of the Display status enable register. The DSPE? i command queries the value (0 or 1) of bit i (0-15).

When a bit becomes set in BOTH the Display status word AND the Display status enable register, bit 1 (DISP) of the Serial Poll status word is set. This causes an SRQ if bit 1 in the Serial Poll enable register is set. To clear a bit in the Display status word, use DSPS?.

---

### DSPS ? {i}

The DSPS? command queries the value of the Display status word. The value is returned as a decimal number from 0 to 65535. The DSPS? i command queries the value (0 or 1) of bit i (0-15).

DSPS? clears the entire word while DSPS? i clears just bit i.

---

### INPE (?) {i} {, j}

The INPE i command sets the Input status enable register to the decimal value i (0-65535). The INPE i, j command sets bit i (0-15) to j (0 or 1).

The INPE? command queries the value of the Input status enable register. The INPE? i command queries the value (0 or 1) of bit i (0-15).

When a bit becomes set in BOTH the Input status word AND the Input status enable register, bit 2 (INPT) of the Serial Poll status word is set. This causes an SRQ if bit 2 in the Serial Poll enable register is set. To clear a bit in the Input status word, use INPS?.

---

### INPS ? {i}

The INPS? command queries the value of the Input status word. The value is returned as a decimal number from 0 to 65535. The INPS? i command queries the value (0 or 1) of bit i (0-15).

INPS? clears the entire word while INPS? i clears just bit i.

---

### INPC ? i

The INPC? command queries the current overload condition of input i. The parameter i selects Ch1 (0) or Ch2 (1).

INPC? returns a value from 0 to 3.

---

<u>return</u>	<u>meaning</u>
0	input is under half scale
1	input is over half scale
2	input is overloaded
3	input is HighV

INPC? always returns the current condition. If an overload occurs and goes away, INPC? will not detect it. Use the status words to detect momentary changes in the overload state.

## Status Word Definitions

The SR785 reports on its status by means of six status words: the Serial Poll, Standard Event, Instrument, Display, Input and Error status words.

Upon power-on, the SR785 may either clear all of its status *enable* registers or maintain them in the state they were in on power-down. The \*PSC command determines which action will be taken.

The status bits are set to 1 when the event or state described in the tables below has occurred or is present.

---

### Serial Poll Status Word

<u>Bit</u>	<u>Name</u>	<u>Set when</u>
0	INST	An enabled bit in the Instrument status word is set
1	DISP	An enabled bit in the Display status word is set
2	INPT	An enabled bit in the Input status word is set
3	IERR	An enabled bit in the Error status word is set
4	MAV	The interface output buffer is not empty
5	ESB	An enabled bit in the Standard status word is set
6	SRQ	An SRQ has occurred (serial poll), or an enabled bit in the this status word is set (*STB?)
7	IFC	No command execution is in progress

The INST, DISP, INPT, IERR and ESB bits are set whenever any bit in both their respective status words AND enable registers is set. Use the INSE, DSPE, INPE, ERRE and \*ESE commands to set the enable register bits. The INST, DISP, INPT, IERR and ESB status bits are not cleared until ALL enabled status bits in the Instrument, Display, Input, Error and Standard Event status words are cleared (by reading the status words or using \*CLS).

Use serial poll or \*STB? to read the Serial Poll status. Use \*SRE to set bits in the Serial Poll enable register.

#### Using \*STB? to Read the Serial Poll Status

The Serial Poll status word is NOT cleared by reading the status using \*STB?. A bit stays set as long as the status condition exists. This is true even for bit 6 (SRQ). Bit 6 will be set whenever the same bit in the Serial Poll status word AND Serial Poll enable register is set. This is independent of whether a serial poll has occurred to clear the service request.

#### Using Serial Poll to Read the Serial Poll Status

Except for SRQ, a bit in the Serial Poll status word is NOT cleared by serial polling the status word. A bit stays set as long as the status condition exists.

When reading the status using a serial poll, the SRQ bit signals that the SR785 is requesting service. The SRQ bit will be set (1) the first time the SR785 is polled

following a service request. The serial poll automatically clears the service request. Subsequent serial polls will return SRQ cleared (0) until another service request occurs. Polling the status word and reading it with \*STB? can return different values for SRQ. When serial polled, SRQ indicates a service request has occurred. When read with \*STB?, SRQ indicates that an enabled status bit is set.

### **Service Requests (SRQ)**

A GPIB service request (SRQ) will be generated whenever a bit in both the Serial Poll status word AND Serial Poll enable register is set. Use \*SRE to set bits in the Serial Poll enable register. A service request is only generated when an enabled Serial Poll status bit becomes set (changes from 0 to 1). An enabled status bit which becomes set and remains set will generate a single SRQ. If another service request from the same status bit is desired, the requesting status bit must first be cleared. In the case of the INST, DISP, INPT, IERR and ESB bits, this means clearing the enabled bits in the Instrument, Display, Input, Error or Standard Event status words (by reading them) or clearing the appropriate bits in the corresponding enable registers. Multiple enabled bits in these status words will generate a single SRQ. Another SRQ (from INST, DISP, INPT, IERR or ESB) can only be generated after clearing the INST, DISP, INPT, IERR or ESB bits in the Serial Poll status word.

The controller should respond to the SRQ by performing a serial poll to read the Serial Poll status word to determine the requesting status bit. Bit 6 (SRQ) will be reset by the serial poll.

For example, to generate a service request when a TRIGGER occurs, bit 0 in the Instrument Status enable register needs to be set (INSE 1 command) and bit 0 in the Serial Poll enable register must be set (\*SRE 1 command). When a trigger occurs, bit 0 in the Instrument status word is set. Since bit 0 in the Instrument status word AND enable register are set, this ALSO sets bit 0 (INST) in the Serial Poll status word. Since bit 0 in the Serial Poll status word AND enable register are set, an SRQ is generated. Bit 6 (SRQ) in the Serial Poll status word is set. Further triggering will not generate another SRQ until the TRIGGER status bit is cleared. The TRIGGER status bit is cleared by reading the Instrument status word (with INST?) or clearing bit 0 in the Instrument status enable register (with INSE). Presumably, the controller is alerted to the trigger via the SRQ, performs a serial poll to clear the SRQ, does something in response to the trigger (read data for example) and then clears the TRIGGER status bit by reading the Instrument status register. A subsequent trigger will then generate another SRQ.

### Standard Event Status Word

---

<u>Bit</u>	<u>Name</u>	<u>Set when</u>
0	unused	
1	unused	
2	QRY	Too many responses are pending
3	DDE	Too many commands received at once
4	EXE	Command cannot execute successfully
5	CME	Command syntax error
6	URQ	A key is pressed or the knob rotated
7	PON	Power is turned on
8-15	unused	reserved by IEEE standard

The Standard Event status word is defined by IEEE-488.2 (1987) and is used primarily to report errors in commands received over the interface. These status bits remain set until read using \*ESR?. They are also cleared by the \*CLS command. Use \*ESE to set bits in the Standard Event enable register.

### Instrument Status Word

---

<u>Bit</u>	<u>Name</u>	<u>Set when</u>
0	TRIG	A measurement has been triggered
1	DISK	A disk operation has been completed
2	OUTP	A hardcopy output operation has completed
3	TACH	A tach pulse has been received
4	CAPT	The capture buffer has filled
5	PAUS	Measurement has been paused
6	STRT	Measurement has been started
7	PLBK	Single Shot Capture Playback has finished
8	PREV	Measurement stopped to wait for average preview
9-15	unused	

The Instrument status word reports on the overall measurement status of the instrument. These status bits remain set until read by INST?. They are also cleared by the \*CLS command. Use INSE to set bits in the Instrument status enable register.

**Display Status Word**

<u>Bit</u>	<u>Name</u>	<u>Set when</u>
0	NEWA	New data is available for DisplayA
1	AVGA	DisplayA linear average has completed
2	STLA	New settled data is available for DisplayA
3	LIMA	DisplayA has failed a limit test
4	SSA	Swept Sine sweep has finished
5	WFA	Display A 1-shot Waterfall has finished
6	WFD	Display A Waterfall has finished drawing
7	unused	
8	NEWB	New data is available for DisplayB
9	AVGB	DisplayB linear average has completed
10	STLB	New settled data is available for DisplayB
11	LIMB	DisplayB has failed a limit test
12	SSB	Swept Sine sweep has finished
13	WFB	Display B 1-shot Waterfall has finished
14	WFB	Display B Waterfall has finished drawing
15	unused	

The Display status bits stay set until read by DSPTS?. They are also cleared by the \*CLS command. Use DSPE to set bits in the Display status enable register.

The AVG, STL, LIM and SS bits are only updated when new data is available for the display (no faster than 8 Hz). Always test for NEWA or NEWB in the status word along with AVG, STL, LIM or SS.

The status updates for the two displays may not be exactly synchronous. Always test the DisplayA status bits separately from the DisplayB status bits. For example, to wait for both displays to finish averaging, do NOT test the Display status word for both AVGA and AVGB simultaneously. The AVGA and AVGB bits may not be set at the same time. Test for AVGA and AVGB separately and wait until both have occurred.

## 5-128 Status Word Definitions

---

### Input Status Word

<u>Bit</u>	<u>Name</u>	<u>Set when</u>
0	LOW1	Input1 has fallen below half full scale
1	HLF1	Input1 has exceeded half full scale
2	OVL1	Input1 has exceeded full scale
3	HIV1	Input1 has exceeded 35 dBV Input range switched to 34 dBV.
4	ARG1	Input1 has AutoRanged
5-7	unused	
8	LOW2	Input2 has fallen below half full scale
9	HLF2	Input2 has exceeded half full scale
10	OVL2	Input2 has exceeded full scale
11	HIV2	Input2 has exceeded 35 dBV. Input range switched to 34 dBV.
12	ARG2	Input2 has AutoRanged
13-15	unused	

These status bits do not necessarily reflect the current input status. The Input status bits stay set until read by INPS?. They are also cleared by the \*CLS command. Use INPE to set bits in the Input status enable register.

To read the current input overload condition, use INPC?.

---

### Error Status Word

<u>Bit</u>	<u>Name</u>	<u>Set when</u>
0	OUTE	An output error has occurred (Print, Plot or Dump)
1	DSKE	A disk error has occurred
2	FLTE	A math error has occurred
3	RAME	RAM memory test fails
4	ROME	ROM memory test fails
5	VIDE	Video memory test fails
6	HELPE	Help memory test fails
7	DSDE	DSP data memory test fails
8	DSPE	DSP program memory test fails
9	DSRE	DSP DRAM memory test fails
10	CAL0	DSP calibration memory fails
11	CAL1	Ch1 calibration memory fails
12	CAL2	Ch2 calibration memory fails
13-15	unused	

The Error status bits stay set until read by ERRS?. They are also cleared by the \*CLS command. Use ERRE to set bits in the Error status enable register.

Bits 7-12 are set by the power on tests. Bits 3-5 are set in the [System]<Diagnostics>.



# Example Program

## Using Microsoft C with the GPIB interface

To successfully interface the SR785 to a PC via the GPIB interface, the instrument, interface card and interface drivers must all be configured properly. To configure the SR785, the GPIB address must be set in the [System]<Remote> menu. The default address is 10; use this address unless a conflict occurs with other instruments in your system. The SR785 will be set to GPIB address 10 whenever a reset is performed (power on with the backspace key down).



Make sure that you follow all of the instructions for installing the GPIB card. The GPIB card cannot be simply unpacked and put into your computer. To configure the card, you may need to set jumpers and switches on the card to set the I/O address and interrupt levels. Refer to your manual for more information.

### Capital Equipment Corp. CEC488 GPIB Card

The CEC488 card contains its low level drivers in ROM. The card address needs to be set so as not to conflict with other devices in your computer. The software interface uses a header file and a link library. There is no initialization required for the drivers.

Use the TEST488 program to test the card installation. Use TRTEST to communicate directly from the keyboard with the SR785. If TRTEST doesn't work, then your programs will not run.

The example is written using the CEC library routines.

### National Instruments GPIB Card

You must run the program "IBCONF" to configure the resident GPIB driver for your GPIB card. Please refer to the National Instruments manual for more information. For example, the following options should be set with IBCONF:

Device Name:	SR785
Device Address:	10
EOS Character:	0Ah (linefeed)
Terminate Read on EOS:	Yes

Once all the hardware and GPIB drivers are configured, use "IBIC". Use "IBWRT" and "IBRD" to send to and receive from the SR785. If you cannot talk to the SR785 via "IBIC", then your programs will not run.

To modify the example to use a National Instruments card, modify the routines where indicated.

### Other GPIB cards

You need to setup and configure your card according to your manual. The example program points out the routines which are interface dependent. Your card should have functions equivalent to those used in the example.

## 5-130 Example Program

---

```
// *****
// Example program using Microsoft C v7.0 and
// the Capital Equipment Corp CEC488 GPIB interface card.
// To use another interface card, modify the GPIB subroutines
// where indicated.

// The SR785 is assumed to be at address 10 (default).

// Link this object file with ieee488.lib (or the library for your GPIB card).

// *****
// Connect the Source output to the Ch1 A and Ch2 A inputs!!
// *****

#include <stdio.h>
#include <stdlib.h>
#include <string.h>
#include <conio.h>
#include <math.h>
#include <ieee-c.h> //This is the CEC header file
                    //Use the .h file for your interface card.
                    //send(), enter(), transmit(), rarray(), tarray() are
                    //CEC routines defined in ieee-c.h

char recv[80]; //global GPIB receive string
char cmd[80]; //global GPIB command string
float rxBuff[401]; //array of IEEE floats to receive binary data
float txBuff[2048]; //array of IEEE floats to send binary data
int status, length; //CEC interface routines use these variables

#define SR785 10 // default GPIB address for SR785

// subroutines at the end of this listing
// you may need to modify these for your interface card!
void TxGpib (int, char *);
void GetGpib (int);
void TxSR785 (char *);
void GetSR785 (char *);
void WaitAvg (void);
double GetData (int, int);
void WaitIFC (void);
void GetSpace (void);

void main (void);

/* ***** */
void main (void)
{
// *****
// You can see the commands received and responses generated by the SR785
// by using [System]<Remote><View Qs> to display the interface buffers
// on the screen.
// *****

int i, mode, type, number, nlen;
long ack;
double dispAVal, dispBVal, octreal, octimag;

printf ("\nSR785 Example Program\n\n");

// ***** initialize your GPIB card here *****
initialize (21,0); //controller
settimeout (7500); //GPIB timeout to 7.5 seconds

// ***** Set the SR785 output interface to GPIB!! *****
TxSR785 ("*CLS"); //clear all status words
TxSR785 ("OUTX0"); //direct SR785 responses to GPIB interface
```

```

// ***** set the Standard Event enable register to catch EXE and CME
// Command errors will set the ESB bit in the Serial Poll status word
TxSR785 ("*ESE 48"); //set bits 4 and 5

// ***** Check communication by querying SR785 IDN string
GetSR785 ("*IDN?");
printf ("[GPIB IDN] = %s\n\n", recv);

printf ("Reset SR785...");
TxSR785 ("*RST"); // reset the SR785 to a default state
printf ("done\n\n");

/* ***** set up FFT averaging for both displays ***** */
TxSR785 ("FAVM 2,1"); //display avg = rms
TxSR785 ("FAVT 2,0"); //type = linear
TxSR785 ("FAVN 2,250"); //number = 250
TxSR785 ("FAVG 2,1"); //avg = On

TxSR785 ("SRCO 1"); // turn on the source

// ***** query these parameters for DisplayA *****
GetSR785 ("FAVM? 0"); mode = atoi(recv);
GetSR785 ("FAVT? 0"); type = atoi(recv);
GetSR785 ("FAVN? 0"); number = atoi(recv);
printf ("Disp Avg = %d : type = %d : number = %d\n\n",mode,type,number);

printf ("Avg start...");
WaitAvg (); //start a linear average and wait until done
printf ("done\n");

TxSR785 ("ASCL 0; ASCL 1"); //autoscale the displays

/* ***** read data points ***** */
dispAVal = GetData (0,40); //read DisplayA bin 40 (peak)
dispBVal = GetData (1,40); //read DisplayB bin 40
printf ("DisplayA = %lf dBVpk\n",dispAVal);
printf ("DisplayB = %lf dBVpk\n",dispBVal);

/* ***** Binary transfer all of DisplayA spectrum ***** */
//We need to send the DSPB?0 command WITHOUT waiting for
//IFC in serial poll status since IFC will not be set until
//AFTER the transfer is complete!
//This section needs to be modified for your GPIB interface

printf ("\nReading entire DisplayA...");
TxGpib (SR785,"DSPB? 0"); //use TxGpib (don't wait for IFC)
transmit ("MLA TALK 10", &status); //make the PC listen, SR785 talk
rarray (rxBuff, 1604, &length, &status); //binary read 1604 bytes
//401 points x 4 bytes/point
WaitIFC (); //serial poll until IFC set, ok to continue

printf ("%d bytes read\n",length); //actual number of bytes read

// ***** print a few points from the spectrum
for (i=37; i<44; i++) {printf ("bin %d = %f dBVpk\n", i, rxBuff[i]);}
printf ("\n");

/* ***** Turn on limit testing ***** */
TxSR785 ("LCLR 0"); //first clear any existing limits for DisplayA
TxSR785 ("LMAX 0,2"); //initialize 2 segments
TxSR785 ("LSEG 0,0,0,12.032E3,-90.0,102.4E3,-90.0"); //upper segment
TxSR785 ("LSEG 0,1,0,1.28E3,-90.0,8.96E3,-90.0"); //upper segment

```

## 5-132 Example Program

---

```
TxSR785 ("LTST 0,1"); //turn on testing

// ***** take a new measurement
printf ("Avg start...");
WaitAvg ();
printf ("done\n");

GetSR785 ("LFAL?0"); //read limit test result
printf ("Limit test = ");
if (atoi(recv)==0) printf ("pass"); else printf ("fail");
printf ("\n\n");

// ***** clean up
TxSR785 ("LCLR 0"); //clear limits, turn off testing
TxSR785 ("FAVG 2,0"); //turn avg off

/* *****
/* ***** Load the arbitrary waveform buffer ***** */
// construct a normalized ramp of 2048 points
// each point is a 4 byte IEEE float between -1.0 and +1.0
for (i=0; i<2048; i++) { txBuff[i] = ((float) i)/2048.F; }

//This section needs to be modified for your GPIB interface

printf ("Loading Arb buffer...");
TxGpib (SR785,"ALOD? 2048"); //use TxGpib so we don't wait for IFC
//SR785 will return a binary 1 to acknowledge
transmit("mla talk 10",&status); //make SR785 a talker, pc a listener
rarray (&ack, 4, &length, &status); //binary read 4 bytes (long int)
//return value should be 1
printf ("%d bytes recvd [val=%ld]...",length,ack);

transmit("mta listen 10",&status); //make SR785 a listener, pc a talker
tarray(txBuff,8192,1,&status); //binary transfer of 8192 bytes to SR785
//set EOI with last byte
WaitIFC (); //serial poll until IFC set, ok to continue
printf ("done\n");

TxSR785 ("STYP 3"); //change source type to Arb
TxSR785 ("I1RG 4; I2RG 4"); //increase the input ranges
TxSR785 ("FSPN 2,3.2E3"); //narrow the span to display the spectrum
GetSpace (); //pause to look

/* *****
/* ***** Switch to Octave Analysis ***** */
TxSR785 ("STYP 2"); //switch to Noise source
TxSR785 ("NTYP 2"); //choose Pink noise
TxSR785 ("MGRP 2,2"); //switch to Octave measurement group

printf ("\n\nWait for octave settle...");
GetSR785 ("DSPS?"); //clear any sticky bits in the Display status word
//and then wait until the settle bits become set
do { GetSR785 ("DSPS?"); } while ( (atoi(recv)&0x0404) != 0x0404 );
printf ("done\n\n");

TxSR785 ("SVTR 0,5"); //save DisplayA to Trace5
TxSR785 ("RCTR 1,5"); //recall Trace5 to DisplayB
GetSR785 ("DSPN? 1"); nlen = atoi(recv); //query the length of DisplayB
printf ("DisplayB length = %d\n",nlen);

/* *****
/* ***** Download an octave trace ***** */
// Make Trace5 a constant -27 dBVpk for all bins.
// Convert -27 dBVpk to Volts (10^(-27/20))=0.04467
// Trace5 is already defined as an Octave measurement
// so the downloaded points are assumed to be Volts squared!
```

```

octreal = 0.04467*0.04467;    //(44.67 mV)^2 real part
octimag = 0.0;                //no imaginary part

//This section needs to be modified for your GPIB interface

printf ("Loading Trace 5...");
TxGpib (SR785,"TASC? 5,34");    //use TxGpib (no wait for IFC RDY)
                                //TASC downloads ASCII data,
                                //use TLOD to download binary data.
transmit("mla talk 10",&status); //make SR785 a talker, pc a listener
rarray (&ack, 4, &length, &status); //binary read 4 bytes (long int)
                                //return value should be 1
printf ("%d bytes recvd [val=%ld]...",length,ack);

// TASC (TLOD) replaces the data in an existing Trace (5).
// The Trace measurement and length are not changed!
// Download 34 bins of data (real, imaginary) with commas to separate,
// and EOI at the end.
// First 33 bins are the spectrum, last bin is the Total Power bin.
// All points are complex.
// For an octave measurement, the imaginary parts are all zero.
transmit("mta listen 10",&status); //make SR785 a listener, pc a talker
for (i=0; i<33; i++) {          //download 33 octave bins
    sprintf (cmd,"DATA '%lf, '", octreal); //real part
    transmit (cmd,&status);
    sprintf (cmd,"DATA '%lf, '", octimag); //imaginary part
    transmit (cmd,&status);
}
// last bin is the total power bin
sprintf (cmd,"DATA '%lf, '", octreal*33); //real part of total power bin
transmit (cmd,&status);
sprintf (cmd,"DATA '%lf' END", 0.0);    //imag part of total power bin,
transmit (cmd,&status);                  //last value, terminate with EOI!
WaitIFC ();                             //serial poll until IFC set, ok to continue

printf ("done\n");

TxSR785 ("RCTR 1,5"); //recall Trace 5 to DisplayB to view the new data

/* *****
/* ***** Define a User Function ***** */
// define Funcl = Oct(1)/Trace5
TxSR785 ("USRO 1,36,-5,5"); //36=Oct(1), -5=divide, 5=Trace5
TxSR785 ("MEAS 0,37");     //DisplayA measures UserFuncl
TxSR785 ("YMAX 0,40");     //Change the top graph scale.
                                //Funcl is Oct(1) normalized by Trace5.
                                //The amplitudes are approximately 0 dB.
                                //The Total Power bin is calculated from the
                                //sum of the octave bins in the function,
                                //(it is not normalized to Trace5).
// ***** end of main program
}

/* *****
void TxGpib (int address,char *command)
// routine to transmit command to a GPIB address
// modify this routine for your GPIB interface

{
send (address,command,&status); // send() is the CEC routine to send
                                // a string to an address.
                                // Sets status = 0 to indicate result OK

if (status != 0) {
// handle transmit errors here
printf ("\nCommand = %s\n",command);
printf ("Error at device %d : status = %d\n",address,status);
GetSpace ();
}
}

```

## 5-134 Example Program

---

```
    }
/* ***** */
void GetGpib (int address)
// routine to get an answer from a GPIB address
// modify this routine for your GPIB interface

{
char temp[80];

enter (temp,80,&length,address,&status); // enter () is the CEC routine
// to enter a string from an address.
// Sets status = 0 to indicate result OK.
// 80 is maxlen, actual received length
// is stored in &length.

if (status != 0) {
// handle receive errors here
printf ("\nError at device %d : status = %d\n",address,status);
GetSpace ();
}

strcpy (recv, temp); //set global recv string
}
/* ***** */
void TxSR785 (char *command)
// Send command to the SR785 and wait until IFC (bit7) is set in
// the Serial Poll status (indicating that the command is finished).

{
TxGpib (SR785,command); // send command to SR785 address
WaitIFC (); //serial poll until IFC set, ok to continue
}

/* ***** */
void WaitIFC (void)
{
//serial poll until IFC (bit7) set (command done)
//modify for your GPIB interface

//stb is serial poll byte
char stb;

do {spoll(SR785,&stb,&status);} while (!(stb&128)); //IFC

if (stb&32) {
// If ESB bit set,
// there must be a command error in the Standard Event status word.
// Handle command errors here.
TxGpib (SR785, "*ESR?"); //clear the Standard Event status word
GetGpib (SR785);
printf ("\nEXE error\n");
GetSpace ();
}
}

/* ***** */
void GetSR785 (char *getcmd)
// query the SR785 for an answer
// getcmd is the query command string

{
TxSR785 (getcmd); //send query command
GetGpib (SR785); //get response into global recv string
}
/* ***** */
void WaitAvg (void)
{
// routine to start a linear avg and wait until both displays are done

int result, avgdone;
```

```

avgdone=0;           //init avgdone status
GetSR785 ("DSPS?"); //clear sticky bits in Display status word first

TxSR785 ("STRT");   //start lin average measurement

do {
    GetSR785 ("DSPS?"); result=atoi(recv); //read display status word
    avgdone = avgdone | (result & 0x0202);
    //Check for AVGA and AVGB status bits until both have occurred.
    //Remember, they may not be set together in the same query!
    //Do not check for both in the same query!
}
while ((avgdone&0x0202)!=0x0202); // both AVGA and AVGB have occurred
}
/* ***** */
double GetData (int disp, int bin)
{
    //routine to move the display marker to a bin and return the data value

    //move the marker in display disp to bin
    sprintf (cmd,"MBIN %d,%d",disp,bin);
    TxSR785 (cmd);

    sprintf (cmd,"DSPY? %d,%d",disp,bin); //read the data value at bin
    GetSR785 (cmd);
    return (atof(recv)); //return the value as a double
}
/* ***** */
void GetSpace (void)      /* Wait for space key */
{
    char ch;

    while ( kbhit() ) getch(); // clear the keyboard buffer
    printf ("\n<Space> to continue, <Q> to quit ");
    do {
        ch = (char) getch ();
        if ((ch == 'q')||(ch =='Q')) exit(0); // exit
    }
    while (ch!=' '); // continue
    printf ("\n");
}
/* ***** */

```





---

## Chapter 6

# File Conversions

---

### In This Chapter

---

<b>Why File Conversion?</b>	<b>6-2</b>
<b>SR785 File Types</b>	<b>6-3</b>
.78D Files	6-3
Buffers: .78C and .78W Files	6-3
<b>Supported External File Types</b>	<b>6-4</b>
ASCII Files	6-4
SDF Files	6-4
Matlab <sup>®</sup> MAT-Files	6-4
Universal File Format	6-4
SR780 .78D Files	6-4
<b>Using the File Conversion Utility</b>	<b>6-5</b>
Using SRTRANS	6-5
Converting Files to ASCII	6-5
Converting Files to SDF	6-7
Converting Files to .78D Format	6-8
Converting Files to MATLAB Format	6-8
Converting Files to Universal File Format	6-9
Converting Files to Capture File Format (.78C)	6-10
Converting Files to Arbitrary Waveform Format (.78C)	6-10

## Why File Conversion?

Data from the SR785 will often need to be exported to external programs such as spreadsheets, plotting programs, or modal-analysis programs. In addition, it is often convenient to be able to import data prepared by external programs into the SR785 for display and comparison. In order to meet both these needs the SR785 is shipped with a file conversion program designed to import and export data as flexibly as possible

Stanford Research Systems is committed to supporting as many file formats as possible- as a result we are constantly adding new formats. Contact Stanford Research Systems to obtain a list of supported programs and file formats and to obtain a free upgrade to the latest file conversion utilities.

## SR785 File Types

Three distinct native SR785 file types are supported by the SR785 file conversion utility.

### **.78D Files**

.78D files, or display files, are generated using the <Trace to Disk> or <Display to Disk> softkeys or by using the TGET?i remote command. The display file contains both the trace data, and a complete record of the instrument's setup when the trace was stored.

### **Buffers: .78C and .78W Files**

The SR785 uses 4 user-allocable buffers to store large amounts of data. These are the Channel 1 capture buffer, the Channel 2 capture buffer, the waterfall buffer, and the arbitrary waveform buffer. Capture and arbitrary waveform buffers have the same format and can be saved to disk as .78C files using the <Buffer to Disk> softkey or the AGET?i and CGET?i,j commands. Waterfall buffers have a different format and can be saved to disk as .78W files using the <Buffer to Disk> softkey or the WGET?i command. The conversion utility can convert ASCII files to .78C files and to convert both .78C and .78W files into ASCII files.

## Supported External File Types

Five different file types are currently supported by the SR785 file conversion utilities.

### ASCII Files

ASCII files consist of plain ASCII text. ASCII files written by the SR785 file conversion utility consist of an alphanumeric header with instrument setup information followed by a number of columns of data values. Columns are separated by spaces, 'TAB' characters, or commas depending on the output option selected. Each line is terminated with a carriage-return line-feed pair. For ASCII input, the conversion utility expects values to be separated by any whitespace character with no leading header data. Characters after the last expected data value are ignored.

### SDF Files

Standard Data Format, or SDF, is a file format supported by Hewlett-Packard instruments and is designed to allow data sharing among various analyzers and external programs. Many modal analysis programs, for example, accept SDF files. The SR785 file conversion utilities can convert version 2 or 3 SDF files to SRS trace (.78D) files and can convert .78D files to SDF files. Note that because of specific differences between HP analyzers and the SR785 it may not be possible to exactly convert all types of measurements..

### Matlab® MAT-Files

MATLAB® is a popular program for doing mathematics and for the analysis of numeric data. Although it is easy to import ASCII data directly into MATLAB®, it is faster and more efficient to use the binary MAT-File format. The conversion utilities include the capability of directly converting a SR785 display file into a binary MAT-File. The MAT-File format used by the conversion utilities is supported by MATLAB version 4 and version 5.

### Universal File Format

The Universal File Format is a specification for an ASCII database capable of describing the geometry, degree-of-freedom, and measurement data used in structural analysis. Most modal analysis programs are capable of importing measurement data formatted in Universal File Format. The SR785 file conversion utilities include the capability of converting SR785 display files into Universal File Type 58 (Function at Nodal DOF) data sets.

### SR780 .78D Files

The file conversion utility will accept SR780 .78D files anywhere a SR785 .78D file can be used. This allows you to convert your old SR780 files to the SR785 format. The file conversion utility, however will not write SR780 files as output nor will it accept SR780 .78C or .78W files as input.

## Using the File Conversion Utility

The disk supplied with the SR785 contains the conversion utility SRTRANS.EXE. There is no installation required, simply copy to the program to your hard disk. The program has a number of options, some of which are generally applicable, and some of which are specific to certain output file types. In the following discussion of the program options, text in [brackets] indicates optional parameters, and text in *italics* indicates that the text should be replaced by the indicated parameters when using the program

### Using SRTRANS

Usage: SRT785 [*options*] *infile* [*outfile*]

*Infile* specifies the file to be converted. Wildcard file characters (\* and ?) are accepted.

*Outfile* specifies the name of the output file. If no name is supplied, the program creates an output file with the same name as the input file and an extension appropriate to the output file type.

Universal Options:

/Oxxx	Output File Type asc = ASCII (default) 78d = .78D SR785 display file sdf = HP SDF (v3) File mat = MATLAB binary file uff = Universal File Format cap = capture buffer (.78C) file arb = arbitrary waveform buffer
/A	Append to output file (uff and ASCII outputs only)
/E <i>filename</i>	Specify error log file. Default is Esrslg.txt
/J	Append to logfile. Default is replace logfile
/?	Help. For Help on options pertaining to a specific output file specify the output file type with no input file, e.g. srtrans /Oasc

Frequently used program options may be kept in a file named "SRT785.INI" in the current directory. The .INI file is scanned for options before the command line options are processed.

### Converting Files to ASCII

Usage: SRT785 /Oasc [*options*] *infile* [*outfile*]

Allowed Input File Types: .78D Files, .78C files, .78W files.

Default Output File Extension: .ASC

ASCII Output Options:

/D	Data only in output, no descriptive header.
----	---

## 6-6 Using the File Conversion Utilities

---

/H1	Header Only, Standard Header (default)
/H2	Header Only, Exhaustive Header. The exhaustive header contains a list of <i>all</i> of the instrument settings.
/Pn	Output n points per line (default is /P1)
/Fstring	Data Format String (default is /F% 14.6e). This is a standard 'C' language printf format specifier which determines the format of numbers in the output file.
/K0	Separate fields with spaces (default)
/K1	Separate fields with commas
/K2	Separate fields with tabs
/A	Append to output file (uff and ASCII outputs only)

Additional Options for .78D and .78W input files:

/Cview1,view2,...	Specify Output Columns (default is /Cx,v) x = x-axis values v = view shown when file was saved mag = magnitude m2 = magnitude squared r = real part i = imaginary part p = phase u = unwrapped phase
/Upk,dB,deg	Specify unit choices (not all choices apply to all views) pk = rms, pk, or pp (rms units, peak units, or peak to peak.) dB = lin, dB, dBm, dBSp1 (lin = linear units) deg= deg, rad (degrees or radians)

Additional Options .78W input files only:

/Wch	Specify display, ch = A or B
/Sn	Specify starting trace in waterfall buffer (0 is first trace)
/Rn	Specify number of traces to extract. -1 means extract to end of buffer (default)

Examples:

Convert the files DATA0.78D through DATA9.78D to ASCII. Show the x-axis values and the real and imaginary parts of the y values in the output files using rms values:

**SRT785 /Oasc /Urms /Cx,r,i DATA?.78D**

Convert the 10th through 15th waterfall traces in WATER.78W to an ASCII file with magnitude squared data values only:

**SRT785 /Oasc /S10 /R6 /Cmag2 WATER.78W**

An example of a .78D file converted with the default program options is shown below:

```

Input File: SRS001.78D
Measure Group: FFT
Measurement: FFT 1
Num of extracted Points: 401
Start Freq: 0 Hz
Span: 102.4 kHz
FFT Lines: 400
Window: BMH
Averaging Mode: Vector
Averaging Type: Exp. / Cont.
Number of Avg: 250
Overload reject: Off

```

X[Hz]	Mag[dBVpk]
0.000000e+00	-5.835680e+01
2.560000e+02	-6.169995e+01
5.120000e+02	-7.247395e+01
7.680000e+02	-9.415229e+01
1.024000e+03	-1.309051e+02
1.280000e+03	-1.354611e+02
1.536000e+03	-1.255198e+02
2.560000e+02	-6.169995e+01
5.120000e+02	-7.247395e+01
7.680000e+02	-9.415229e+01
1.024000e+03	-1.309051e+02
1.280000e+03	-1.354611e+02
1.536000e+03	-1.255198e+02

etc.,

## Converting Files to SDF

Usage: SRT785 /Osd *[options] infile [outfile]*

Allowed Input File Types: .78D Files.

Default Output File Extension: .DAT

For input files which contain nodal degree-of-freedom informat, the following SDF Output Options can be used to change the DOF parameters:

<i>/Mrnn</i>	Reference node number
<i>/Mrdk</i>	Reference node direction (see table below)
<i>/Msn</i>	Response node number
<i>/Msdk</i>	Response node direction
	Node Direction id's
0	- $\theta_z$
1	- $\theta_y$
2	- $\theta_x$
3	-Z
4	-Y
5	-X
6	Scalar
7	+Z
8	+Y
9	+X
10	+ $\theta_x$
11	+ $\theta_y$





$p$  = phase  
 $u$  = unwrapped phase  
*/Upk,dB,deg* Specify unit choices (not all choices apply to all views)  
 $pk$  = rms, pk, or pp (rms units, peak units, or peak to peak)  
 $dB$  = lin, dB, dBm , dBSpl (lin = linear units)  
 $deg$  = deg, rad (degrees or radians)

Example:

Convert an SR780 file named MYFILE.78D into a MAT file named M1.MAT containing a MATLAB array named data1 with 2 columns containing the magnitude and phase of the input data:

**SRT785 /Omat /Gdata1 /Cmag,p MYFILE.78D M1.MAT**

### Converting Files to Universal File Format

Usage: SRT785 /Ouff [*options*] *infile* [*outfile*]

Allowed Input File Types: .78D Files.

Default Output File Extension: .UFF

UFF output options:

<i>/Mrinodename</i>	Reference node name
<i>/Mrnn</i>	Reference node number
<i>/Mrdk</i>	Reference node direction (see table below)
<i>/Msinodename</i>	Response node name
<i>/Msnr</i>	Response node number
<i>/Msdk</i>	Response node direction
	Node Direction id's
0	$-\theta_z$
1	$-\theta_y$
2	$-\theta_x$
3	$-Z$
4	$-Y$
5	$-X$
6	Scalar
7	$+Z$
8	$+Y$
9	$+X$
10	$+\theta_x$
11	$+\theta_y$
12	$+\theta_z$

### Converting Files to Capture File Format (.78C)

Usage: SRT785 /Ocap [*options*] *infile* [*outfile*]

Allowed Input File Types: ASCII

Default Output File Extension: ..78C

Capture output options:

/Bn	Base Frequency n=0 for 100kHz, n=1 for 102.4 kHz (default)
/Ln	Length of Buffer. n is number of kPoints. (default= /L2)
/Nn	Sampling Frequency. Sampling frequency = $256 \text{ kHz}/2^n$ or $262\text{kHz}/2^n$ depending on the base frequency. Default is n=0.
/In	Full Scale Input Range in dBV. Default is n=0.

The input values are ASCII floating point numbers. The numbers represent values in volts and must not exceed the full scale input range set with the /In option. If the number of values found in the input file is less than the number of points specified with the /Ln option, the output file is zero padded.

Example:

Convert the ASCII floating point values found in the file C1.TXT into a 10 kpoint capture file with an input range of 20 dBV (10V) :

**SRT785 /Ocap /L10 /I20 C1.TXT**

### Converting Files to Arbitrary Waveform Format (.78C)

Usage: SRT785 /Oarb [*options*] *infile* [*outfile*]

Allowed Input File Types: ASCII

Default Output File Extension: .78C

Arbitrary Waveform output options:

/Bn	Base Frequency n=0 for 100kHz, n=1 for 102.4 kHz (default)
/Ln	Length of Buffer. n is number of kPoints. (default= /L2)
/Nn	Sampling Frequency. Sampling frequency = $256 \text{ kHz}/2^n$ or $262\text{kHz}/2^n$ depending on the base frequency. Default is n=0.

The input values are ASCII floating point numbers with a range of  $\pm 1.0$ . These values represent the source voltage when the arbitrary source is on with an amplitude of 100%. If the number of values found in the input file is less than the number of points specified with the /Ln option, the output file is zero padded.