


CUBE POLARIZING BEAMSPLITTERS

Features


- Durable and convenient
- Optimised for popular laser wavelengths

Polarizing film is coated on the internal face of a cube beamsplitter. Thin film polarizers utilize the polarization which occurs on reflection deviated by 90° angle. Cube polarizing beamsplitters can be optimized for a particular wavelength to give superior performance for laser application can be optimised for a particular wavelength to give superior performance for laser applications.

The four outer faces are all anti-reflection coated.

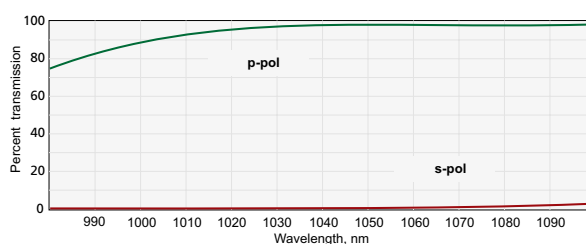

Typical transmission curve @ 1064 nm

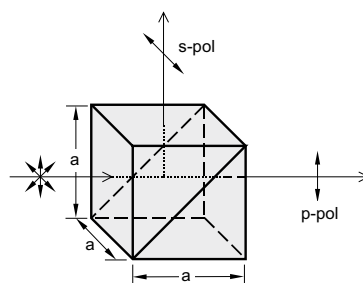

Specifications

Material	BK7, UV FS
Surface quality	40 – 20 scratch & dig (MIL-PRF-13830B)
Surface flatness	$\lambda/4$ @ 633 nm
Clear aperture	90% of the face size
Beam deviation	<3 arcmin
Dimension tolerance	± 0.3 mm
Laser damage threshold	0.3 J/cm ² 10 ns pulses at 1064 nm

Please contact us for polarizing cubes
with an extinction ratio of up to $T_p/T_s > 500:1$

BK7 GLASS, $T_p/T_s > 200:1$

Wavelength range, nm	Reflection s-pol, %	Transmission p-pol, %	Side axa, mm	Unmounted	Mounted
				Catalogue number	Catalogue number
532	$R_s > 99.5$	$T_p > 95$	10 × 10	430-0101	430-0101-M2Ps
532	$R_s > 99.5$	$T_p > 95$	15 × 15	430-0151	430-0151-M2P
532	$R_s > 99.5$	$T_p > 95$	20 × 20	430-0201	430-0201-M2P
532	$R_s > 99.5$	$T_p > 95$	25 × 25	430-0251	430-0251-M2P
633	$R_s > 99.5$	$T_p > 95$	10 × 10	430-0104	430-0104-M2Ps
633	$R_s > 99.5$	$T_p > 95$	15 × 15	430-0154	430-0154-M2P
633	$R_s > 99.5$	$T_p > 95$	20 × 20	430-0204	430-0204-M2P
633	$R_s > 99.5$	$T_p > 95$	25 × 25	430-0254	430-0254-M2P
780	$R_s > 99.5$	$T_p > 95$	10 × 10	430-0102	430-0102-M2Ps
780	$R_s > 99.5$	$T_p > 95$	15 × 15	430-0152	430-0152-M2P
780	$R_s > 99.5$	$T_p > 95$	20 × 20	430-0202	430-0202-M2P
780	$R_s > 99.5$	$T_p > 95$	25 × 25	430-0252	430-0252-M2P
1064	$R_s > 99.5$	$T_p > 95$	10 × 10	430-0103	430-0103-M2Ps
1064	$R_s > 99.5$	$T_p > 95$	15 × 15	430-0153	430-0153-M2P
1064	$R_s > 99.5$	$T_p > 95$	20 × 20	430-0203	430-0203-M2P
1064	$R_s > 99.5$	$T_p > 95$	25 × 25	430-0253	430-0253-M2P

UV FS, $T_p/T_s > 100:1$

Wavelength range, nm	Reflection s-pol, %	Transmission p-pol, %	Side axa, mm	Unmounted	Mounted
				Catalogue number	Catalogue number
308	$R_s > 99$	$T_p > 90$	10 × 10	430-1105	430-1105-M2Ps
308	$R_s > 99$	$T_p > 90$	15 × 15	430-1155	430-1155-M2P
308	$R_s > 99$	$T_p > 90$	20 × 20	430-1205	430-1205-M2P
355	$R_s > 99$	$T_p > 90$	10 × 10	430-1107	430-1107-M2Ps
355	$R_s > 99$	$T_p > 90$	15 × 15	430-1157	430-1157-M2P
355	$R_s > 99$	$T_p > 90$	20 × 20	430-1207	430-1207-M2P

Please contact us if you need polarizing beamsplitters of other wavelengths, other sizes or other configurations.

Housing accessories

for unmounted cube polarizing beamsplitters


Prism Holders

840-0160, 840-0170

Find more at
EksmaOptics.com


Mounting Suggestion

Adjustable Polarizer Holder of Side Driver 840-0195

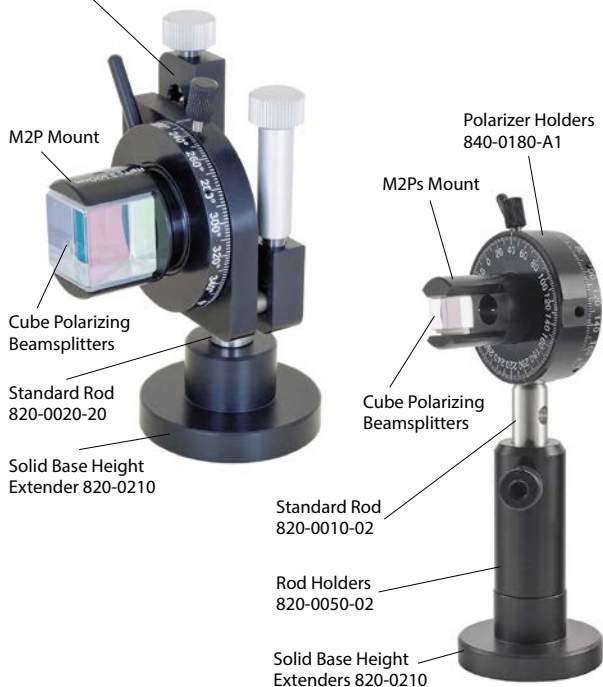


Mounting

Example of **M2P** mount for Cube Polarizing Beamsplitters of 20×20 mm standard dimensions.

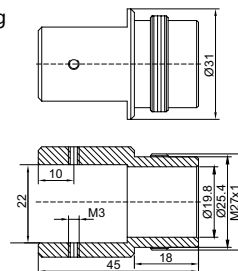

M2P mount can be used with 840-0180, 840-0195, 840-0020.

Example of **M2Ps** mount for Cube Polarizing Beamsplitters of 10×10 mm dimensions.

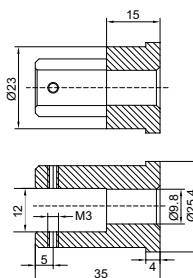

M2Ps mount can be used with 840-0180, 840-0020.


Drawings of M2P and M2Ps for prisms of other dimensions are available on request.

HIGH ENERGY POLARIZING CUBE BEAMSPLITTERS

Features

- Optically contacted for high power applications
- Typical damage threshold: >15 J/cm², 10 ns, 10 Hz at 1064 nm
- Precision surface quality
- High extinction ratio

The four outer faces are all anti-reflection coated.

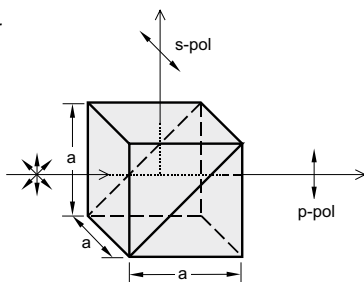

Specifications

Material	UV FS (laser line) N-SF2 (broadband) N-SF15 (dual wavelength)
Size tolerance	±0.2mm
Extinction ratio	>1:500
Surface quality	20 – 10 scratch & dig (MIL-PRF-13830B)
Wavefront distortion	λ/10 at 633 nm
Beam deviation	<3 arcmin
Clear aperture	>85% of size
Reflection s-polarization	>99.5%
Transmission p-polarization	>97% for laser line
	>90% for broadband

Housing accessories

for unmounted cube polarizing beamsplitters

Prism Holders 840-0160, 840-0170

Find more at EksmaOptics.com



Operating wavelength, nm	Material	Laser Damage threshold, 10 ns, 10 Hz	Side a × a, mm	Unmounted	Mounted
				Catalogue number	Catalogue number
343	UV FS	>3 J/cm ² at 343 nm	12.7 × 12.7	435-1128	435-1128-M2Ps
355	UV FS	>3 J/cm ² at 355 nm	12.7 × 12.7	435-1127	435-1127-M2Ps
515	UV FS	>6 J/cm ² at 515 nm	12.7 × 12.7	435-1128	435-1128-M2Ps
532	UV FS	>6 J/cm ² at 532 nm	12.7 × 12.7	435-1121	435-1121-M2Ps
800	UV FS	>8 J/cm ² at 800 nm	12.7 × 12.7	435-1122	435-1122-M2Ps
800	UV FS	>8 J/cm ² at 800 nm	25.4 × 25.4	435-1222	435-1222-M2P
1030	UV FS	>15 J/cm ² at 1030 nm	12.7 × 12.7	435-1124	435-1124-M2Ps
515 + 1030	N-SF15	>10 J/cm ² at 1030 nm	25.4 × 25.4	435-1224D	435-1224D-M2P
1064	UV FS	>15 J/cm ² at 1064 nm	12.7 × 12.7	435-1123	435-1123-M2Ps
420 – 680	N-SF2	>1 J/cm ² at 532 nm	12.7 × 12.7	436-1121	436-1121-M2Ps
700 – 1080	N-SF2	>2 J/cm ² at 1064 nm	12.7 × 12.7	436-1123	436-1123-M2Ps



グローバルにネットワークを広げ、最先端の科学をお客様に提供

本社：〒134-0088 東京都江戸川区西葛西6-18-14 T.I.ビル

Tel. 03-3686-4711

営業所：〒532-0003 大阪府大阪市淀川区宮原4-1-46 新大阪北ビル

Tel. 06-6393-7411

URL: <https://www.tokyoinst.co.jp> Mail: sales@tokyoinst.co.jp



超高真空・極低温走査型プローブ顕微鏡
高速分光測定装置、クライオスタット



Nd:YAGレーザー、Ti:Sレーザー
OPOLレーザー

- 本カタログに記載されている内容は、改良のため予告無く変更する場合があります。（製品の仕様、性能、価格などはカタログ発行当時のものです）
- 本カタログに記載されている内容の一部または全部を無断で転載することは禁止されております。
- 本カタログに記載されているメーカー名、製品名などは各社の商標または登録商標です。