



*9103 shown with HV option

description

The 9103 is a versatile, general-purpose picoammeter designed to accurately measure DC current from the nA to mA range. The easy-to-read user interface is designed to provide optimal control and quick access to all functions. Features such as an ASCII programming language and data recording/graphing are standard.


The 9103 is available with high-speed and high-voltage options, for sample rates up to 500 reads/sec, and the ability to float up to 5 kV DC.

applications

- Electron & ion beam current measurements
- Mass spectrometer current measurements
- Photodiode current & leakage measurements
- Beam and particle monitoring
- Spectrometer measurements
- Data logging current comparisons



benefits

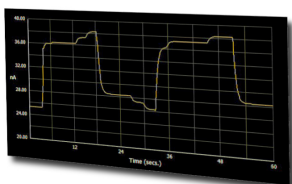

USB Interface: PC software included

Data Recording: Sample data can be recorded and saved

Data Graphing: Real-time graphing with polarity options

Programmable: ASCII command interface

Sync: Multiple 9103's can be synchronized using RBD's Actuel software to create a multichannel current data logger

Bias (Optional): Increases accuracy of electron or ion beam current measurements by reducing secondary electron emission (*not available on High-Voltage models*)

configurations

The 9103 is available in 4 configurations:

- **Standard:** sample rate 40 reads/sec.
- **High-Speed:** sample rate 500 reads/sec.
- **High-Voltage:** sample rate 40 reads/sec., floats up to 5 kV DC
- **High-Speed / High-Voltage:** sample rate 500 reads/sec., floats up to 5 kV DC

specifications

Resolution & Accuracy per Range

Range	Resolution	Accuracy / Offset (18° – 28°C, 0 - 70% RH)*	
2 nA	100 fA	0.5 ± %	+500 fA
20 nA	1 pA	0.5 ± %	+3 pA
200 nA	10 pA	0.4 ± %	+20 pA
2 µA	100 pA	0.4 ± %	+200 pA
20 µA	1 nA	0.4 ± %	+2 nA
200 µA	10 nA	0.4 ± %	+20 nA
2 mA	100 nA	0.4 ± %	+200 nA

*Temperature Coefficient: 0°–18°C & 28°–50°C. For each °C, add $0.1 \times (\% \text{ rdg} + \text{offset})$ to accuracy spec.

General:

USB Interface	160 mW power consumption
2 nA Min. DAC Resolution	2 nA Min. : 0.1 pA, 2 mA Min. 100 nA
Sampling Rate	Standard: 40 reads/sec. High-Speed option: 500 reads/sec.
Range	2 nA to 2 mA with 100 fA resolution
Input Protection	50 mA (fused)
Recorder Output Voltage	± 2.000 V per range <i>not available with High-Voltage option</i>
Warm-up	1 hour to rated accuracy
Environment	Operating 0° C to 50 ° C
Mechanical Dimensions	55H x 170W x 165L (mm) 2.11H x 6.68W x 6.30L (inches)
Net. Weight	0.816 kg / 1.5 lbs
Connections <i>Standard-Voltage</i>	Input: Isolated BNC. Recorder Option Analog Output: Banana jacks
Connections <i>High-Voltage</i>	HV Input: SHV HV Reference: SHV
Bias Options	No Bias; Internal (fixed +90 V DC); External (Dual BNC) <i>not available with High-Voltage option</i>
Safety	Conforms to USB Power Specification - for use by qualified personnel who are trained in the use of test and measurement instrumentation
High Voltage Option	5 kV DC Isolation
Accessories Included	Installation and operation manuals, USB driver,, Actuel software, USB power cable, low-noise BNC cable (<i>not available with HV option</i>)

All specifications are subject to change without notice.





グローバルにネットワークを広げ、最先端の科学をお客様に提供

本 社：〒134-0088 東京都江戸川区西葛西6-18-14 T. ビル

Tel. 03-3686-4711

営業所：〒532-0003 大阪府大阪市淀川区宮原4-1-46 新大阪北ビル

Tel. 06-6393-7411

URL: <https://www.tokyoinst.co.jp> Mail: sales@tokyoinst.co.jp



超高真空・極低温走査型プローブ顕微鏡
高速分光測定装置、クライオスタット



Nd:YAGレーザー、Ti:Sレーザー
OPOLレーザー

- 本カタログに記載されている内容は、改良のため予告無く変更する場合があります。（製品の仕様、性能、価格などはカタログ発行当時のものです）
- 本カタログに記載されている内容の一部または全部を無断で転載することは禁止されております。
- 本カタログに記載されているメーカー名、製品名などは各社の商標または登録商標です。